H 1000

| INSTALLATION MANUAL


INSTALLATIEHANDLEIDING

MANUAL D'INSTALLATION

INSTALLATIONANLEITUNG

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MANUAL DE INSTALACION

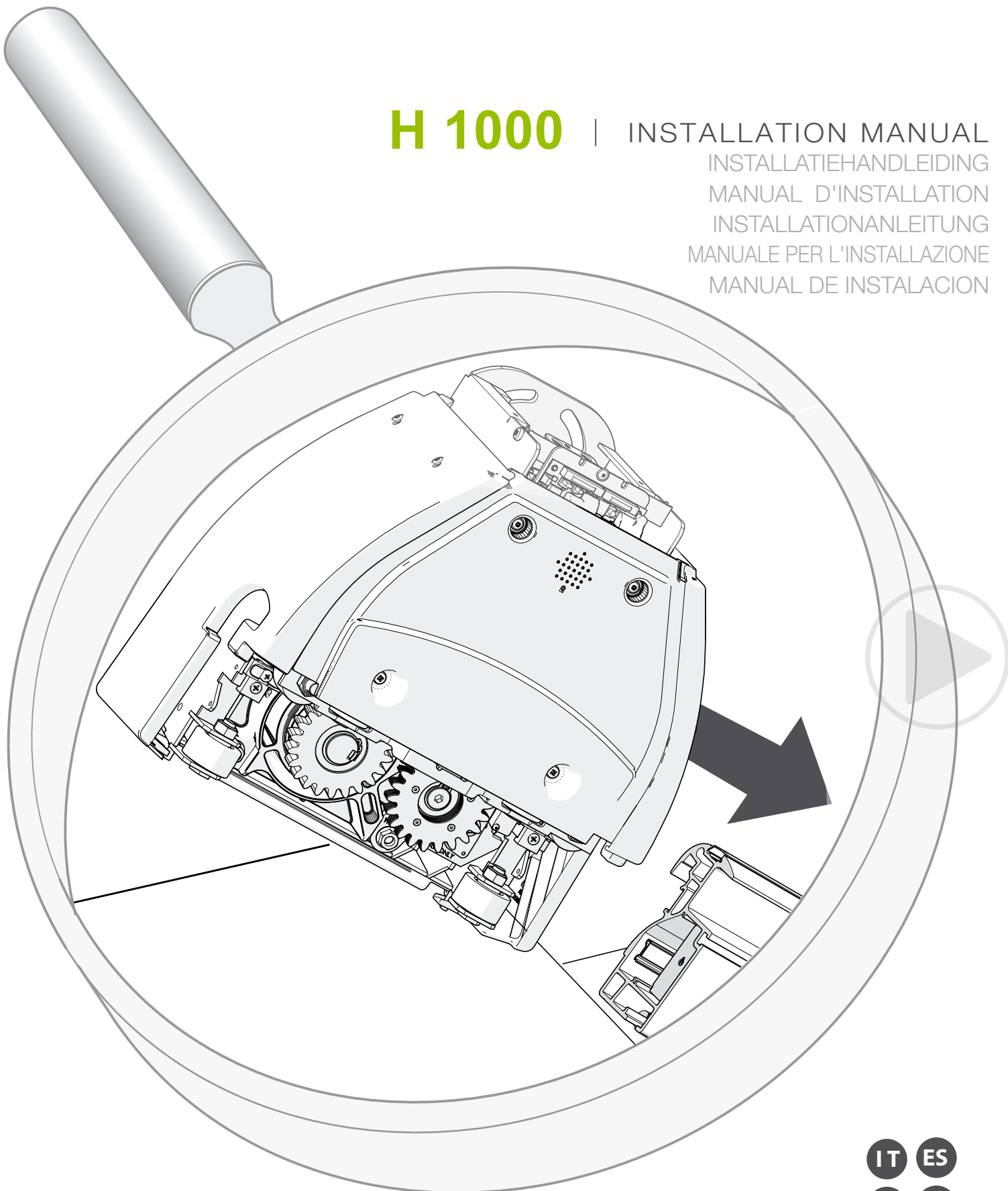Manual Smart
Powered Smart

 **savaria**®

EN All rights reserved. No part of this publication may be reproduced and/or published by means of print, photocopy, microfilm or in any other way or form, without prior written permission from the manufacturer. This also applies for the accompanying drawings and diagrams. The manufacturer reserves the right to change parts at any time, without prior or direct notification to the customer. The contents of this manual may also be changed without prior warning. This manual is intended for the standard models of the lift. The manufacturer therefore cannot be held liable for any damage resulting from specifications that differ from those of the standard model of the lift supplied to you.

This manual has been compiled with the best intentions and every possible care but the manufacturer cannot accept responsibility for any errors in this manual or for the consequences thereof.

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EN
The English manual is the original language

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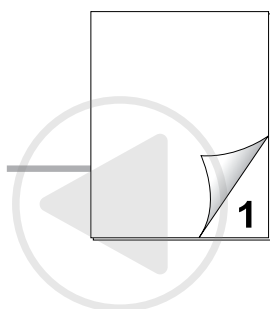
1 - USED TOOLS AND PARTS


2

Thank you for choosing a Handicare stairlift.
This manual is for the Handicare1000 lift with a manual seat or powered seat type. The manual includes instructions for the standard and the optional features.

2 - RAIL

6

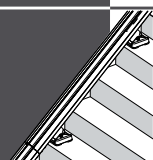
Cut [p.06]
Clean [p.08]
Joint [p.10]
Foot [p.12]
Charge wire [p.14]

Earth wire [p.16]
Wiring upstairs [p.18]
Fixation [p.20]
Power supply [p.22]

3 - POWERPACK


24

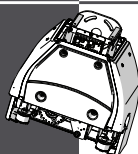Preparing [p.24]
Placing [p.26]
OSG test [p.28]
Disconnection shorting links [p.30]

4 - FINISH RAIL


32

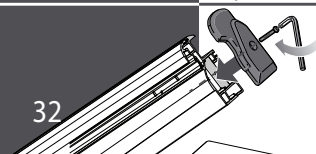Charger upstairs [p.32]
Endstop [p.34]

5-CHASSIS


36

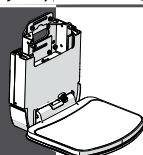Placing chassis [p.36]
Placing seat base [p.38]
Fixation seat base [p.40]
Fixation chassis [p.42]

6- SEAT


44

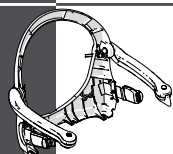Fixation seat back [p.44]

7- WIRING


46

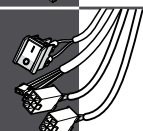Seat loom [p.46]
Communications [p.48]
Seat base wiring [p.50]

8- FINISH

54

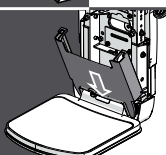
Bottom front cover [p.54]
Connections top front cover [p.56]
Fixation top front cover [p.58]
Upholstery [p.60]
Track [p.62]

Vertical alignment [p.64]
Adjustments [p.66]
Adjustments [p.68]
Hinged rail [p.70]

9 - REMOTE CONTROLS

72

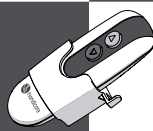Mounting [p.72]
Check linking [p.74]
Changing channels [p.76]
Linking [p.78]

10 - CHECK

80

Measurements [p.80]
Swivel - OSG safety [p.82]
Option functions [p.84]
Charge points [p.86]
Safety stops [p.88]

Armrests [p.90]
Keyswitch - stops [p.92]
Instruct user [p.94]
Declaration [p.96]
Placing label [p.98]

In this manual, the paperclip symbol refers to a relevant attachment.

ATTACHMENTS

100

A - Rail cuttings [p.101]
B - 3 rail parts [p.103]
C - Emergency operation(track) [p.104]
D - Emergency operation(seat) [p.106]
E - Foot conversion[p.108]
F - PERCH seat bracket[p.110]
G - Outdoor cover [p.112]

H - OSG test [p.116]
I - Reinforcement plate [p.122]
J - TorqueWrenchTest [p.124]
K - Protection equipment [p.125]
L - Competence [p.126]
M - Installation forms [p.129]

STAIRLIFT TIGHTENING TORQUE

UNLESS OTHERWISE STATED



Size (mm)	Material	Max. Torque (Nm)					
		Bolt Grade					
		4.6	4.8	5.8	8.8	10.9	12.9
M6	Self Colour	3	3.8	4.8	8.6	12	14.5
M6	Zinc	3.4	4.3	5.6	9.9	14	16.5
M8	Self Colour	7.2	9.1	11.7	21	29	35
M8	Zinc	8.3	10.5	13.6	24	34	40
M10	Self Colour	14.4	18	23.4	42	58	70
M10	Zinc	16.6	21	27	48	67	81

Nm

HAND TIGHT

Not loose, but not really tight either.
Somewhere around 1 or 2 Nm.

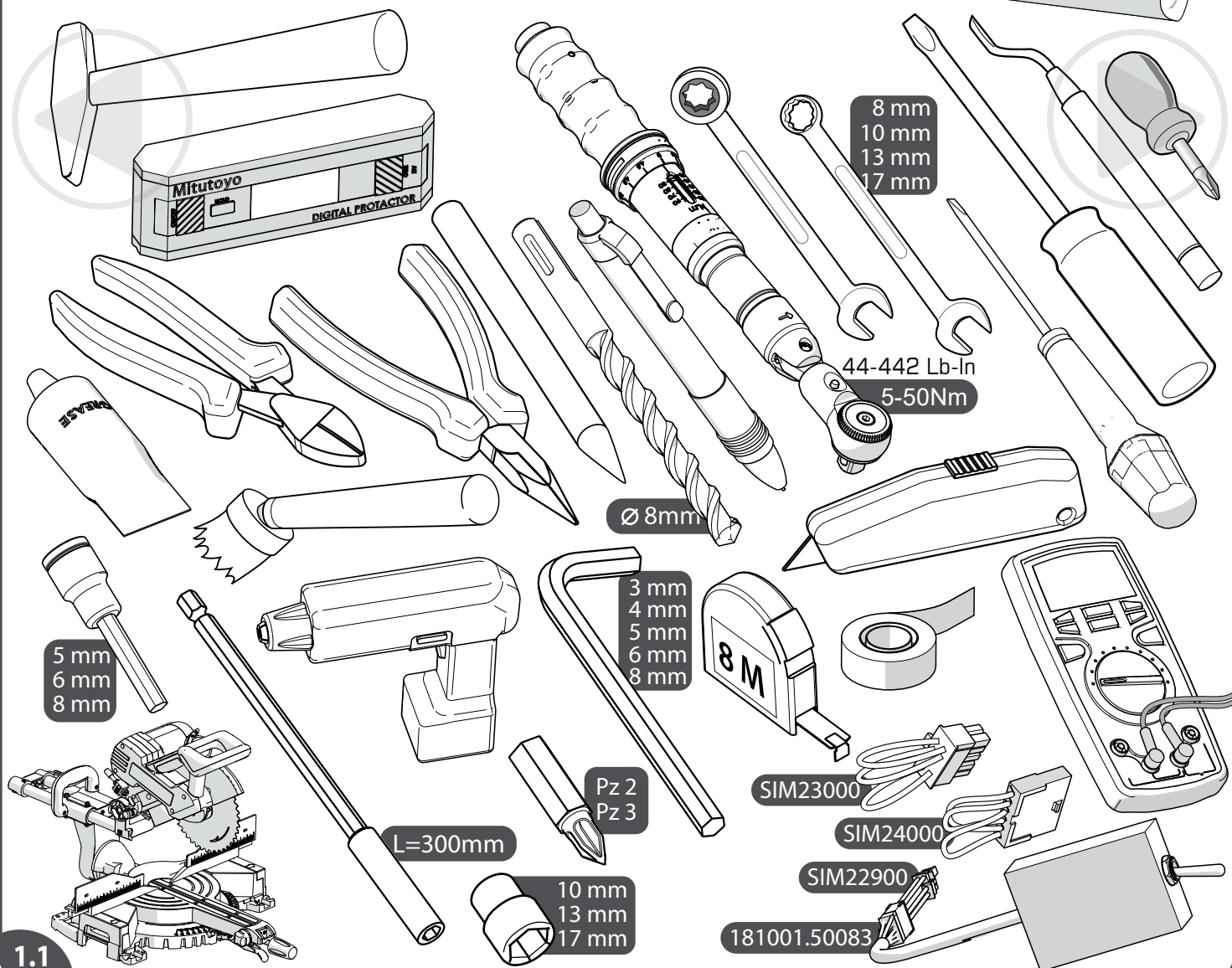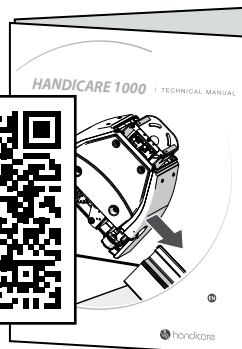Only a technician certified for this lift should install this lift (for more detailed information about the required competence see attachment H of this manual). **General required knowledge can be found in the 1000 Technical Manual. The latest instructions for the lift and instructions for options are available for download at www.handicarepartners.com.**

Ensure good and sufficient tools; minimum required, in addition to installation kit, to install the stairlift:

File	Allen keys {2,5/3/4/5/6/8 mm}	Pencil	Crimper	Voltmeter
Hacksaw	Tape-measure {min. 8 m}	Sockets {10/13/17 mm}	Pliers	Hammer
Drill + drills {8 mm}	Spanners {8/10/13/17 mm}	Shorting connectors	Brush	Marker pen
Insulation tape	Screwdrivers with fine flat head	Stanley knife	Loading toggle	Spirit level
EN Pozi-drive bits {Pz 2/3}	Bit extenders {300 mm}	Side cutters	Torque wrench	Grease

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
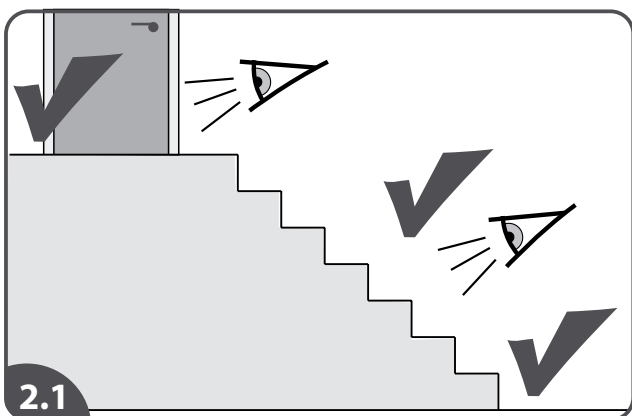
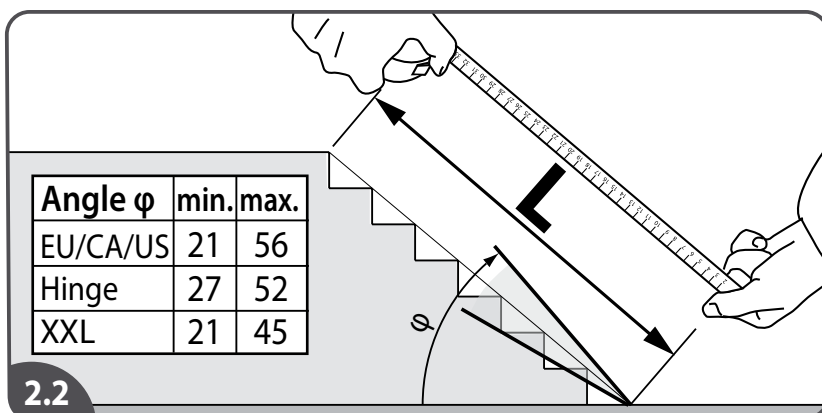
Diagram B is a circular inset showing a close-up of the saw blade cutting through a workpiece. A curved arrow indicates the direction of the cut, which is angled upwards. The workpiece is held in place by a clamp.

..... **mm**

Waste

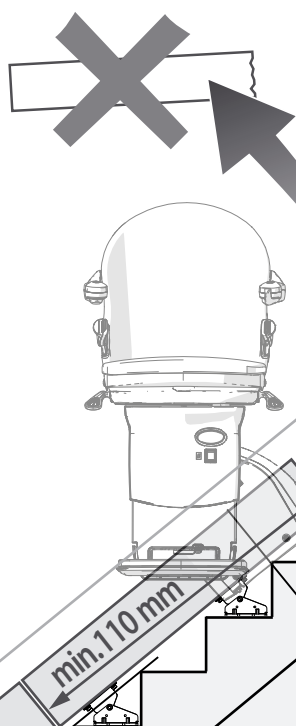2.3

R

REMOVE CONTACT AND STRIP

**REMOVE
MAGNET**

ALUMINIUM CUTTING

$$R - (C + D) - L =$$

Standard finish:
footplate horizontally
aligned to the top nose

*** FOR MORE ACCURATE C&D**

[2.1] Before you start, inspect the stairs. Clarify with your customer regarding any floor heating, alarm wires, pads or pipework that may affect/hinder installation prior to starting.

[2.2] Measure the distance between the upper floor and the ground floor (diagonally). This is the length **L**.

[2.3] Check that there are no obstacles on the first floor, close to the stairs. The length to be sawn is rail length(**R**) - **L** - (**C+D**). Saw the rail part to this length, at right angles. You can use an average value for **C+D**. For a more accurate length **C+D**, please refer to the table in Attachment A. Remove the charge contact, cover and magnet before cutting the rail.

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





IT

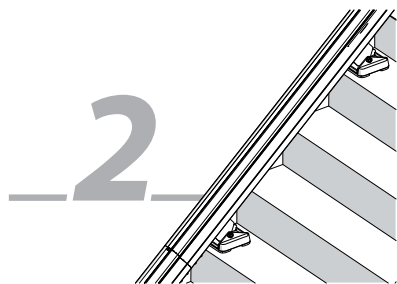RAIL
CLEAN

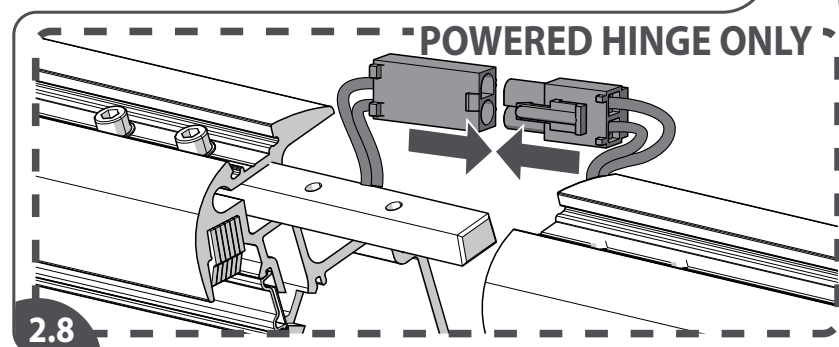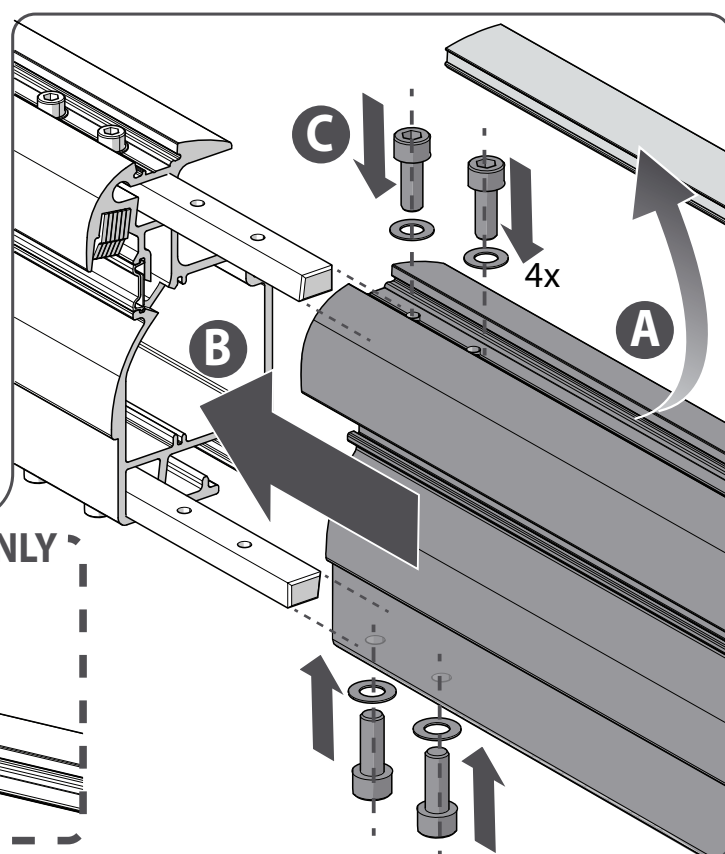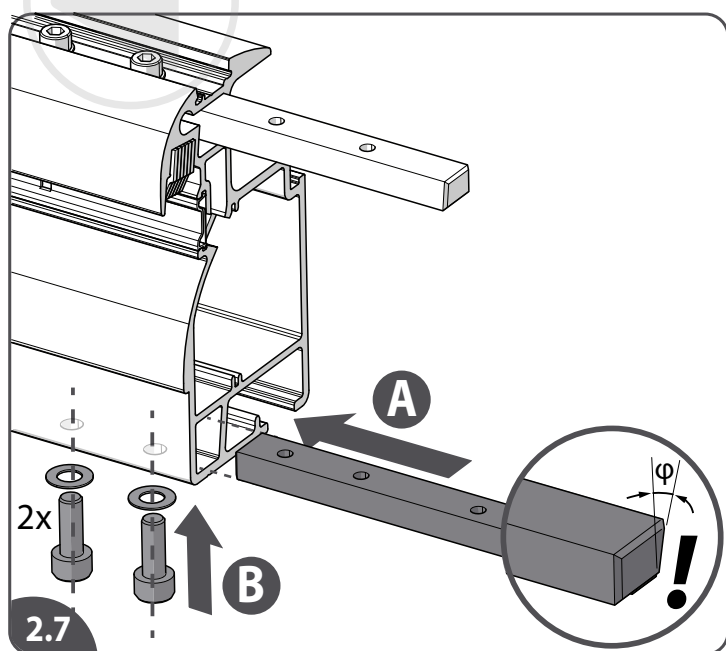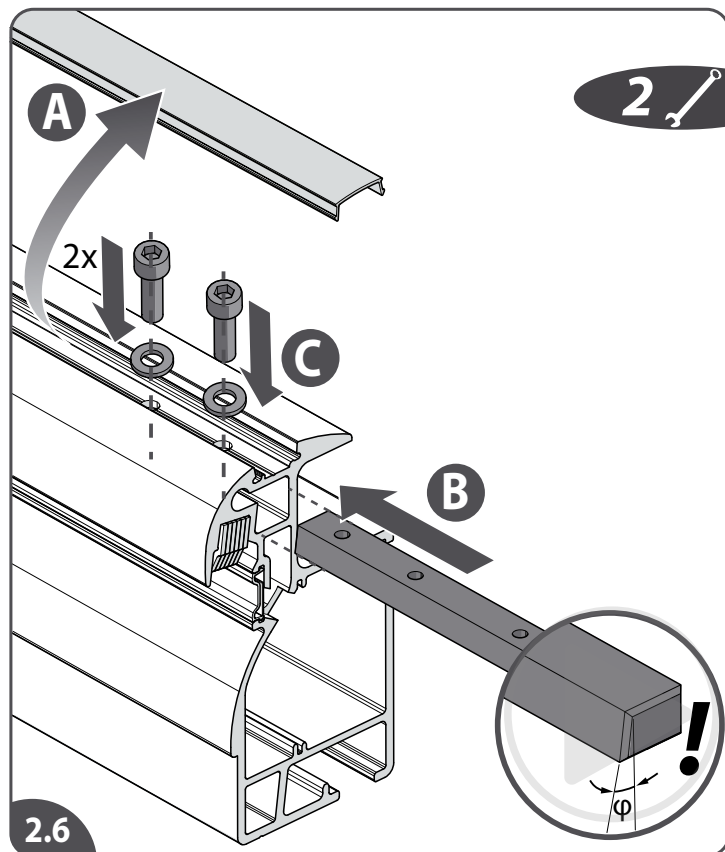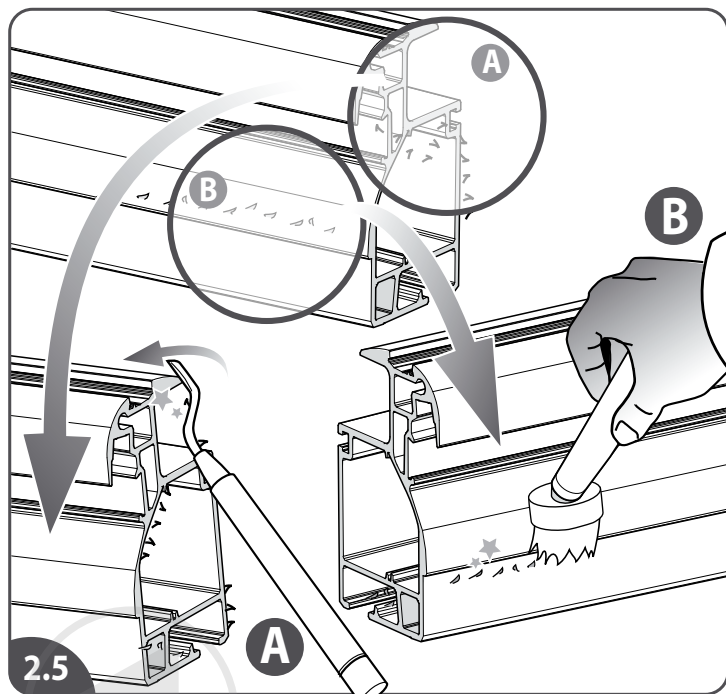- [2.5] Deburr the sharp edges on the sawn end {A}. Remove as much debris as possible from the rail {B}.
- [2.6] Remove the track cover {A}. Fit the first splice bar by sliding the bar into the bottom rail, observing the correct orientation {B}. Fit the bolts in the bottom rail {C}. Do not fit the top bolts in the splice bar: these should be fitted later.
- [2.7] Fit the second splice bar by sliding the bar into the bottom rail, observing the correct orientation {A}. Fit the bolts in the bottom rail {B}.
- [2.8] Remove the track cover {A}. Fit the top rail part tight up against the lower rail part {B}. Fit the bolts (4x) in the top rail {C}.

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





IT

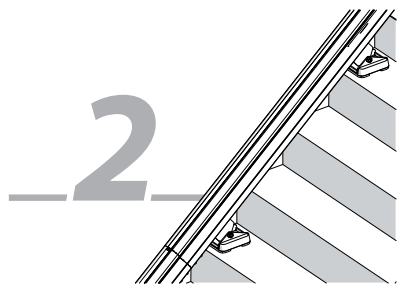RAIL JOINT

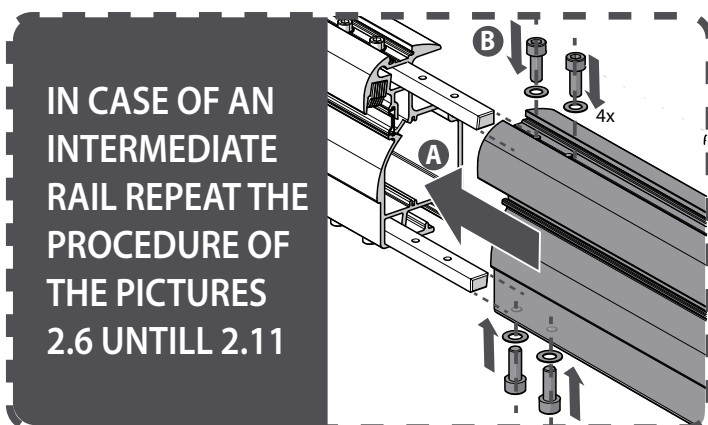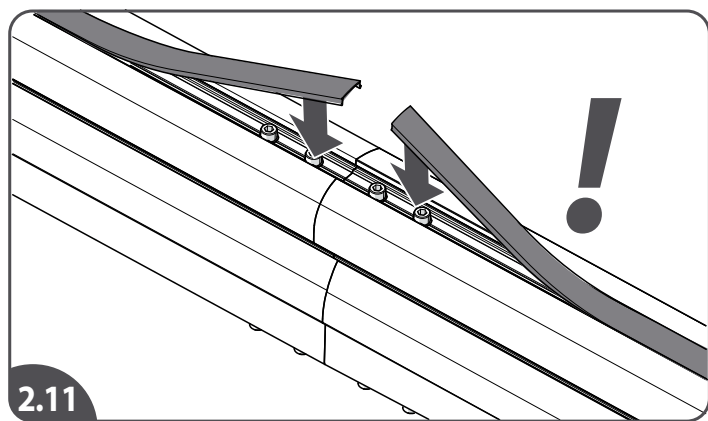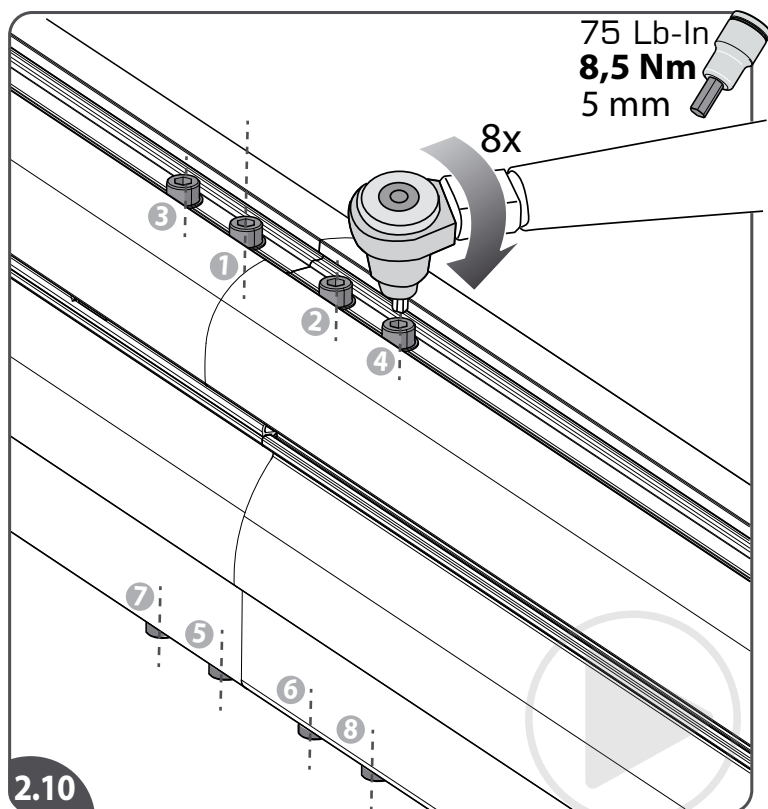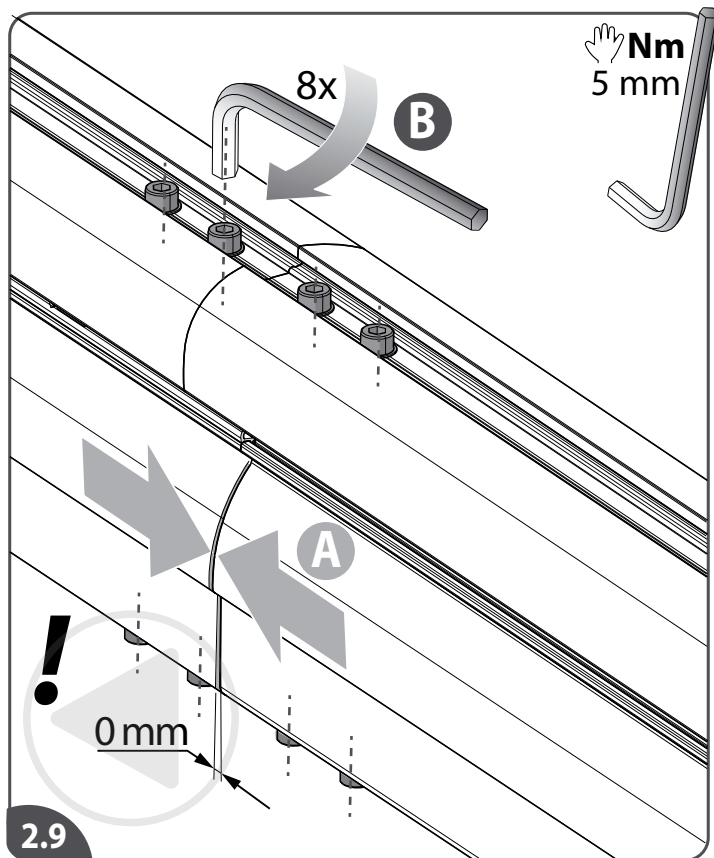IN CASE OF AN
INTERMEDIATE
RAIL REPEAT THE
PROCEDURE OF
THE PICTURES
2.6 UNTILL 2.11

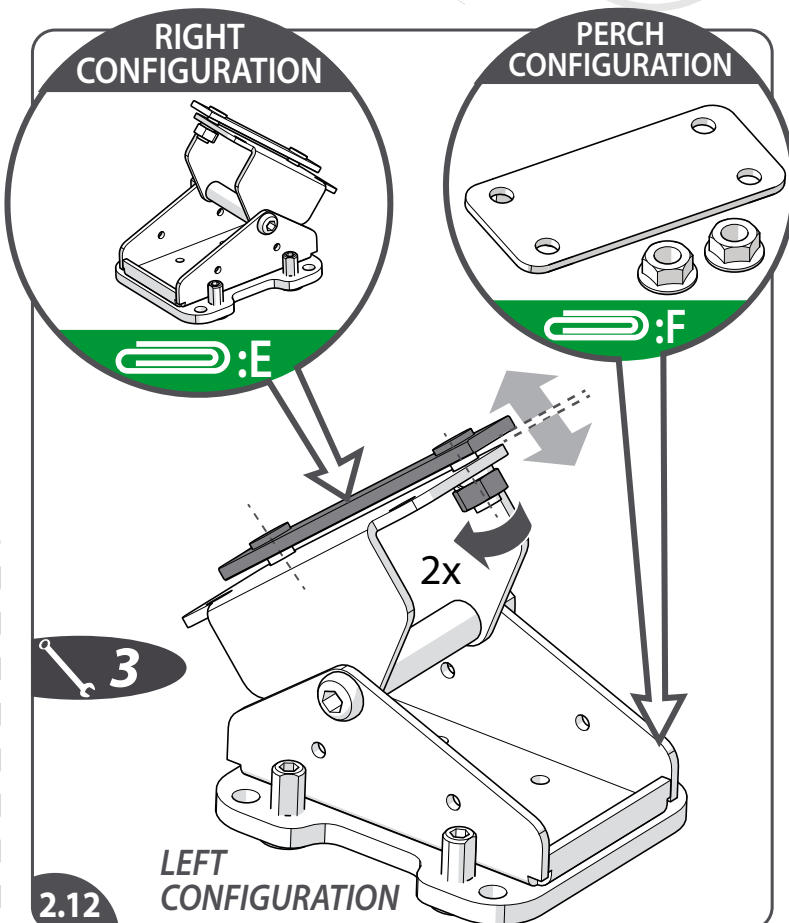- [2.9] Fit the other rail part tight up against the lower rail part {A}. Tighten (8x) the bolts securely {handtight} {B}.
 [2.10] Fasten the top and bottom splice bar by alternately tightening the bolts {8,5 Nm | 75Lb-In}.
 [2.11] Reinstall the track covers.

Note: be sure track covers are flush. It's easy to leave them slightly proud and the micro switches on the pack can rub and get damaged.

In case of an intermediate rail repeat the procedure of the pictures 2.6 until 2.11.

- [2.12] Adjust the brackets (loosen the bolts a bit). If needed change the configuration of the brackets (see attachment E (right configuration) or attachment F (PERCH seat configuration)).

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





IT

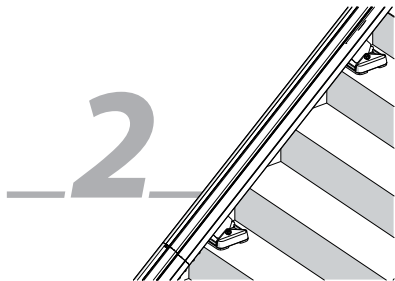RAIL FOOT

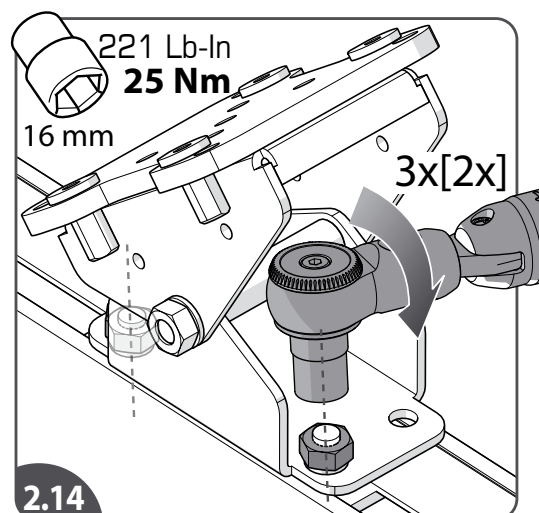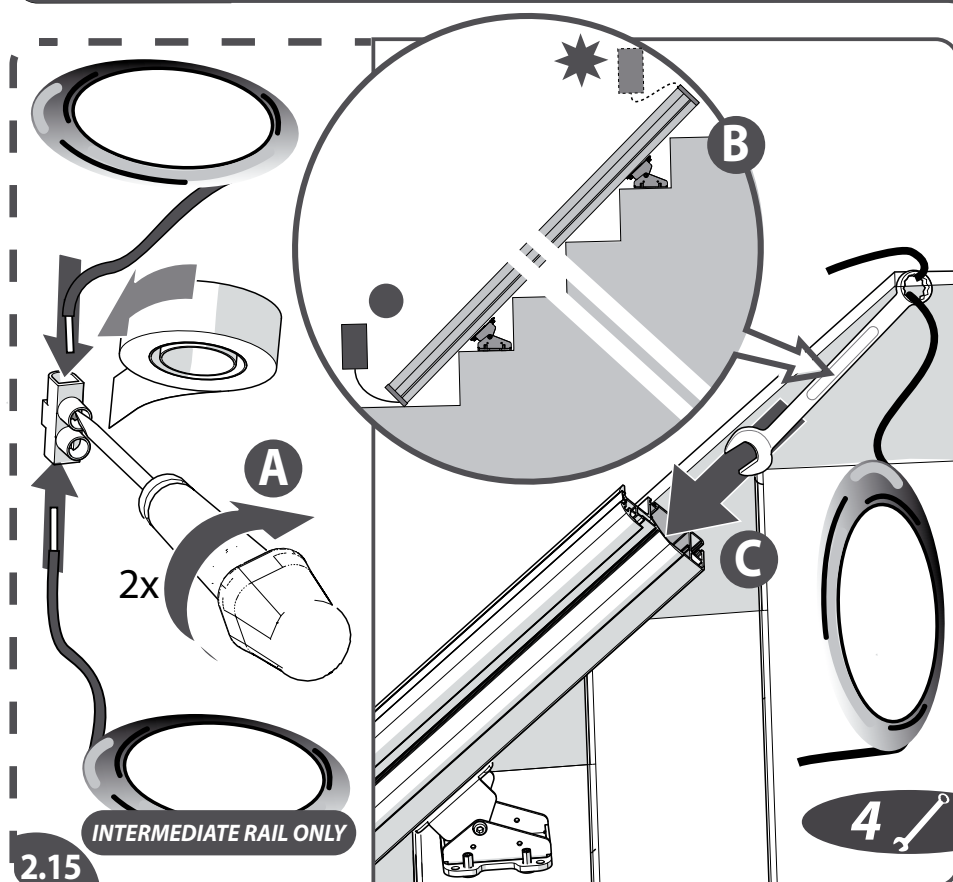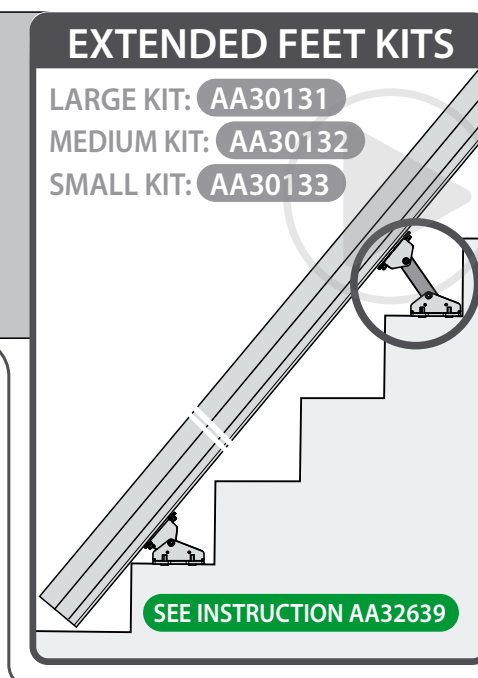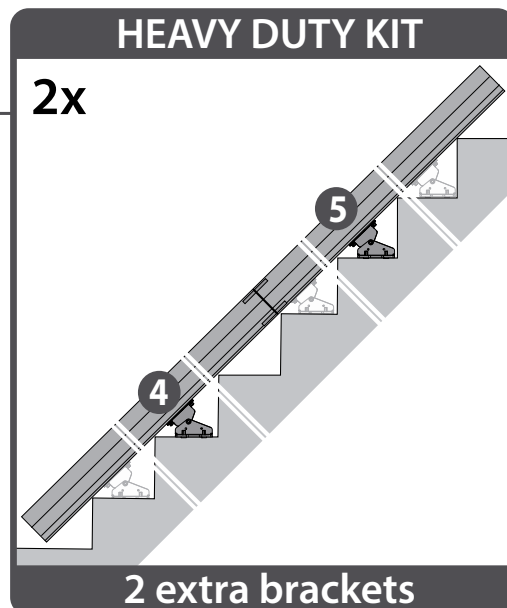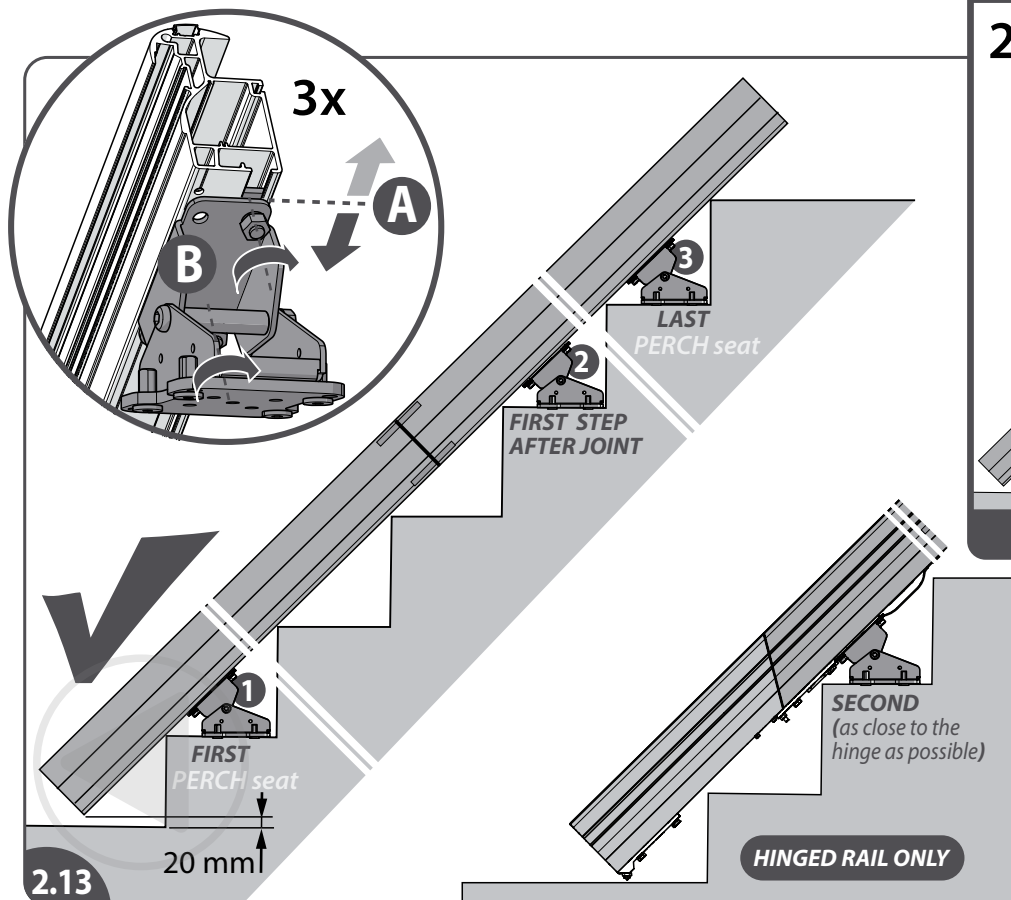[2.13] Place the bracket on the rail{A}. Position the bracket on the rail{B}. Position the brackets as shown (brackets on the first/last step of the stairs and a bracket on the first step after the joint). In case of uneven steps or desired larger rail angle, use the extended feet (see instruction AA32639).

[2.14] Tighten the bolts of the brackets securely (25 Nm/221 Lb-In).

[2.15] Choose the best location for the charger, either at the top or bottom of the stairs{B}. Connect the supplied internal charge wire with a spanner and drop it into the main chamber of the rail{C}. In case of a hinged rail the charger needs to be installed at the top of the stairs. In case of an intermediate rail the two supplied charge wires needs to be connected together and the terminal block taped{A}.

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IT

[2.16] Remove the insulation of the charge wires (15 mm).

[2.17] Split the duo connector block in two{A}. Cut off the connectors of the charger{B}. Remove the insulation of the wires of the charger (15 mm).

[2.18] Twist the two ends bare wires together{A}. Fasten these wires in a connector block{B}. In case the charger is mounted downstairs fasten also the +wire (white) of the charger in this block{C}.

[2.19] Insulate the connector block.

[2.20] In case the charger is mounted upstairs: put the wires in the main chamber of the rail{A}. Slide the end cap bracket into the rail{B}.

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




IT

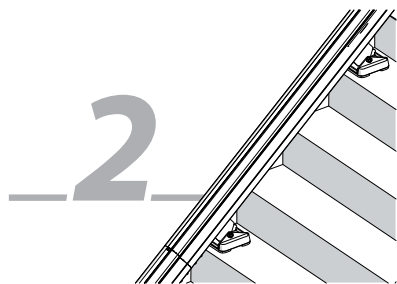RAIL EARTH WIRE

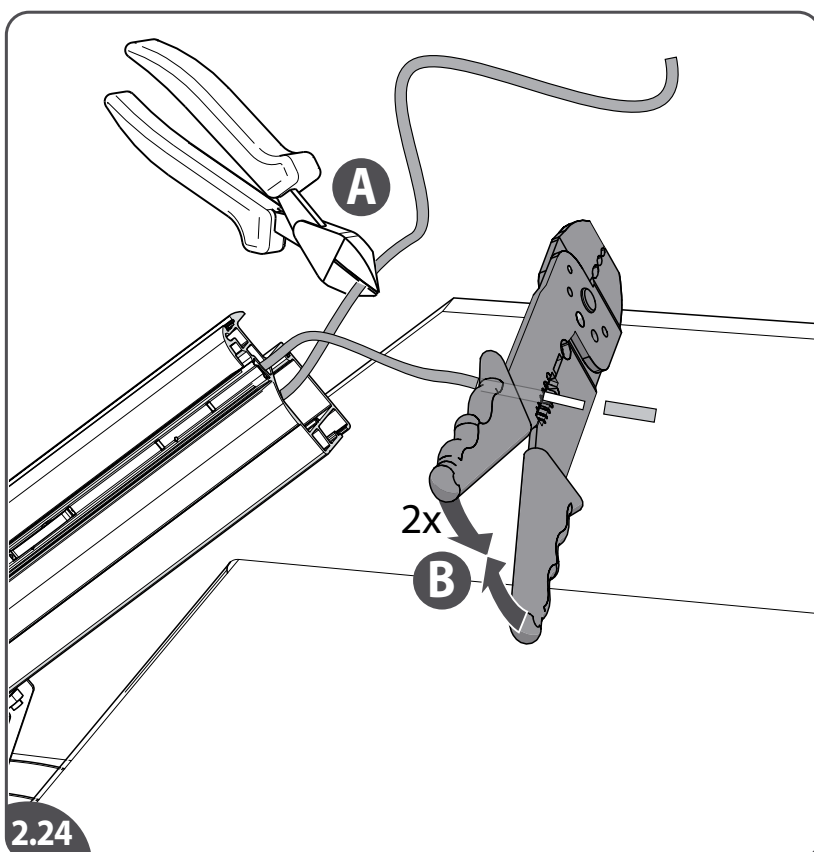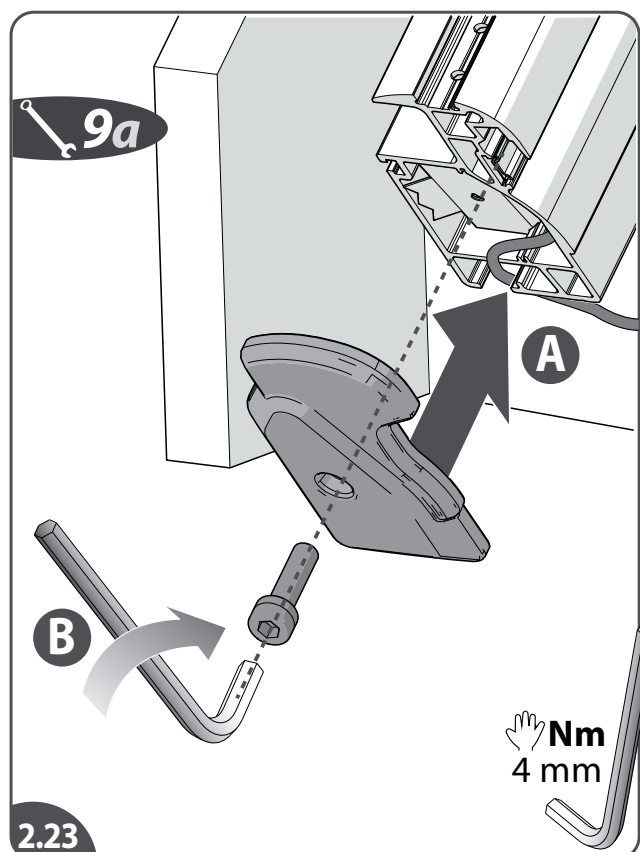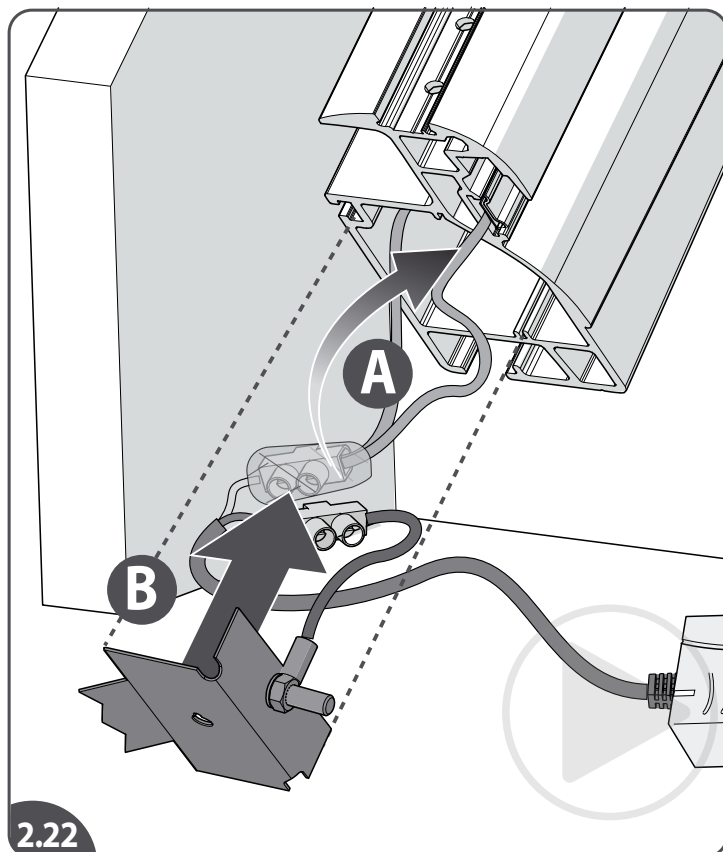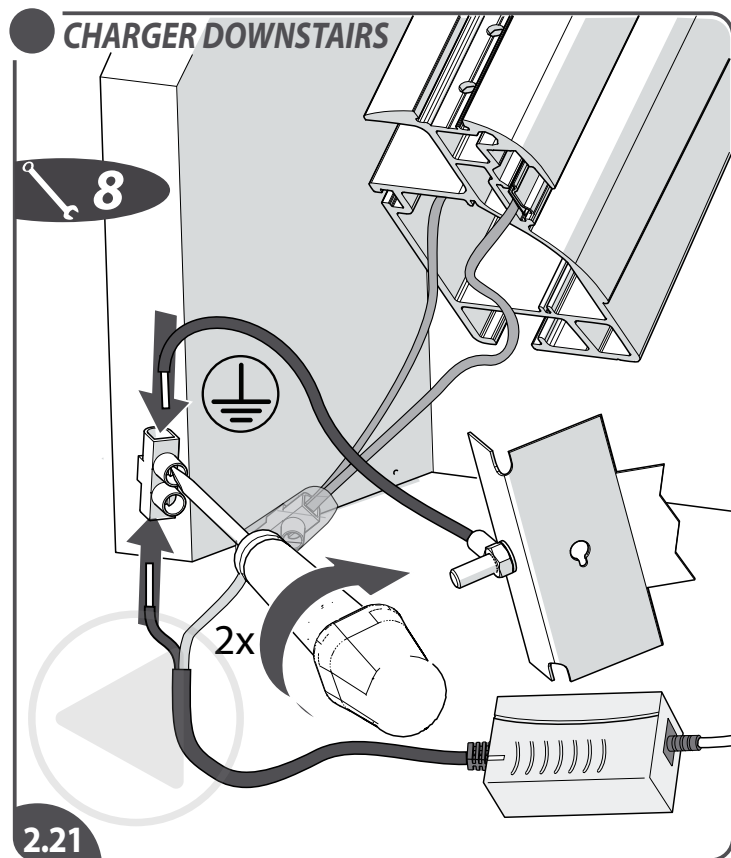- [2.21]The charger downstairs: fasten the -wire (black) of the charger and the earthwire in another connector block.
 [2.22]Put the wires in the main chamber of the rail{A}. Slide the end cap bracket(earth) into the rail{B}.
 [2.23]Install{A} and carefully fasten (handtight) the end cap{B}.
 [2.24]Cut the excess wire from the internal charging wire {A}. Remove the insulation of both wires (15 mm){B}.

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




IT

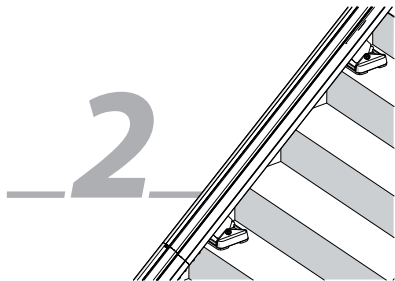RAIL WIRING UPSTAIRS

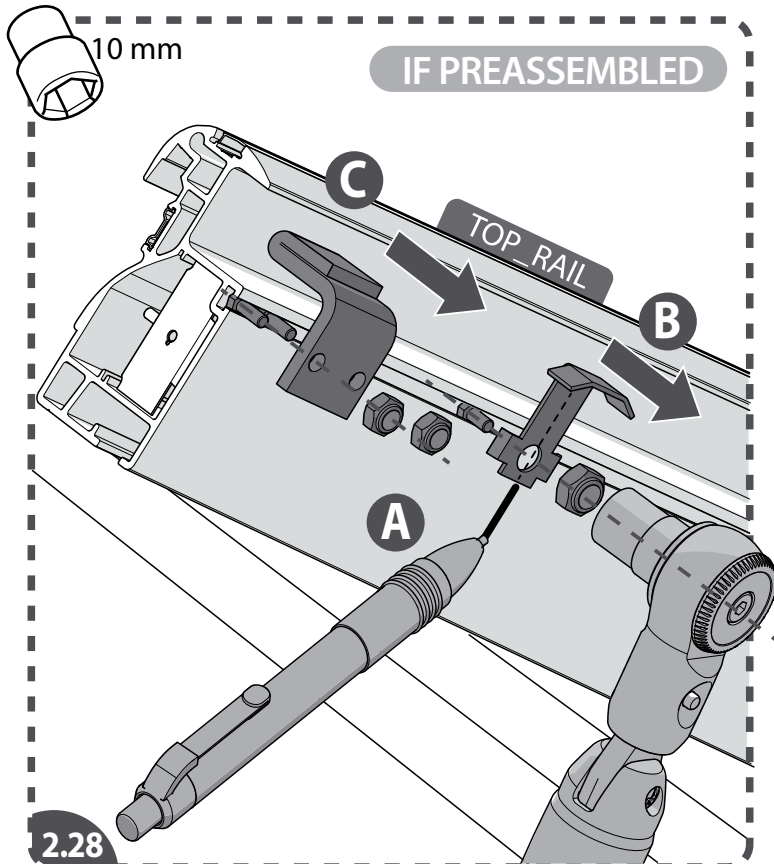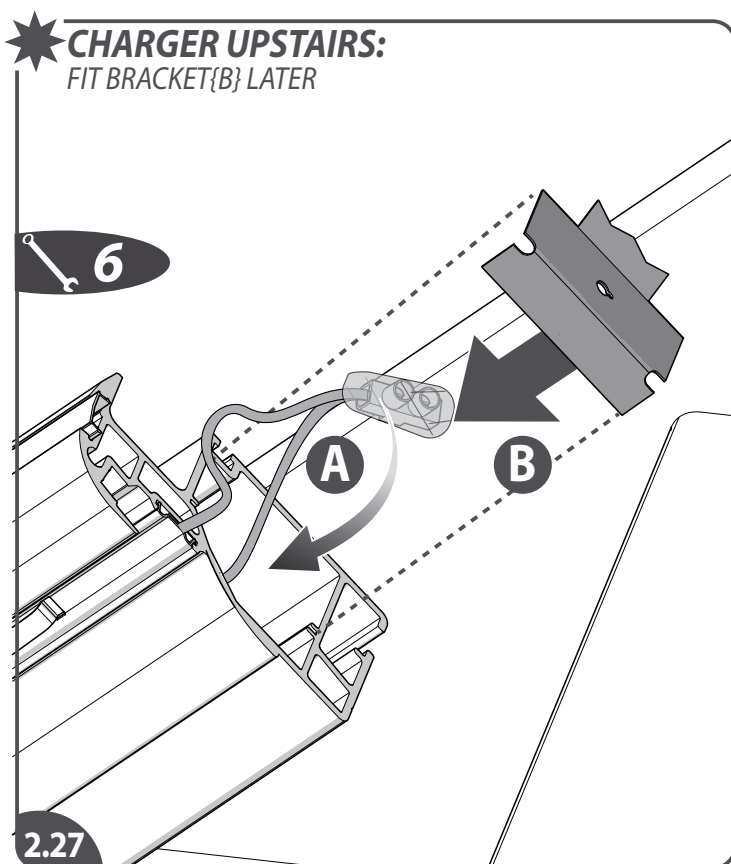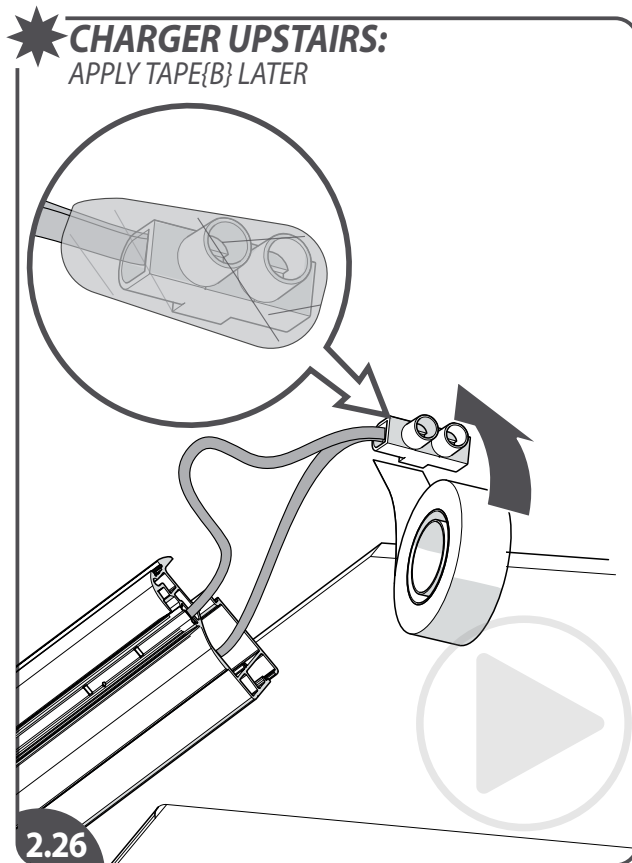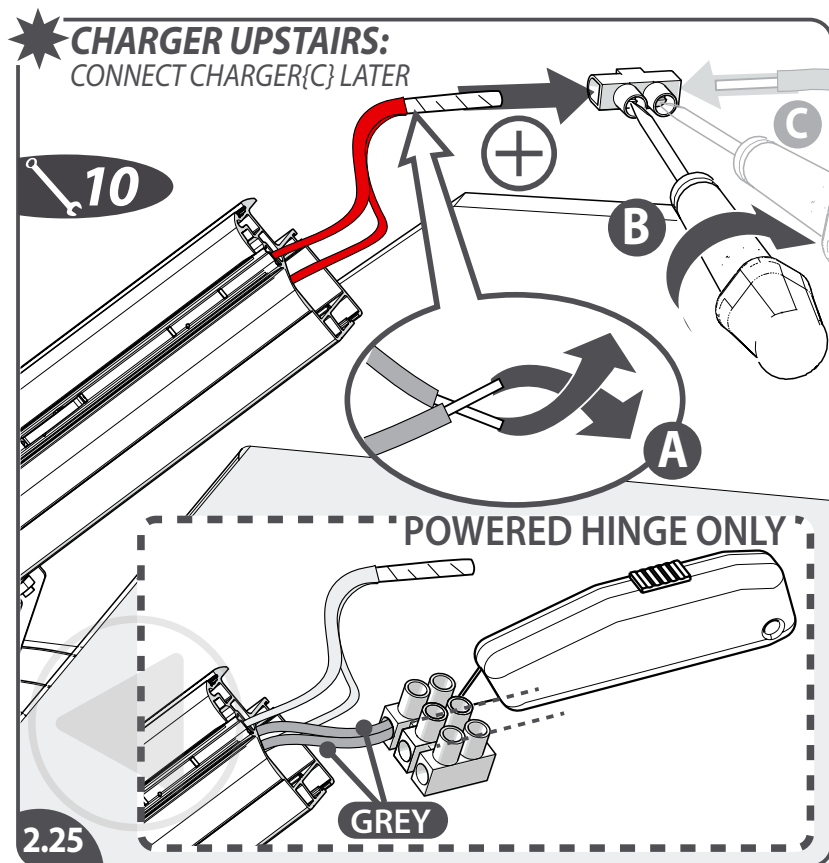[2.25] Only in case of a powered hinge: split the 3-way connector block in three parts. Twist the two ends bare wires together{A}. Fasten these wires in a connector block{B}. In case the charger is mounted upstairs: connect the charger later{C}.
 [2.26] Insulate the connector block.
 [2.27] Put the wires in the main chamber of the rail{A}. Slide the endkap bracket into the rail{B}.
 [2.28] In case endstop/endlimit are preassembled in the top rail part, mark the position of the endlimit{A}. Remove the endlimit{B} and endstop{C}.

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





IT

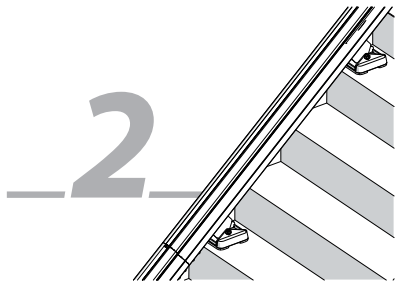RAIL FIXATION

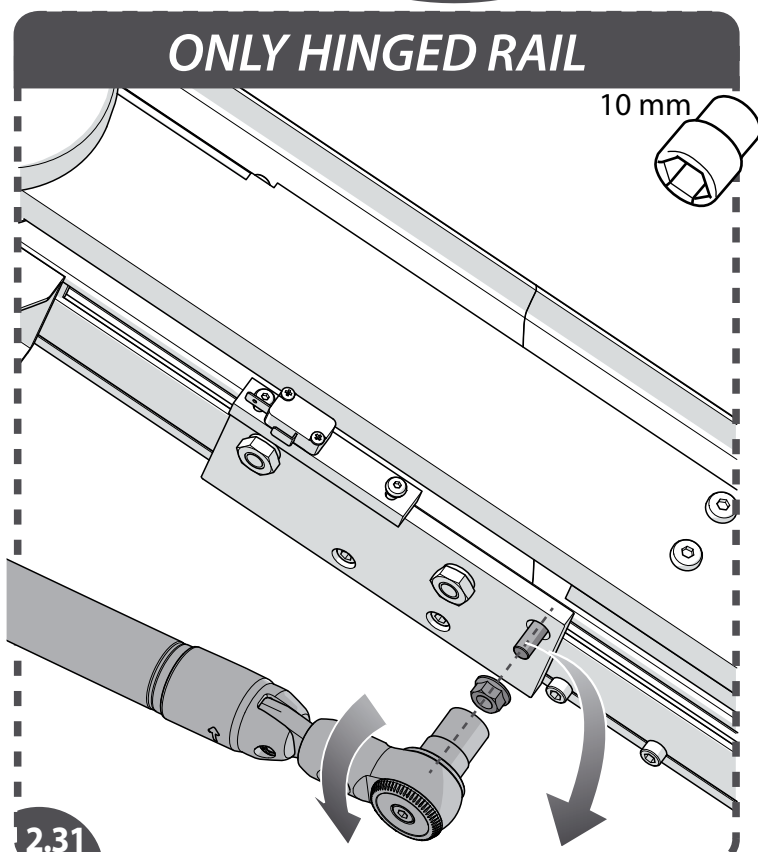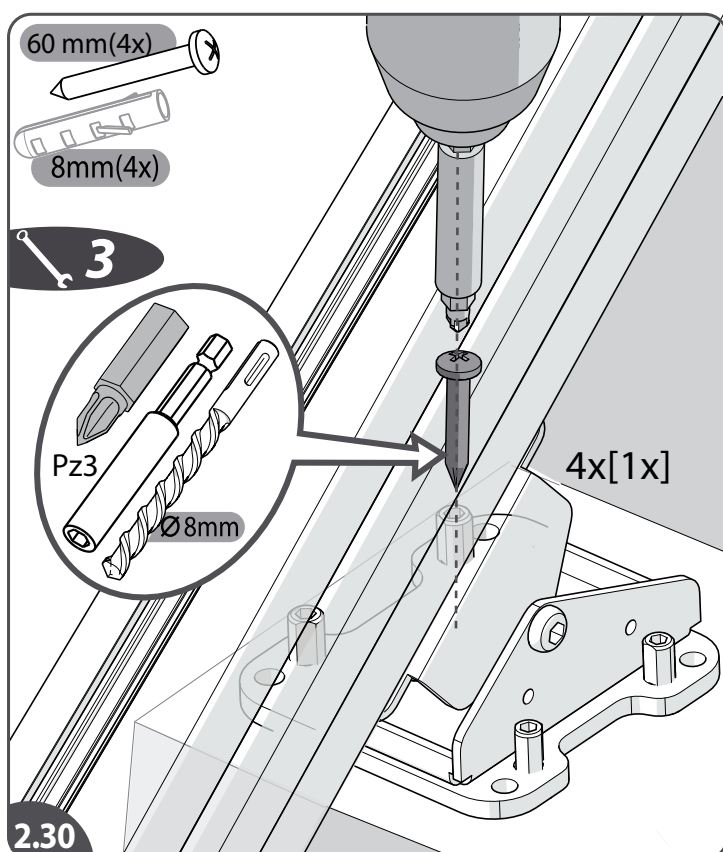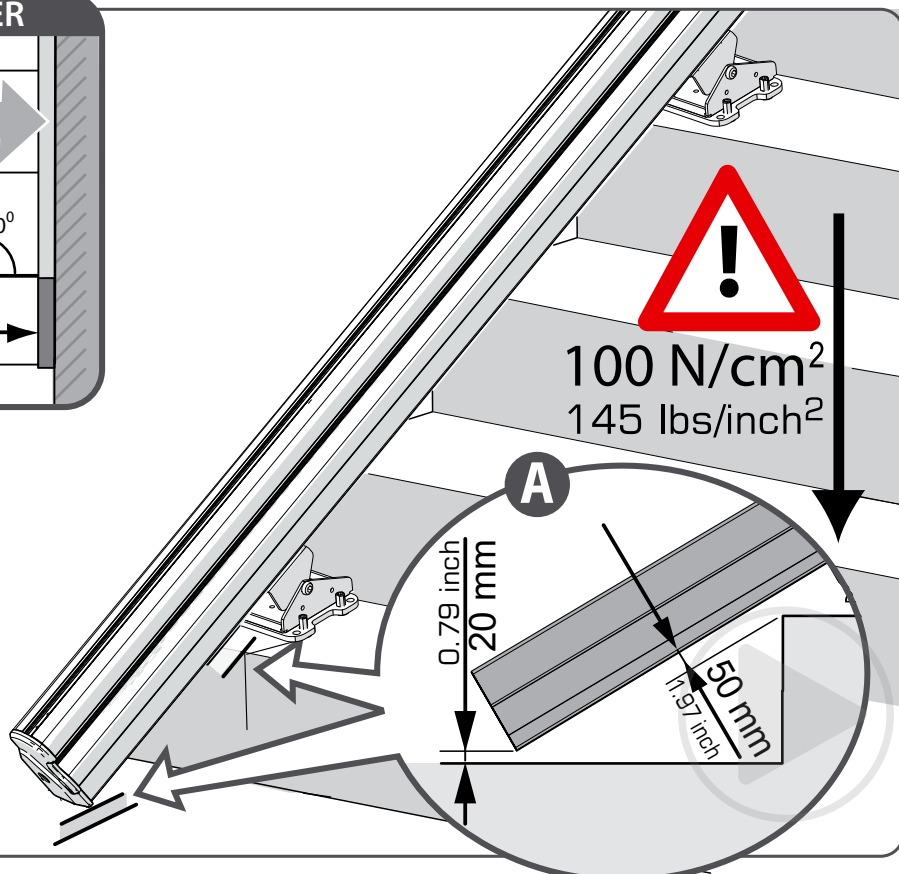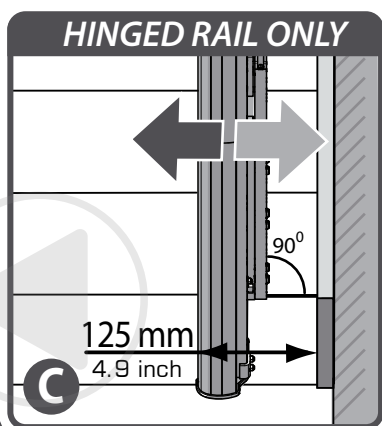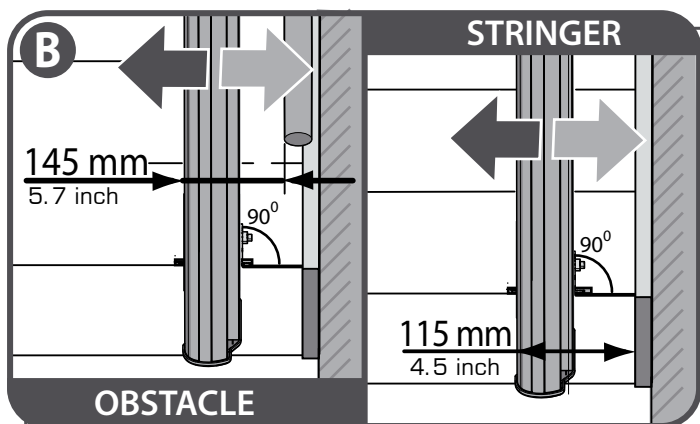[2.29] Position the track and achieve a minimum 50 mm (1.97 inch) clearance between the rail and the stairs{A}. Next, position the rail in relation to the wall/stringer/obstacle{B}. The minimum allowable clearance from the wall or obstacle is 145 mm(5.7 inch). From stringer it is 115 mm(4.5 inch). In case of a hinged rail a clearance of 125 mm (4.9 inch) is needed. Heavy duty units may also necessitate the track being positioned further into the staircase (up to an additional 50 mm).

[2.30] Secure the brackets in position using 1 screw (for wooden stairs) or rawl plug/screw (for concrete/stone stairs). Ensure that a minimum clearance of 50 mm (1.97 inch) is maintained between the rail and the stairs, and a minimum of 20 mm (0.79 inch) between the rail end and the floor.

EN [2.31] *Only hinged rail:* remove the transport bolt+nut underneath the hinge.

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




POWER PACK PREPARING

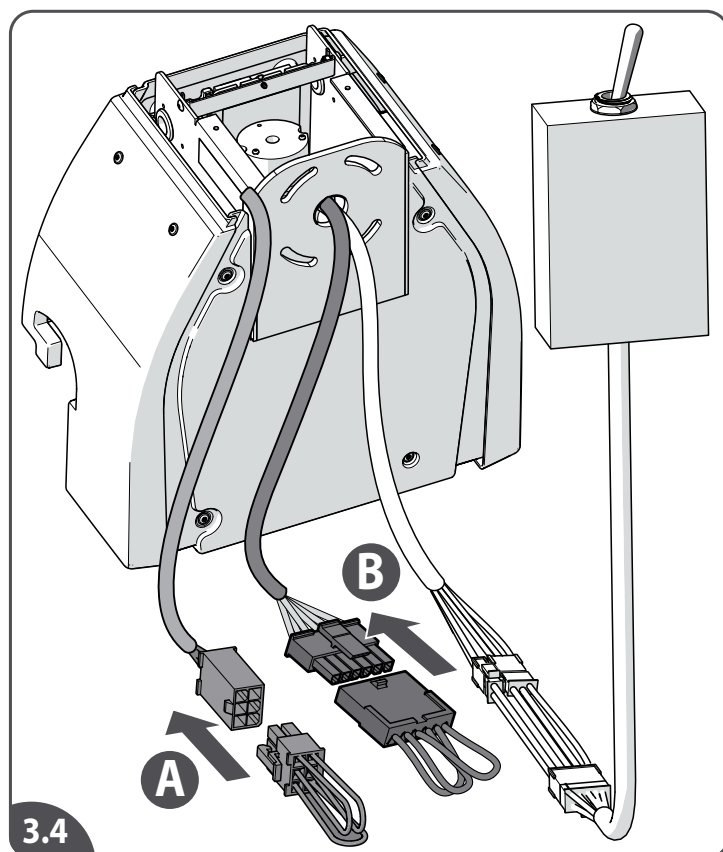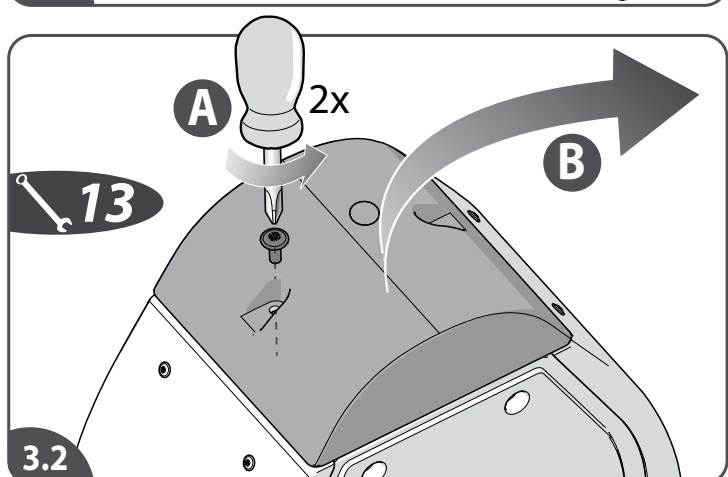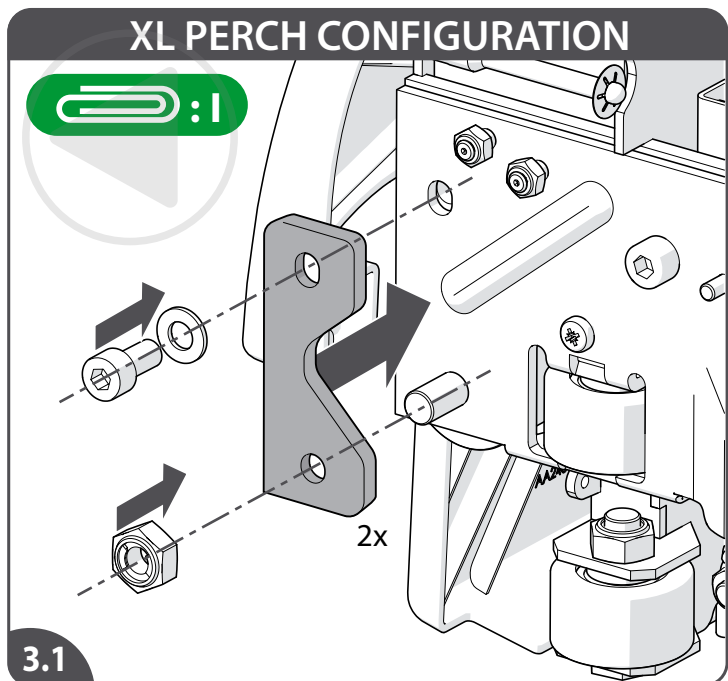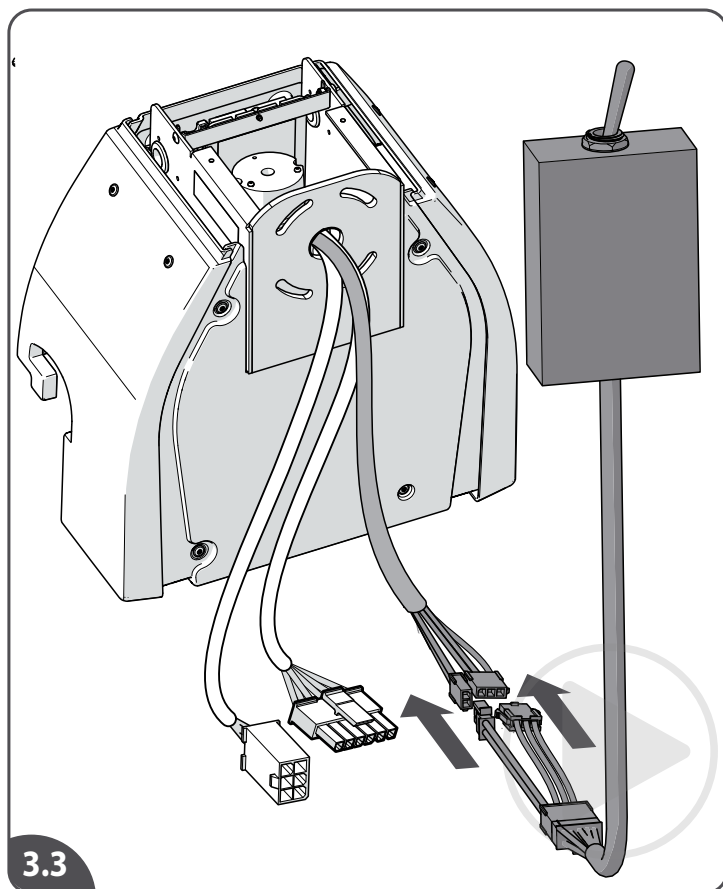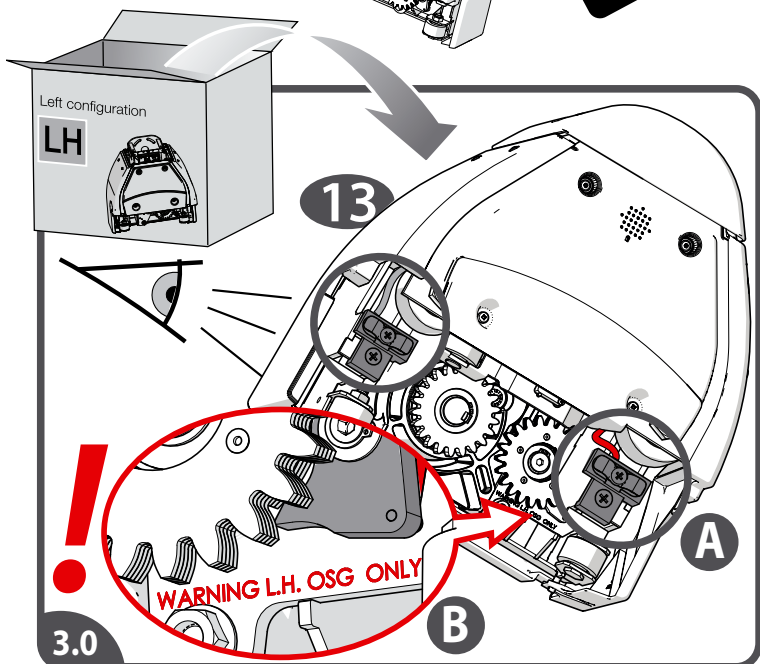[3.0] Remove the power pack from the box and put it on the floor at the top of the stairs. Visually inspect the contacts on the power pack to ensure that they are not damaged in any way. The powerpack can have two contacts: a chargecontact (the lower one) and a communication contact (the upper one) for the hinge option{A}. Check if the power pack has the right configuration (LH on a LH lift, RH on a RH lift). Otherwise the OSG is not functioning!{B}. Check all rollers are free running by hand.

[3.1] In case of a XL PERCH configuration, install extra plates on the power pack (see attachment I).

[3.2] Loosen the top cap of the power pack and remove it.

[3.3] Connect the loading toggle (SIM22900+181001.50083).

EN [3.4] Connect the shorting link (switches footrest){A}. Connect the shorting link (switches underpan){B}.

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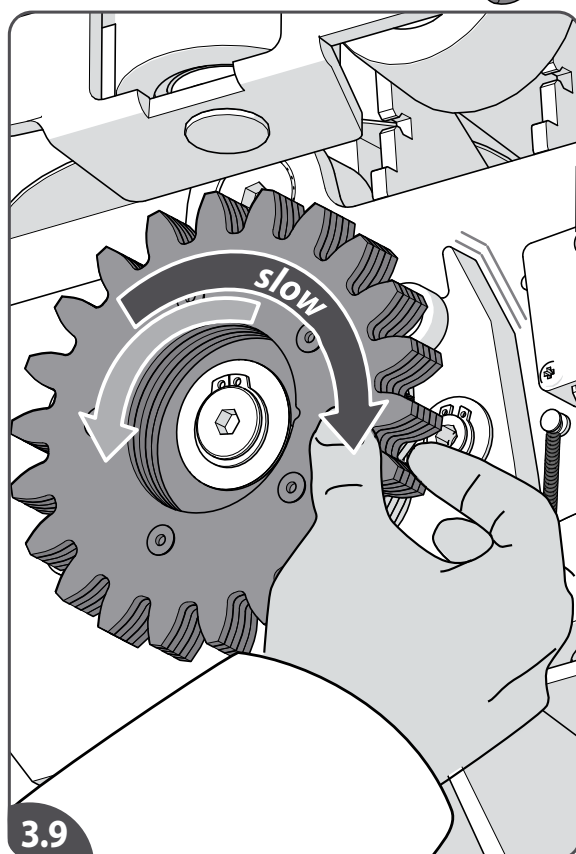
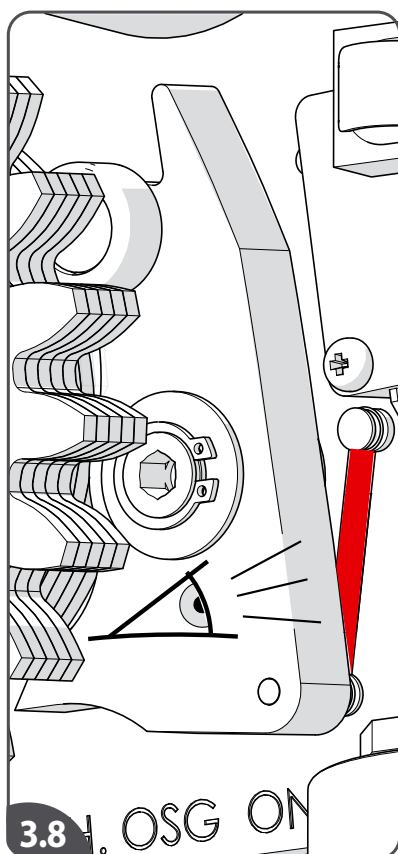
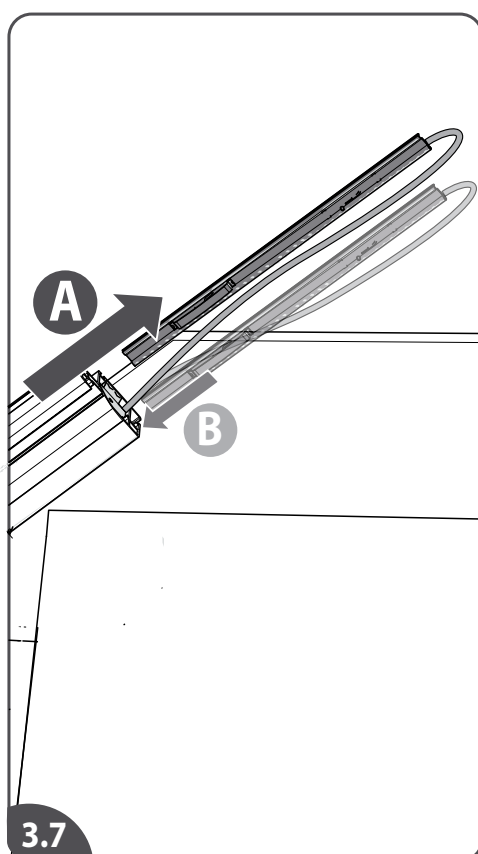
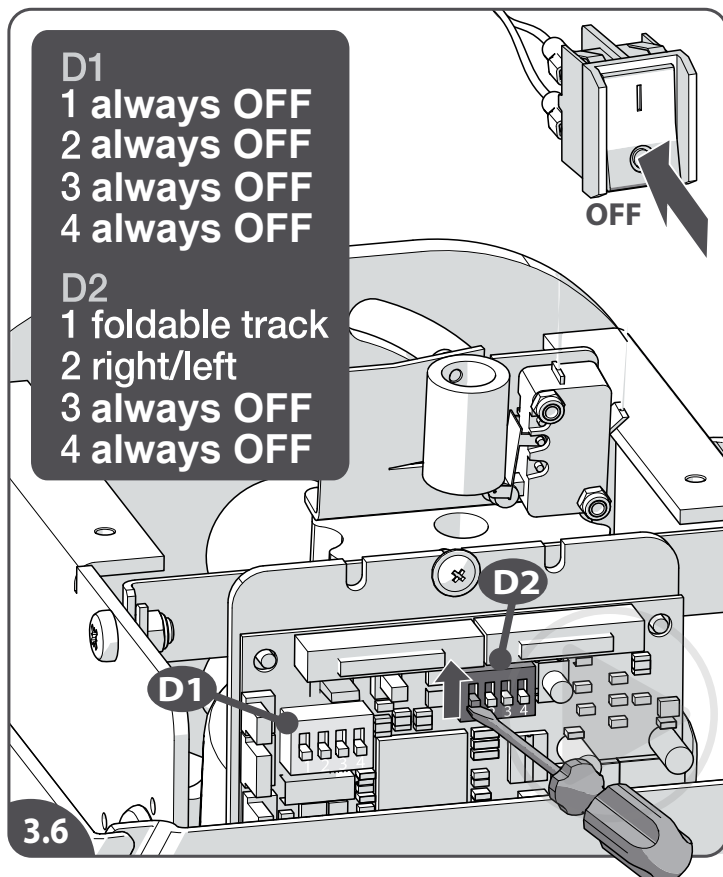
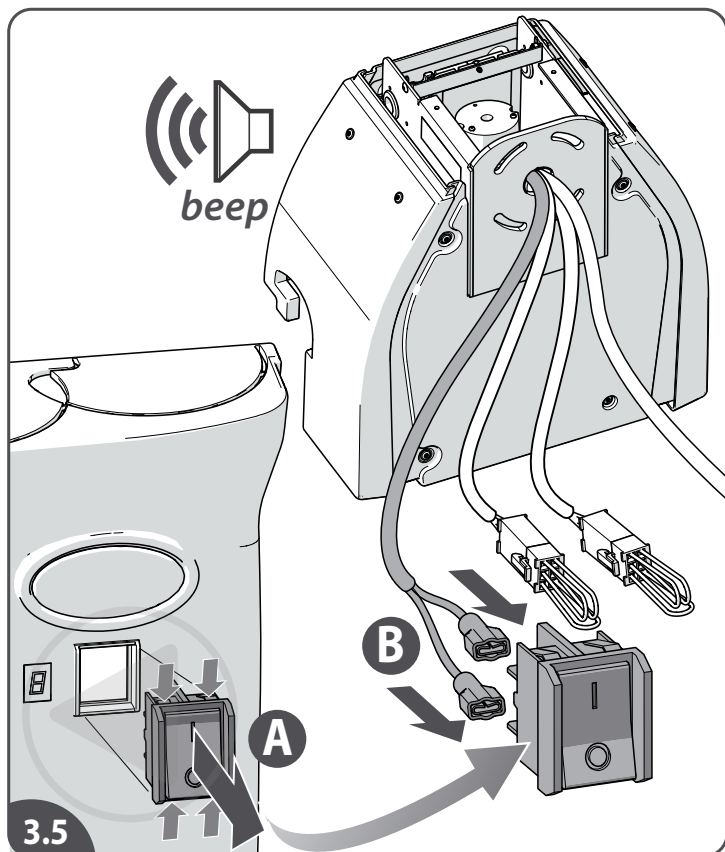
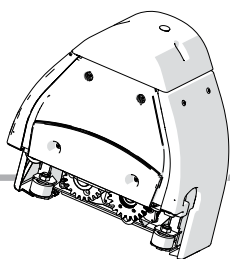
[3.5] Remove the on/off-switch from the front cover{A} and connect the switch to the battery looms{B}.
 [3.6] If needed, change the settings on the dip switches on the main PCB (power supply/power switch should be off before changing): e.g. in case of a hinged rail (D2-1) or in case of a right hand configuration (D2-2). Look at the inside of the front cover for the settings.

Note: Always remember after any changes to main PCB to switch off for 15 seconds and switch back on again.

[3.7] Remove the charge contact{A}. Put this contact temporarily in the bottom chamber of the rail{B}. The OSG should be inspected before loading to the rail and tested when the power pack is on the rail:
 [3.8] Visually inspect the oscillator spring is located over the anchor capstan (do not probe the spring; this is an instrument spring and can be damaged).

EN [3.9] Rotate the gear slowly in both directions with your forefinger on the OSG gear; the oscillator should move in a controlled manor pulled by the tension springs.

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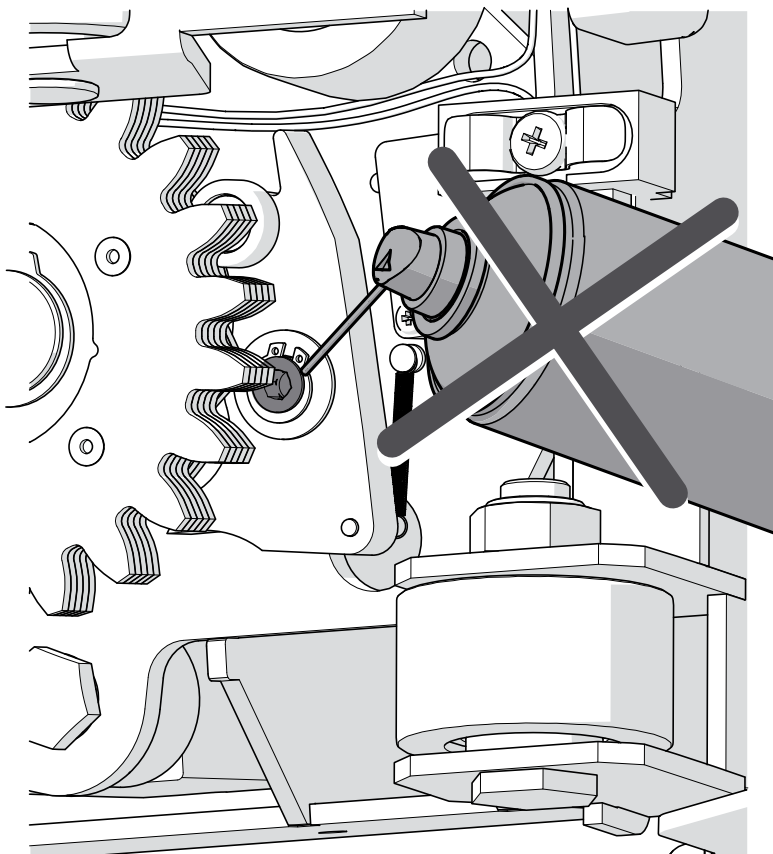
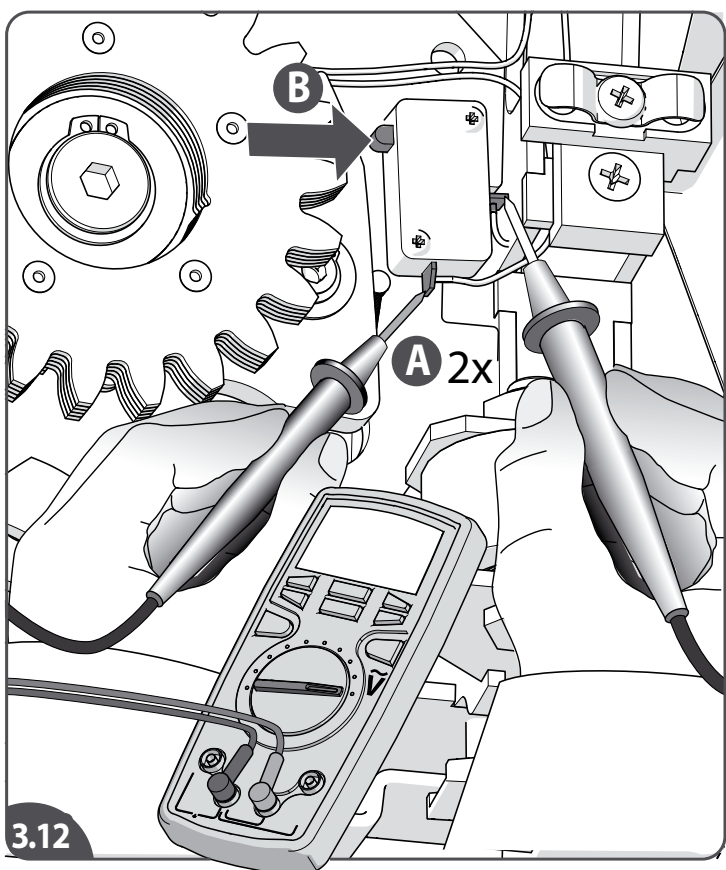
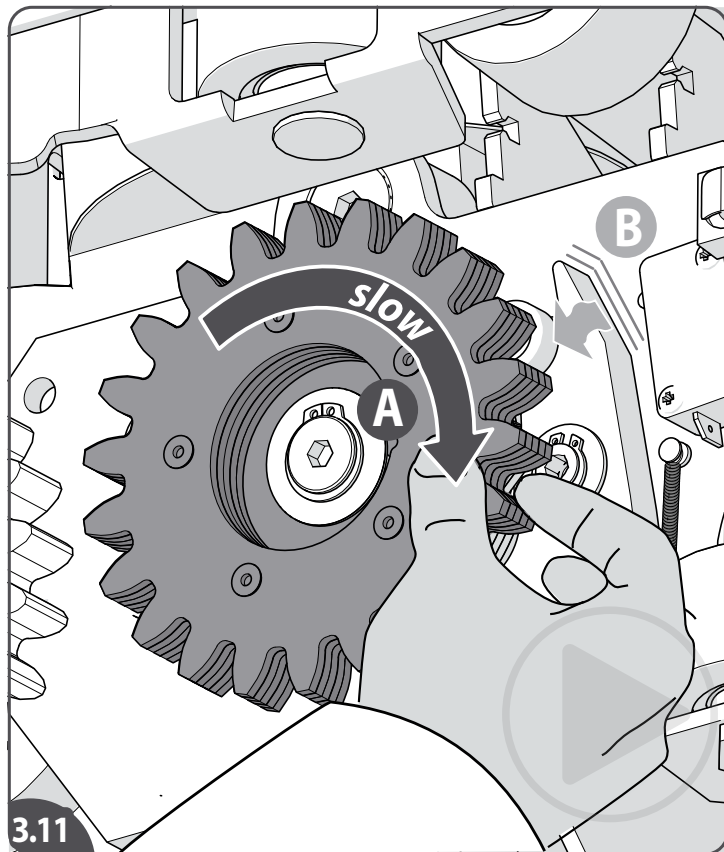
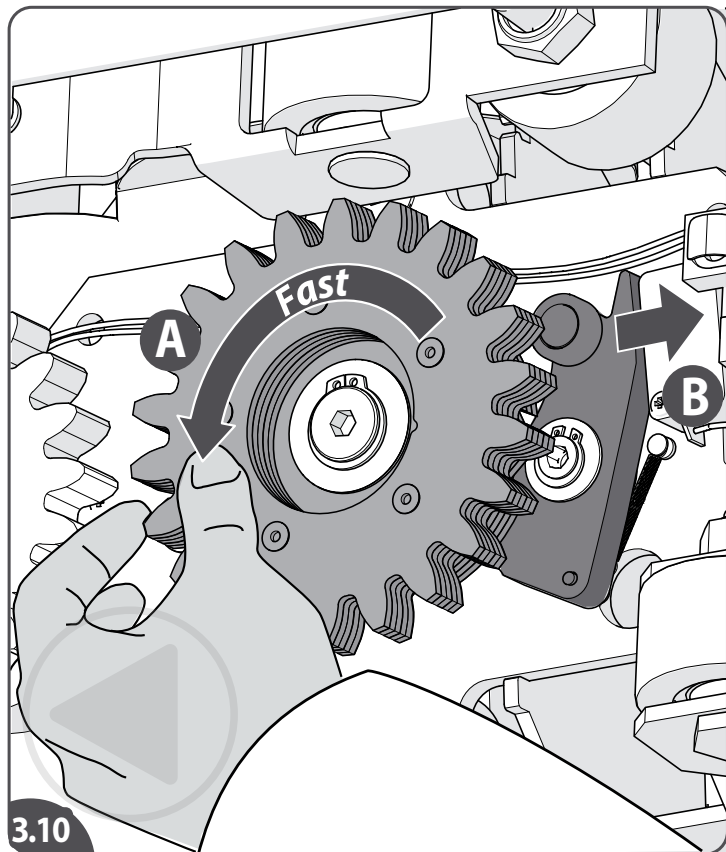
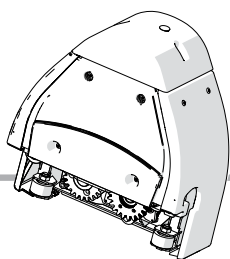
[3.10] Accelerate the rotation of the gear in a counterclockwise direction with your forefinger{A}. The oscillator will engage and lock the OSG gear as shown{B}.

[3.11] Rotate slowly in the opposite direction to free the oscillator.

[3.12] Test with a multi meter both microswitch terminals for continuity{A}. Depress the microswitch while probing continuity and listen for interruption of continuity{B}.

Note: no lubrication should be used on the oscillator pivots: they are dry life bearings. Lubrication will compromise free rotation and cause false tripping.

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



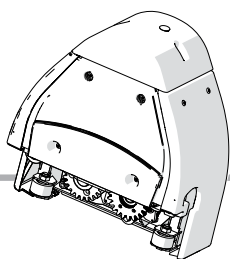
POWER PACK

DISCONNECT SHORTING LINKS

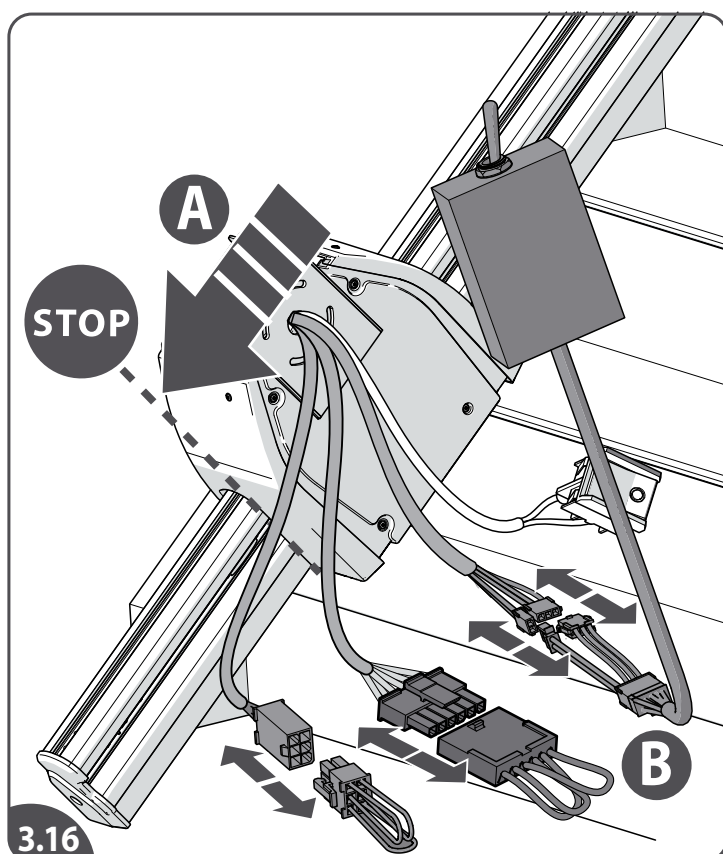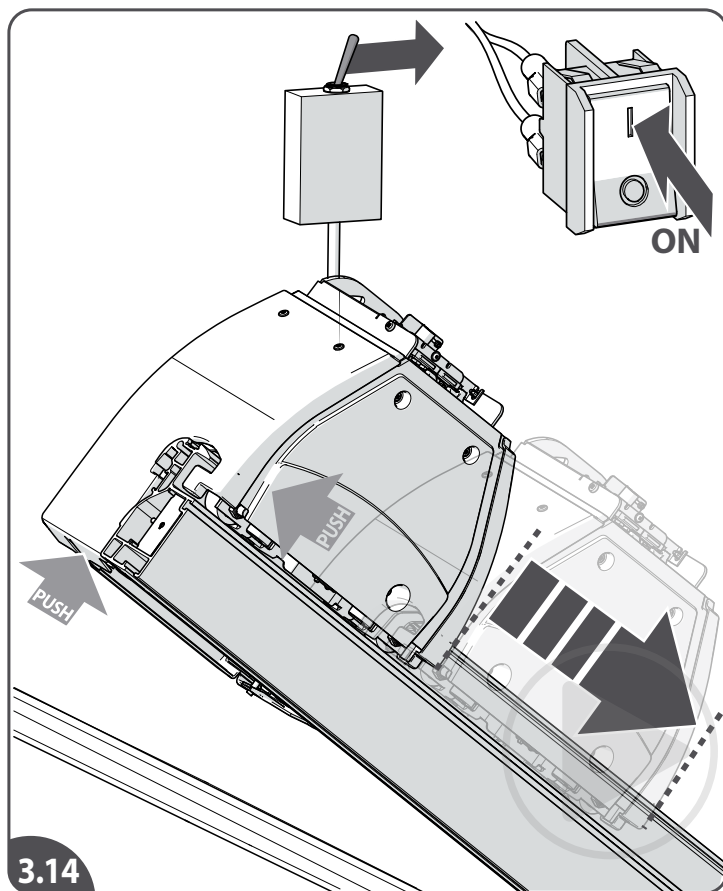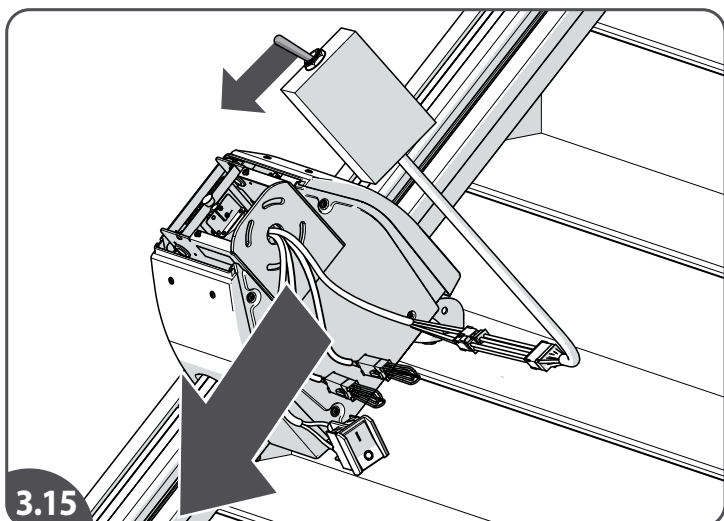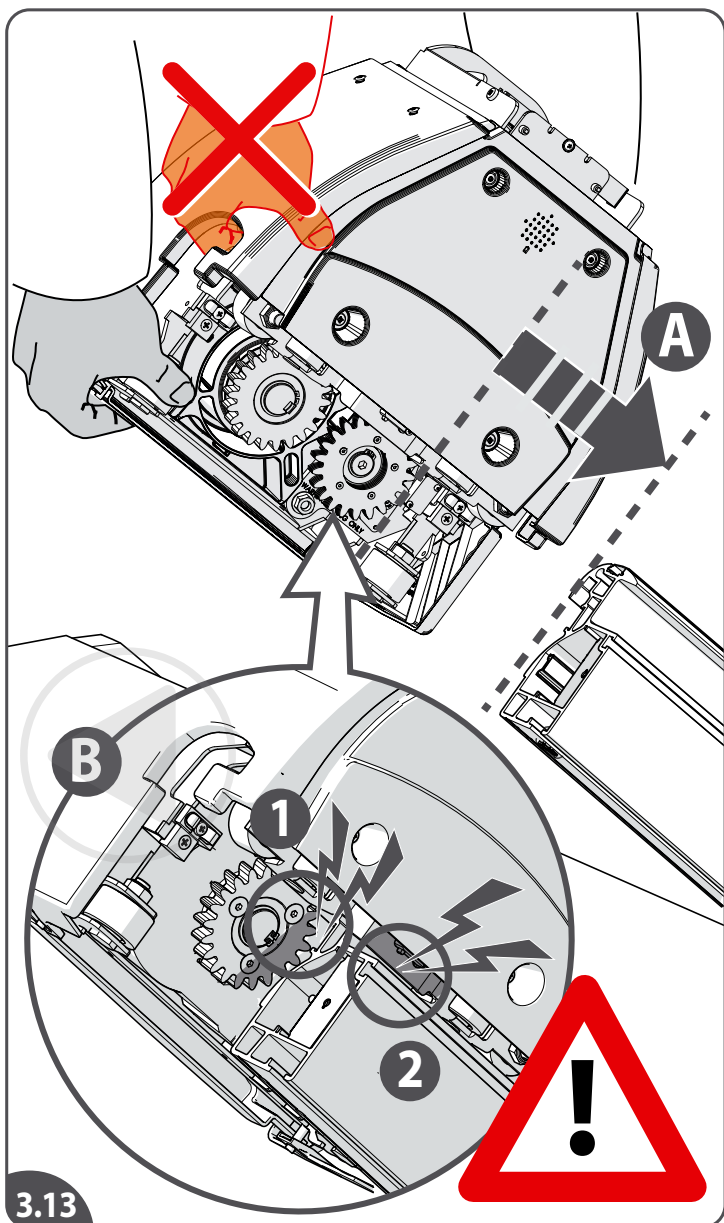- [3.13] Carefully place the power pack on the rail{A}. [B] Be careful that the gearwheel is not hitting the rail{1} and the frame of the power pack does not scrape the surface of the rail{2}.
- [3.14] To prevent this, push the unit up and to the front. Drive the power pack into the rail by operating the joystick.
- [3.15] Drive the power pack down until just above the charge contact.
- [3.16] Disconnect the loading toggle and the shorting links{B}.

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



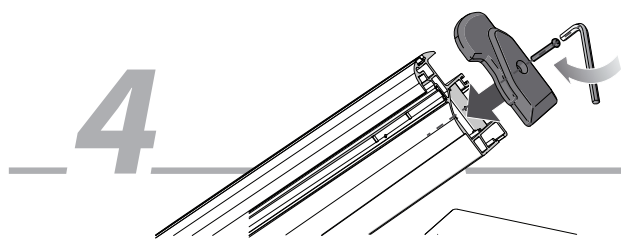
FINISH RAIL

CHARGER UPSTAIRS

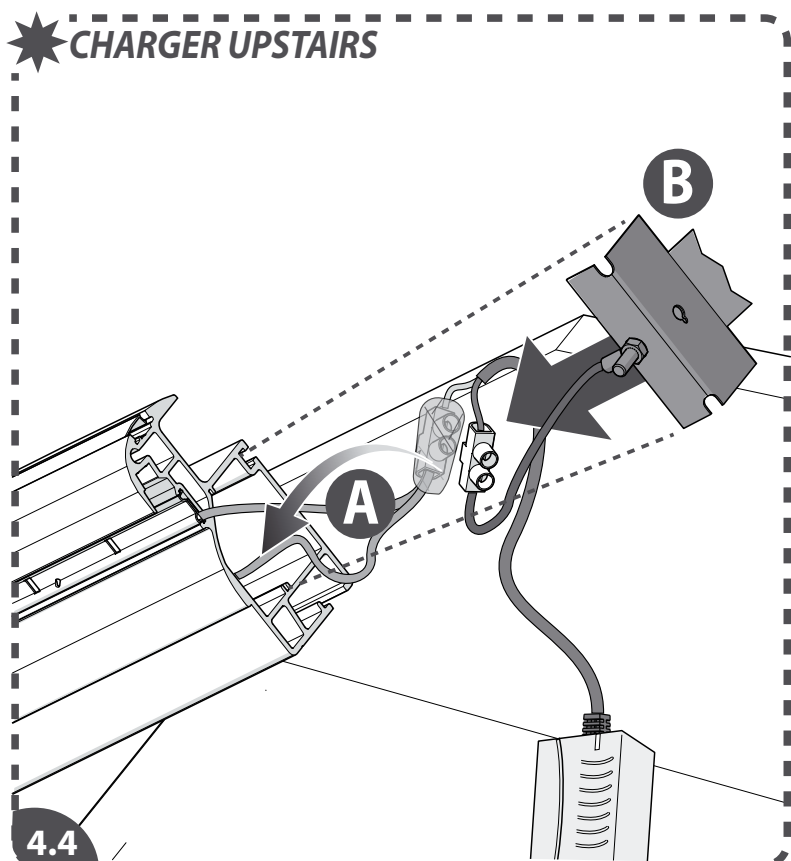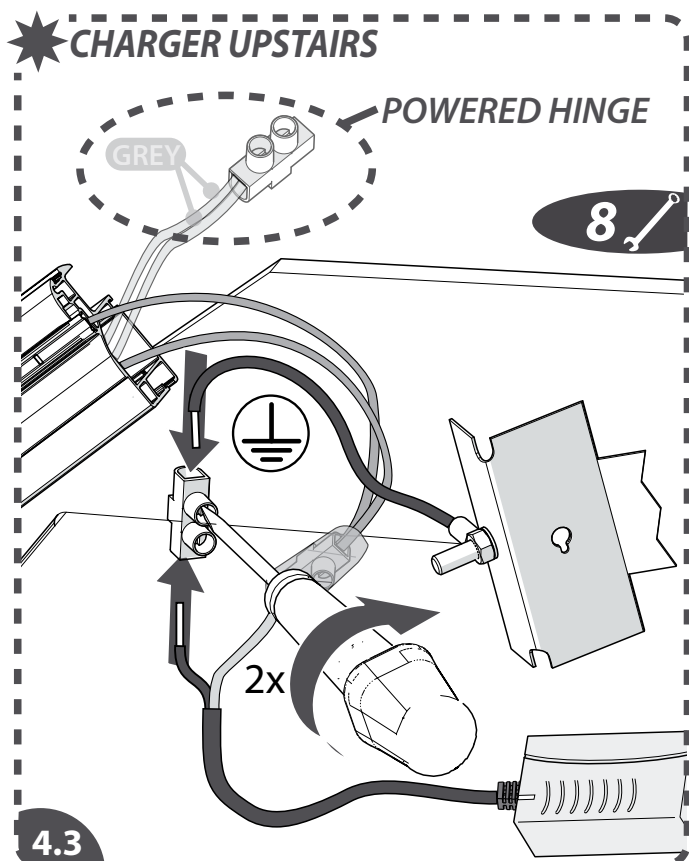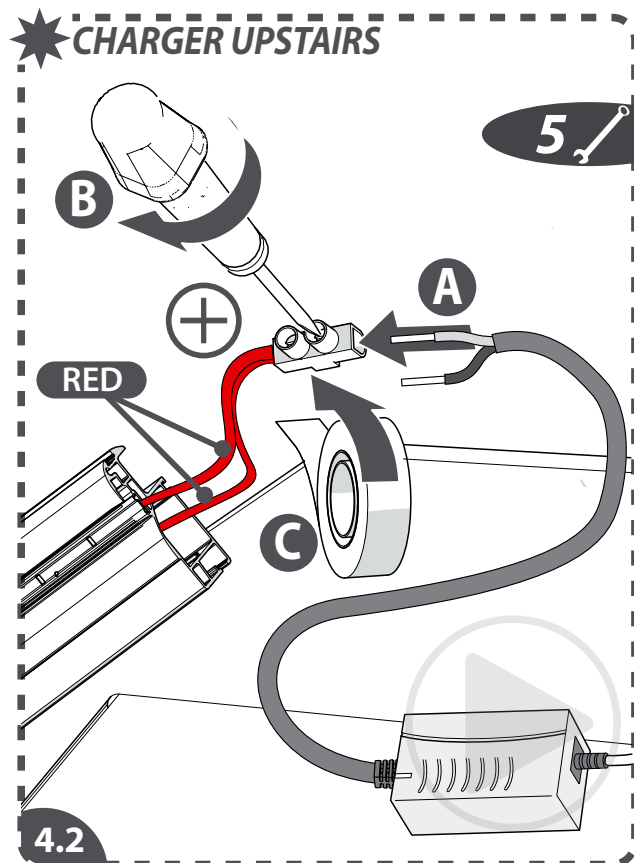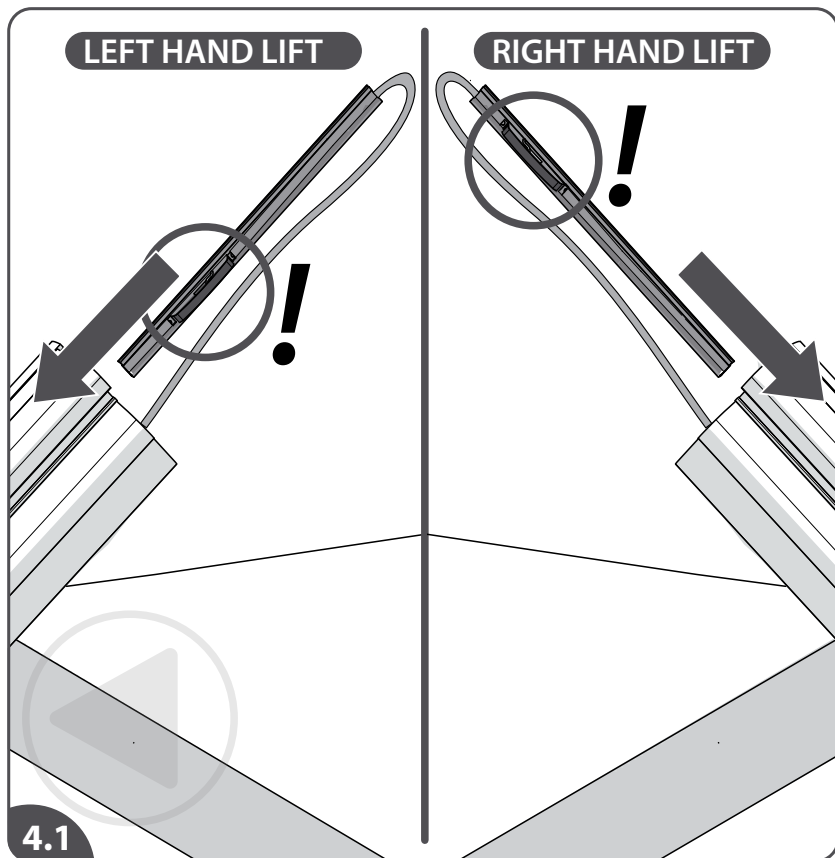- [4.1] Refit the charge strip (contact on down side (left configuration) !).
- [4.2] In case the charger is mounted upstairs: fasten the +wire (white) of the charger in the connector block{A+B}. Insulate the connector block{C}.
- [4.3] Fasten the -wire (black) of the charger and the earthwire in another connector block. In case of a powered hinge: be carefull do not connect the -wire (black) of the charger to the grey communication wires !
- [4.4] Put the wires in the main chamber of the rail{A}. Slide the earth bracket into the rail{B}.

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




4

FINISH RAIL ENDSTOP

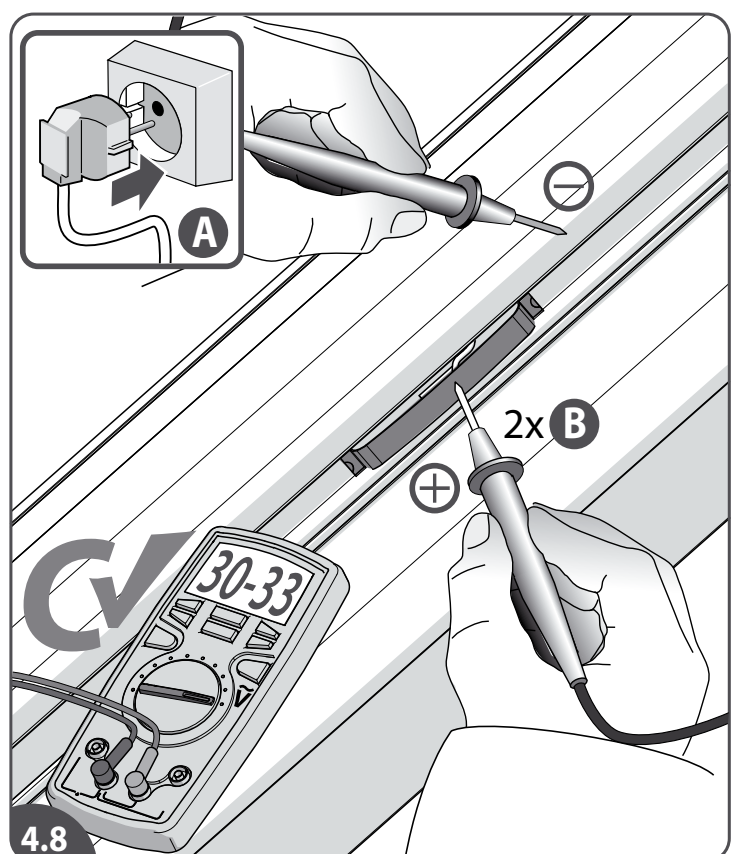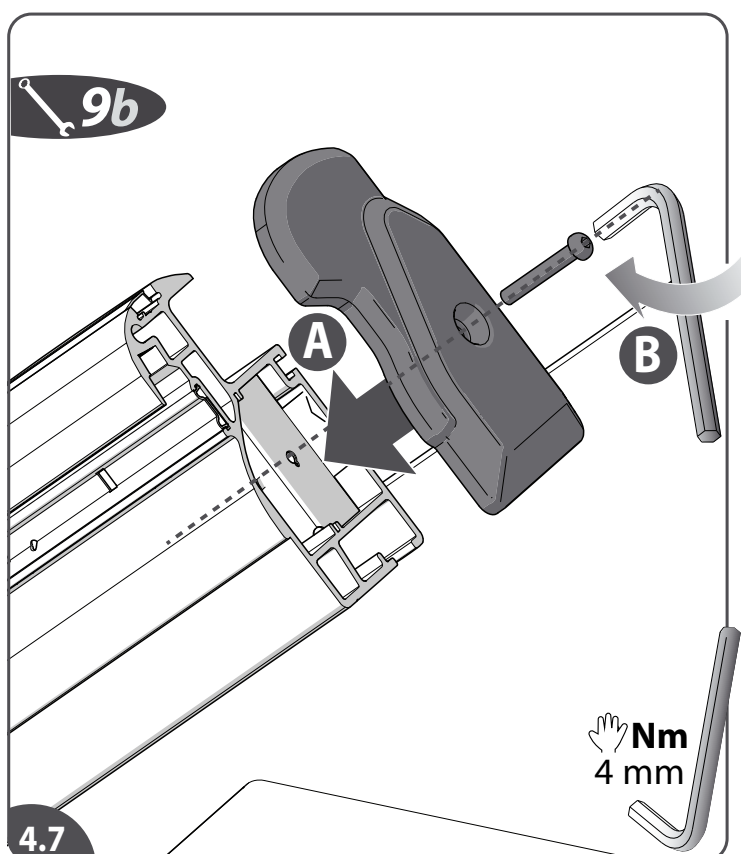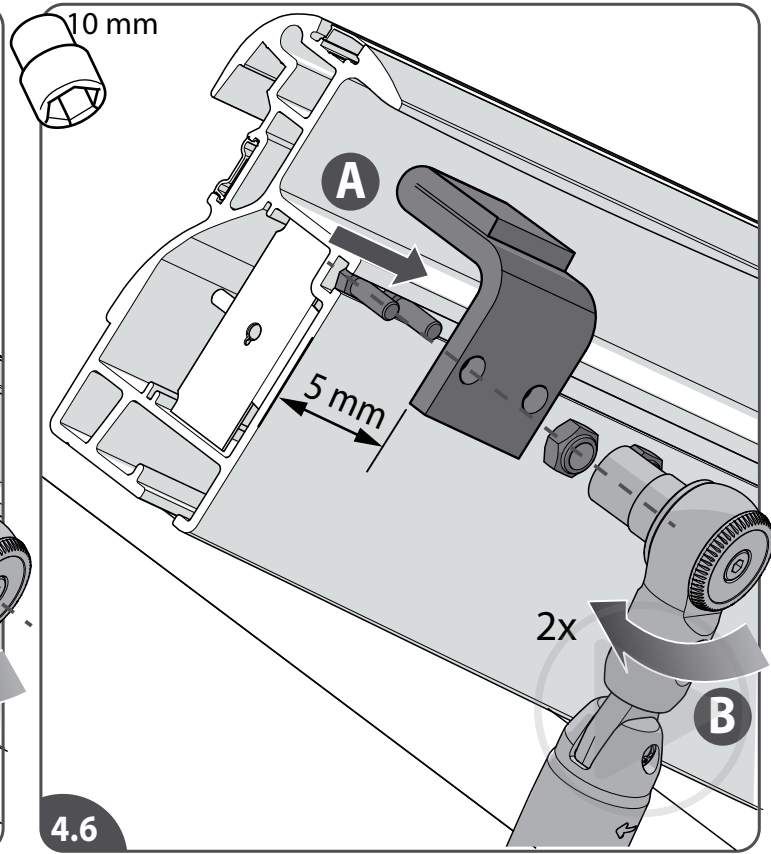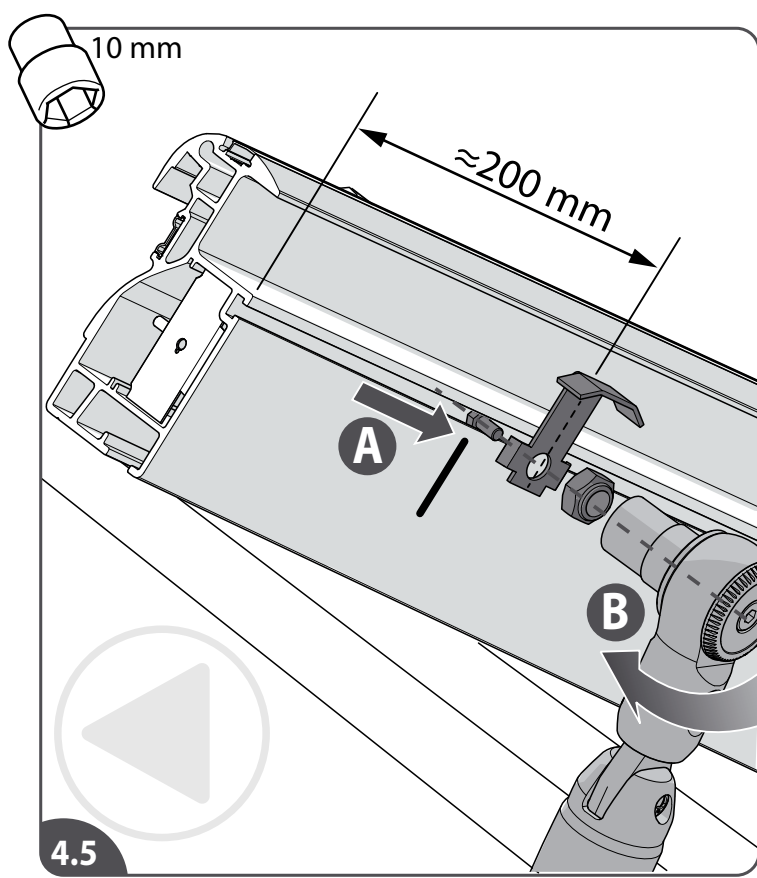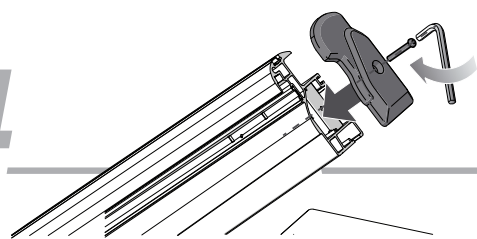- [4.5] Slide the bolt of the limit actuator in position into the rail (approx 200 mm from the end of the rail or on the marked position) {A}. Install the limit actuator and tighten the bolt securely {B}.
- [4.6] Slide the bolts of the endstop into the rail {A}. Install the endstop and tighten the bolts securely {B}.
- [4.7] Install {A} and carefully fasten (handtight) the end cap {B}.
- [4.8] Put the plug in the socket {A}. Check the voltage of the charge contacts (2x). It should be between 30 V - 33V {B}.

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
IT


5

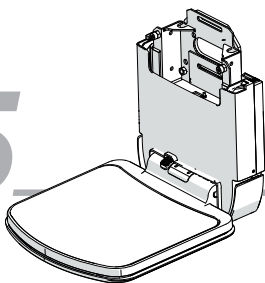
YOU MUST BE LOGGED IN TO
THE HANDICARE PARTNERSITE!

CHASSIS


PLACING CHASSIS

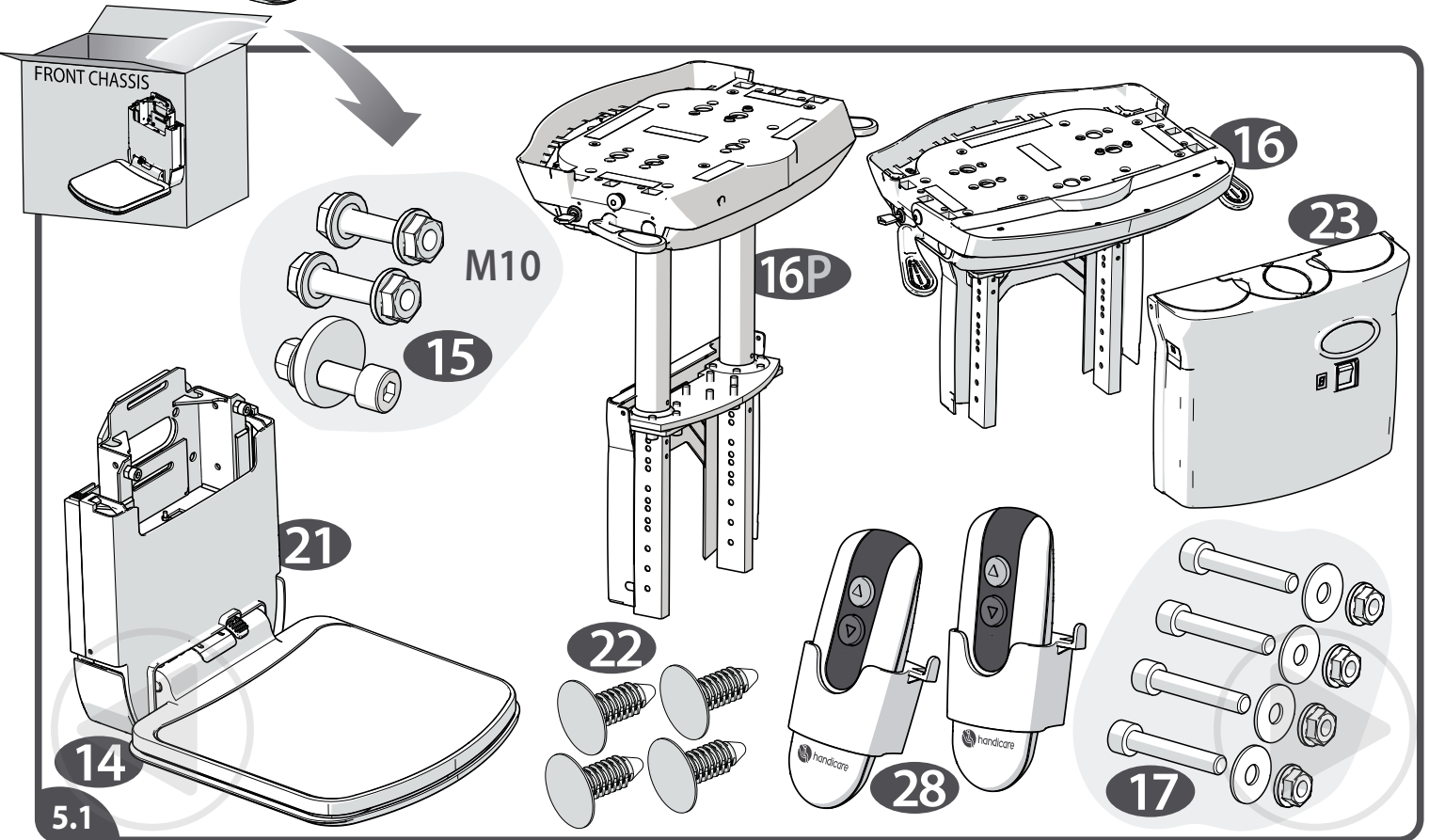5.1

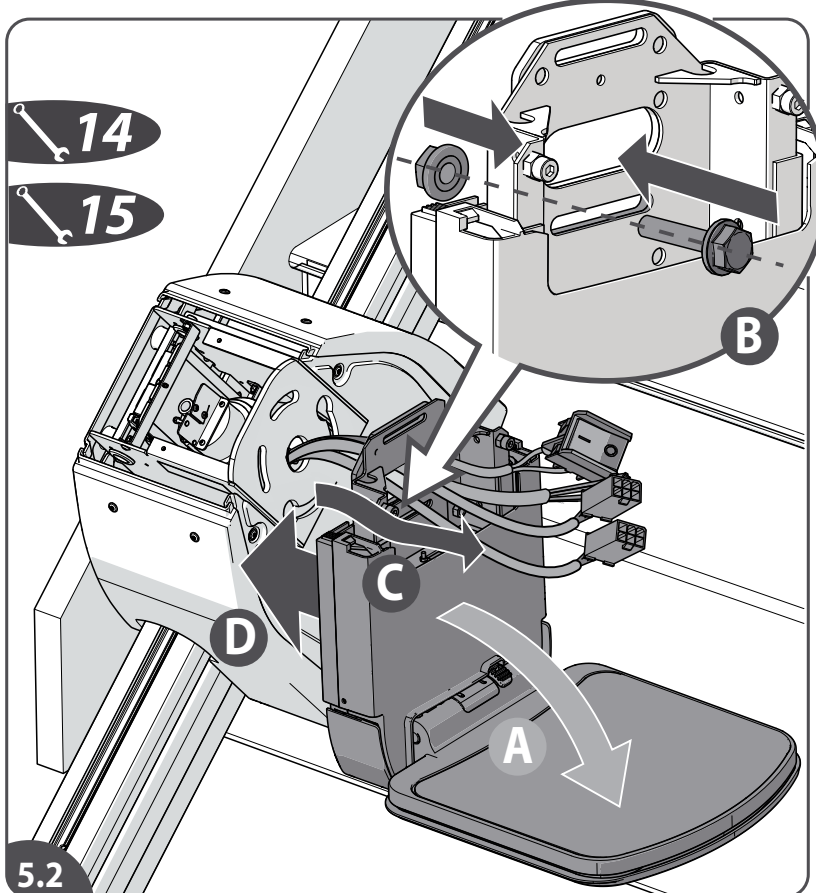5.2

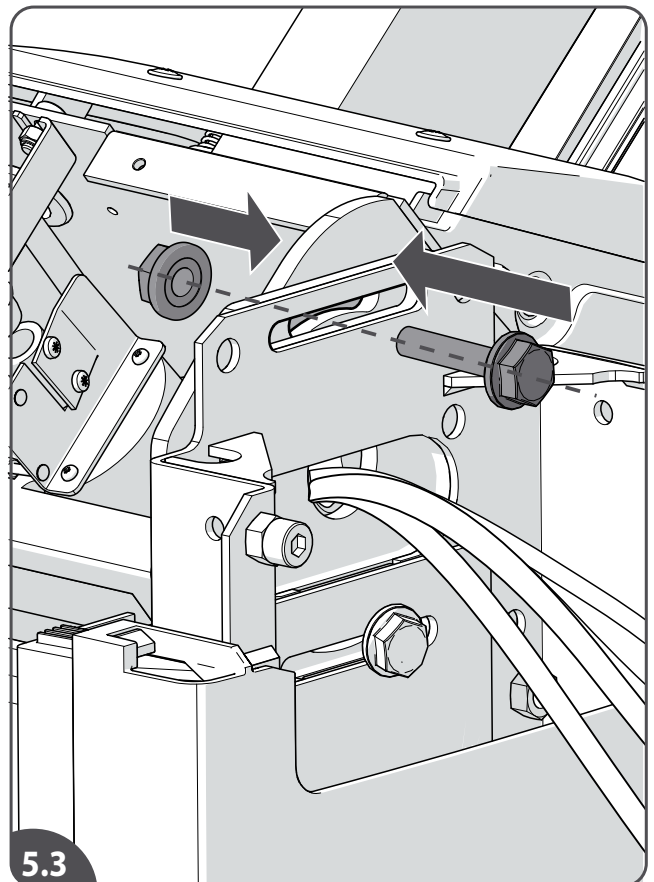5.3

- [5.1] Remove the front chassis components from the box. Remove the seat base components from the box.
 [5.2] Carefully pull down the footplate{A}. Install the first bolt/nut on the chassis. Do not tighten this bolt yet; only a few turns{B}. Guide the power pack looms through the hole of the chassis{C}. Position the chassis against the power pack{D}.
 [5.3] Install the second chassis bolt into the frame. Do not tighten this bolt yet; only a few turns.

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




5

CHASSIS PLACING SEAT BASE

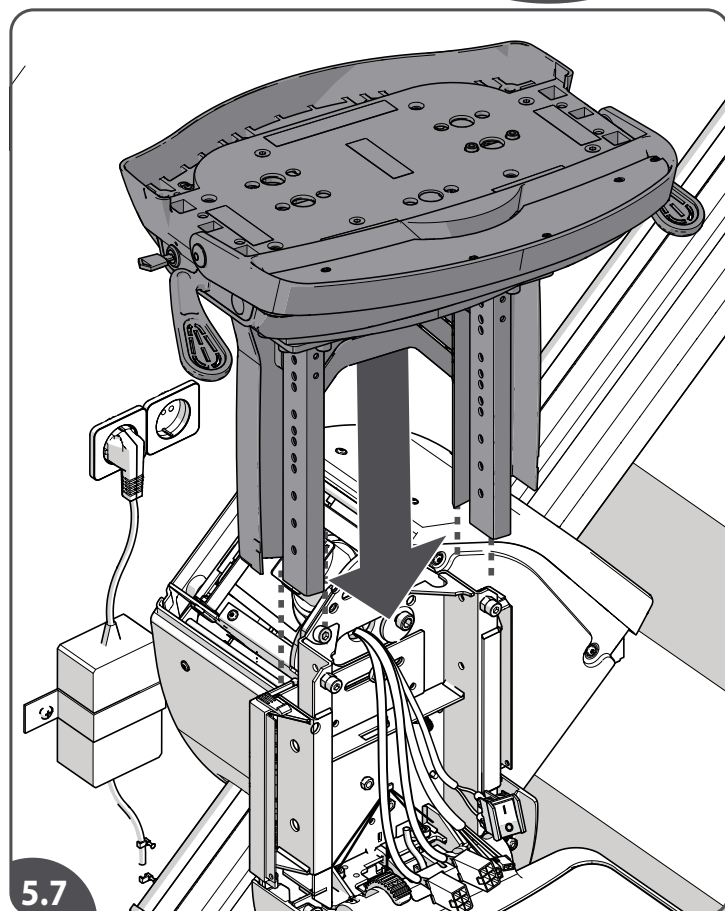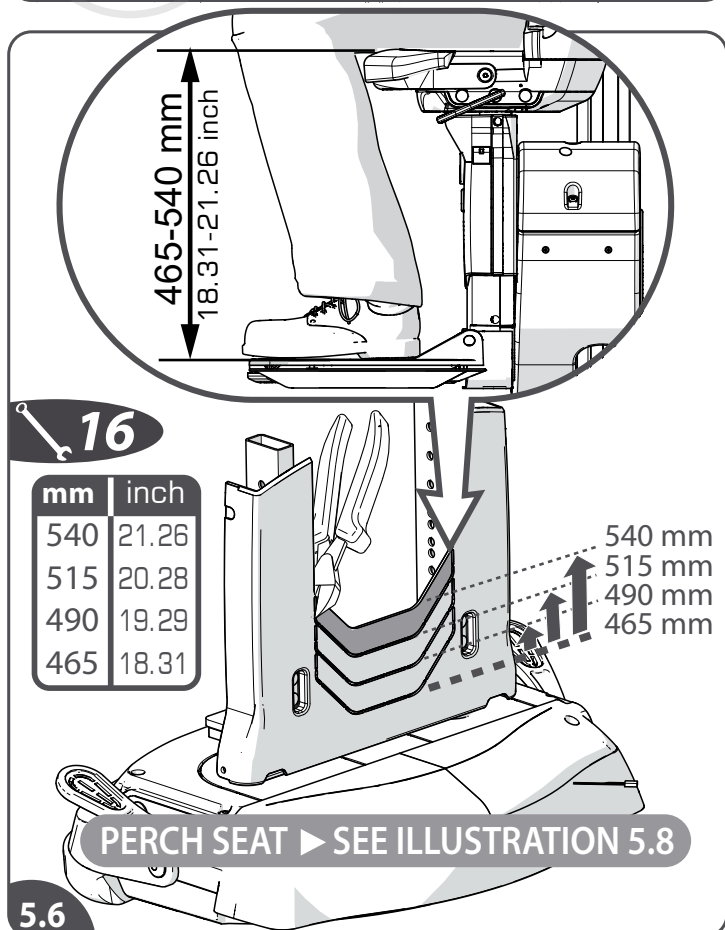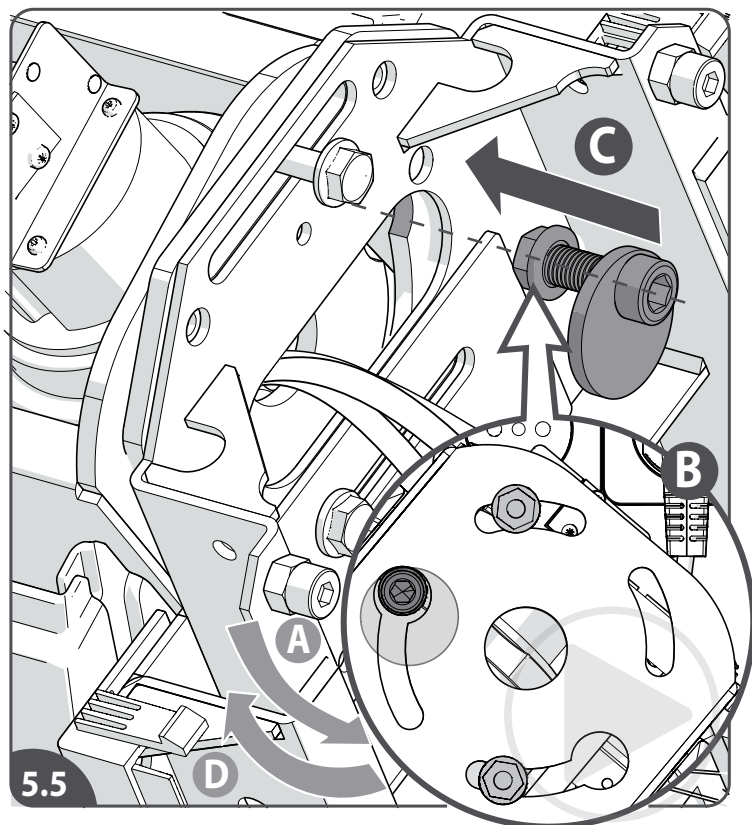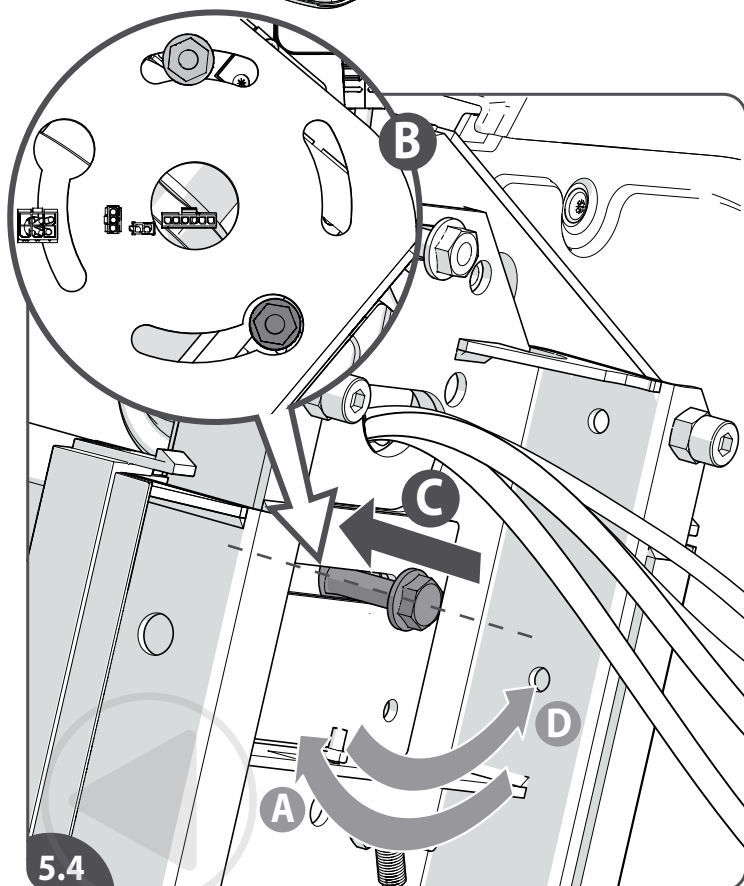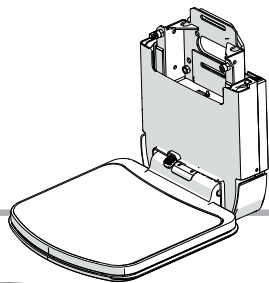[5.4] Turn the chassis{A} until the first bolt is exactly opposite the round hole in the mounting plate{B}. Hook the bolt into the power pack{C}.

[5.5] Turn the chassis{A} until the third chassis bolt, together with the fixing washer, is exactly opposite the round hole in the mounting plate{B}. Hook the bolt into the power pack; this bolt is always placed on top{C}.

[5.6] Depending on the main user height remove one or more cutouts.

[5.7] Slide the seatbase into the frame of the chassis.

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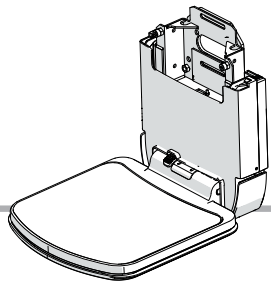
17

5.10

PERCH seat: [5.8] Depending on the main user height remove one or more cutouts.
 [5.9] Slide the seatbase into the frame of the chassis.
 [5.10] Keep the seat base in the correct height position and install the four frame bolts.

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




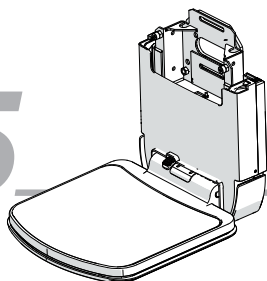
FIXATION CHASSIS

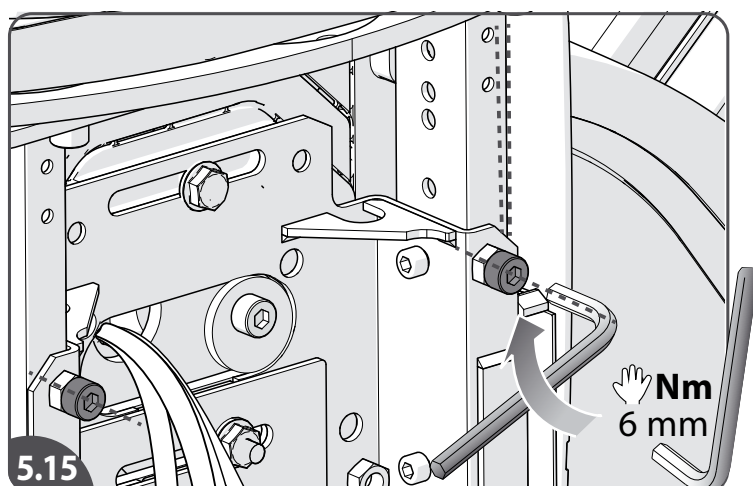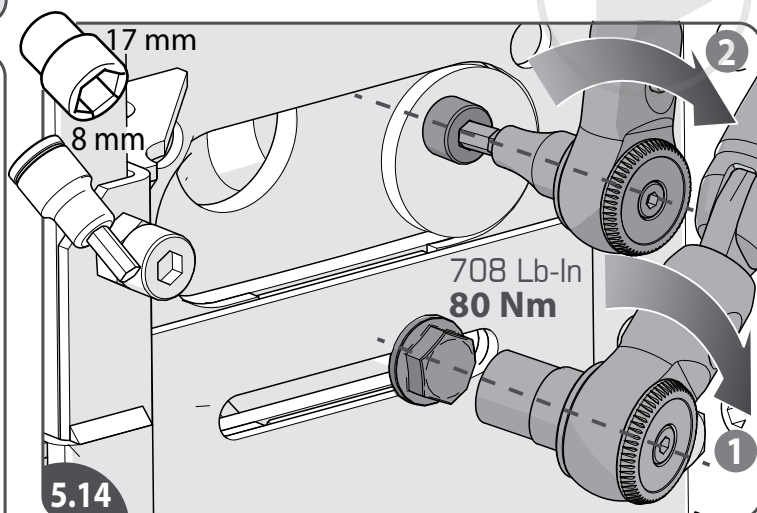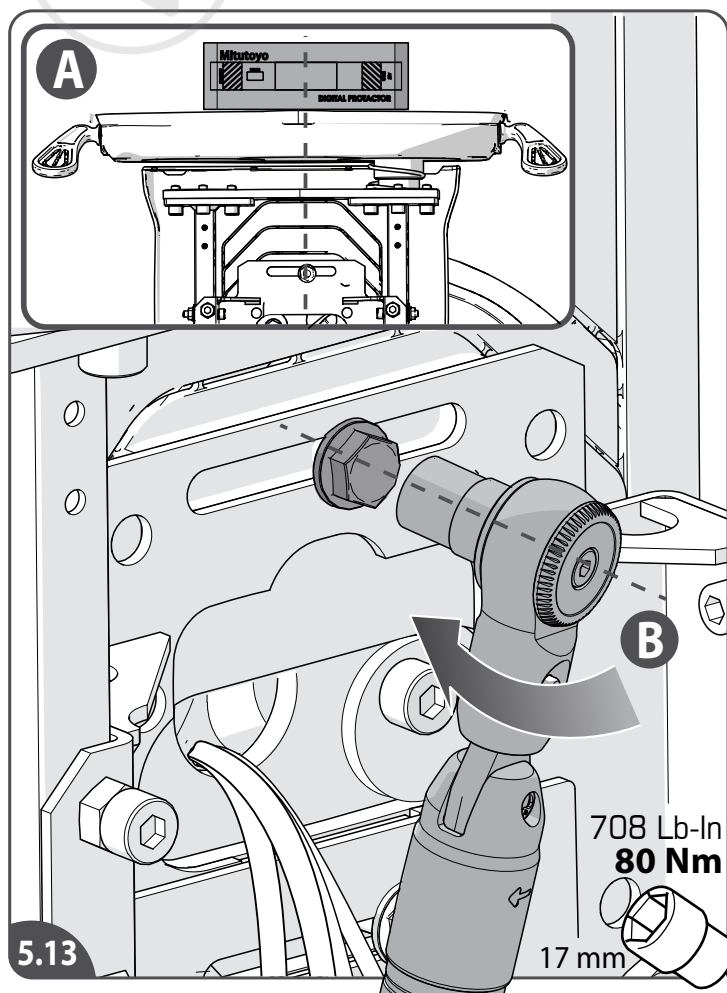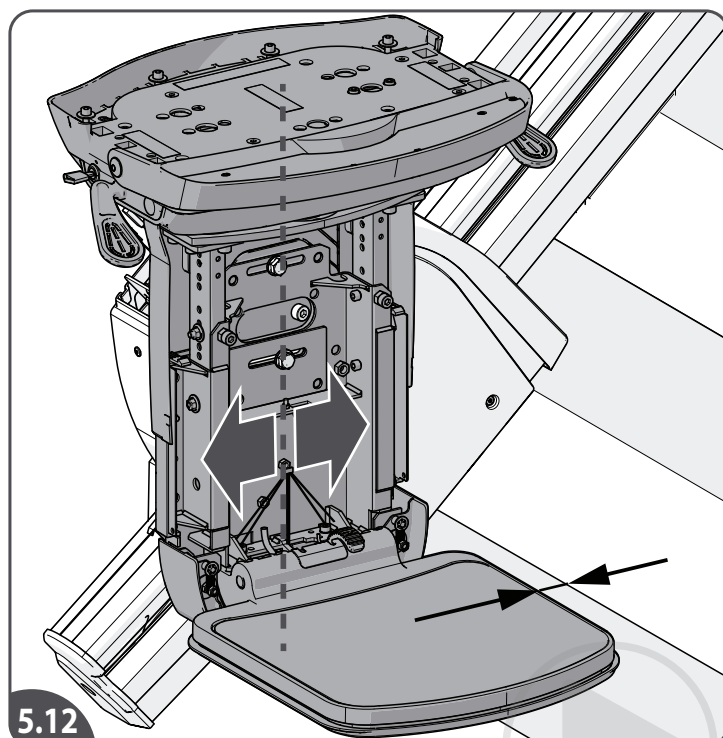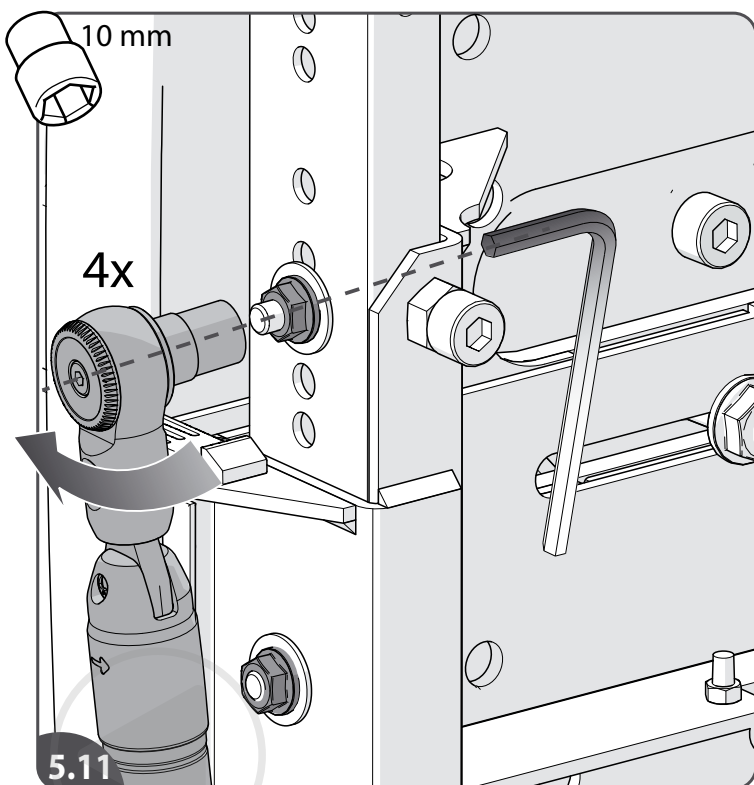[5.11] Tighten the frame bolts securely.

[5.12] Check if the distance between the footrest and stairs is correct. If needed slightly change the position of the chassis in the lateral direction.

[5.13] Use a spirit level to level the lift {A} and fasten the upper chassis bolt securely (80 Nm/708 Lb-In) {B}.

[5.14] Fasten the other chassis bolts securely (80 Nm/708 Lb-In).

[5.15] Tighten the two allen bolts in front of the frame (hand tight).

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



SEAT FIXATION SEAT BACK

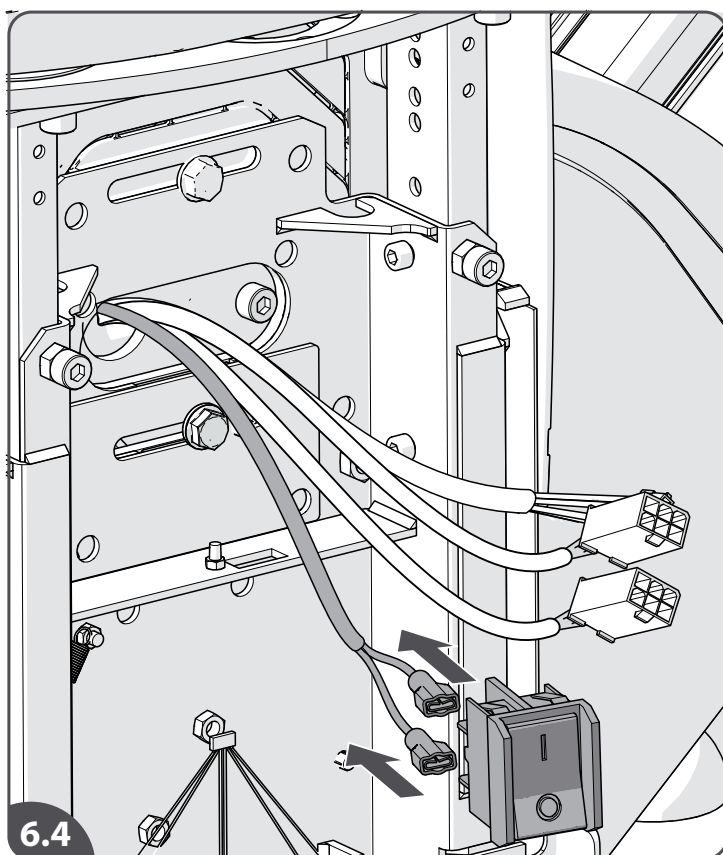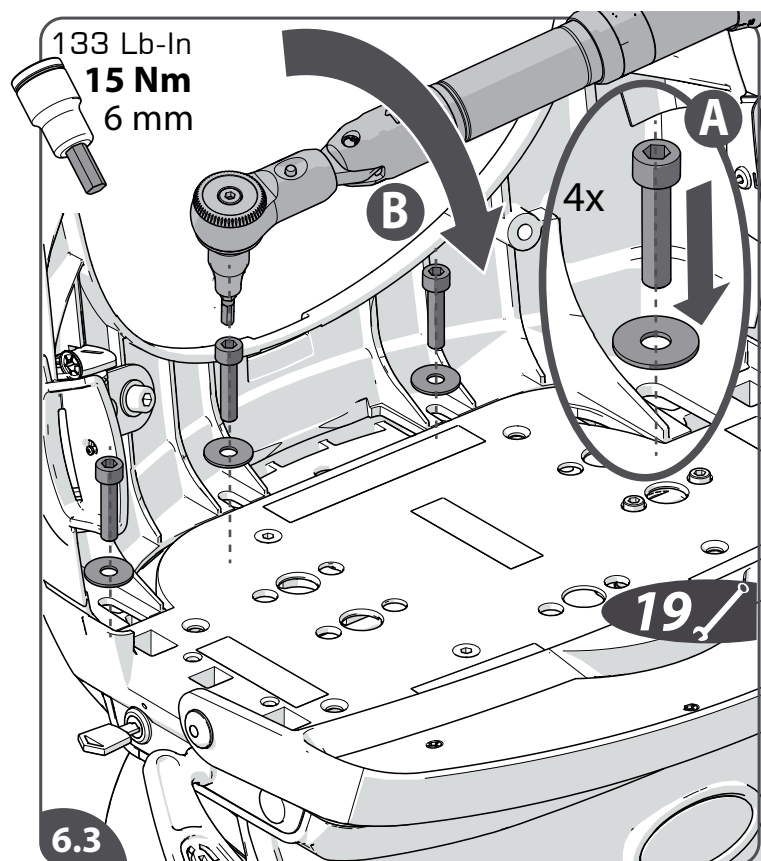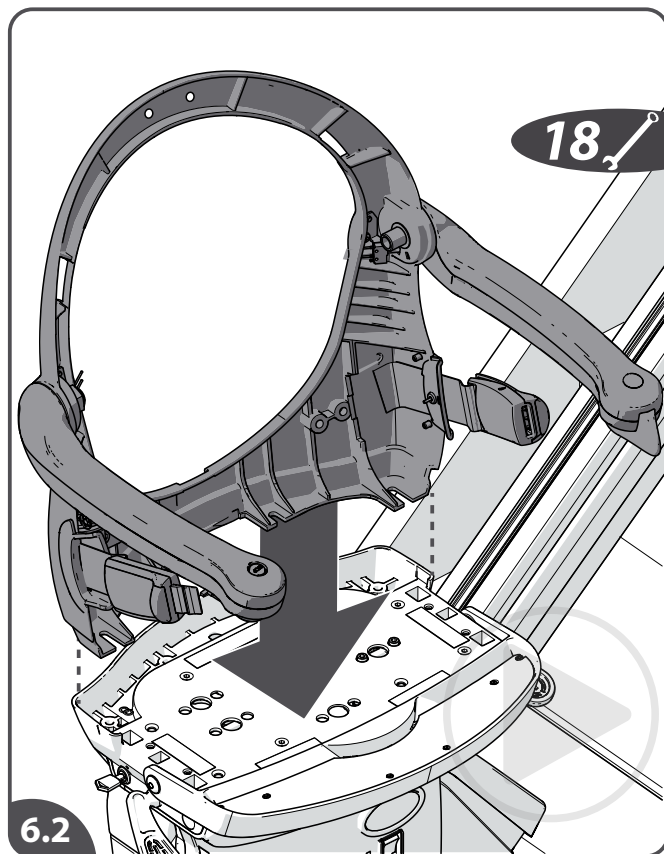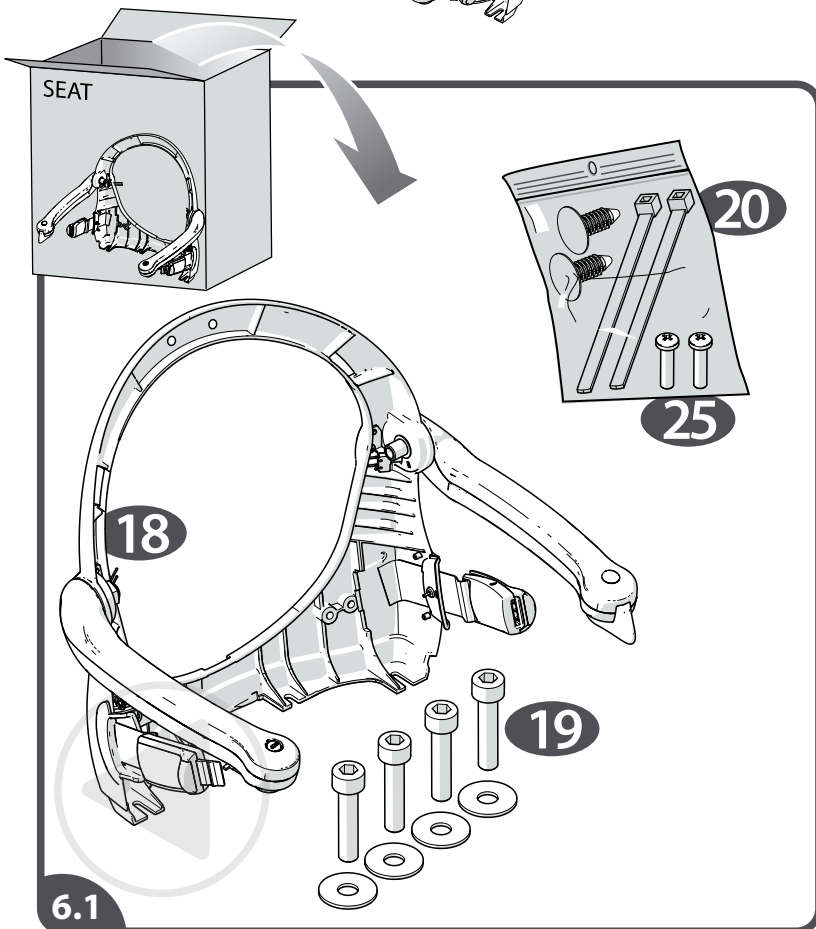- [6.1] Remove the seat back components from the box.
- [6.2] Position the seat back on the seat base.
- [6.3] Install the washers+bolts(4x) {A} and fasten the bolts [15 Nm | 233 Lb-In]{B}.
- [6.4] Remove the on/off-switch.

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



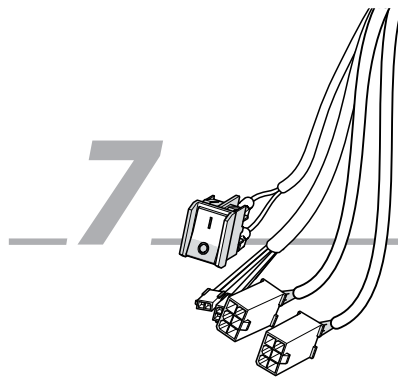
SEAT WIRING

SEAT LOOM

IN THE COLOURED SCHEME IN THE BACK OF THIS MANUAL YOU CAN FIND THE MARKED CONNECTIONS.

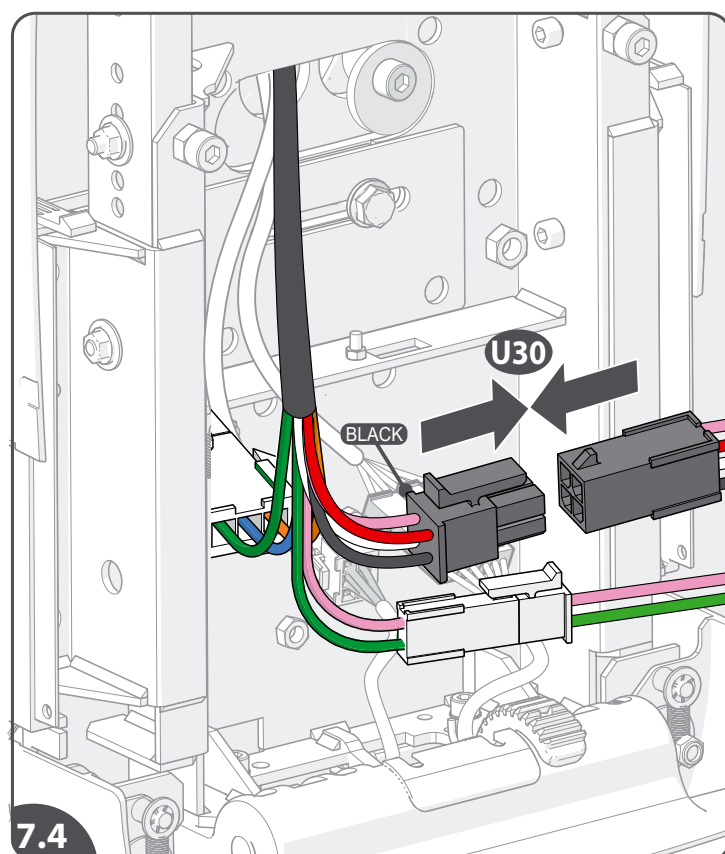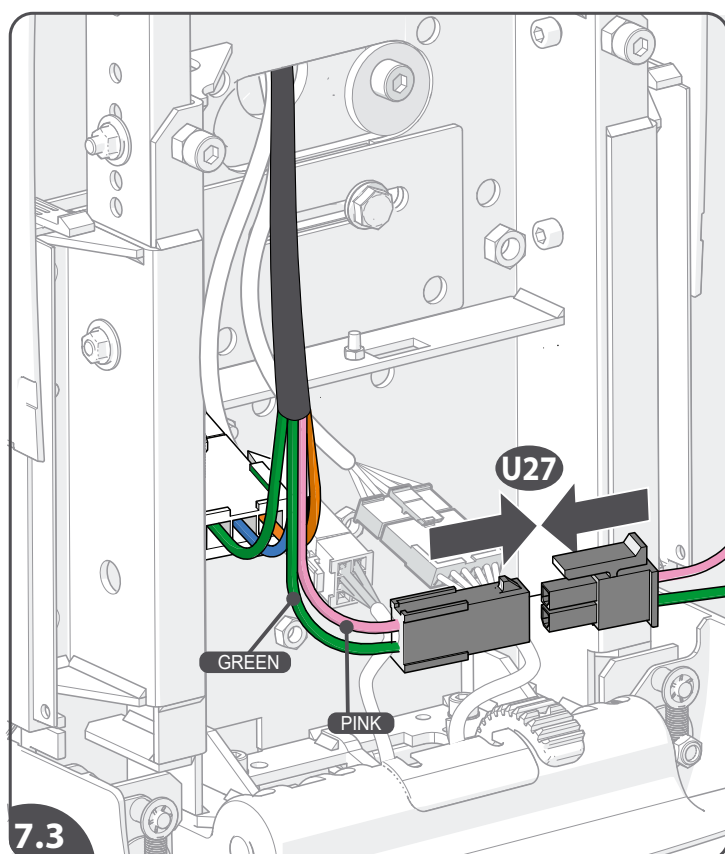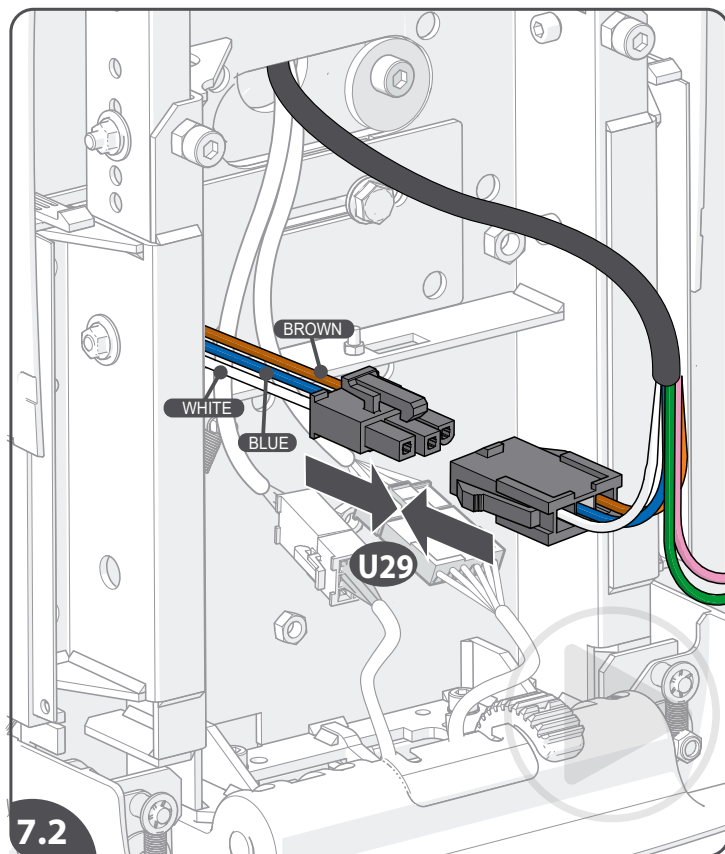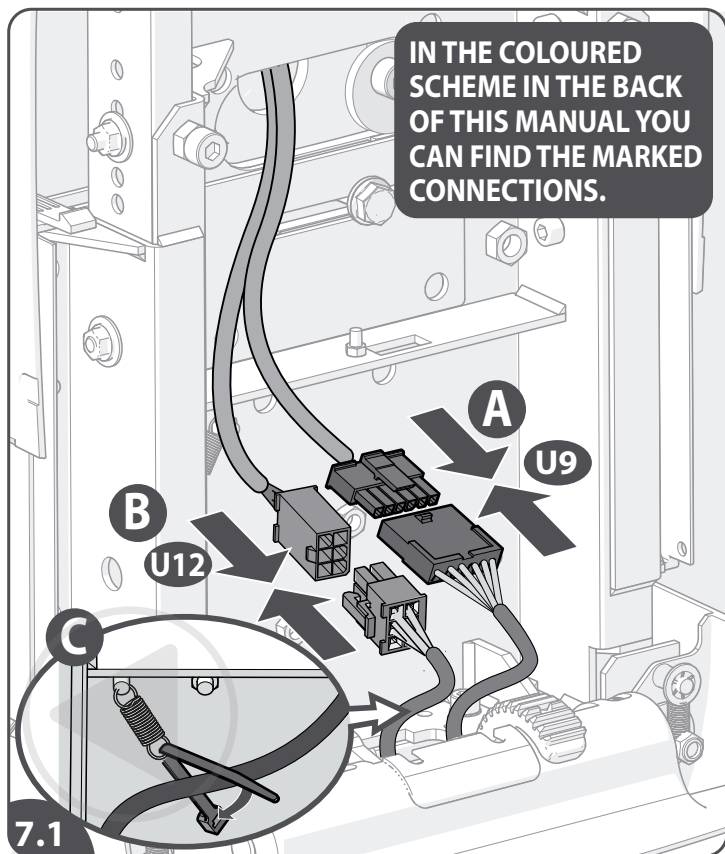- [7.1] Connect the underpan loom (U9){A}. Connect the footrest loom (U12){B}. Fasten this loom with a tie-wrap to the chassis frame {C}.
- [7.2] Connect the toggle loom (U29).
- [7.3] Connect the keyswitch loom (U27).
- [7.4] Connect the black connector (U30).

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
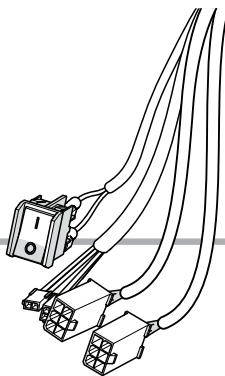
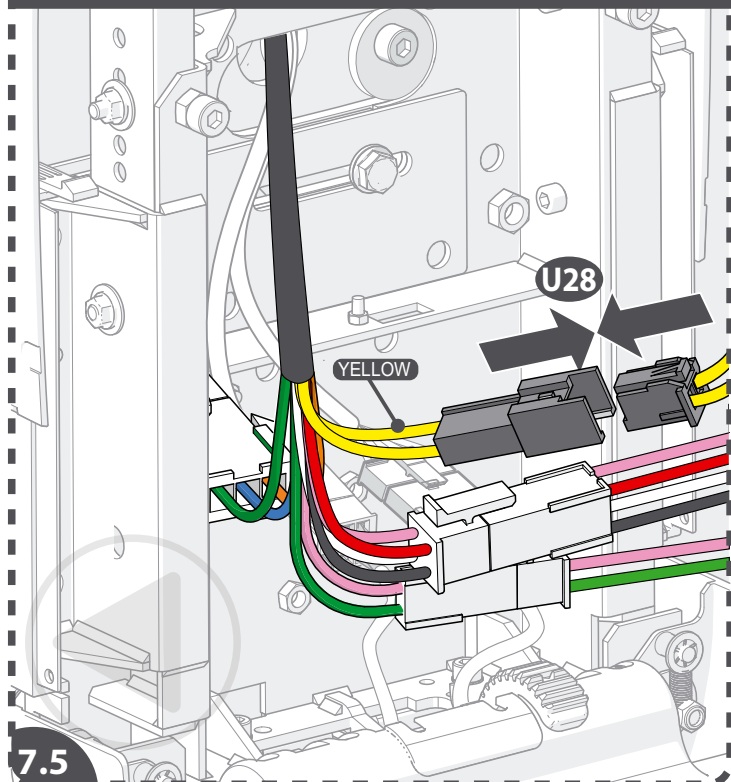
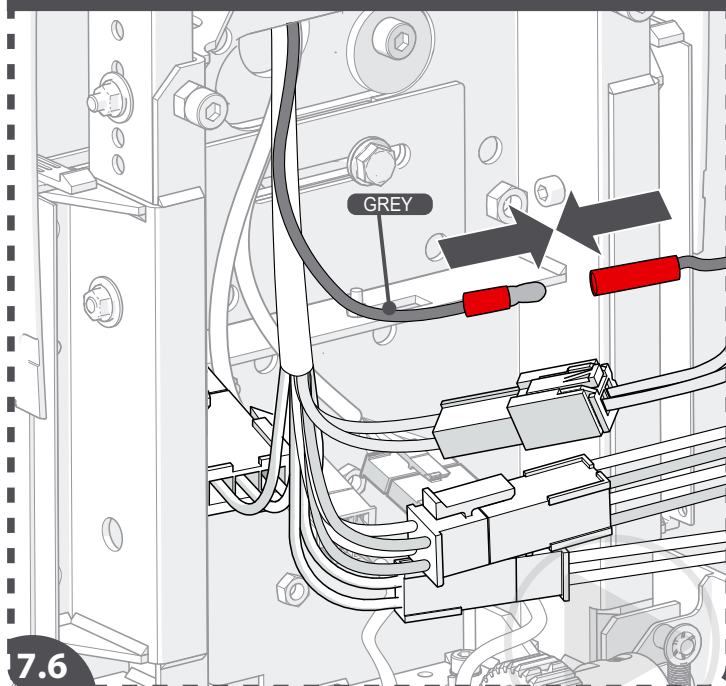
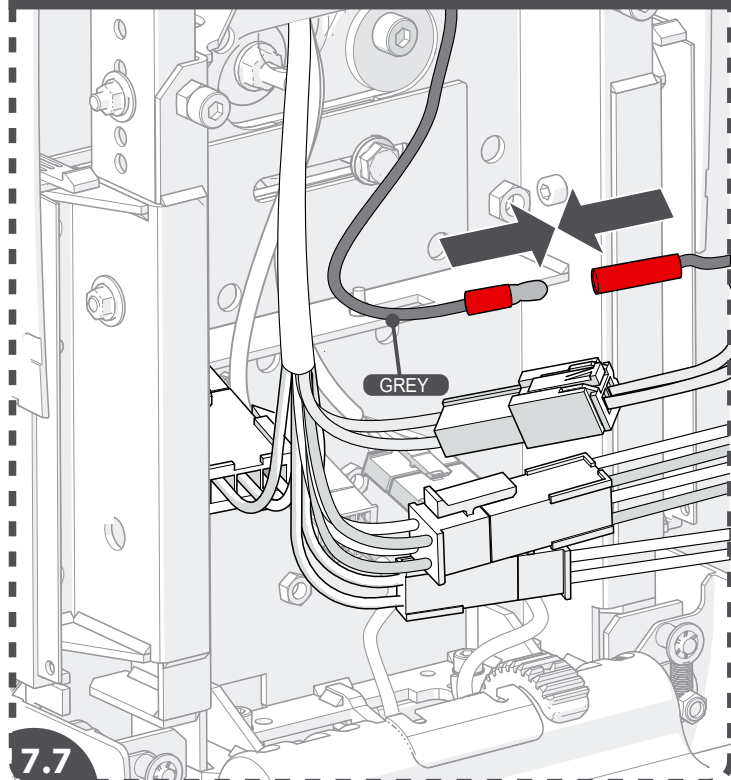
ONETOUCHALERT

▶ AA30046

See the specific manual for more information

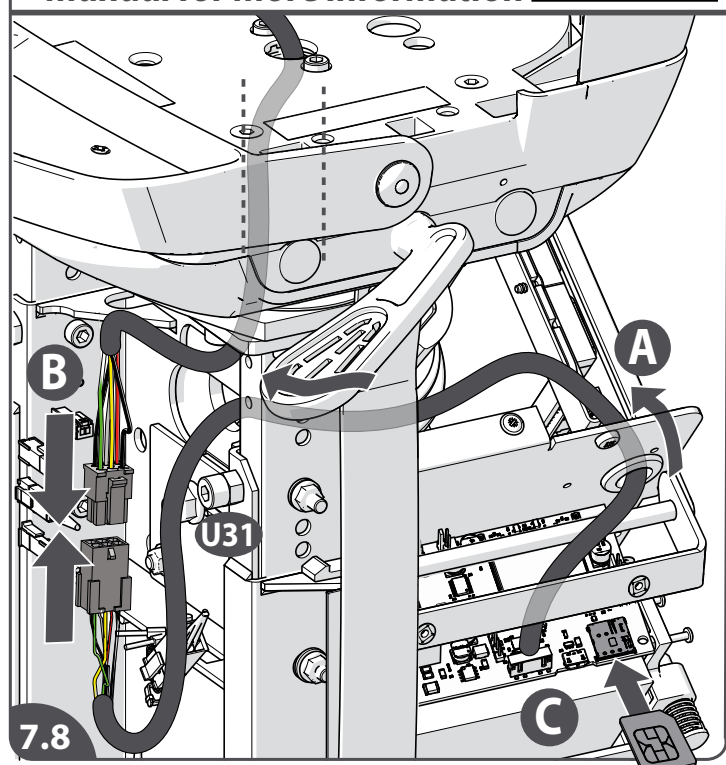- [7.5] *Only the powered footrest*: connect the powered footrest loom (U28).
 [7.6] *Only the powered footrest*: connect the communications cable (grey wire).
 [7.7] *Only the powered swivel*: connect the communications cable (grey wire). In case of a dual swivel (HAN126-172001.52088) change the dipswitches settings (1=OFF, 2=ON, 3=ON, 4=ON).
 [7.8] If present on the lift, guide the OTA loom {print part} through the power pack and chassis {A}. Connect the OTA loom parts (U31) {B}. Please refer to instruction AA30046 for more information. Install the SIM card {C}.

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

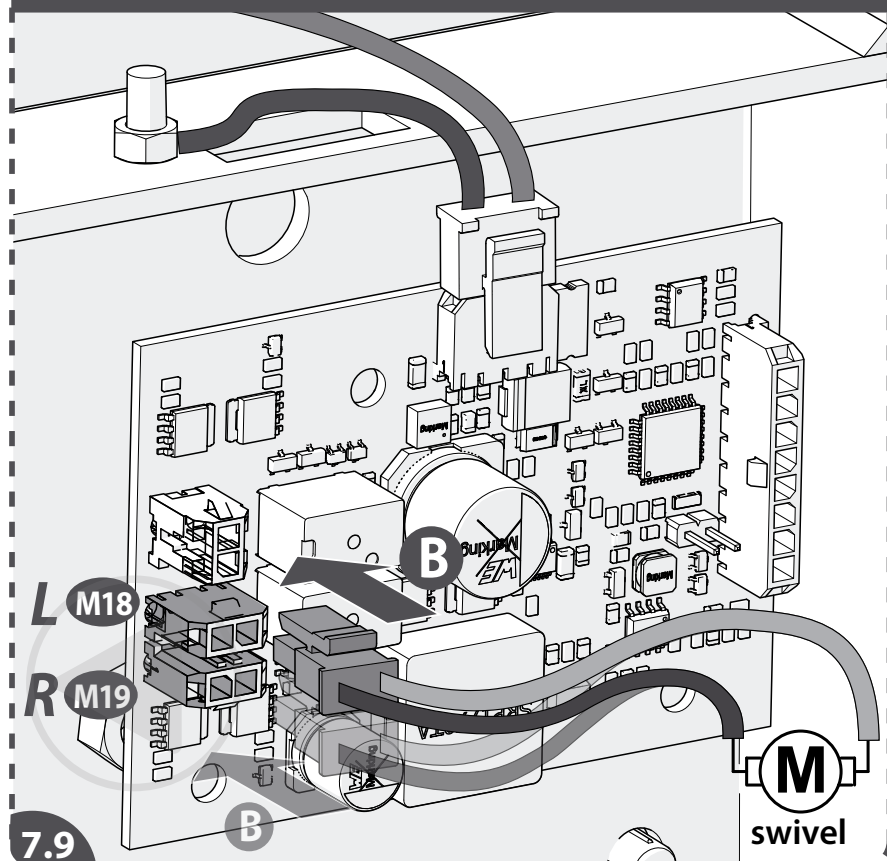
ONLY DUAL SWIVEL

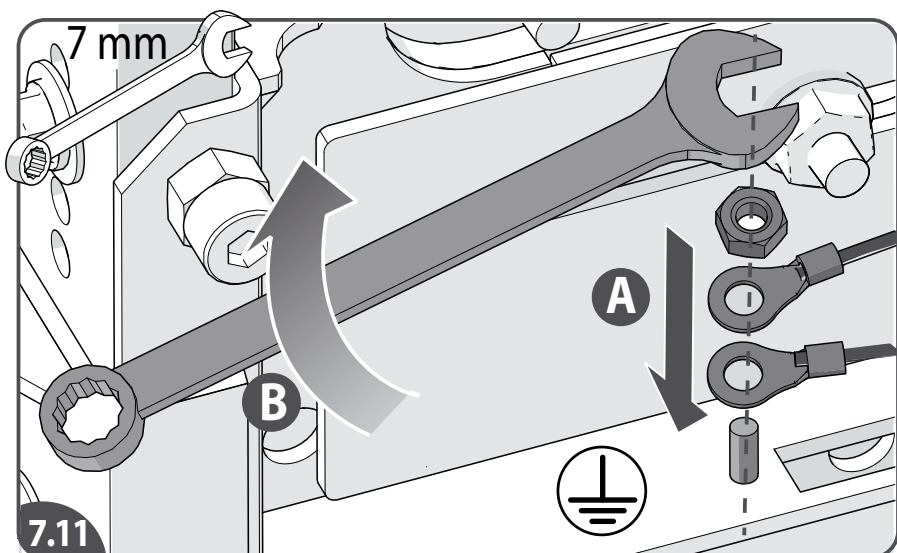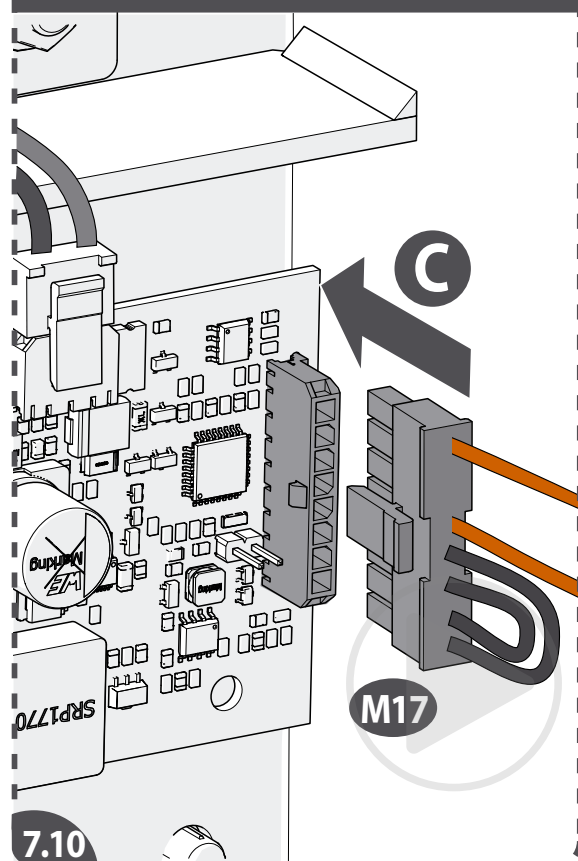Powered options:

[7.9] Connect the swivel motor loom with the PCB (LH Installation to M18, RH Installation to M19).

[7.10] Connect the dual swivel loom (M17).

[7.11] Install {A} and fasten the earth wires {B}.

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



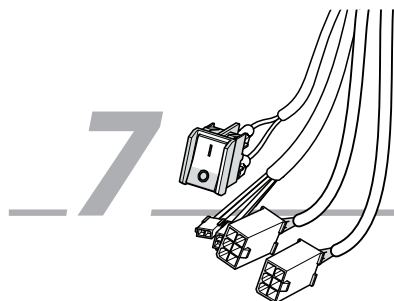
SEAT BASE WIRING

ONLY POWERED FOOTREST

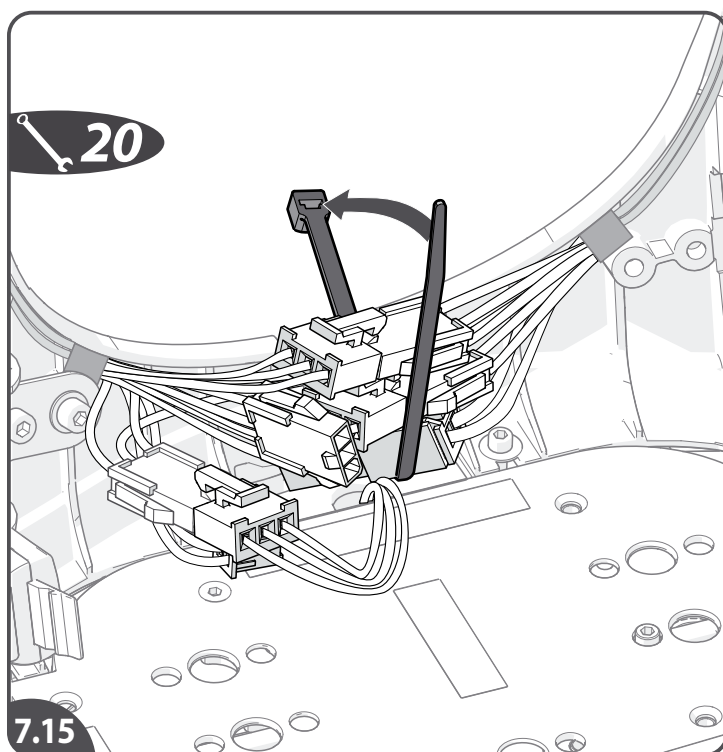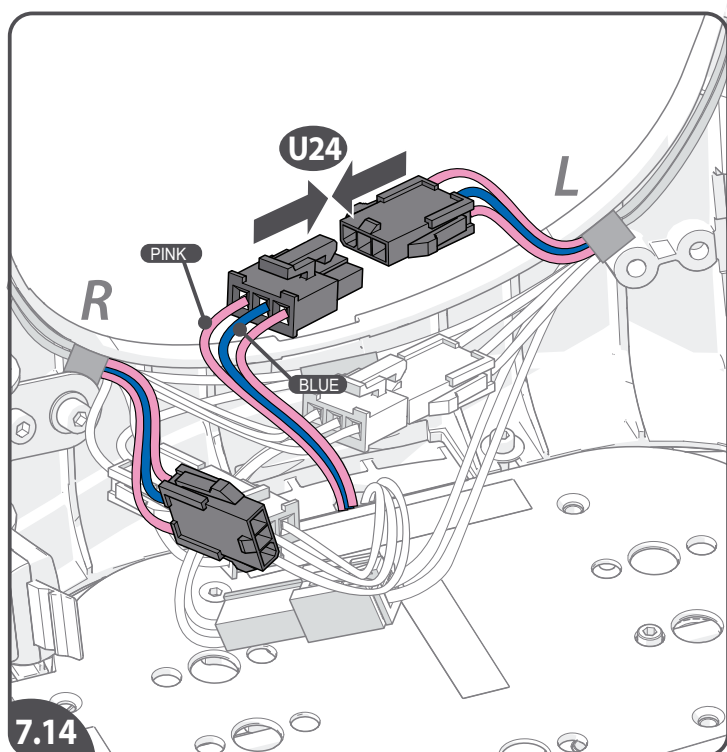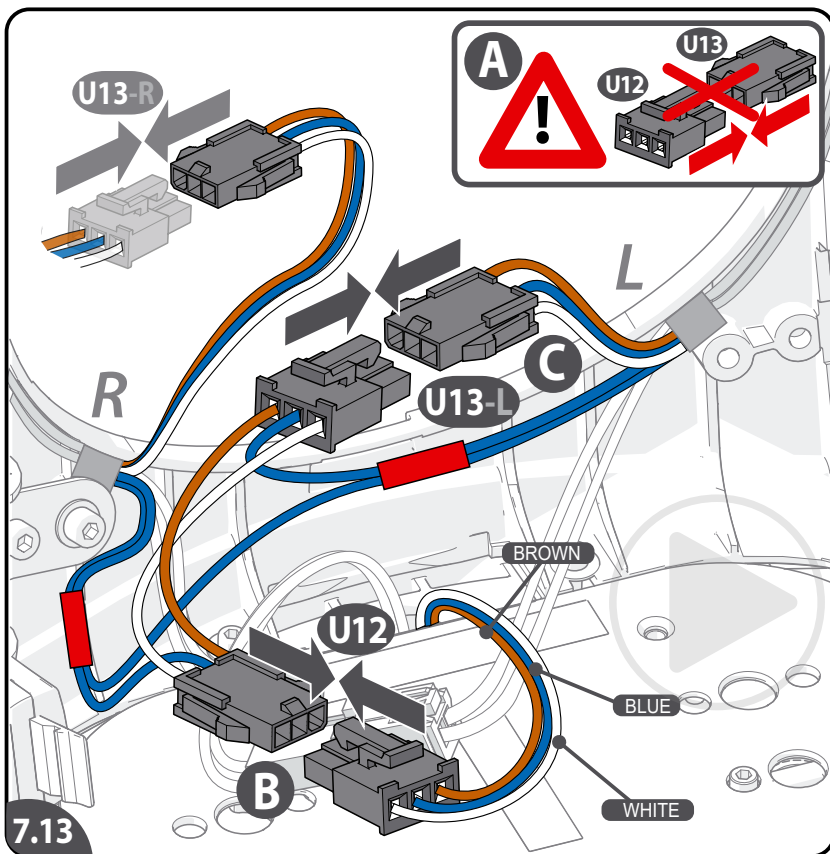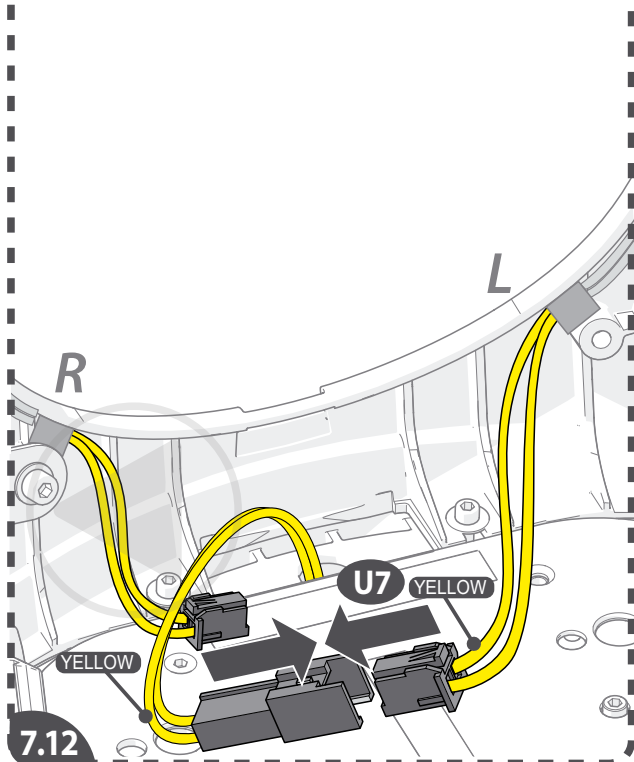- [7.12] *Only the powered footrest*: connect the footrest loom to the footrest switch {U7} (pay attention to the correct armrest!).
- [7.13] Connect the toggle loom to the armrest safety loom {U12}{B}. Do not bypass the armrest safety loom by connecting the toggle loom directly to the toggle {A}. Connect the armrest safety loom to the toggle {U13} (pay attention to the correct armrest!) {C}.
- [7.14] Connect the e-stop loom to the stopbutton {U24} (look for the correct arm rest!).
- [7.15] Tie the looms/connectors together with a tie-wrap.

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
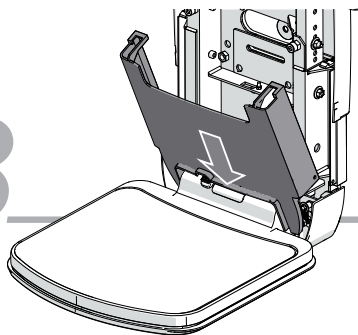
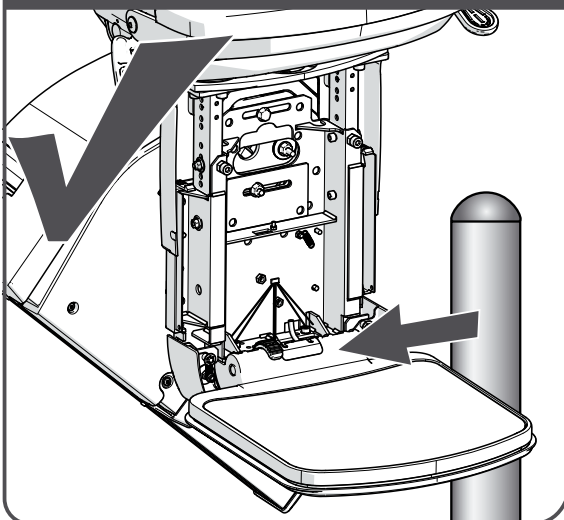
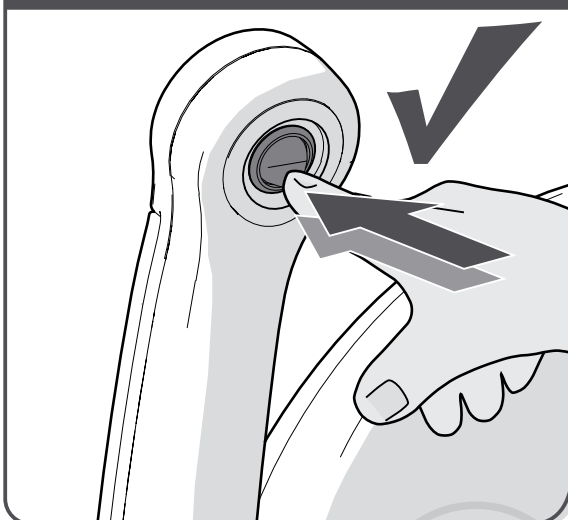
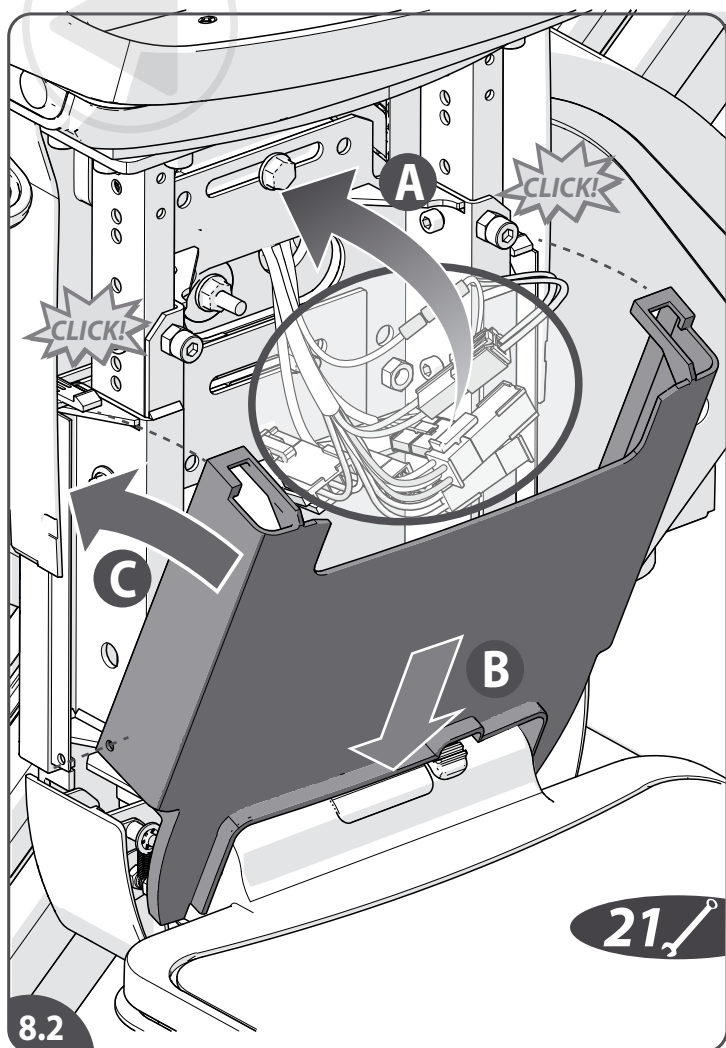
8.2

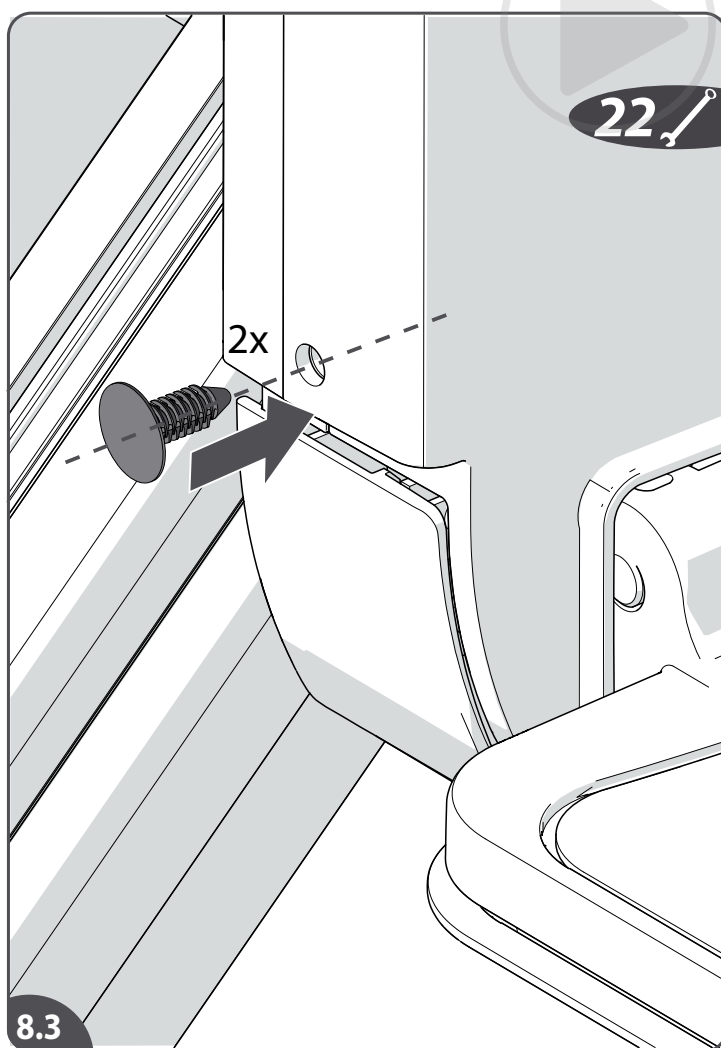8.3

- [8.1] Connect the on-/off switch{A}. Check the functions and safeties before installing the covers. Look at chapter 10 how to do this.
- [8.2] Place the wiring neatly in the upper part of the chassis frame{A}. Slide the bottom front cover in the bottom of the chassis{B}. Turn the cover in the right position{C}.
- [8.3] Fit the two blanking plugs in the holes of the cover.

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



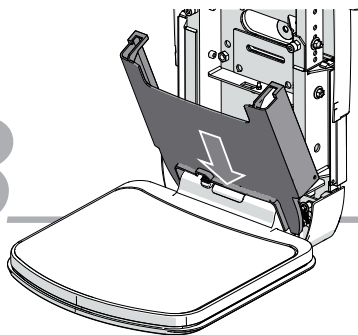
FINISH

CONNECTIONS TOP FRONT COVER

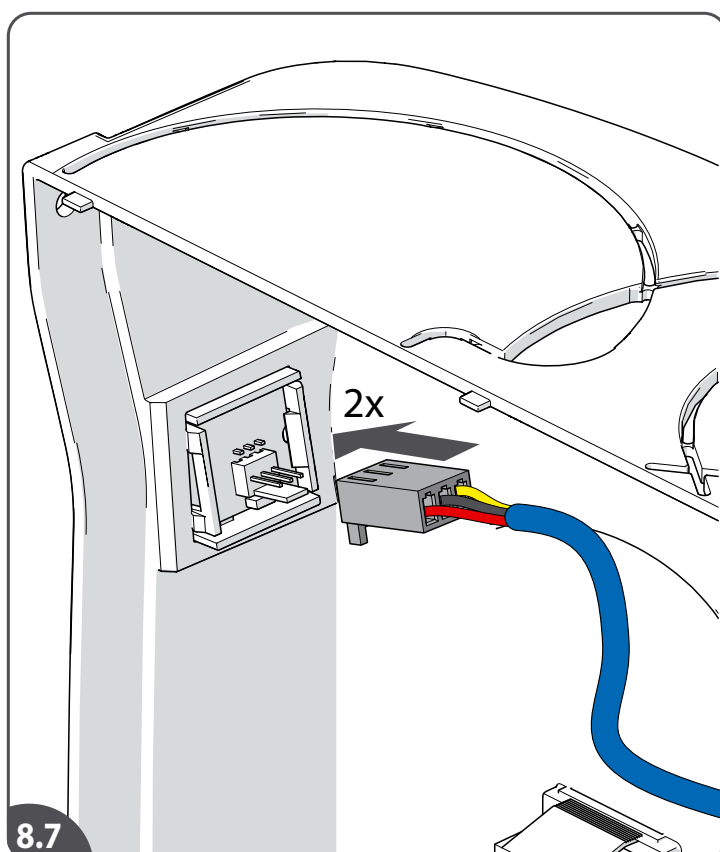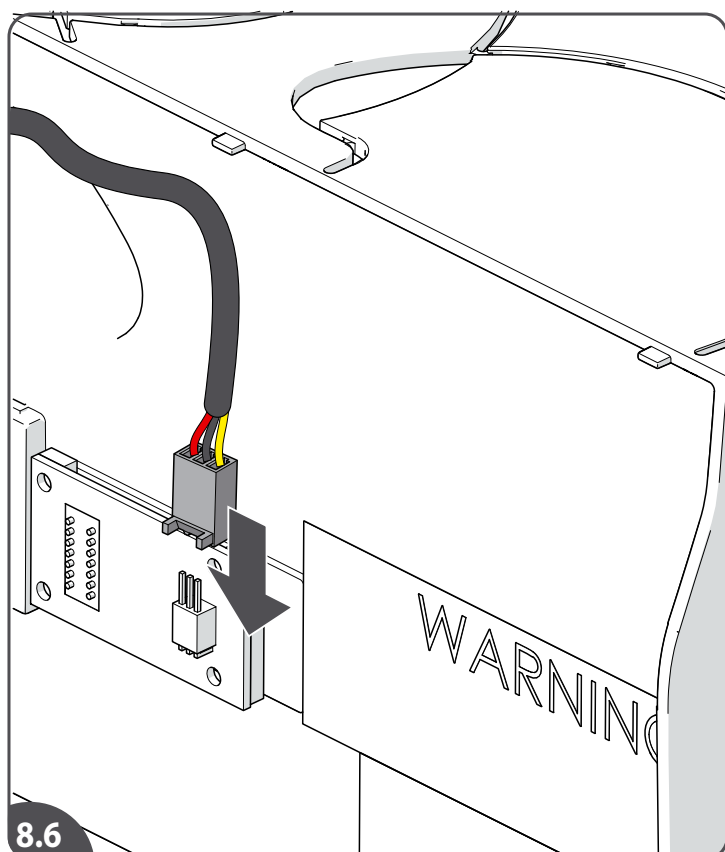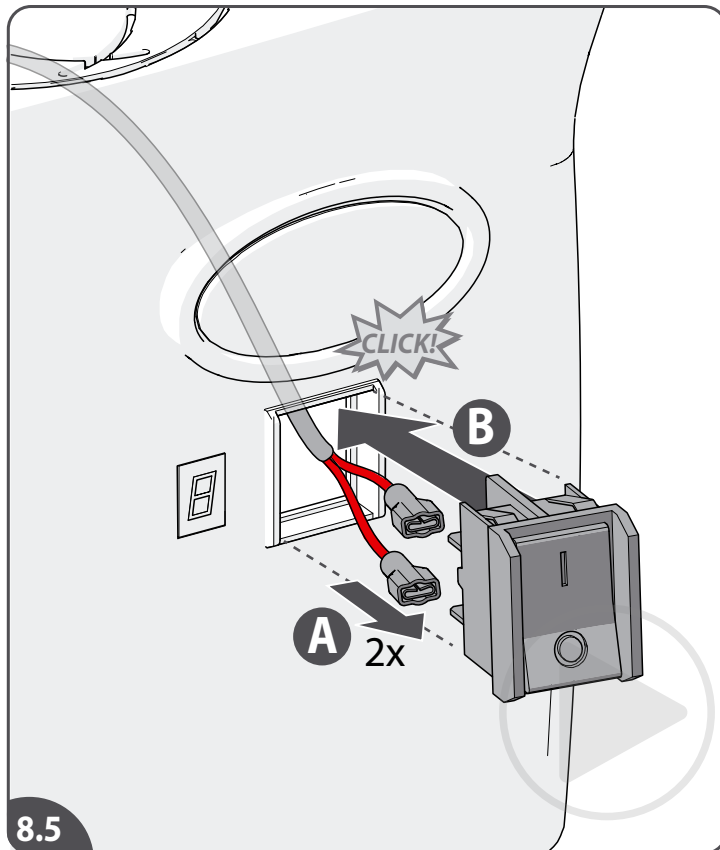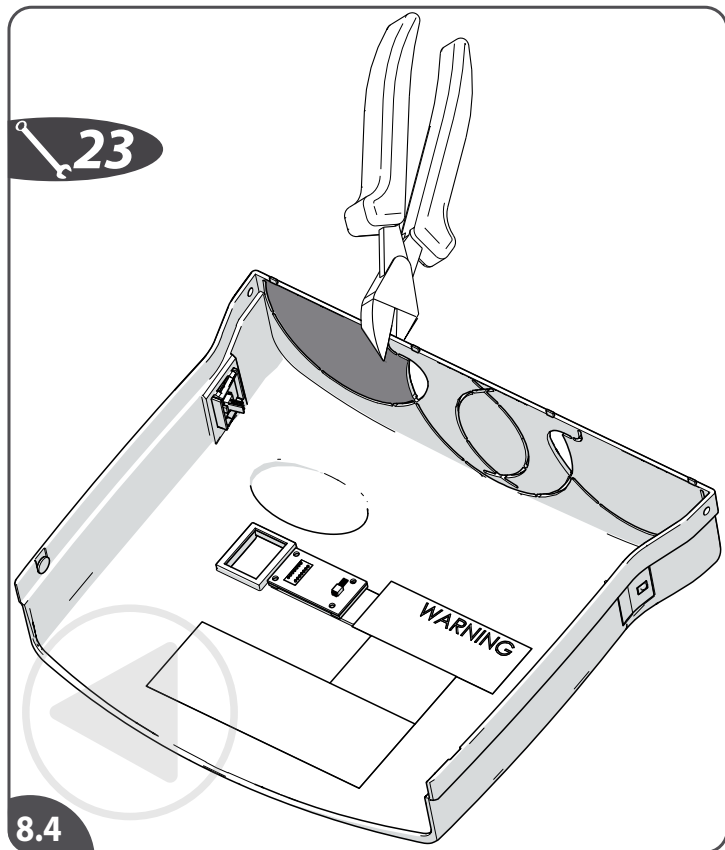- [8.4]If not done yet: cut out the appropriate groove in the upper front cover.
- [8.5]Reconnect the on/off-switch{A}. Click the on-/off switch in the front cover{B}.
- [8.6]Connect the diagnostic display.
- [8.7]Connect the two infrared sensors (indoor lifts only).

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



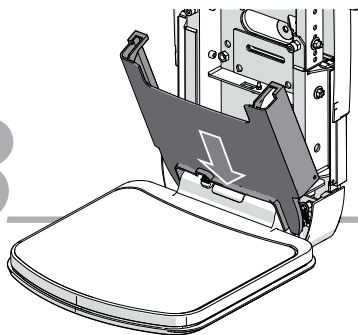
FIXATION TOP FRONT COVER

RF(outdoor) ONLY

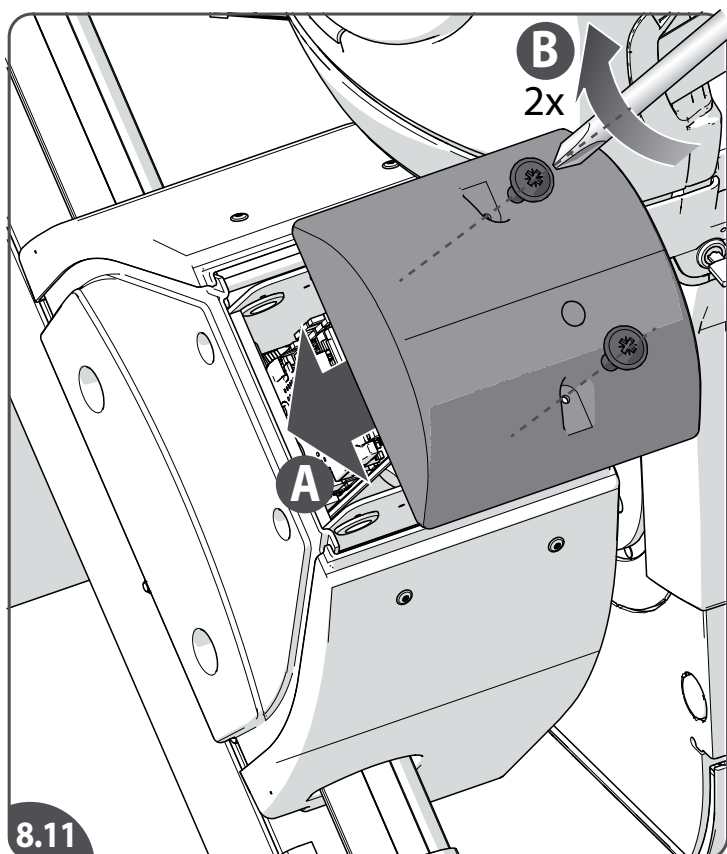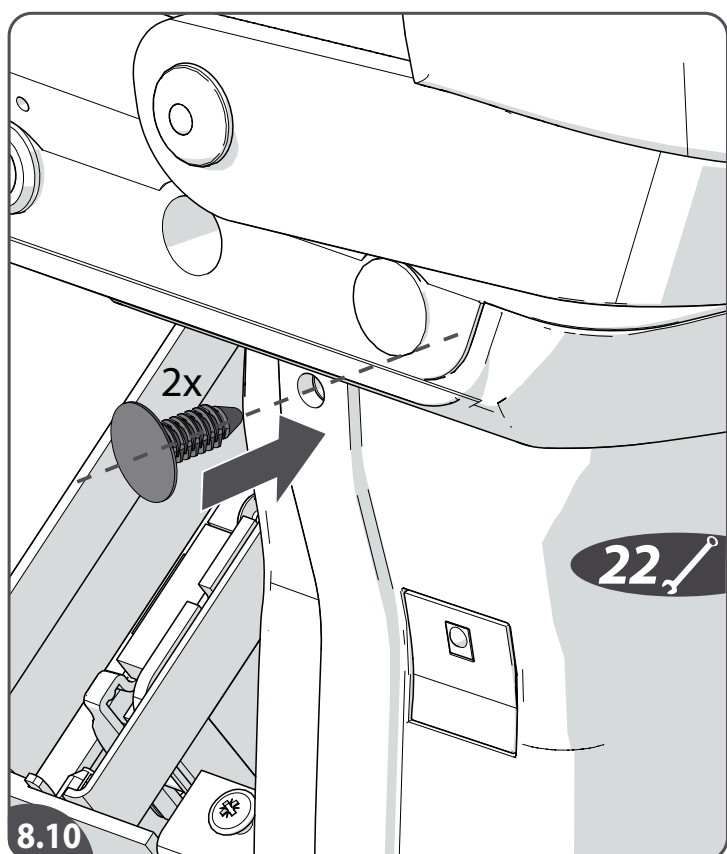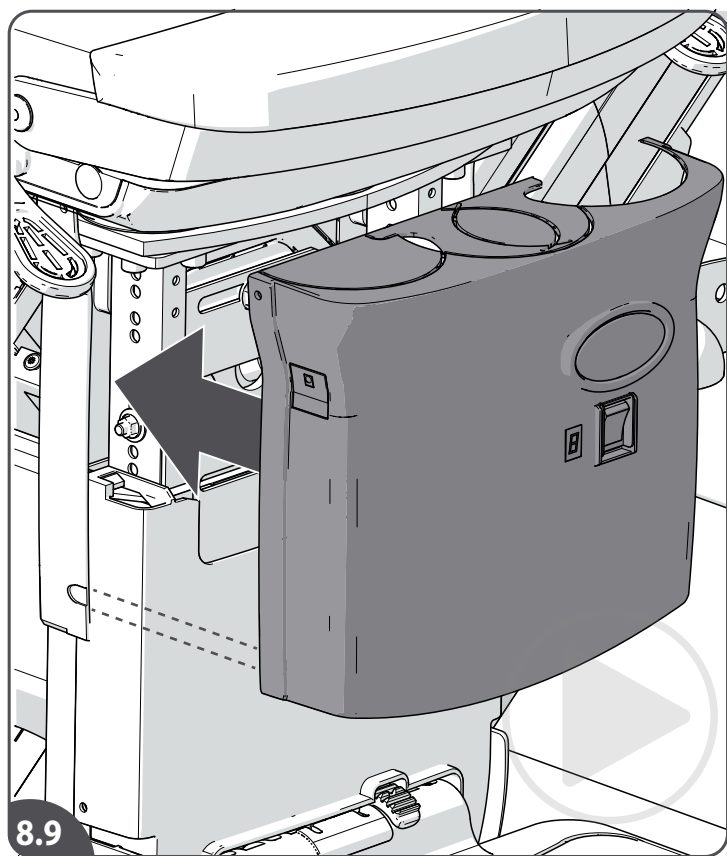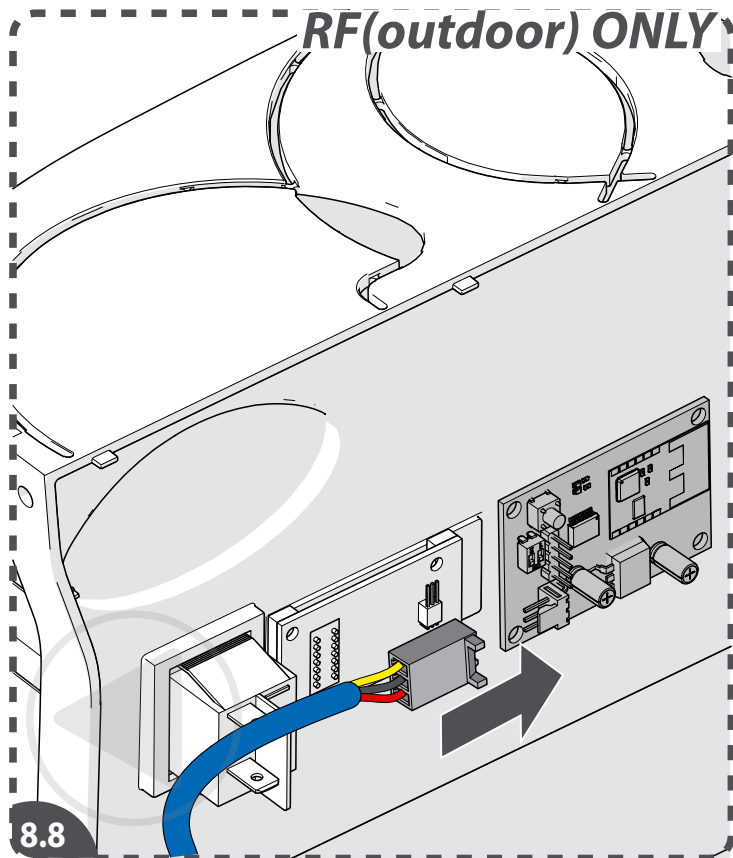Only RF remote controls (outdoor): [8.8] Connect the RF board.
 [8.9] Slide the upper front cover in position.
 [8.10] Fit the two blanking plugs in the holes of the cover.
 [8.11] Position the power pack cover {A} and fasten the screws {B}.

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





8

FINISH UPHOLSTERY

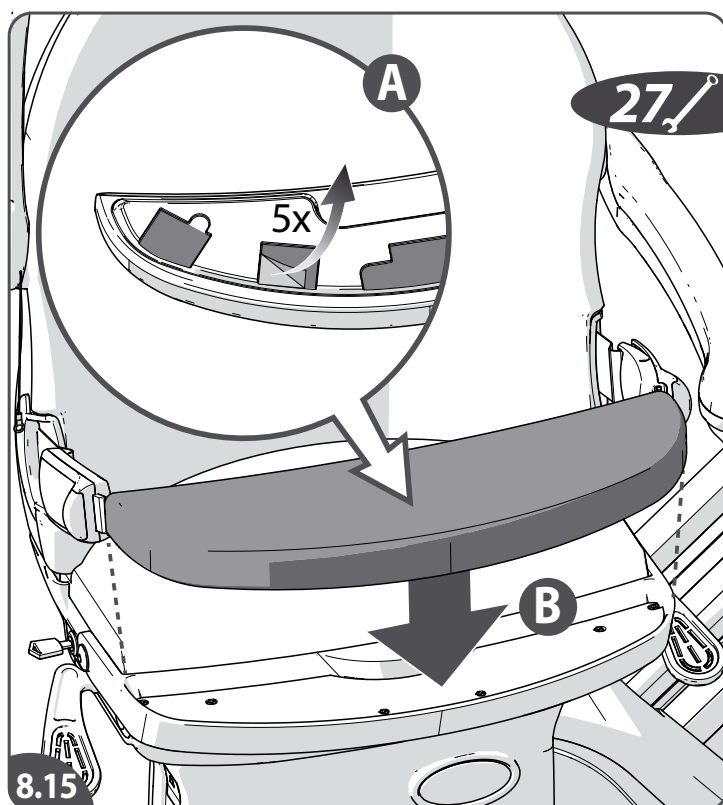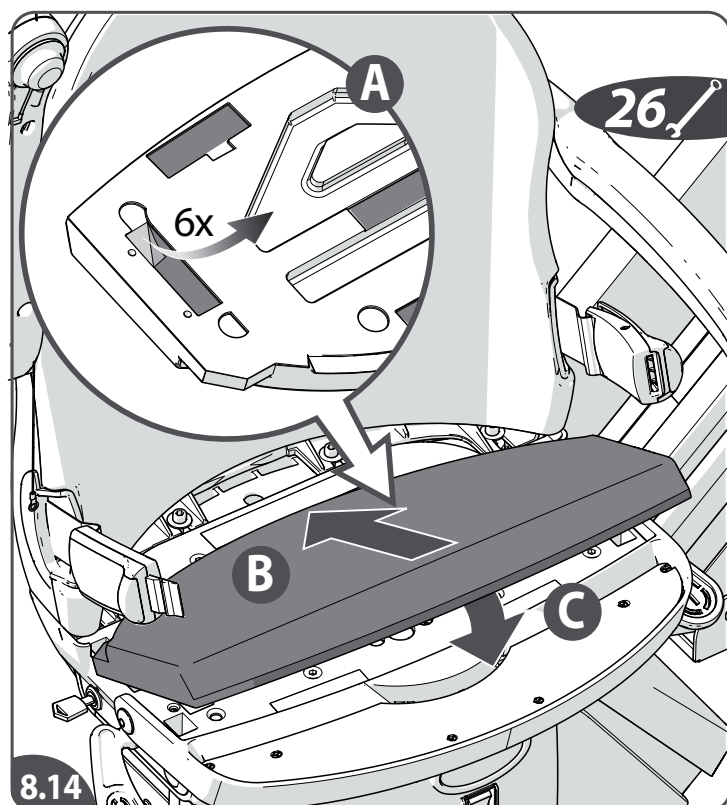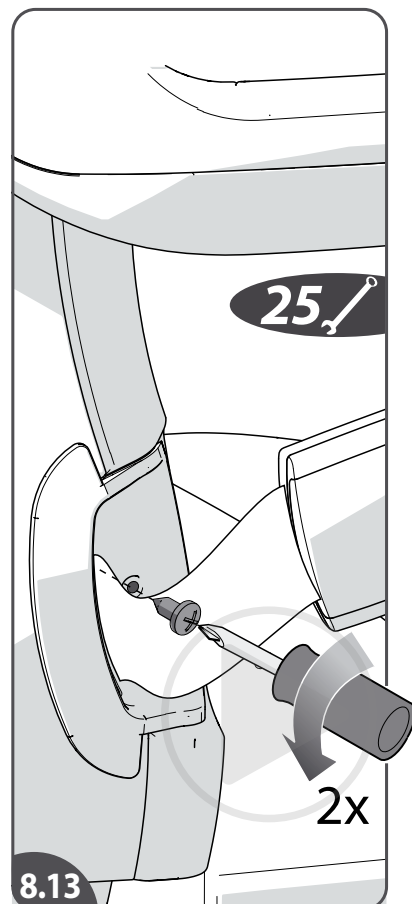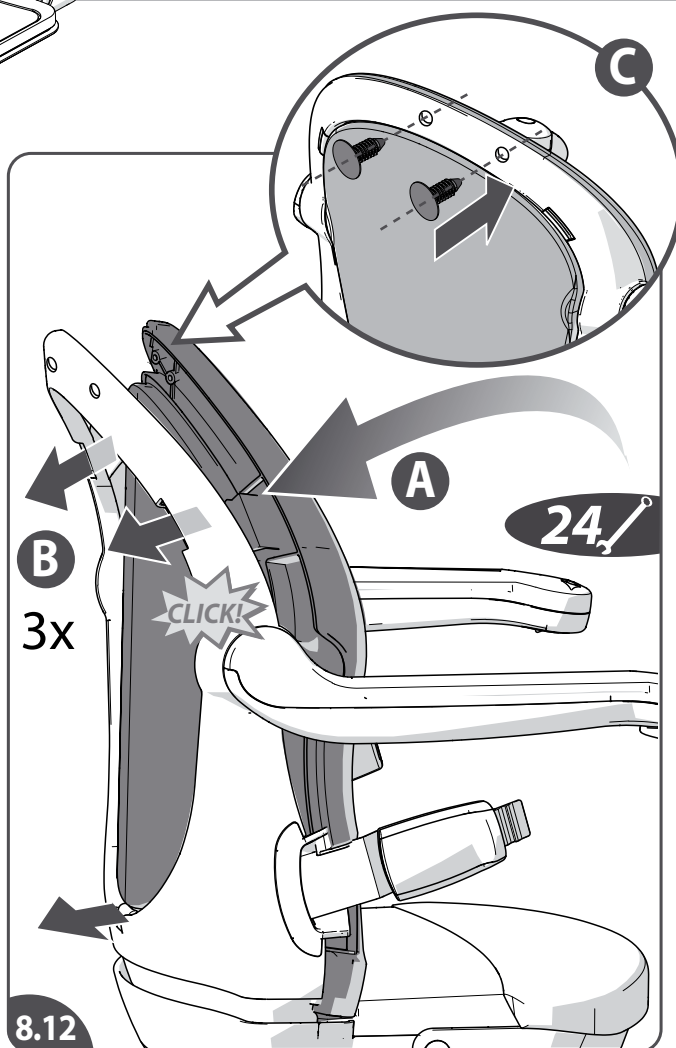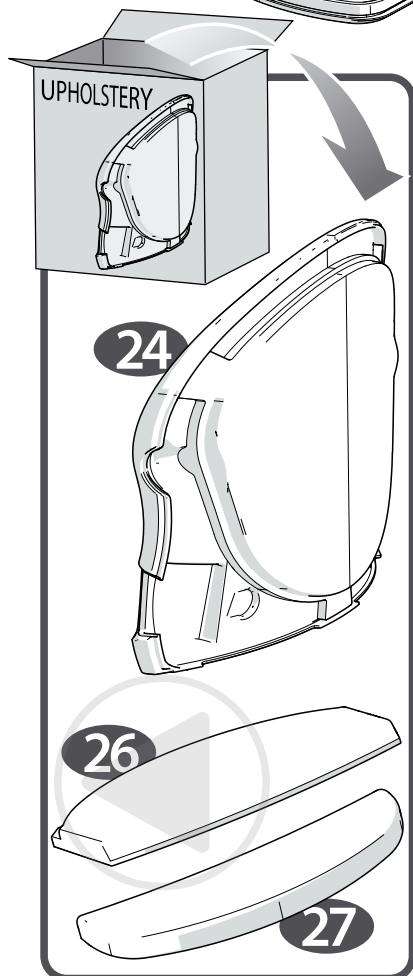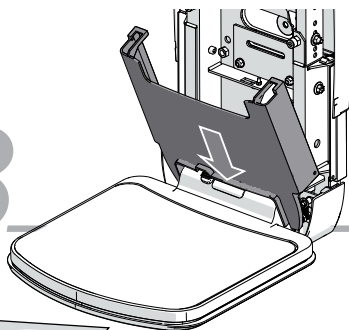Remove the upholstery components from the box.

[8.12] Install the back seat upholstery {A}. Press firmly on the locks to fix the upholstery {B}.

[8.13] Fasten the upholstery with the screws (2x).

[8.14] Remove the backings of the velcro pads (6x) {A}. Fit the main seat cushion {B+C}.

[8.15] Remove the backings of the velcro pads (5x) {A}. Fit the seat flap cushion {B}.

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




IT

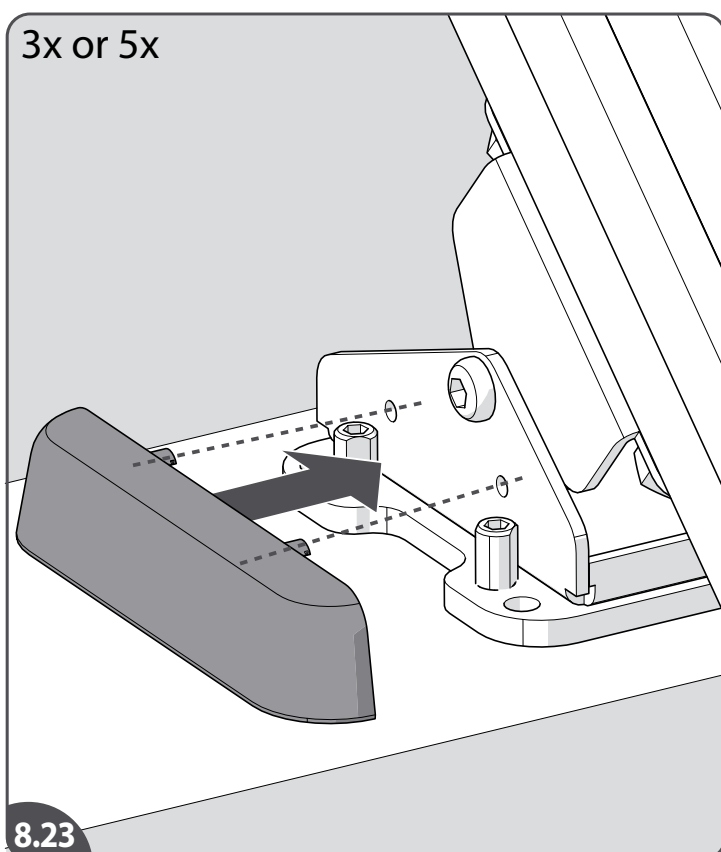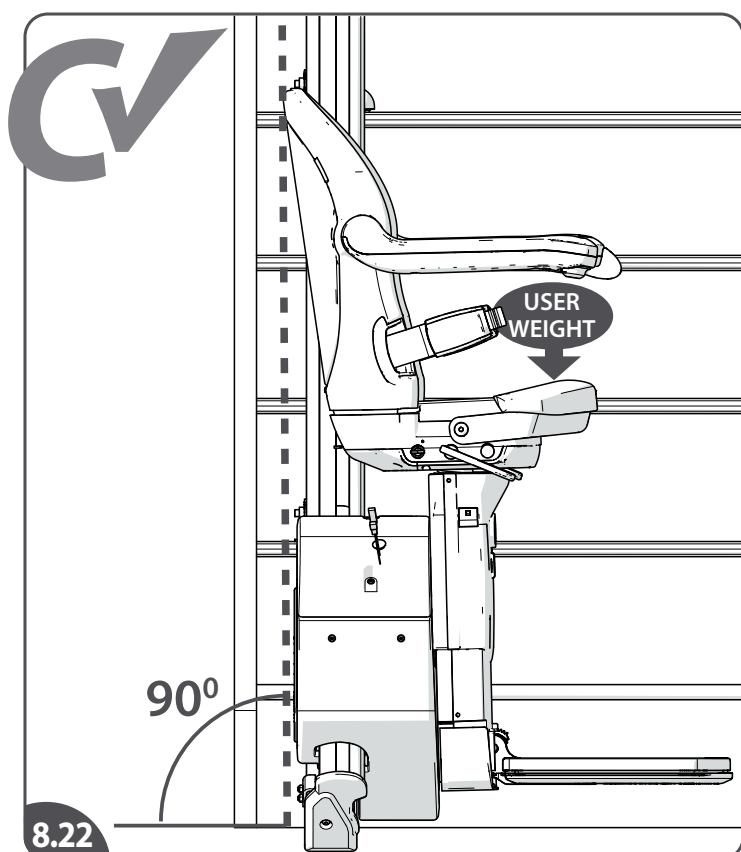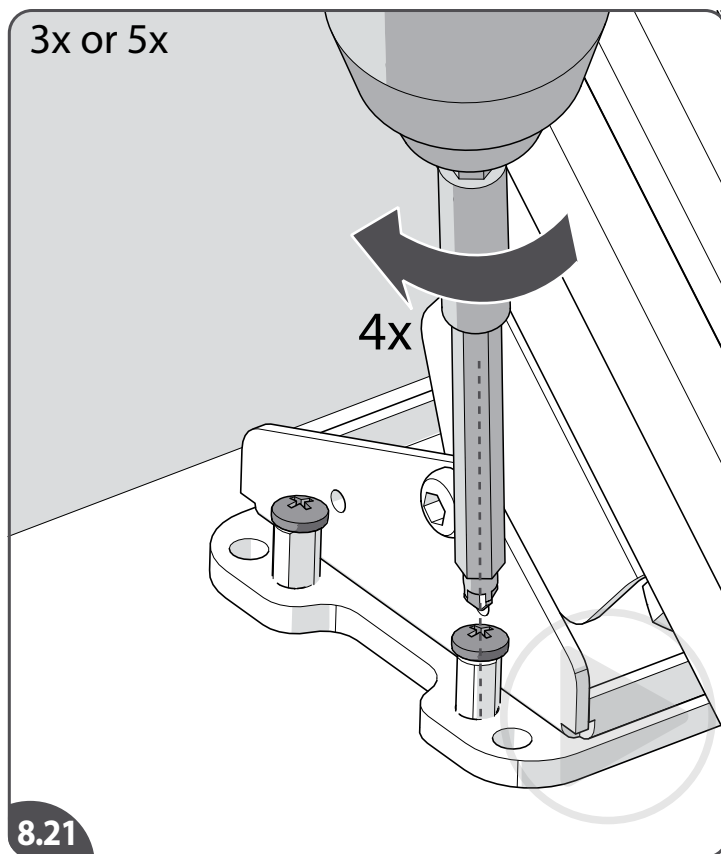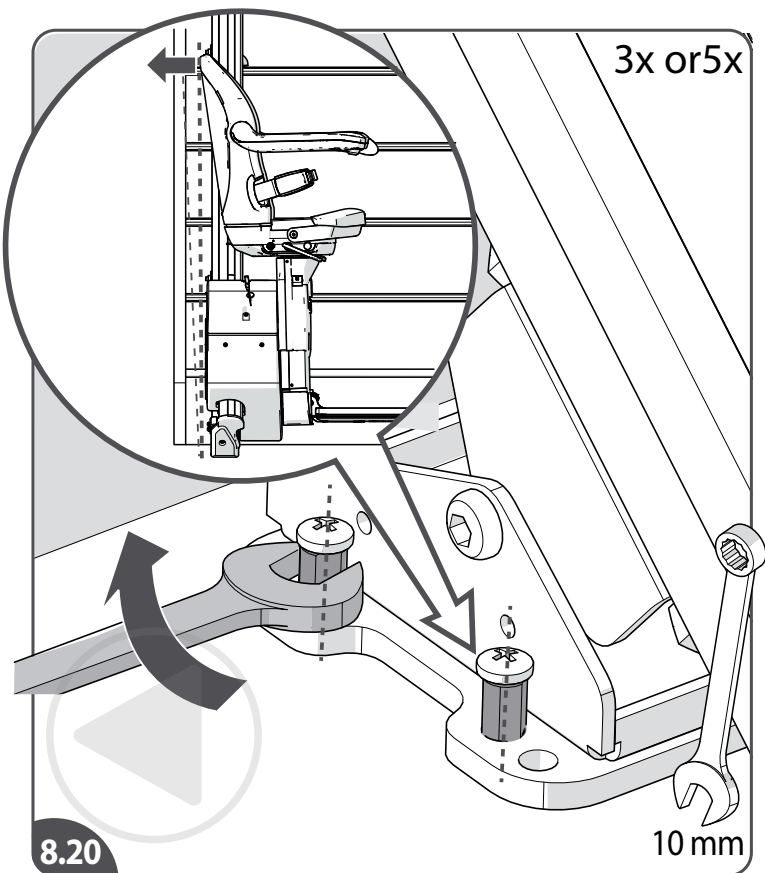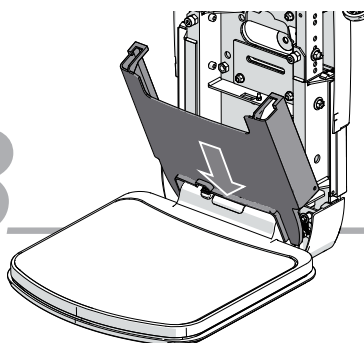[8.20]Vertically align the track by screwing the adjusters up or down on each foot. Adjust the vertical alignment of the track such that it is slightly angled towards the wall. The top of the seat should clear the wall by approximately 20 mm when the alignment is correct.

Note: You must ensure that the alignment of the track is consistent along the complete length.

[8.21]Fully tighten the three screws at each foot.

[8.22]Recheck the vertical alignment of the track and if necessary correct using the front foot adjuster.

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




8

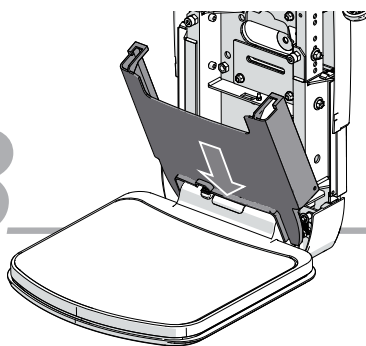FINISH ADJUSTMENTS

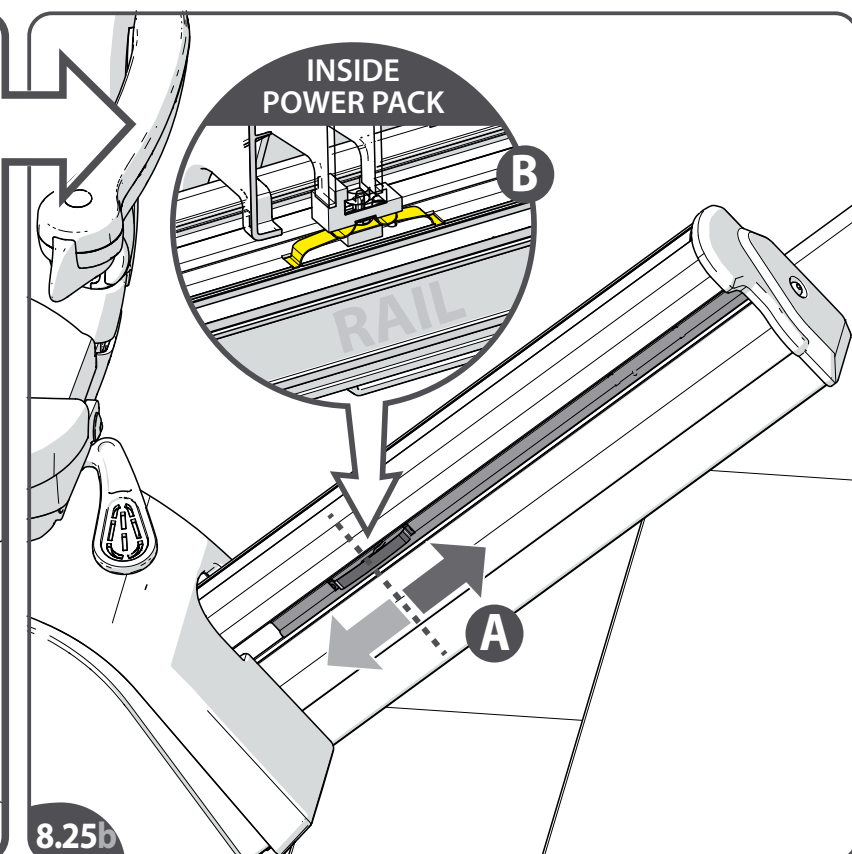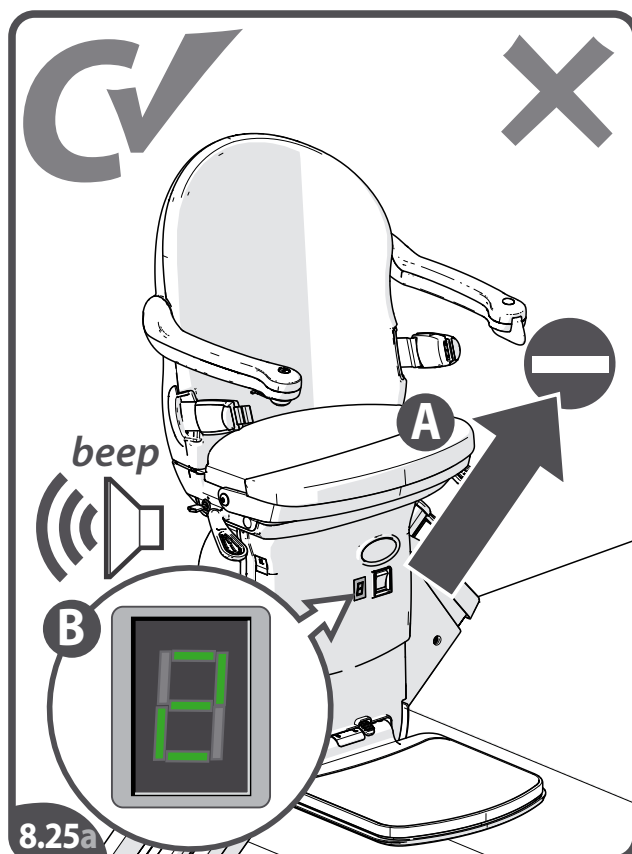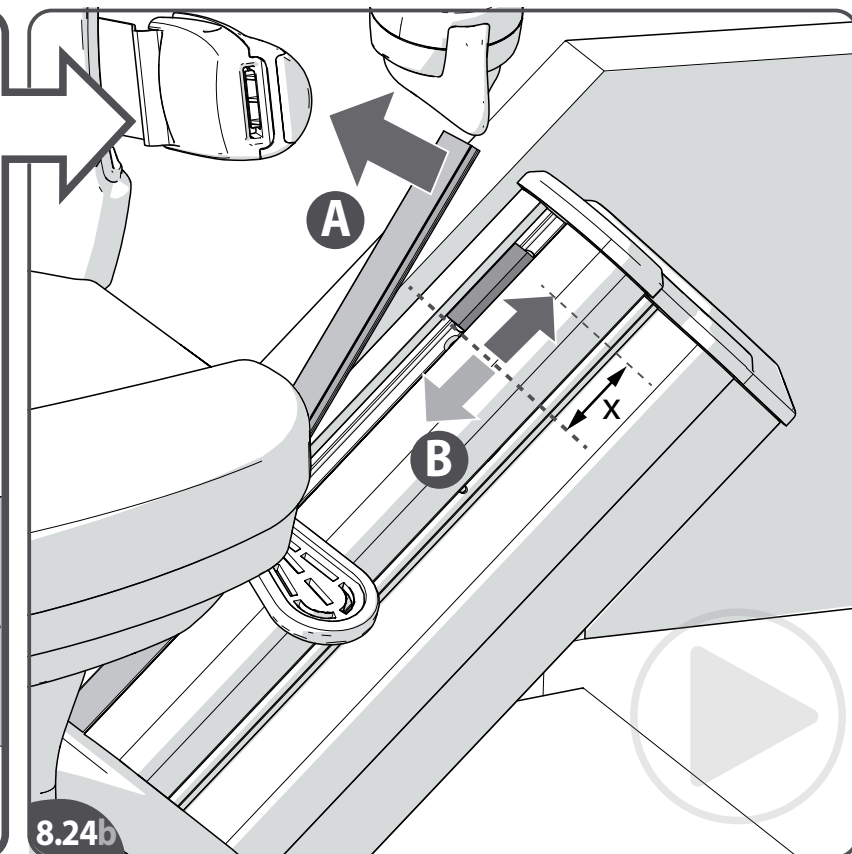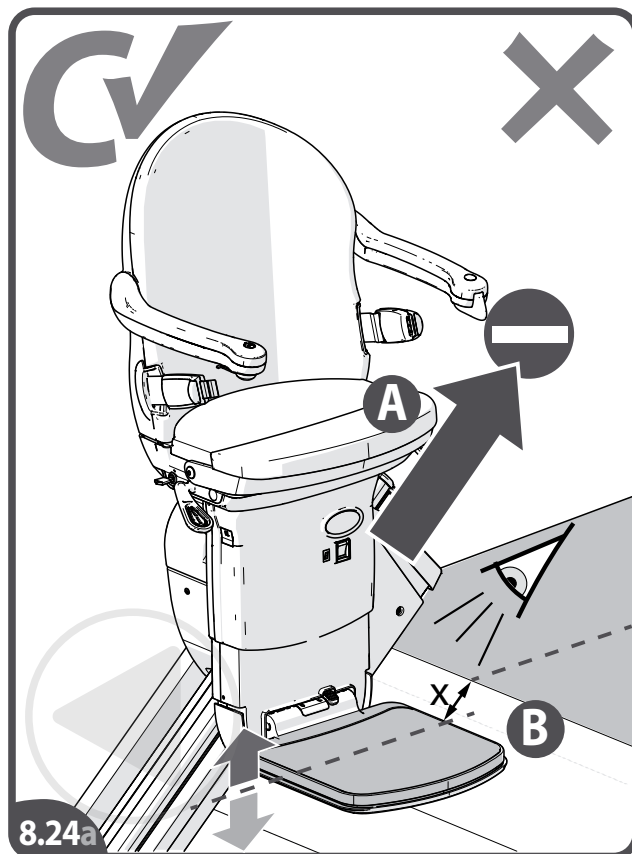[8.24a] Drive the stairlift to the top of the track{A}. Check if the footrest is levelling with the landing{B}.
 [8.24b] If the footrest is not levelling with the landing, drive the lift halfway down and remove the plastic cover on top of the rail{A}. Slide the magnet for the end limit up or down (with the distance the footrest is too high or too low){B}. Drive the lift back to the top of the stairs and check that the footplate now finishes level with the landing.
 [8.25a] Again drive the stairlift to the top of the track{A}. Check if the lift is charging{B}.
 [8.25b] If the lift is not charging (the lift beeps and the display shows a "2") adjust the charge contact; slide the contact slightly up or down{A}, until the contacts (power pack and rail) are aligned properly{B}.

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




IT

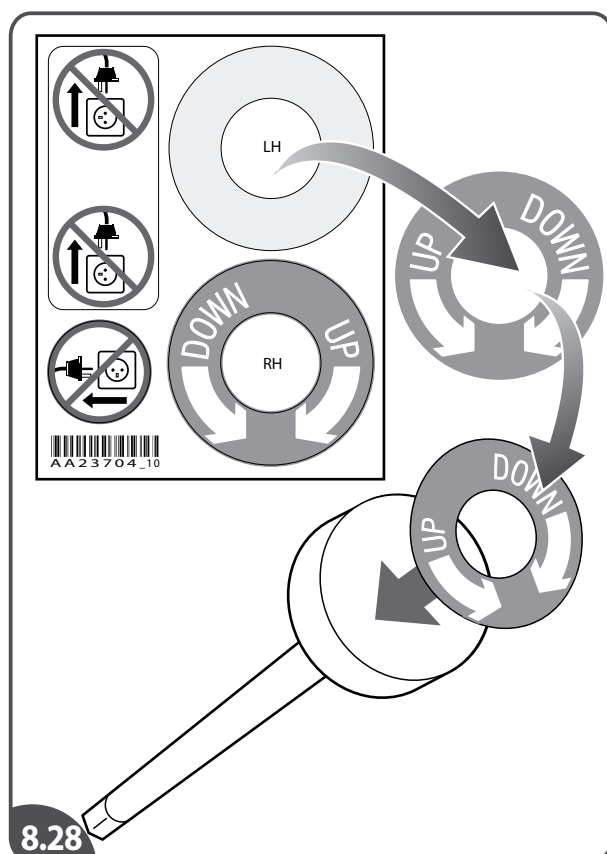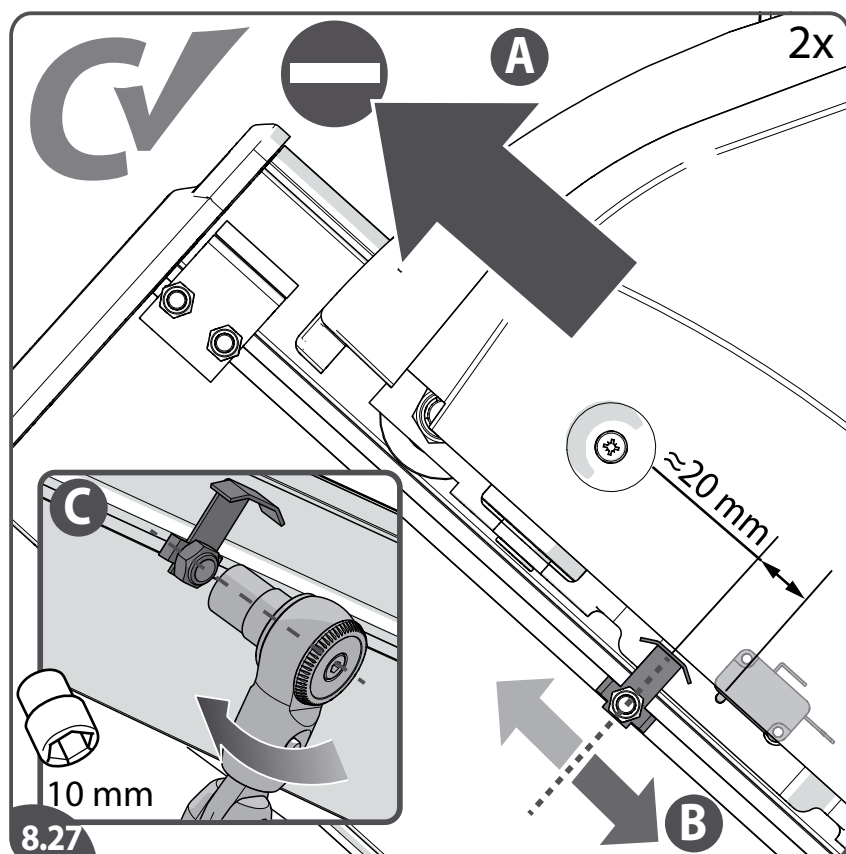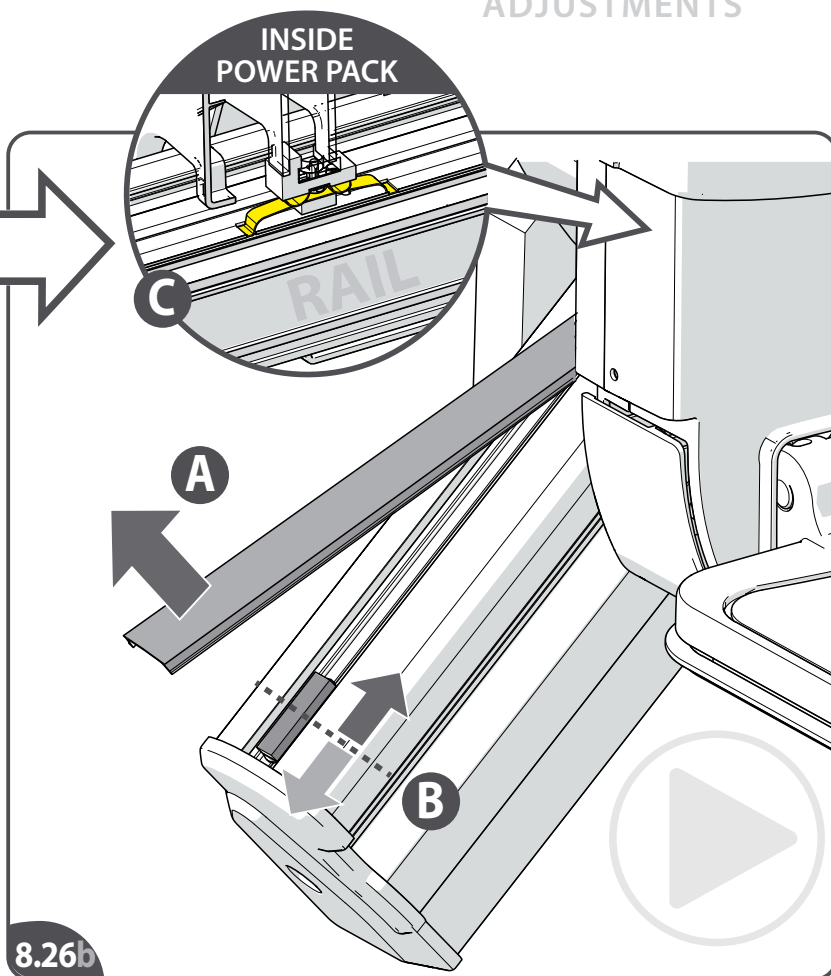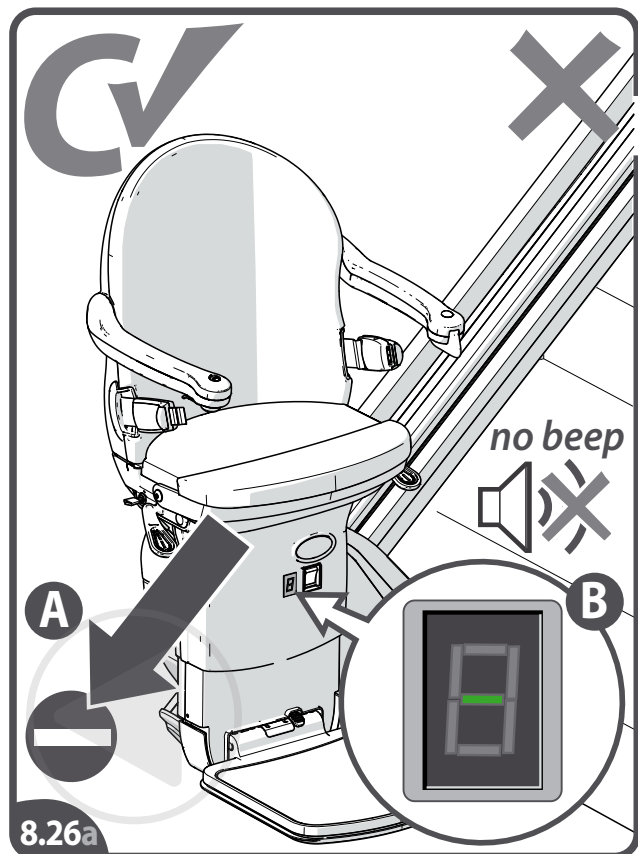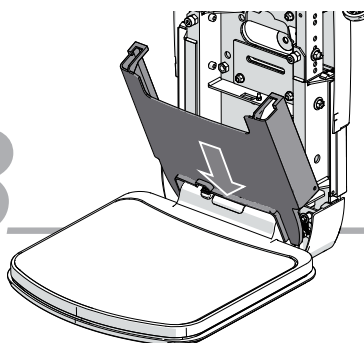[8.26a] Drive the stairlift to the bottom of the track {A}. Check if the lift is charging (the lift is not beeping and the display shows a "-") {B}. The end limit is factory set and it should work without adjustment.
 [8.26b] If the lift is not charging adjust the magnet for the end limit: drive the lift halfway down and remove the plastic cover on top of the rail {A}; slide the magnet slightly up or down {B}, until the contacts (power pack and rail) are aligned properly {C}.
 [8.27] Drive the lift to the top and bottom of the stairs and check the positions of the final limit cams (top and bottom) {A}. The cams should be positioned approx. 20 mm further than the end position of the final limit switch.
 If needed adjust {B} and refasten {C} the cams. The bottom final limit cam is factory set and should not need any adjustment. [8.28] Fit the correct label on the hand winder.

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


8

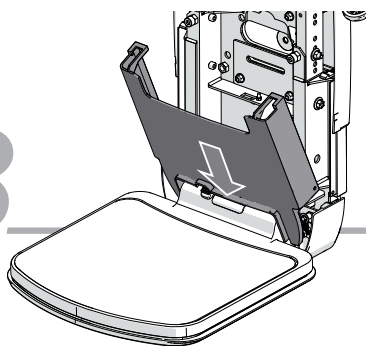FINISH
GREASE RAIL

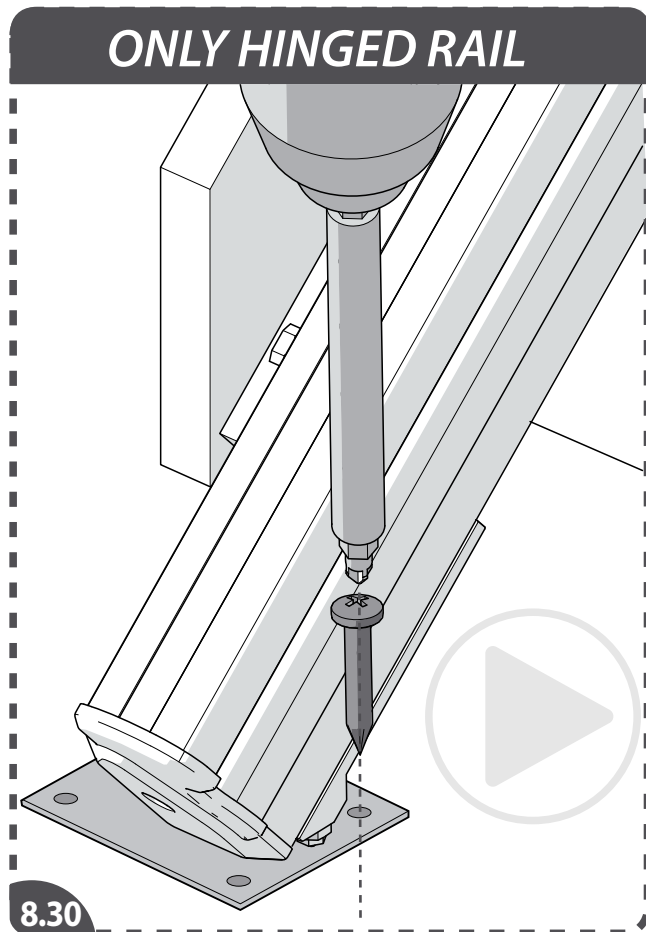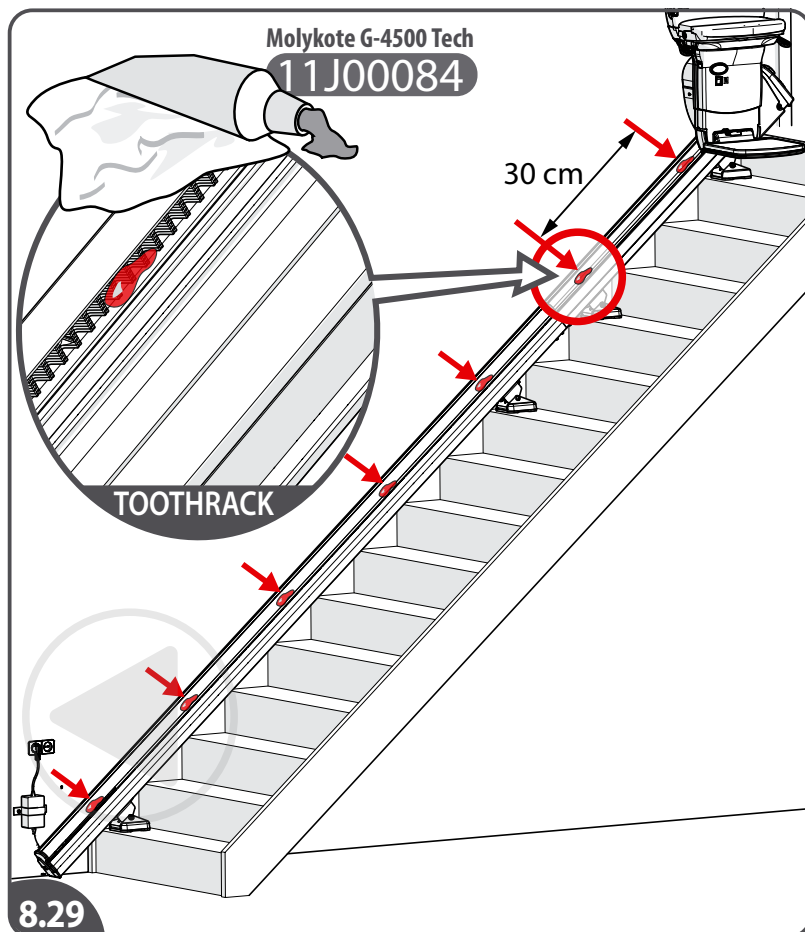[8.29] Lightly lubricate (like 1 pea drop every 30 cm) the teeth of the rack to allow the carriage to travel smoothly.
 [8.30] Fit the floor plate underneath the hinge on the ground floor. Lower the hinge manually and position the plate. Lift the hinge and screw the floor plate down using the fixing provided. Adjust the feet on the hinge so that when the hinge is fully closed the feet locate into the floor plate and support the end of the hinge section.

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




9

REMOTE CONTROL MOUNTING

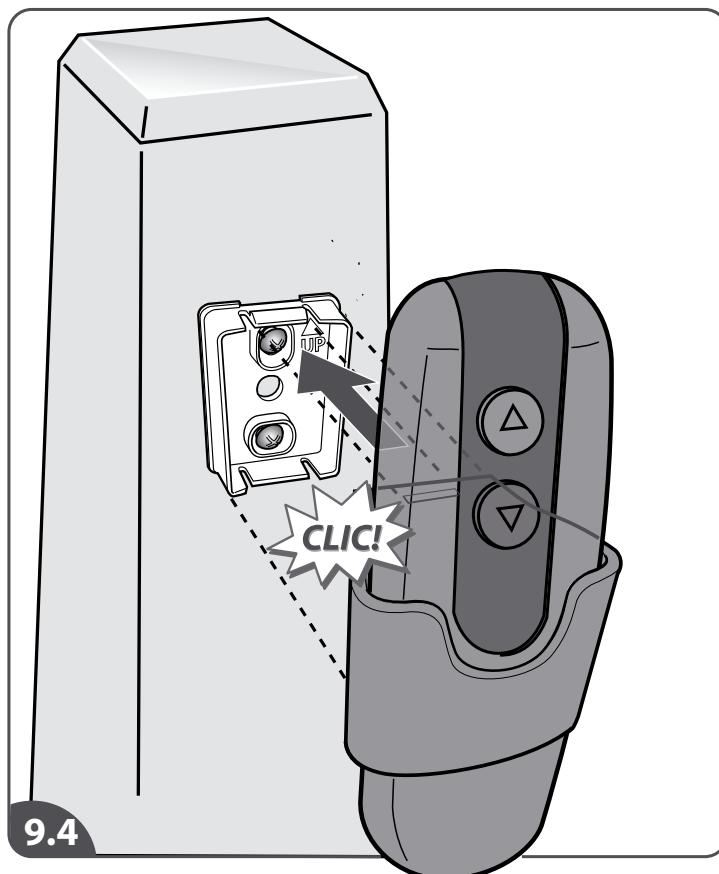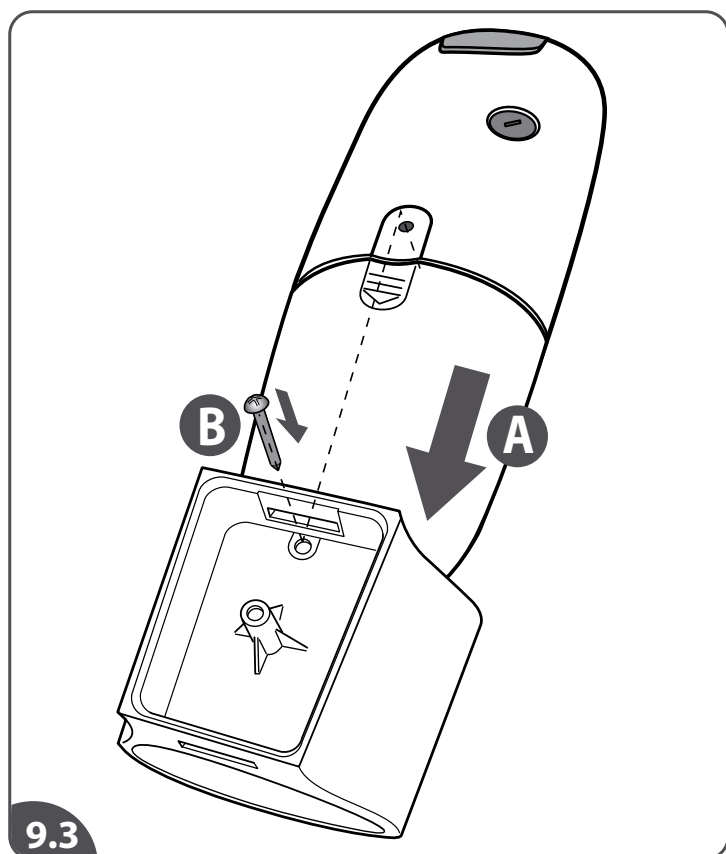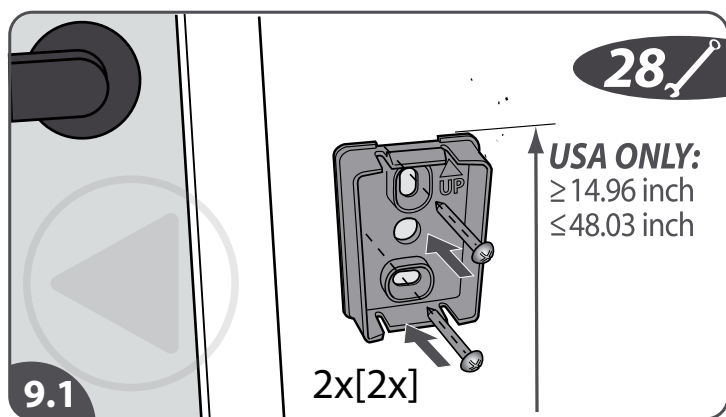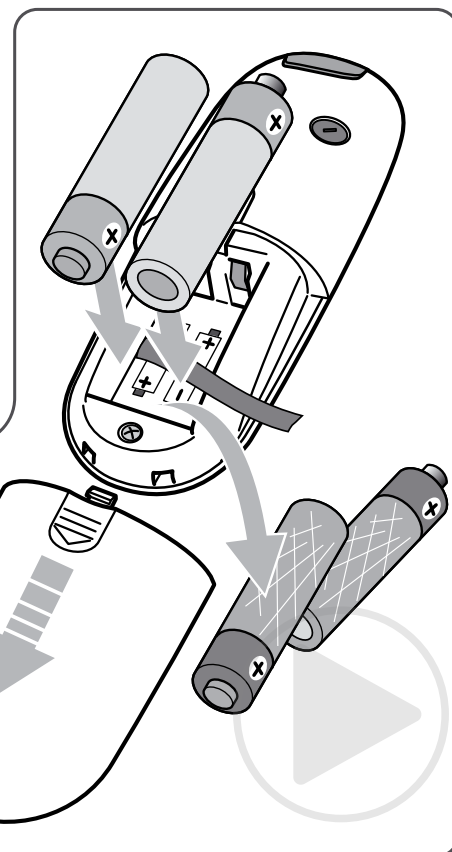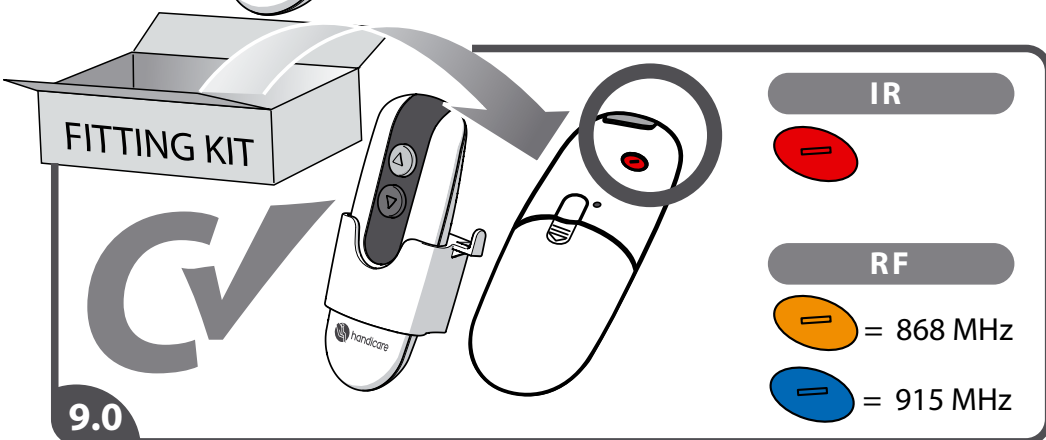[9.0] Check the amount and type of remote controls. The cap on the back of the remote control indicates the type: **Red=IR, Orange=RF**{868MHz}, **Blue=RF**{915MHz}. The supplied remote controls should work out of the box. If not or in case of multiple lift, they can be reprogrammed (see manual, page 70).

[9.1] Measure the position for the remote control bracket. Mount the bracket of the remote control with two screws on a solid background. Mount a remote control at the top and bottom of the stairs.

[9.2] Insert the batteries into the remote control.

[9.3] Mount the remote control into the holster{A} and fix it with a screw{B}.

EN [9.4] Click the holster/remote control over the wall plate. Leave the remote control in the holster so that it does not get lost.

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
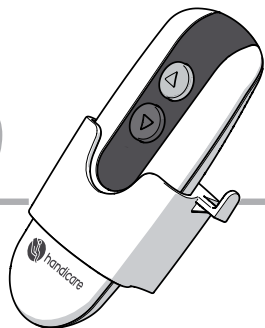
CHECK LINKING


HINGED RAIL ONLY

- = Lift moves **UP**
- = Lift moves **DOWN**
- = Lift moves to **PARKING** position


9.5

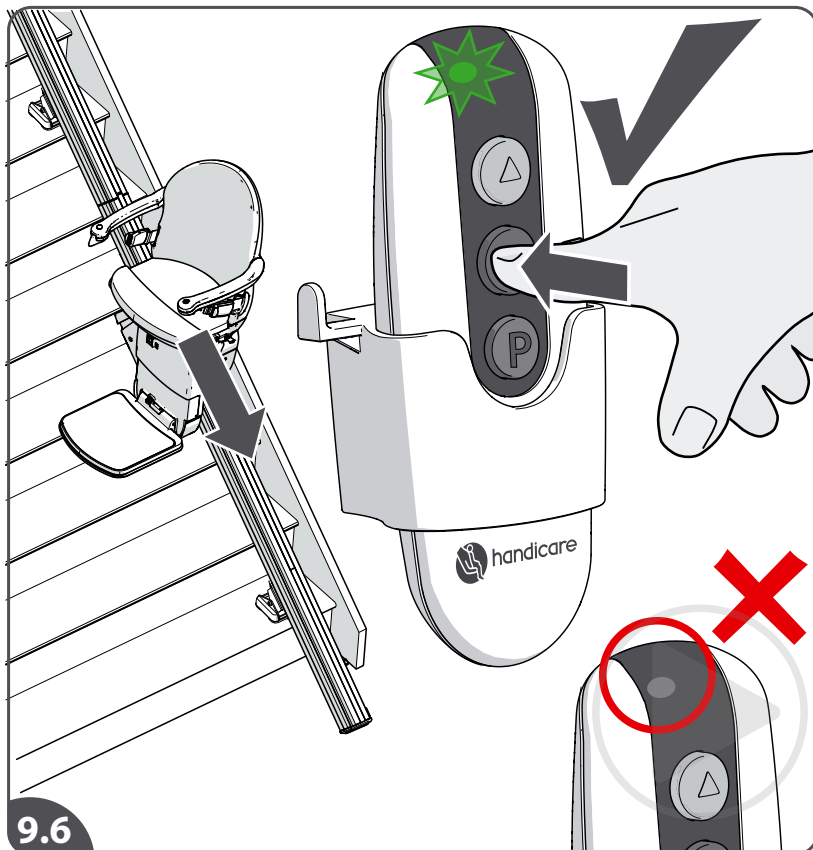9.6

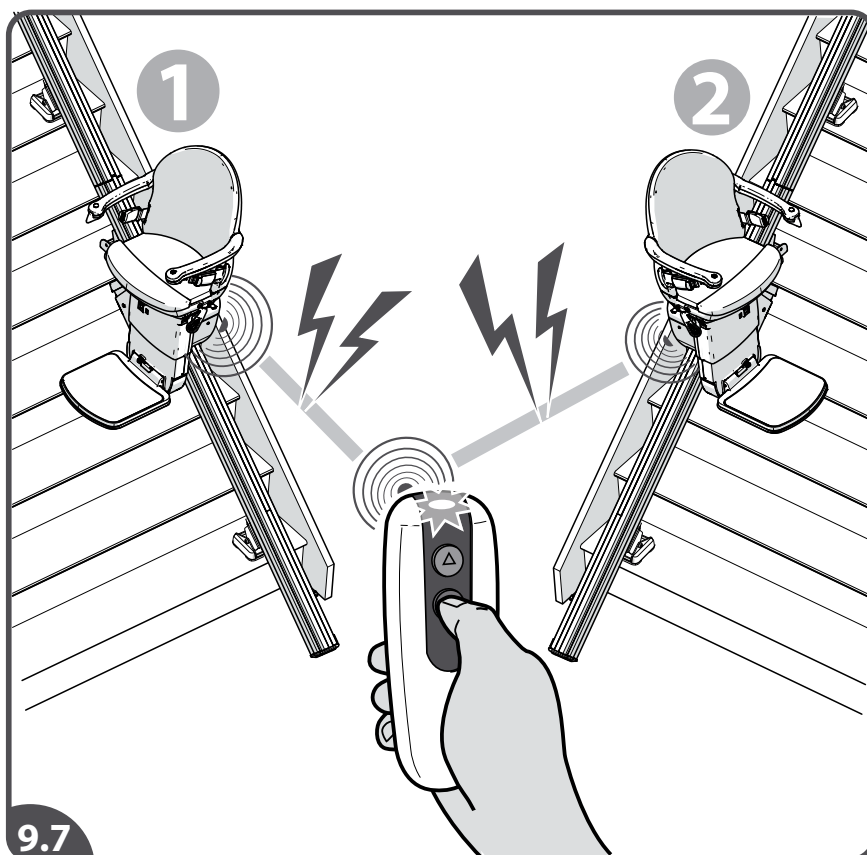9.7

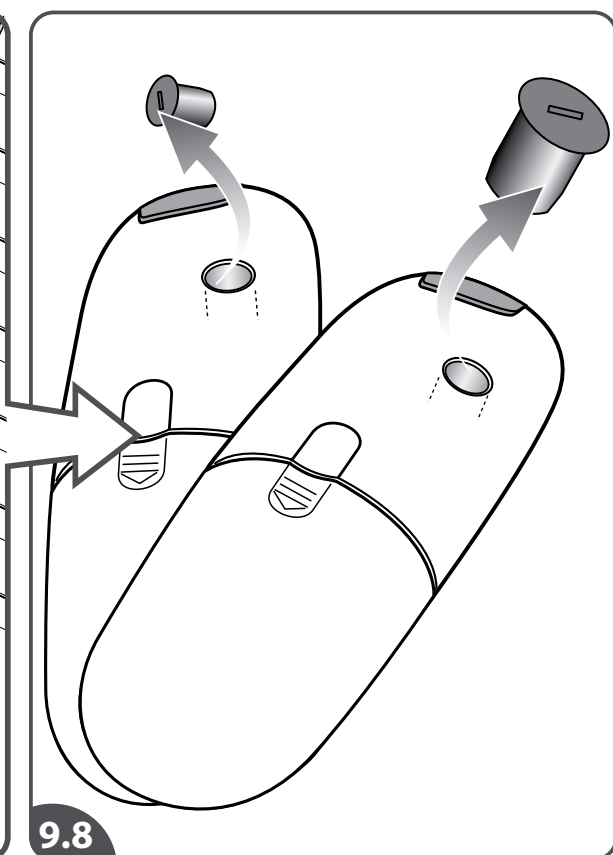9.8

[9.5]Function of the buttons: green = **UP**, red = **DOWN**, blue = **PARK** (only hinged rail).

[9.6]Press and hold the button (only press one button at a time) of the remote control to check if the remote control is linked with the lift. The remote control will show a flashing green light (searching). Once this signal has connected with your lift the green light will stop flashing and stay lit while your lift is moving.

[9.7]In case there is interference in the radio signal, change the channel of the radio signal of the lift.

Note: another lift nearby can interfere with the signal .

[9.8]**IF NEEDED FOLLOW STEPS 9.8 – 9.15** >To change the channel of the radio signal it is necessary to change the dip switches in the remote controls. In case of **RF** also the dip switches in the unit need to be changed. Remove the caps on the backside of the remote controls.

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





10

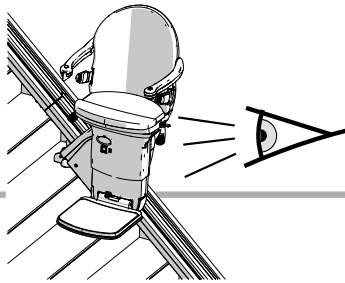CHECK MEASUREMENTS

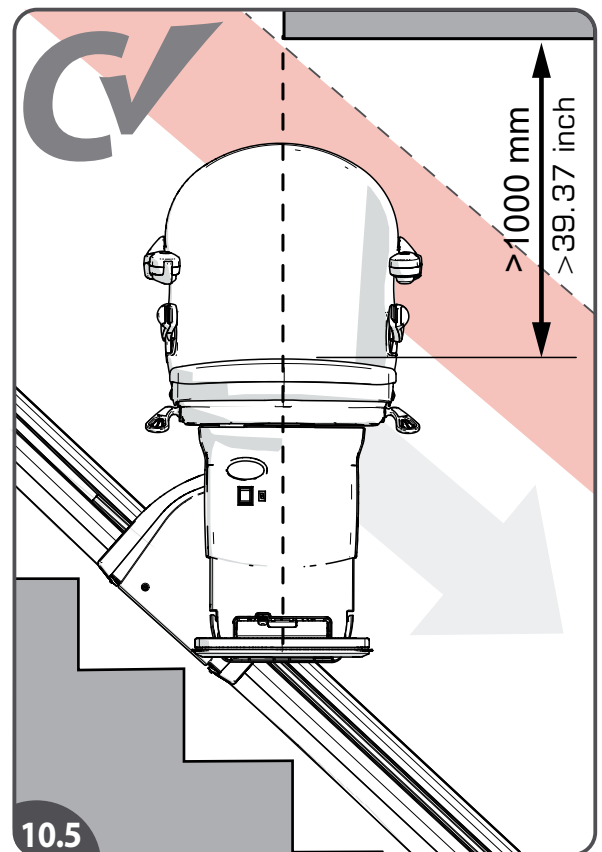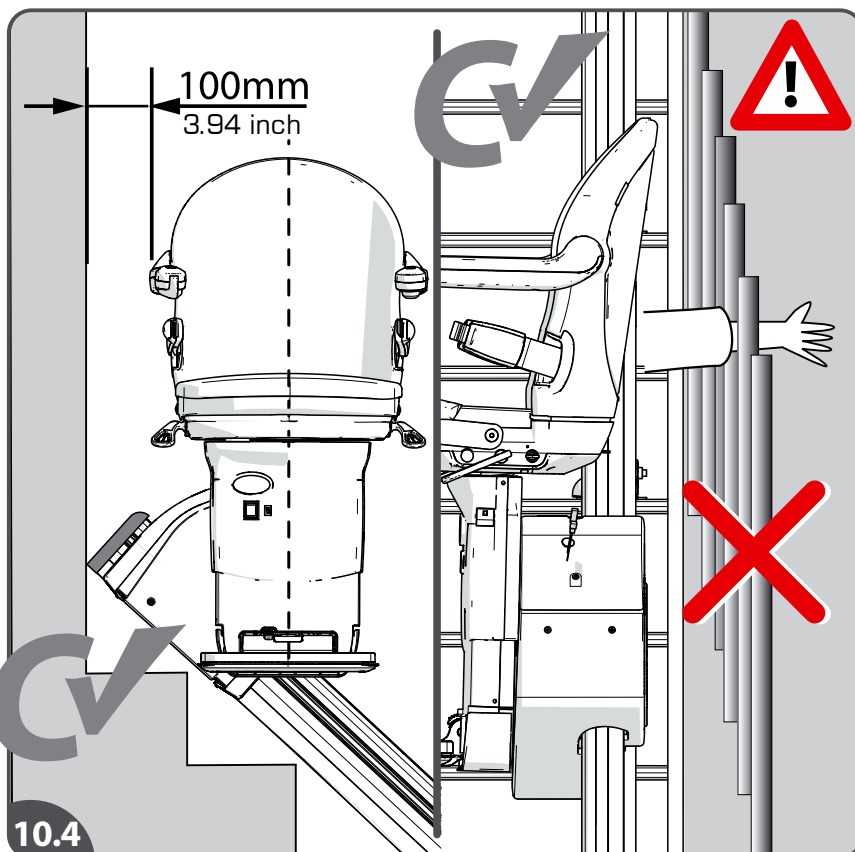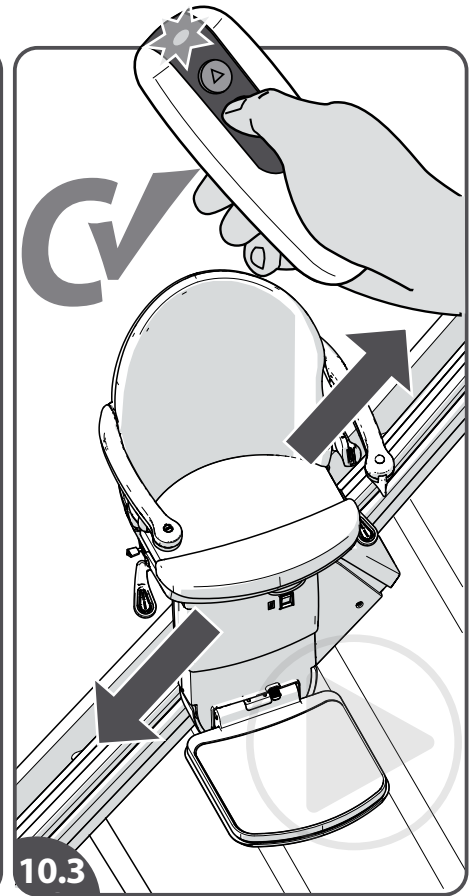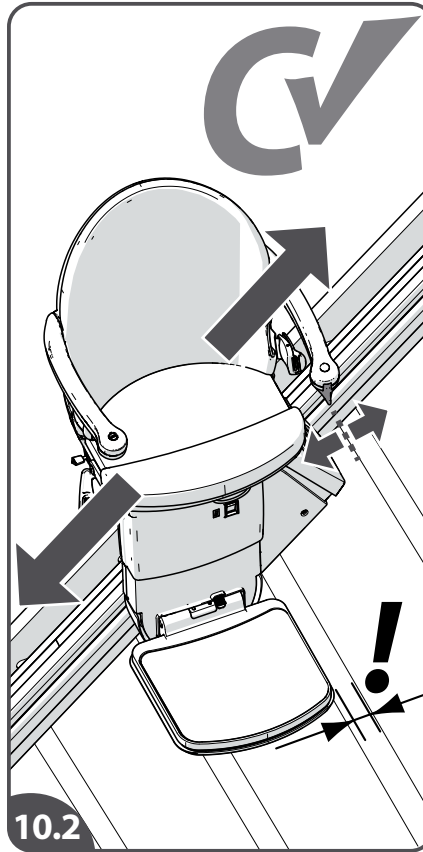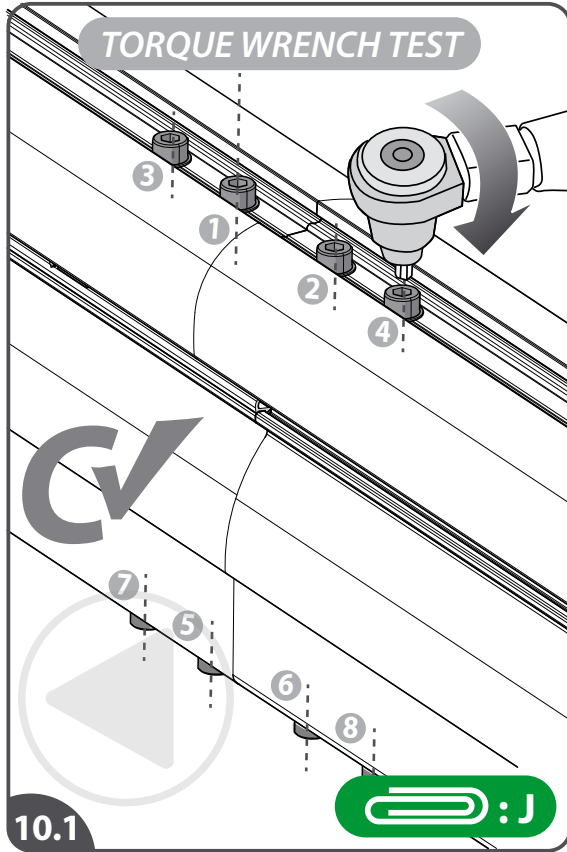- [10.1] Check the torque of the load-carrying fasteners. See attachment I how to do this.
- [10.2] Move the joystick up and down. Check whether the stairlift responds correctly.
- [10.3] Operate the stairlift with the remote control. Check whether the lift responds correctly.
- [10.4] Make sure that there is nothing behind the stairlift that body parts could get trapped in or on.
- [10.5] Check the position of the bulkhead.

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




IT

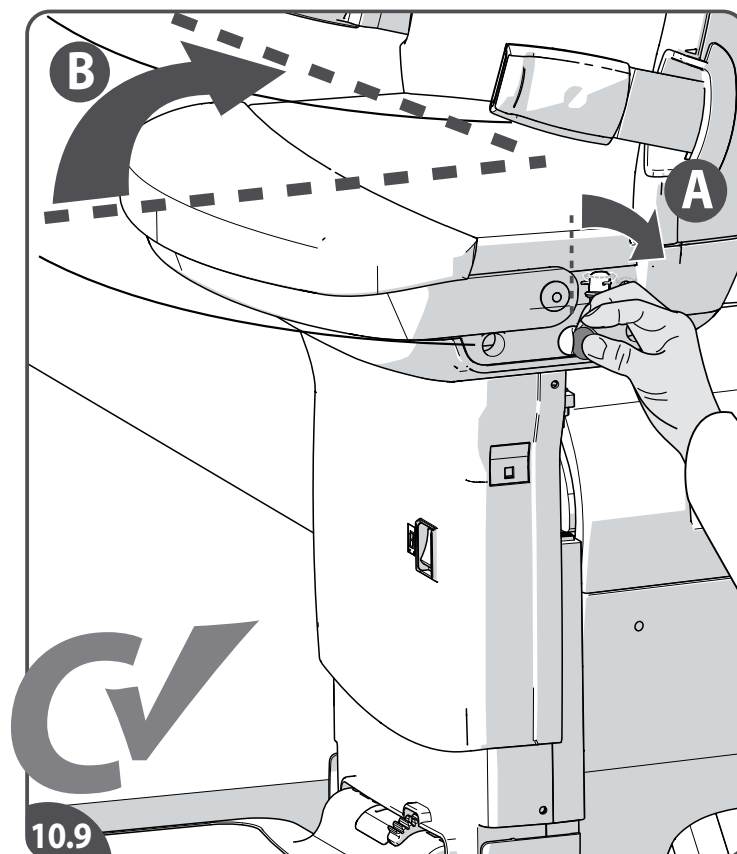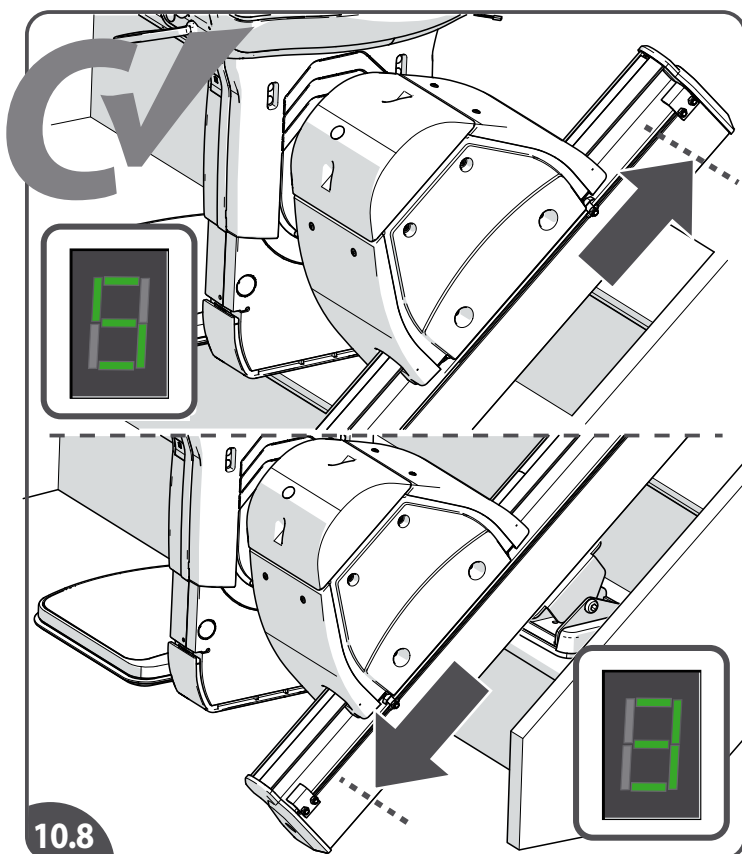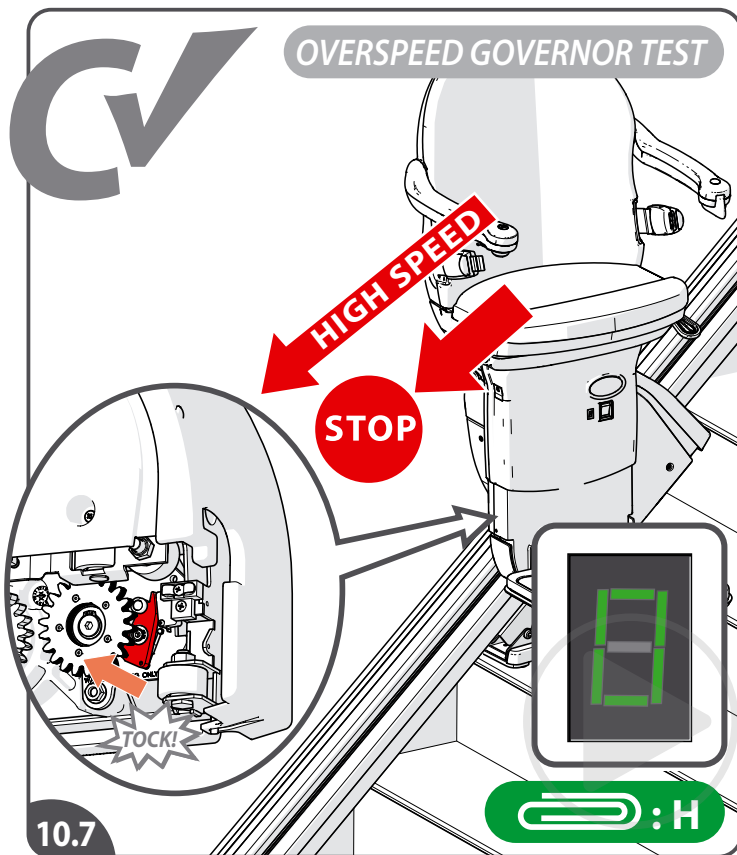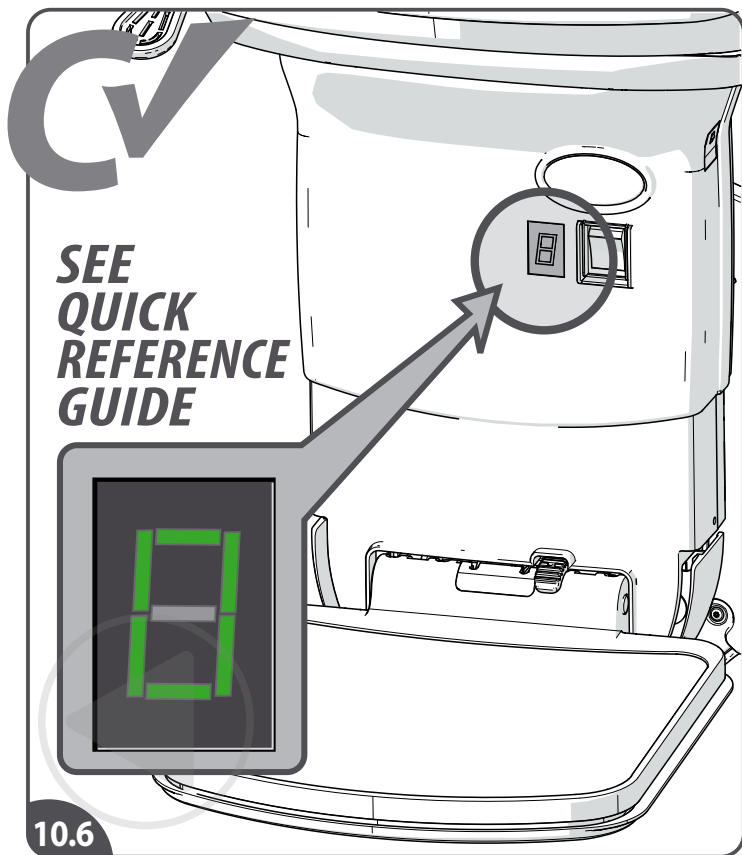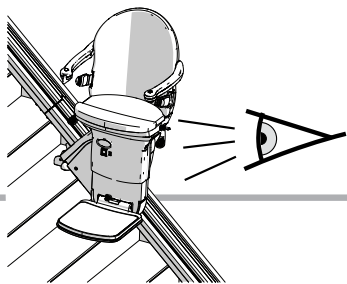- [10.6] Check the diagnostic display; make sure there are no fault messages. See the user manual for the error codes.
- [10.7] Check the OverSpeed Governor(OSG) safety (see attachment H).
- [10.8] Check that the stairlift makes a correct landing at the top and bottom of the rail.
- [10.9] Check the swivel release: turn the release knob and swivel the seat manually.

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



10

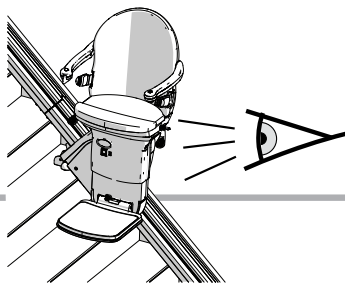CHECK.
OPTION FUNCTIONS

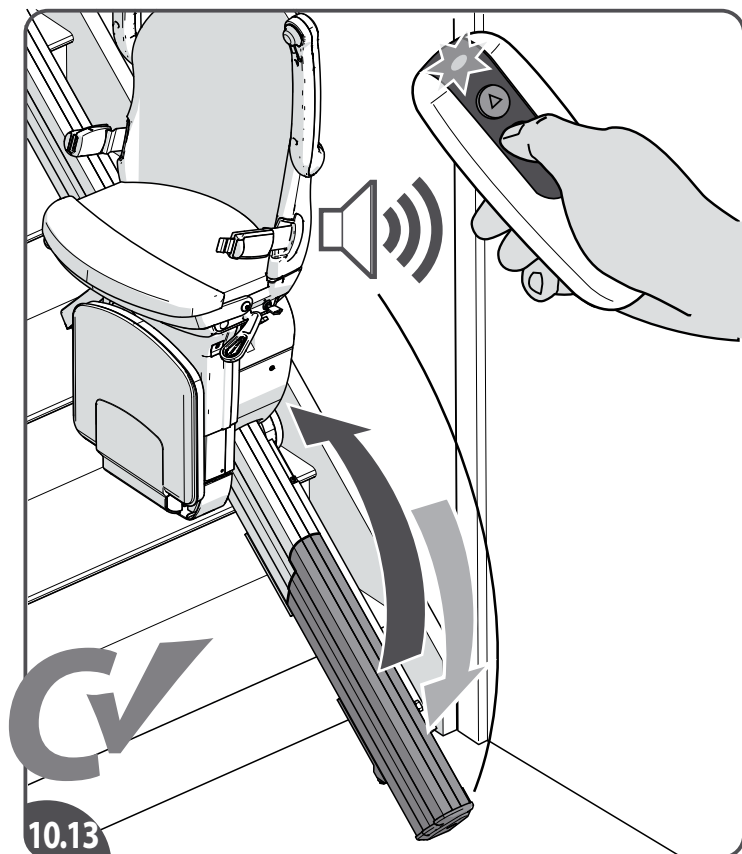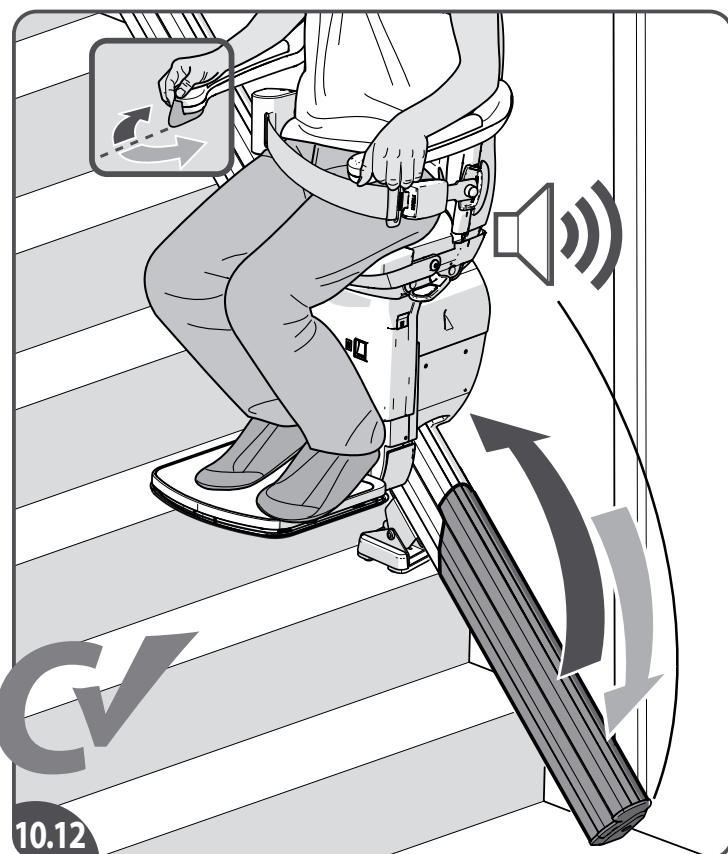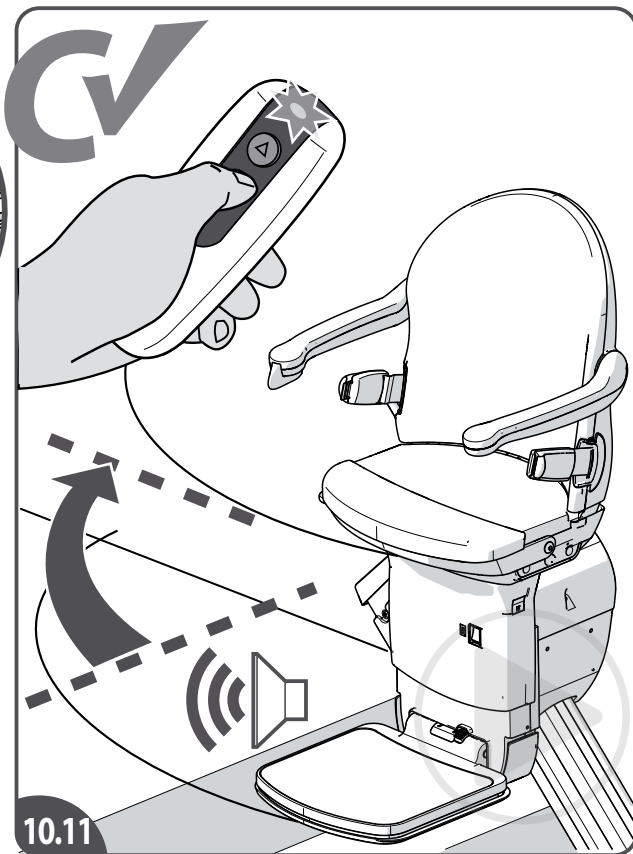- [10.10] Check the powered swivel function. A lift (with an automatic swivel configuration) is able to swivel when it reaches the charge point (33 volt). The lift is not able to swivel when the keyswitch is off.
- [10.11] Check the powered swivel function again while using the remote control.
- [10.12] Check the powered hinge track function.
- [10.13] Check the powered hinge track function again while using the remote control.

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




10

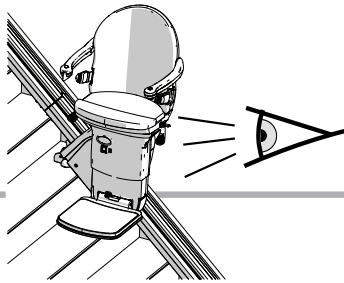CHECK.
CHARGE POINTS

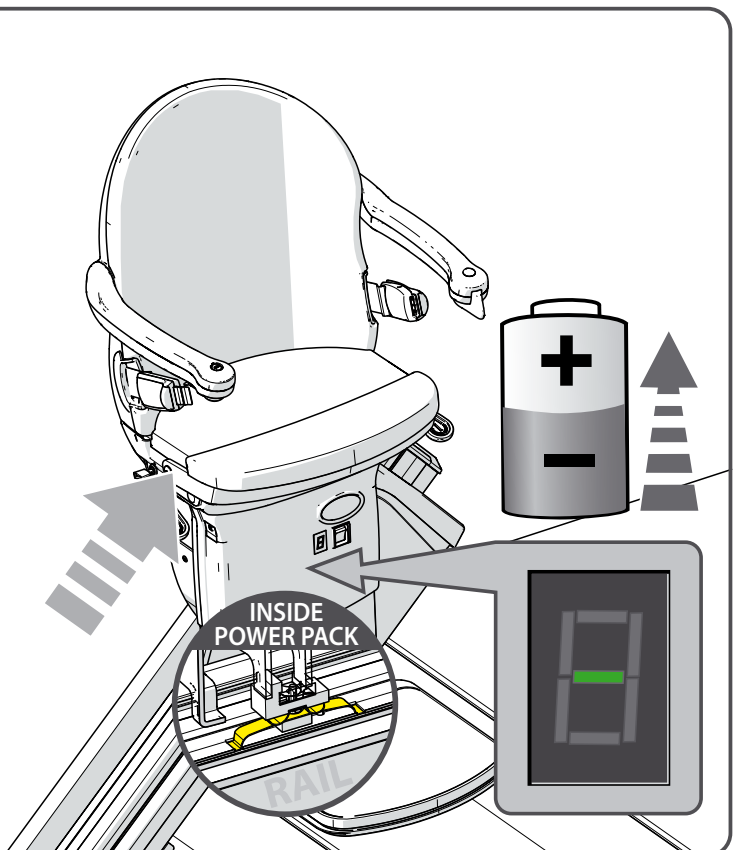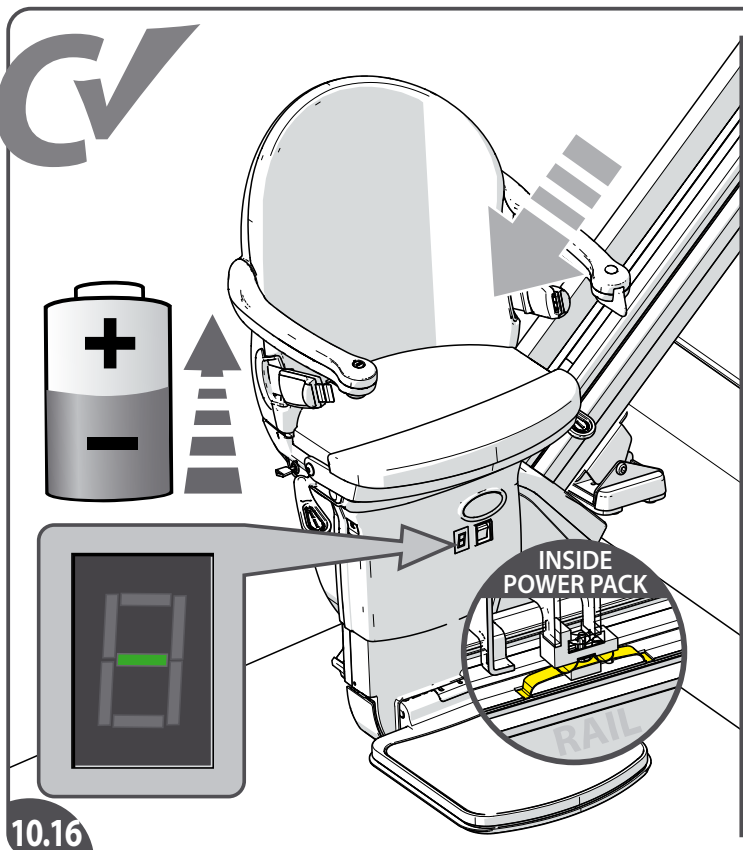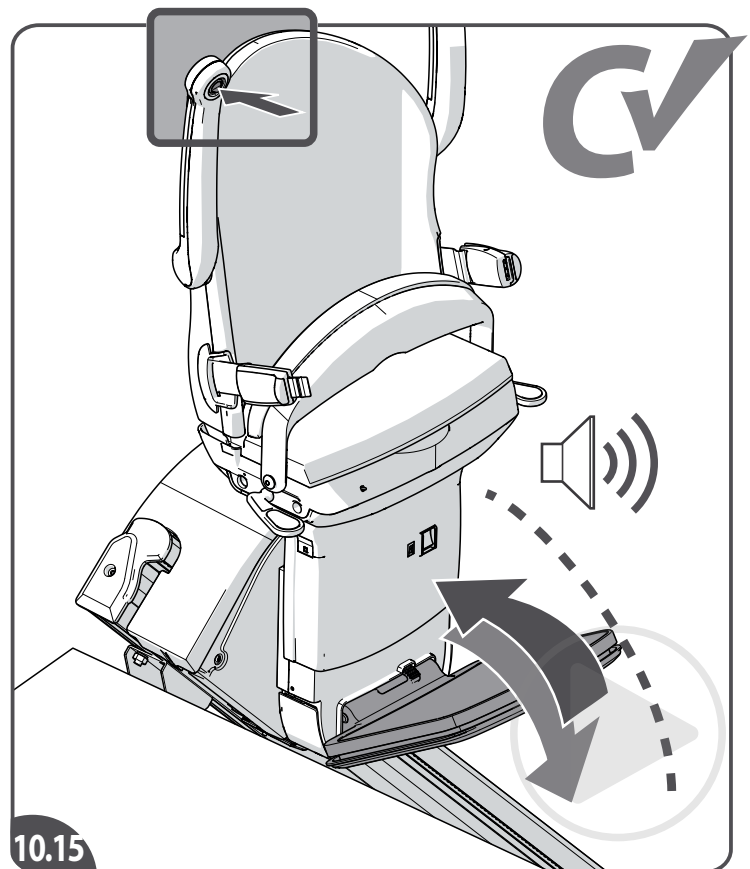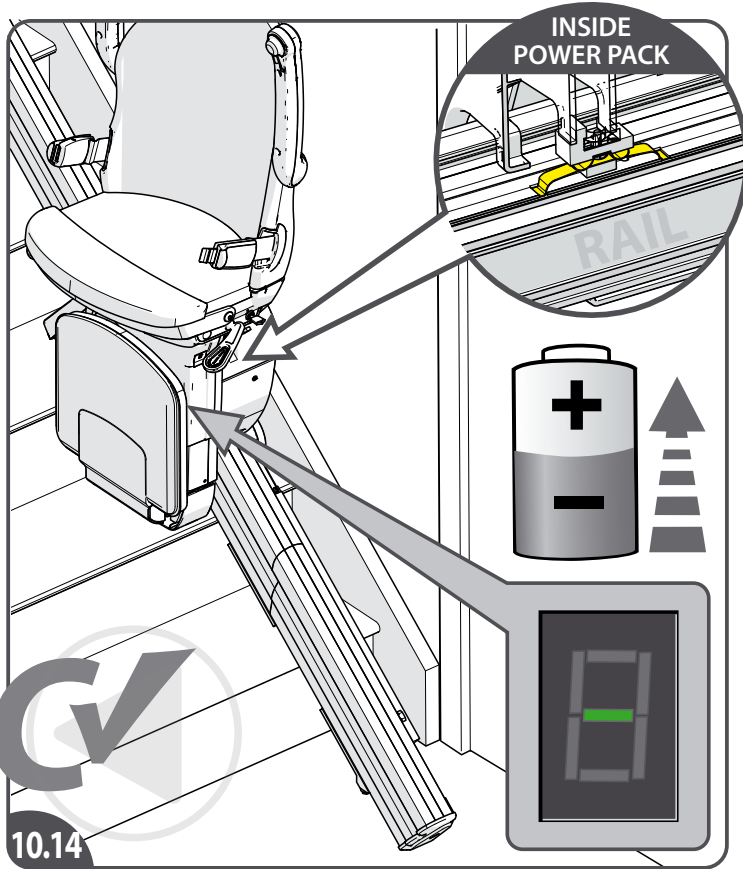[10.14] If the lift has an intermediate charge point (hinged rail), test the park function. Ensure that the diagnostic display shows a “-” symbol at the charging point.

[10.15] Check the powered footrest function. The footrest is not able to move when the keyswitch is off.

[10.16] Check that the lift charges correctly: drive the lift to the bottom charge contacts and ensure that “-” shows in the diagnostic display. Drive the lift to the top charge contacts and ensure that “-” shows in the diagnostic display.

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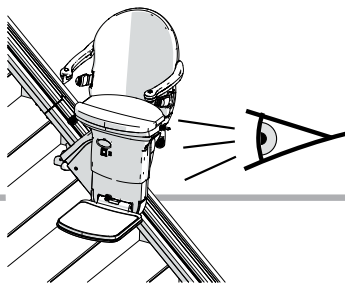
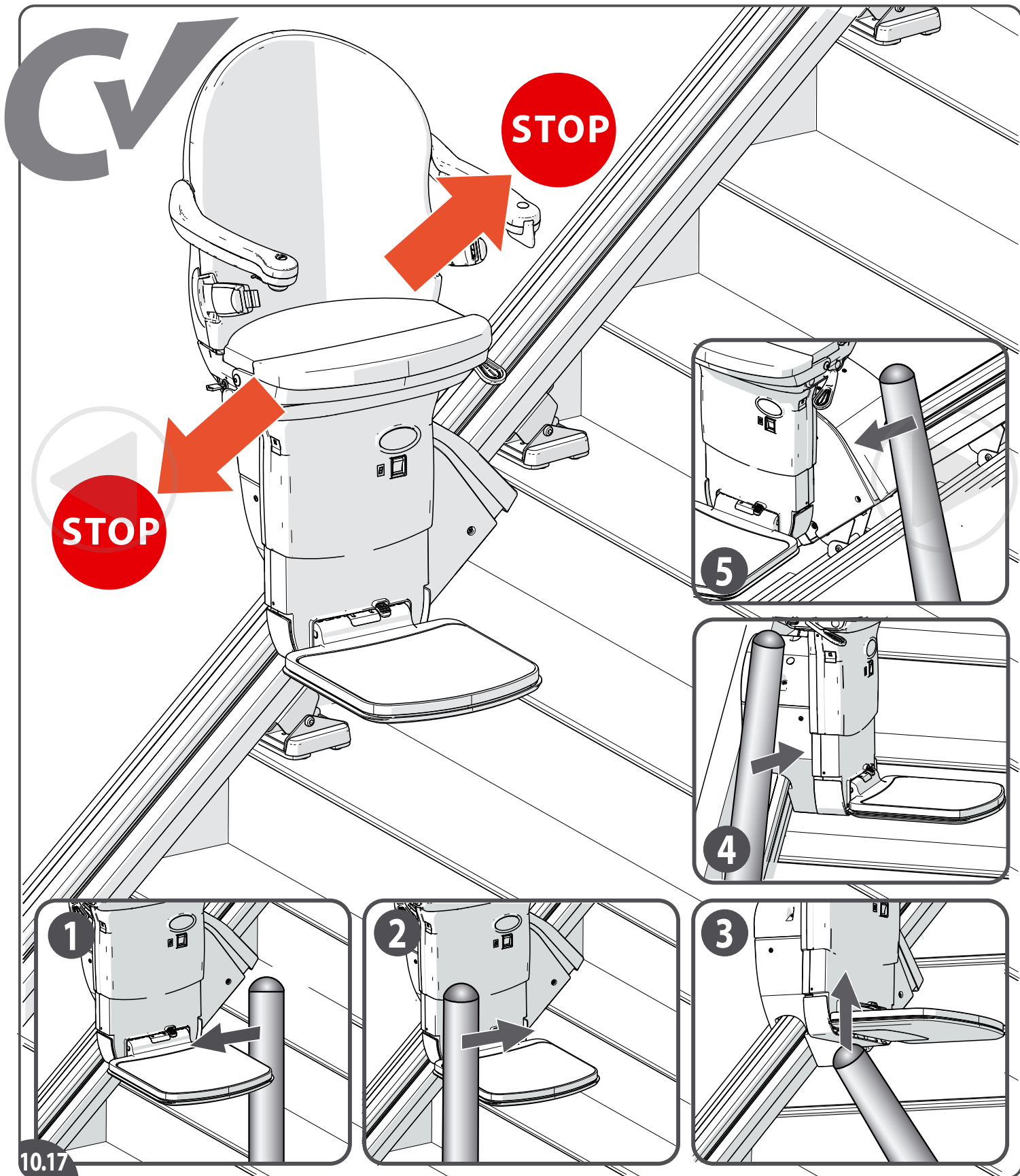
10.17

[10.17] Check the safety stops.

Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to the requirements of the latest adopted ASME A18.1 code.

Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.

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




10

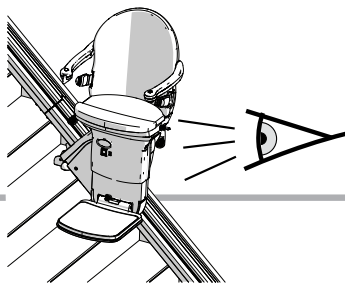CHECK ARMRESTS

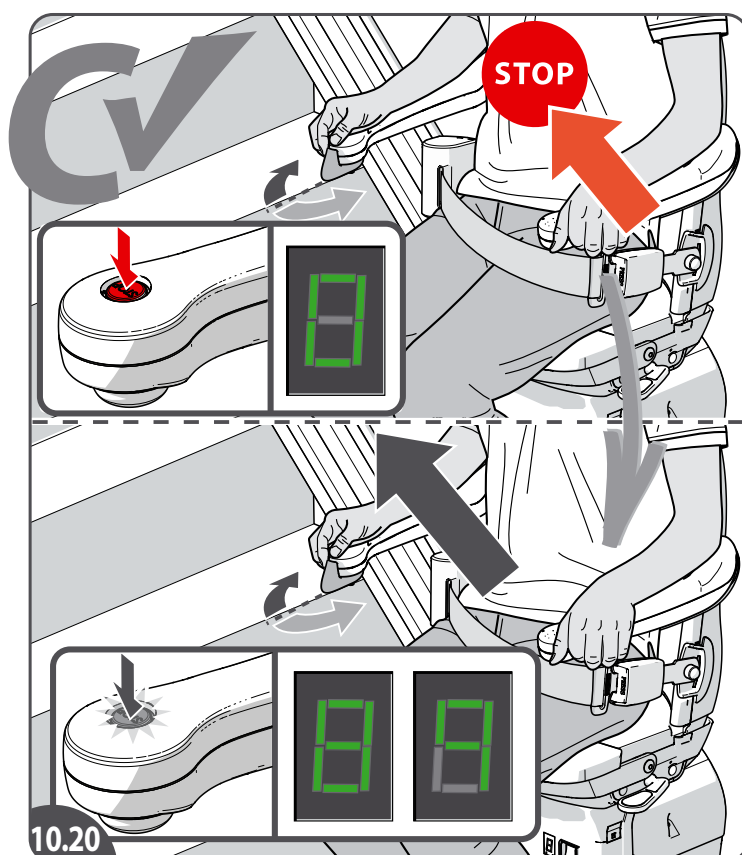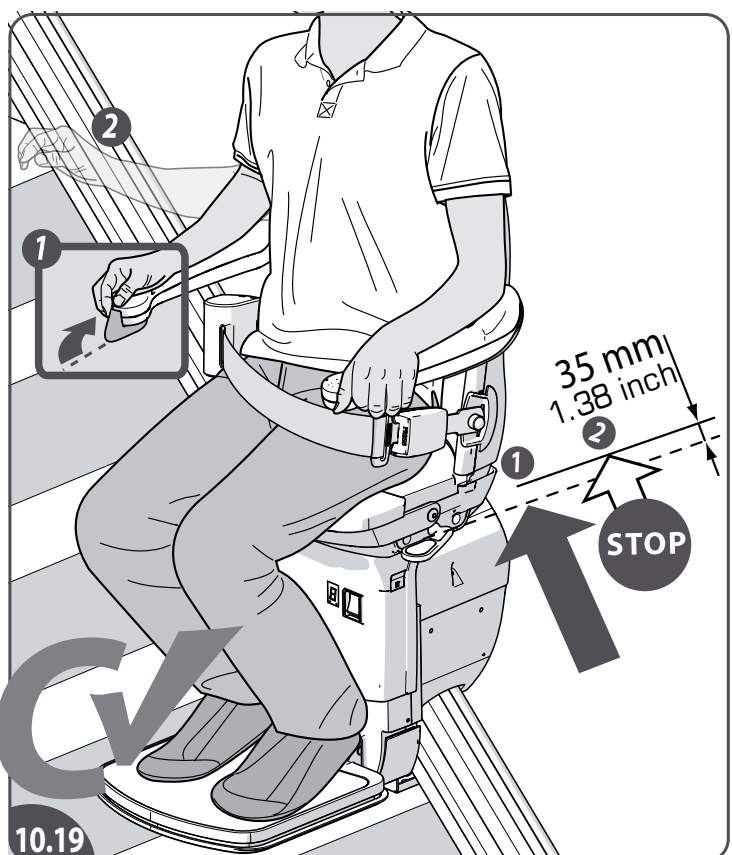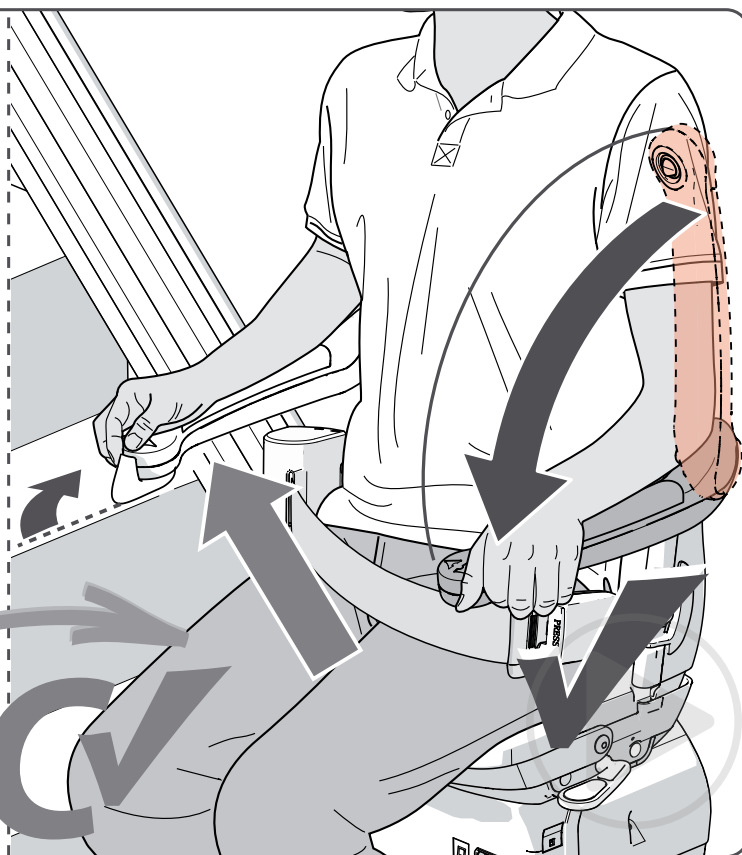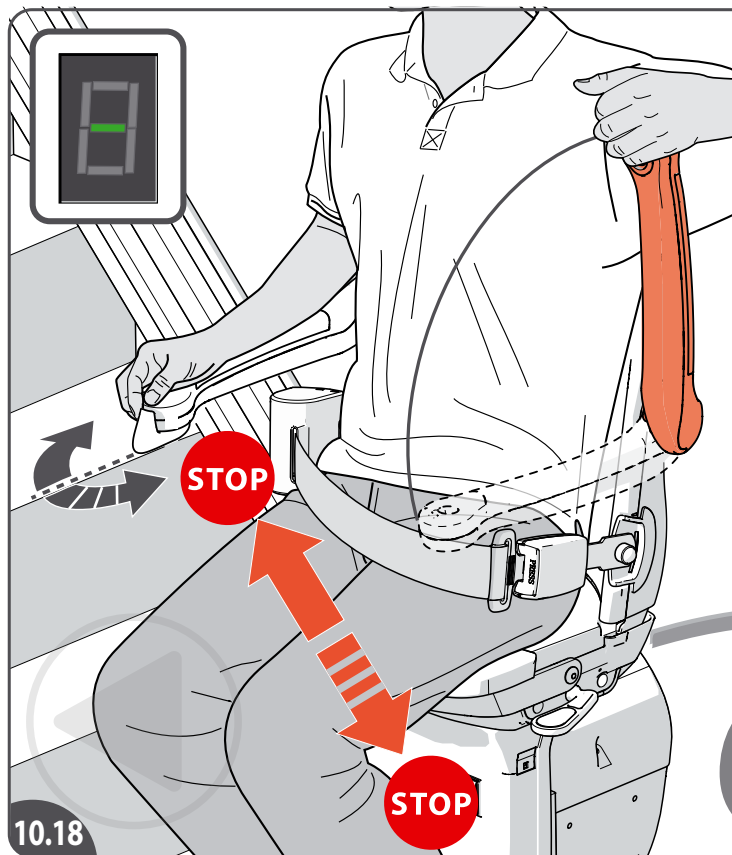[10.18]Fold the armrest up. The stairlift must stop immediately. Put the armrest back down and check that the stairlift operates normally again.

[10.19]Release the joystick toggle switch. The stairlift should stop within 35 mm (1.38 inch).

[10.20]Press the stop button. The lift must stop immediately. Press the stop button again; check that the stairlift operates normally.

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




IT

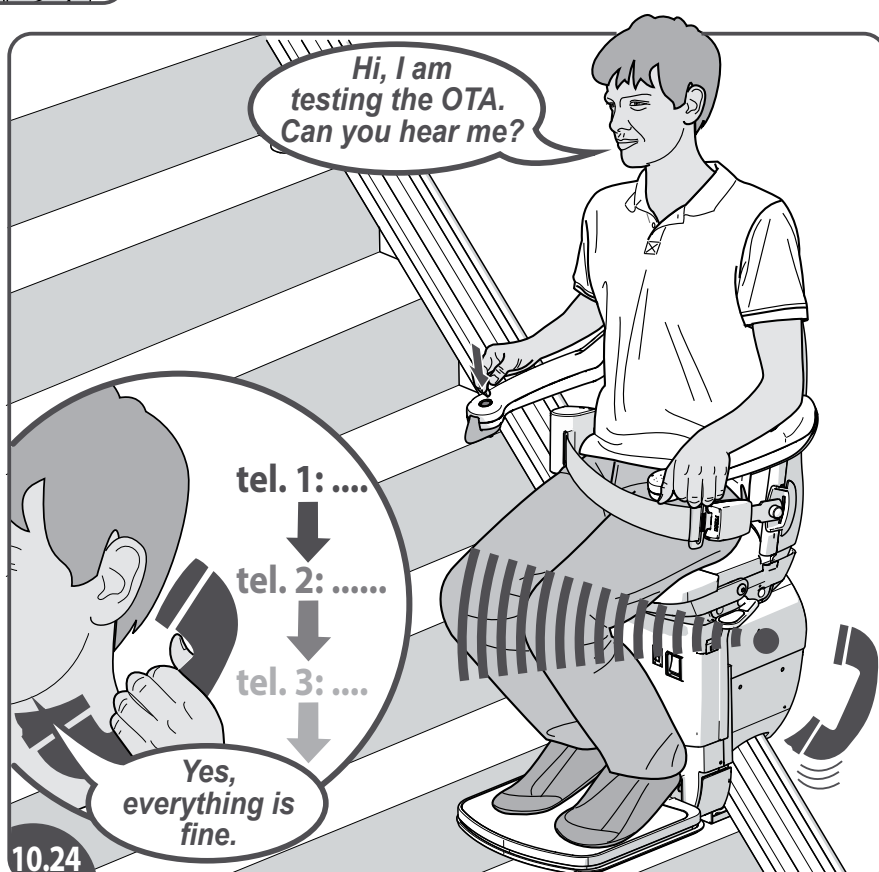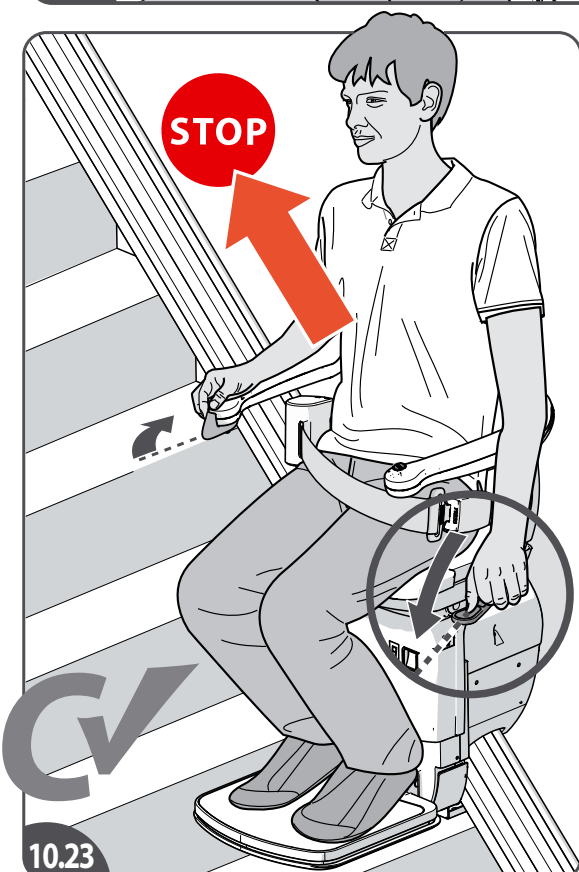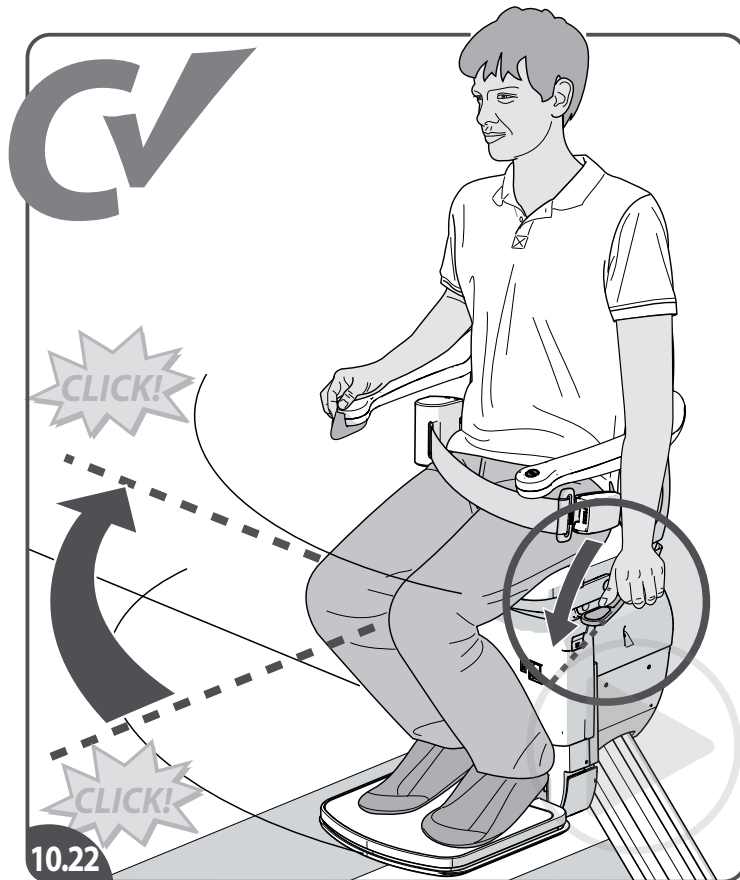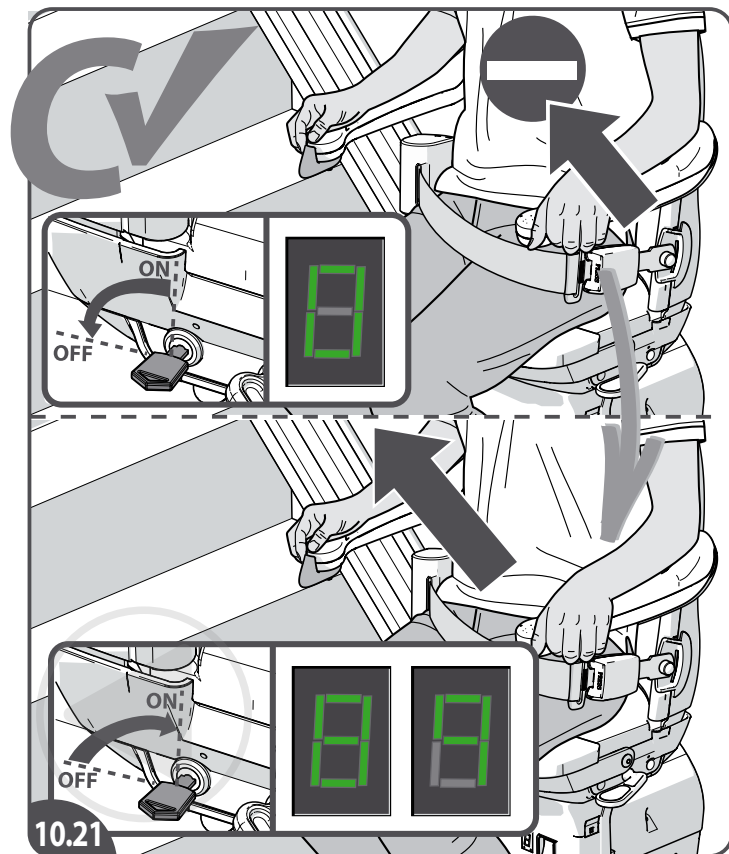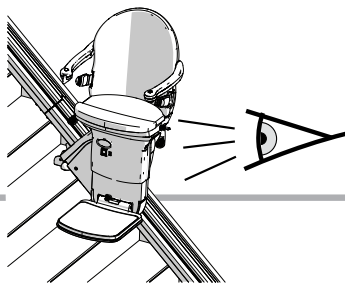- [10.21] Turn the key to the OFF position and ensure that the stairlift will not move. The diagnostic display should show a "0". Turn the key to the ON position. The stairlift should be able to move again.
- [10.22] Check the manual swivel function and that the end stops work satisfactorily (it clicks at the end).
- [10.23] While moving the stairlift, swivel the seat; the stairlift must stop immediately.
- [10.24] It is the responsibility of the Dealer to ensure that the Integrated TWVC System is checked and test calls are performed at installation and normal stairlift service intervals. As it is the responsibility of the customer (end-user) to provide and manage the Sim Card used within the Integrated TWVC System.

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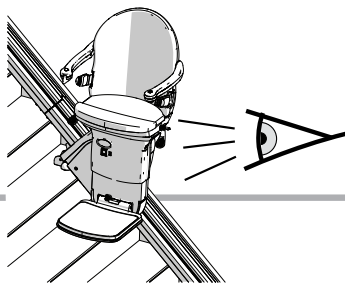
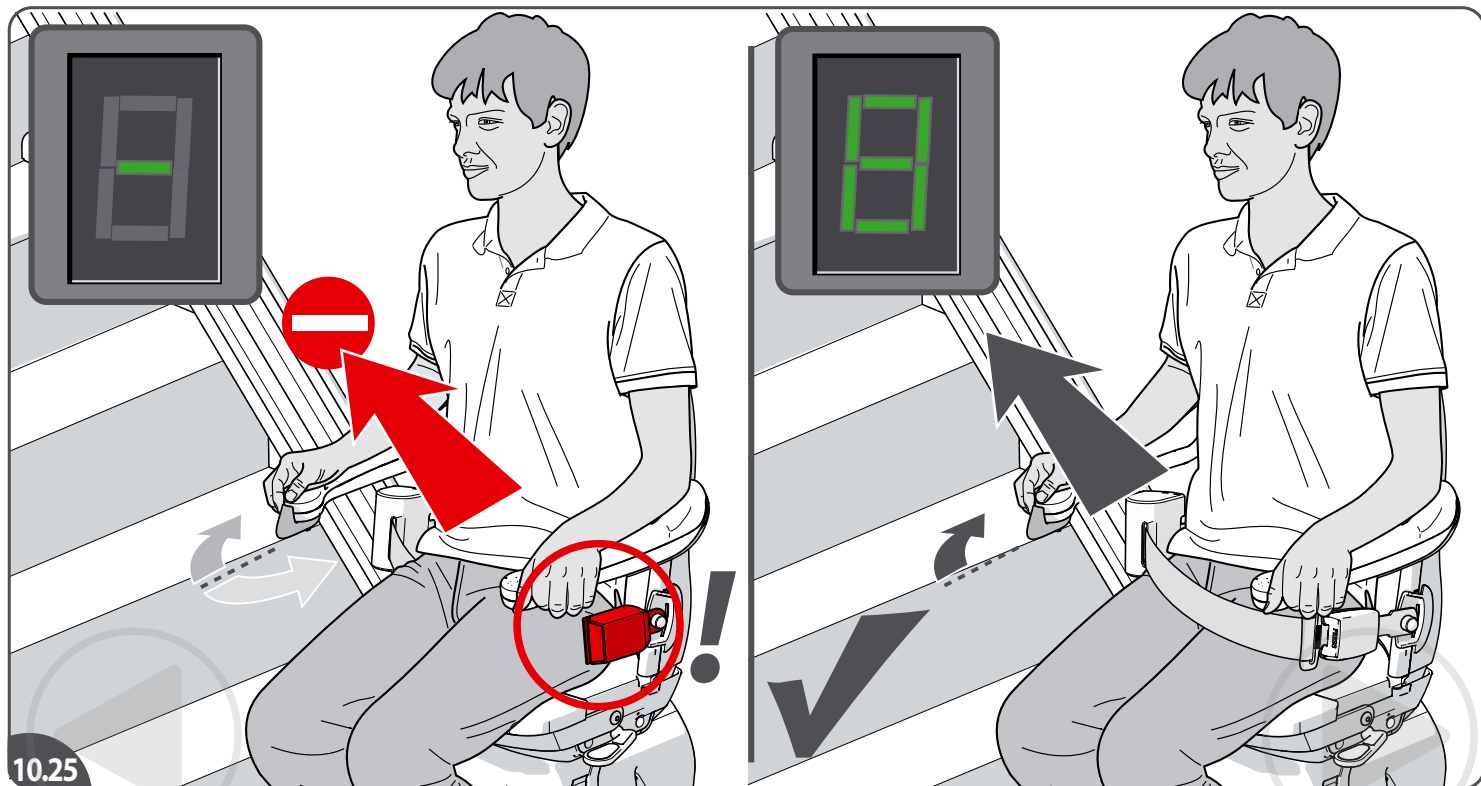
10.26

[10.25] Check the seat belt safety.

[10.26] Instruct the user by telling him/her all users operating instructions, who are explained in the user manual.

Note: Ensure that the site is cleaned and tidied before demonstrating the lift to the customer.

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IT

[10.27] A full ride is made along the complete travel, include hinged rail guide sections and swivelling of the seat where applicable; fill in the installation checklist. Corrective actions required for boxes that are ticked as "No".

[10.28/29/30] The declaration of conformity is legal after the declaration engineer and the declaration instructed person have been signed off by the engineer and the user. The declaration engineer and the declaration instructed person shall be saved by the dealer and made available to Handicare on request.

EN *Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to the requirements of the latest adopted ASME A18.1 code. Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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NL *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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ES *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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DE *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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
FR *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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
IT *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

10

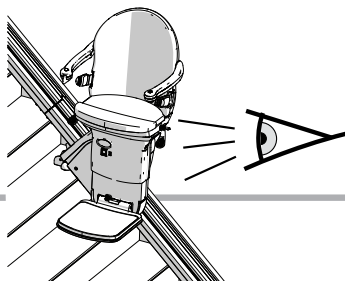CHECK PLACING LABEL

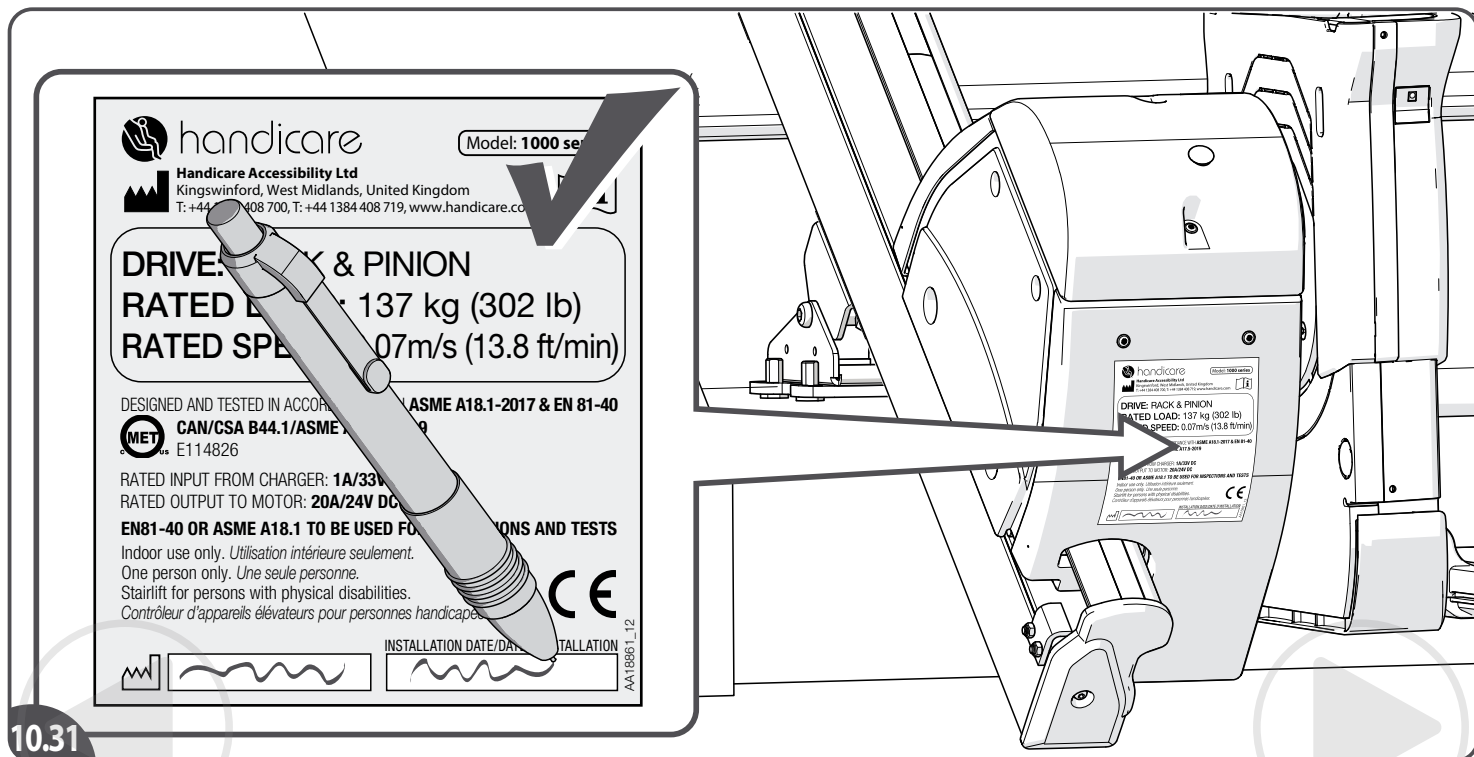10.31

[10.31] Fill in the installation dates on the label. The label is located on the side of the power pack cover.

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APPENDICE

RAIL CUTTING (mm)

EN [1] Measure the distance between the landing floor and the ground floor (diagonal) = **L**

[2] Verify the related cut-out dimensions **C** and **D** in the table

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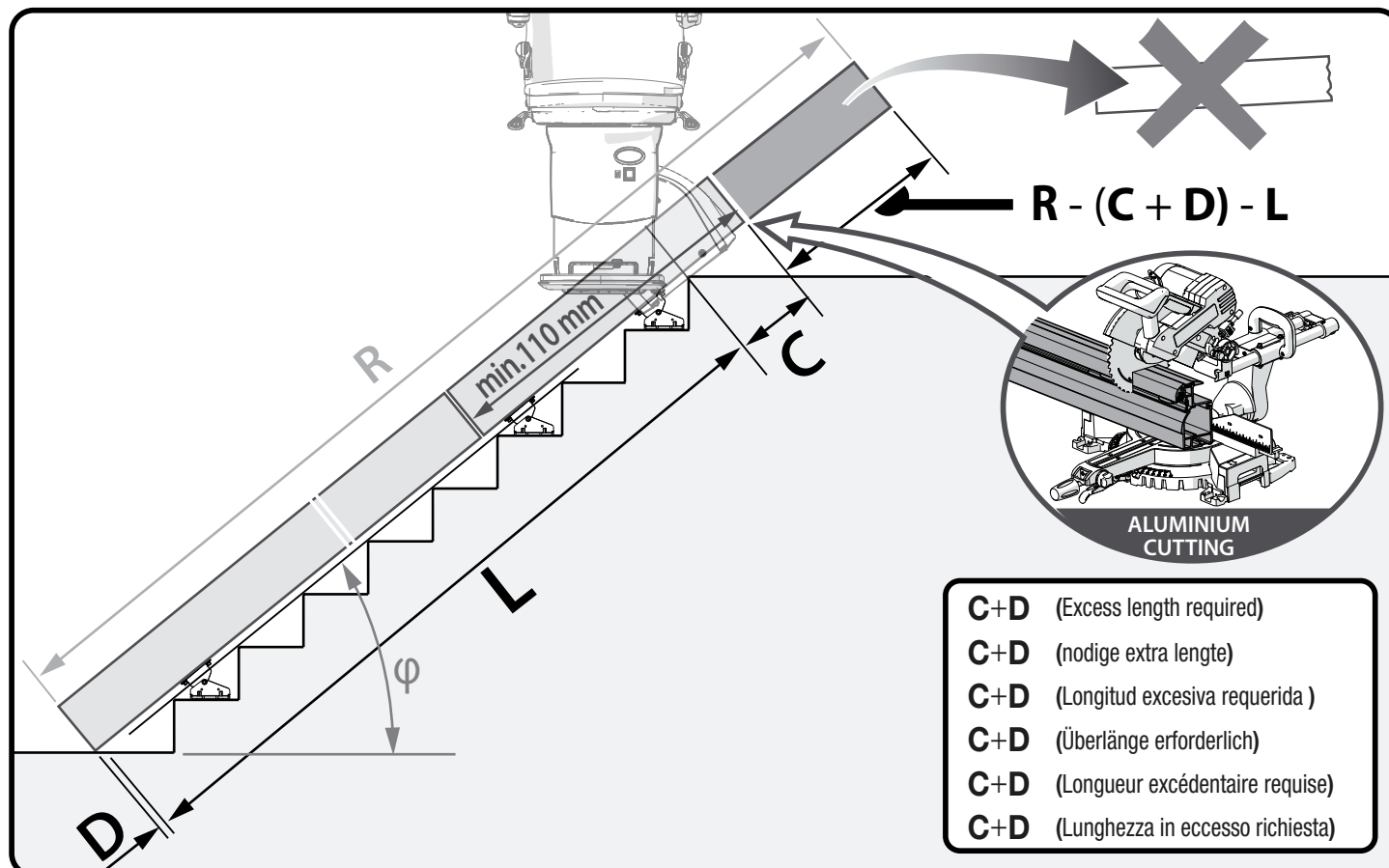
A **RAIL CUTTING (inch)**


[1] At First measure the distance between the landing floor and the ground floor (diagonal) = **L**

[2] Verify the related cut-out dimension **C** in the table


[3] Saw the rail part: **R - L - (C+D)**

Φ (Angle)	C+D (Excess length required)
(38°)	8.03
(39°)	8.27
(40°)	8.50
(41°)	8.70
(42°)	8.94
(43°)	9.13
(44°)	9.30
(45°)	9.48
(46°)	9.65
(47°)	9.92

Φ (Angle)	C+D (Excess length required)
(48°)	10.0
(49°)	10.2
(50°)	10.3
(51°)	10.4
(52°)	10.5
(53°)	10.6
(54°)	10.7
(55°)	10.8
(56°)	10.9




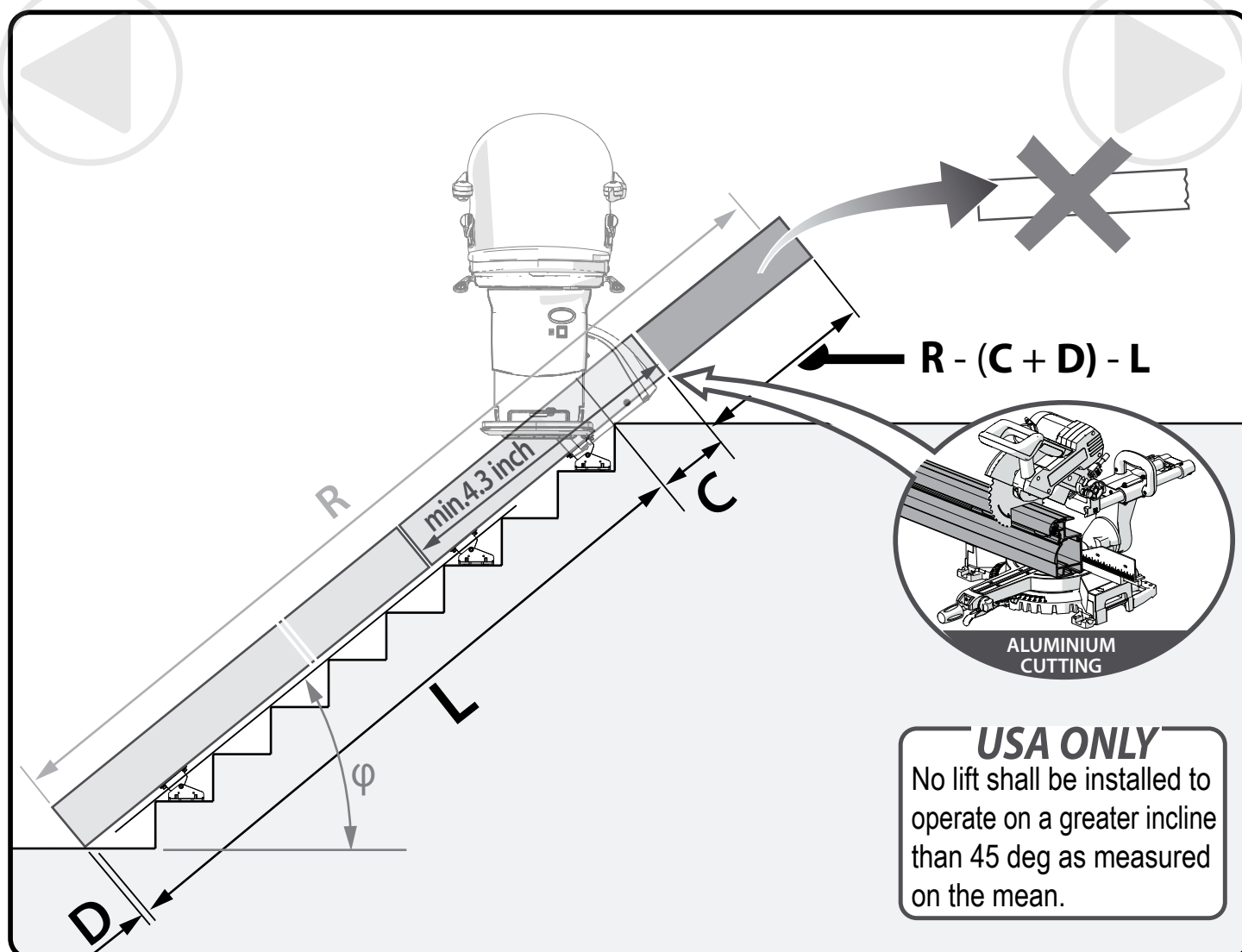
B **3 RAILPARTS (INTERMEDIATE KIT)**

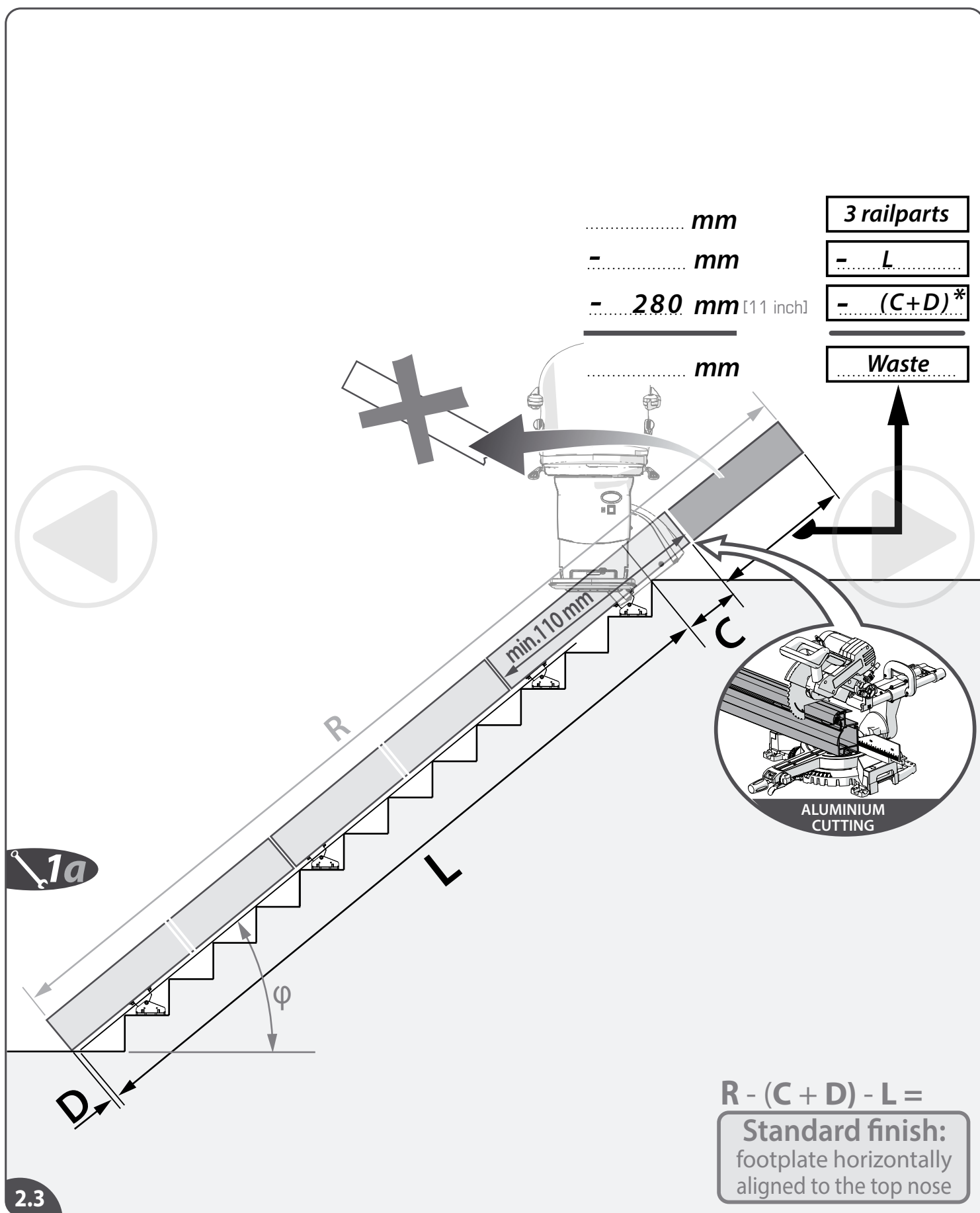2.3

C *EMERGENCY OPERATION > TRACK*

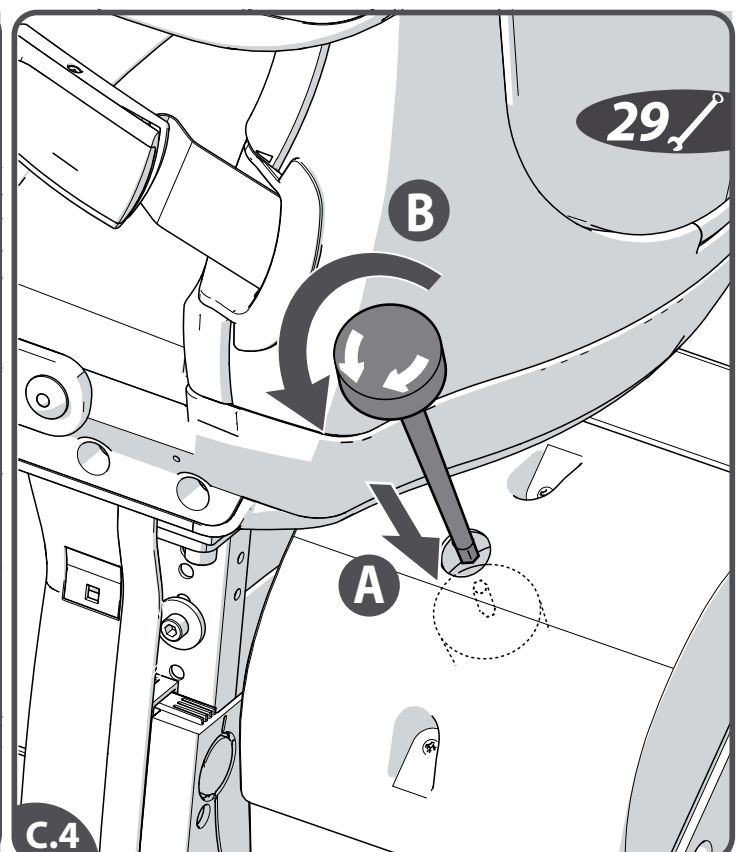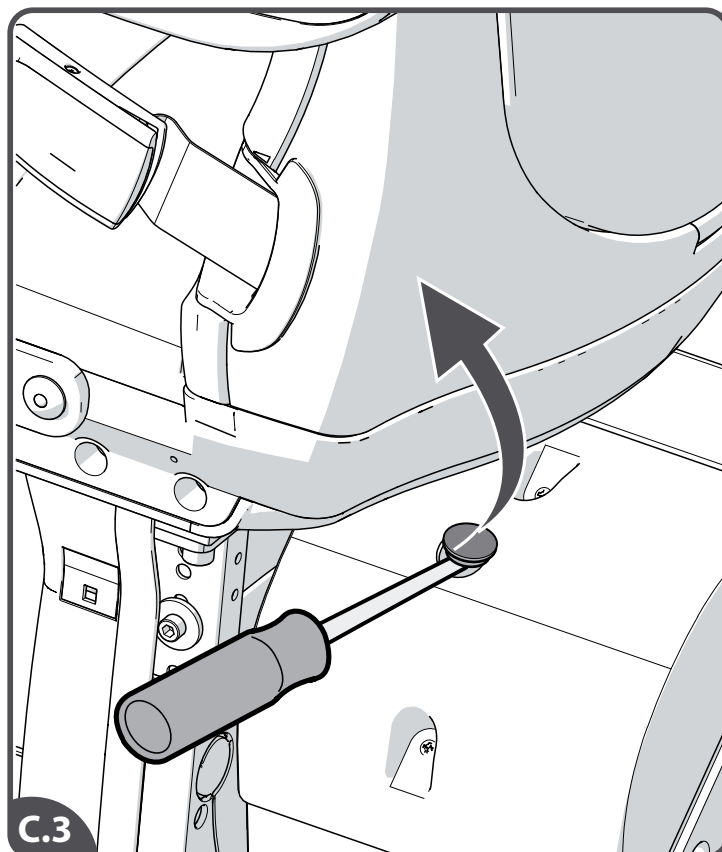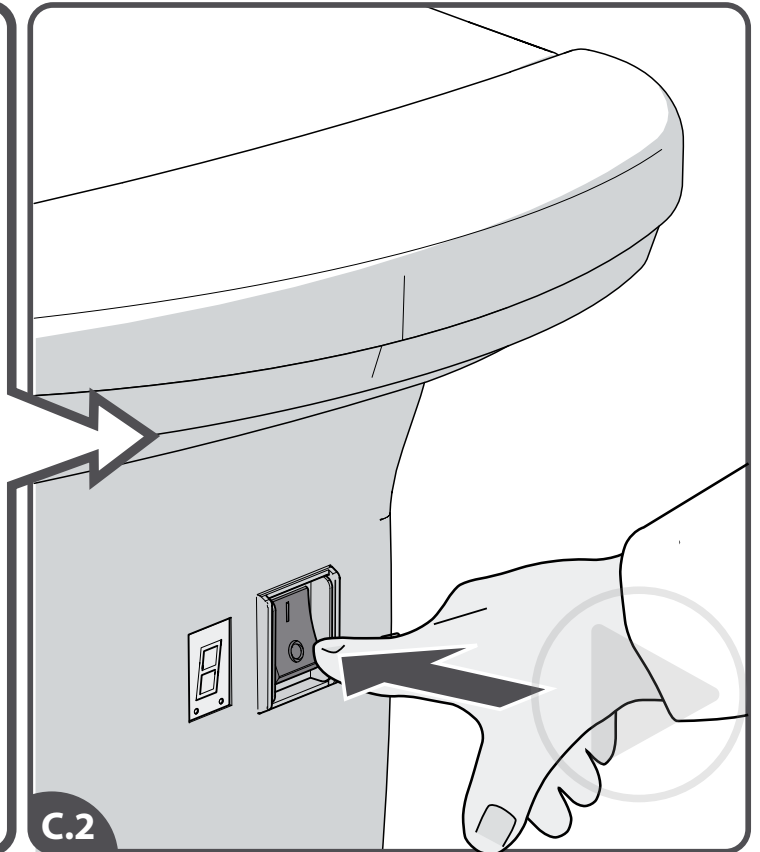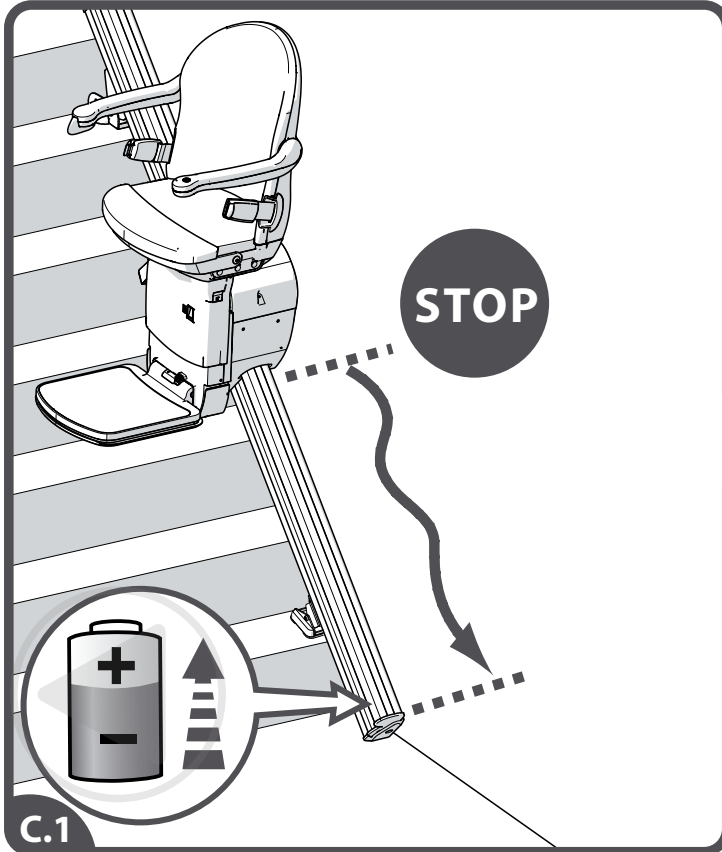- [C.1] When the stairlift stops unexpectedly and no longer moves any further, the hand winder can be used to move the lift to a charging point.
- [C.2] Turn off the lift.
- [C.3] Remove the cap from the unit with a screwdriver.
- [C.4] Turn the hand winder as indicated on the winder. Then wind the lift in the required direction.

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



IT

D EMERGENCY OPERATION > SEAT

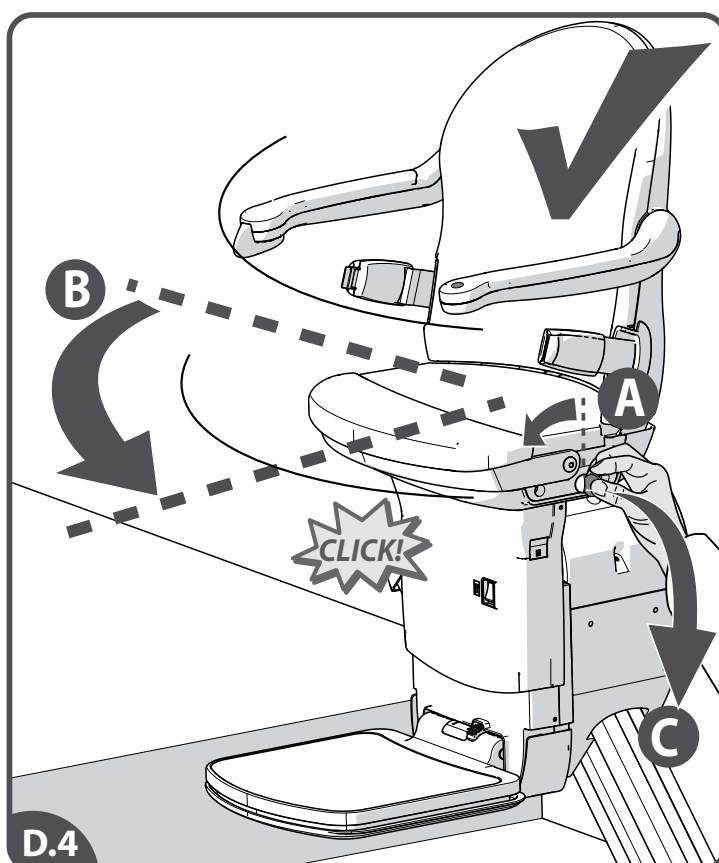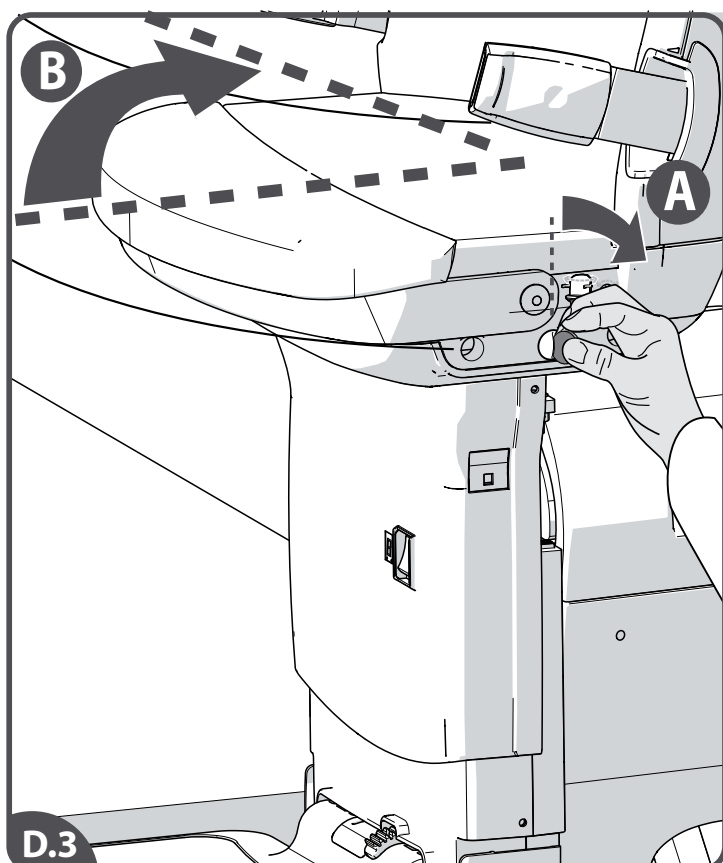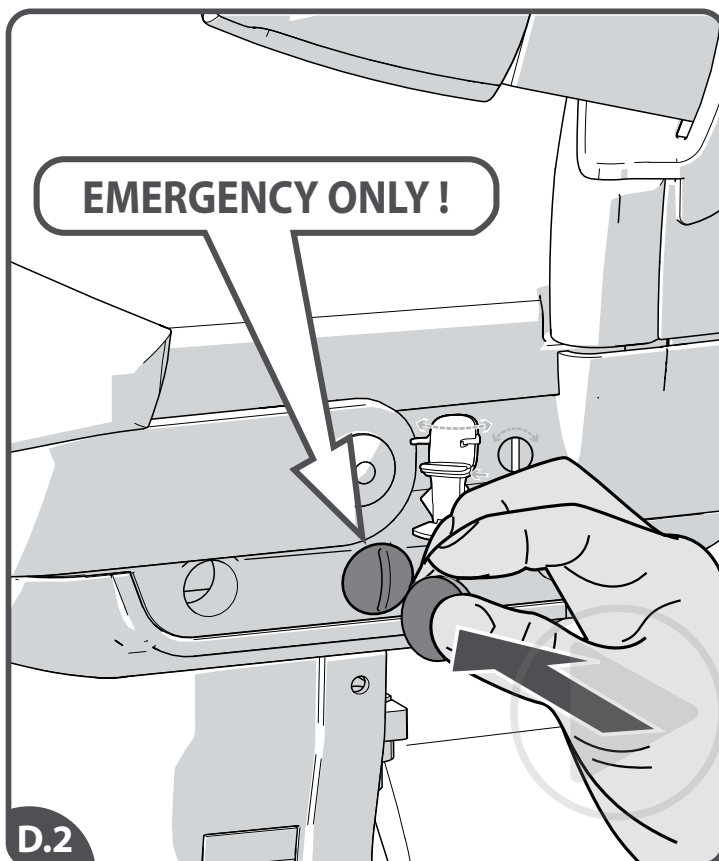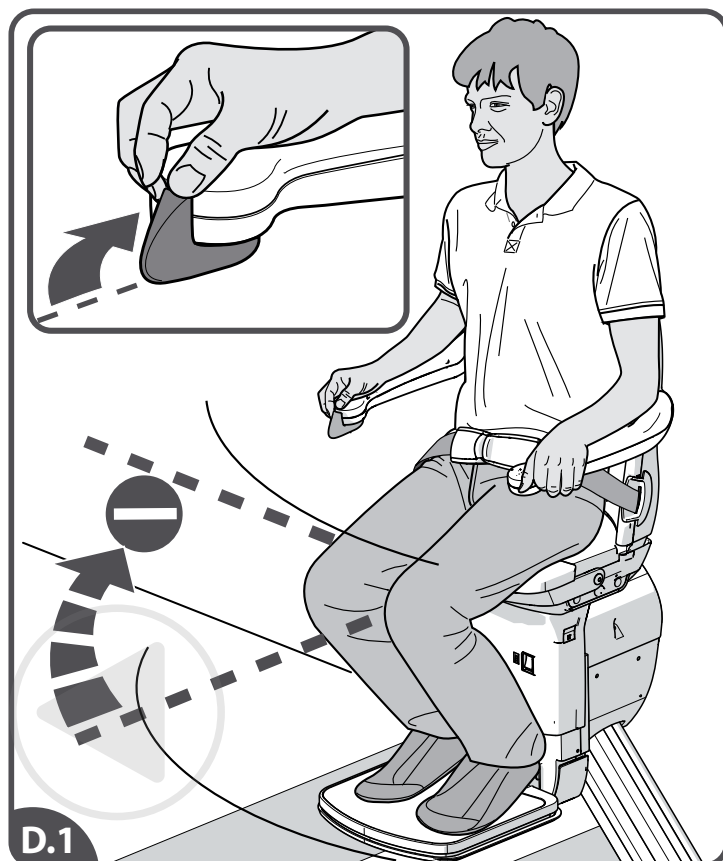- [D.1] If the power swivel is not functioning it can be swivelled manually. Do this only in an emergency.
- [D.2] Insert a coin in the button on the side of the seat.
- [D.3] Swivel the seat manually.
- [D.4] The seat must be returned to the start position. Remove the coin. Do this before you use the lift again.

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
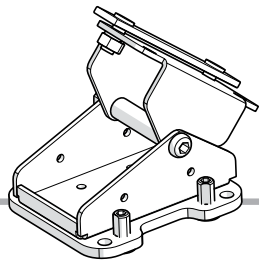
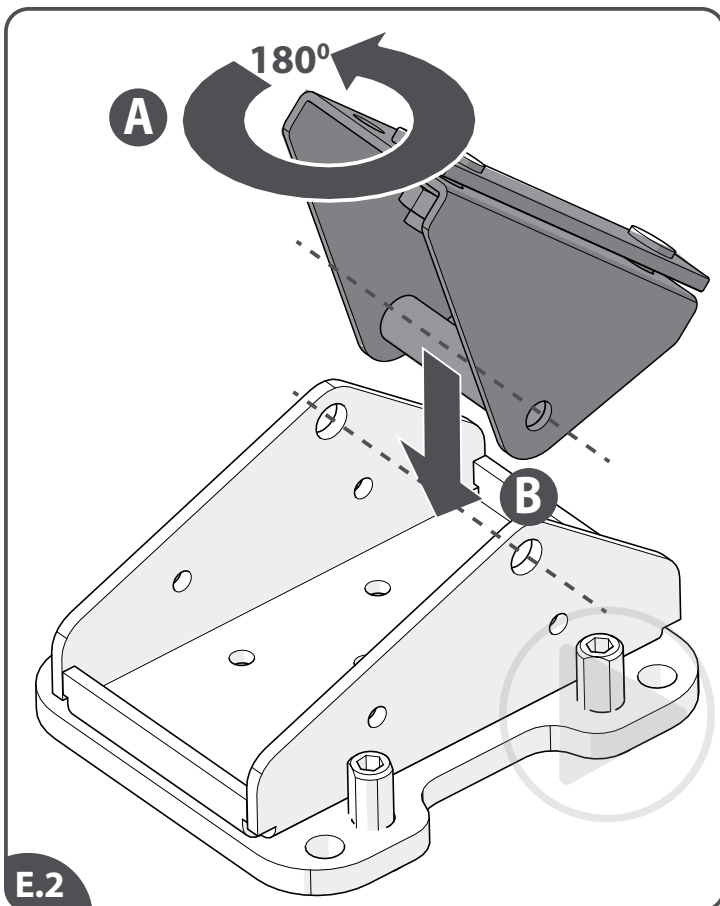
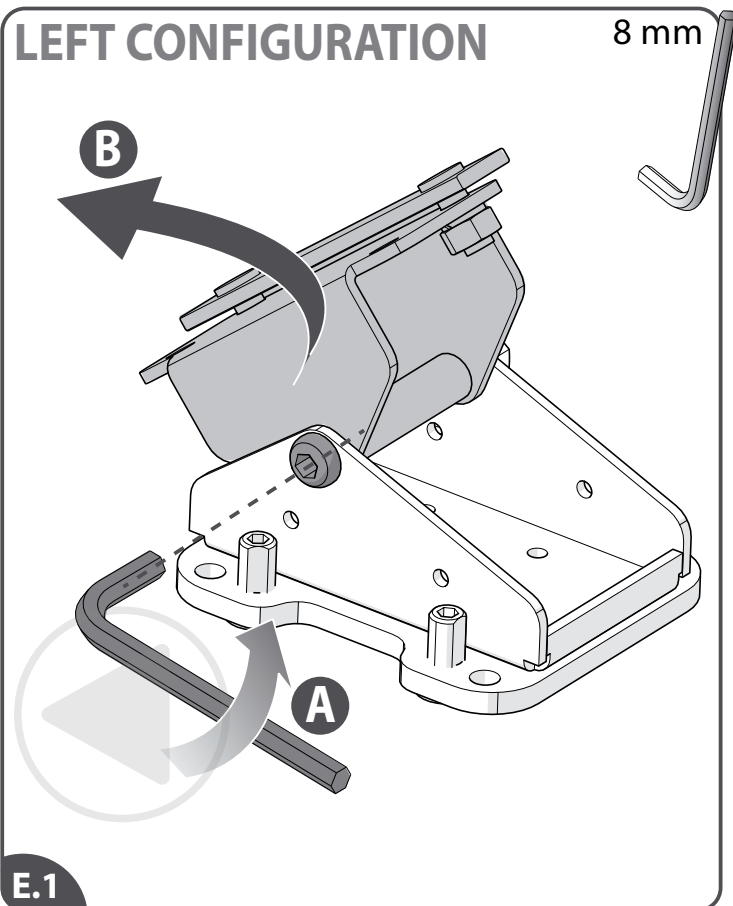
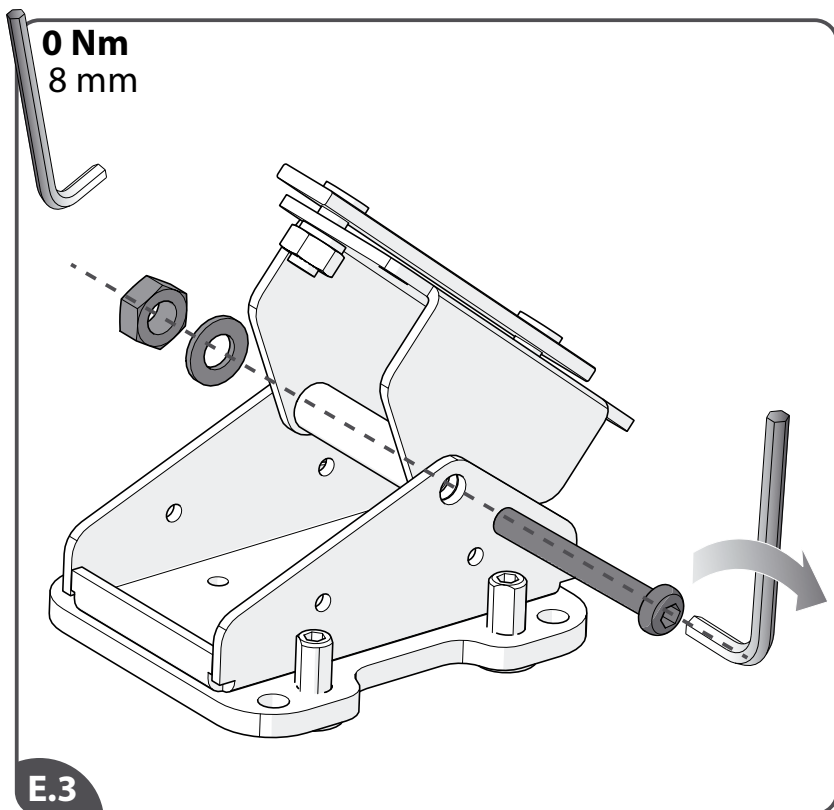
RIGHT CONFIGURATION

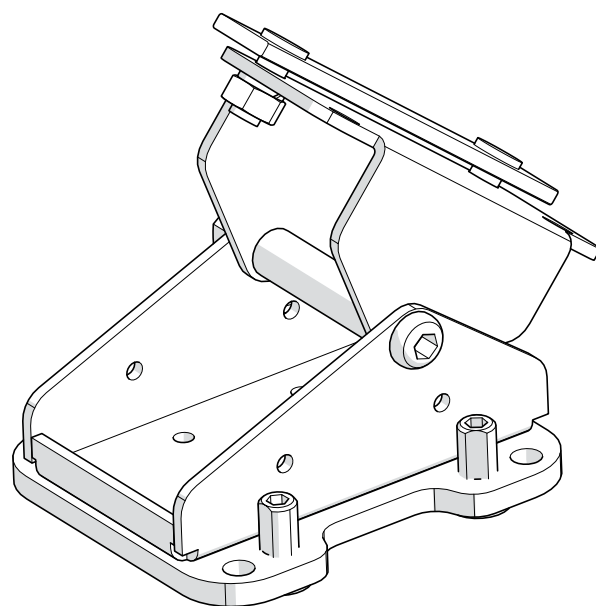3x or 5x

- [E.1] Loosen the pivot bolt{A} and remove the upper foot part{B}.
- [E.2] Rotate the upper foot part 180°{A} and position this part{B}.
- [E.3] Fix the upper foot part to the alternative pivot point. Do not tighten the bolt yet.
- [E.4] Repeat this procedure for the other feet.

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
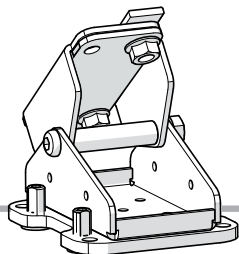
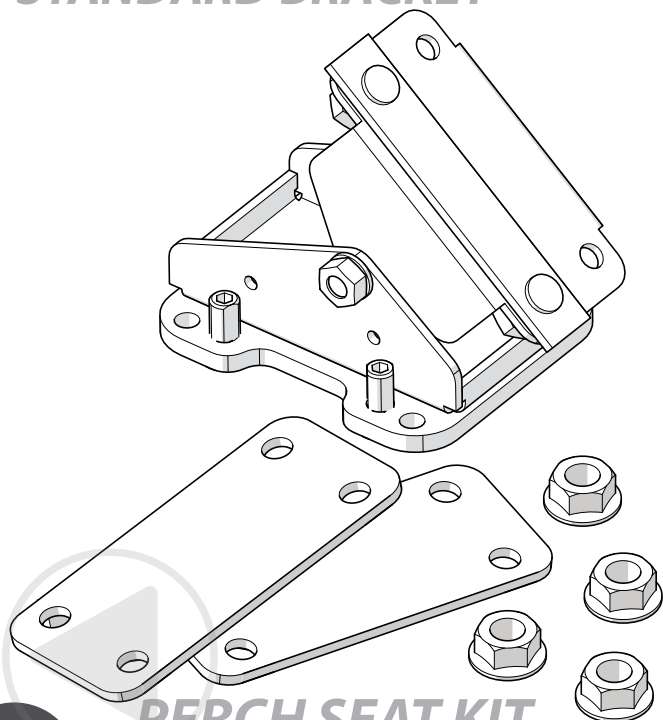
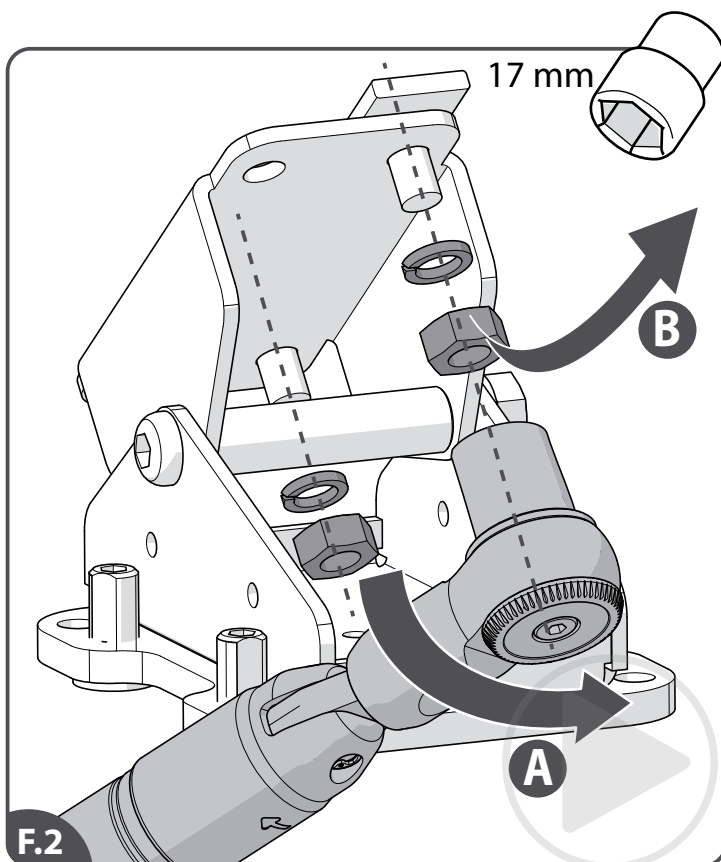
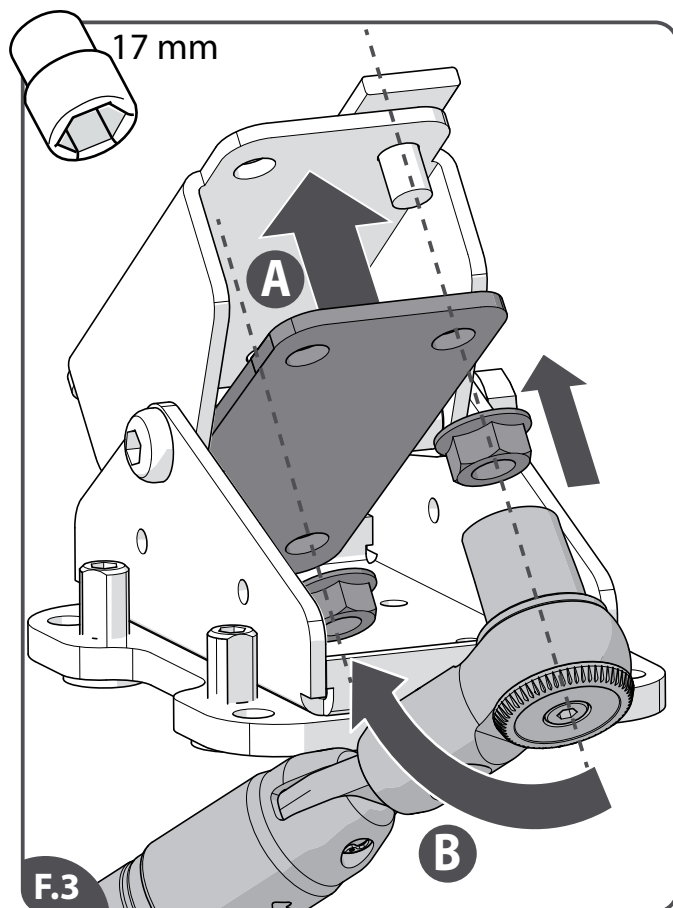
PERCH SEAT CONFIGURATION

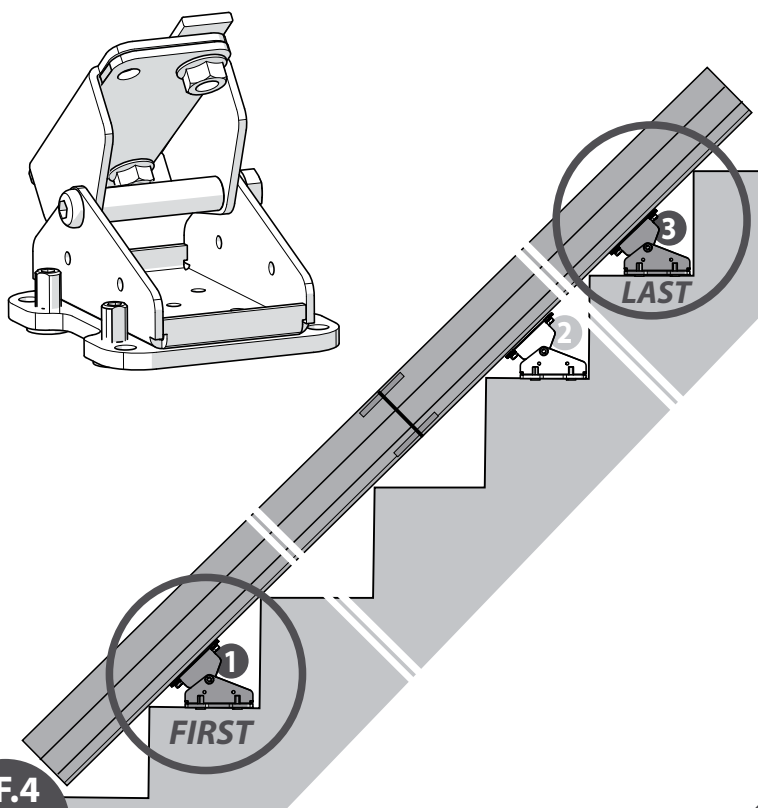F.4

- [F.1] Use the PERCH seat kit to configure the standard foot.
- [F.2] Remove the bolts/washers(2x).
- [F.3] Install the extra plate and flange nuts(2x).
- [F.4] Repeat this procedure for the second foot. Position the two PERCH seat feet at the top and bottom of the rail.

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






G

OUTDOOR COVER

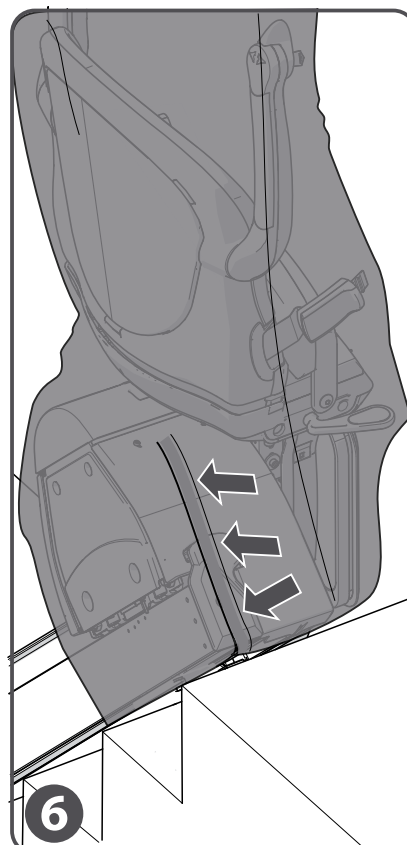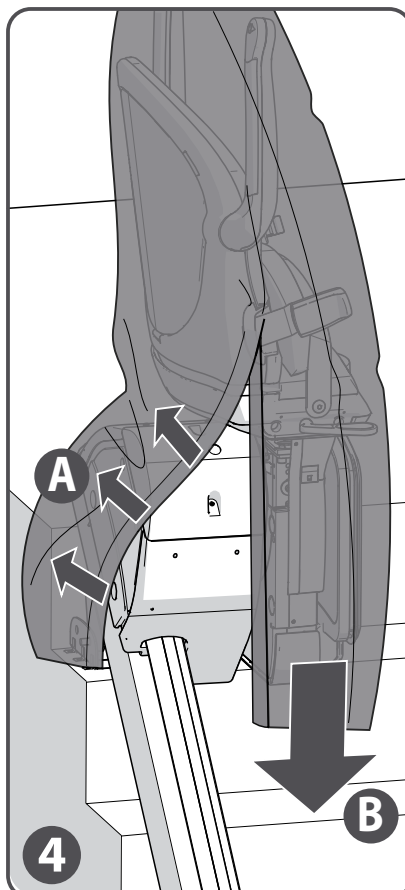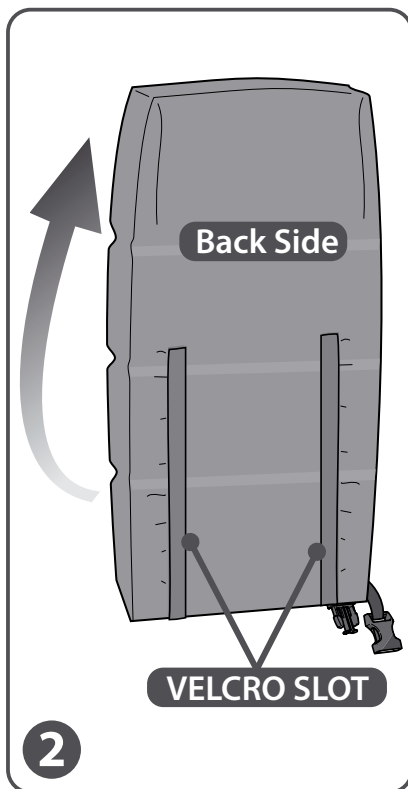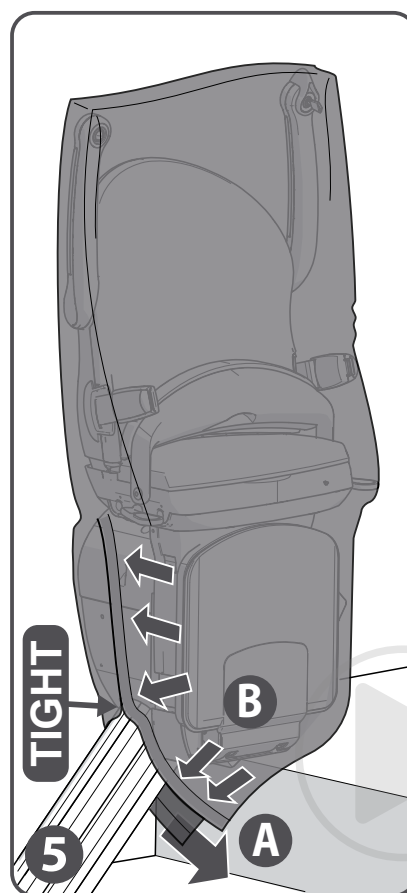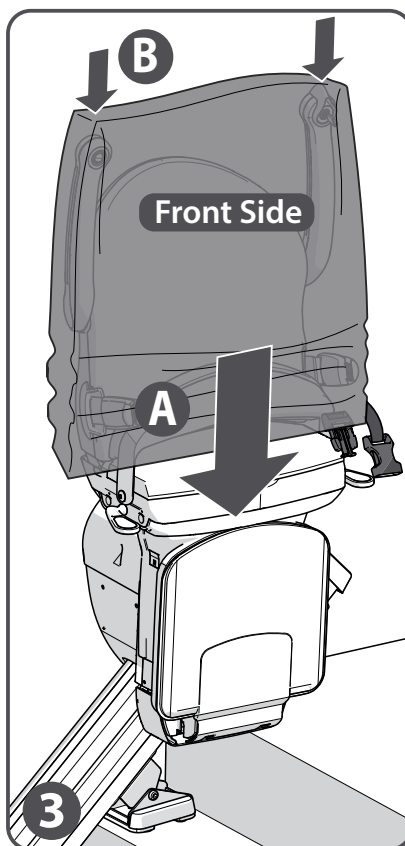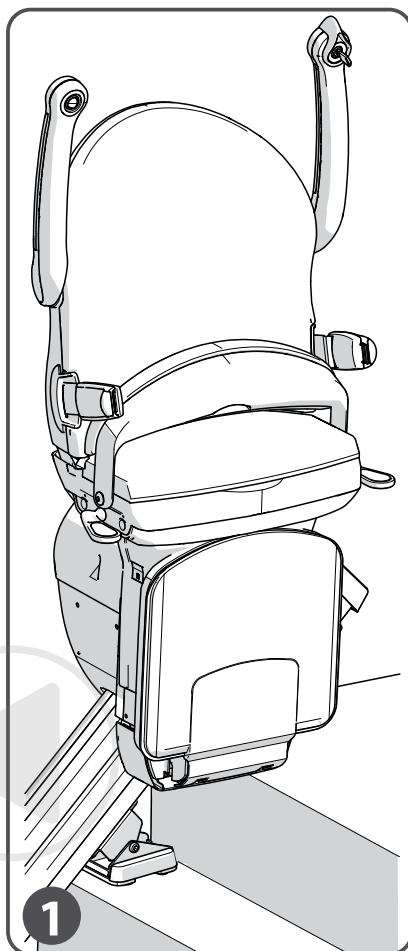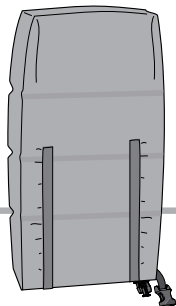It is recommended that the cover is always fitted and removed by an able-bodied person.

- [1] The seat must be in the upright position (arms up, seat base up, and footplate up).
- [2] Unfold the cover, the front of the cover is plain, the rear of the cover has Velcro slots on each side.
- [3] Place the cover over the seat with the front of the cover facing towards you{A}. Feed the arms of the seat into the corner seams of the cover{B}.
- [4] Unseal the Velcro sides{A} and pull the bottom of the cover down as far as possible{B}.
- [5] Pull the bottom of the rear side of the cover underneath the track to the front{A}. Seal the Velcro sides around the track where needed as tightly as possible, making sure it forms a continuous seal with no gaps{B}.
- [6] On the other side of the chair: depending on the configuration of the lift, one of the Velcro sides may cover the end of the track.

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





IT

G **OUTDOOR COVER**

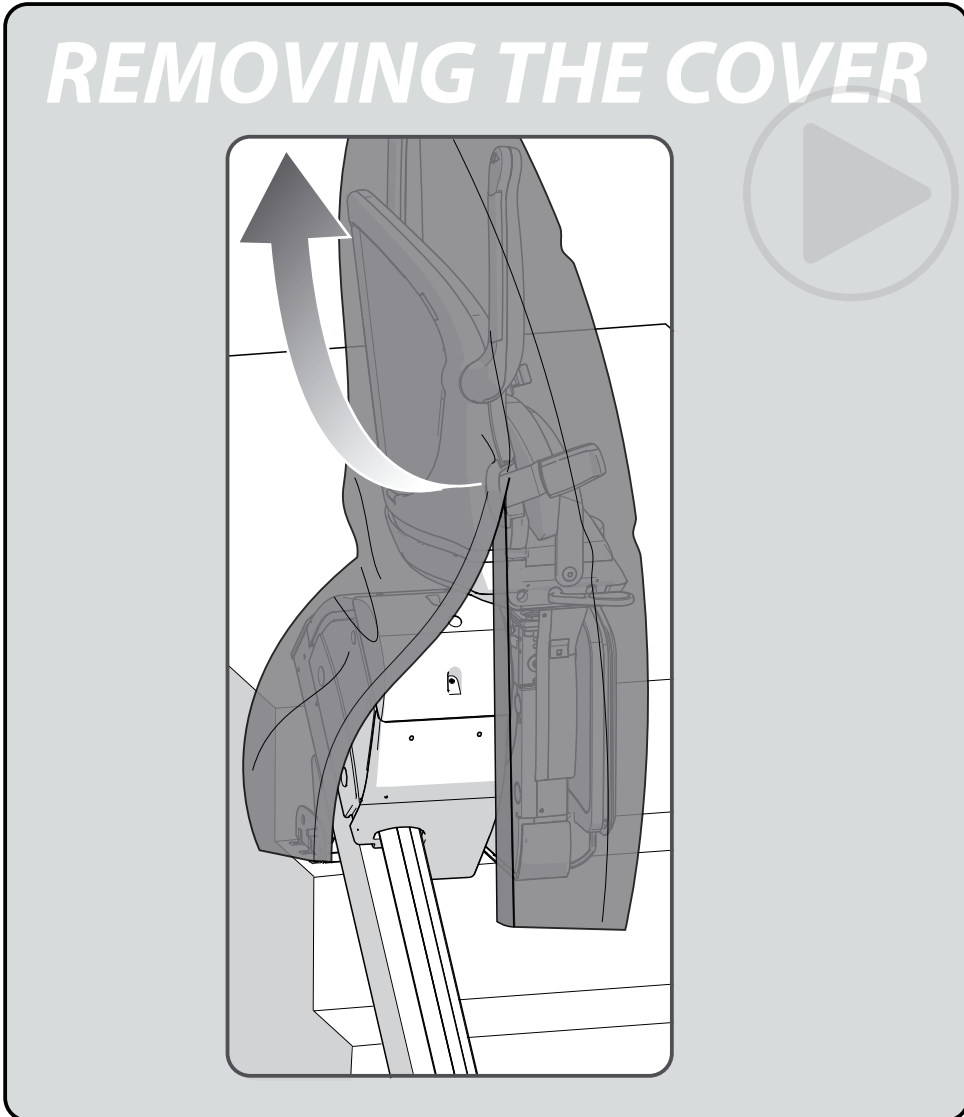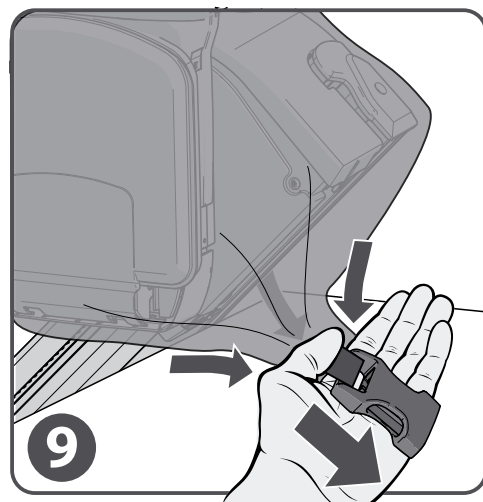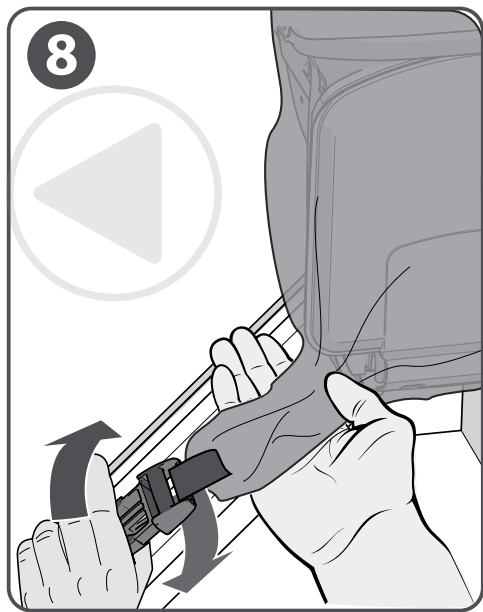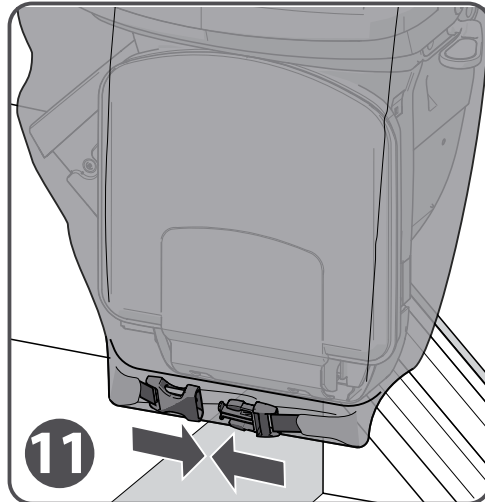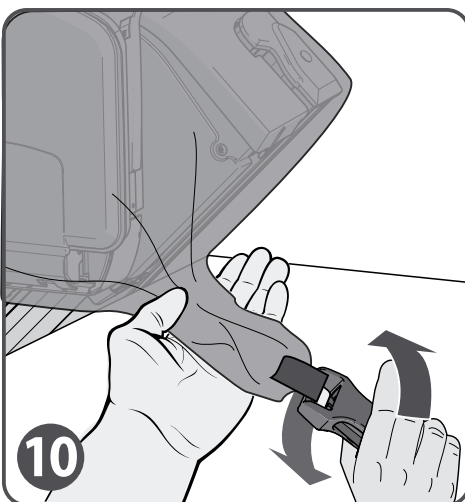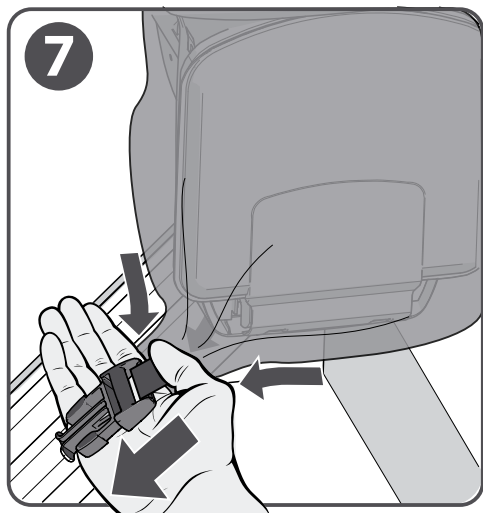[7]At the bottom of the cover, hold the left hand side clip in your left hand. Tuck/fold any excess of the cover to the left of the clip towards the centre of the cover (Make sure you are holding some of the cover from the rear of the lift). [8]Rotate the clip in your left hand and turn in a clockwise direction, gathering up the material on the left hand side. Stop when the material is in a tightly compressed roll. [9+10]Repeat this procedure on the other side of the lift. [11]Connect the two clips together.

Removing the Cover:

[1]Unclip the clip located at the underside of the lift. [2]Unravel the material and unseal the Velcro sides. [3]Lift the cover off the seat, fold the cover and store it in a safe location or take it with you to be refitted at the other end of the staircase. [4]Unfold the armrests, seat base and footplate (arms down, seat base down, and footplate down). [5]The lift is now ready to be used.

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



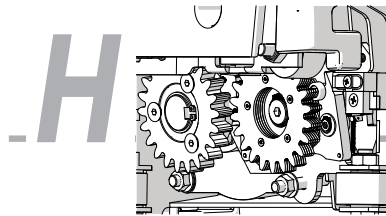
IND_OSG TEST

LIFTS AFTER 07-2024

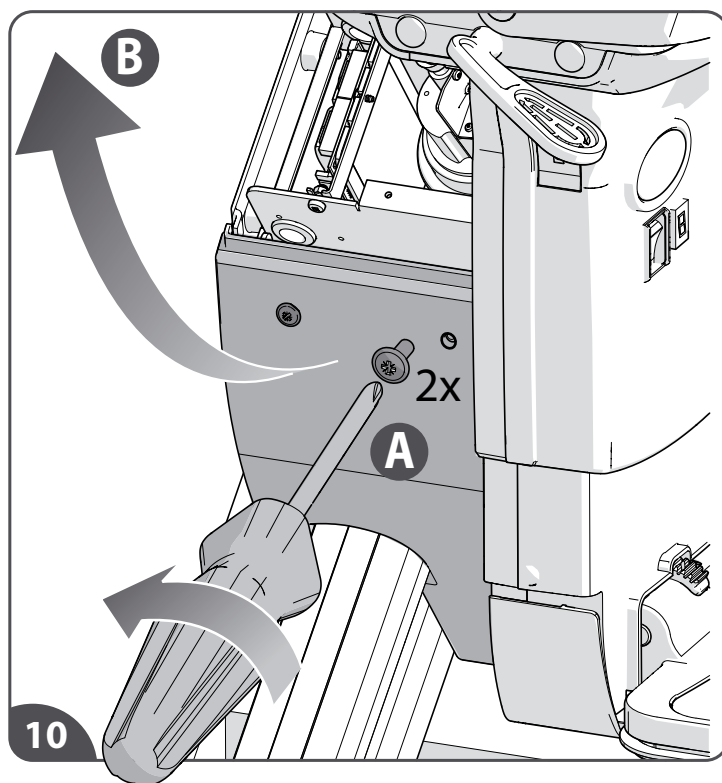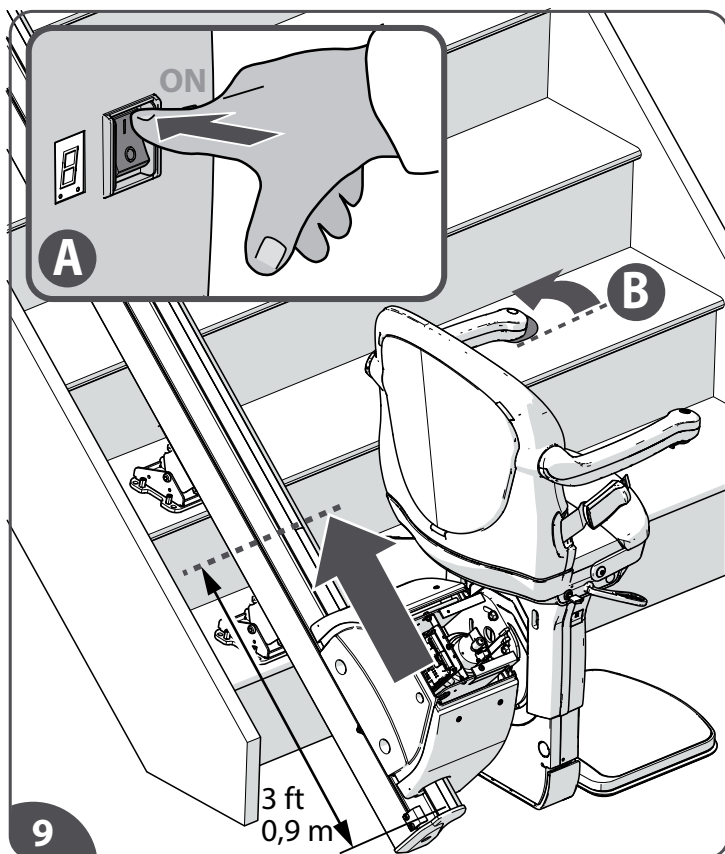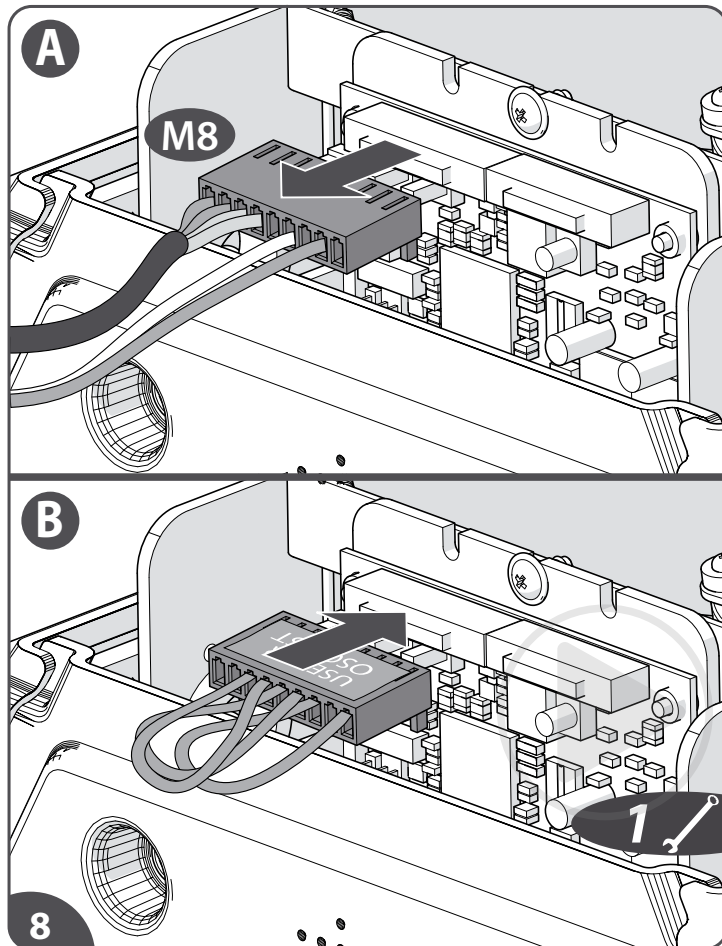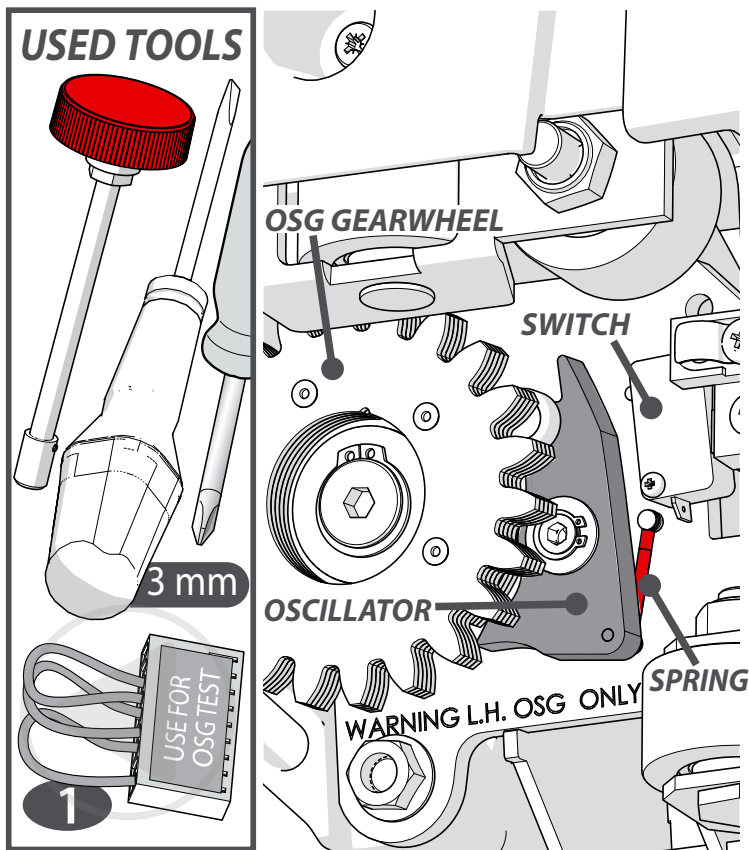The OSG assembly is a component of gears, cams, and catchment pockets. An oscillator arm track the cams and is pulled into engagement by an instrument spring.

[8] Disconnect the 9-way seatloom on the HAN125 PCB{A}. Connect the supplied OSG test connector{B}.

[9] Turn on the lift{A}. Drive the lift so that it is at least 3 ft(0,9 m) from a boarding point{B}.

[10] Remove the lower side cover (for a RH lift: upper side cover) of the power pack.

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



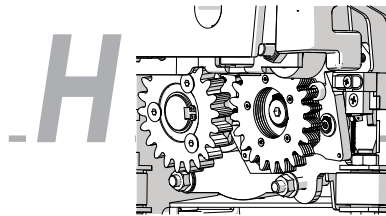
IND_OSG TEST

LIFTS AFTER 07-2024

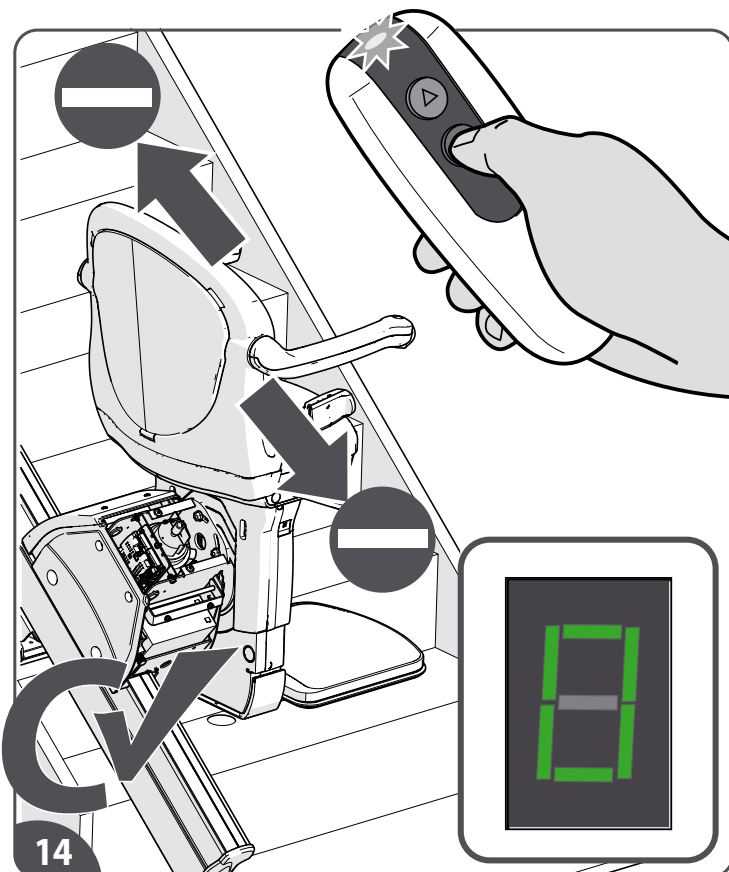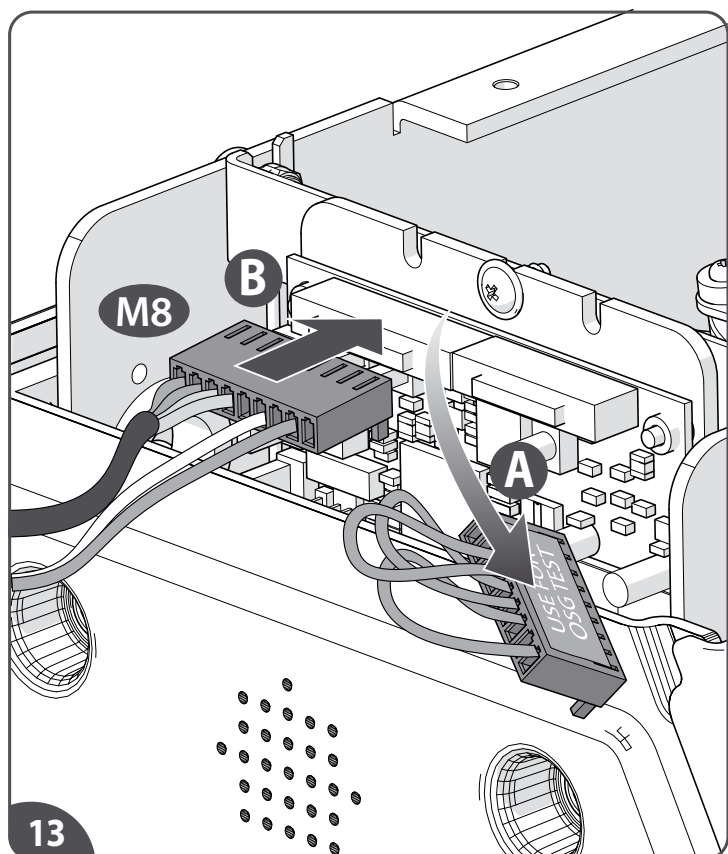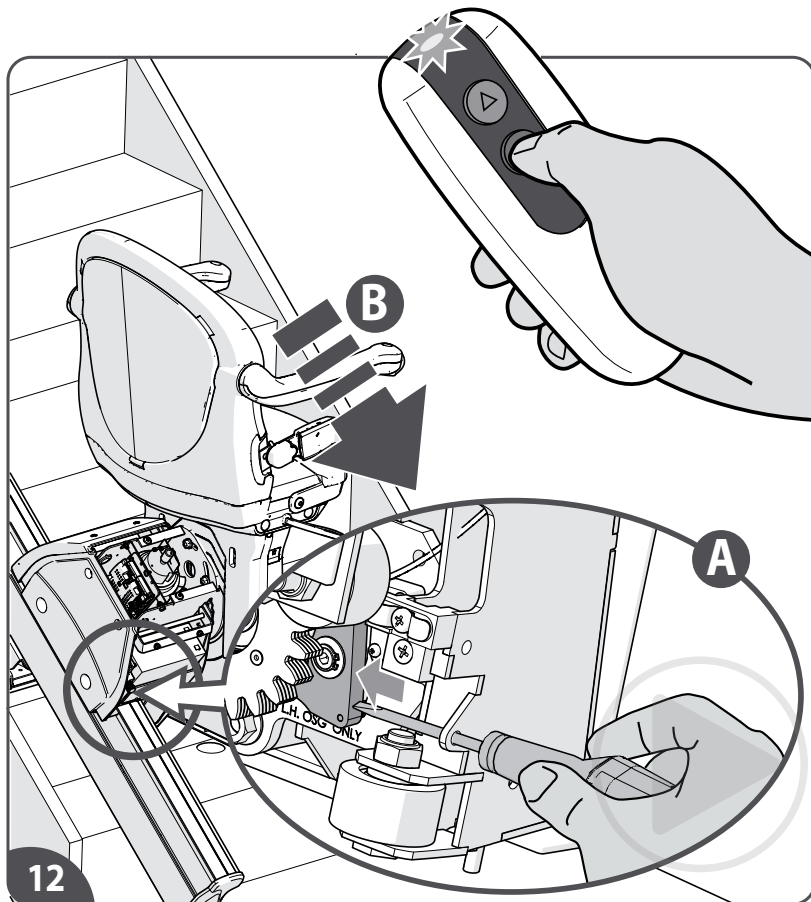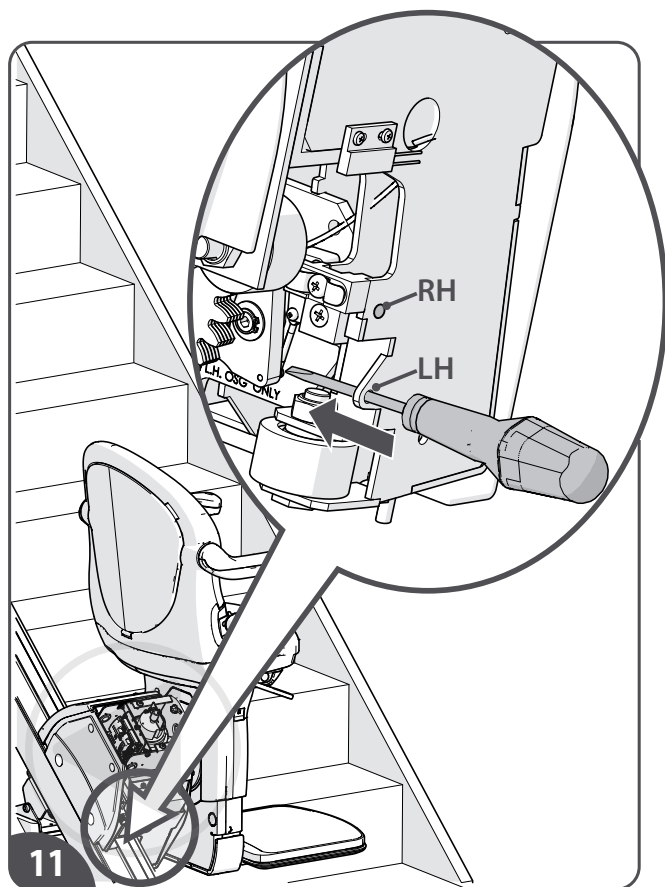- [11] Fit a 3 mm flat screwdriver in the access hole of the OSG (access RH on top, access LH on bottom of power pack frame).
- [12] Maintaining a light pressure on the OSG activation lever from underneath the power pack{A}, drive the lift DOWN using the remote control until the OSG is activated{B}.
- [13] Remove the OSG test connector{A}. Reconnect the original 9-way seatloom on the HAN 125 PCB{B}.
- [14] Attempt to drive the lift. Operation of the OSG is confirmed by the lift not driving and Code "0" on the display on the frontcover (Safeties / OSG active).

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



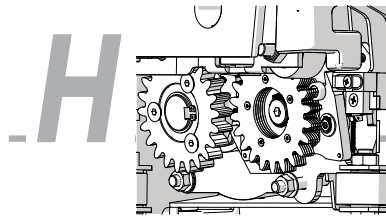
IND_OSG TEST

LIFTS AFTER 07-2024

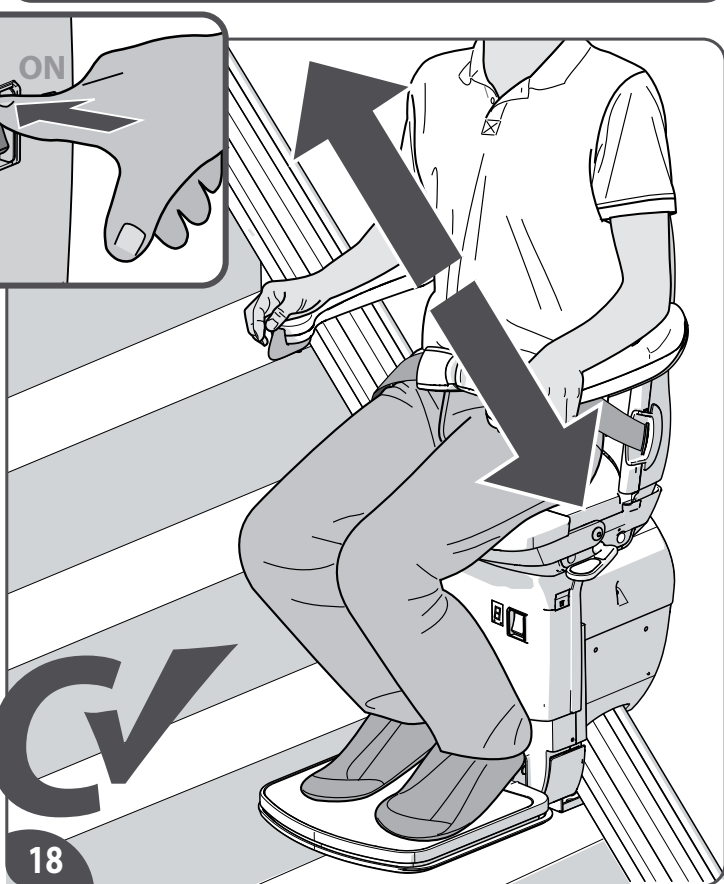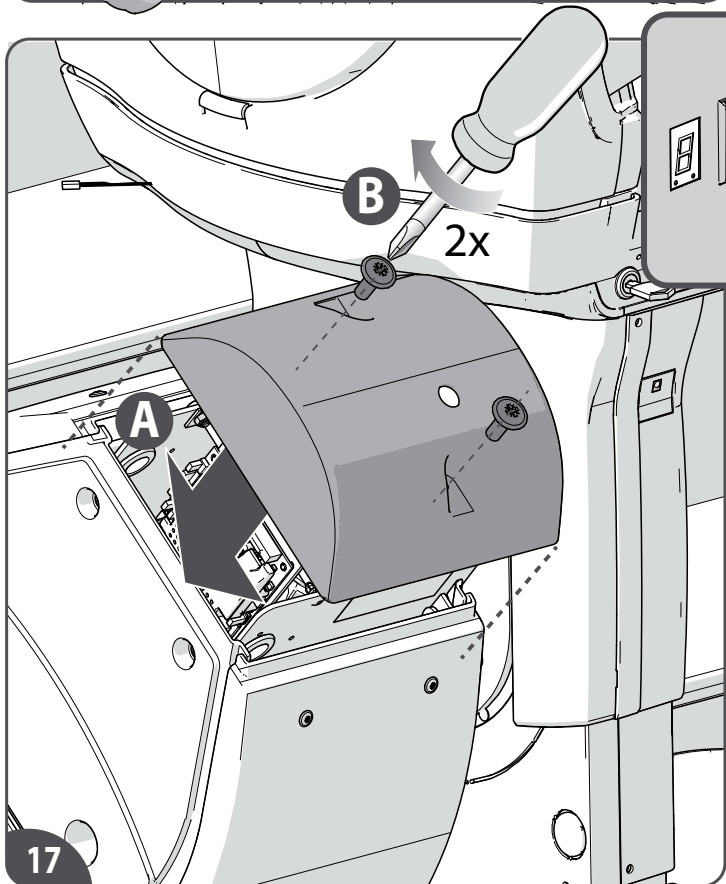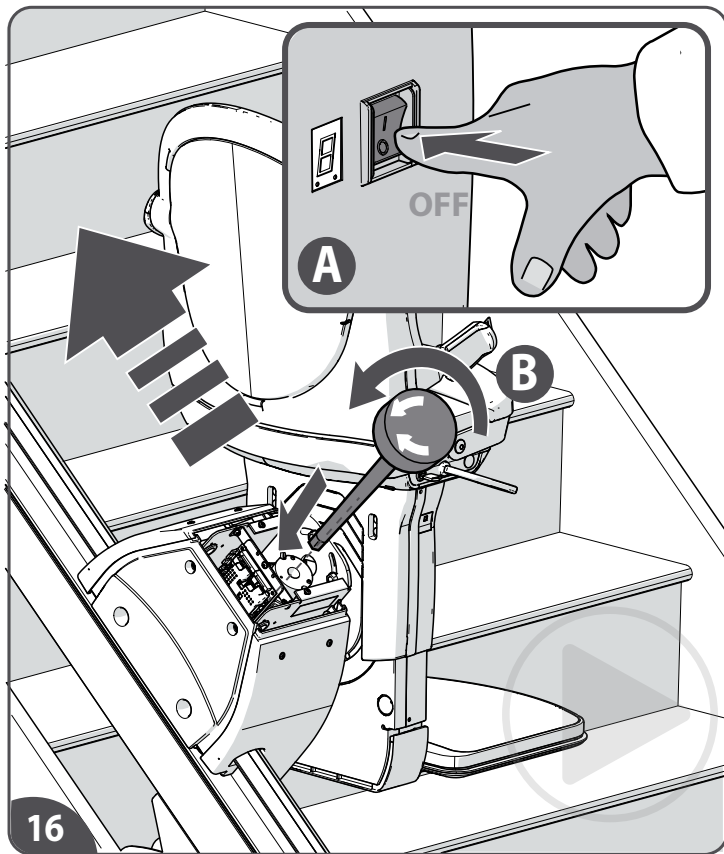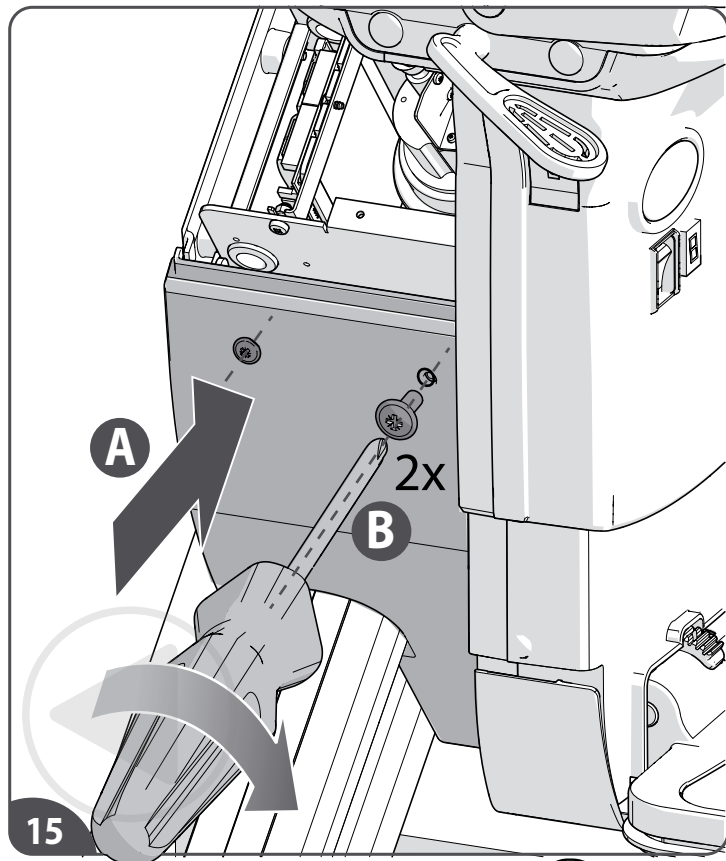- [15] Remove the screwdriver from the OSG assembly and reinstall the side cover.
- [16] Turn off the lift {A}. Use the hand winder to move the lift up a short distance to release the OSG {B}.
- [17] Refit {A} and fasten {B} the upper cover from the power pack.
- [18] Turn on the lift. Drive the lift UP / DOWN to confirm that the OSG has been released and the lift drives normally.

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





IT

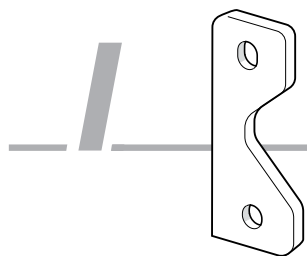REINFORCEMENT PLATE

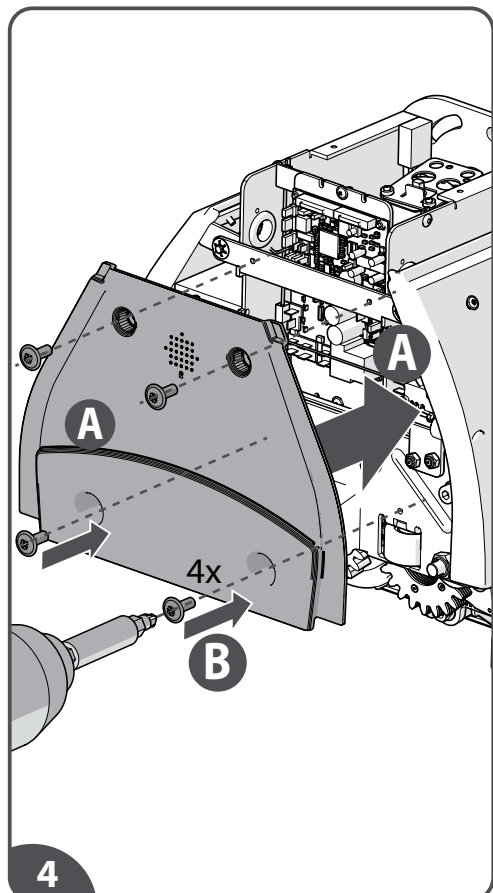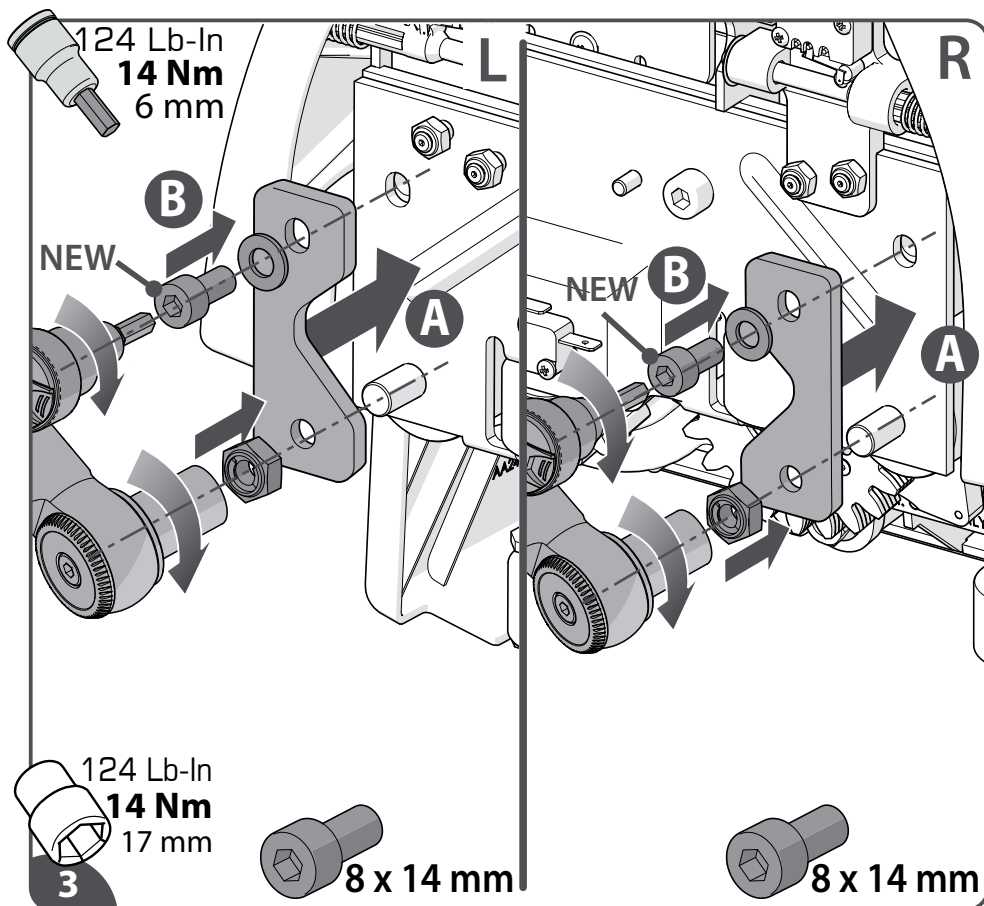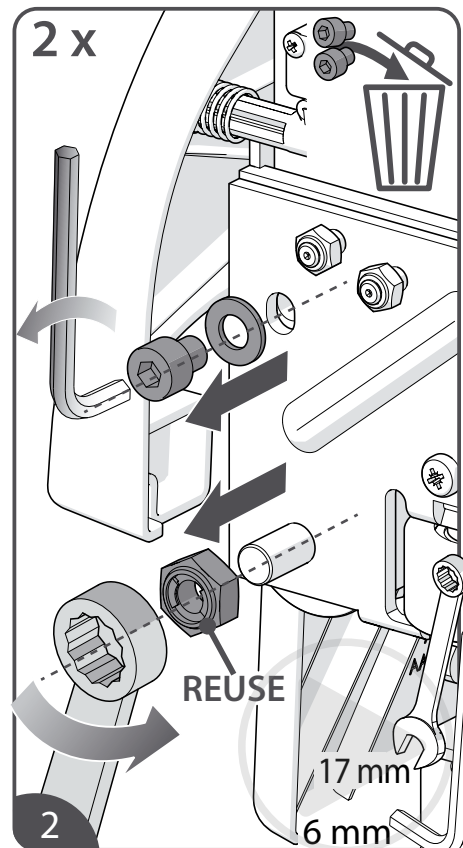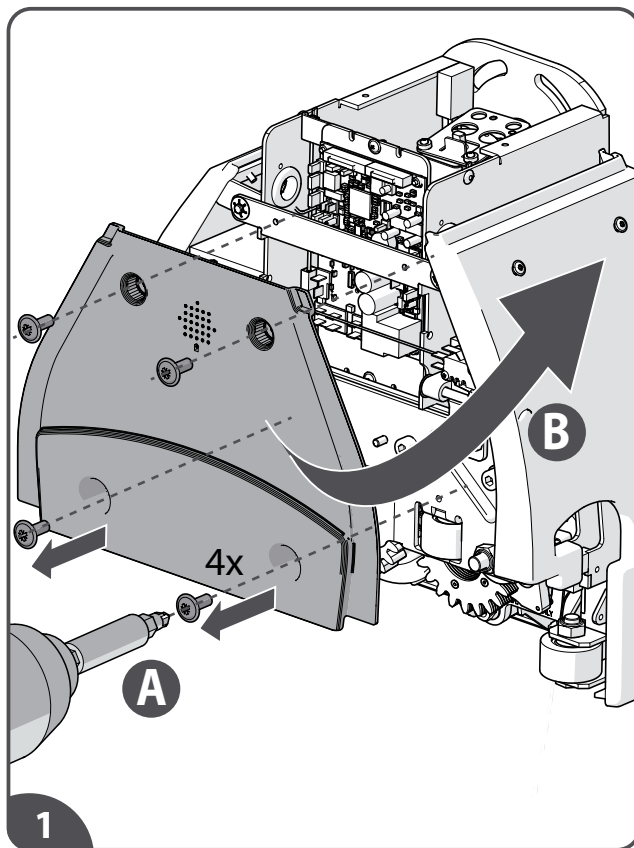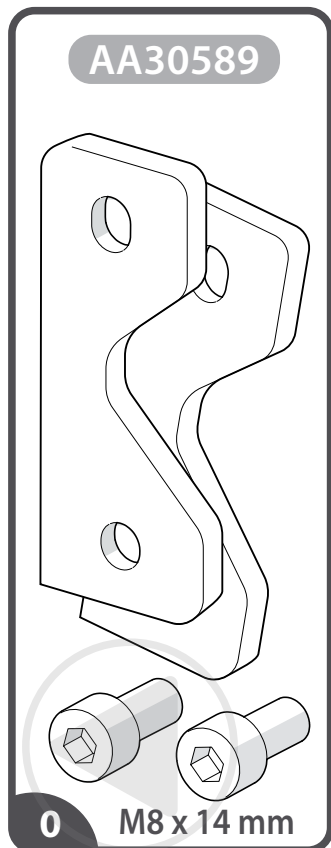- [0] Check the kit.
- [1] Remove the back cover.
- [2] Remove the bolt and nut from the roller plate on the back side. Reuse the nut.
- [3] Install on the left and right side the reinforcement plate.
- [4] Reinstall the back cover.

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




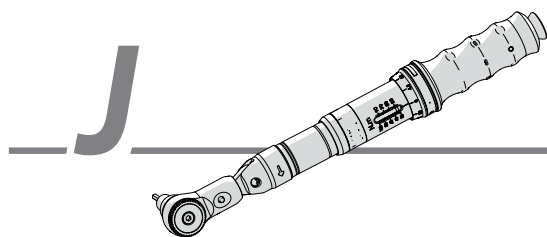
TORQUE WRENCH TEST

CHECK FASTENERS

After the test drive, check the torque of the load-carrying fasteners with a torque wrench

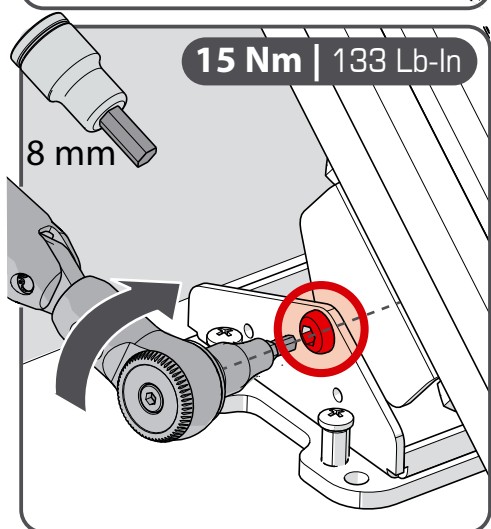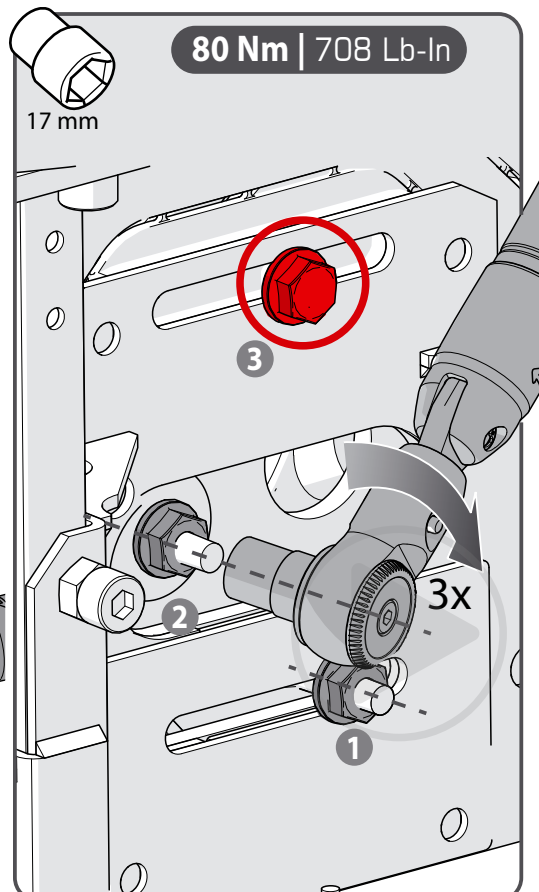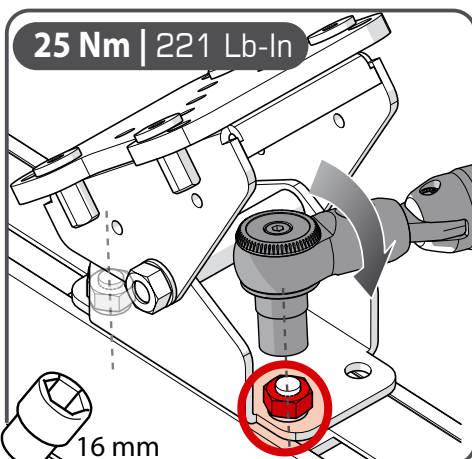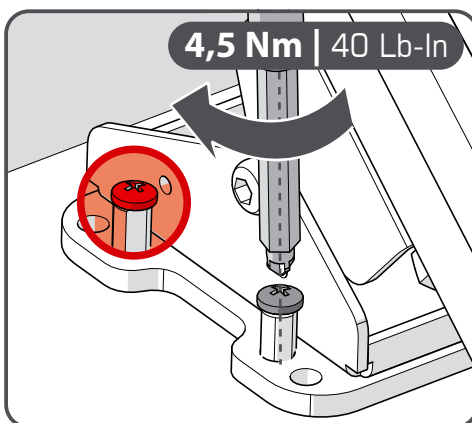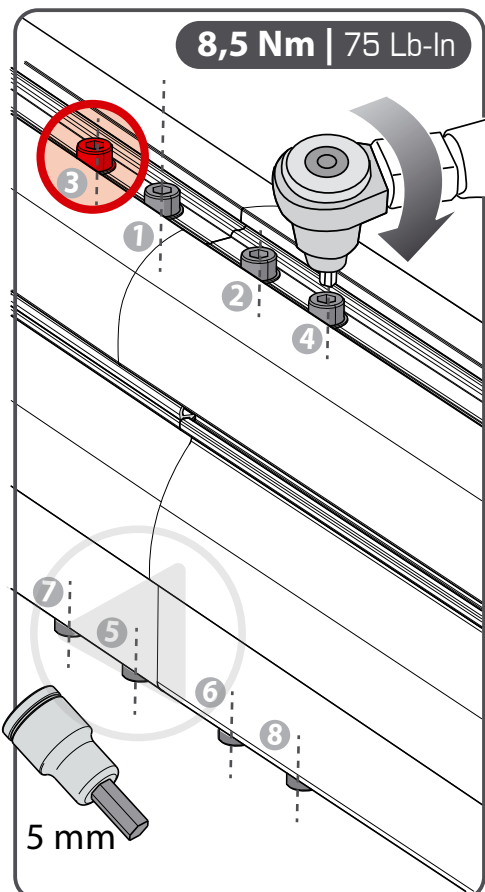K PROTECTION EQUIPMENT EXPLAINED

When performing any work on the stairlift, take the following precautions:

- Wear safety shoes type S3 or equivalent protection (S3 is an European classification; outside of Europe there are other classifications).
- Wear safety glasses for eye protection, i.e. to prevent injury because of parts suddenly striking into the direction of your face such as springs or when drilling in floors and walls.
- Use smooth fitting gloves to move and manage heavy / sharp parts of the unit and the rail.
- Wear a safety cap for installation with low bulkheads.
- Wear ear plugs when drilling and cutting. Exposure to sounds above 85 decibels can cause hearing loss.
- To prevent getting caught in moving parts of the stairlift: do not wear loose clothing or jewelry, such as necklace, bracelets and tie long hair.

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VERY IMPORTANT

EN

This product must ONLY be installed by properly qualified and trained stairlift installation engineers, who have satisfactorily completed a Handicare-approved training course on how to install and/or commission the Handicare 1000-stairlift. Grave consequences regarding the safety of the stairlift in question and of its users may result from any violation of this requirement.

Intended Audience

This manual is intended for use by such qualified and trained stairlift installation engineers ONLY; it is NOT intended to show an unskilled person how to install and/or commission the stairlift, nor does it allow an unskilled person to install and/or commission the stairlift in an appropriate/safe manner.

This manual may NOT be used to replace the necessity for training.

The procedures described in this manual are 100% mandatory and MUST be followed.

Required Knowledge

All users of this manual MUST have the following knowledge:

Sufficient understanding of the key features, functions and capabilities of the Handicare 1000 stairlift;

Sufficient mechanical and electrical skills to be able to understand the function and purpose of each component of the stairlift, and sufficient mechanical and electrical skills as required for the correct installation and commissioning of the stairlift;

Sufficient Health and Safety awareness.

Disclaimer

It is the responsibility of the stairlift engineer to check that the manual is applicable to the type and version of the stairlift being installed; he/she MUST not fit a stairlift by reference to the manual if it is not the correct manual.

The stairlift engineer must read the manual to understand all the information to enable safe installation and/or commissioning of the stairlift.

In the absence of any negligence or breach of duty on the part of Handicare Accessibility Limited, no liability or responsibility will be accepted by it for any loss or damage arising from anyone's failure to follow and adhere to the guidance in the manual in connection with the installing and/or commissioning of the stairlift.

The latest and any previous versions of this manual can be obtained from www.handicarepartners.com

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Unit serial no:

Seat serial no:

- 1) Confirm that the installation dates have been completed on the EC Declaration of Conformity
- 2) Confirm that the remote controls are mounted on the wall and are functioning properly
- 3) Confirm that the safety belt is secure and functions correctly
- 4) Confirm that all warnings relating to the maximum working load on the lift have been affixed

2 Rail

- 1) Confirm that the rail joints are properly secured
- 2) Confirm that the rail supports are properly secured
- 3) Confirm that the folding rail is functioning correctly

3 Electrical drive unit tests

- 1) Confirm that the batteries are securely fitted in the carriage structure
- 2) Confirm that the following have been effectively indicated to the user:
 - a) location of the charge contacts.
 - b) "OFF CHARGE" alarm.
 - c) required action when stairlift is not charging

4 Clearance dimensions

- 1) Confirm the correct clearance dimensions as specified in the installation manual are maintained throughout the full travel of the stairlift
- 2) Confirm that there are no open bars/openings at the rear of the lift
- 3) Confirm that other users can safely use the staircase (minimum clearance of 400 mm between the folded seat and opposing obstacles). Install handrails if required.

5 Functional tests

- 1) Confirm that the final limit switches (upper and lower) operate correctly (position and operation)
- 2) Confirm that the control lever, key switch, emergency stop button, and the detection switches in the armrests operate correctly
- 3) Confirm that the safety detection devices operate correctly (footrest/unit/armrest)
- 4) Confirm that the swivel seat interlock and switch operate correctly
- 5) Confirm that the stopping distance of the lift is a maximum of 2 cm
- 6) Confirm that the restart time of the lift is at least 1 second
- 7) Confirm that the alarm device operates correctly when activated
- 8) Confirm that the lift has successfully passed the static (torque wrench test) and dynamic test ride in accordance with the instructions in the installation manual

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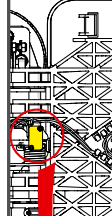
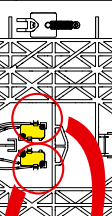
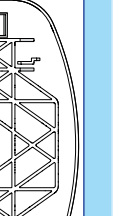
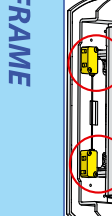
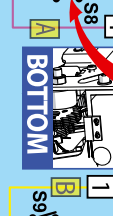
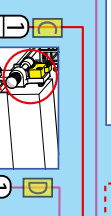
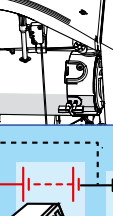
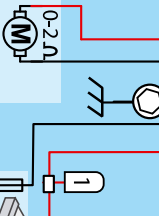
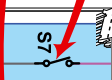
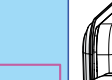
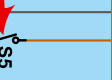
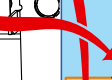
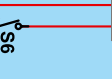
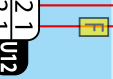
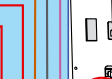
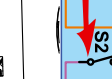
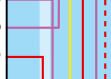
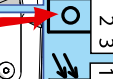
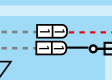
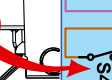
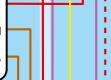
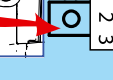
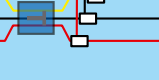
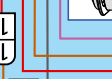
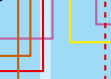
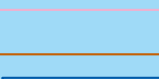
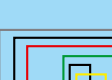
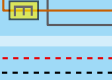
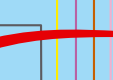
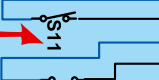
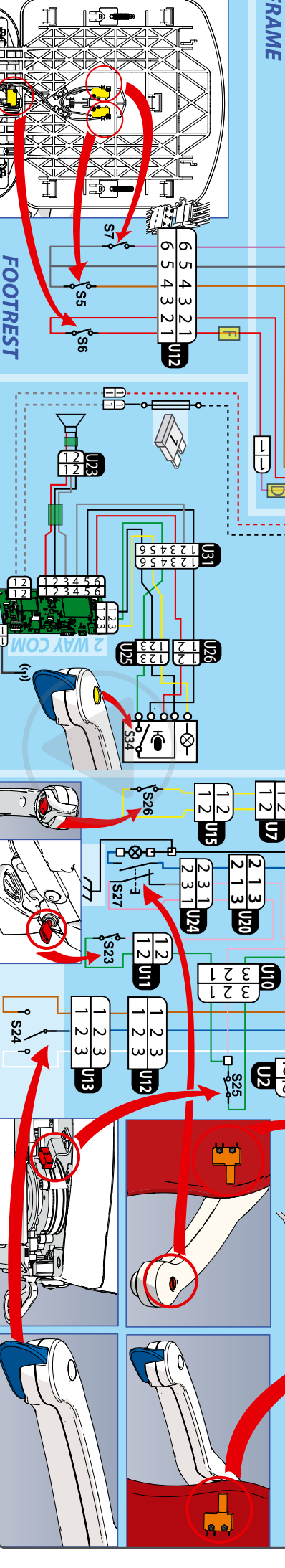
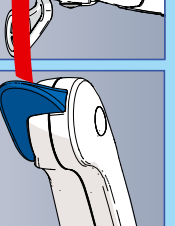
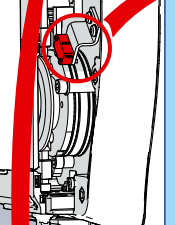
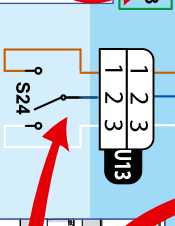
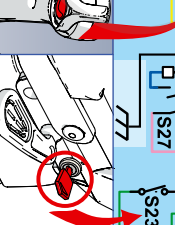
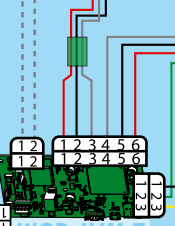
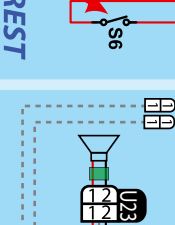
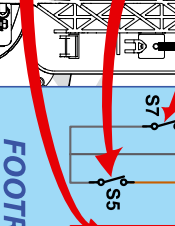
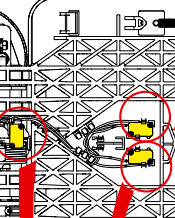
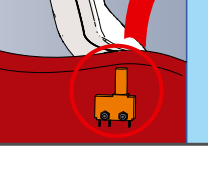
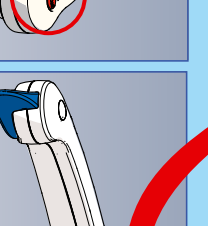
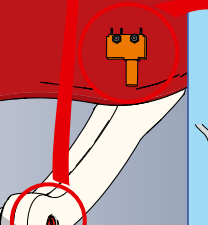
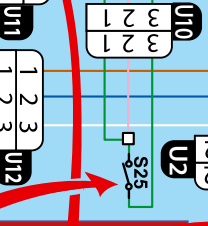
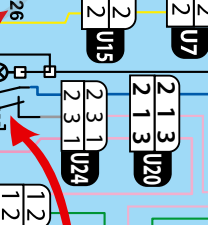
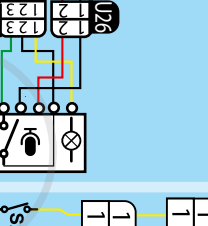
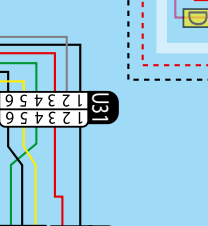
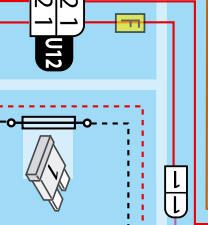
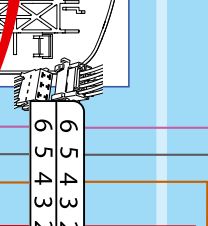
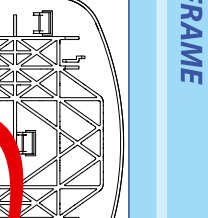
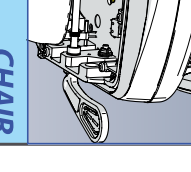
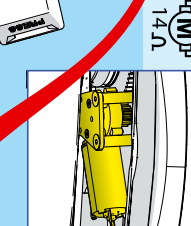
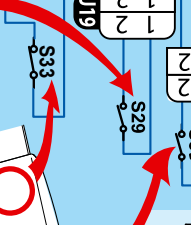
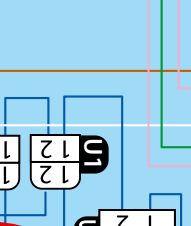
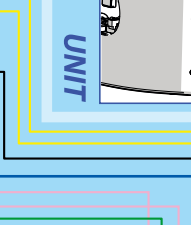
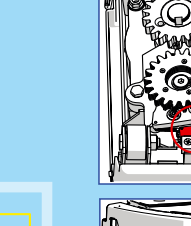
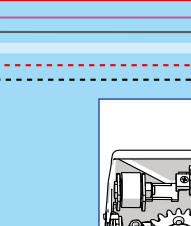
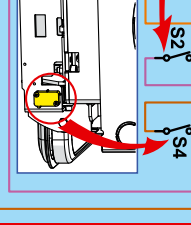
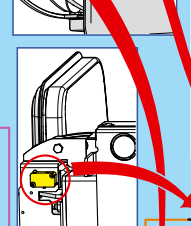
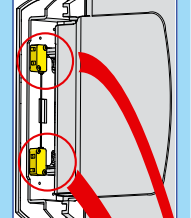
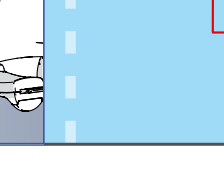
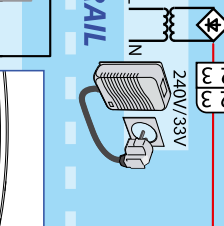
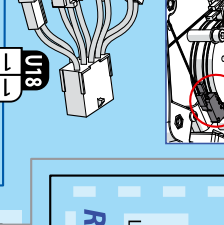
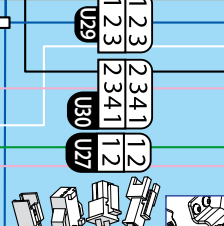
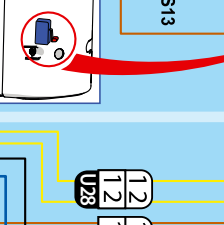
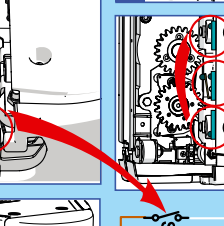
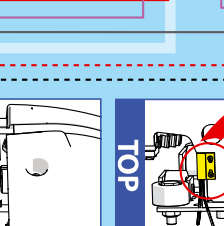
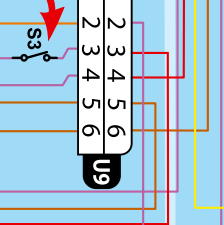
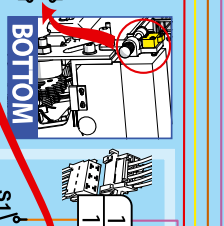
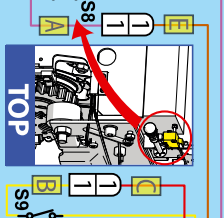
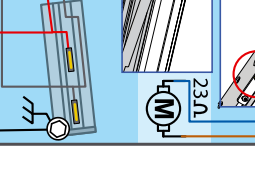
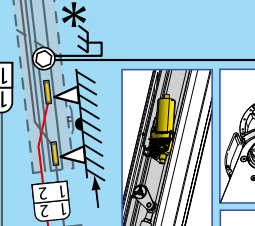
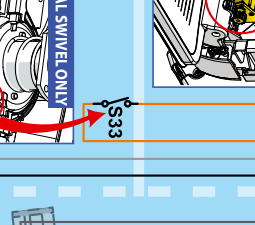
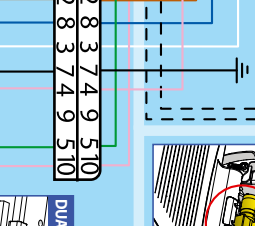
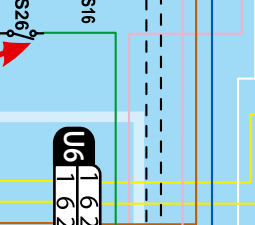
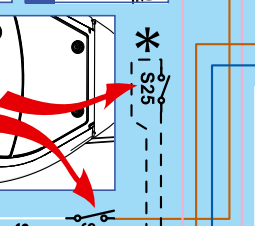
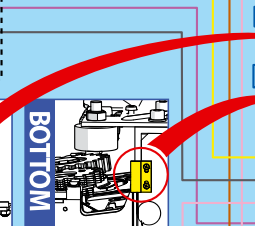
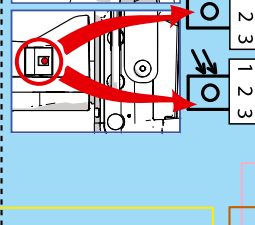
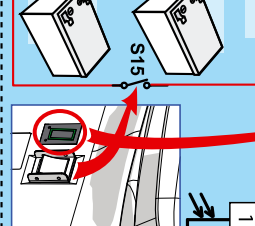
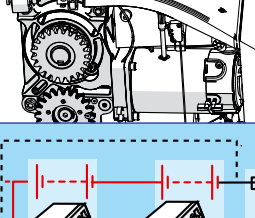
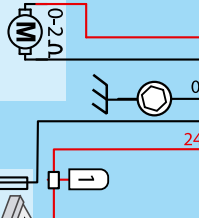
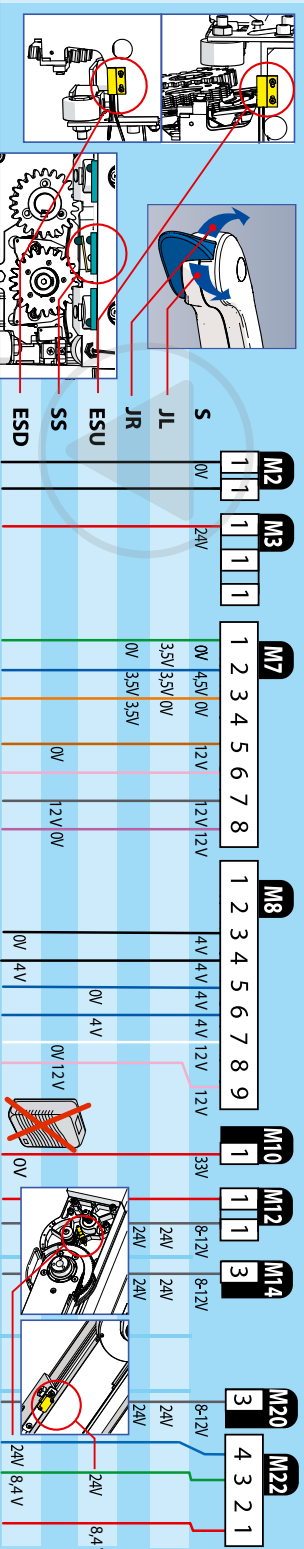
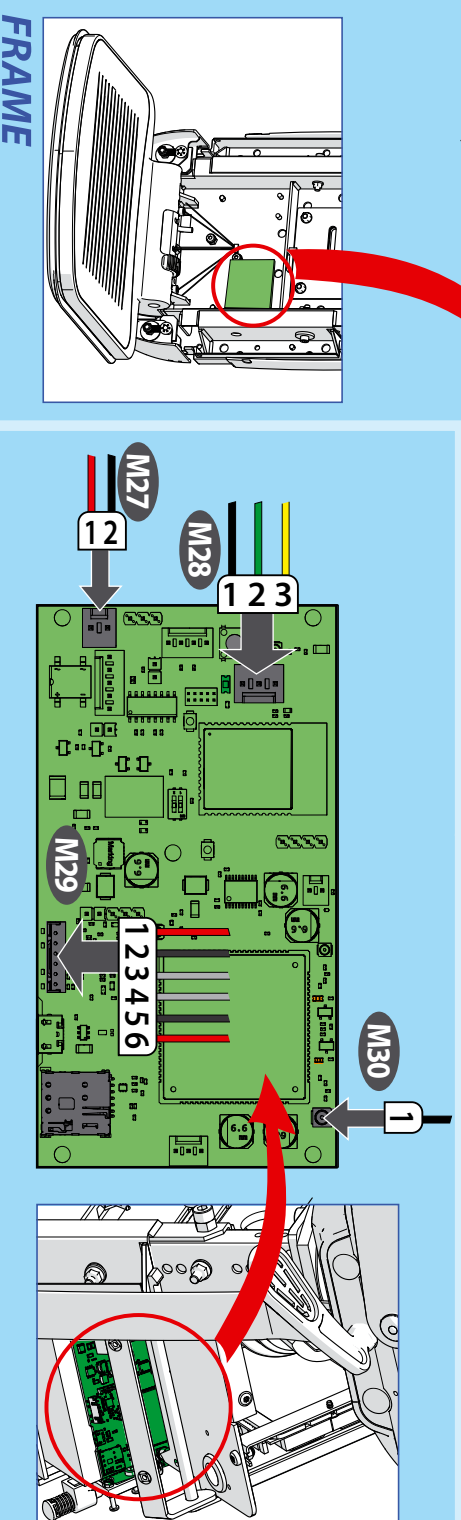
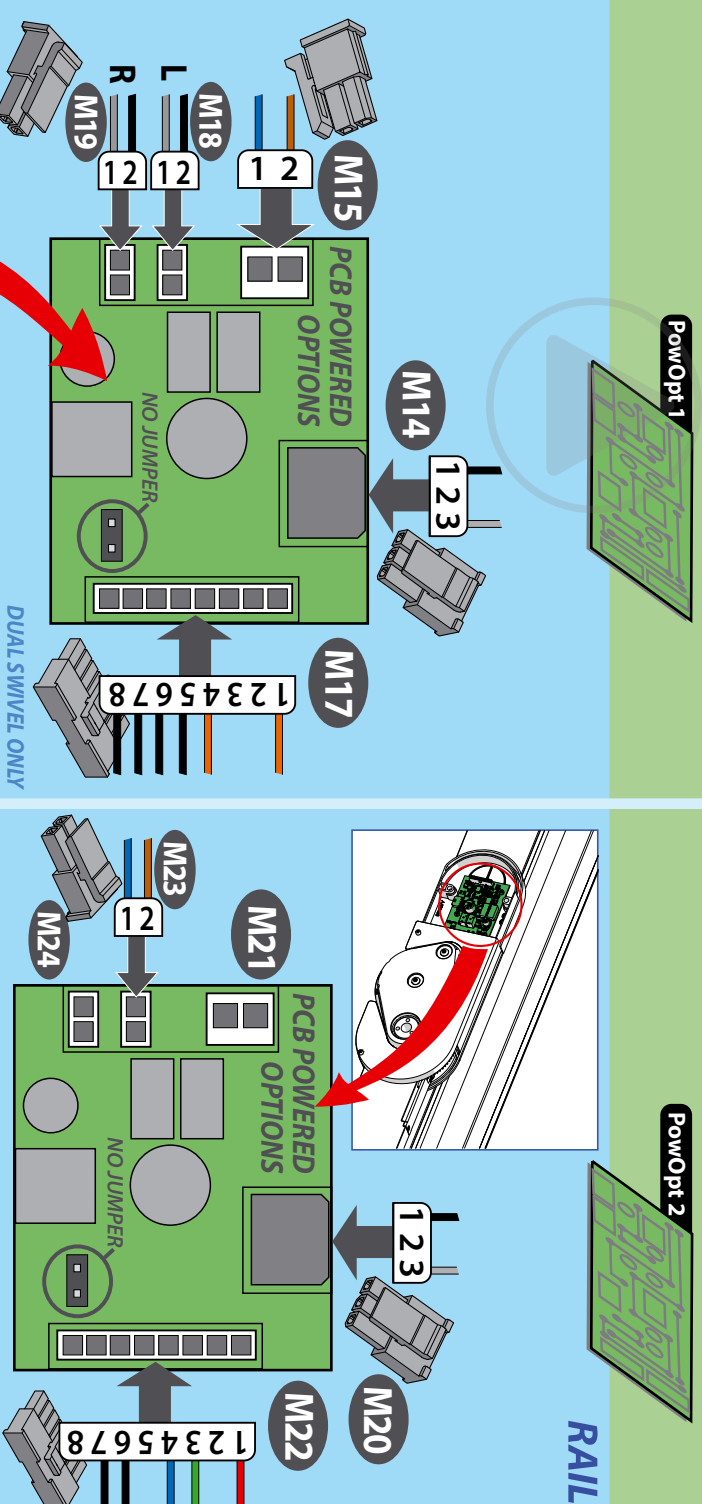
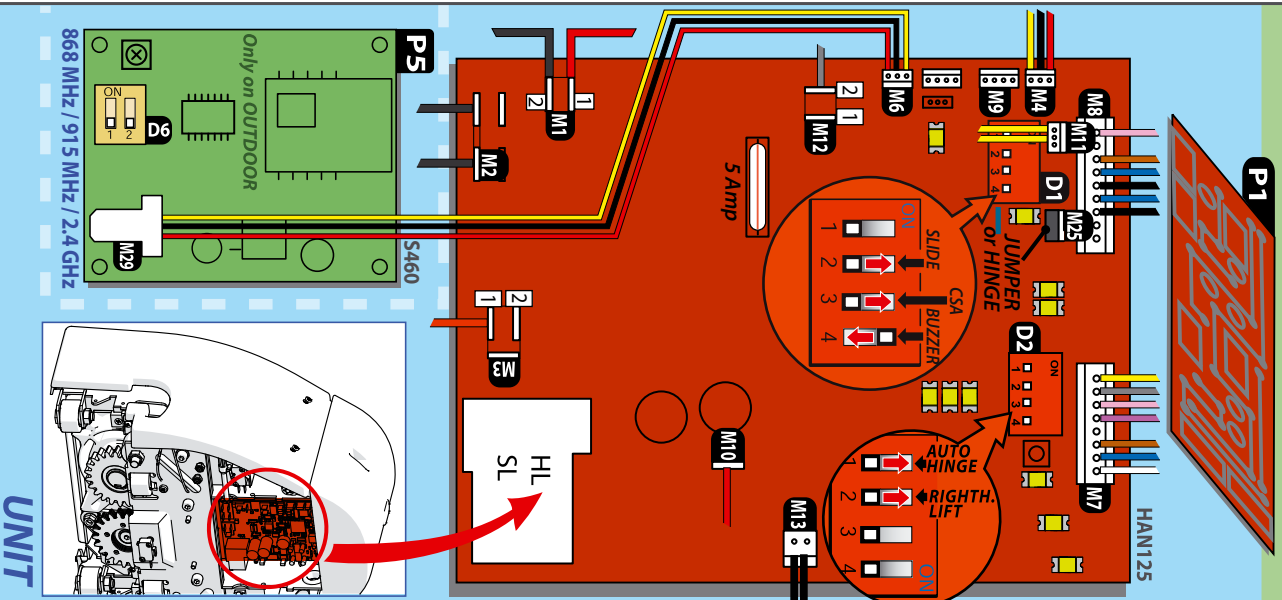
Savaria Lifts (UK) Ltd

82 First Avenue,
Pensnett Estate Kingswinford,
West Midlands DY6 7FJ,
England T +44 (0)1384 408700
enquiries@handicare.co.uk

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Savaria Canada

2 Walker Drive Brampton,
ON L6T 5E1 Canada
T 1.905.791.5555
T 1.855.728.2742
info@savaria.com

www.handicare .com

Date	Version	Change
17-11-24	0N	Website links and QR codes
15-12-24	0O	Chassisbolts
30-01-25	0P	Position excentric washer
06-02-25	0Q	Hinged and H125 PCB updated
26-07-25	0R	Add extended foot/add XL plates/ change attachment symbols
17-11-25	0S	Stair inclination change and Inspection item 3.3 changed
26-11-25	0T	Text correction page 41/43
29-01-26	0U	Installation forms
10-06-26	0V	Preinstalled connections PowerOpt PCB
19-08-26	10	Grease toothrack every 30 cm

INSTALLATION MANUAL
INSTALLATIEHANDLEIDING
MANUAL D'INSTALLATION
INSTALLATIONANLEITUNG
MANUALE PER L'INSTALLAZIONE
MANUAL DE INSTALACION
MANUAL DE INSTALACAO
INSTALLASJONSMANUALEN
INSTALLATIONS MANUAL
INSTRUKCJA

AA30068_10_190826

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Inc. 10888 Metro Court