
FREECURVE 4* | INSTALLATION MANUAL

**early european adoption program*

Classic
Elegance
Alliance

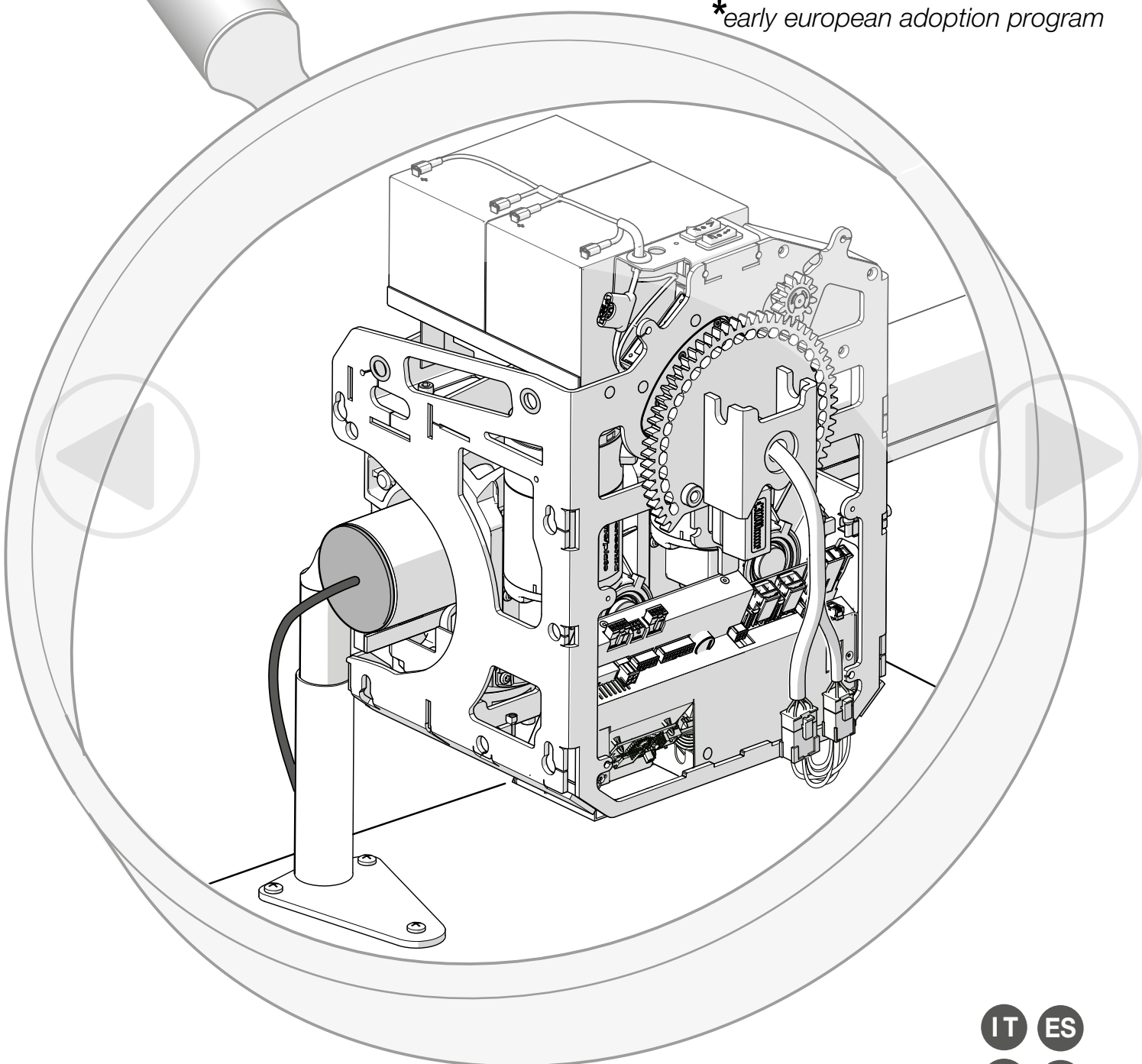 **savaria**®

EN All rights reserved. No part of this publication may be reproduced and/or published by means of print, photocopy, microfilm or in any other way or form, without prior written permission from the manufacturer. This also applies for the accompanying drawings and diagrams. The manufacturer reserves the right to change parts at any time, without prior or direct notification to the customer. The contents of this manual may also be changed without prior warning. This manual is intended for the standard models of the lift. The manufacturer therefore cannot be held liable for any damage resulting from specifications that differ from those of the standard model of the lift supplied to you.

This manual has been compiled with the best intentions and every possible care but the manufacturer cannot accept responsibility for any errors in this manual or for the consequences thereof.

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EN

The Dutch manual is the original language

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
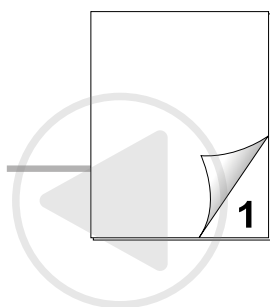
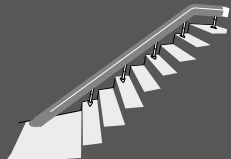
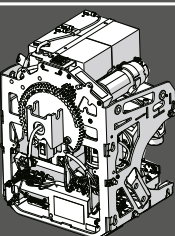
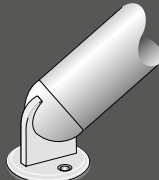
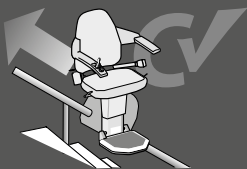
1  2	Thank you for choosing a Handicare stairlift. This manual is for the Handicare FC4 lift with a Classic or Elegance/Alliance seat type. The manual includes instructions for the standard and the optional features.
2  6	Inspect the stairs [p.06] First charge contact [p.08] Positioning charger [p.10] Stanchions [p.12] Bottom rail section [p.14] Other railparts [p.16] Joint [p.18] Fasten joint [p.20] Other charge contacts [p.22] Terminate power supply [p.24] Charger [p.26] Labels [p.28]
3  30	Locking plate [p.30] Frontcover [p.32] Batteries [p.34] Shortening links [p.36] Data input [p.38] Store USB [p.40] LED colours [p.42] Mapping on computer [p.44] Calibrate [p.46] Batteries remote control [p.48] Linking remote controls [p.50] Charge points [p.52] Installing covers [p.54] Fasten and check covers [p.56]
4  56	Mounting frame [p.56] Footrest loom [p.58] Front cover [p.60] Seatbelt [p.62] Manual swivel [p.64] Powered swivel [p.66]
5  68	Mounting remote controls [p.68] Testride [p.70] Fix stanchions [p.72] Grease toothrack [p.74]
6  76	Dimensions [p.76] Safety belt [p.78] Safety detection [p.80] Emergency stop [p.82] Armrest safety [p.84] Instruct customer [p.86] Fill in forms [p.88] Load label [p.90]
In this manual, the paperclip symbol refers to a relevant attachment.  : A  A  92	A - Hinged rail [p.94] B - A-rails [p.104] C - Remote controls [p.108] D - Wall support [p.110] E - Emergency operation [p.112] F - Footrest adjustment [p.118] G - Lockable switch [p.120] H - Active Seat [p.126] I - Torque Wrench Test [p.132] J - Upper stanchion [p.133] K - Protection Equipment [p.134] L - Competence [p.135] M - Installation forms [p.138] N - Electrical diagrams [p.154]

 : K{page 134}



Nominal	Torque	Torque
M 5	5,5 Nm	48.7 Lb-In
M 6	9,5 Nm	84.1 Lb-In
M 8	20 Nm	177 Lb-In
M 10	46 Nm	407.1 Lb-In

 **Nm** *Hand tight*
Not loose, but not really tight either.
Somewhere around 1 or 2 Nm.

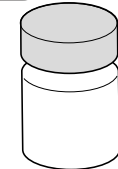Art.nr. 11J00114
 (Interflon Fin Grease)

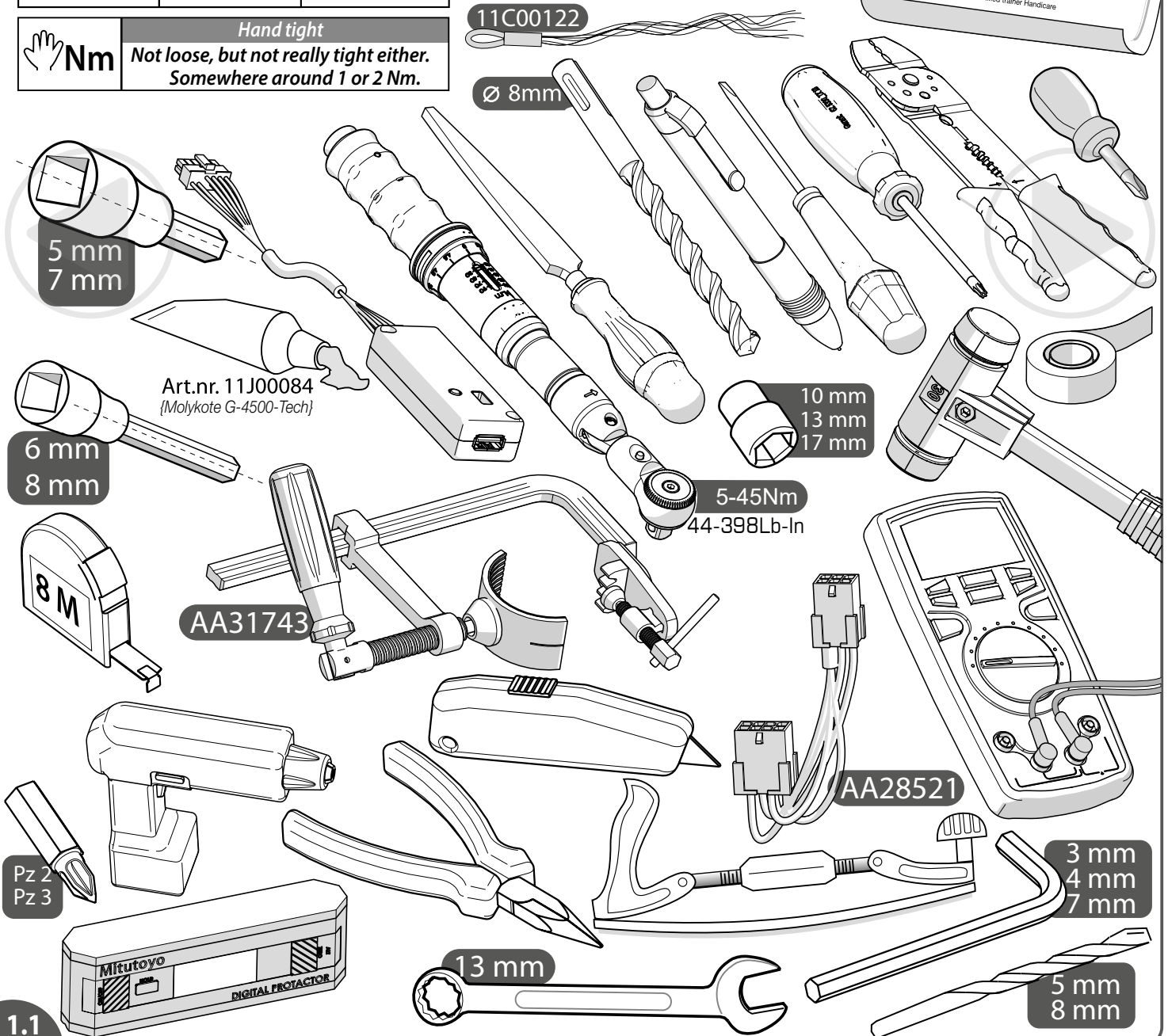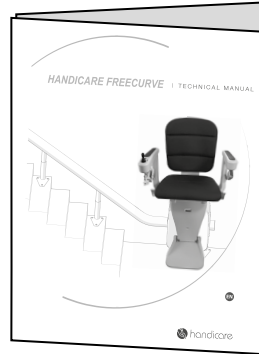Only a **certified** technician should install this lift (for more detailed information about the required competence see attachment of this manual). General required knowledge can be found in the Freecurve Technical Manual. The latest instructions of the lift and instructions for options are available for download at www.handicarepartners.com. Ensure good and sufficient tools; the shown tools are minimum required to install the stairlift. In case of installation a new part on a existing lift always replace the part by a part with the same article number.

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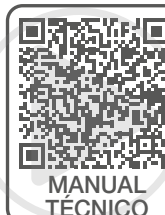
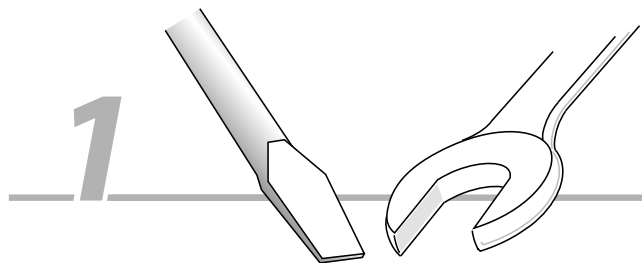
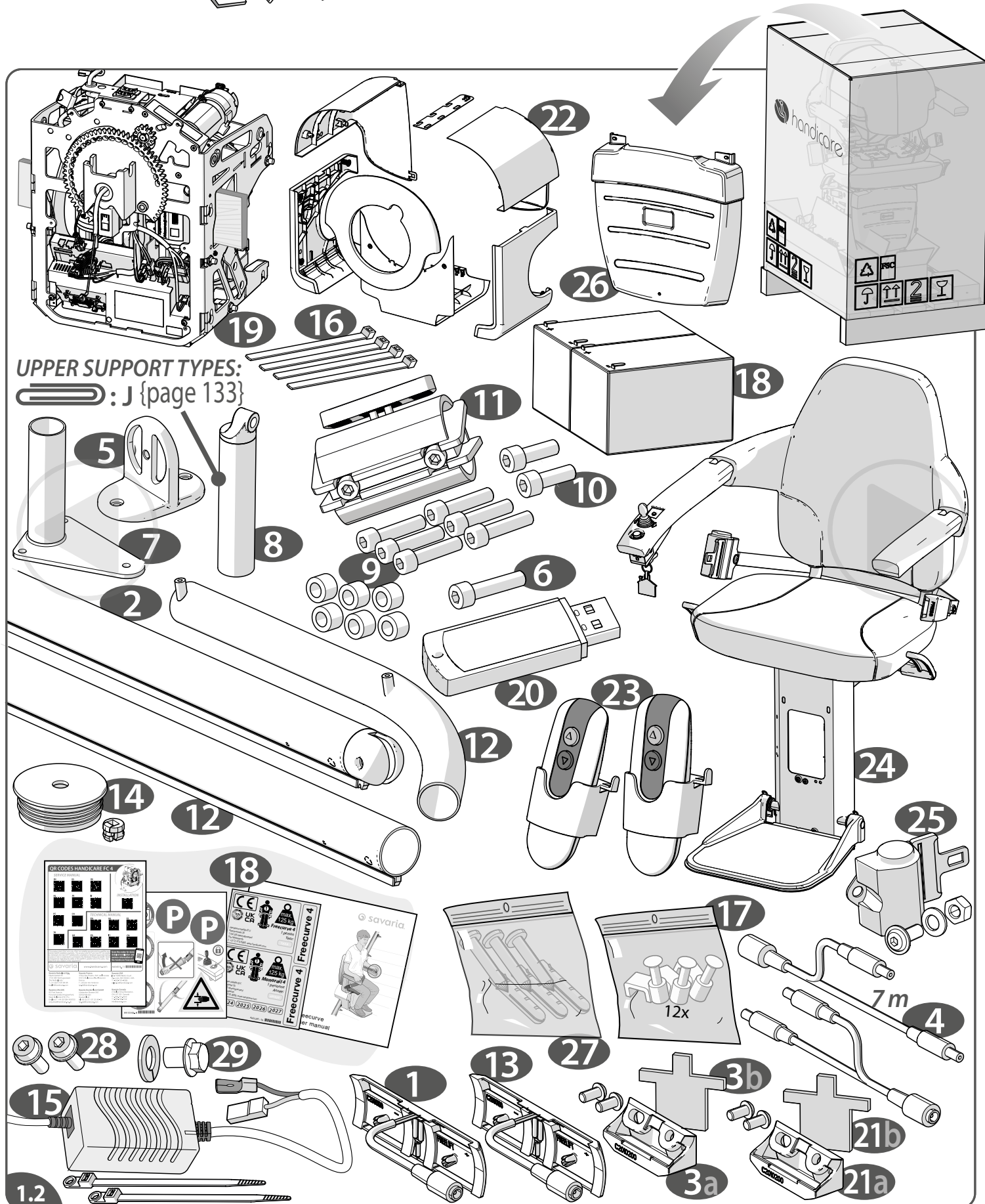
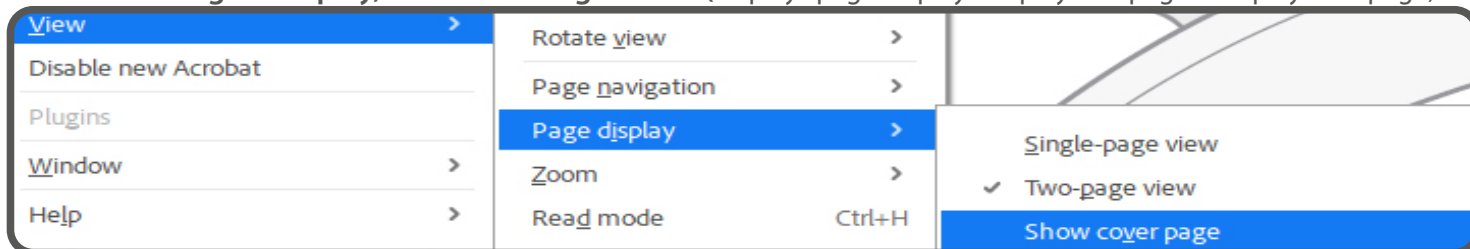
[1.2] Check that all the parts are present. The part numbers shown next to the parts are given in the installation manual at the time of installation.

When replacing parts on an existing stairlift, always use parts with identical part numbers.

EN **Hand tight** = not loose, but not really tight either. Somewhere around 1 or 2 Nm.

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


TIP

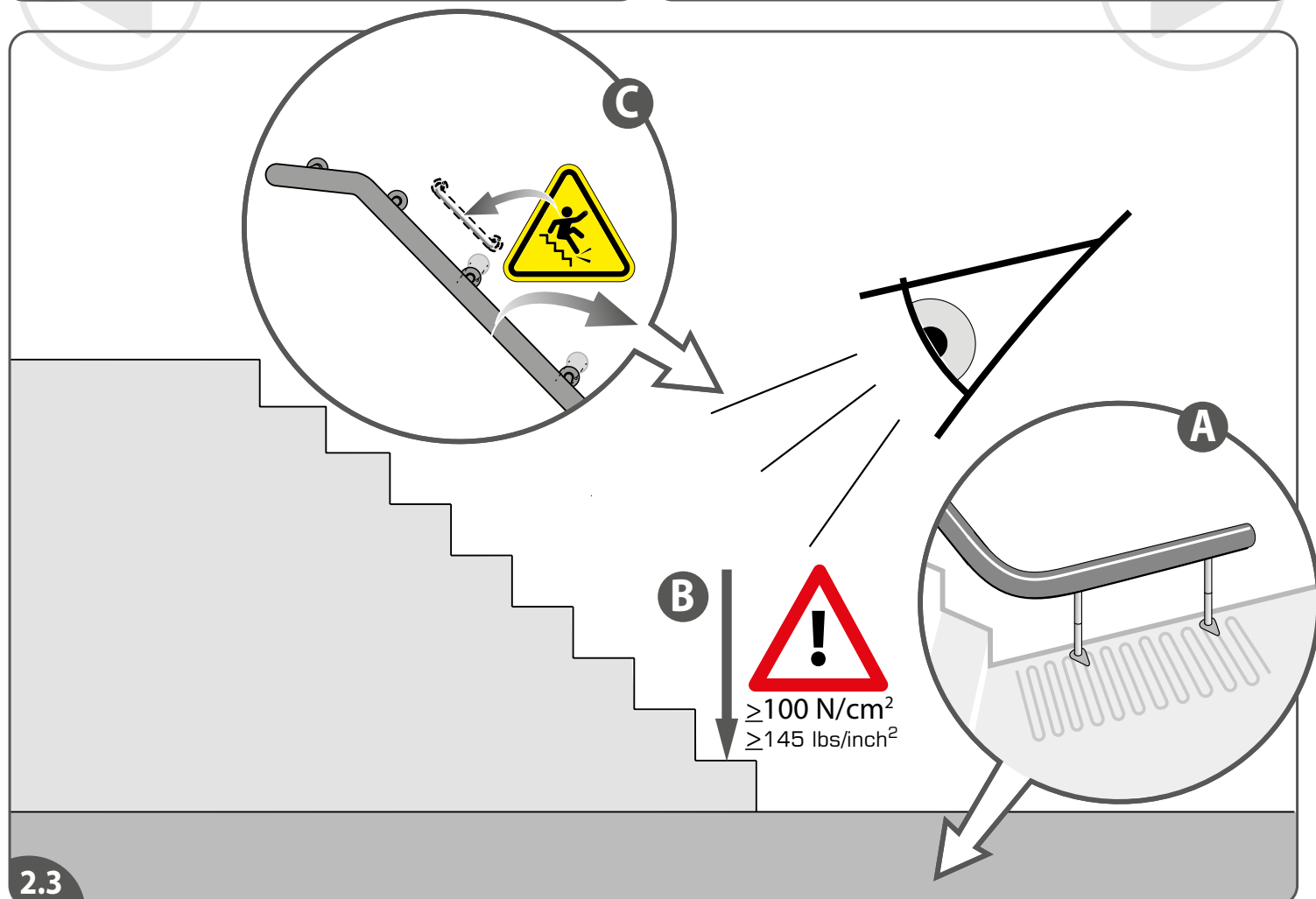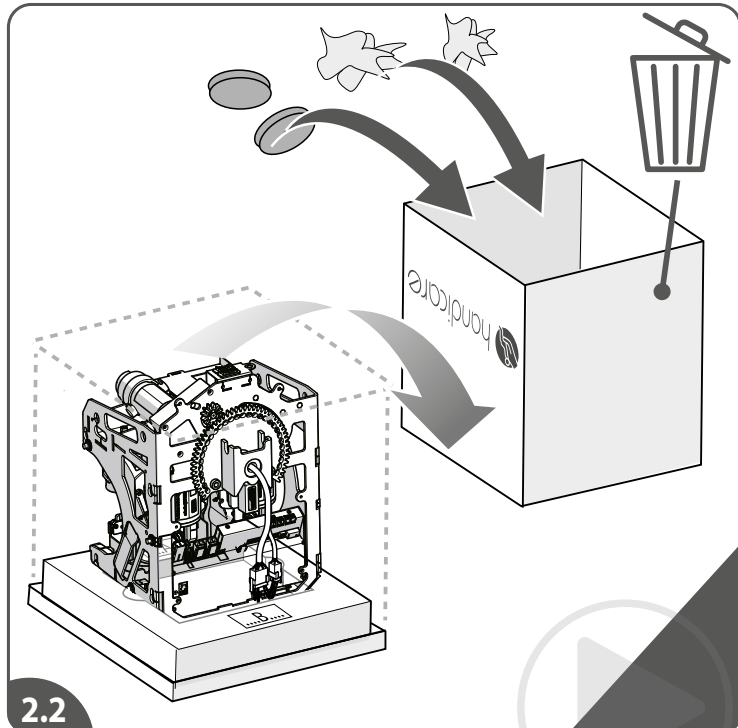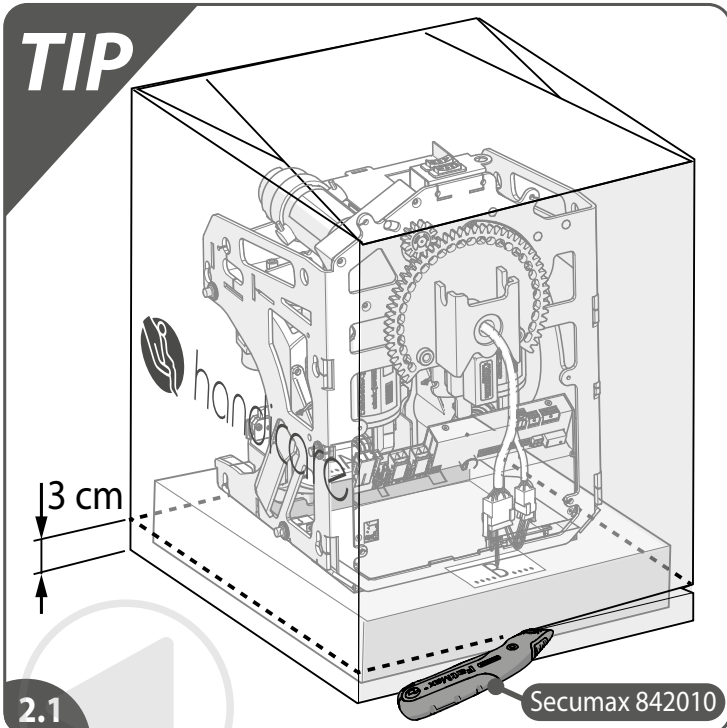[2.1/2.2] Cut the box at the bottom with a Secumax knife to remove the power pack.

[2.3] Before you start, inspect the stairs. Clarify with your customer regarding any floor heating, alarm wires, pads or pipework that may affect/hinder installation prior to starting{A}.

Check if the stairs has sufficient strenght{B}.

If it is necessary to remove the banister, place substitute equipment to ensure the customer will not fall{C}.

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




C0052050
C0052030

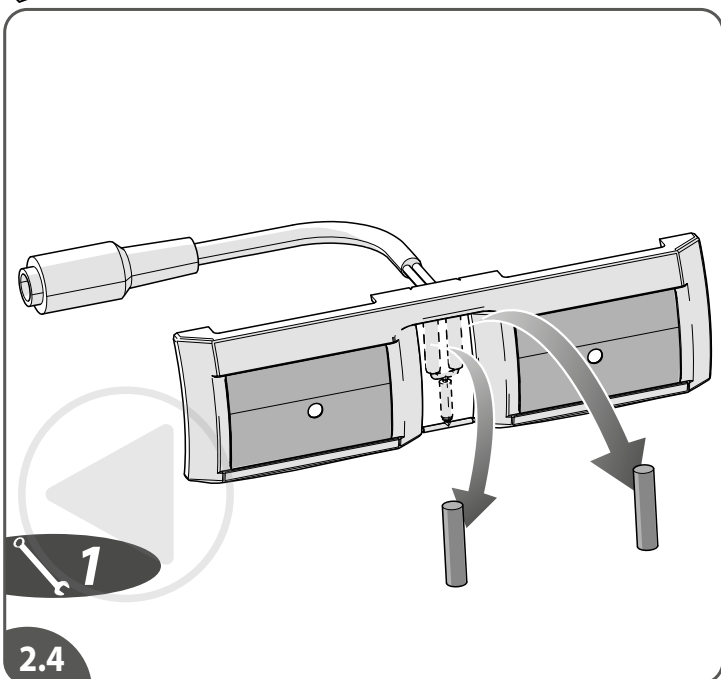INTERNAL CURVES ONLY

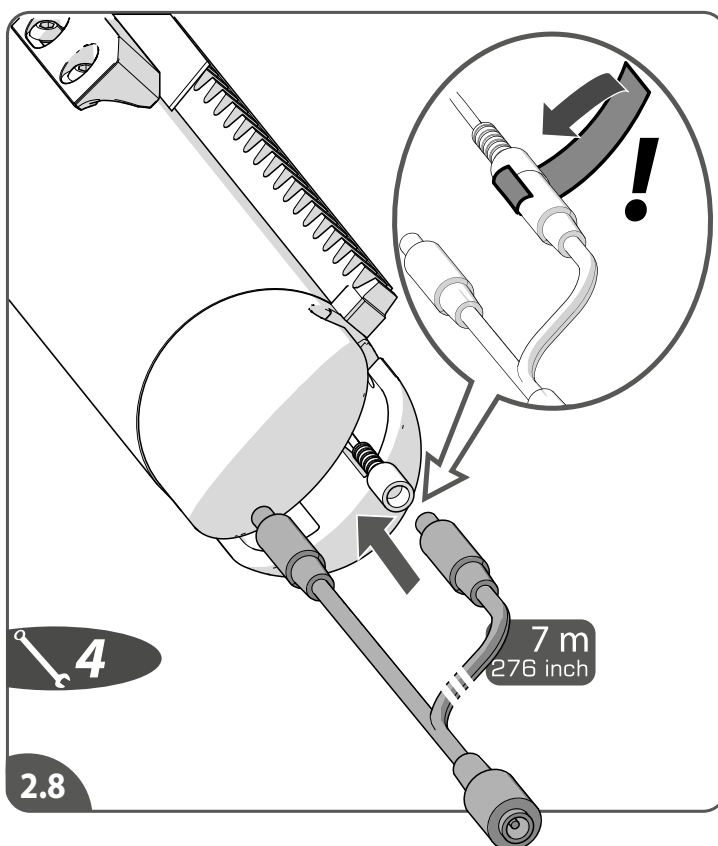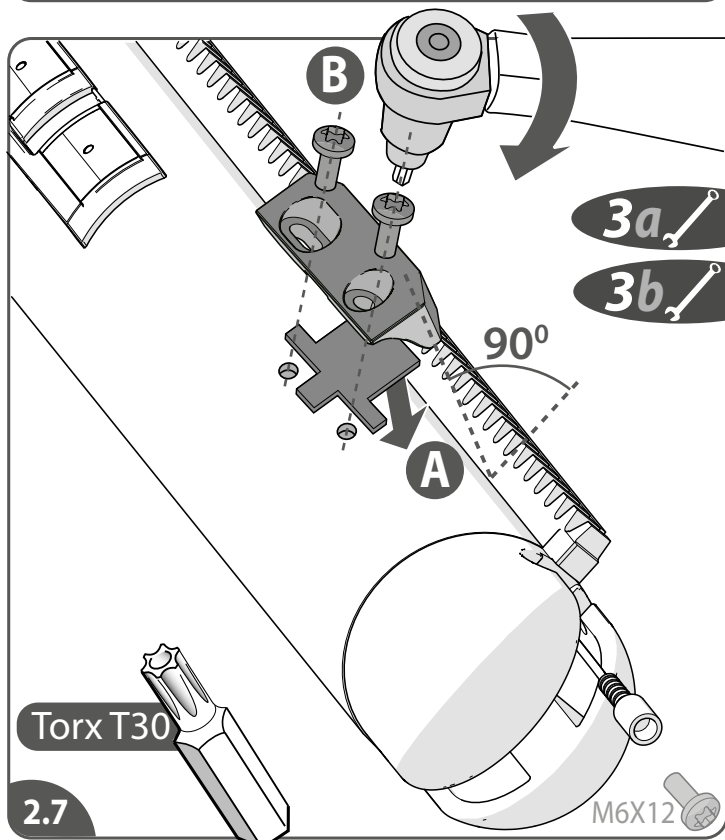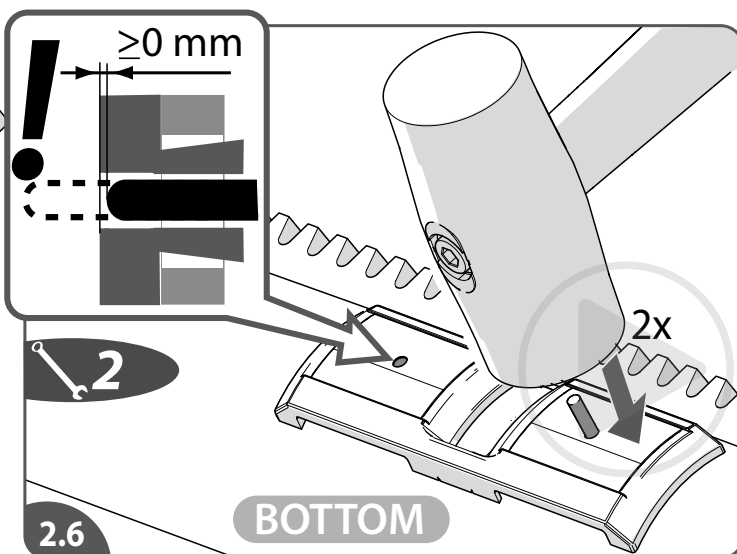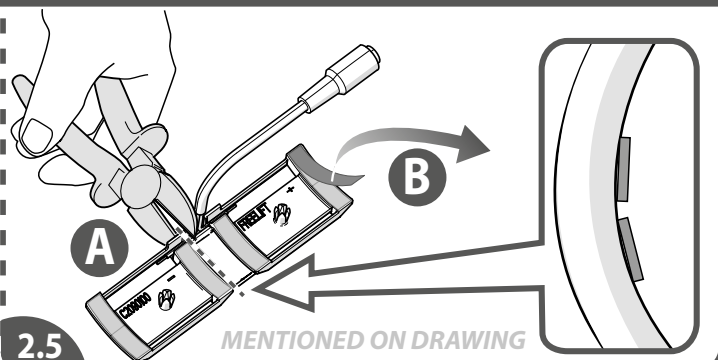If present, first install the hinged rail parts (see attachment A, page 94).

[2.4] Install the bottom charge contact: remove the pins.

[2.5] If the charge contact is in an internal curve, cut it in two parts{A}. Remove the tape{B}.

[2.6] Run the cable into the hole of the rail. Place the charge contact on the rail and fix it: gently hit the pins in the charge contact

[2.7] Position the bracket{A}. Install the endstop (90°) and fasten the endstop{B}.

Note: There is no thread on the endstop holes; bolts are selftapping.

EN [2.8] Connect the charging wire. Secure the connection with tape.

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
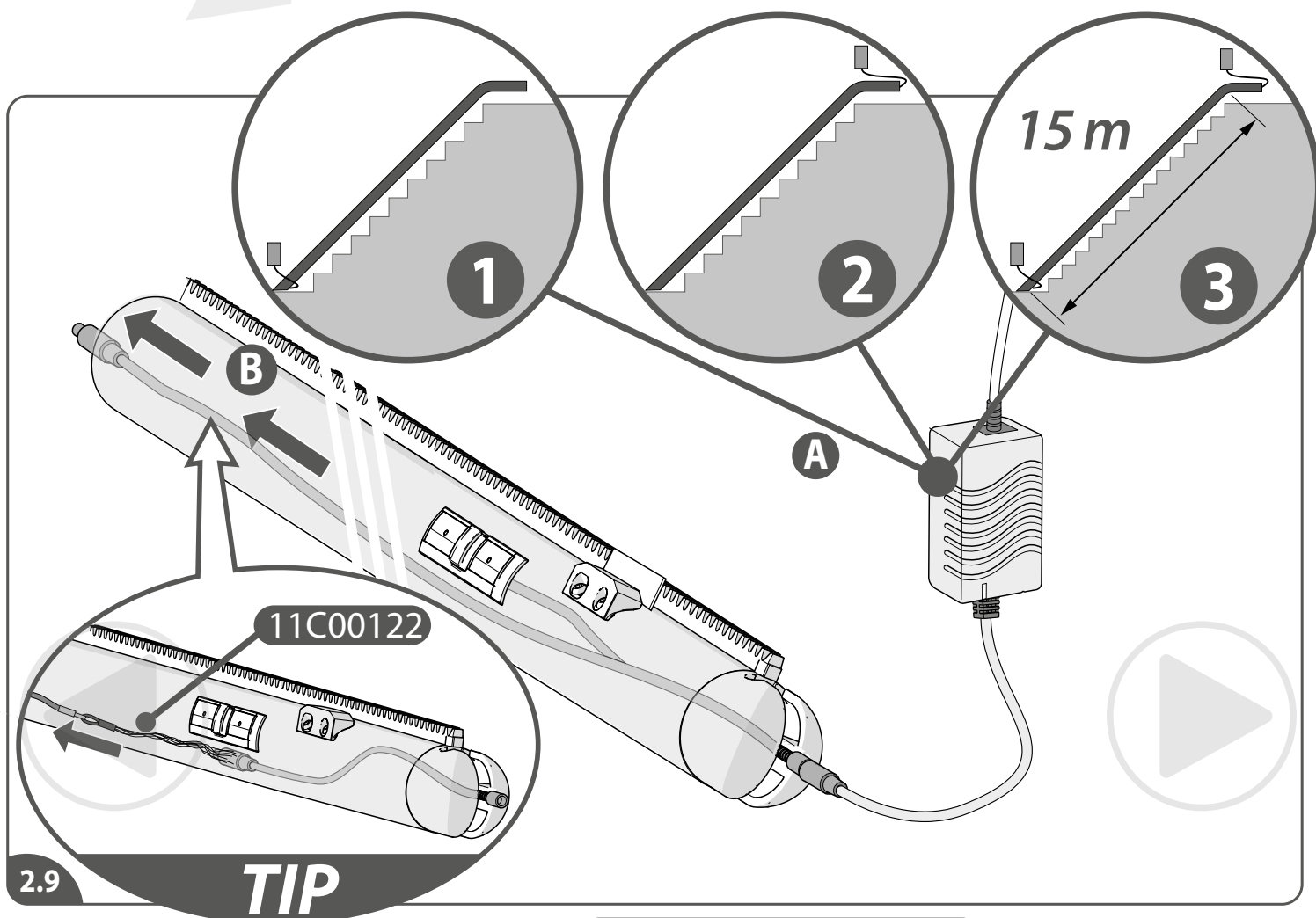
short

158/[168]

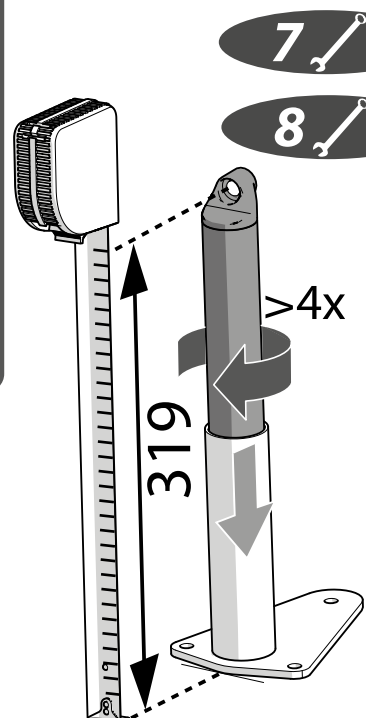[2.9] Decide where to install the power supply: depending on application this could be at the top {2} or bottom {1} of the stairs. Or both {3} in case the stairs is longer than 15 m. Run the charging wire through the rail {B}.
 [2.10] Check the drawing for the right stanchions.
 [2.11] Install the stanchion on the foot: turn around the stanchion at least 4 times. Keep turning till the right height is reached. Check the right height (can be find on the drawing).

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
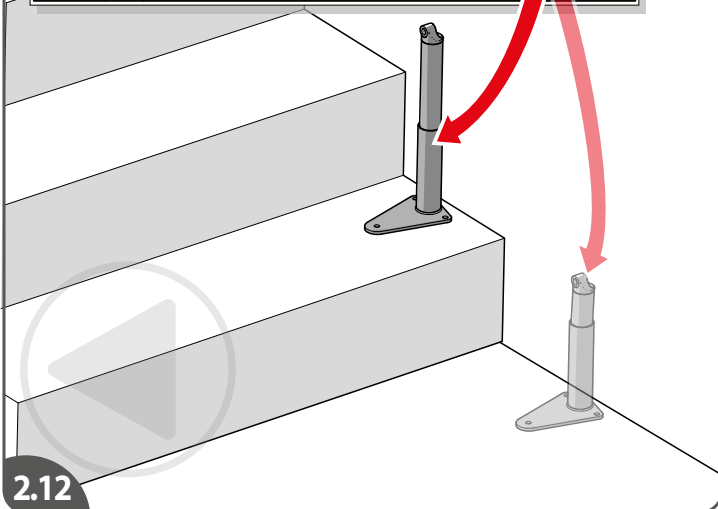
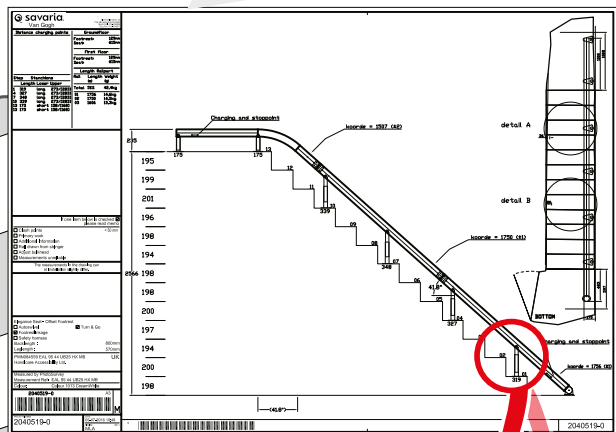
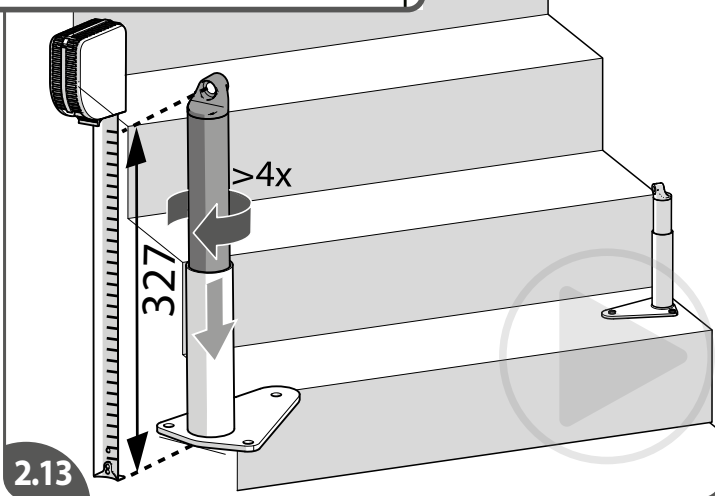
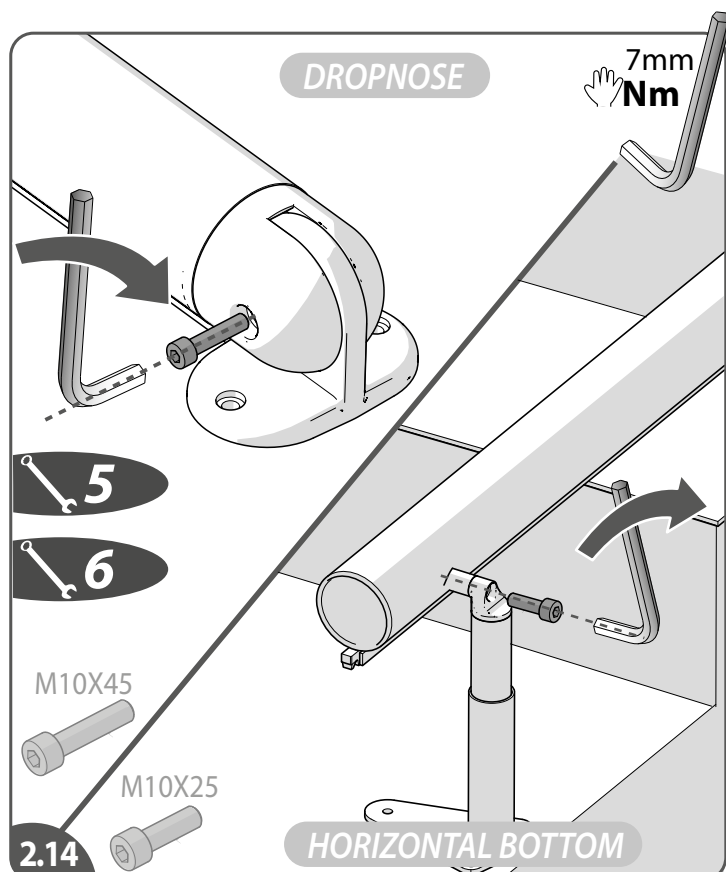
2.14

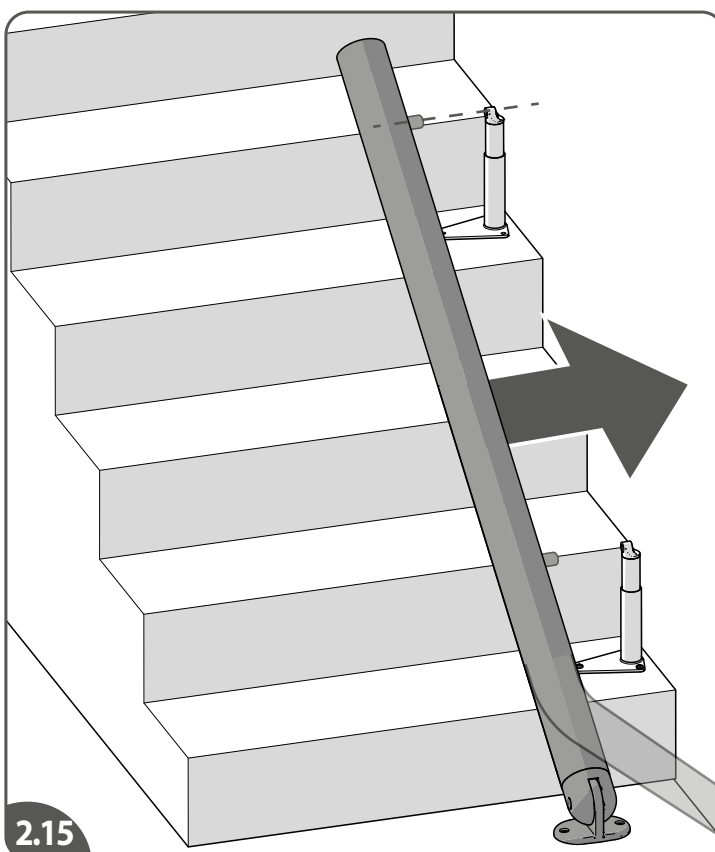2.15

- [2.12]Position the stanchion as shown in the drawing.
- [2.13]Install the next stanchion in an other foot: follow the same procedure and position the stanchion as shown in the drawing.
- [2.14]Depending on the railending, install the bracket or foot at the bottom rail section.
- [2.15]Position the bottom rail section.

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






IT

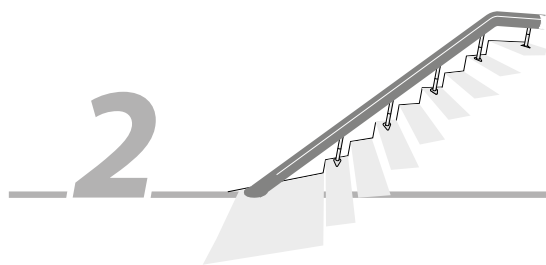RAIL BOTTOM RAIL SECTION

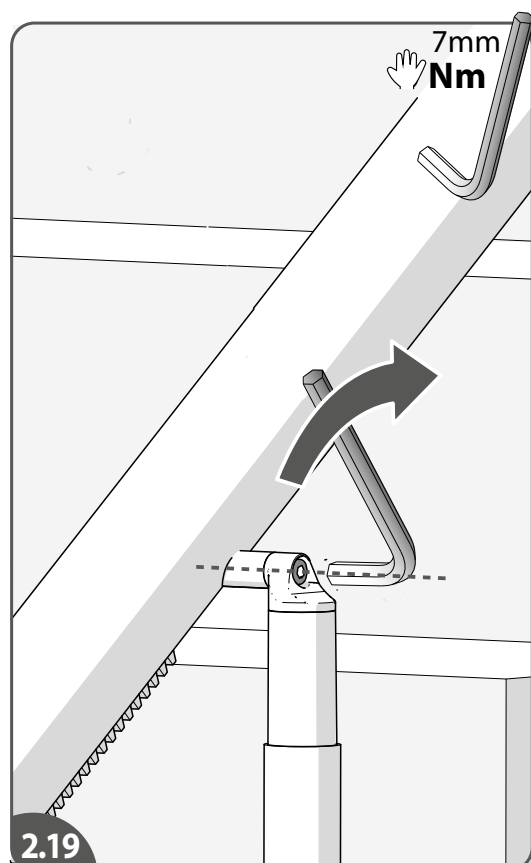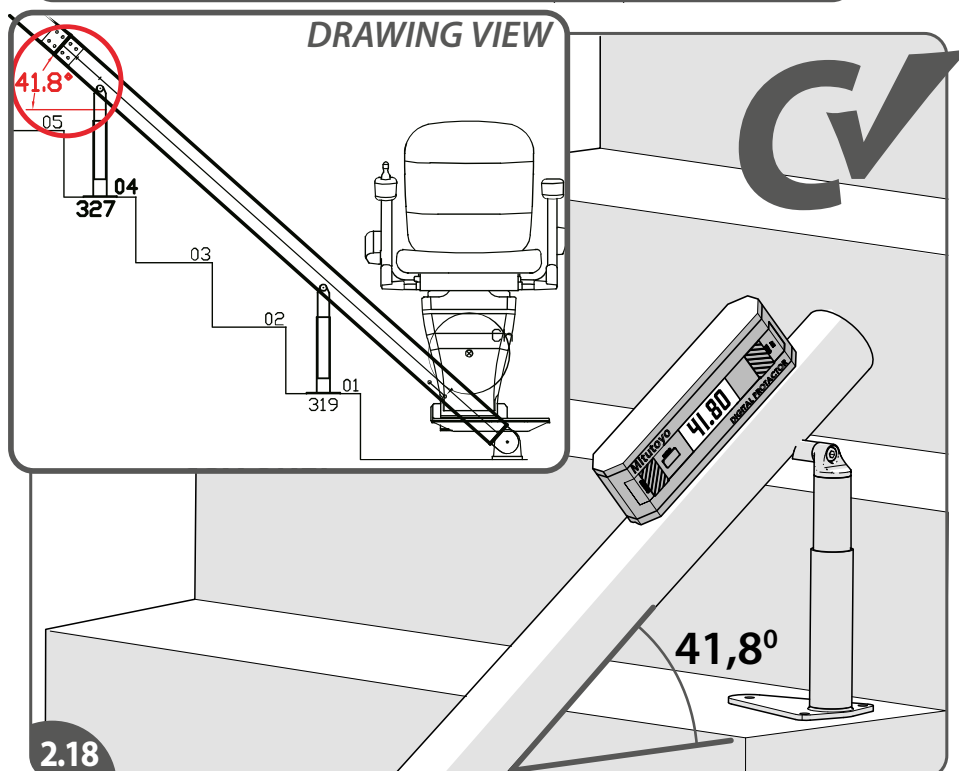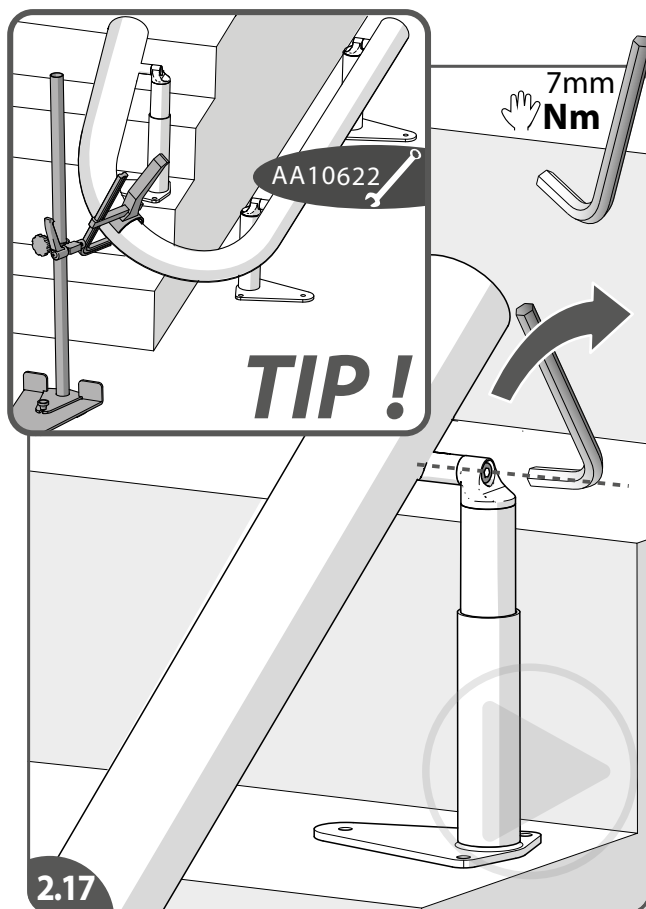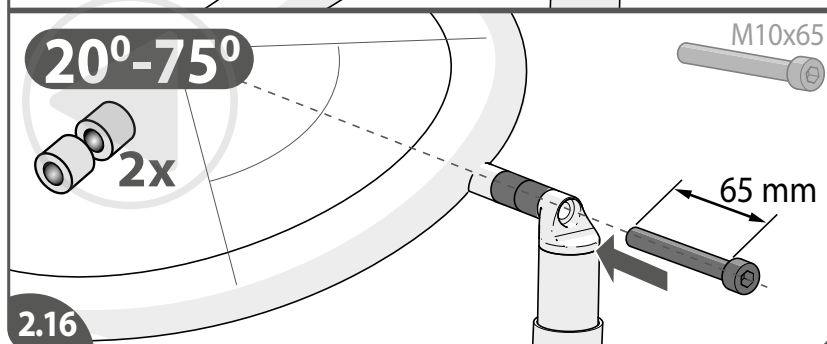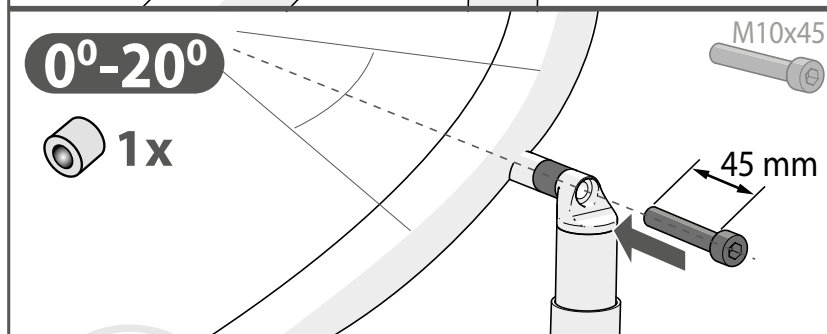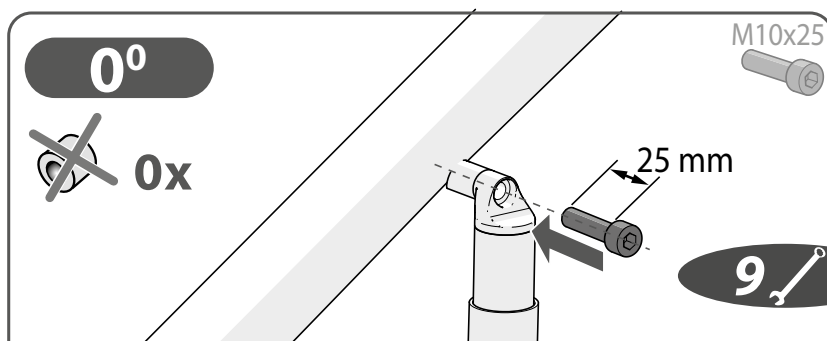[2.16] If needed install extra spacers between rail and stanchion (external curves only).

[2.17] Fasten the last stanchion to the first rail section [hand tight]. After the complete track installation the torque of the stanchion-rail bolts need to be set on 10 Nm (88 Lb-In).

TIP: when a rail part does not remain stable during assembly, use the **installation assistant**, available on www.handicarepartners.com

[2.18] Compare the rail angle with the mentioned angle in the drawing. If the actual rail angle is different to the rail drawing then adjust the heights of the stanchions accordingly, the angle of the rail is the critical measurement.

EN

[2.19] Fasten the lower stanchion(s) to the rail [hand tight].

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



IT

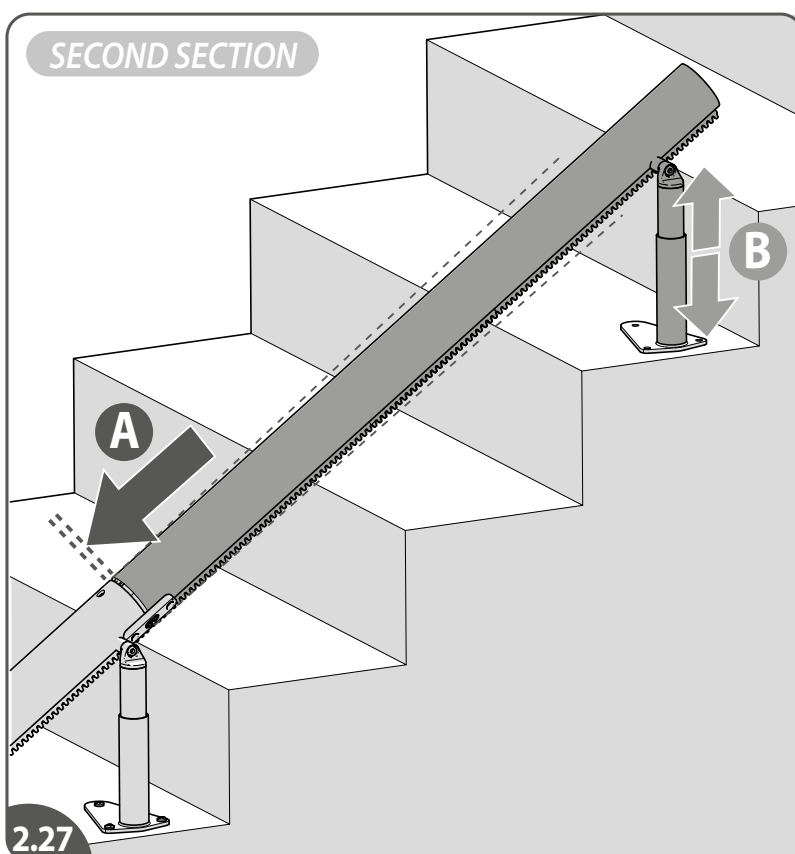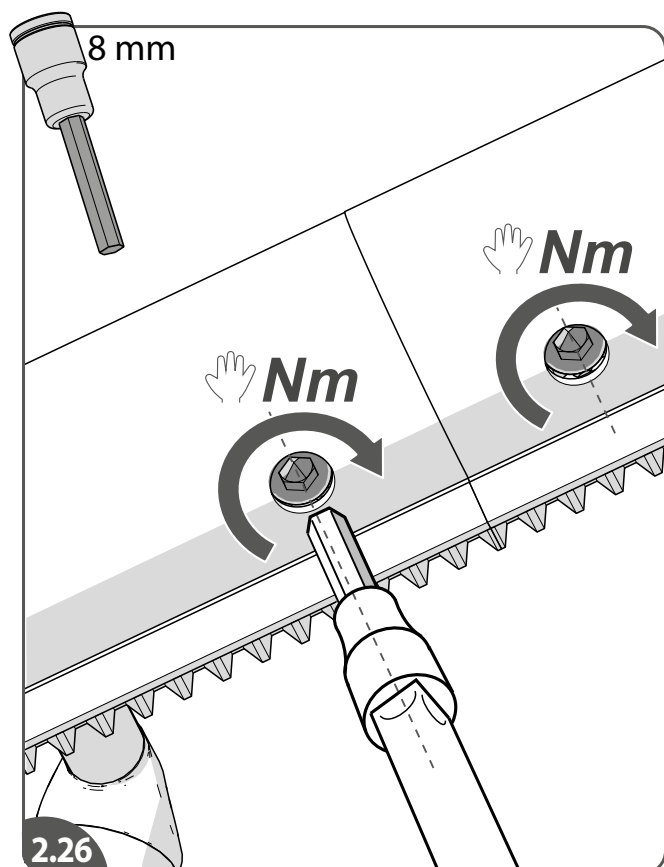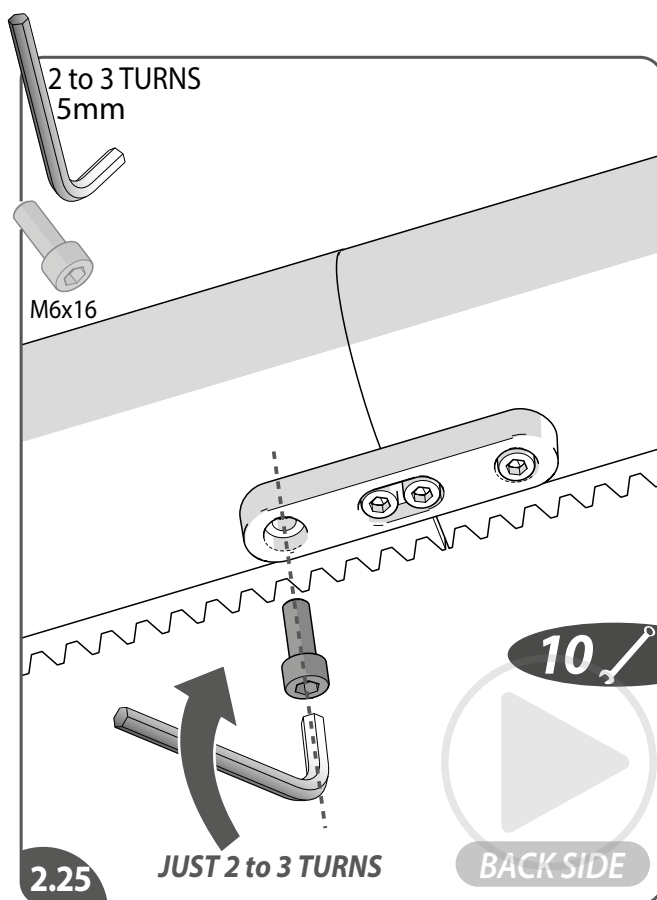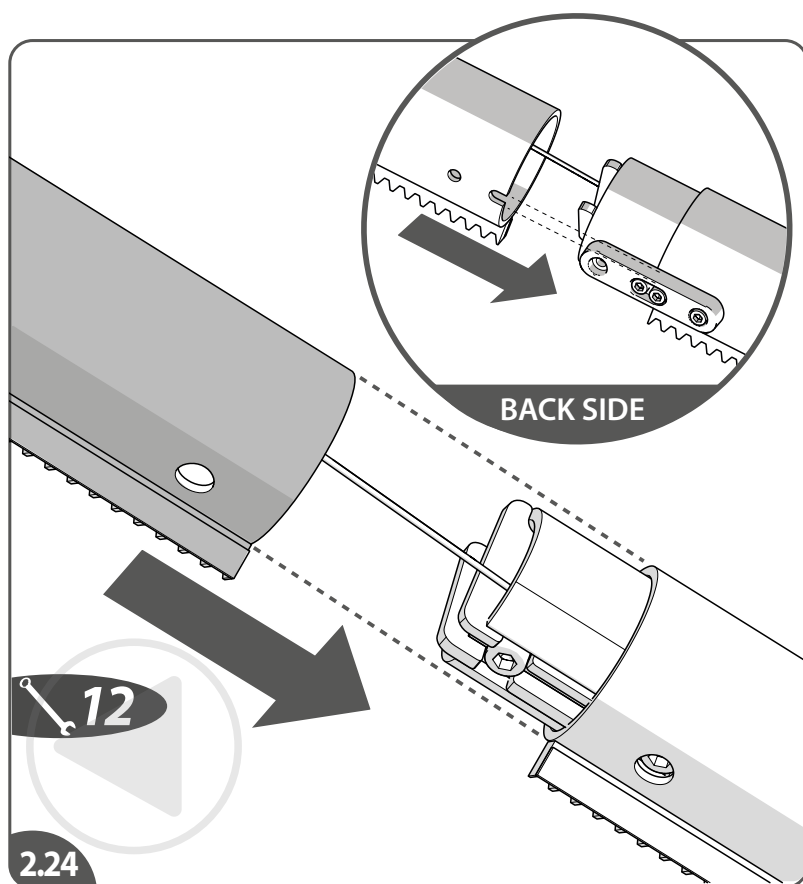- [2.24] Guide the charge wire and place the next rail part.
- [2.25] Insert the last fixation bolt (just 2 to 3 turns) in the joint.
- [2.26] Start to fix the joint.
- [2.27] Close the cap between the railparts{A} by adjusting the stanchions{B}.

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





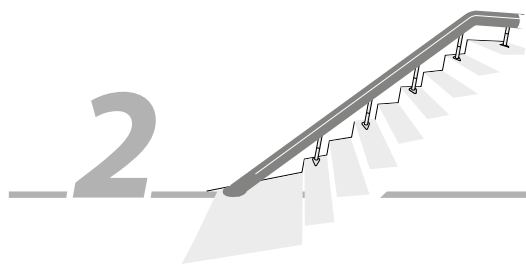
RAIL

FASTEN JOINT

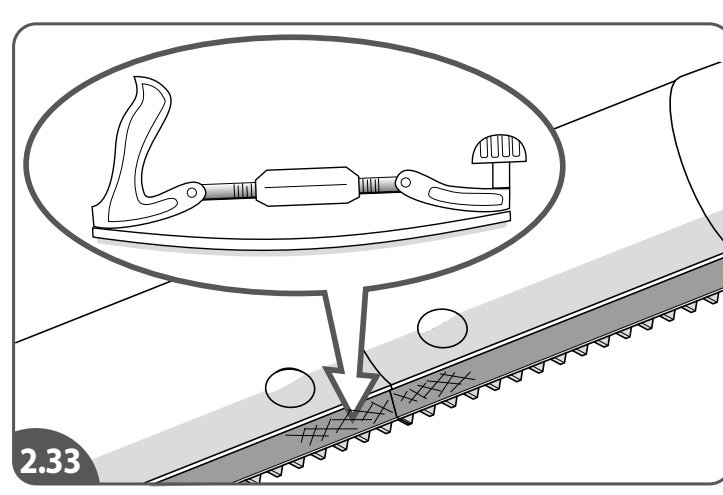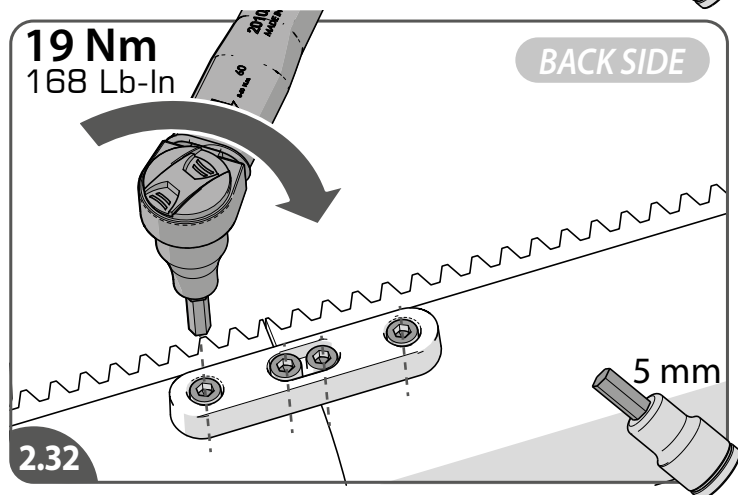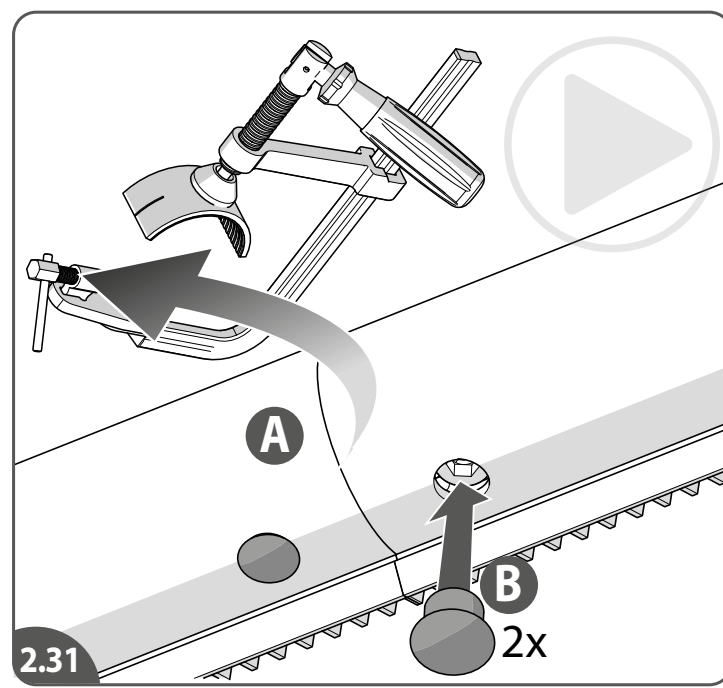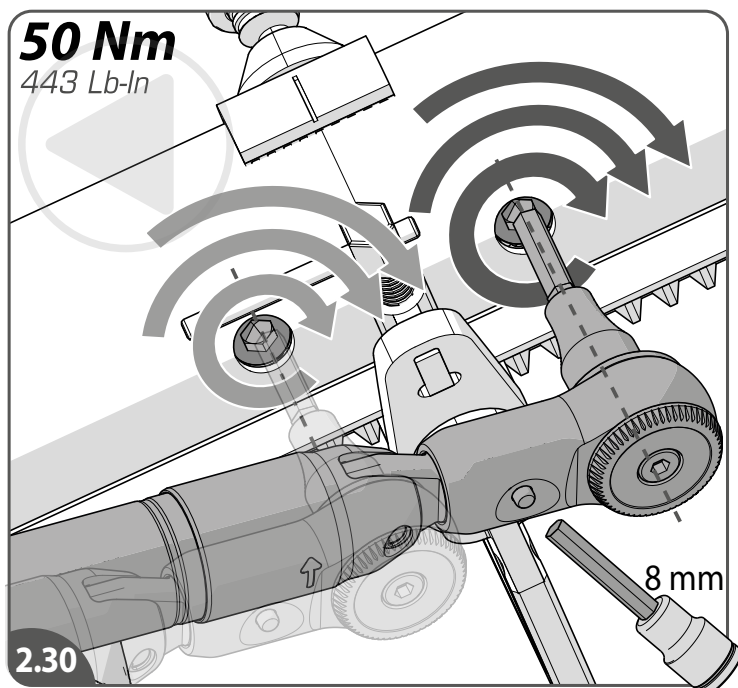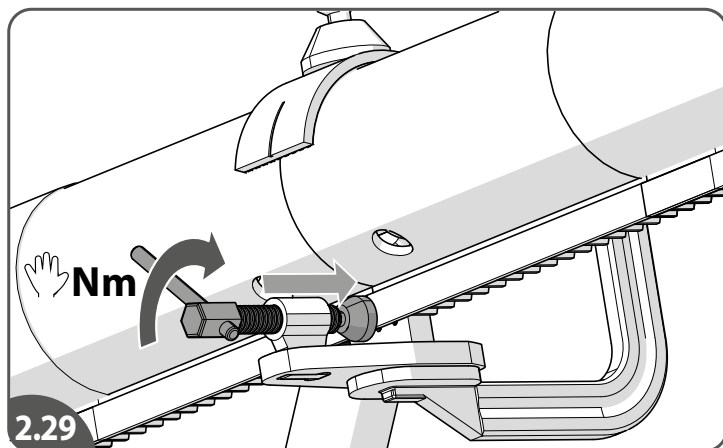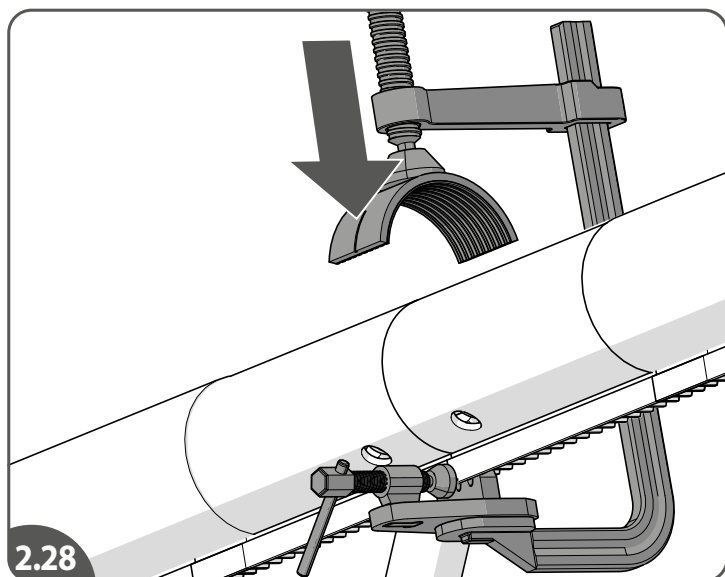- [2.28] Install the tool (*Handiclamp*) to align the toothtracks of the rail parts.
- [2.29] Align the toothtracks of the rail parts.
- [2.30] Torque it on alternating sides to 50 Nm (442 Lb-In).
- [2.31] Remove the toothrack tool{A}. Put the caps in the holes{B}.
- [2.32] Tighten the bolts on the backside.
- [2.33] Smooth the domex strip with the body file.

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





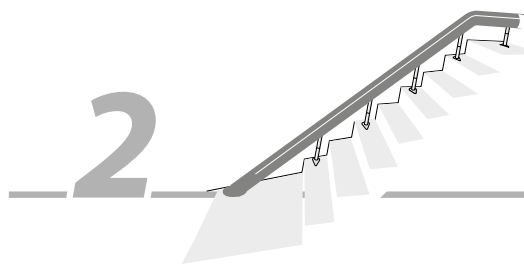
RAIL

OTHER CHARGE CONTACTS

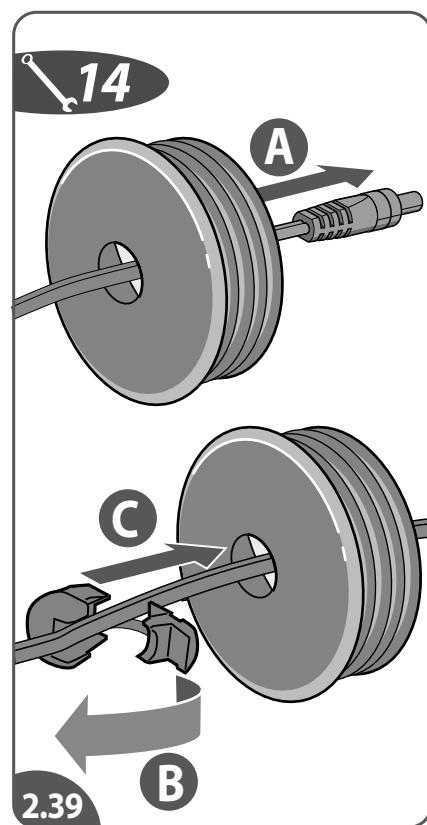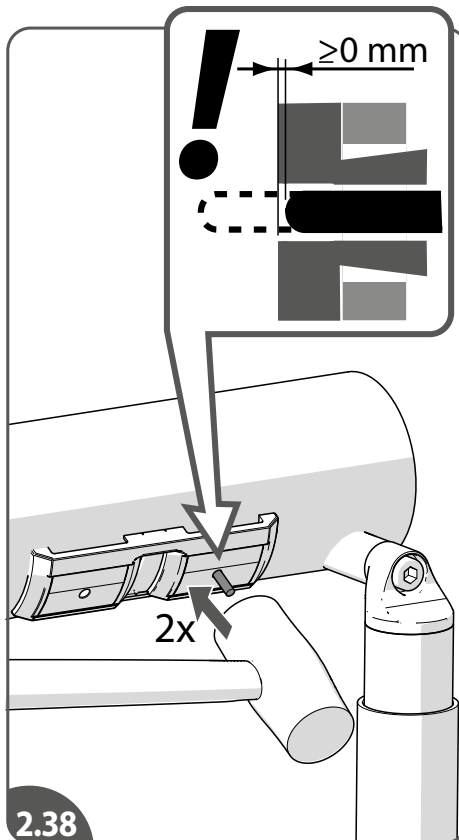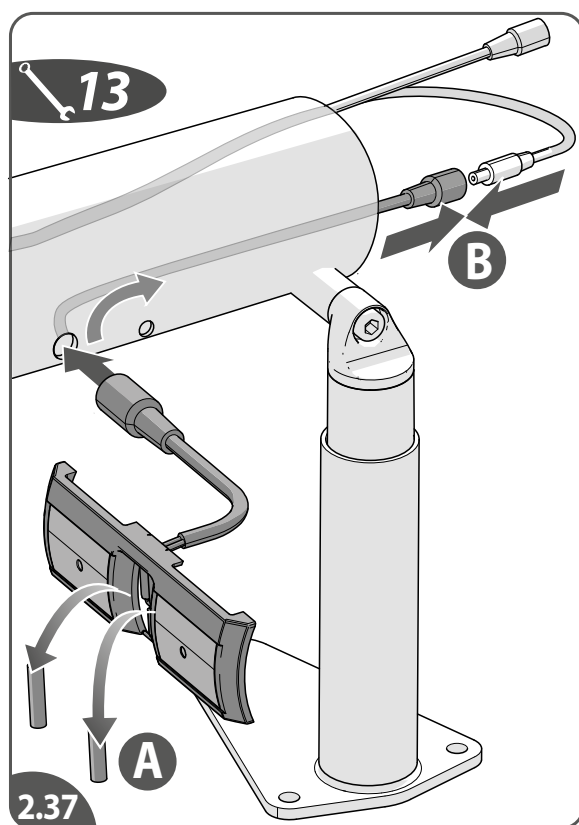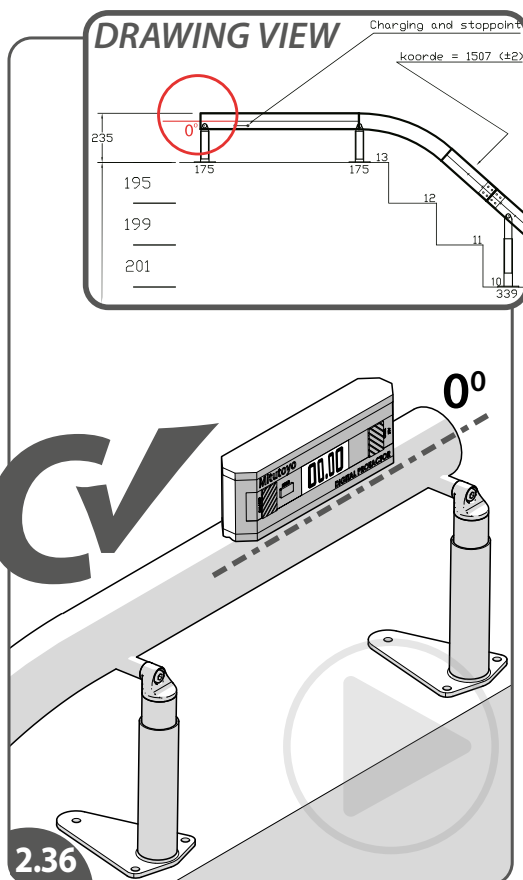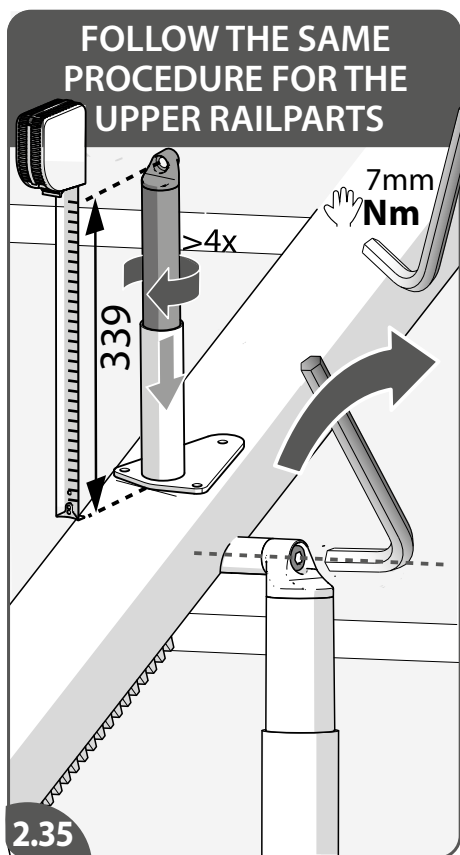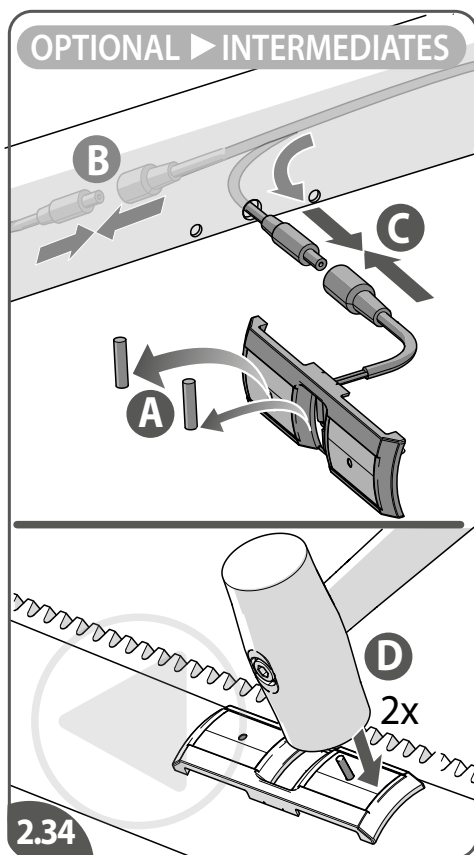[2.34] In case of optional intermediate charge contacts, connect the extra charging cable{B}. Connect the intermediate charge contact{C} and install charge contact on the rail{D}.
 [2.35] Follow the same procedure to install the remaining rail parts.
 [2.36] Compare the rail angle with the mentioned angle.
 [2.37] Remove the pins from the charge contact{A}. Connect the cable of the charge contact{B}.
 [2.38] Place the charge contact on the rail and fix it. Gently hit the pins in the charge contact.
 [2.39] Guide the cable of the charger through the hole of the endcover{A}. Close the clip{B} and press it in the endcover{C}.

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
2 *RAIL*

TERMINATE THE POWER SUPPLY

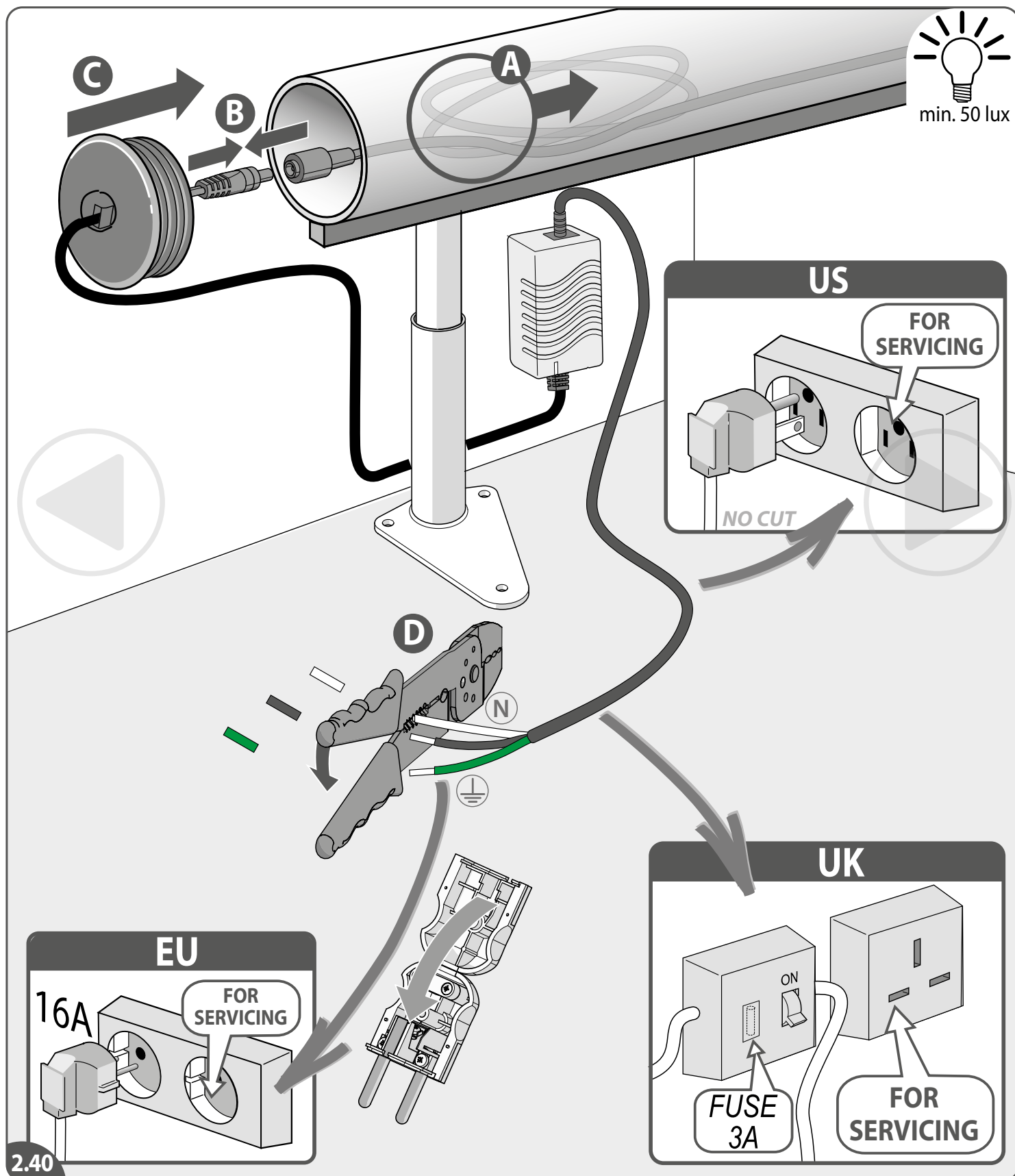[2.40]Fold up the excess cable and insert it into the rail {A}. Connect the charger to the charging wire{B}. Press the endcover in de rail{C}. Terminate the power supply to the ac supply{D}. All mains wiring must comply with the local standards in force at time of installation and with section 4.3 of EN60204/1 and also with the requirements of NEC and ANSI/NFPA70. 3-pin plug not supplied.

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




IT

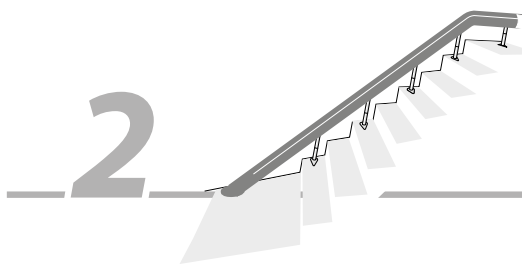RAIL CHARGER

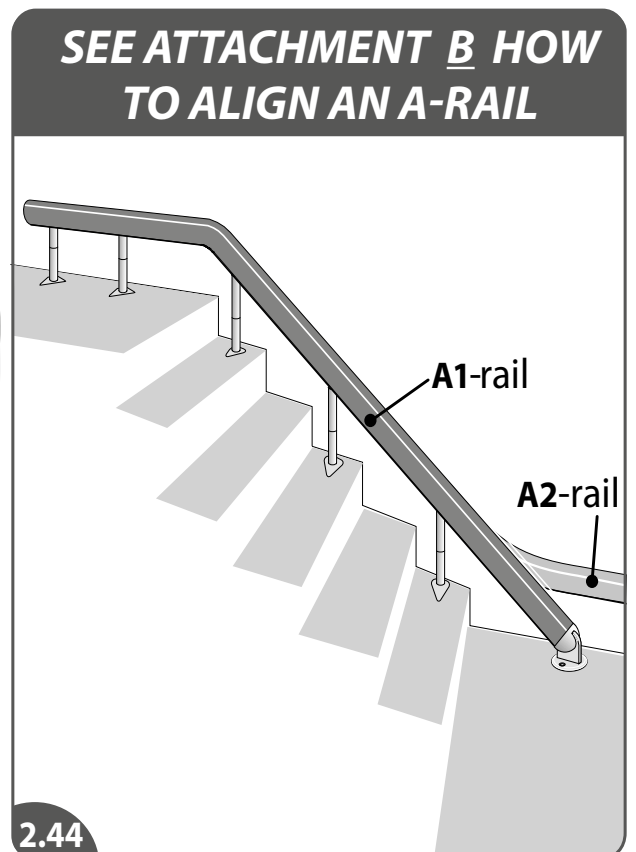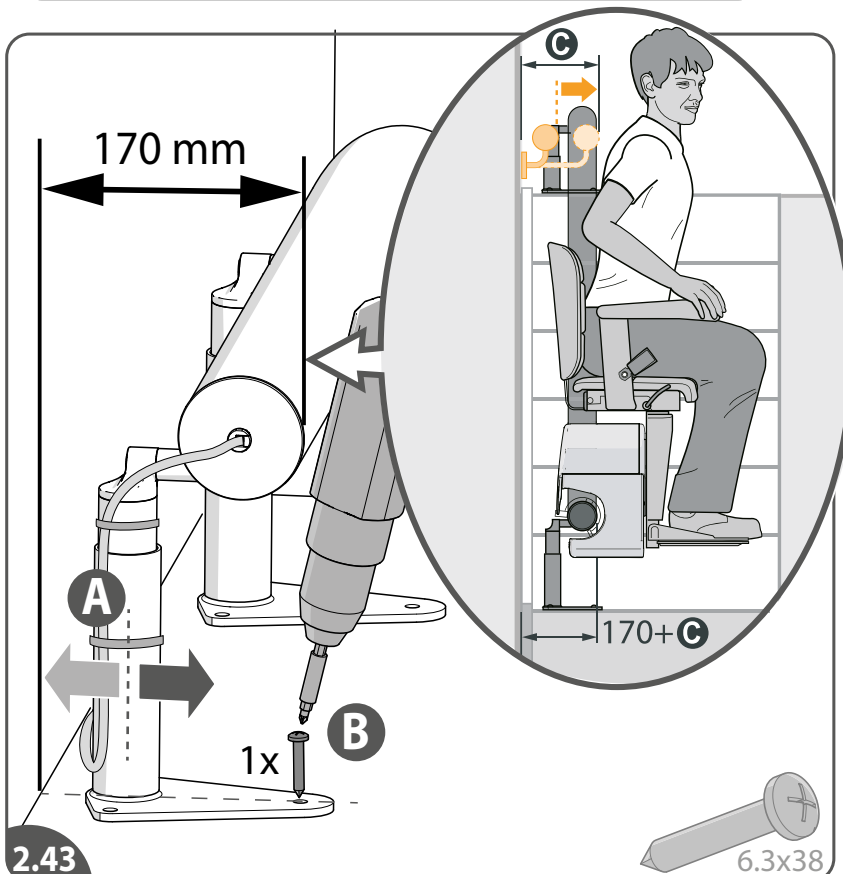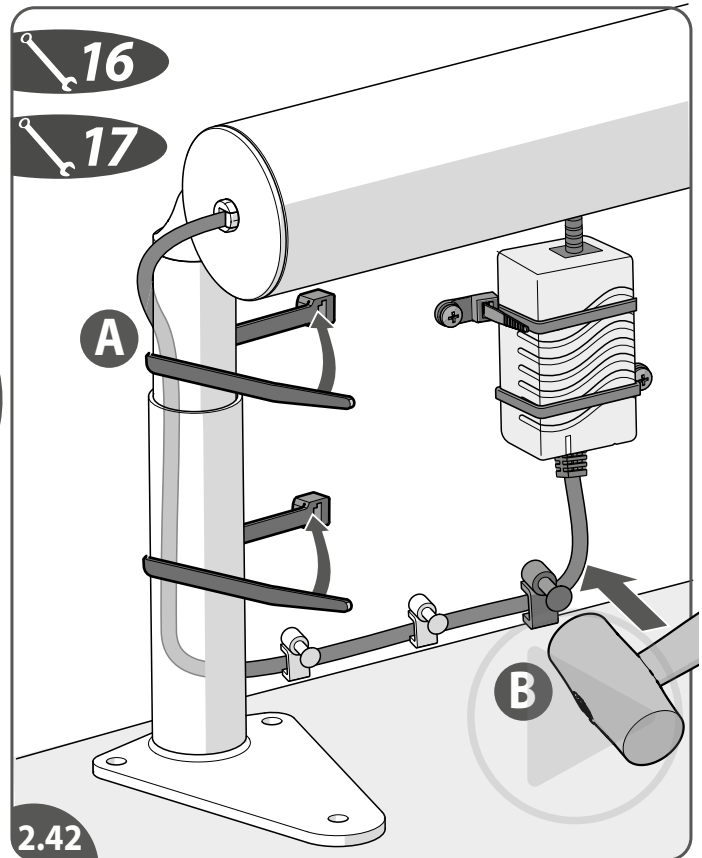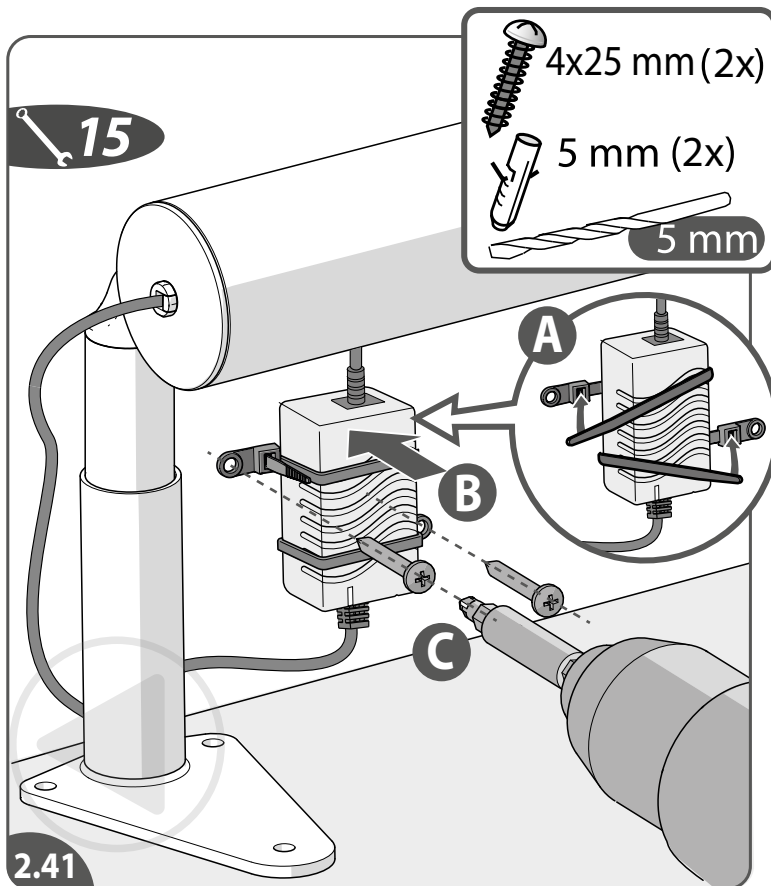- [2.41] Position the charger on the wall and fasten it with tiewraps and screws.
 [2.42] Fasten the cable of the charger to the stanchion{A}. Fasten the cable of the charger to the wall{B}.
 [2.43] Position the rail to the wall{A}. If there is an obstacle increase the distance to the wall. If the rail is not stable, fix some stanchions with a few screws{B}.
 [2.44] In case of an **A**-rail align the track. See attachment **B** how to do this.

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



IT

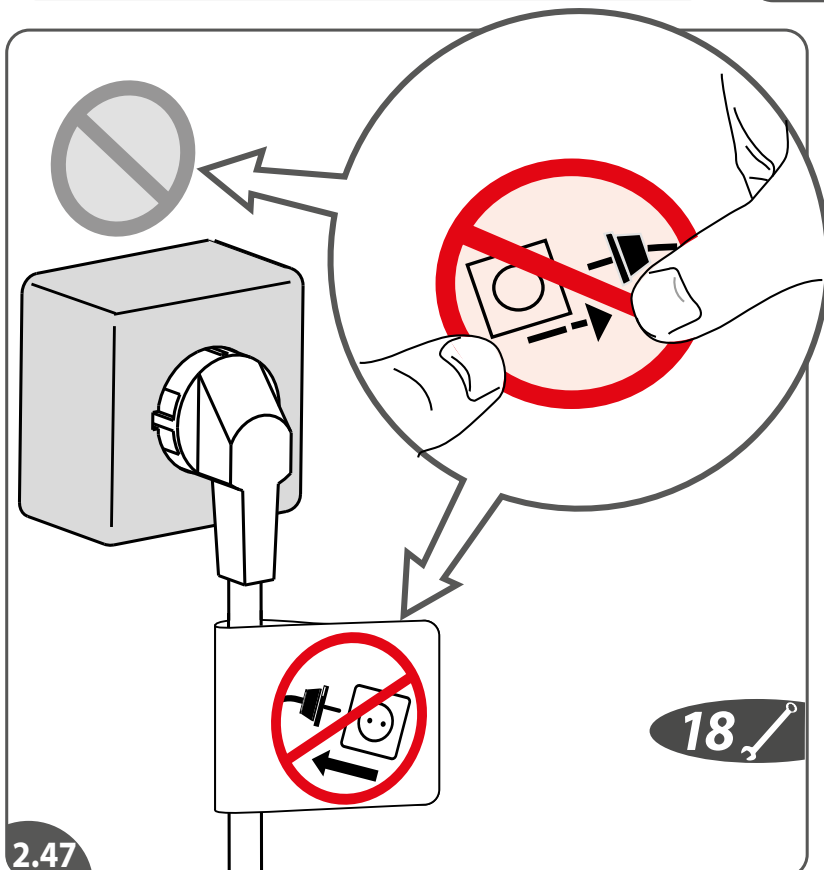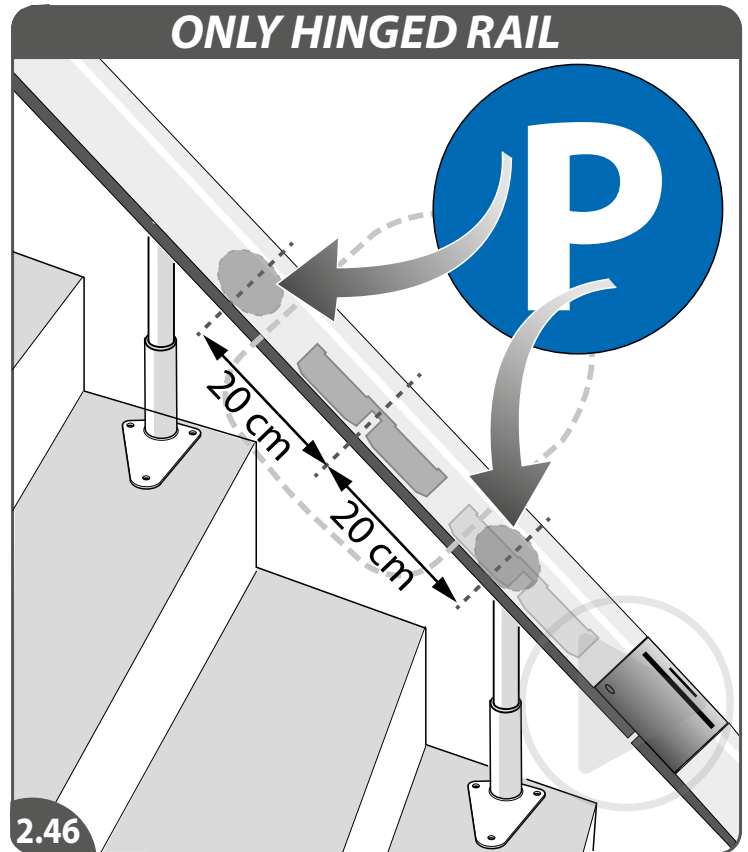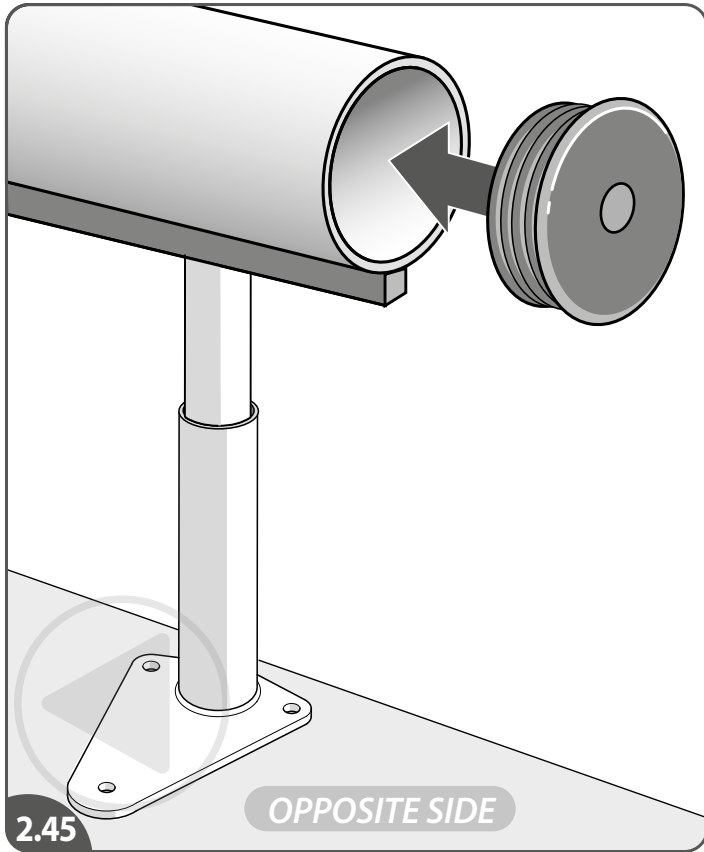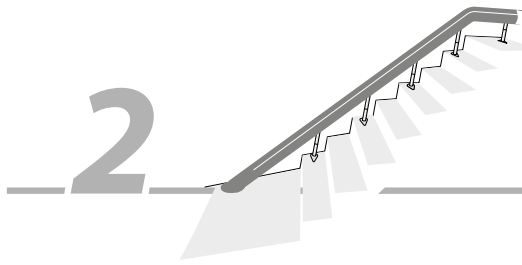[2.45] Mount the endcover on the opposite side of the rail.
[2.46/2.47] Place the warning labels.

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



3

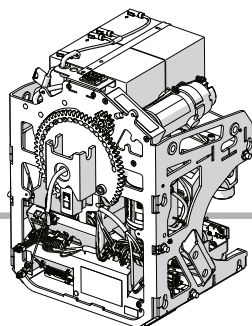POWER PACK LOCKING PLATE

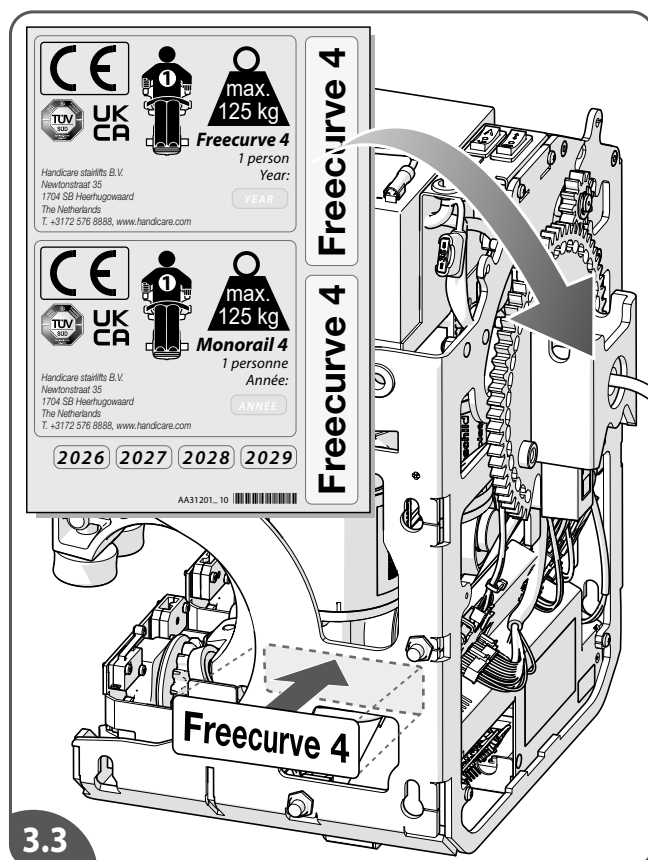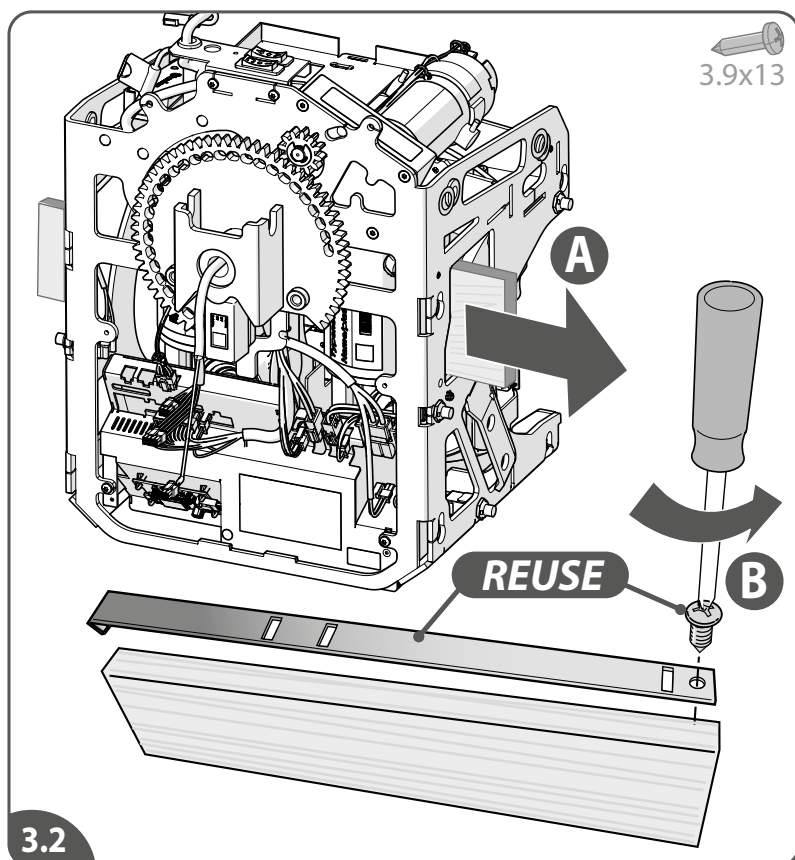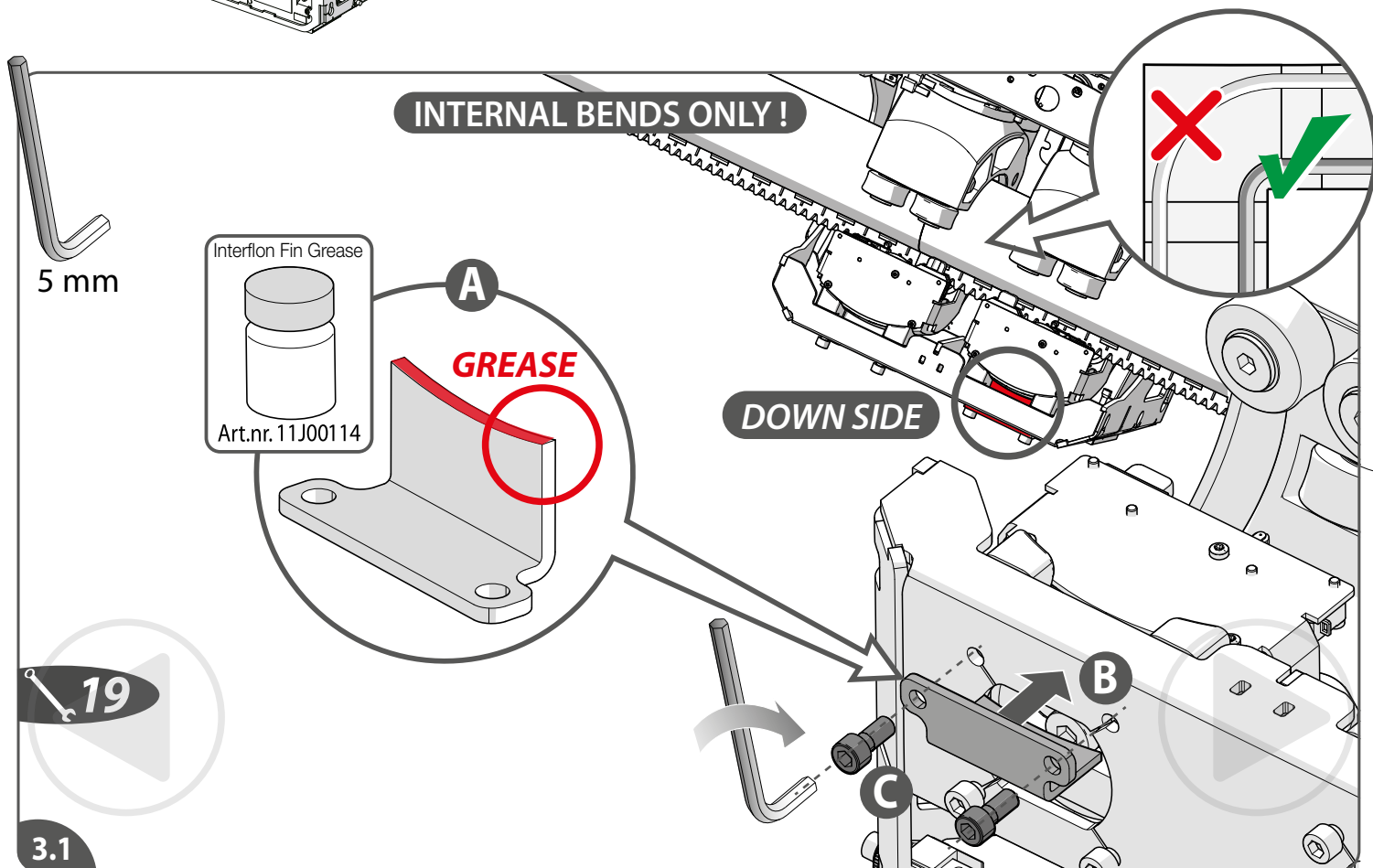- [3.1] In case of only internal curves, install the bracket (on the motor, at the down side). Grease the curved surface of the motorframe.
- [3.2] Remove the transport bar{A}. Remove the strip from the transport bar{B}. Please keep this screw to tighten the battery strip later in the installation.
- [3.3] Put a **Freecurve 4**-label on the side of the powerpack.

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




3

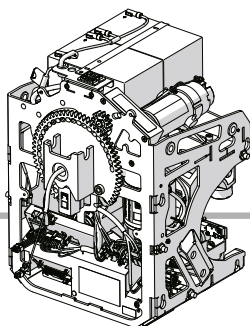POWER PACK FRONTPLATE

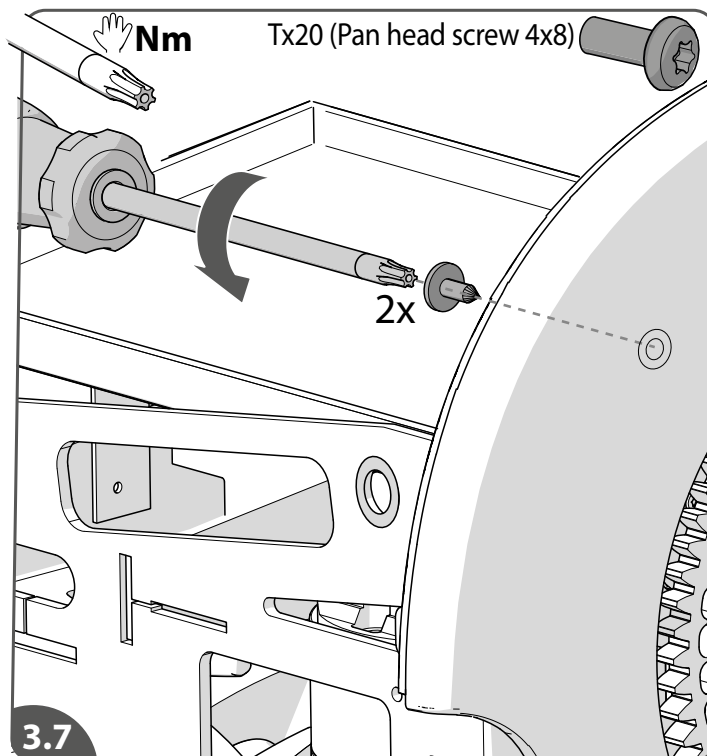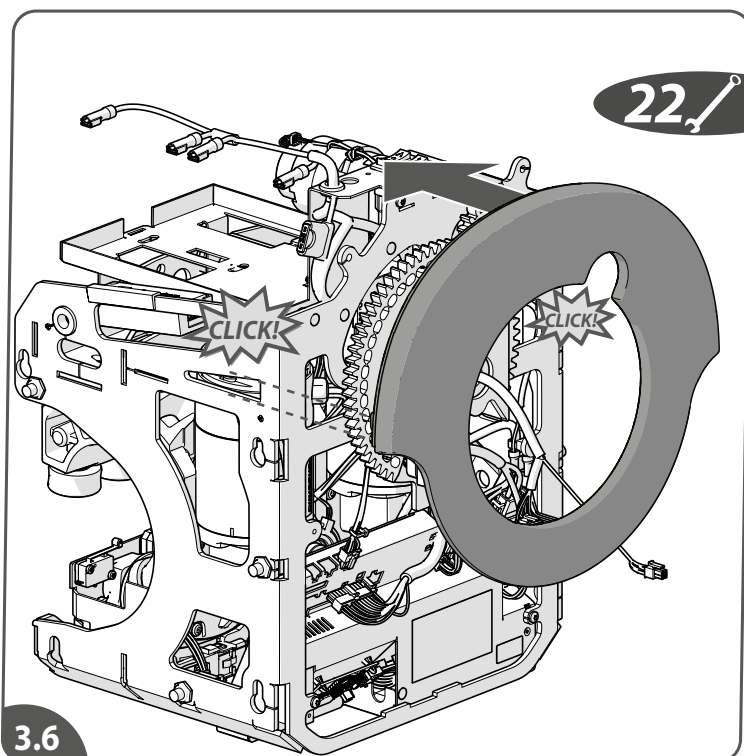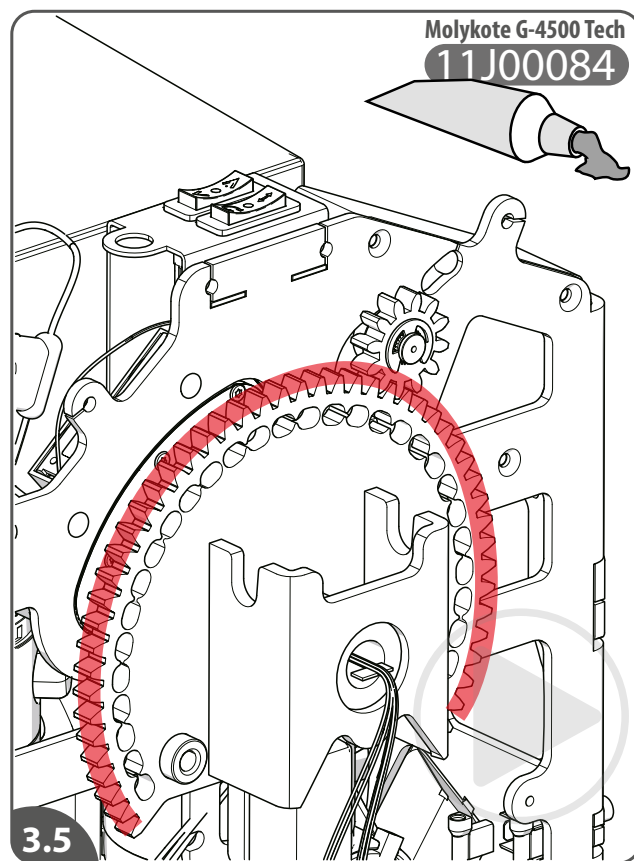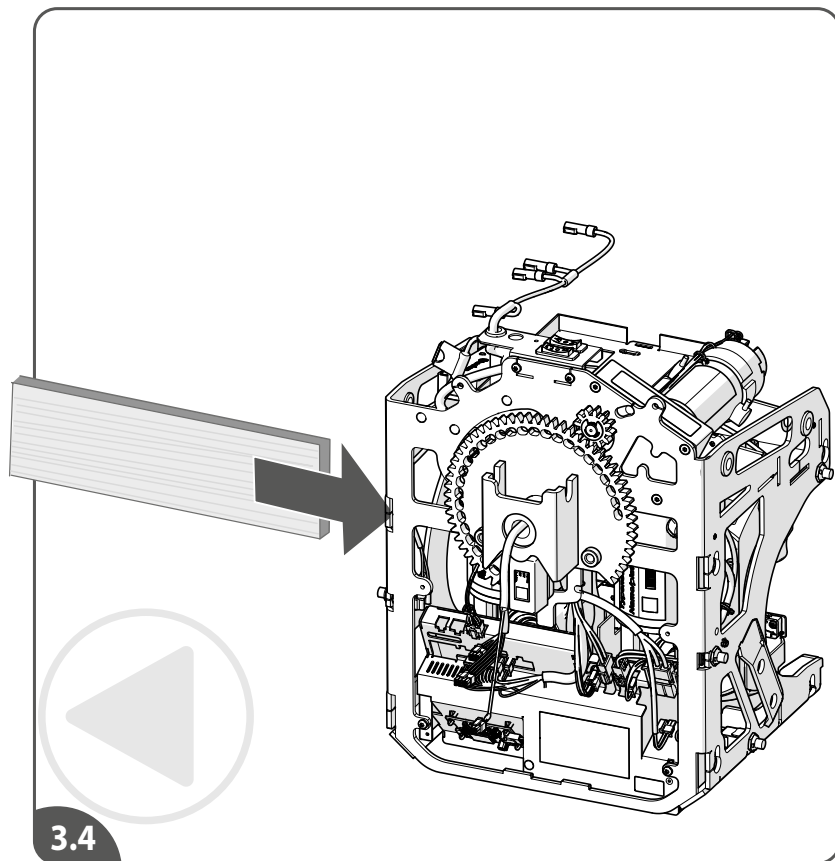- [3.4] Replace the transport bar in the powerpack.
- [3.5] Grease (11J00084 {*Molykote G-4500-Tech*}) the gearsegment of the powerpack.
- [3.6] Install the frontcover.
- [3.7] Fasten the frontcover (handtight).

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
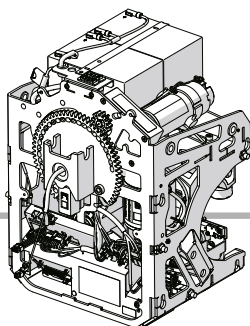
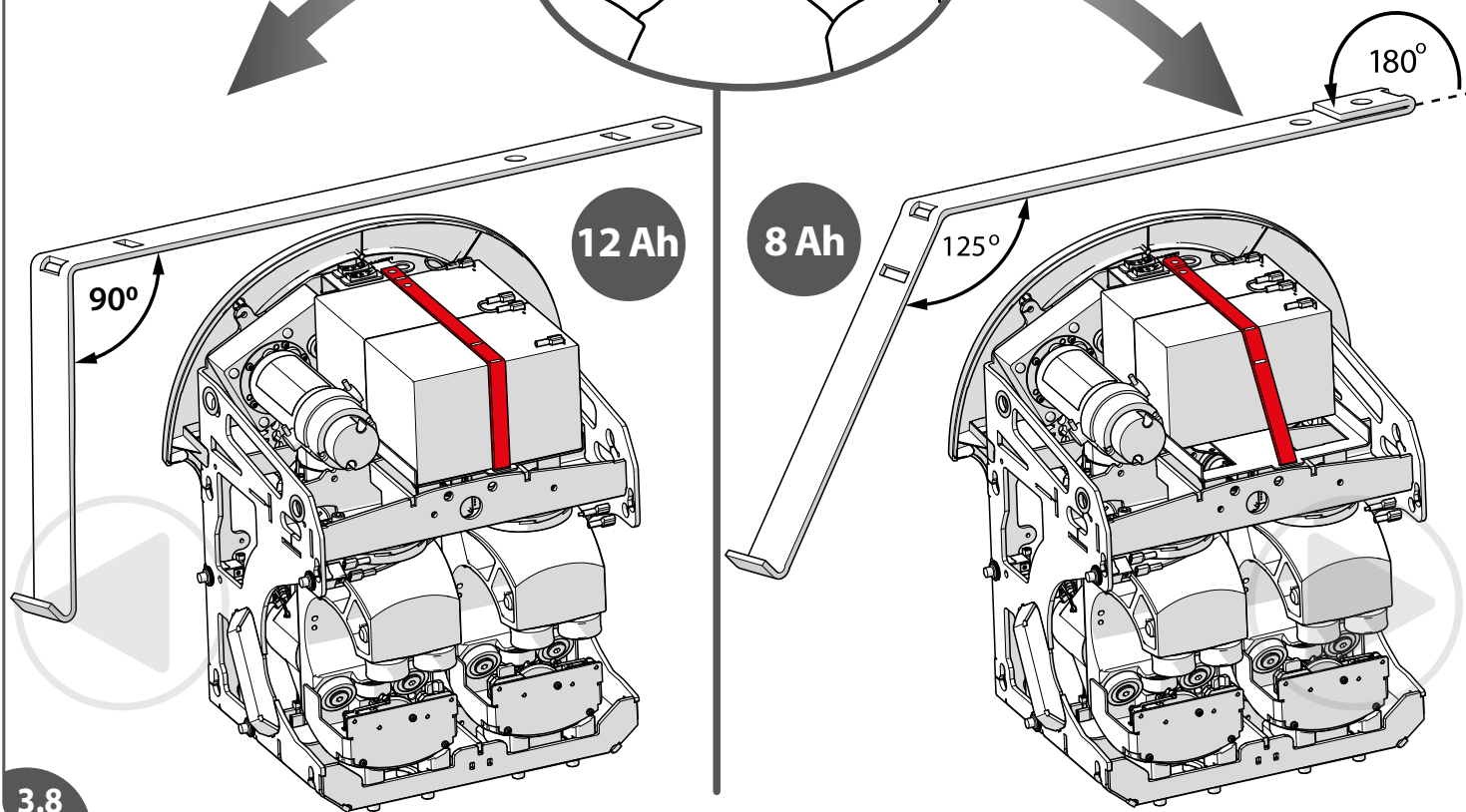
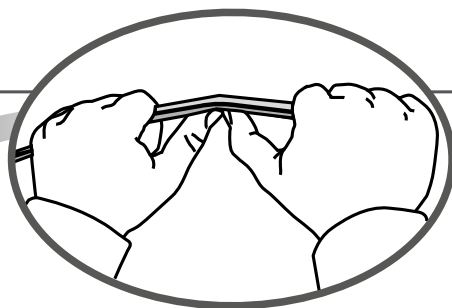
3.8

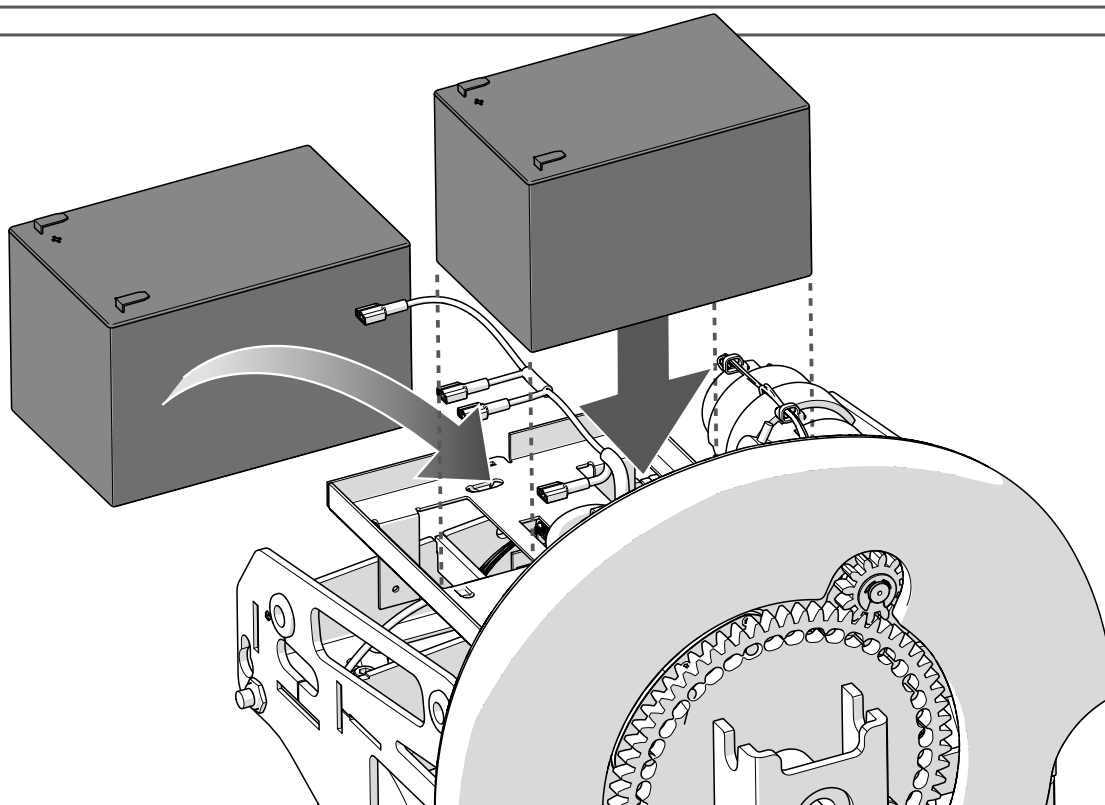3.9

[3.8] Depending on the capacity of the batteries, the strip should be bent on different positions.
[3.9] Install the batteries in the powerpack.

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




3

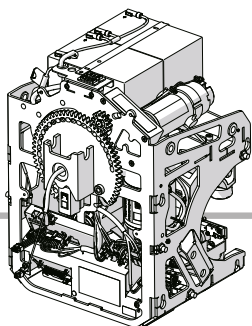POWER PACK SHORTENING LINKS

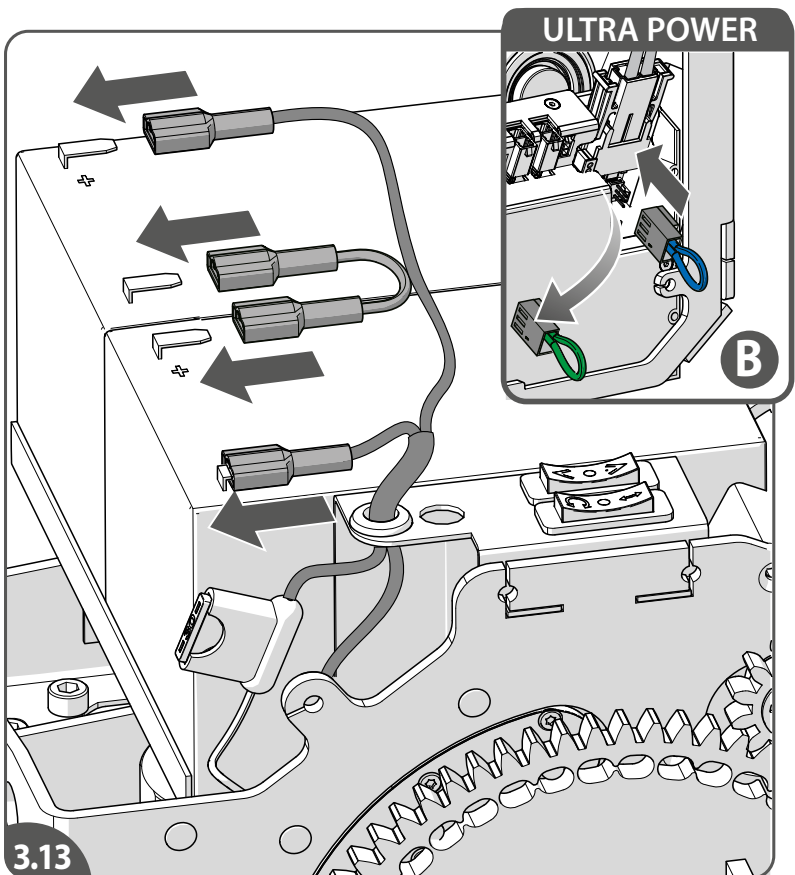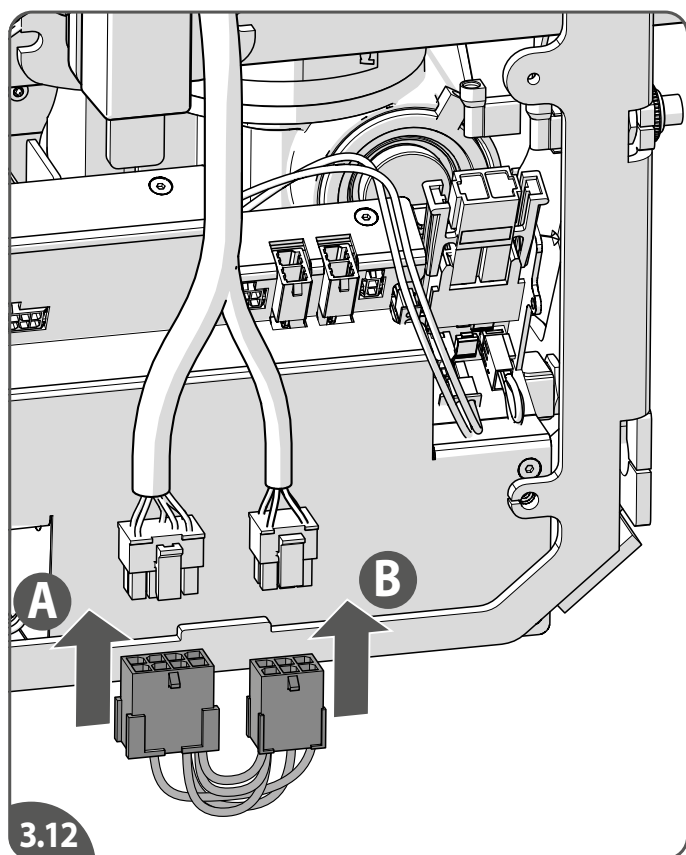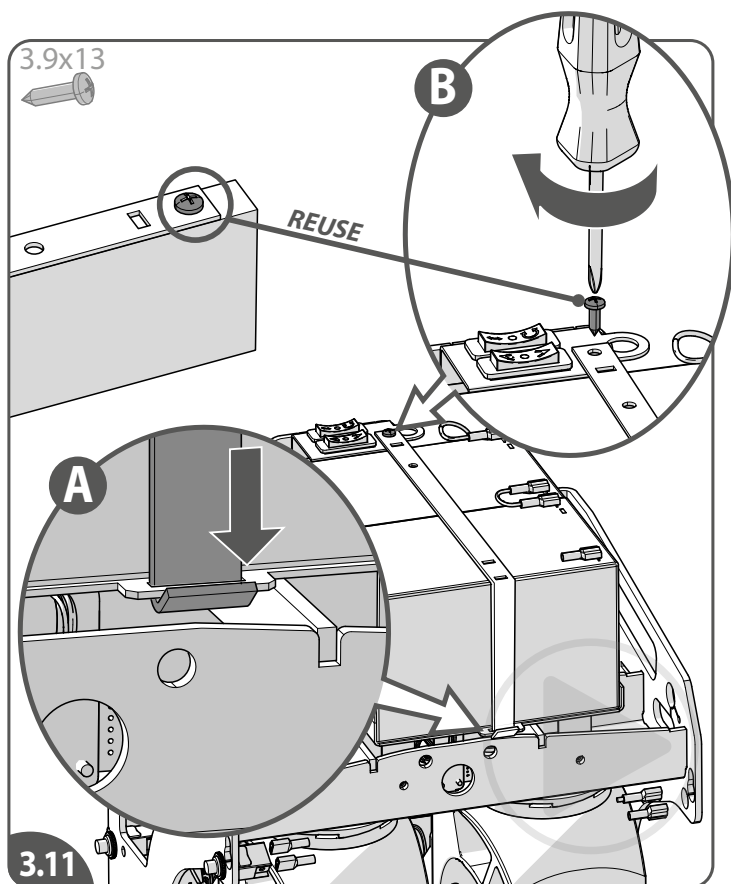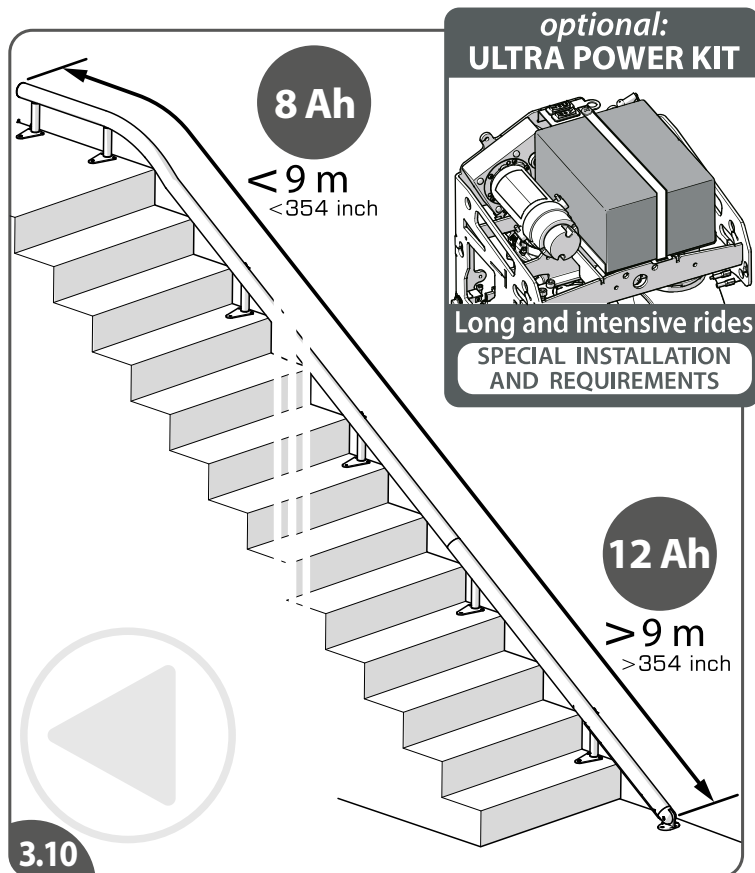[3.10] The required capacity of the batteries depends on the length of the rail and the intensity of the usage of the lift.

[3.11] Hook the strip in the cutout of the battery frame{A}. Fasten the strip to the battery frame{B}.

[3.12] Connect the shunts.

[3.13] Connect the batterie wires. In case of a **ultra power kit** remove the *green* jumper connector and install the *blue* jumper connector{B}.

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




IT

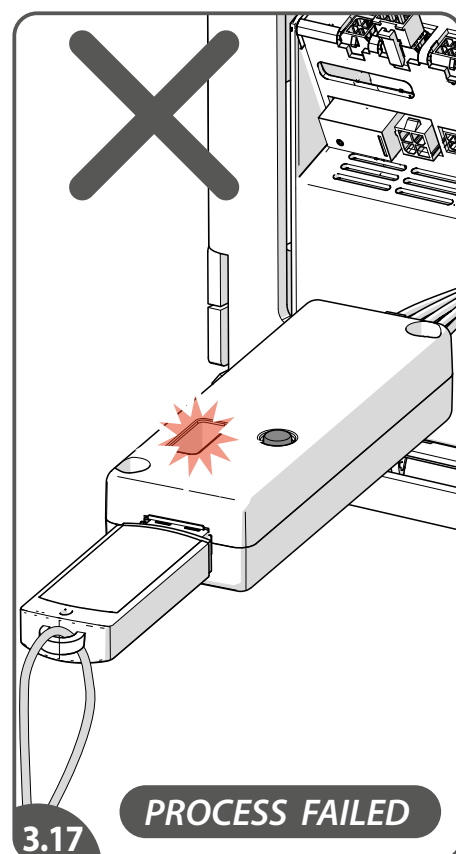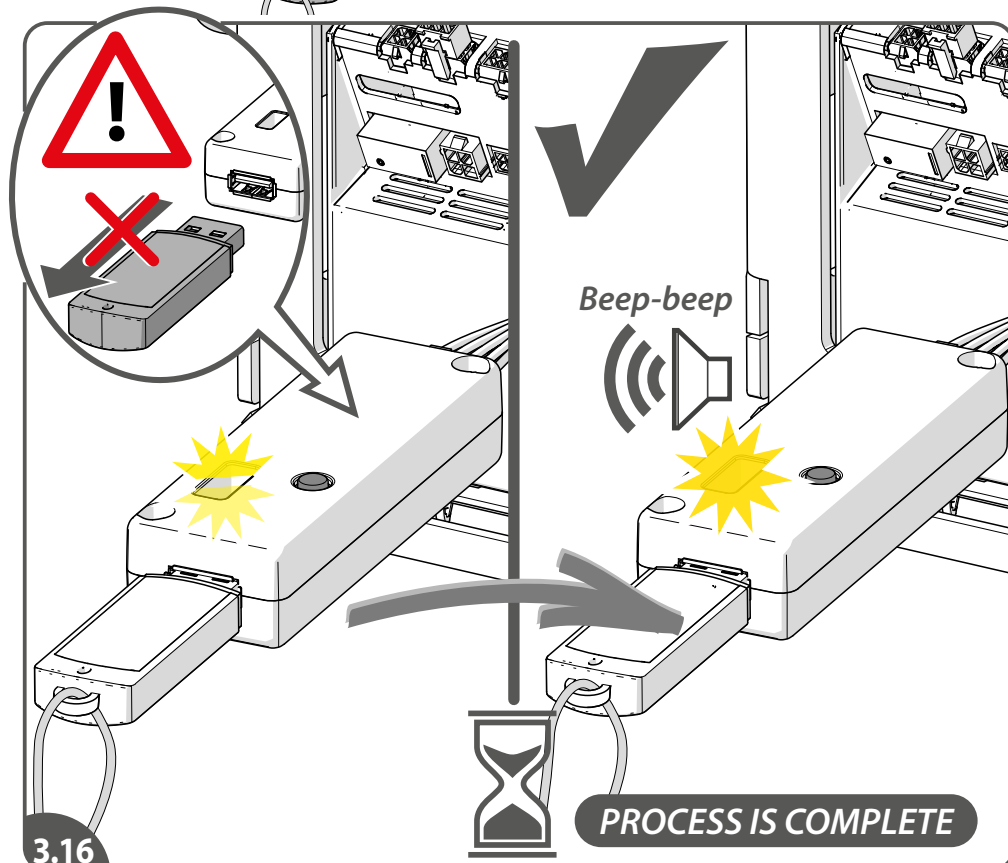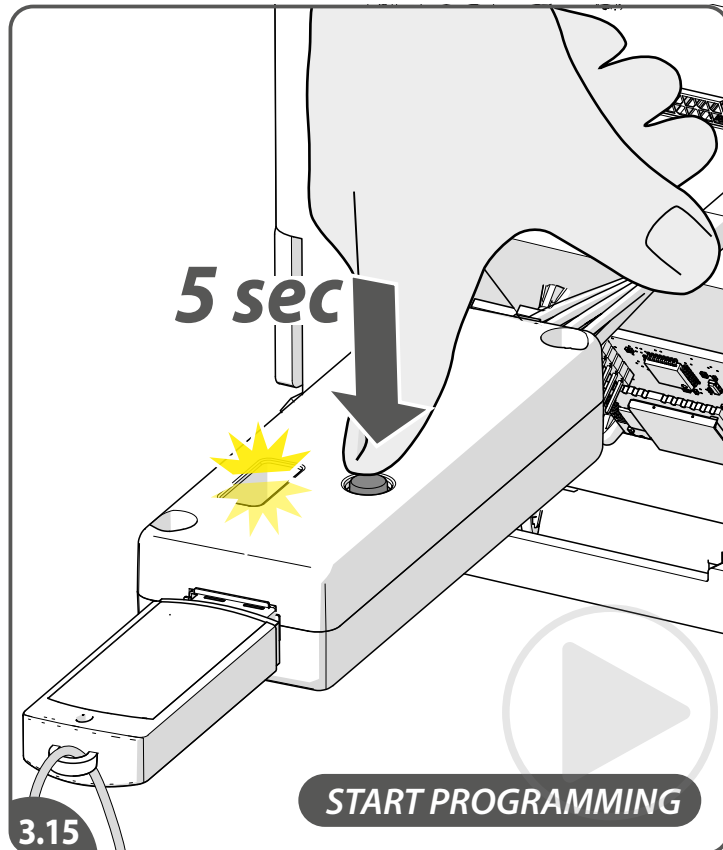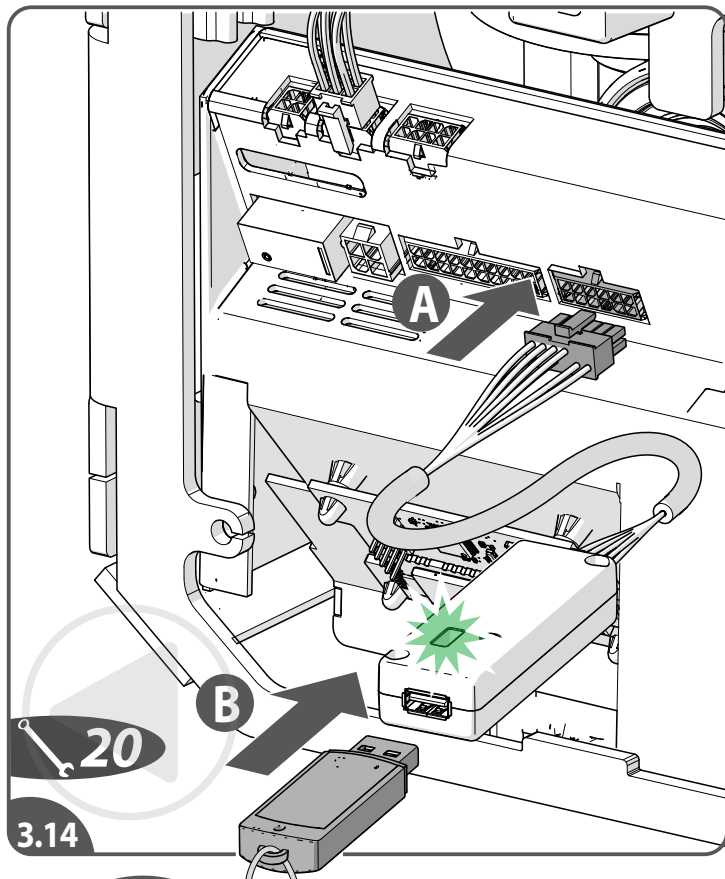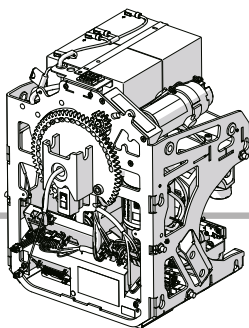[3.14] Connect the connector of the rail data tool to the PCB{A}. Insert the USB stick into the rail data tool{B}. When the LED is solid green it is okay to proceed.
 [3.15] Press and hold the button for approx. 5 sec until yellow LED starts blinking (copy data from USB to PCB has started). Release the button.
 [3.16] Wait until the yellow LED becomes solid and the pcb gives 2 beeps: the process is complete. Do not remove the tool or USB stick during this process.
 [3.17] The red LED will activate when the process has failed. On the next page are the possible statuses of the LED of the data tool explained.

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

IT

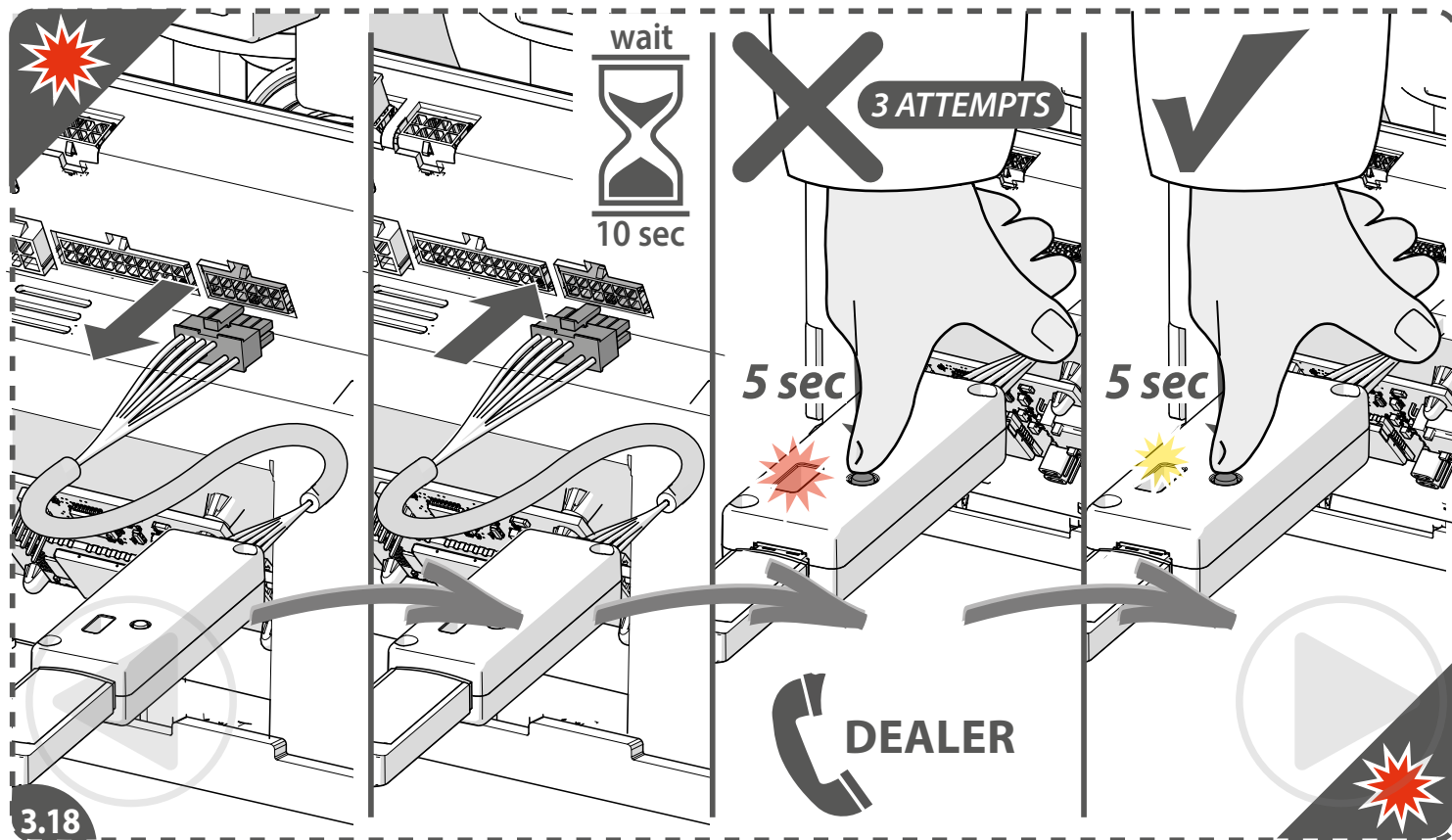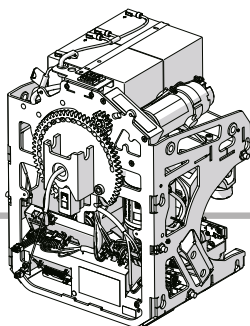GREEN	Off = No Power Flashing = indicate Power on On Solid = USB stick detected
YELLOW	Off = No Communication Flashing = indicate communication from USB to Processor On Solid = transfer of data to Processor completed ok
BLUE	Off = No Communication Flashing = indicate communication from Processor to USB On Solid = transfer of data from Processor to USB (incl. error log) completed ok
RED	Off = All good Flashing = no trackdata file found on usb stick or file corrupt On Solid = Error Flashing = indicate track data file on usb stick is corrupt or usb stick is of not-compatible format

[3.18]If the process fails, disconnect the rail data tool. Wait 10 seconds and reconnect the datatool/USB stick. Press the button for another 5 seconds. If the process still fails after 3 attempts, call the dealer.

[3.19]Explanation of the LED colours on the USB stick.

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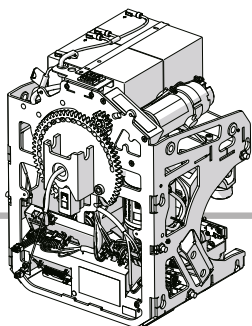
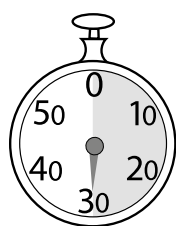
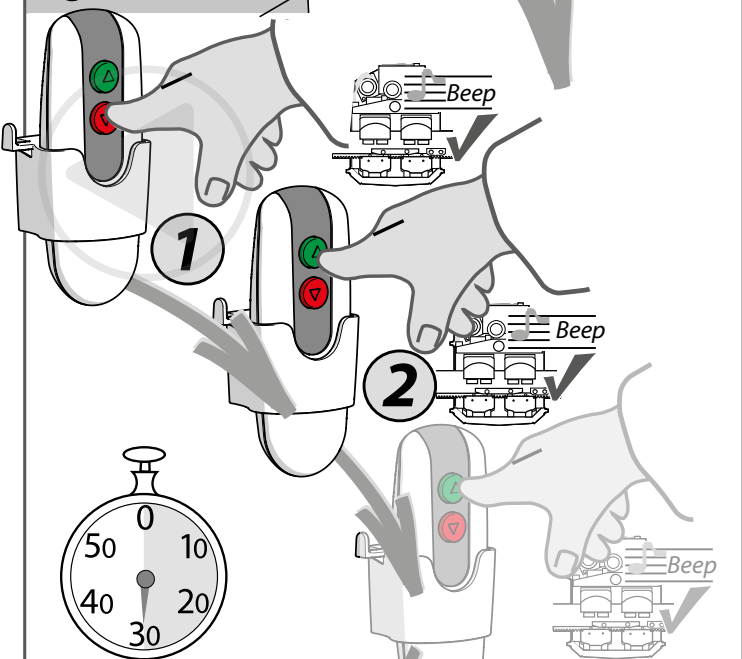
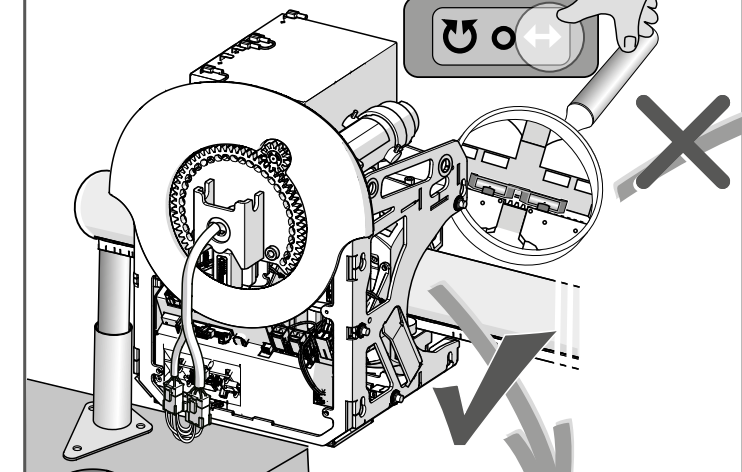
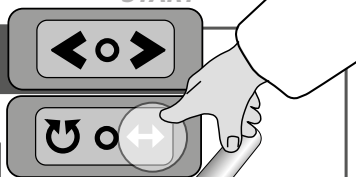
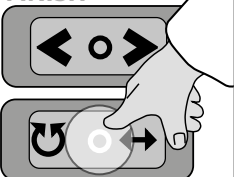
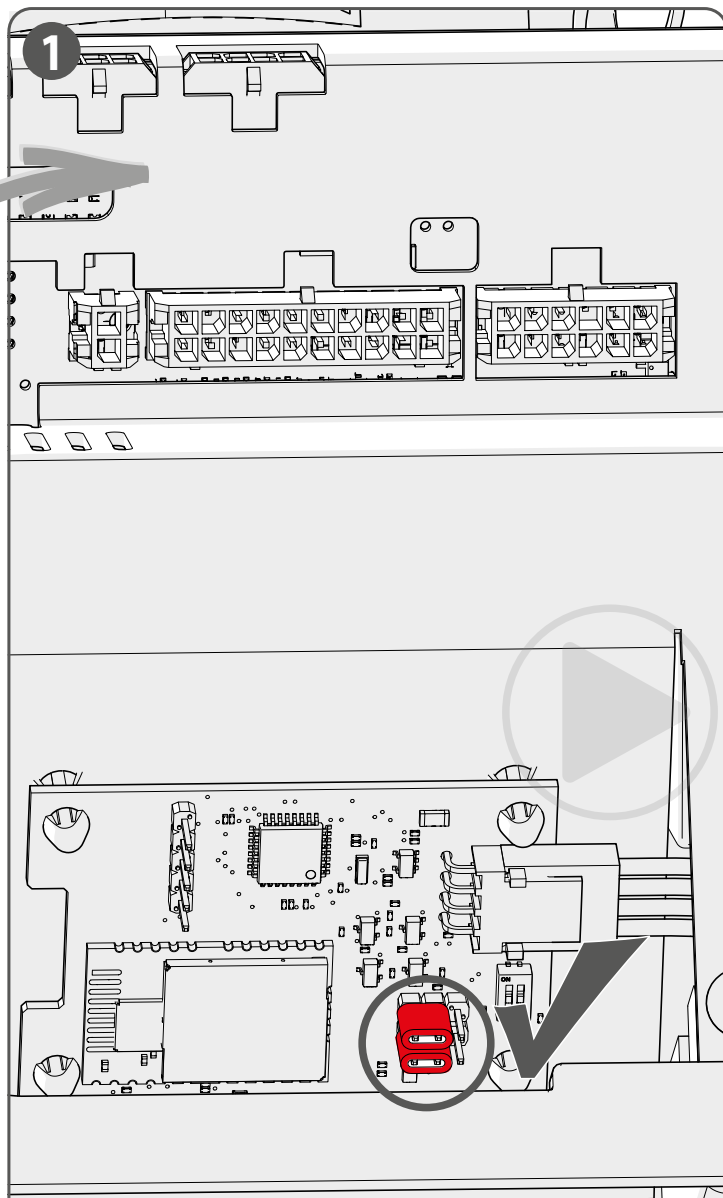
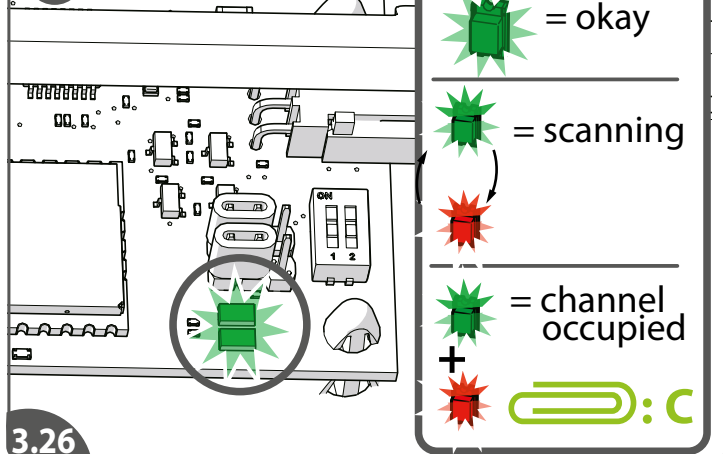
-  = okay
-  = scanning
-  = channel occupied
-  +  : C

3.26

[3.25]Link the remote controls: press the calibration switch to traction mode. Link the remotes one by one by pressing a button on every remote (beep: linking is okay). If all the remotes are linked confirm the procedure by putting the calibration switch in neutral position. Linking of the remotes have to take place within 30 seconds. If not, repeat the procedure.

[3.26]If no pairing of the remote controls is possible check the jumper setting{1}.

If the remote controls are still not connecting (LEDs on the board), see attachment **C** to solve the problem: the dipswitch settings may not be correct {2}.

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

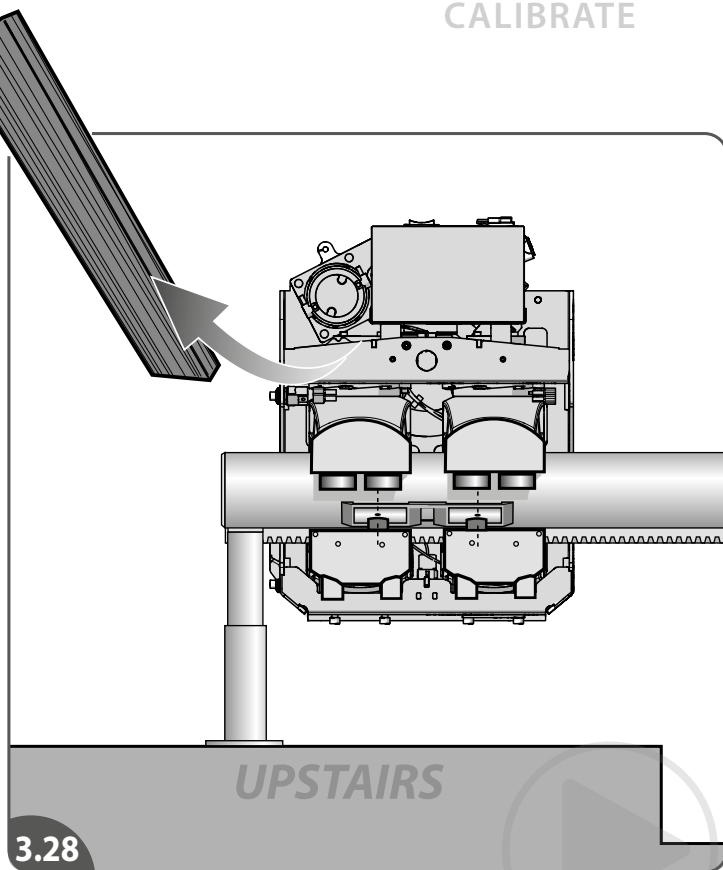
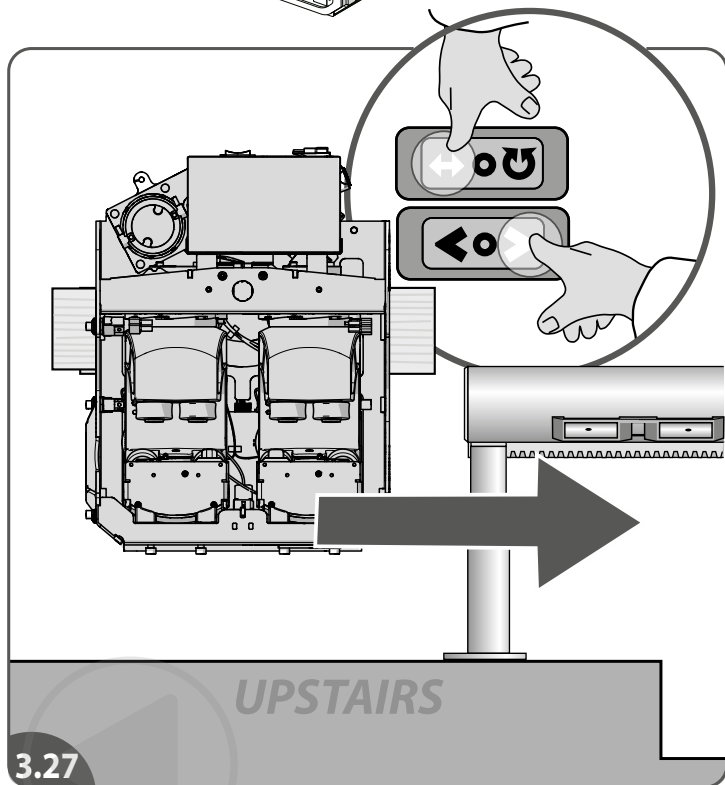
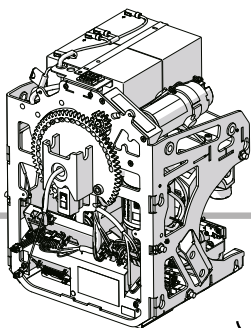
CALIBRATE TRANSLATION

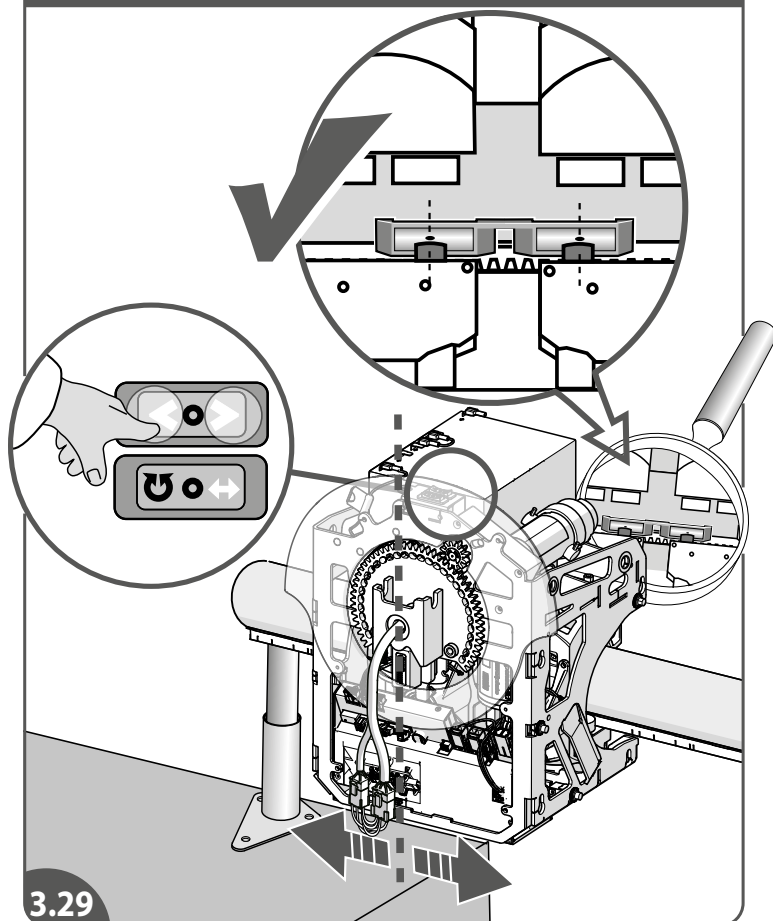CALIBRATE ROTATION

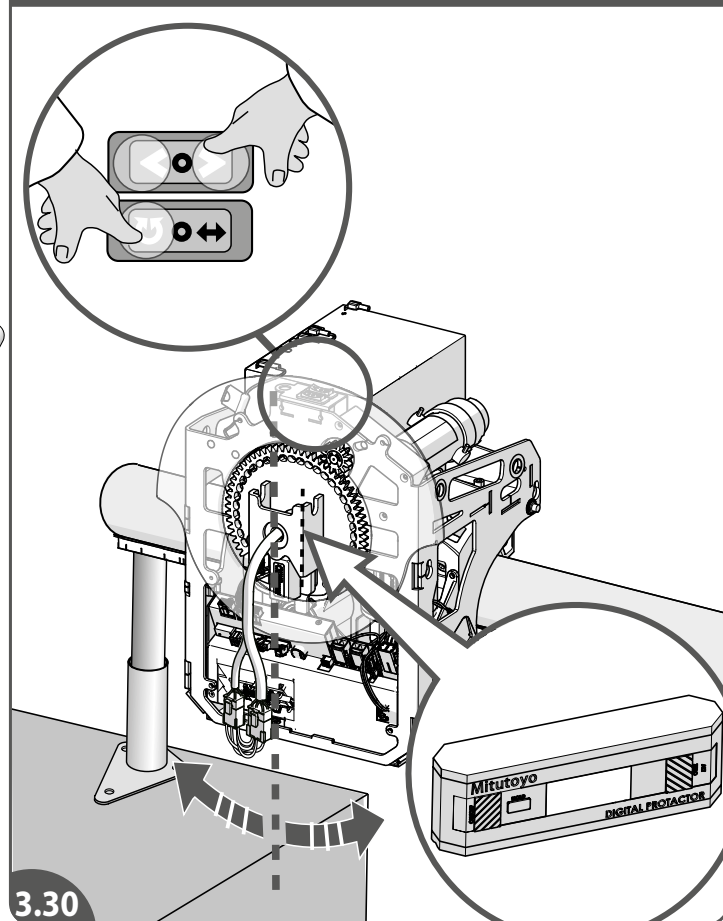- [3.27] Drive the powerpack onto the track.
- [3.28] Remove transport bar.
- [3.29] Ride the powerpack to the calibration point, this is the center of the upper charge contact
- [3.30] Calibrate the rotation of the power pack.

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



FINISH CALIBRATION

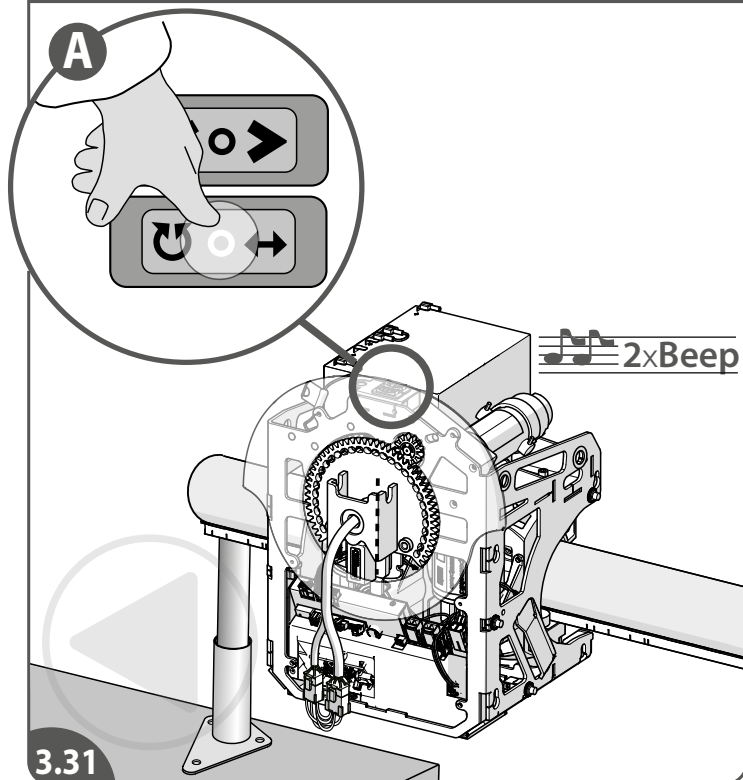TESTRIDE

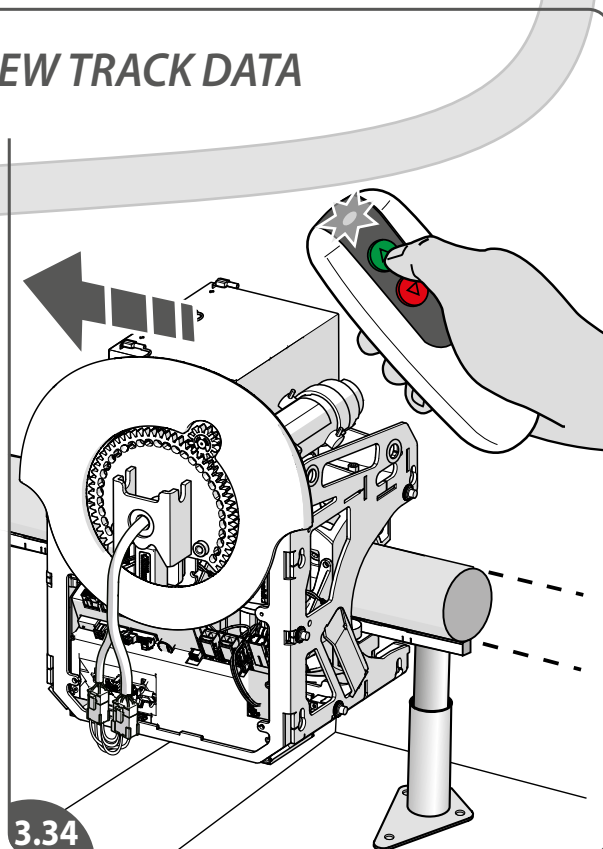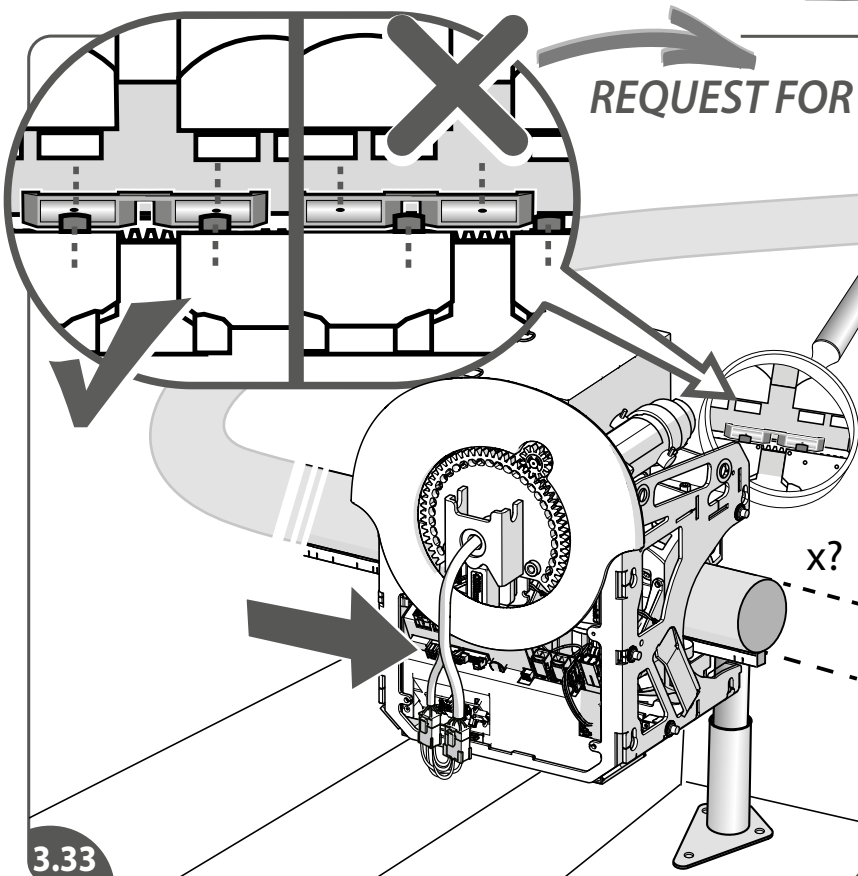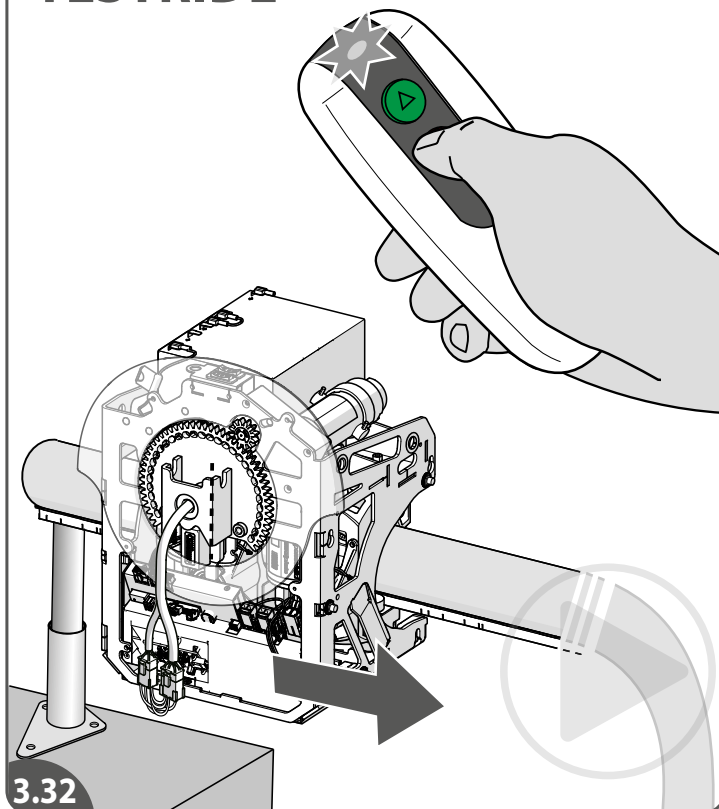[3.31]Close the calibration procedure by putting the calibration switch on the unit to the neutral position{A}. By closing the calibration the lift beeps twice{B}.
 [3.32]Make a testride.
 [3.33/3.34]Check the stop on the chargestoppoints. Request Savaria for new track data in case the locations of the chargestoppoints needs to be changed.

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




IT

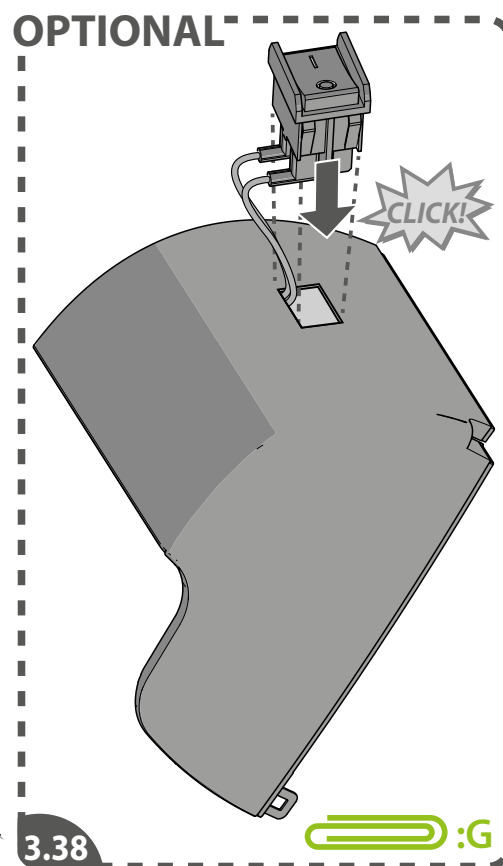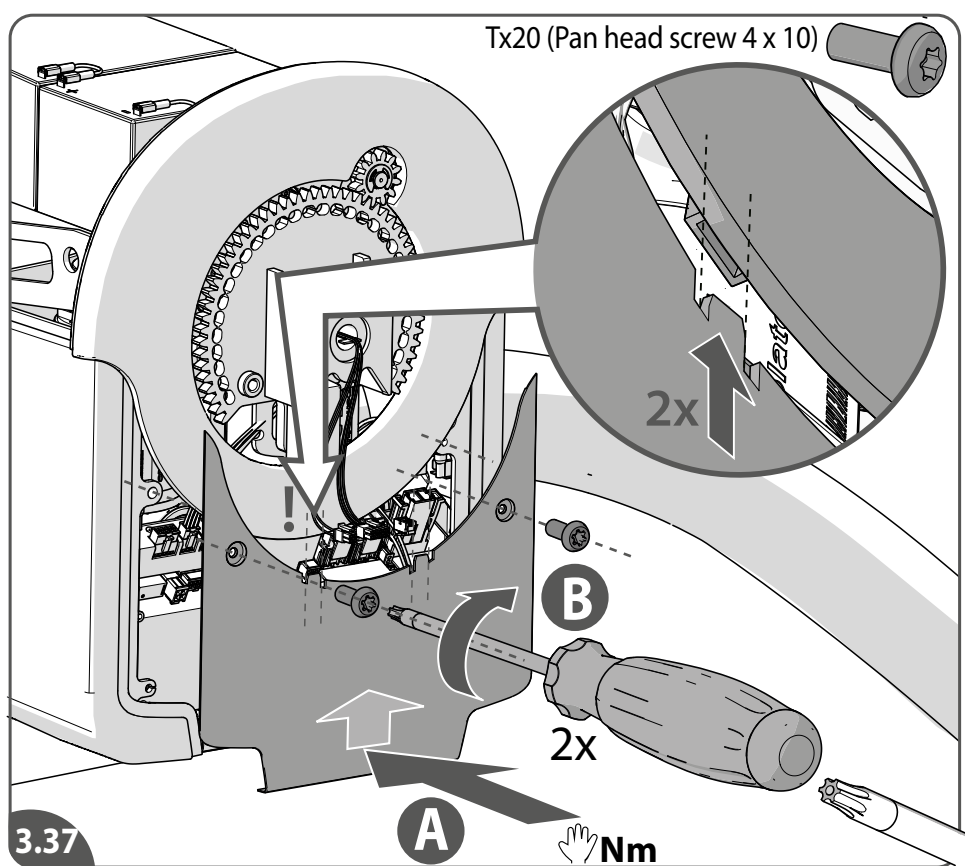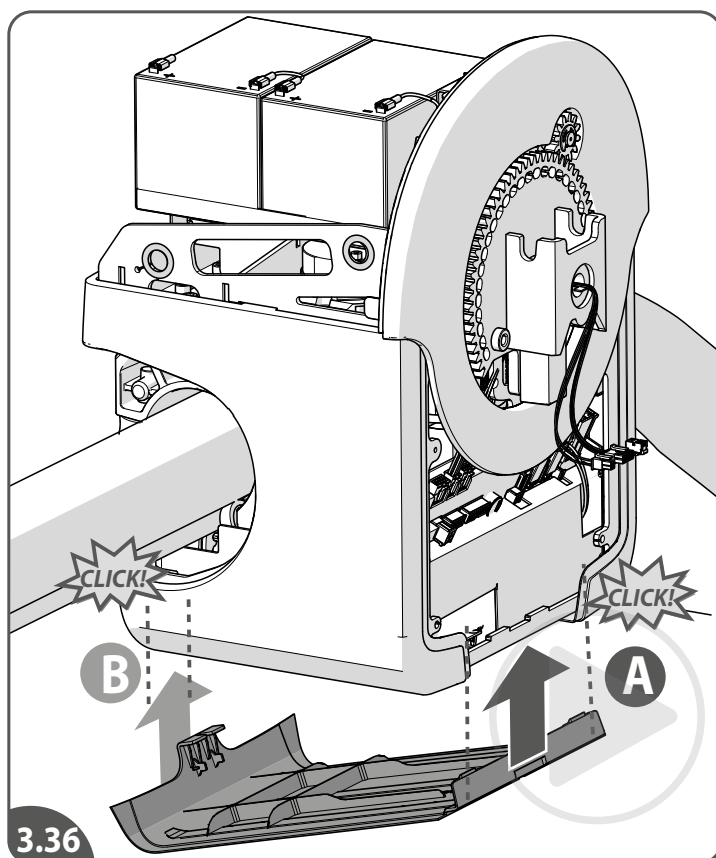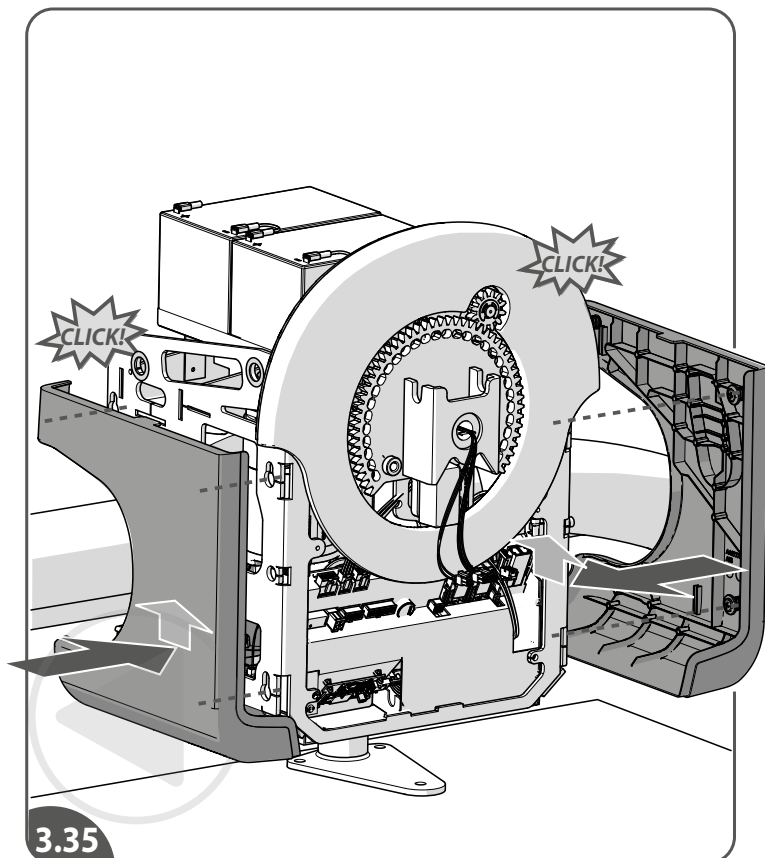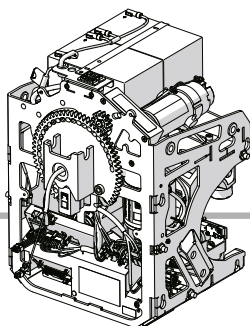[3.35] Install the side covers. Check the function of the side covers.

[3.36] Install the bottom plate.

[3.37] Install the other frontplate (handtight). In doing so, slide the tabs of the bottom front plate into the slots of the top cover.

[3.38] *Optional*: see attachment **G** (page 120) how to install a lockable switch on the top cover.

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



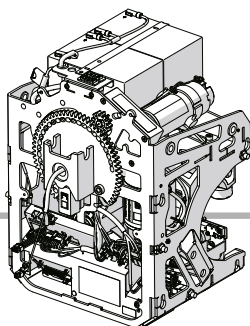
POWER PACK

FASTEN AND CHECK COVERS

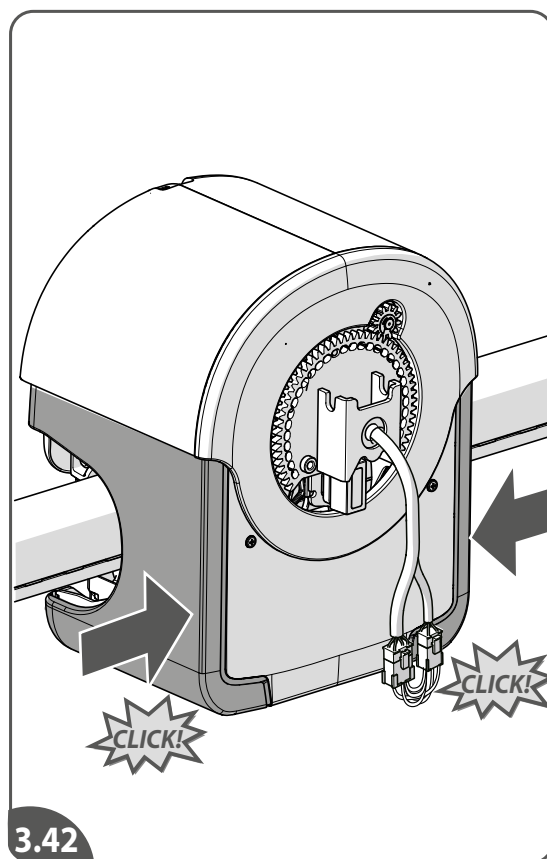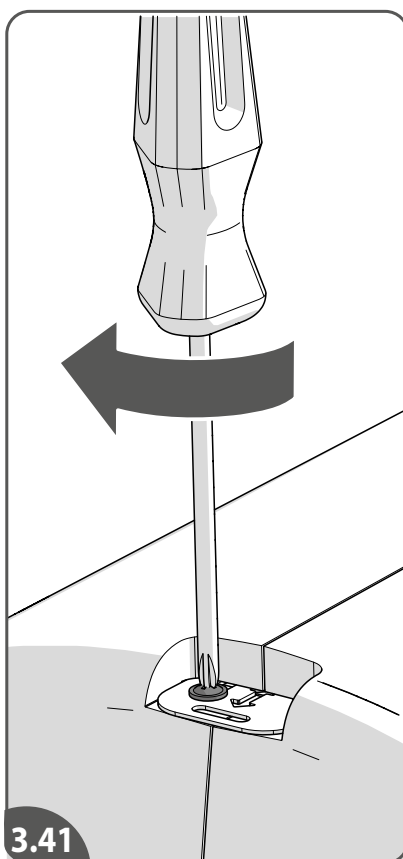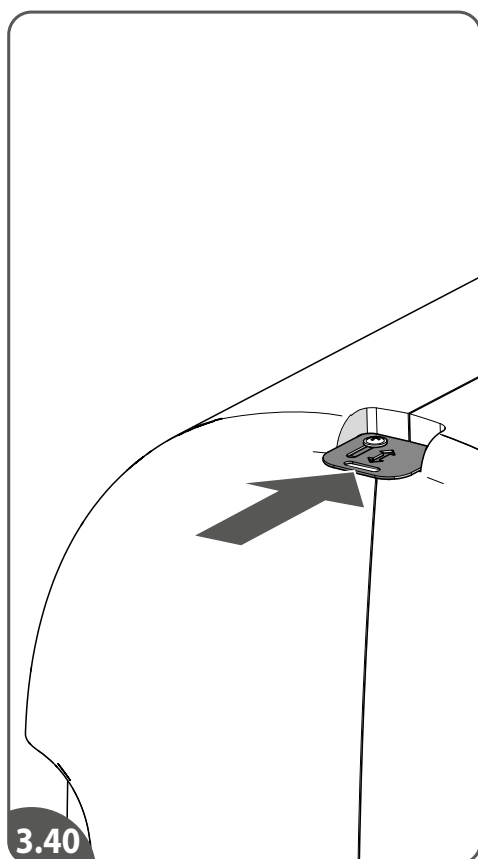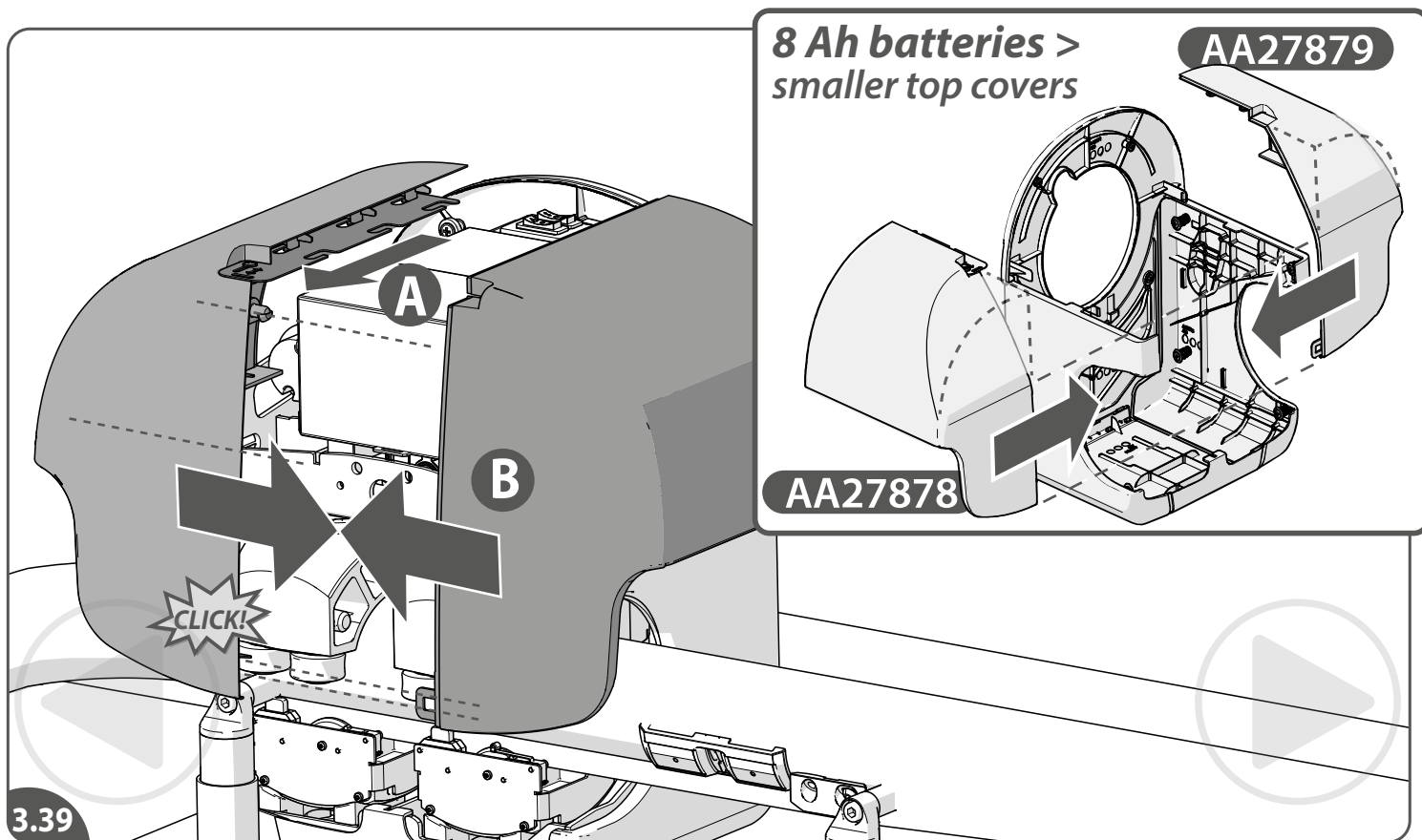- [3.39] Slide the lock{A}. Install the top covers{B}.
 [3.40] Lock the top covers. In case of applying the low capacity batteries use the smaller top covers.
 [3.41] Fasten the lock.
 [3.42] Check the function of the side covers.

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





IT

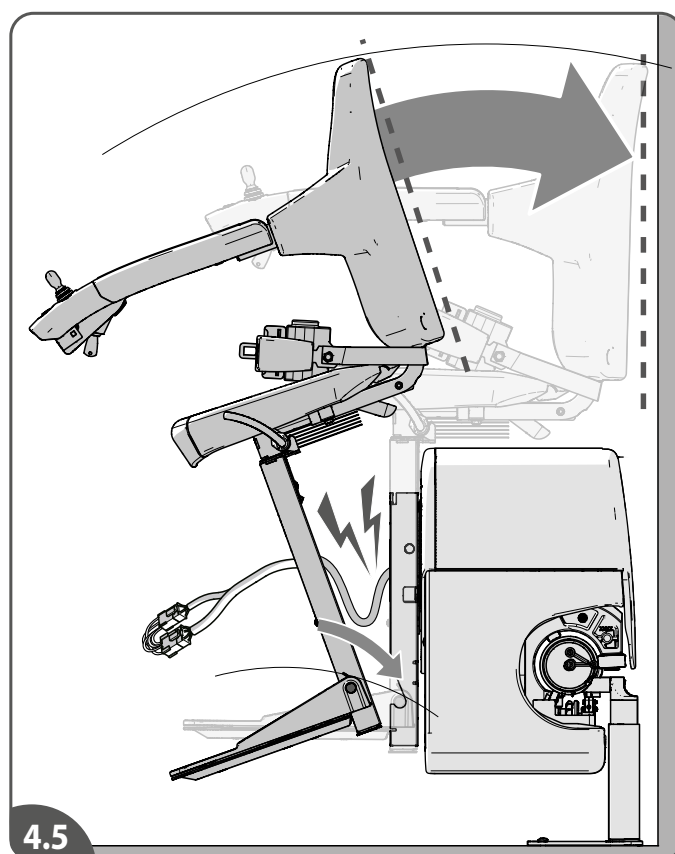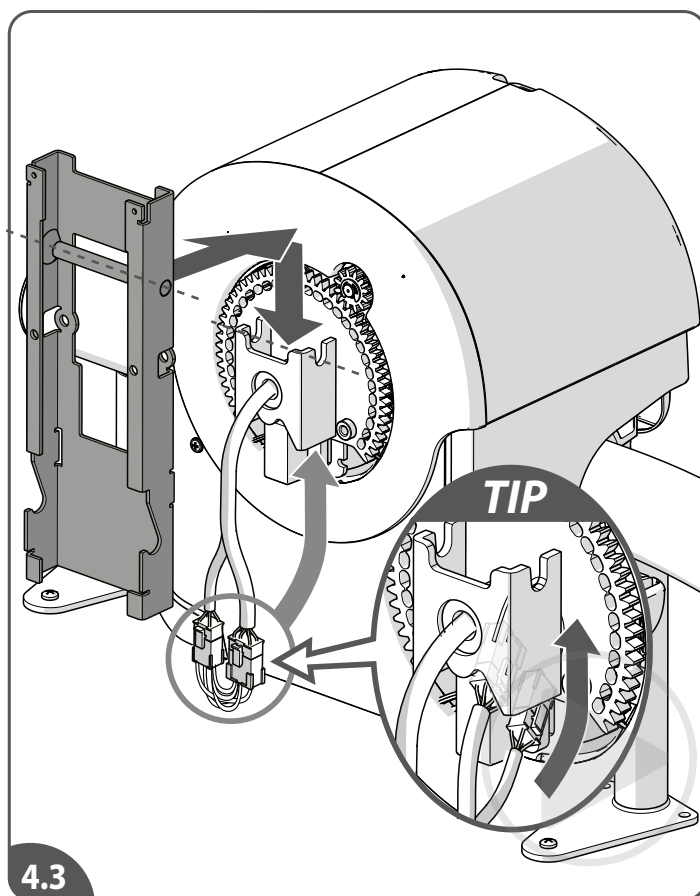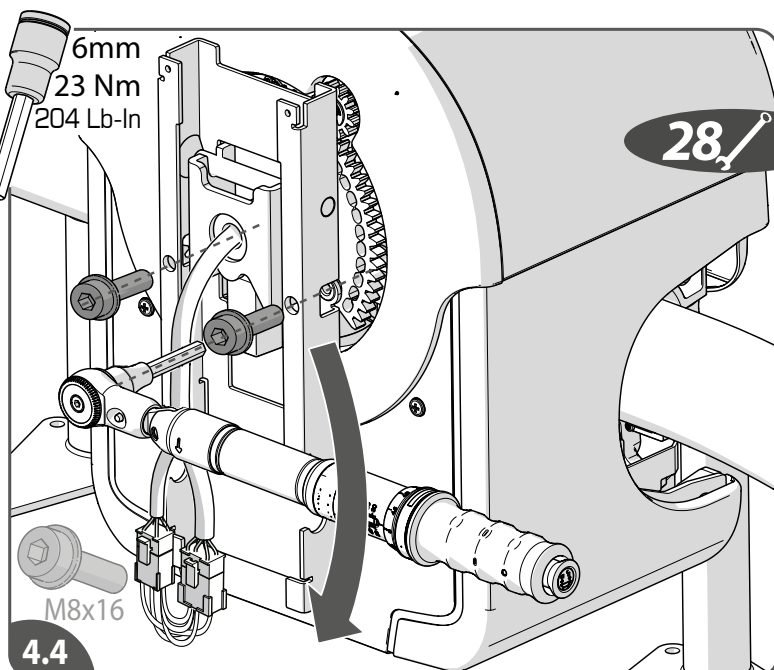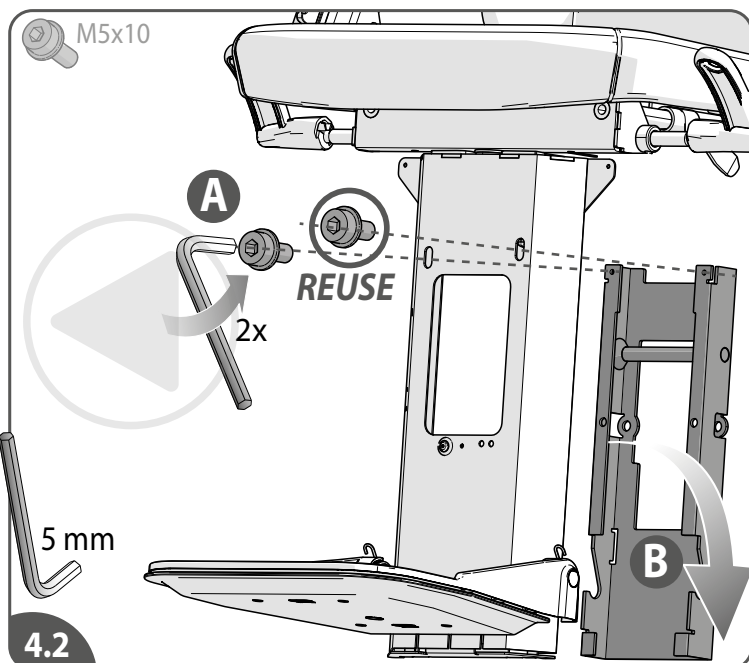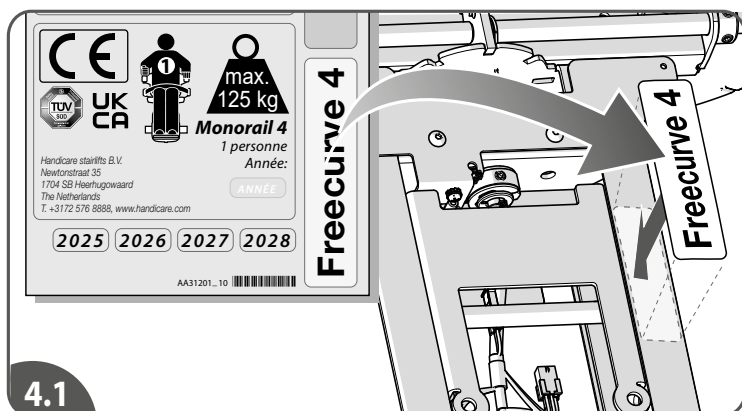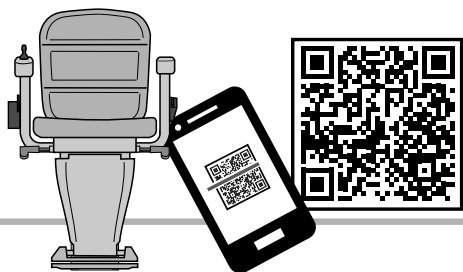- [4.1] Put a **Freecurve 4**-label on the backside of the seat frame.
 [4.2] Install the chair at the bottom of the stairs. Remove the frame from the seat. Reuse the bolts to fasten the seat to the frame later on.
 [4.3] Install the frame on the powerpack.
TIP: tuck the seat cable away before installing the frame/seat.
 [4.4] Fasten the frame [23 Nm | 204 Lb-In].
 [4.5] Install the seat on the frame. Don not pinch the cables!

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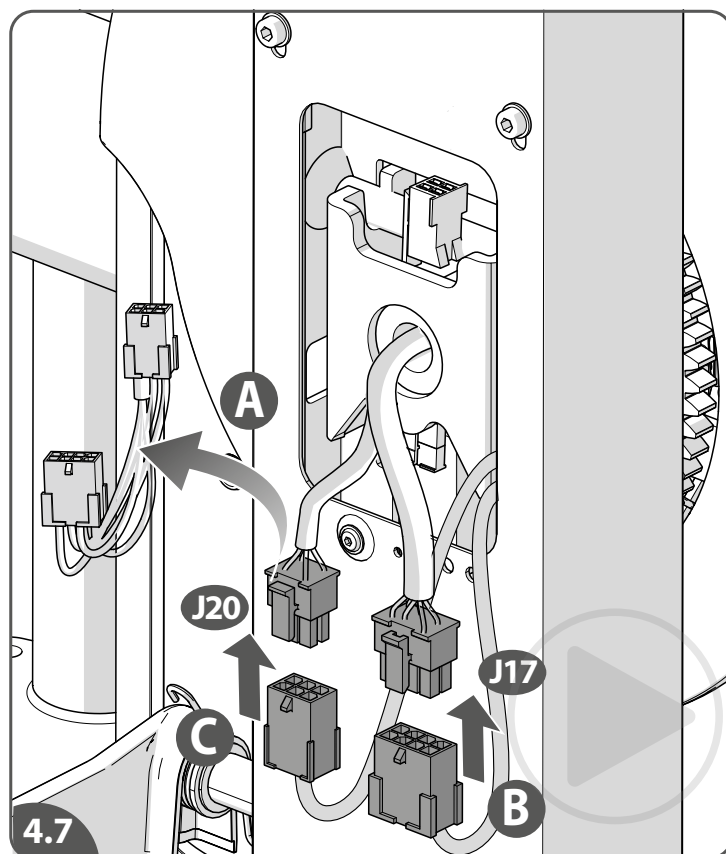
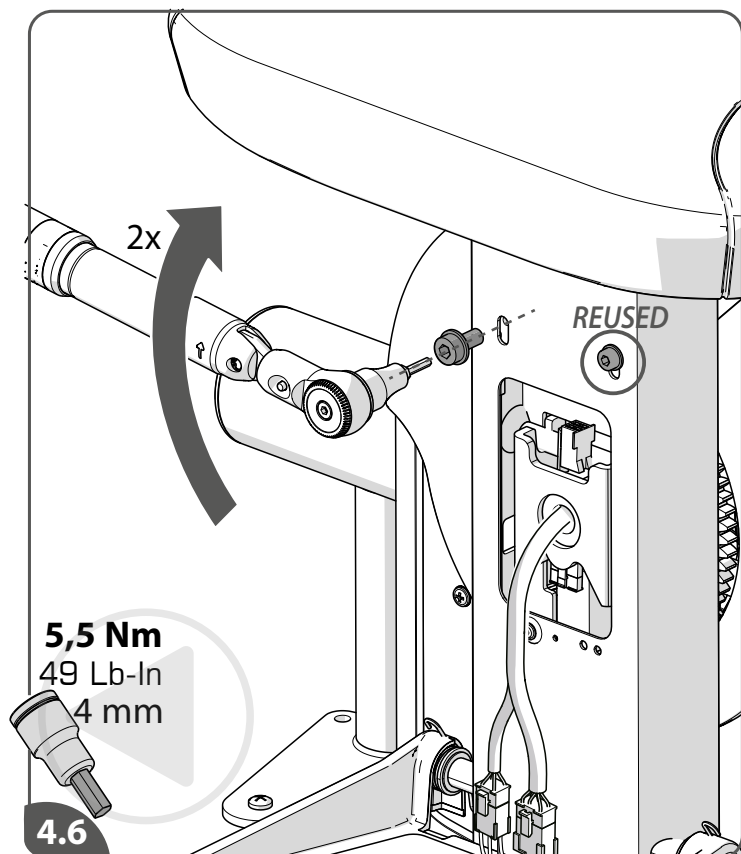
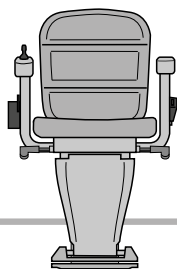
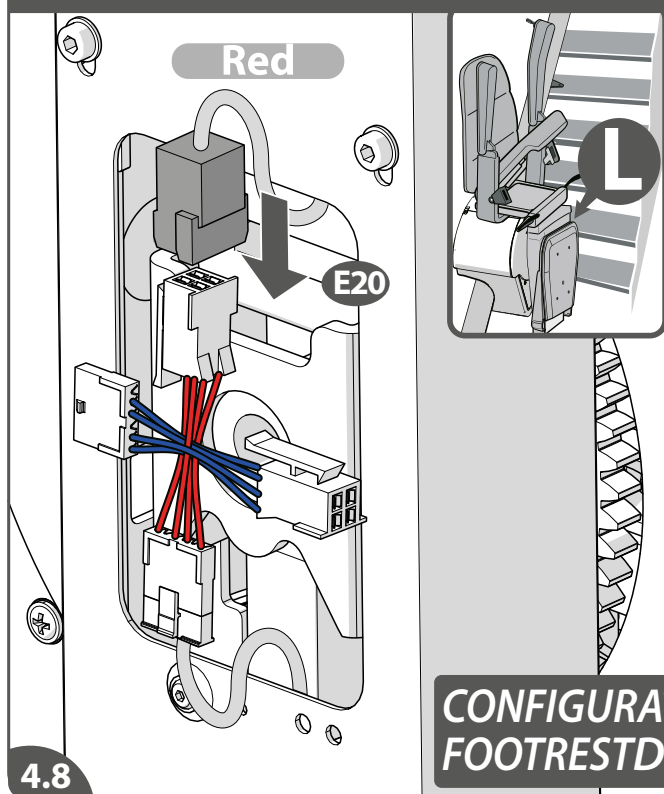
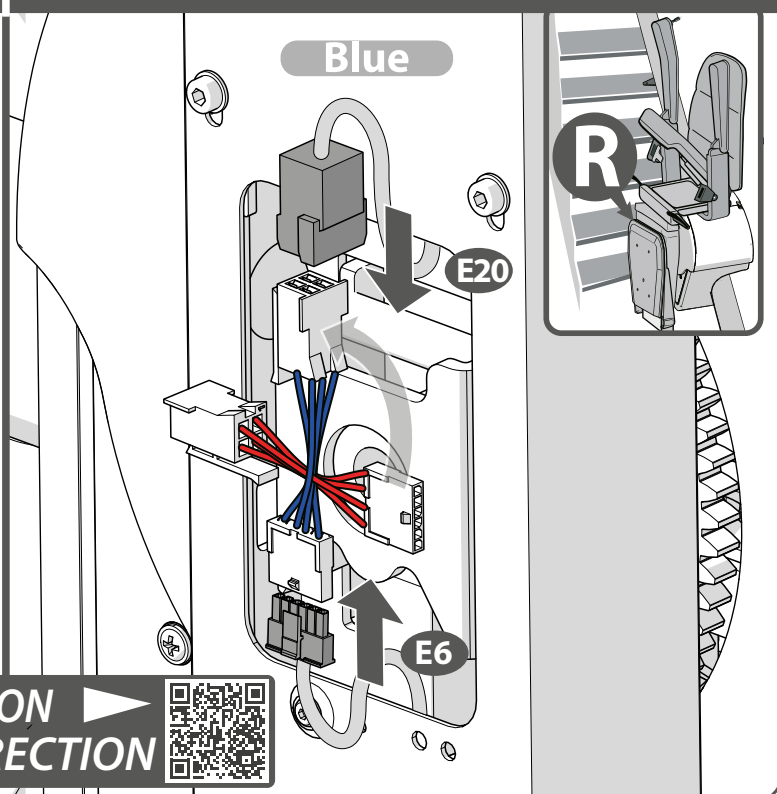
CONFIGURATION
FOOTRESTDIRECTION

[4.6] Fasten the seat on the frame [5,5 Nm | 49 Lb-In].

[4.7] Remove the shunts {A}. Connect the seat loom to the powerpack **[J17]** {B}. Connect safety loom to the powerpack **[J20]** {C}.

[4.8] Connect the footrest loom: in case of a left configuration with the red wires **[E20]**. In case of a right configuration, switch footrest connectors and connect the footrest loom with the blue wires **[E20+E6]**.

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




FINISH LIFT

MOUNTING REMOTE CONTROLS

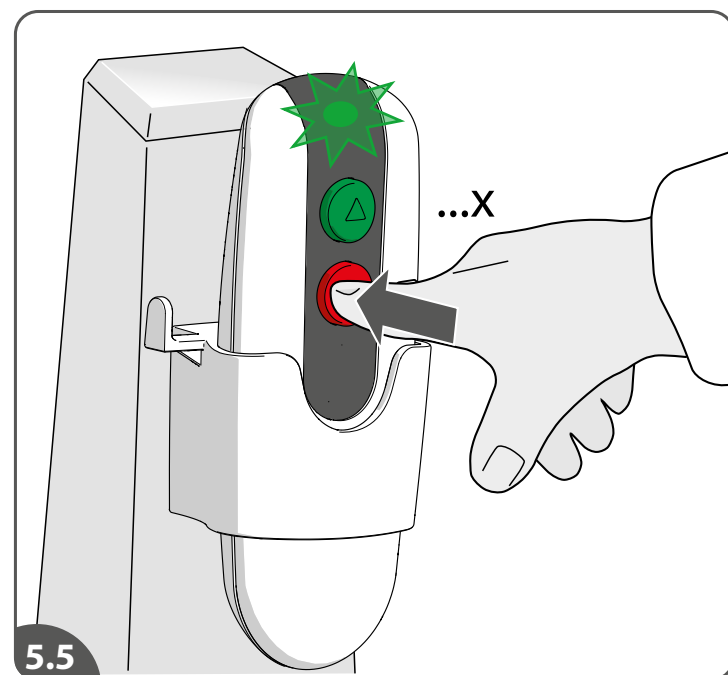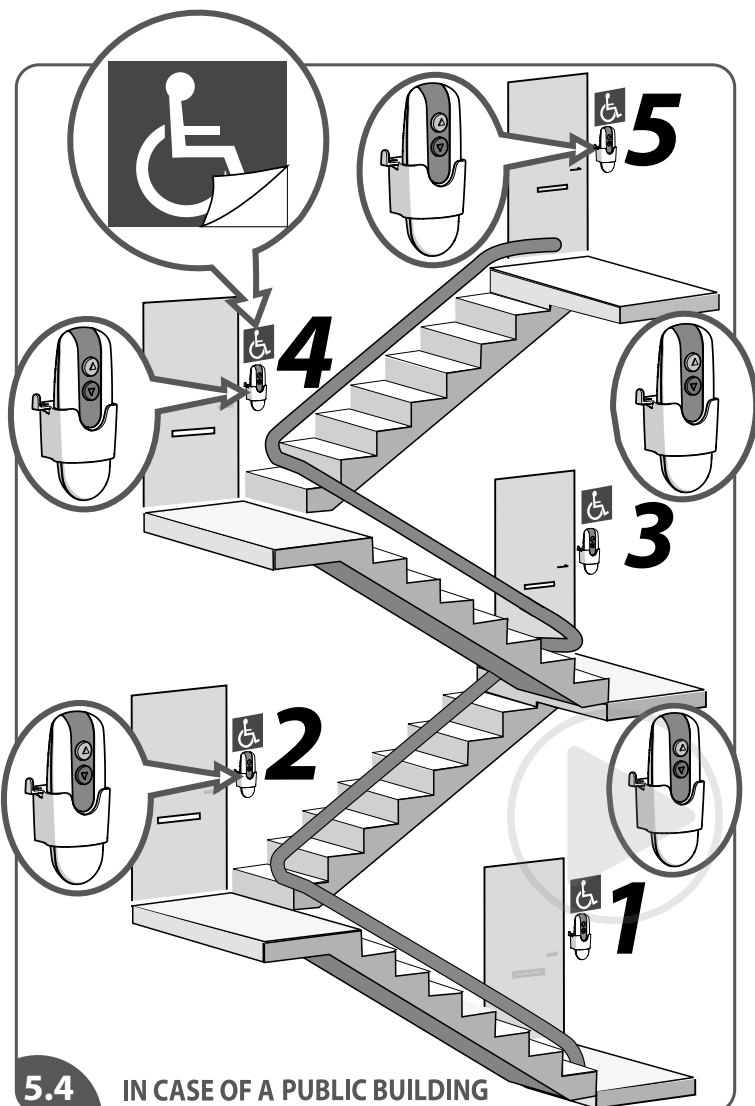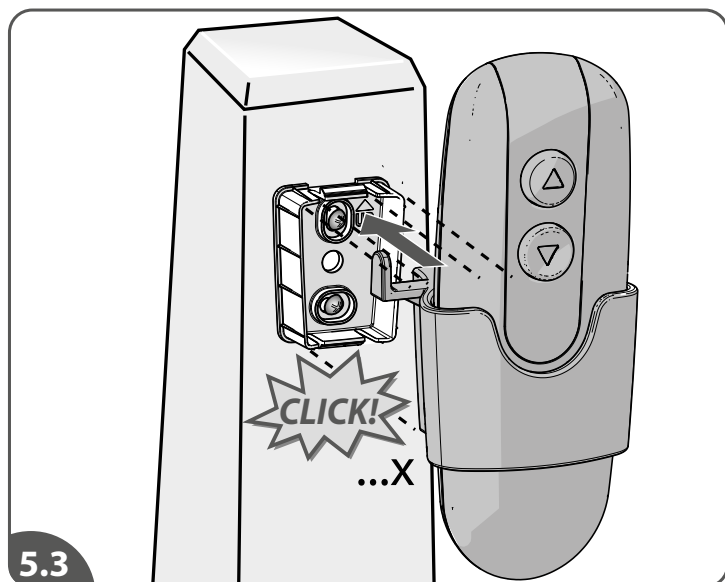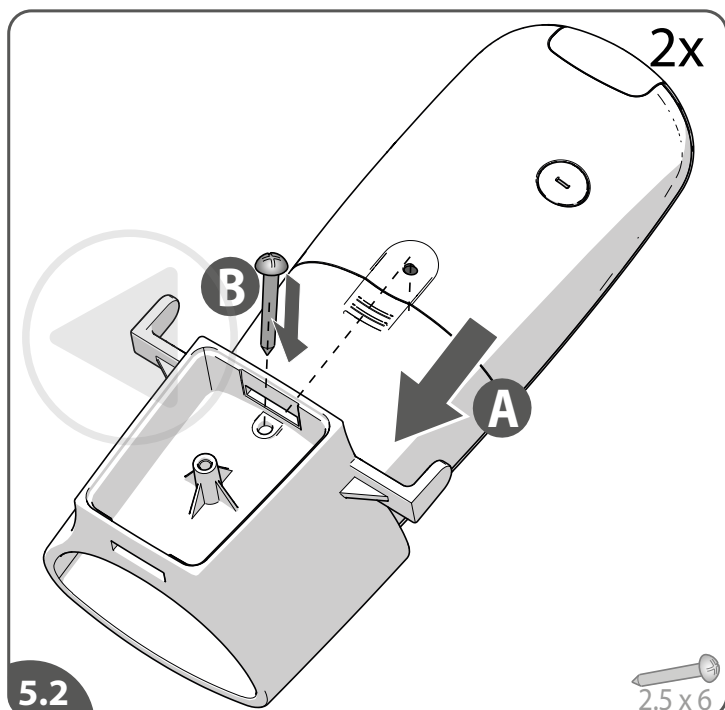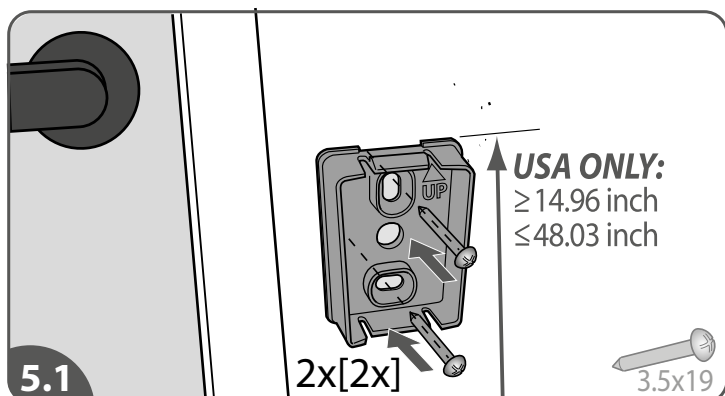- [5.1] Measure the position for the remote control bracket. Mount the bracket of the remote control with two screws on a solid background. There must be a minimum of 2 remote controls, with at least one at the top and bottom of the installation.
- [5.2] Mount the remote control into the holster{A} and fix it with a screw{B}.
- [5.3] Click the holster/remote control over the wall plate. Leave the remote control in the holster so that it does not get lost.
- [5.4] In case of a public building, place the disabled person-labels (supplied by dealer) near the remote controls.
- [5.5] Check if the remote controls are still connected.

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

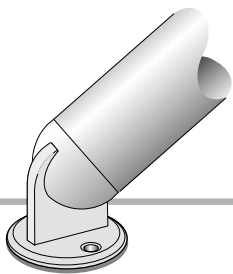
FINISH LIFT

TESTRIDE

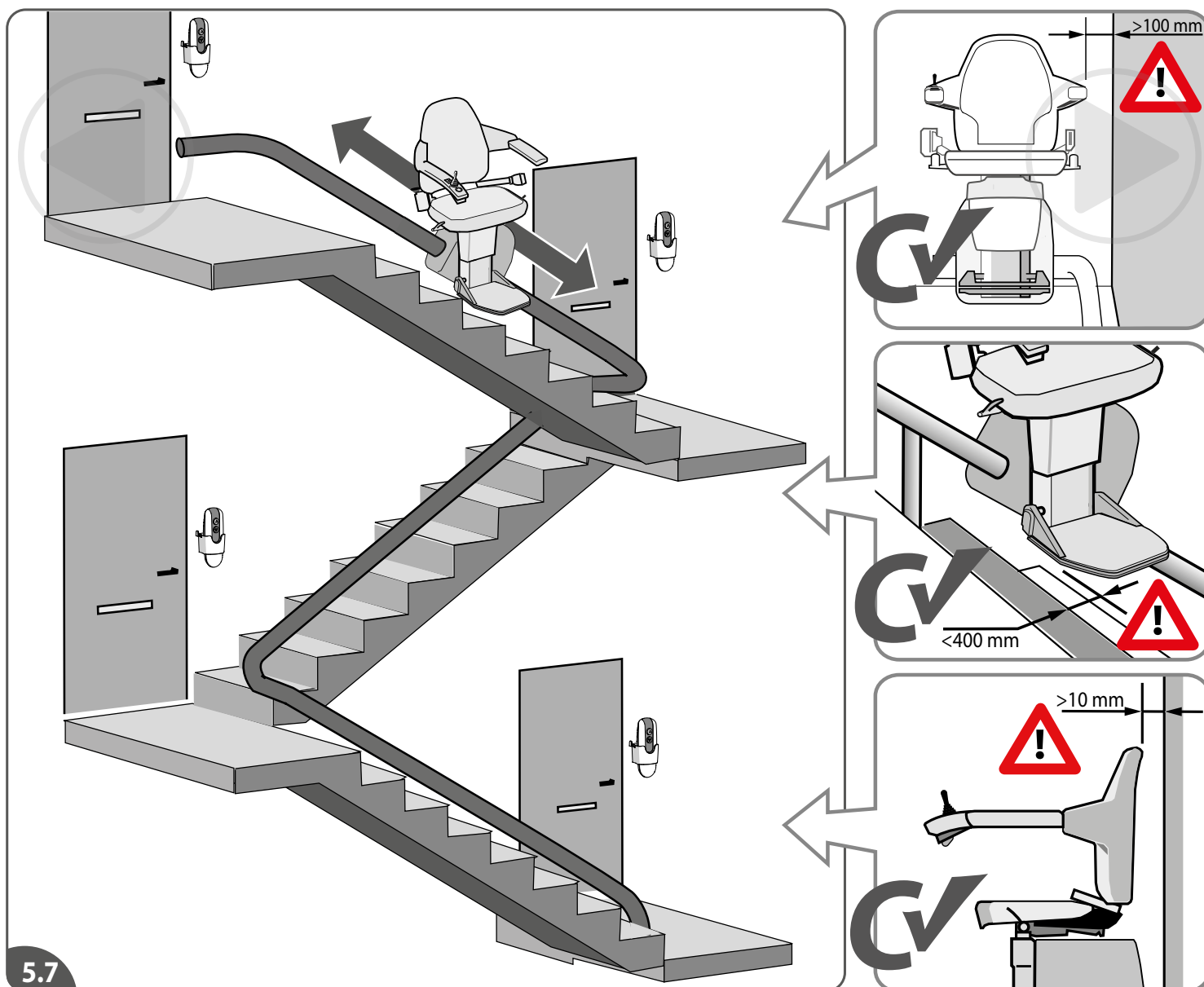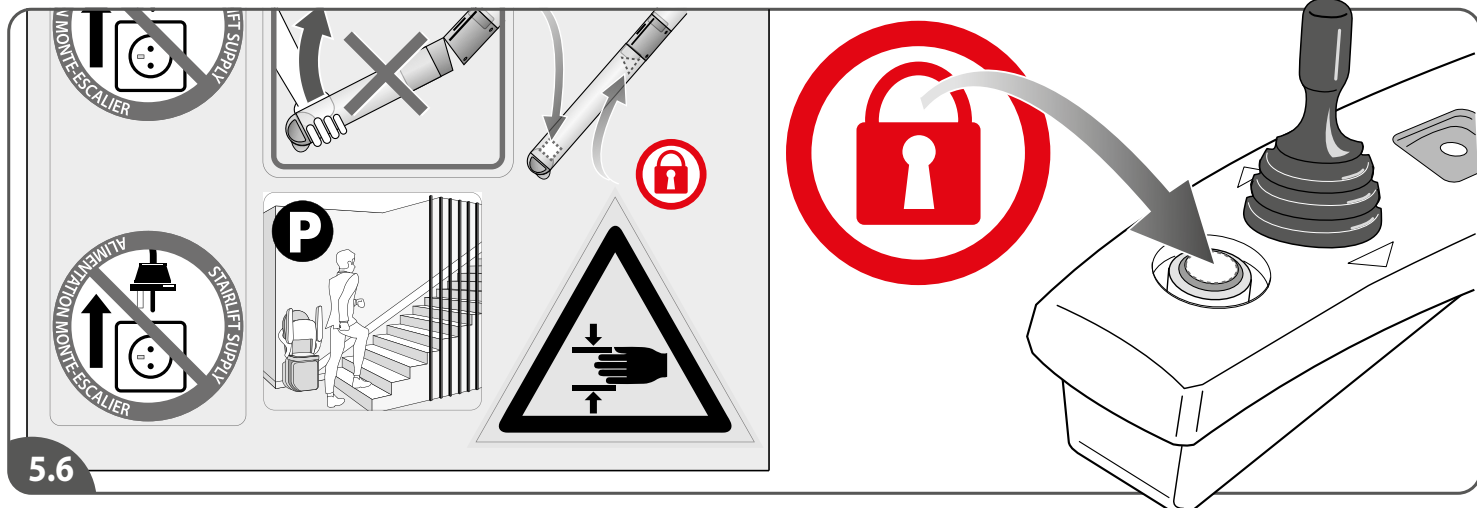- [5.6] Place the label of the emergency stop.
 [5.7] Make final check of the track position: ride the lift down and up.

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




IT

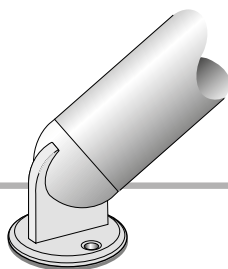TIP

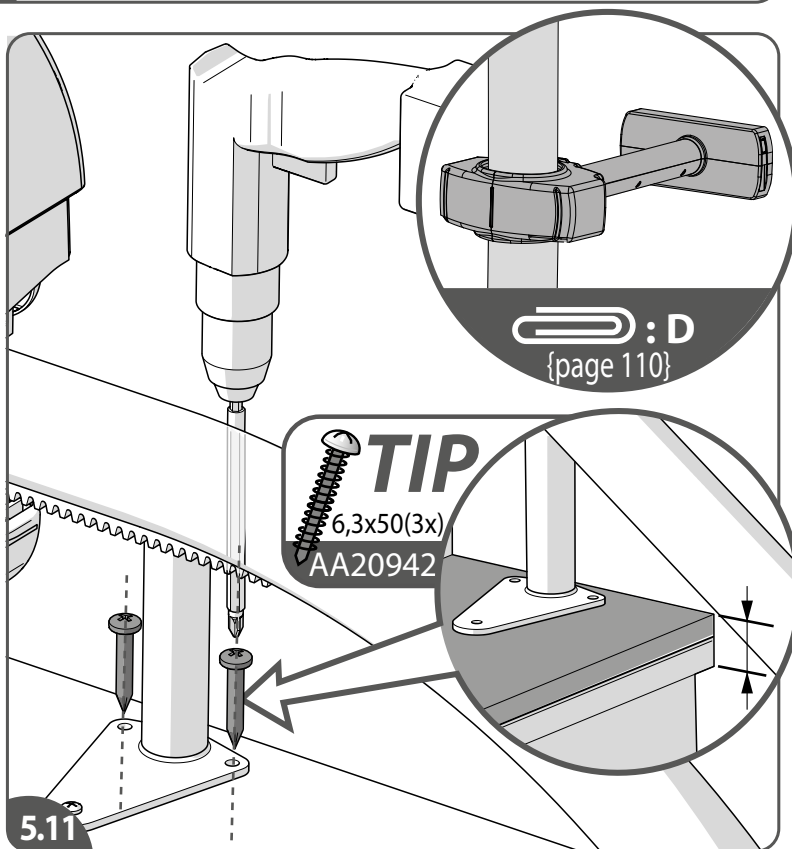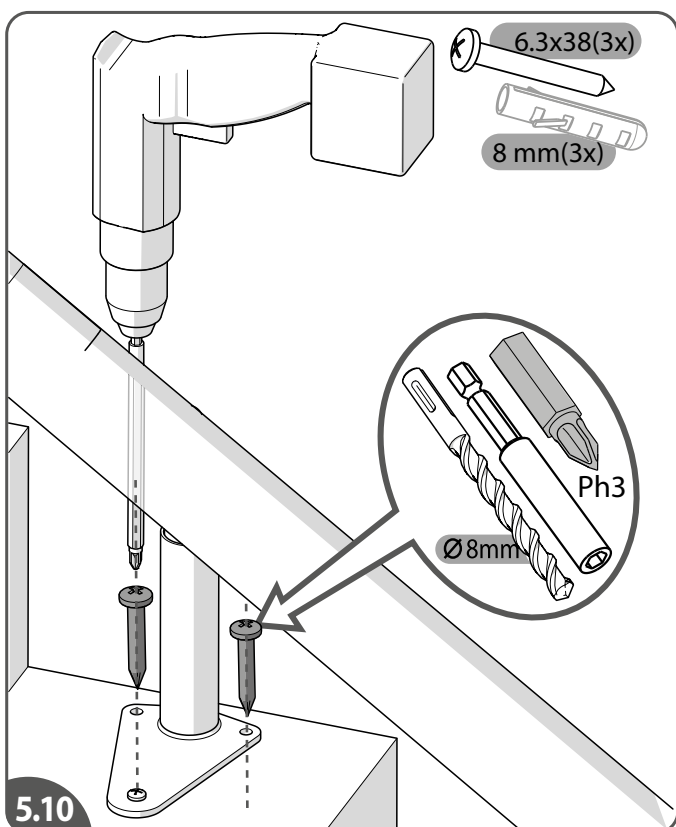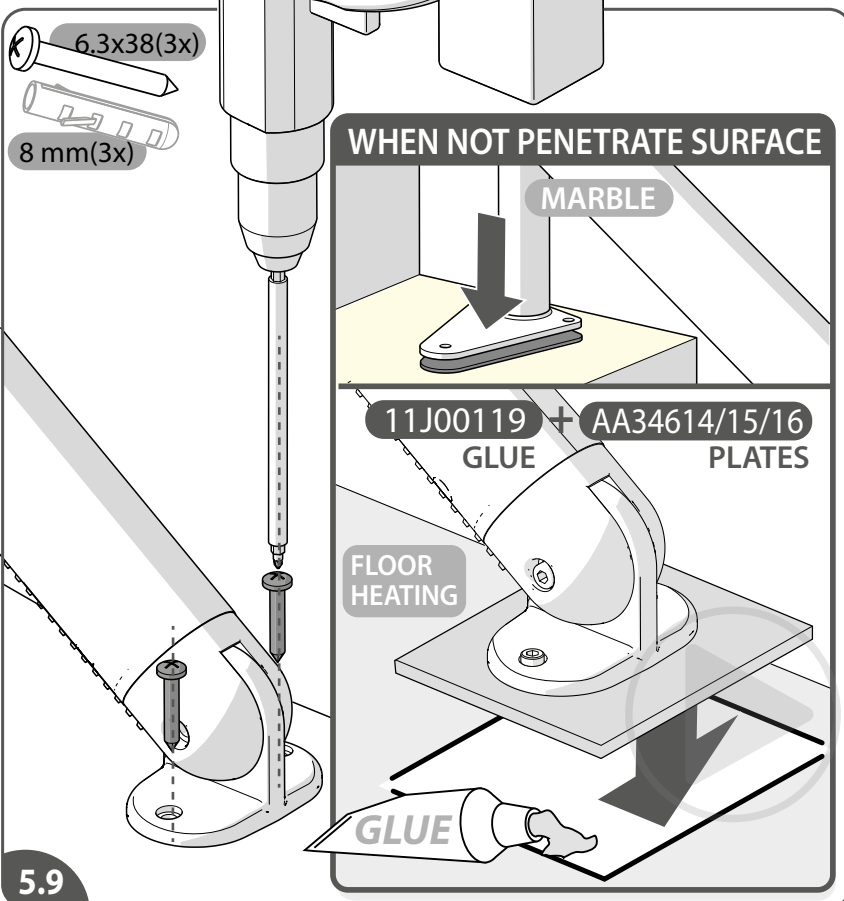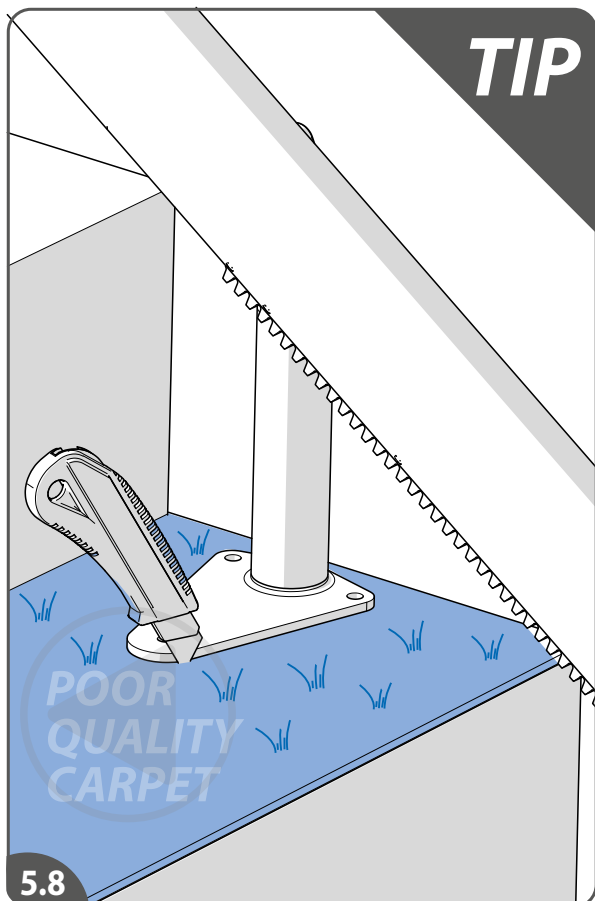[5.8] **TIP:** cut away bad quality carpet underneath the stanchion.

Make sure there is no space between the foot and the step before fixing the rail.

[5.9] Fasten all the feet to the steps/floor. In some cases (like a marble surface) it is necessary to glue the stanchion. Kits containing plates in various sizes/colours are available to increase the bonding surface.

[5.10] If needed, drill 8 mm holes in the wall/floor/stairs to fix the bracket or stanchion and use the plugs.

[5.11] In some cases (like **A**-rails) it is necessary to use 1 or 2 wall brackets to fix the stanchion (zie attachment **D**, p 110).

Longer screws will be necessary in case of thick carpet on the stairs.

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
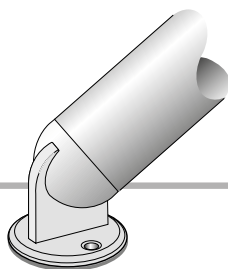
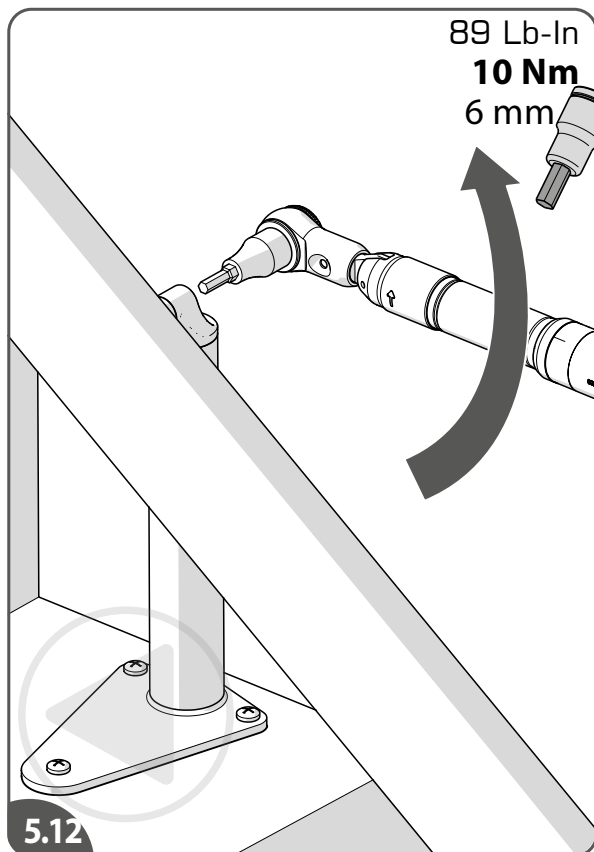
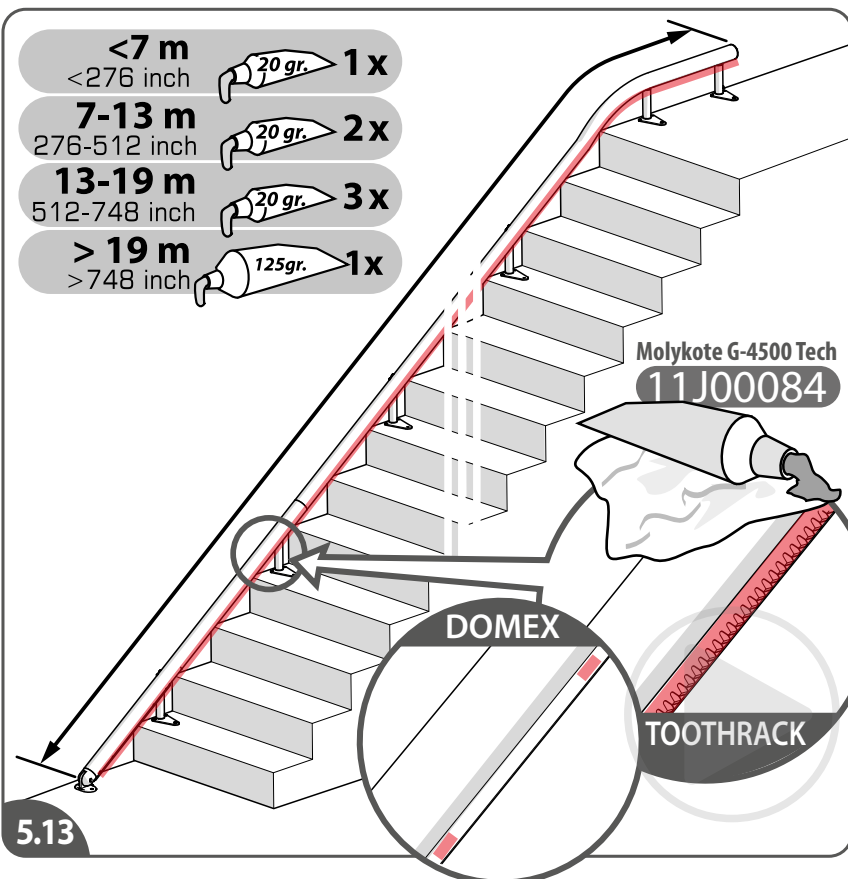
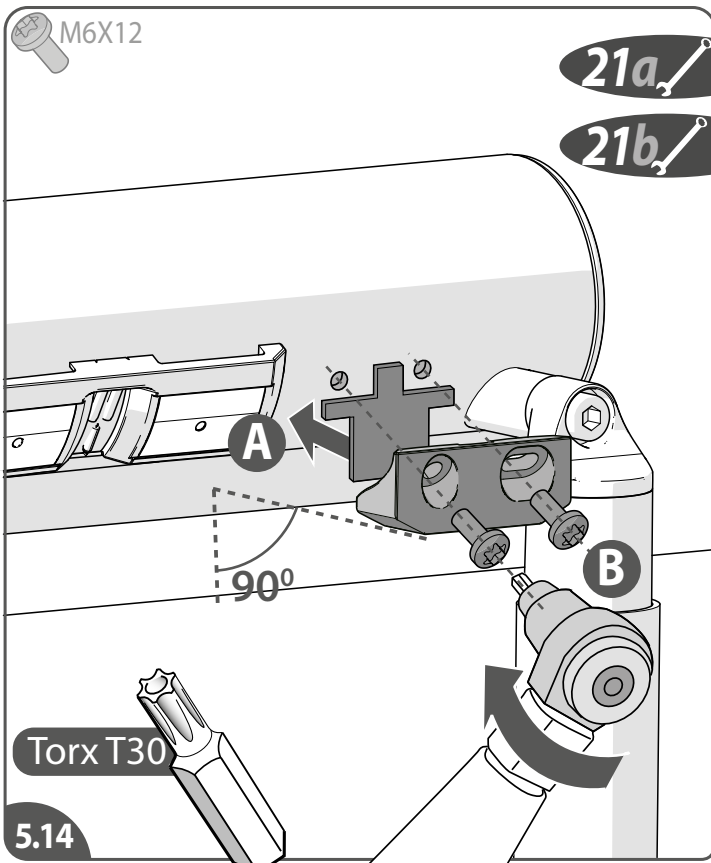
5.14

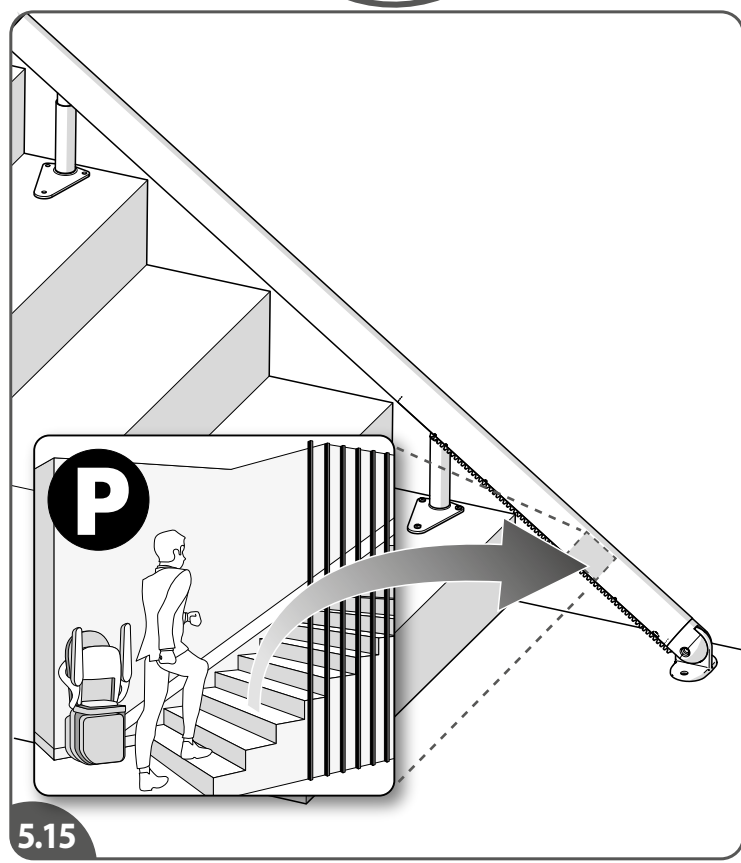5.15

- [5.12]Torque up all the rail bolts [10 Nm | 89 Lb-In].
- [5.13]Grease the toothrack over the entire length of the rail. Grease the domex on a few spots a little bit. For the quantity to be applied, see the illustration.
- [5.14]Position the bracket{A}. Install the endstop (90°) and fasten the endstop{B}.
- [5.15]Position the label on the parking/loading place where there is minimal 50 cm free walking space.

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



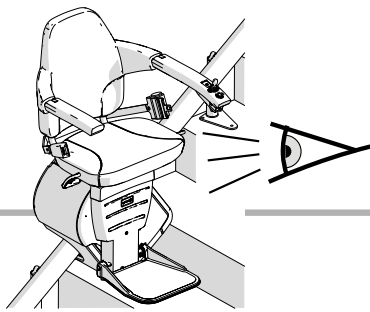
TORQUE WRENCH TEST

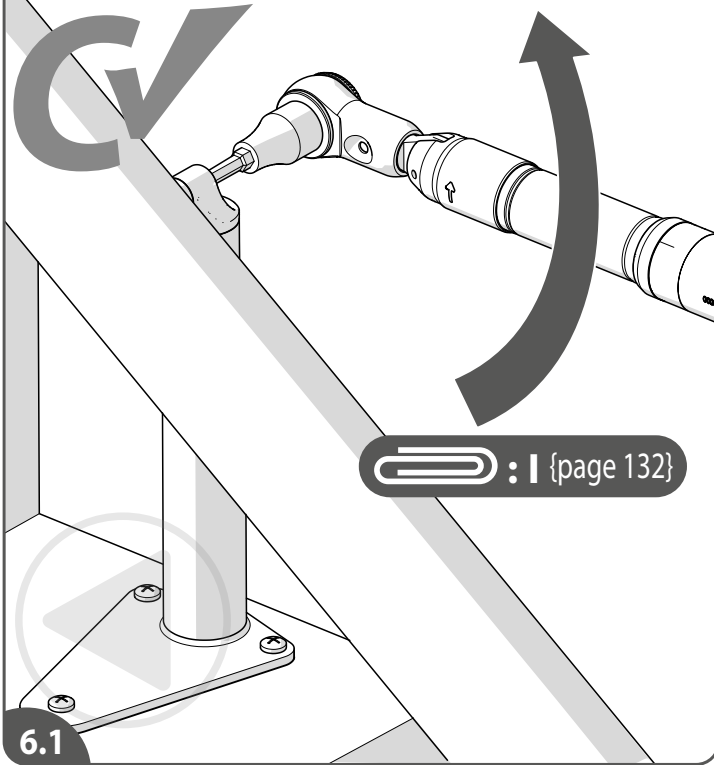: I {page 132}

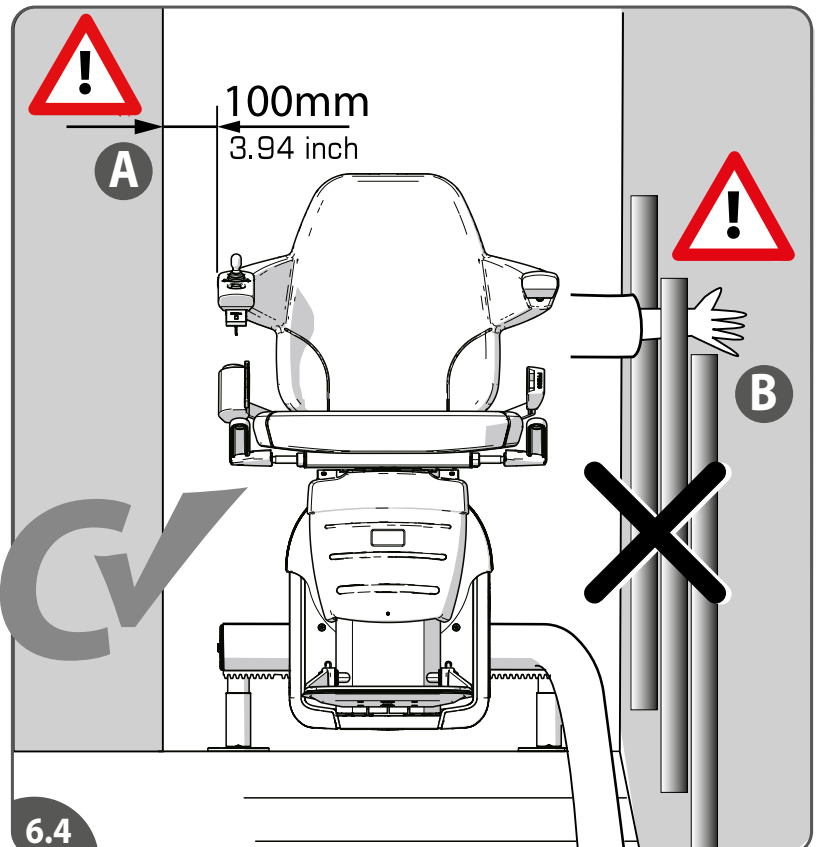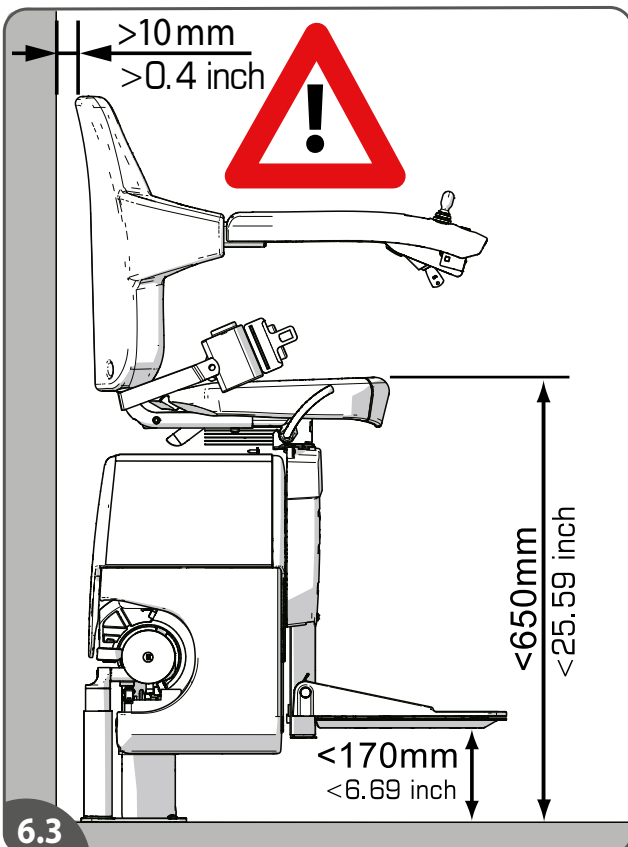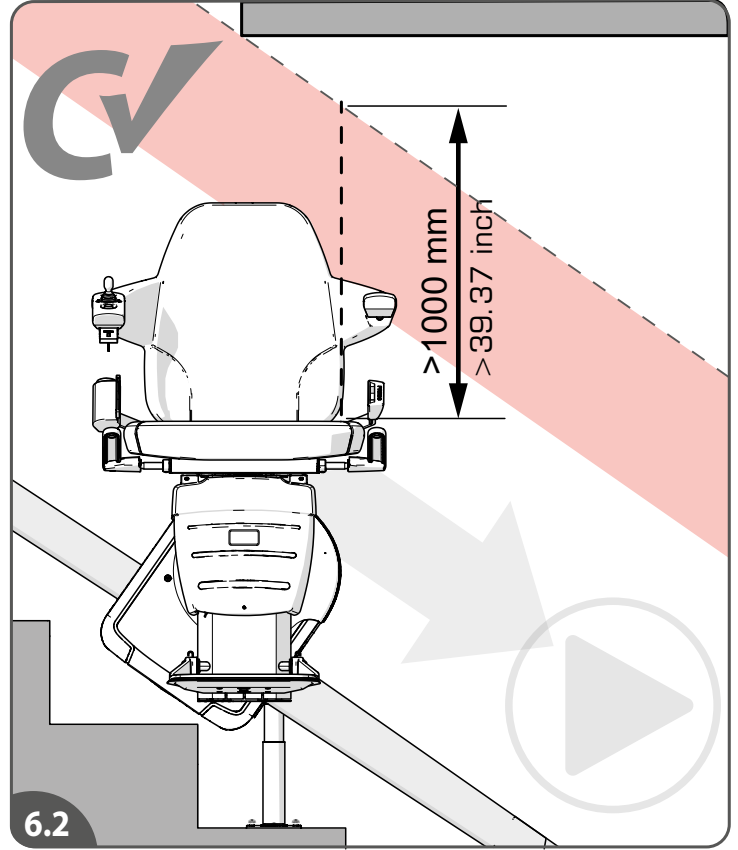- [6.1] Check the torque of the load-carrying fasteners. See attachment **I** (page 132) how to do this.
- [6.2] Check the position of the bulkhead: check that there is no obstacle within 1000 mm (39.97 inch) above the seat in the direction of movement of the lift.
- [6.3] Check the chair dimensions: the backseat at least 10 mm (0.4 inch) from the wall, the seatbase maximum 650 mm (25.59 inch) above the floor and the footrest maximum 170 mm (6.69 inch) above the floor.
- [6.4] Make sure that there is nothing behind the stairlift that body parts could get trapped in or on.

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





6

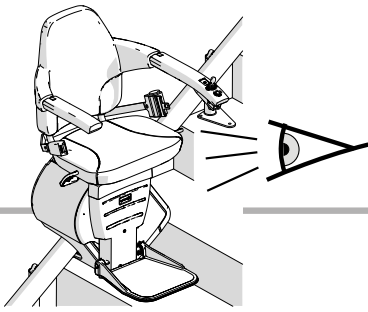CHECK LIFT SAFETY BELT

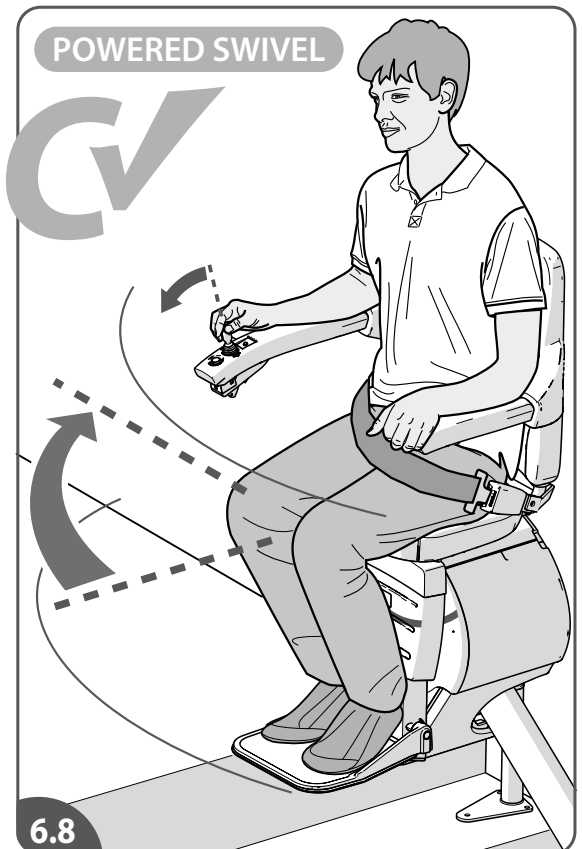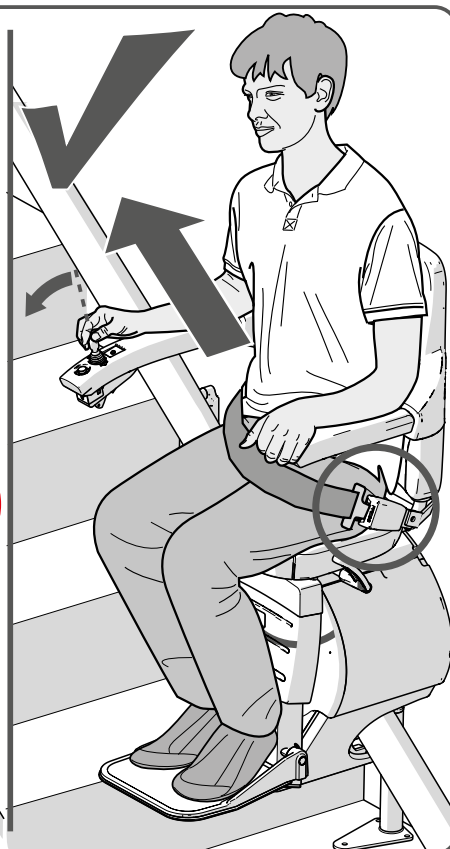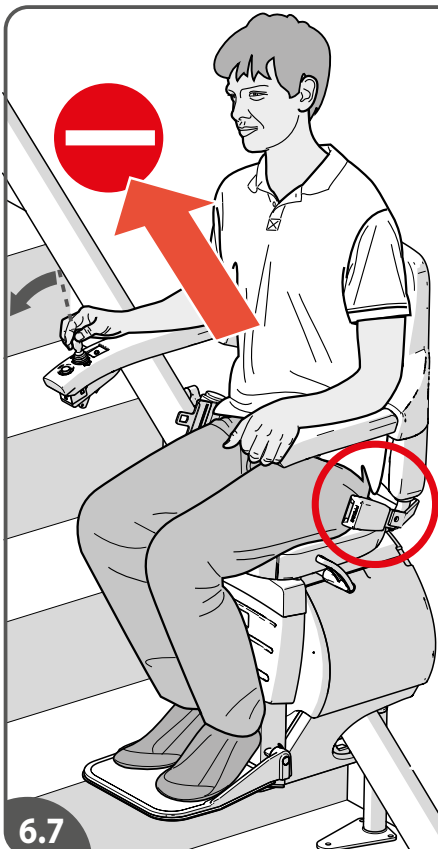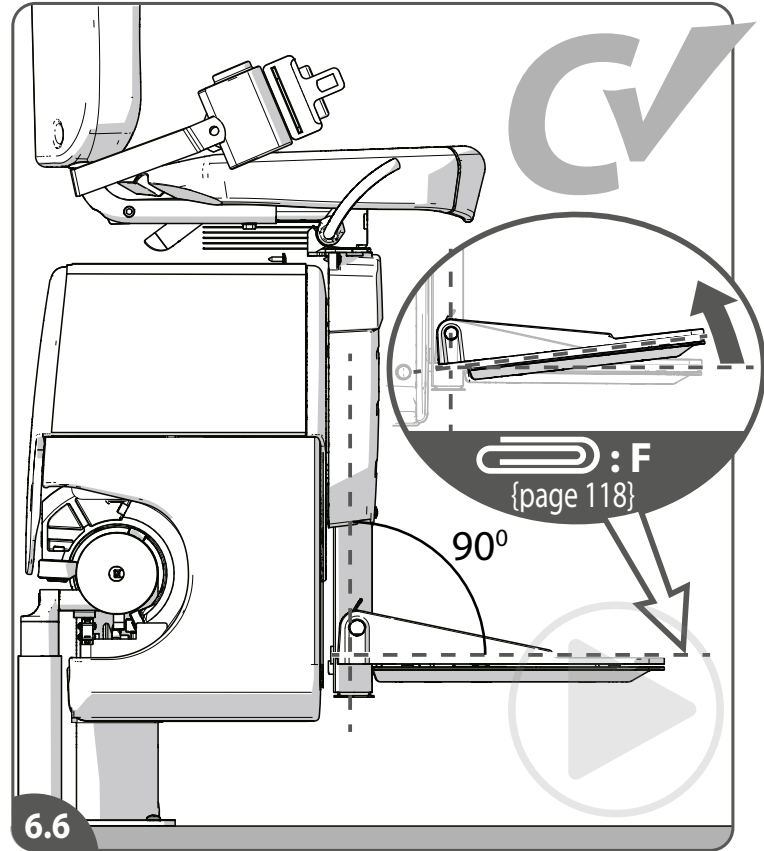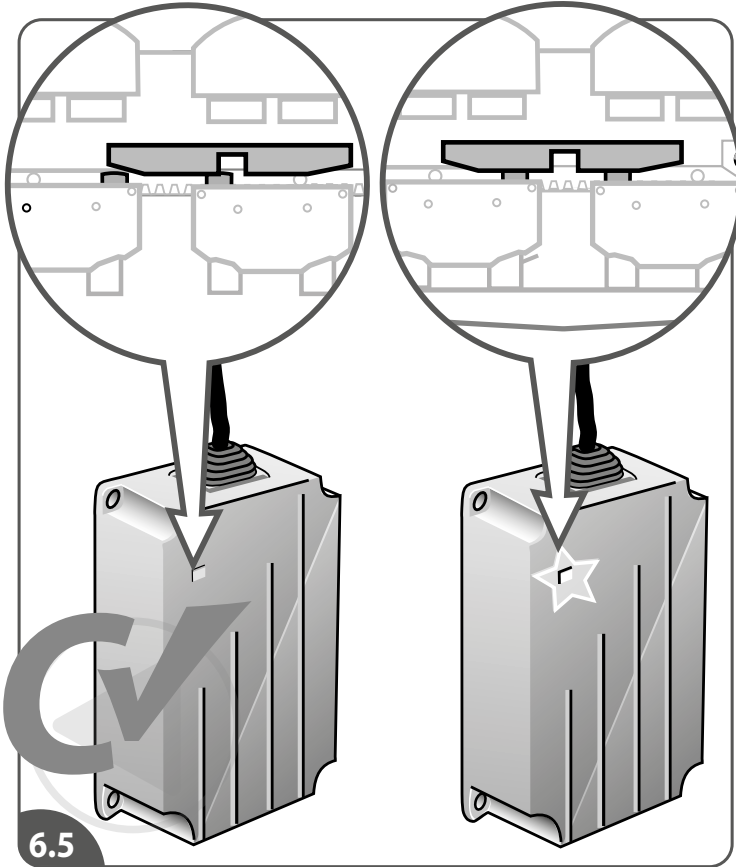- [6.5] Check that the lift is charging.
- [6.6] Check the position of the footrest.
- [6.7] Check the seatbelt.
- [6.8] Check the powered swivel function.

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


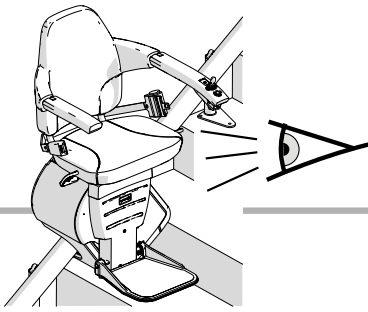
CHECK LIFT SAFETY DETECTION

MANUAL SWIVEL

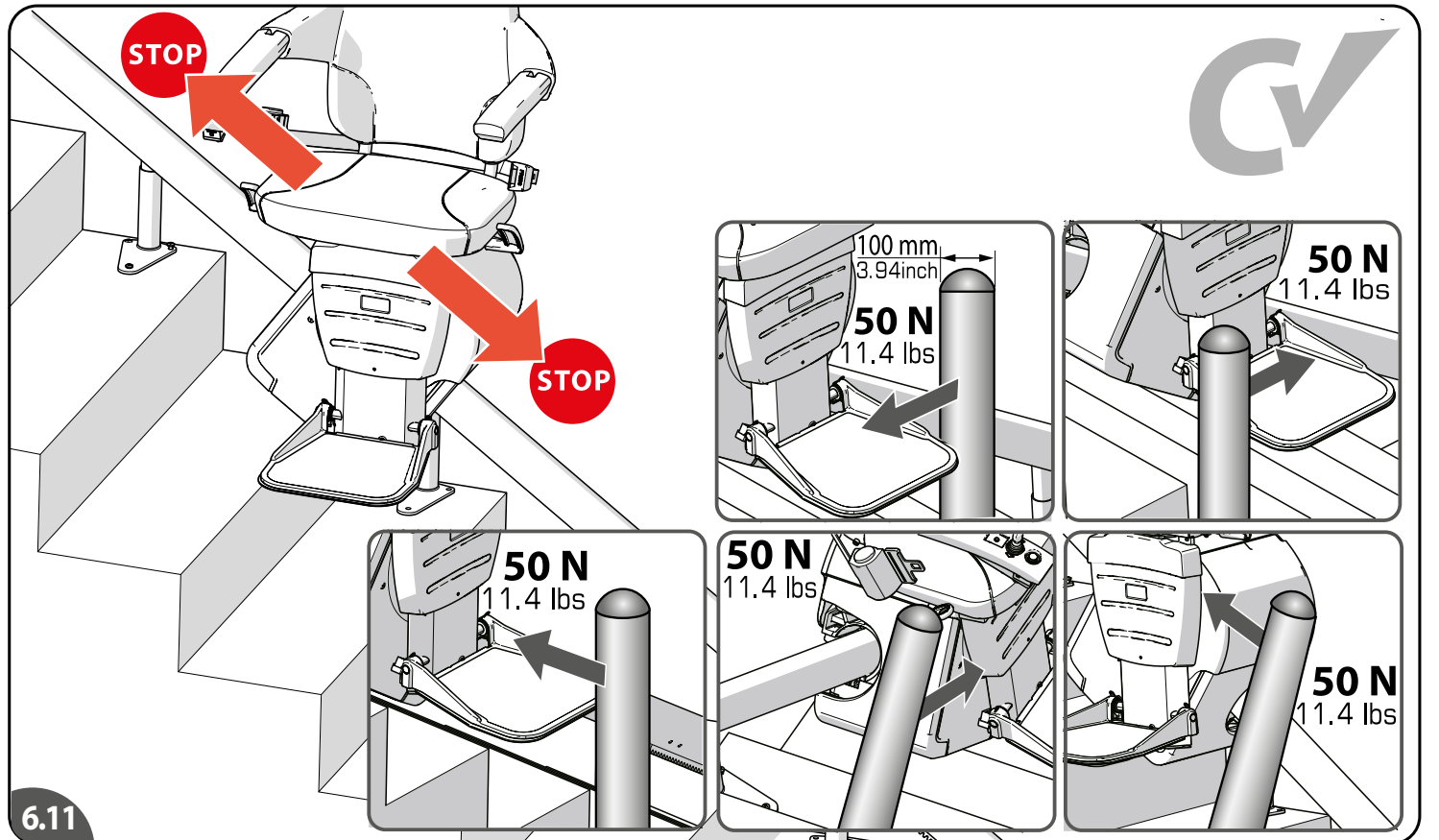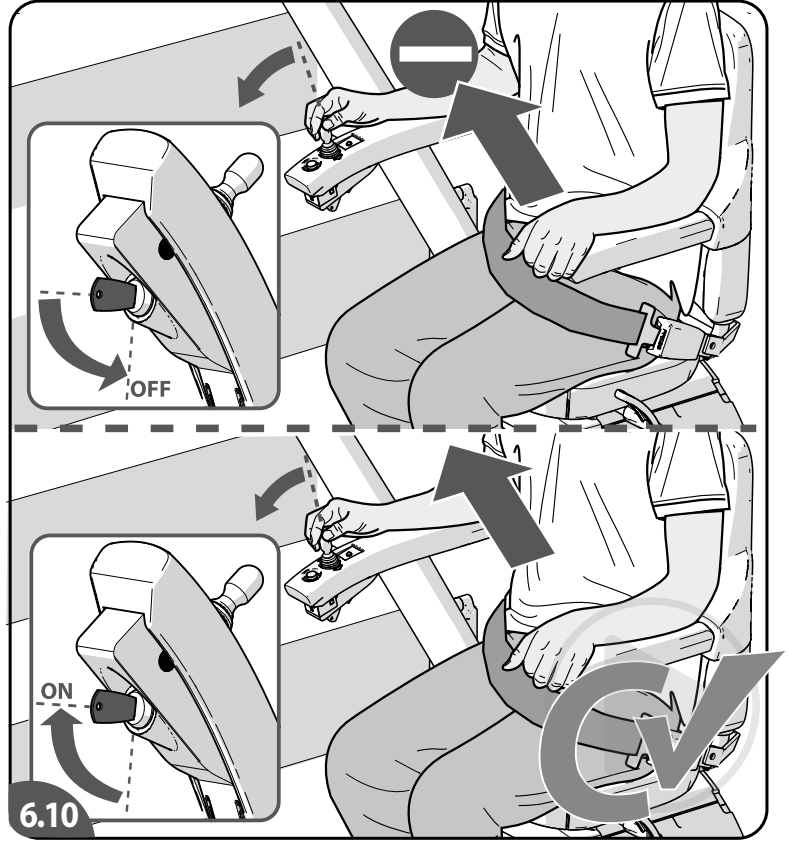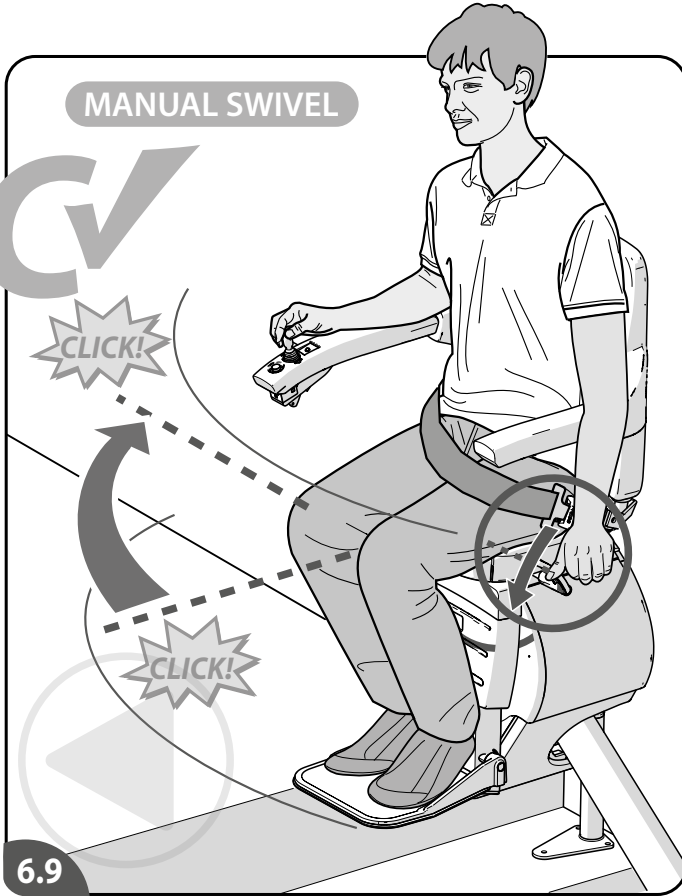[6.9] Check the manual swivel function and that the end stops work satisfactorily (it clicks at the end).
 [6.10] Turn the key to the OFF position and ensure that the stairlift will not move. Turn the key to the ON position. The stairlift should be able to move again.
 [6.11] Check the safety stops.

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


6

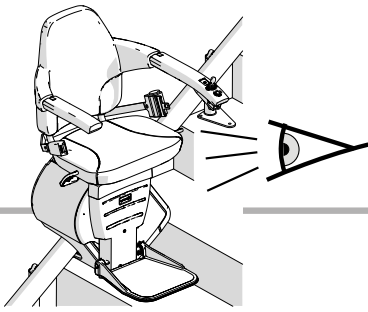CHECK LIFT EMERGENCY STOP

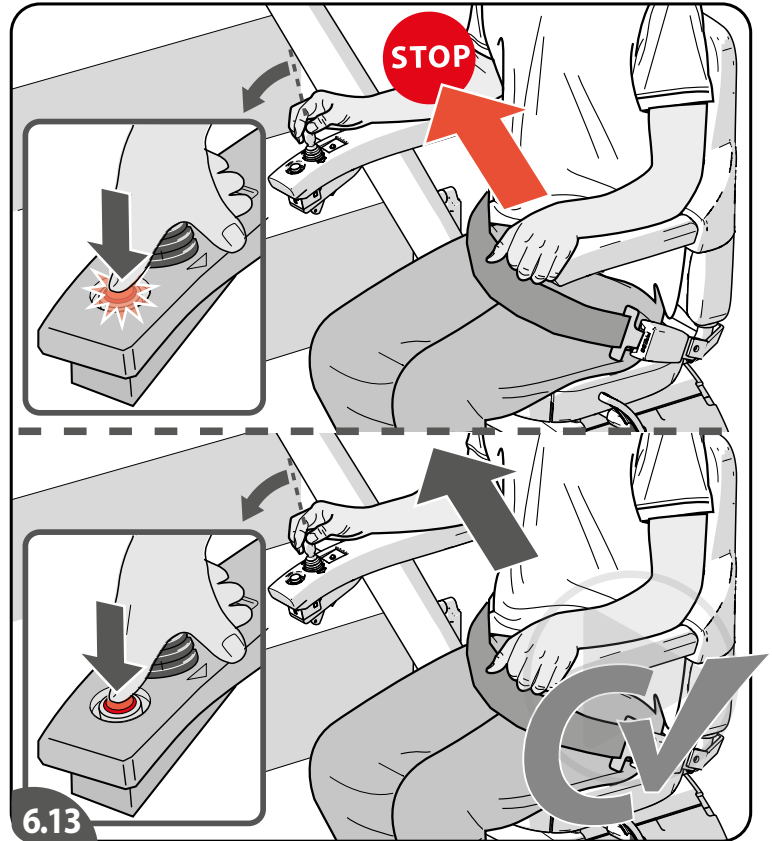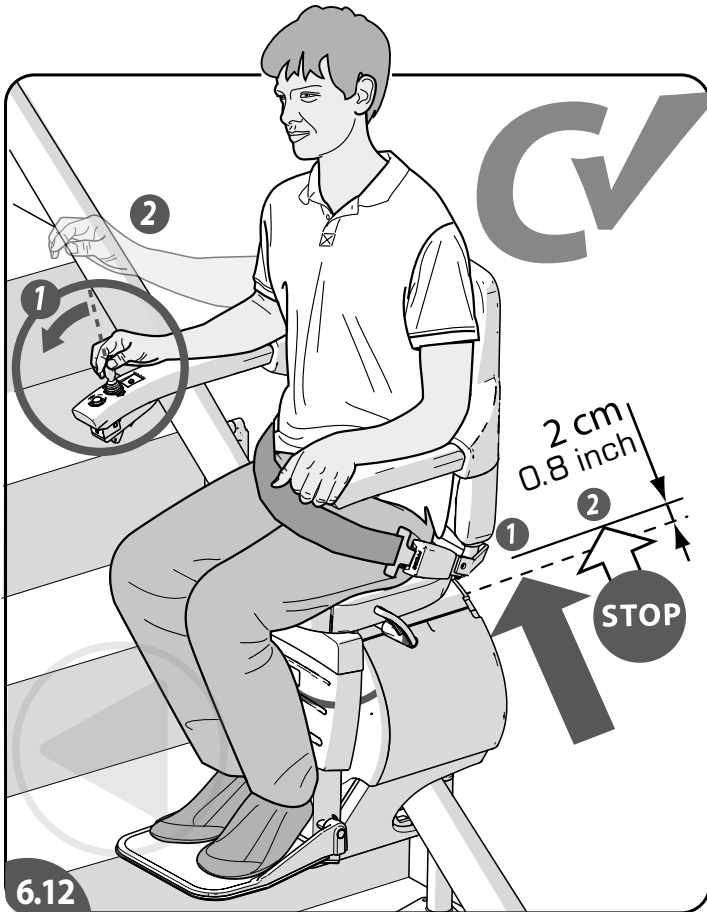- [6.12]Release the joystick. The stairlift should stop within 2 cm (0.8 inch).
- [6.13]Press the stop button. The lift must stop immediately. Press the stop button again; check that the stairlift operates normally.
- [6.14]Check that the stairlift makes a “soft” landing at the top and bottom of the rail.
- [6.15]Check the powered footrest function.

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


6

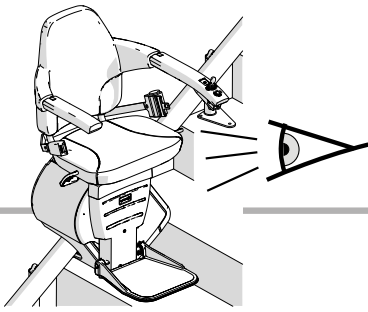CHECK LIFT ARMREST SAFETY

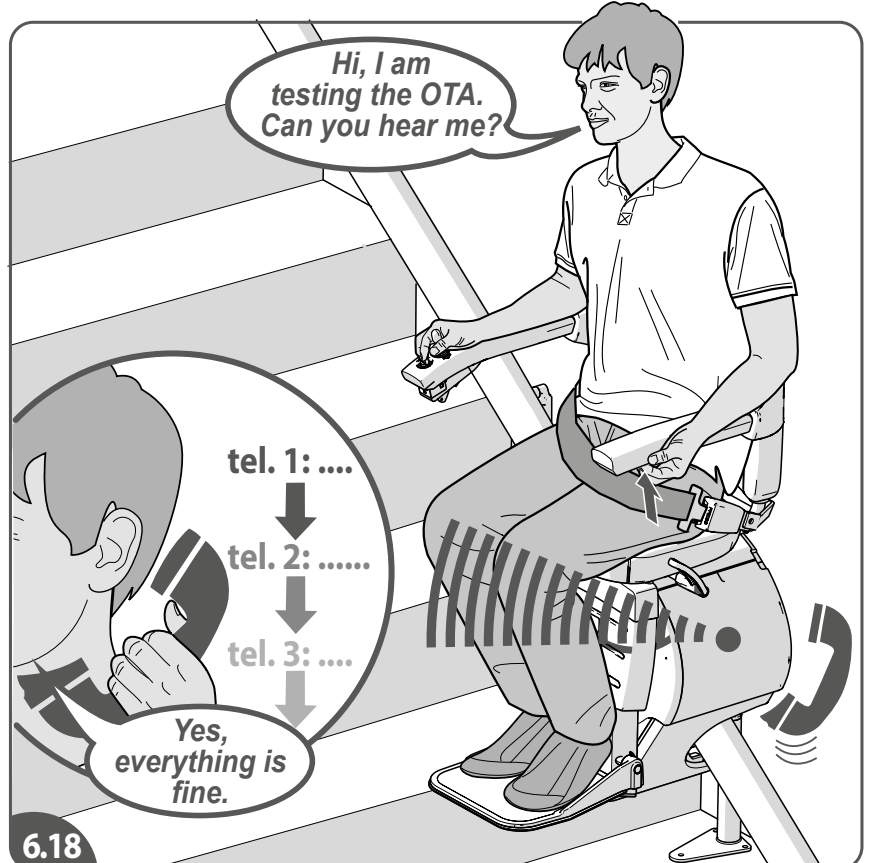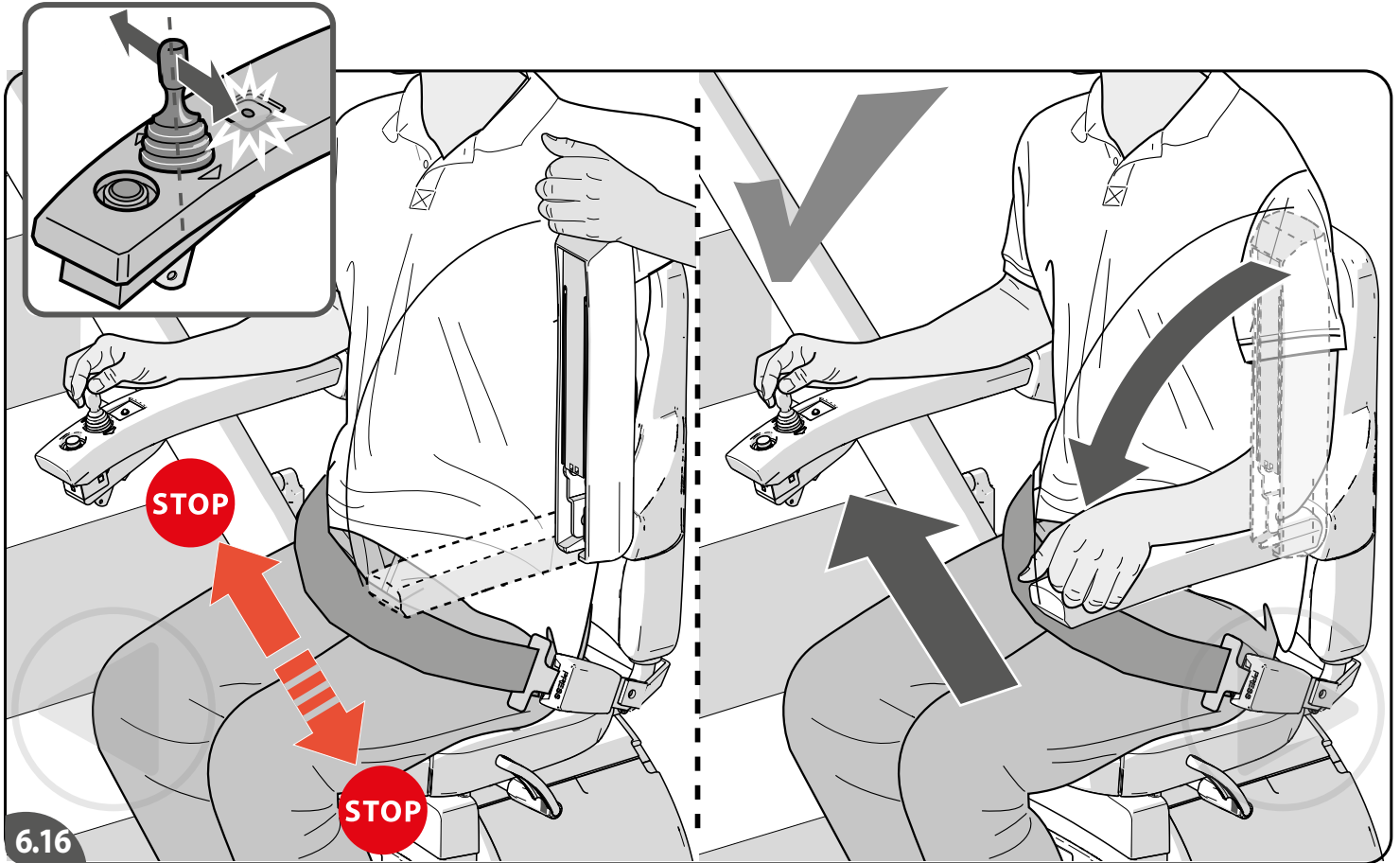[6.16] Fold the armrest up. The stairlift must stop immediately. Put the armrest back down and check that the stairlift operates normally again.

[6.17] While moving the stairlift, swivel the seat; the stairlift must stop immediately.

[6.18] It is the responsibility of the Dealer to ensure that the Integrated OTA System is checked and test calls are performed at installation and normal stairlift service intervals. As it is the responsibility of the customer (end-user) to provide and manage the Sim Card used within the Integrated OTA System.

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
6

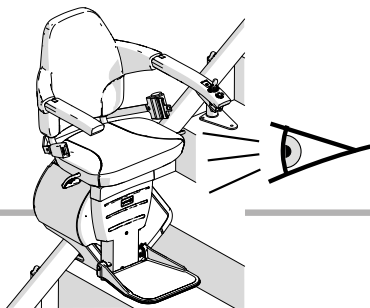CHECK LIFT
INSTRUCT CUSTOMER

[6.19] Instruct the user by telling him/her all users practice.

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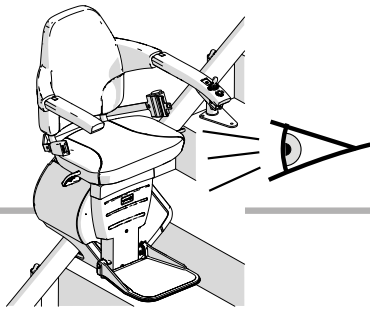
Location: Manufacturer: Type:

Unit serial no: Serial no:

1) Confirm that the installation data have been completed on the EC Declaration of Conformity

2) Confirm that the remote control is mounted on the wall and is functioning properly

3) Confirm that the safety belt is secured and functions correctly

4) Confirm that all warnings according to the maximum working load on the lift have been affixed

2 Rail

1) Confirm that the rail joints are properly secured

2) Confirm that the rail supports are properly secured

3) Confirm that the folding rail is functioning correctly

3 Electrical drive unit tests

1) Confirm that the batteries are securely fitted in the carriage (figure)

2) Confirm that the following have been effectively indicated to the user:

a) location of the charge contacts

b) "OFF CHARGE" alarm

c) required action when stairlift is not charging

4 Clearance dimensions

1) Confirm the correct clearance dimensions as specified in the installation manual are maintained throughout the full travel of the stairlift

2) Confirm that there are no open bars/openings at the rear of the lift

3) Confirm that other users can safely use the staircase (minimum clearance of 400 mm between the folded seat and opposing obstacles). Install handrails if required

5 Functional tests

1) Confirm that the final limit switches (upper and lower) operate correctly (position and operation)

2) Confirm that the control lever, key switch, emergency stop button, and the detection switches in the armrests operate correctly

3) Confirm that the safety detection devices operate correctly (footrest/unit/armrest)

4) Confirm that the safety seat interlock and switch operate correctly

5) Confirm that the stopping distance of the lift is a maximum of 2 cm

6) Confirm that the restart time of the lift is at least 1 second

7) Confirm that the alarm device operates correctly when activated

8) Confirm that the lift has successfully passed the static (torque wrench test) and dynamic test ride in accordance with the instructions in the installation manual

1) Confirm that the installation data

Location:

Unit serial no:

Signature: M

{page 138}

6.21

DECLARATION ENGINEER

I have thoroughly examined this lift and that the foregoing is a correct report of the lift compound with approved components and installed according to the instructions that I/we have instructed the user by telling him/her all users practice without further instruction. The user manual has been handed over to the user.

SEE USERMANUAL

Name Engineer:

Qualifications:

Signing date:

Signature:

If employed by a company/association > name/address:

DECLARATION INSTRUCTED PERSON

I declare that the above mentioned engineer(s) has declared all with the complete operation as instructed person.

Address:

Phone number:

Signature:

6.22

DECLARATION ENGINEER

I have thoroughly examined this lift and that the foregoing is a correct report of the lift compound with approved components and installed according to the instructions that I/we have instructed the user by telling him/her all users practice without further instruction. The user manual has been handed over to the user.

Name Engineer:

Qualifications:

Signing date:

Signature:

If employed by a company/association > name/address:

DECLARATION INSTRUCTED PERSON

I declare that the above mentioned engineer(s) has declared all with the complete operation as instructed person.

Address:

Phone number:

Signature:

DECLARATION OF CONFORMITY

I declare that the Freecurve is in conformity with the provisions of the Machinery Directive 2006/42/EC and EMC Directive 2014/30/EU.

in accordance with the following standards:

EN 81-40, EN ISO 14121, EN 60204-1, EN 60555

Year of build: 2023

Date put into service:

6.23

DECLARATION ENGINEER

I/We certify that on I/we thoroughly examined this lift and that the foregoing is a correct report of the lift compound with approved components and installed according to the instructions that I/we have instructed the user by telling him/her all users practice without further instruction. The user manual has been handed over to the user.

Name Engineer:

Qualifications:

Signing date:

Signature:

If employed by a company/association > name/address:

DECLARATION INSTRUCTED PERSON

I declare that the above mentioned engineer(s) has declared all with the complete operation as instructed person.

Address:

Phone number:

Signature:

EC DECLARATION OF CONFORMITY

I declare that the Freecurve is in conformity with the provisions of the Machinery Directive 2006/42/EC and EMC Directive 2014/30/EU.

in accordance with the following standards:

EN 81-40, ISO 9386-2, EN ISO 12100-1, EN ISO 12100-2, EN ISO 14121, EN 60204-1, EN 60555 (including the EMC standards) ASME A18.1, CSA B-355

Heerhugowaard, The Netherlands, 15 april 2023.

Handtekening project:

Year of build: 2023

The declaration of conformity is legal after the installation checklist and the declaration instructed person. The installation checklist and the declaration instructed person shall be saved by the dealer.

[6.20] A full ride is made along the complete travel, include hinged rail guide sections and swivelling of the seat where applicable; fill in the installation checklist.

[6.21/22/23] The declaration of conformity is legal after the declaration engineer and the EC-declaration of conformity have been signed off by the engineer and the user. The declaration engineer and the EC-declaration of conformity shall be saved by the dealer and made available to Savaria on request.

EN ***Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to the requirements of the latest adopted ASME A18.1 code. Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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NL ***Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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ES ***Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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DE ***Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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
FR ***Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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IT ***Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.***

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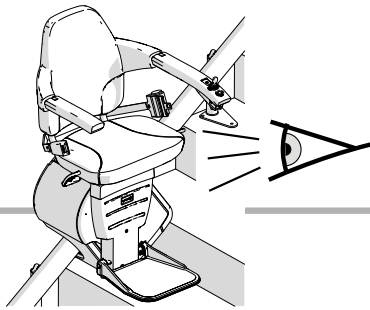DRIVE: RACK & PINION
RATED LOAD: 276lbs (125kg)
RATED SPEED:
21.7 ft/min (0.11m/s)

CONFORMS TO ASME A18.1 & ASME A17.5
CERTIFIED TO CSA B44.1 & CSA B355
DESIGNED AND TESTED IN ACCORDANCE WITH EN 81-40
Indoor use only. Utilisation intérieure seulement.
One person only. Une seule personne.
Stairlift for persons with physical disabilities.
Monte-escalier pour personnes handicapées.
ASME A18.1 TO BE USED FOR INSPECTIONS AND TESTS

RATED INPUT FROM CHARGER: 1.5A/24V DC
OUTPUT TO MOTORS: 30A/24V DC
DUTY CYCLE: 50 50%
IP RATING: IP12
DATE OF INSTALLATION: _____
NAME OF INSTALLER: _____

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2026 **2027** **2028** **2029**

AA31201

[6.24] Take the right warning label out of the manual and put the correct stickers for year and weight on the label {A}. Mount the label on the chassis {B}.

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



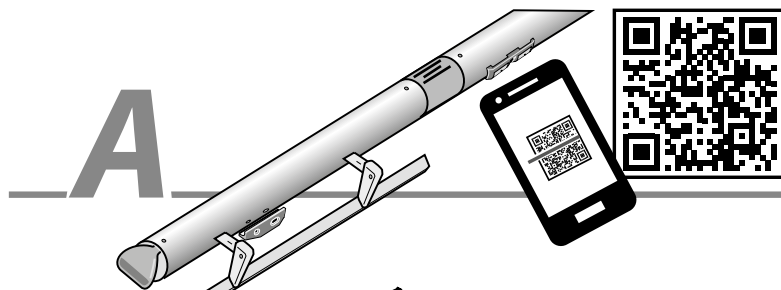
HINGE

CHARGE POINT/HINGEPOINT

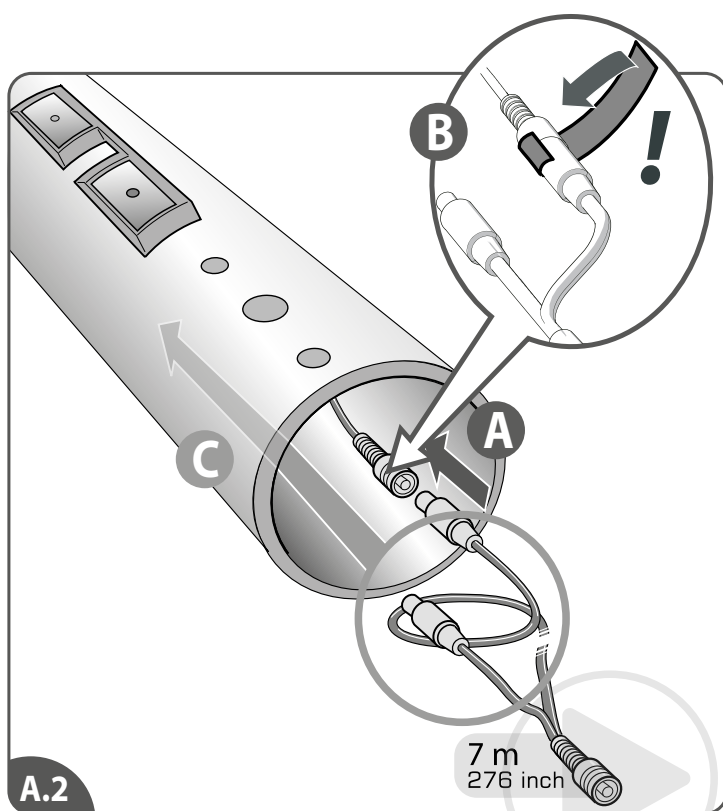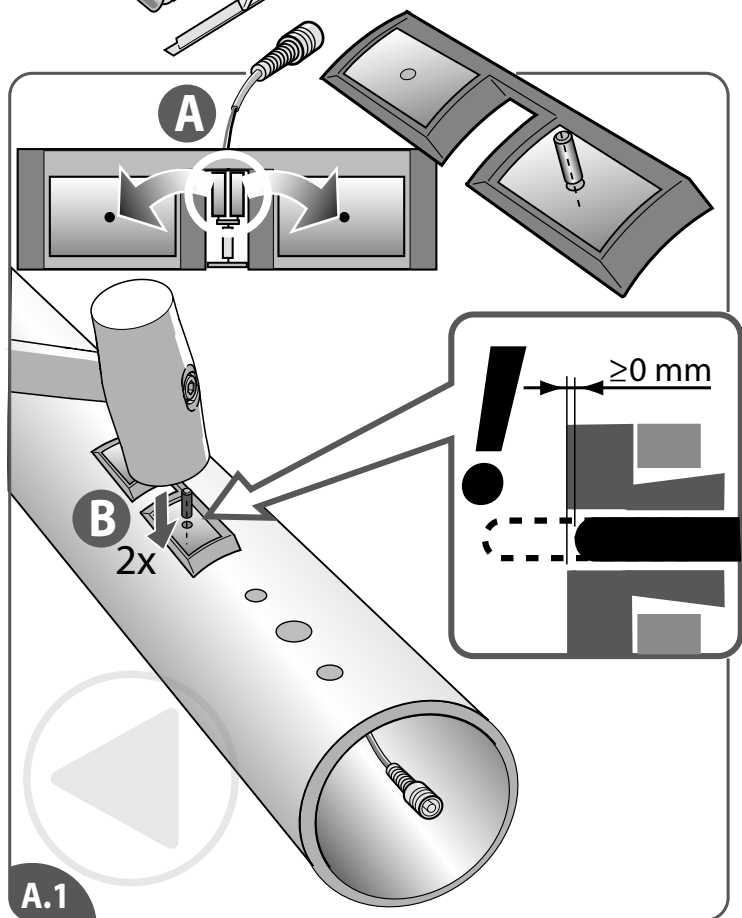INSTALL THE HINGED RAIL PART IN THE SAME WAY AS AN ORDINARY RAIL PART. SEE CHAPTER 2 FOR A DETAILED INSTALLATION.

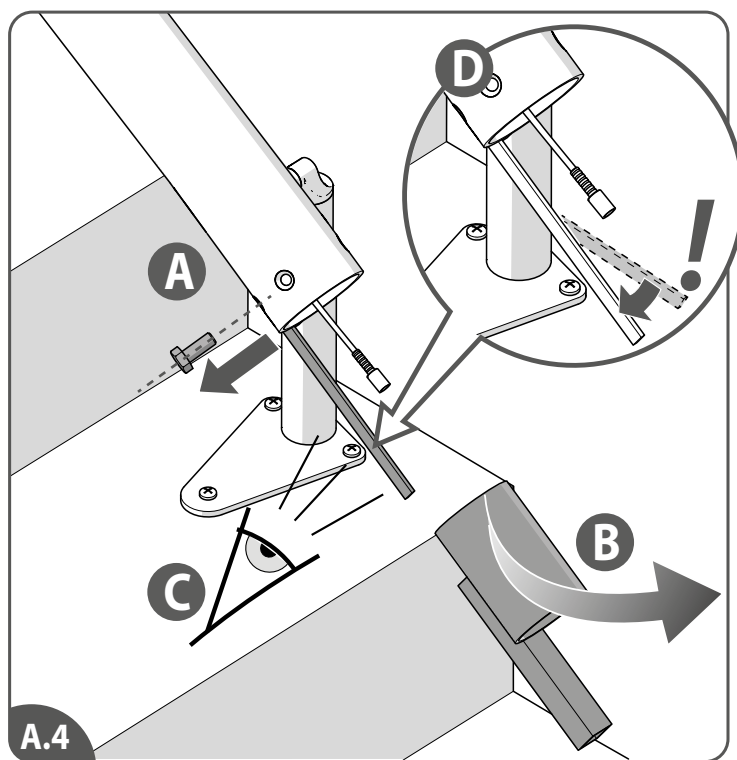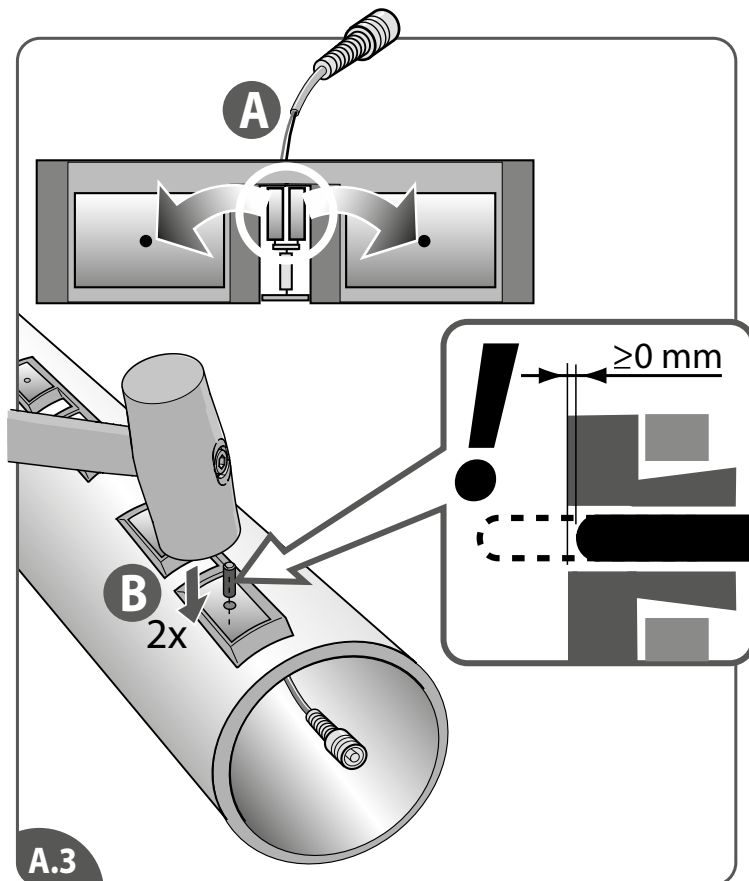[A1]Install the bottom charge contact: remove the pins{A}. Run the cable into the hole of the rail, position the charge contact on the rail and fix it. Gently hit the pins in the charge contact {B}.

[A2]Connect the charging wire{A}. Secure the connection with tape{B}. Run the charging wire through the rail{C}.

[A3]Install the hinged rail point: remove the pins{A}. Run the cable into the hole of the rail, position the hinged rail point on the rail and fix it. Gently hit the pins in the hinged rail point {B}.

Install the hinged rail part in the same way as an ordinary rail part. See chapter 2 for a detailed installation.

EN [A4]Remove the toothrack protection{A+B}. Inspect the domexstrip{C}. If needed, remove the bending in the strip{D}.

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



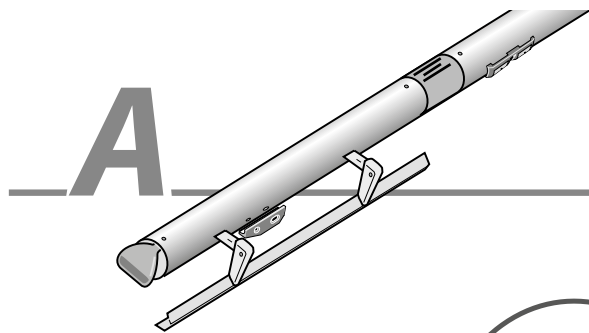
HINGE MOTOR

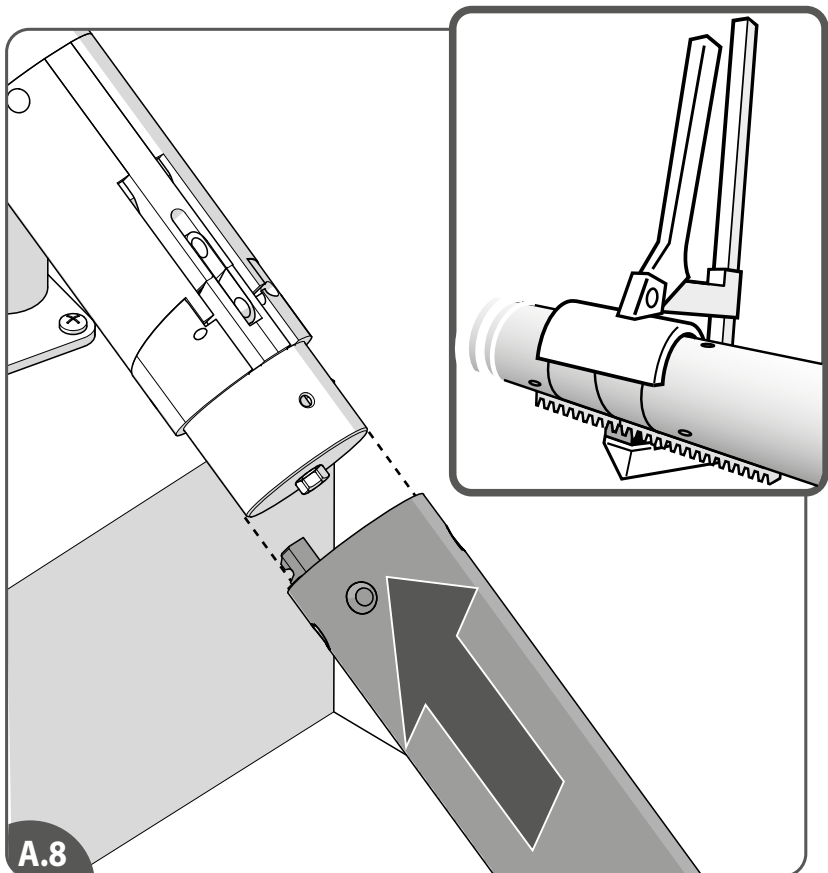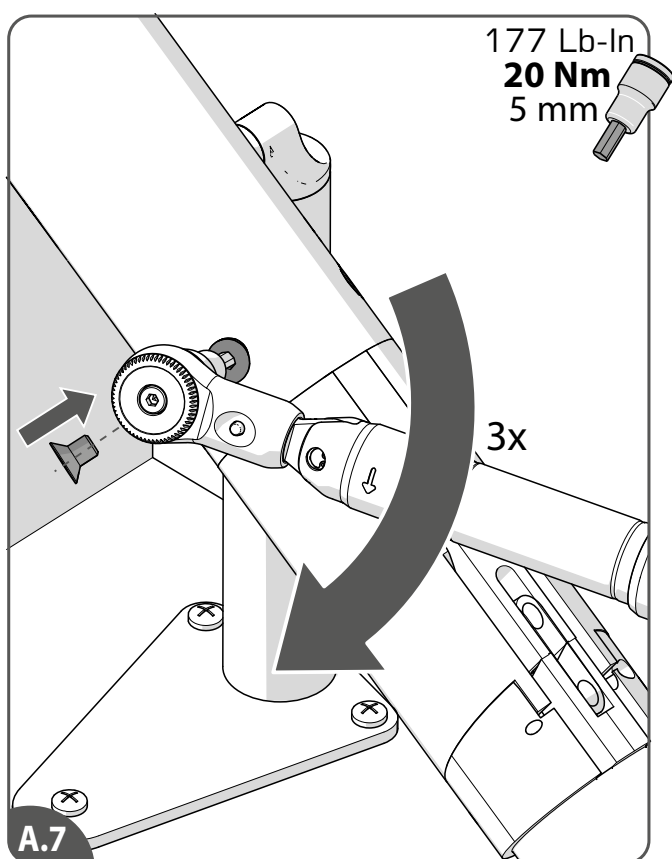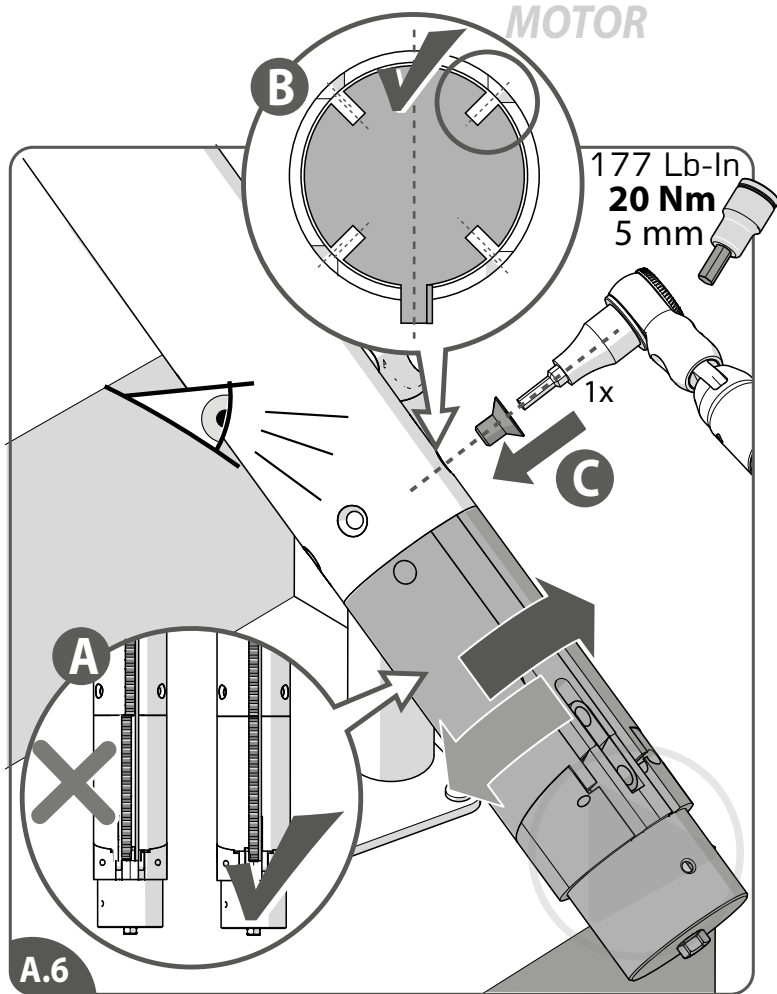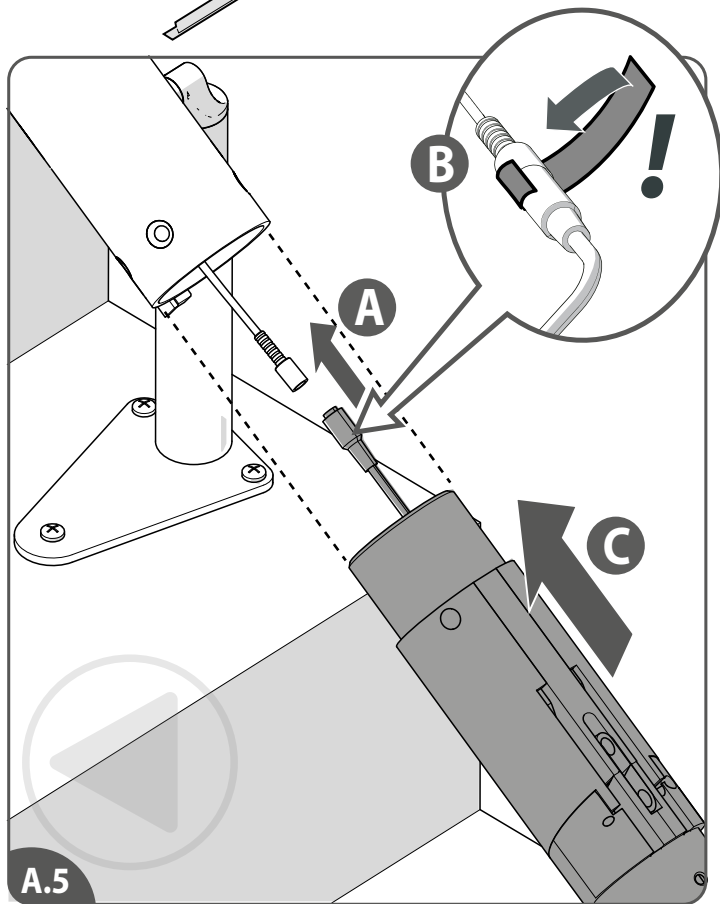- [A5] Connect the hinged rail motor cable{A}. Secure the connection with tape{B}. Install the hinge motor{C}.
 [A6] Align the toothrack of the motor with the toothrack of the rail{A} and install the first motor bolt{B}. Fasten this bolt [20 Nm | 177 Lb-Inch]{C}.
 [A7] Install the other motor bolts. Fasten these bolts [20 Nm | 177 Lb-Inch].
 [A8] Install the lower hinged rail part. Align the toothracks.

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IT

[A9]Align the toothrack of the motor with the toothrack of the rail{A} and install the first motor bolt{B}. Fasten this bolt [20 Nm | 177 Lb-Inch]{C}.

[A10]Install the other motor bolts. Fasten these bolts [20 Nm | 177 Lb-Inch].

[A11]If necessary, shorten the domex strip a little for a proper connection of the two domex strip parts.

[A12]Install the endcover and fasten the two screws of the endcover. In case of plush carpet change the position of the endcover{B}.

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





IT

B A-RAILS [ADDITIONAL] A1 vs A2

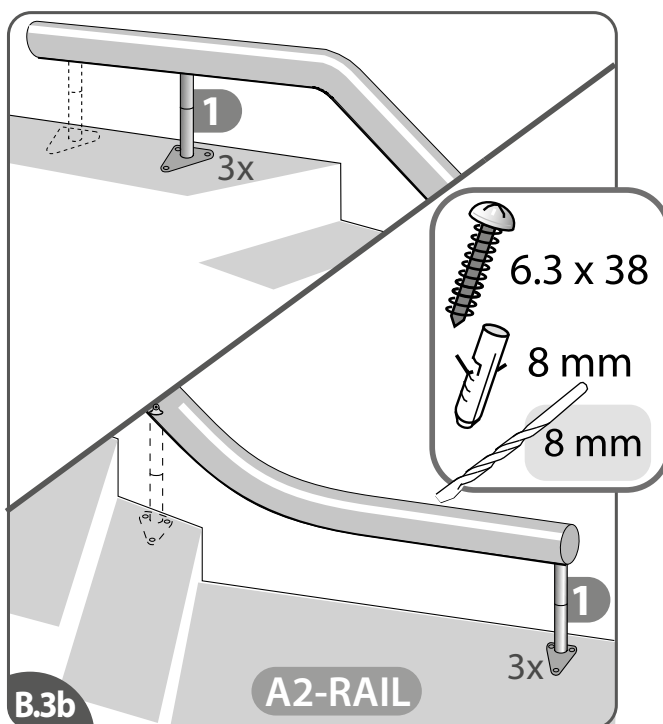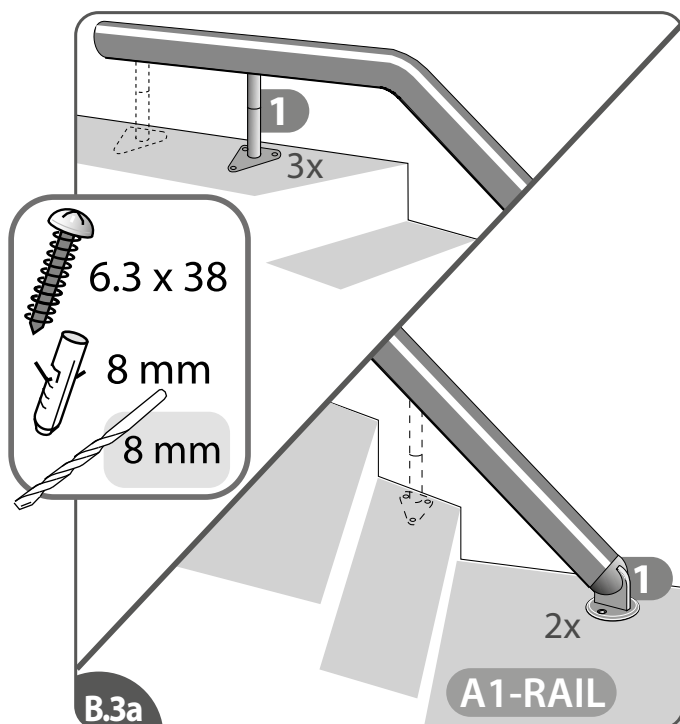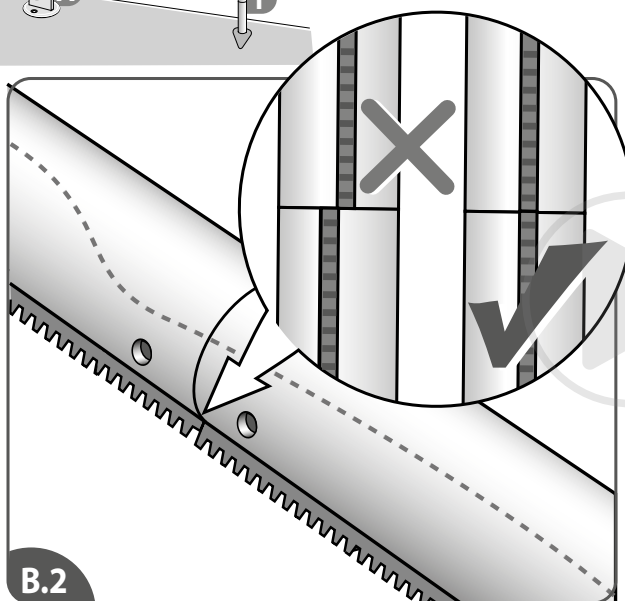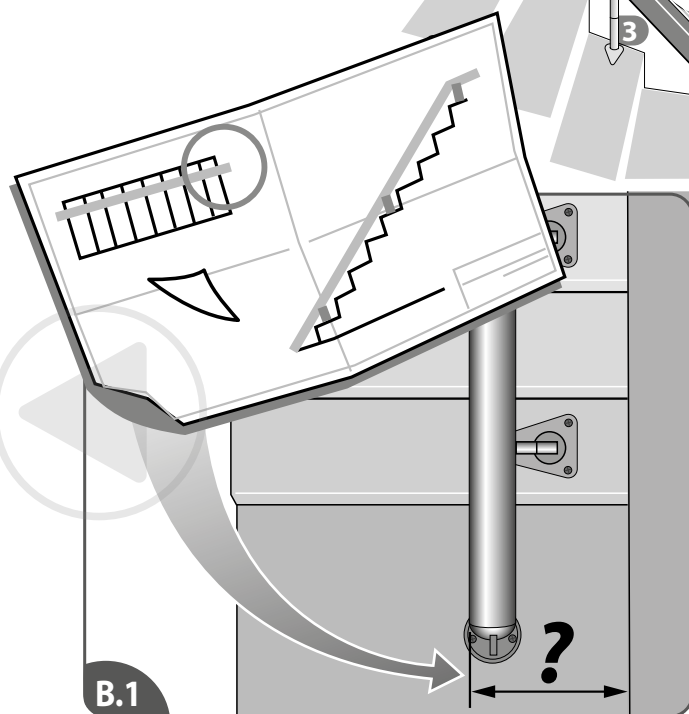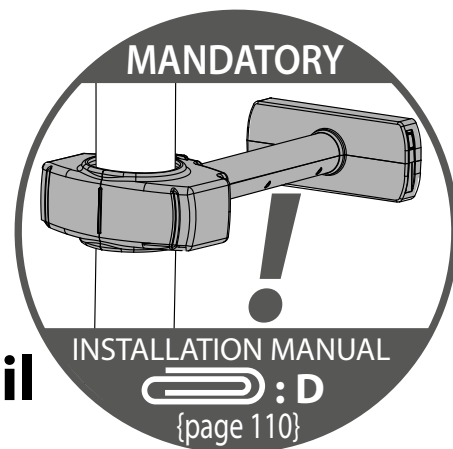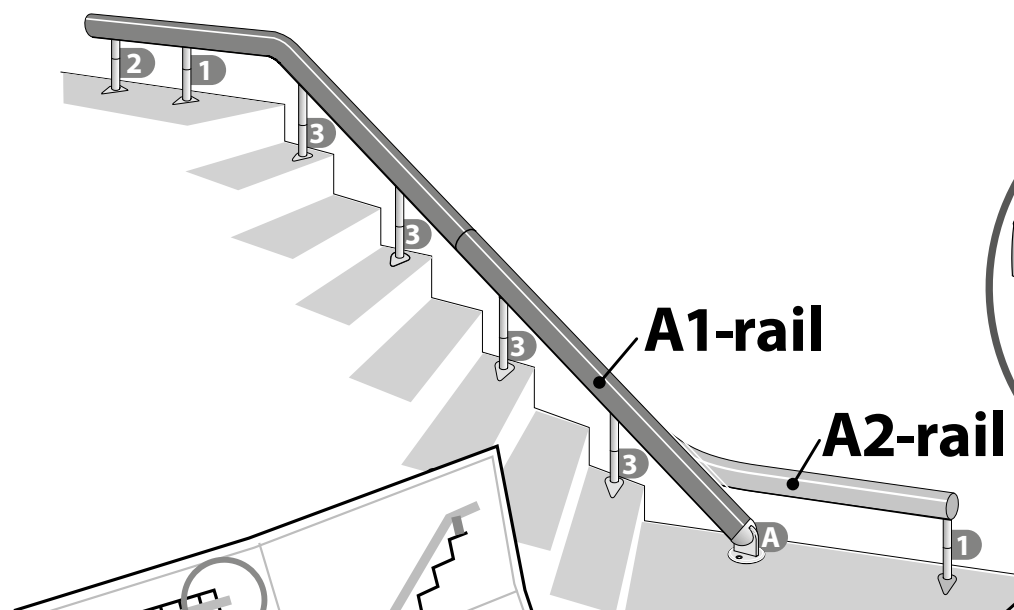Note: wall brackets are mandatory on **A**-rails.

[B1] Position the rail on the stairs. The distance to the wall is mentioned on the enclosed drawing.

[B2] After installing of the joint according to this manual, do an extra check on the alignment of domexstrip / toothrail.

[B3a] Fix two stanchions{1} on the bottom/topfloor according to this manual.

[B3b] In case of rail **A2** you start with the lower stanchion at the bottomfloor and the first stanchion at the top landing.

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




IT

B A-RAILS [ADDITIONAL] PERPENDICULAR

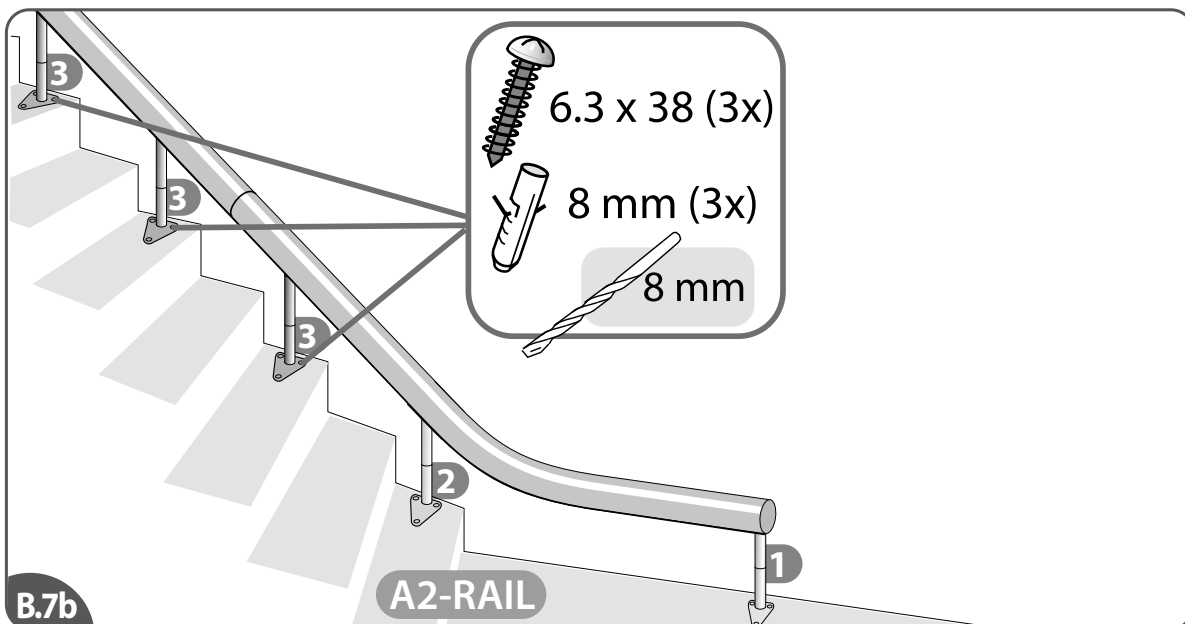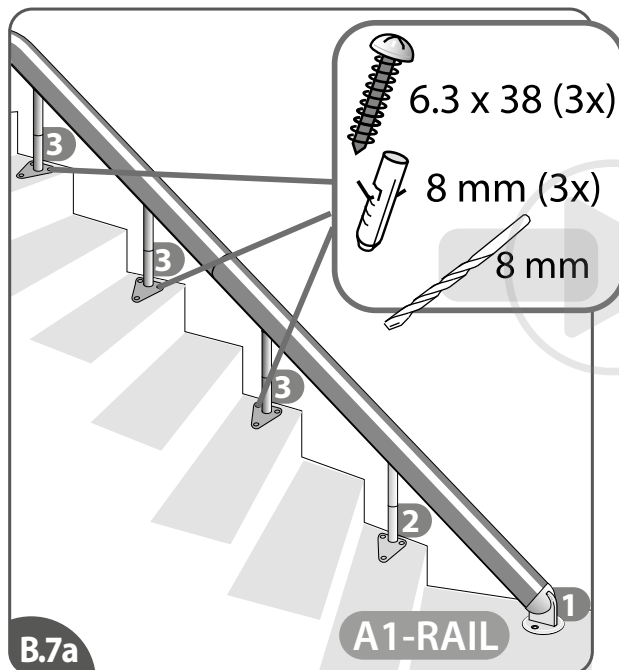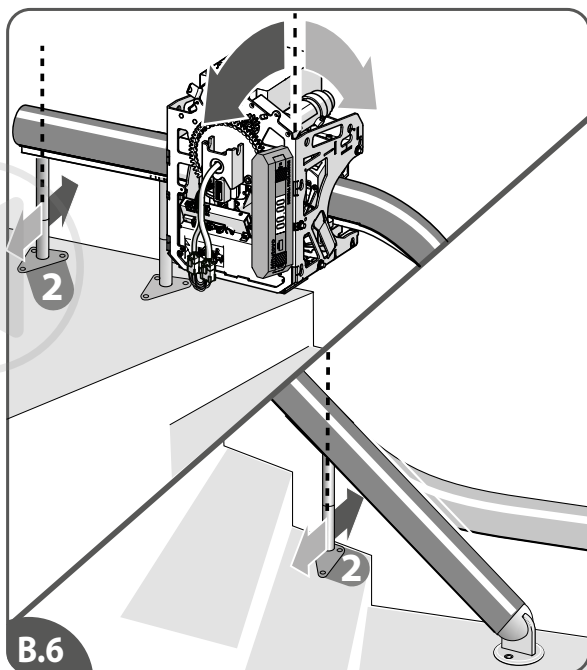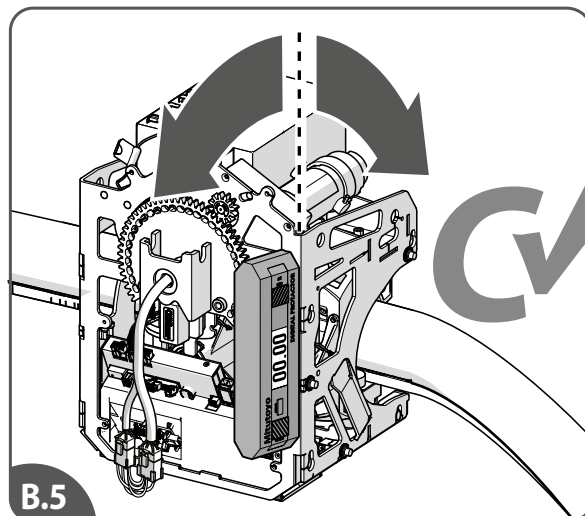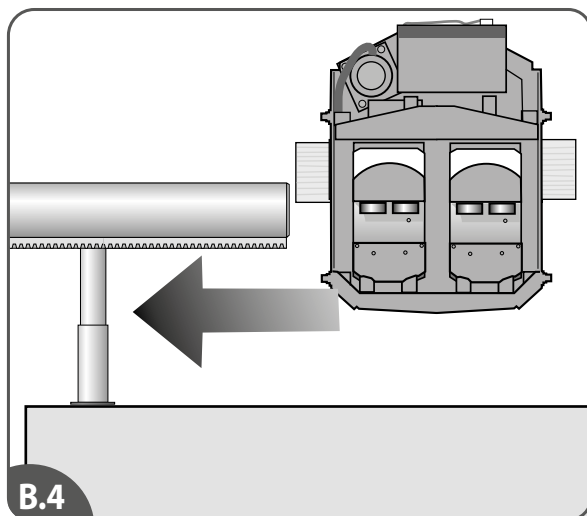- [B4] Ride the unit on the rail and position it in the curve of the rail.
 [B5] Check with a spirit level if the unit is exactly vertical.
 [B6] If the unit leans forwards or backwards, the rail shall be adjusted by repositioning of the upper and lower stanchion{2}. This will till the rail forwards or backwards. In case of rail **A1** the stanchions on the staircase first.
 [B7a] Fix the other stanchions to the stairs{3}.
 [B7b] In case of rail **A2** also fix the other stanchions{3}.

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



CHANGING CHANNELS

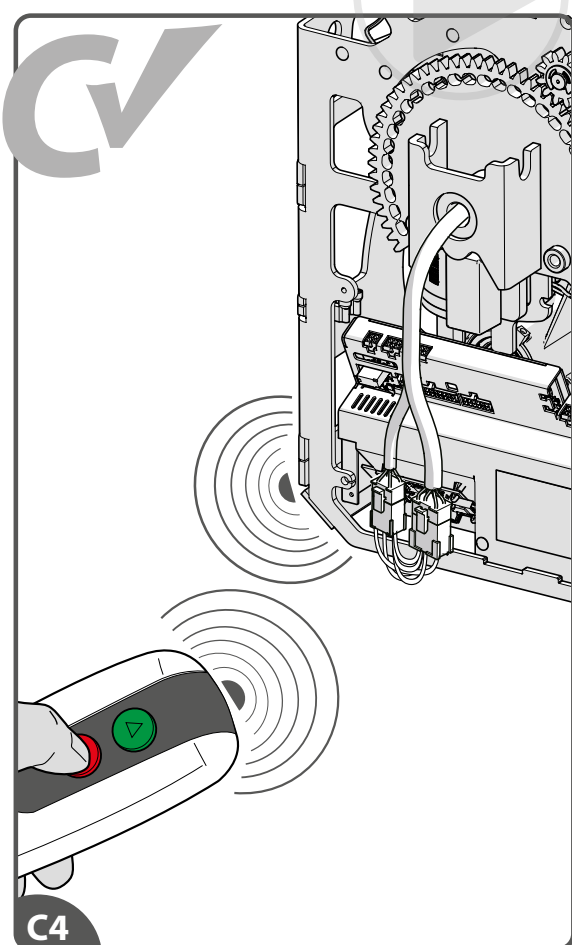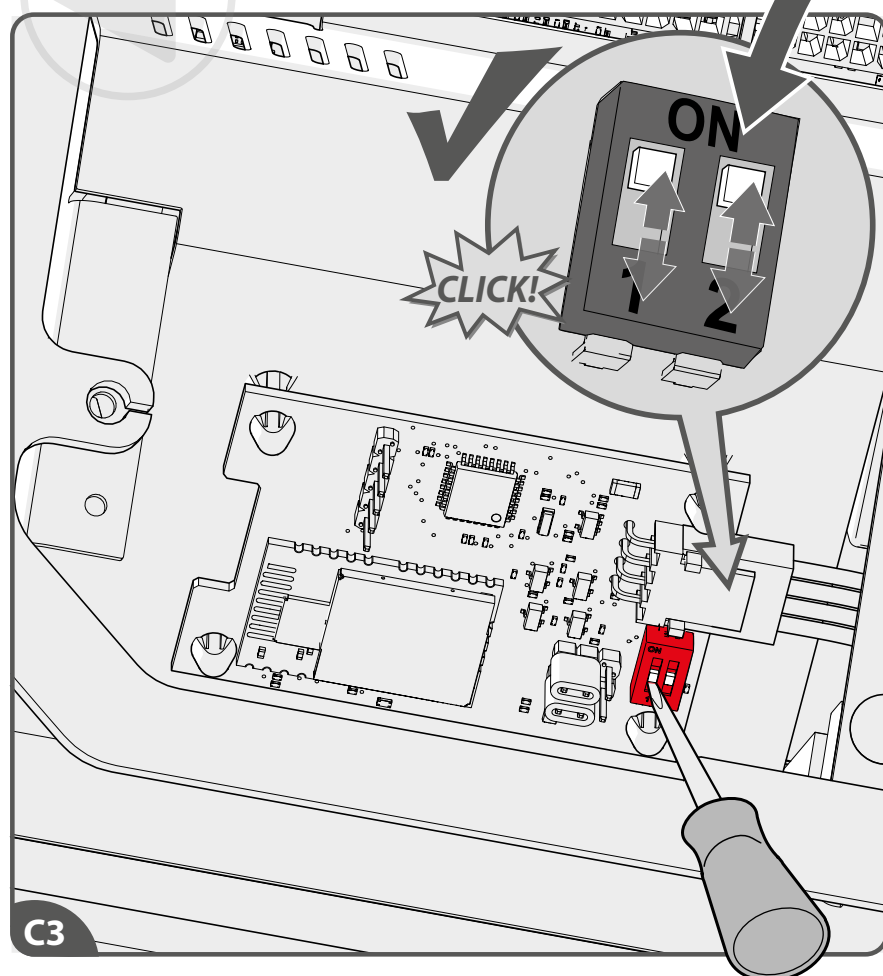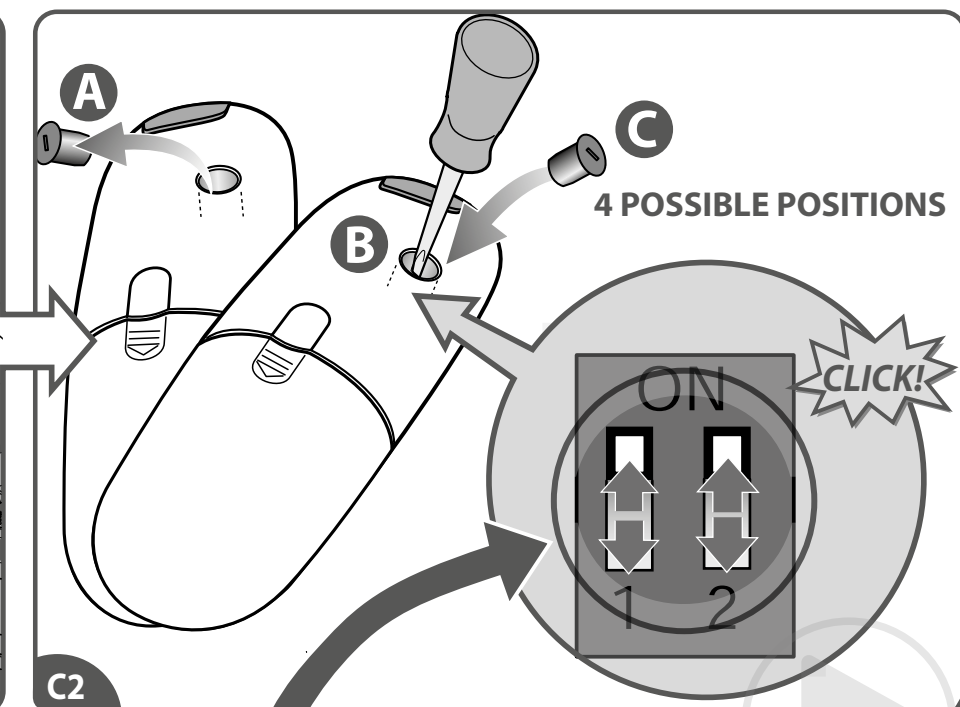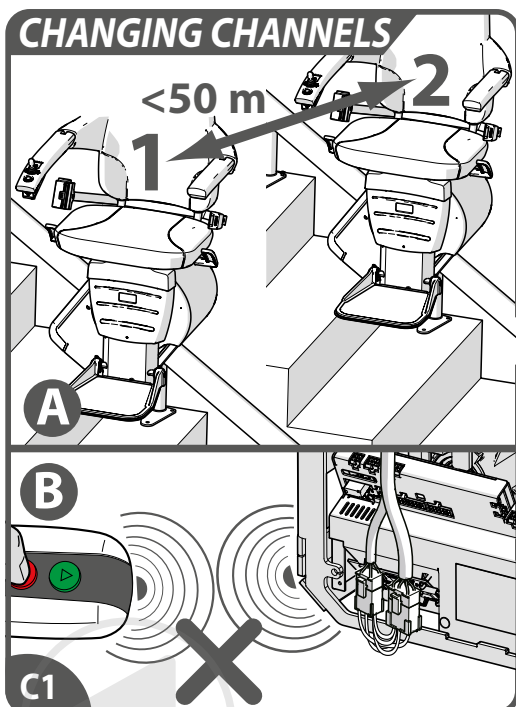Changing channels

[C1] If there is more than one lift in the range of the remotes {A} or otherwise interference {B}, you can set the remotes on a different channel. There are four channels.

[C2] Remove the caps on the backside of the remotes {A}. Set the dipswitches {B} (**note:** pay attention to the position of the word “**on**” as this is not always on the same position). Reinstall the caps {C}.

[C3] Set the dipswitch on the board.

[C4] Check if the remote controles are still connected. If not, repeat the normal linking procedure on page 48.

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




IT

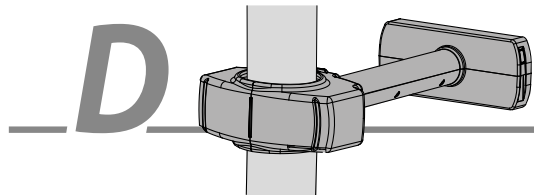WALL SUPPORT FIXATION

SEE SALESMANUAL (AA28974) IN WHICH SITUATIONS WALL SUPPORTS ARE MANDATORY

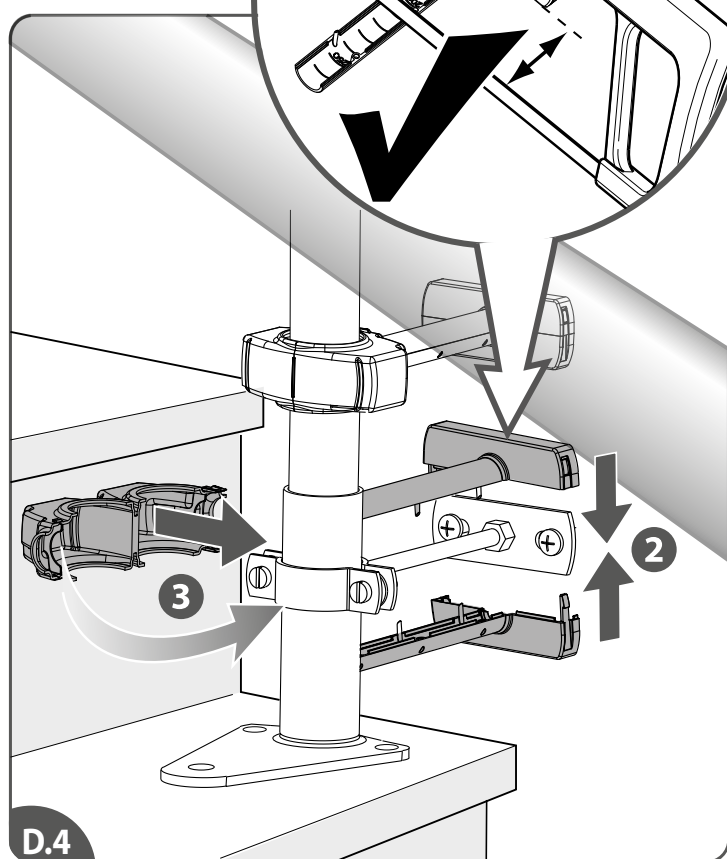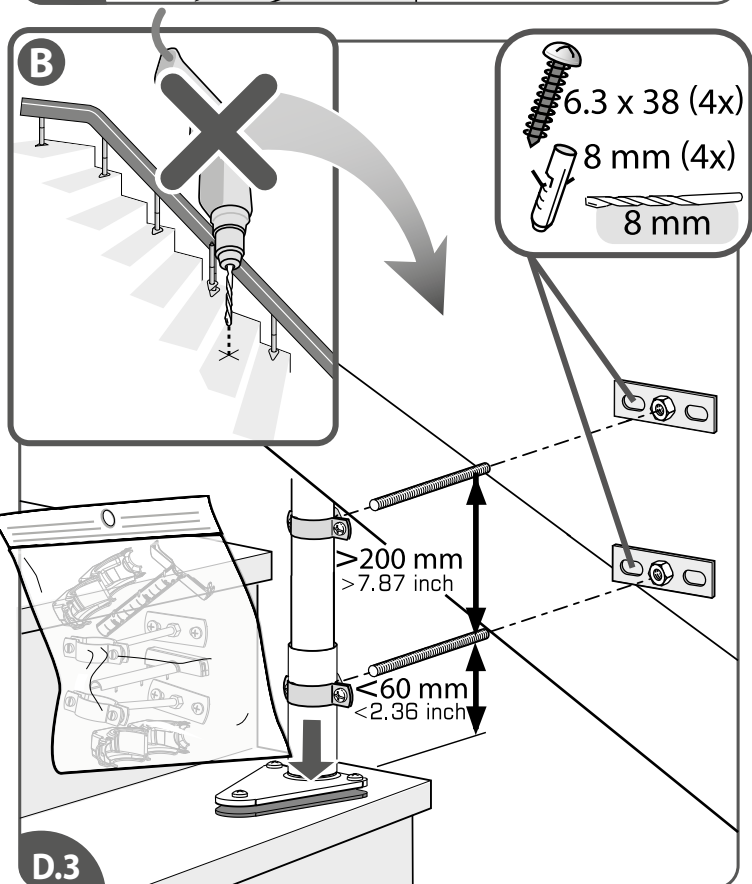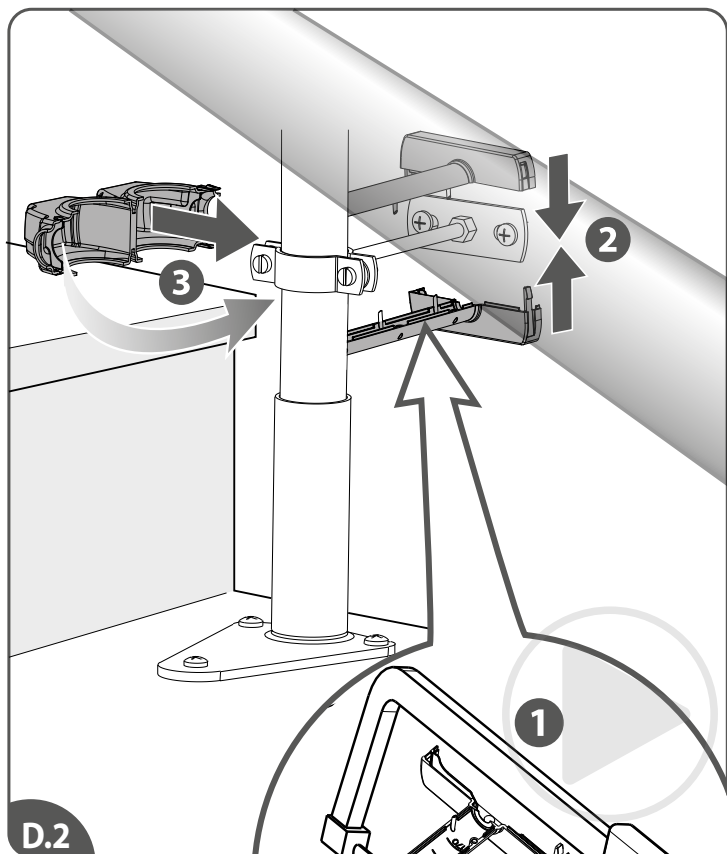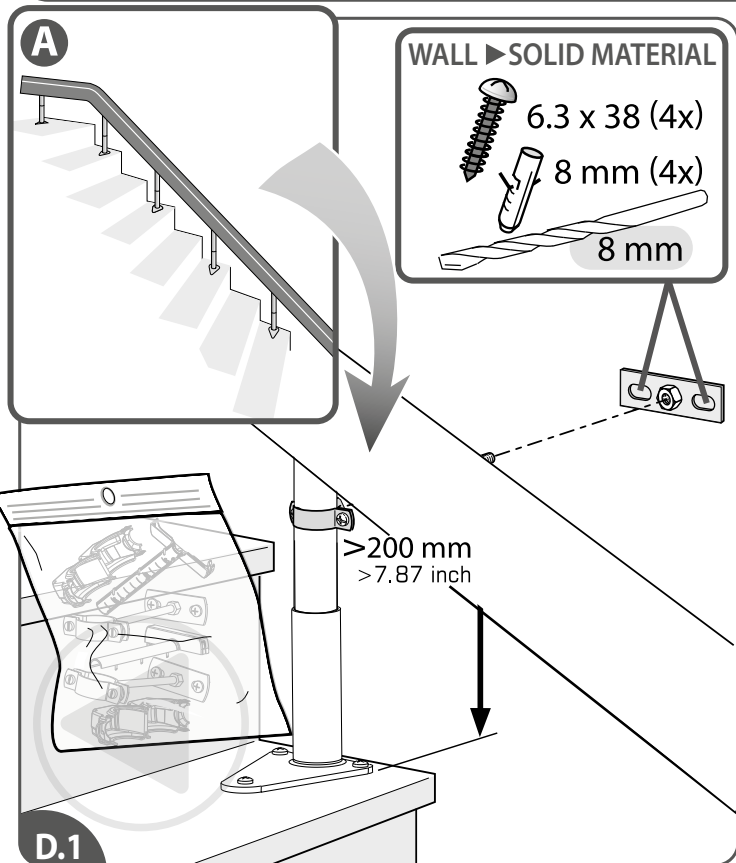Note: see salesmanual (AA28974) in which situations wall supports are mandatory.

When using wall brackets the wall shall be of solid material e.g. not a gypse/dry wall. The wall shall be minimum half brick. Bricks shall not be hollow.

[D1]In case it is allowed to fix the stanchion with screws to the stairs, just install one wall bracket. Install the bracket at least 200 mm [7.87 inch] above the stairs/floor. Bring the threaded rod in position and fix the wallplate with screws.

[D2]Cut the bracket covers to the right length{1}. Close the covers around the rod{2}. Install the front bracket cover{3}.

[D3]In case it is not allowed to fix the stanchion with screws to the stairs, install two wall brackets. Install the bracket less than 60 mm [2.36 inch] above the stairs/floor.

EN [D4]Repeat the above described procedure to install the lower wall bracket.

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
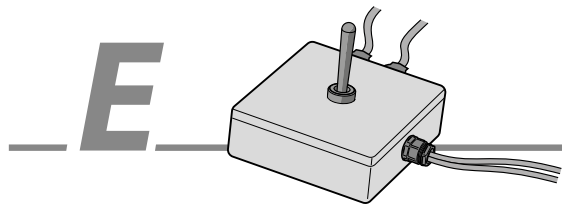
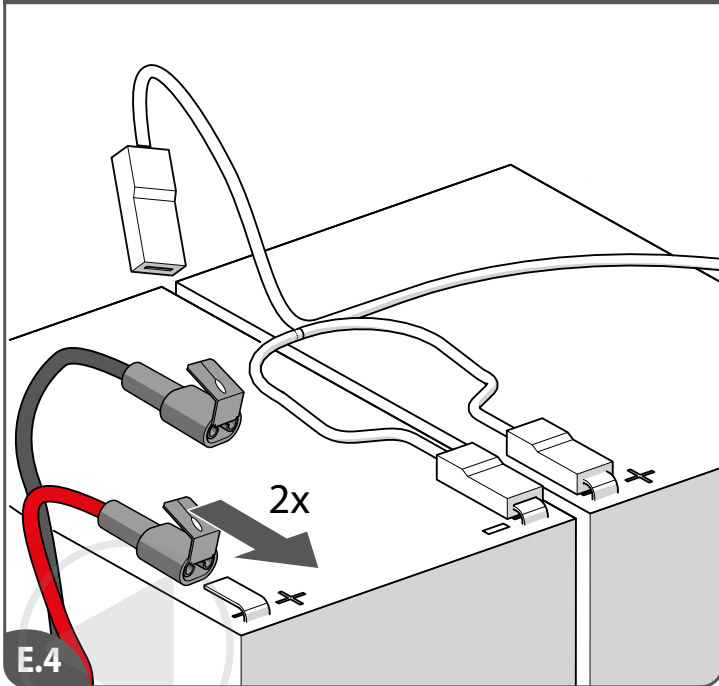
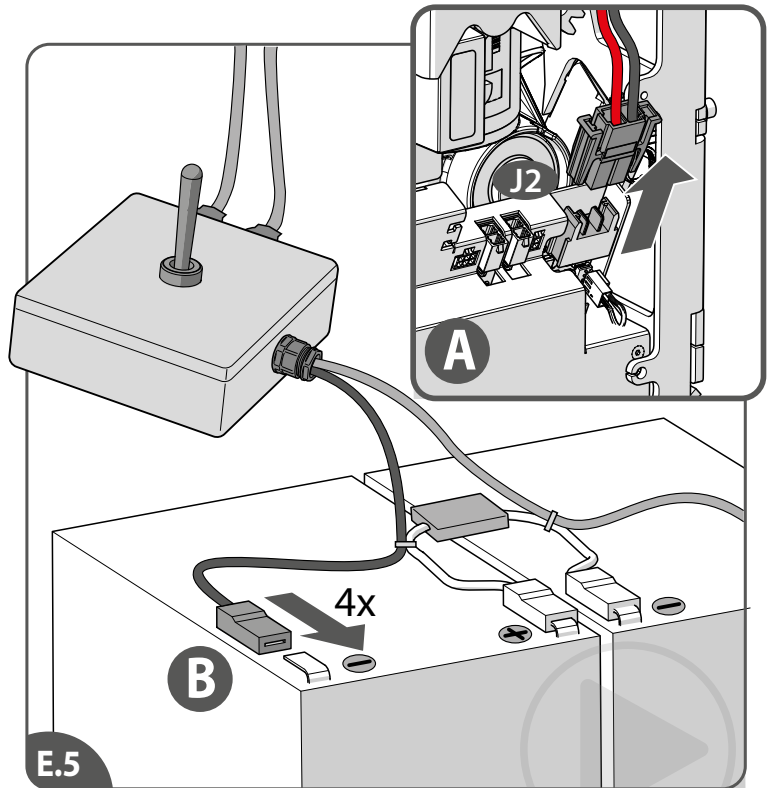
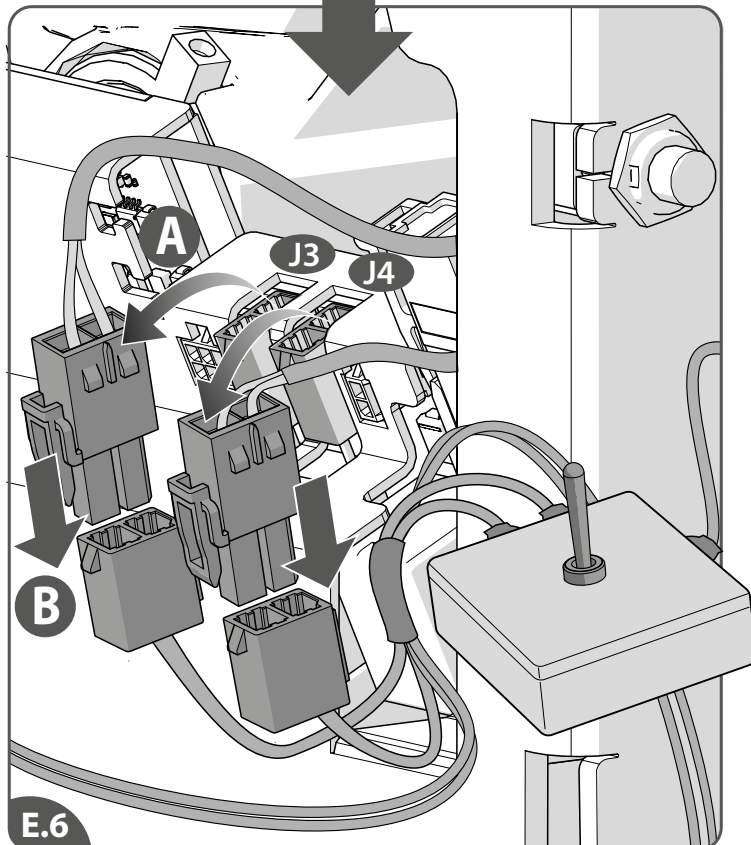
E.6

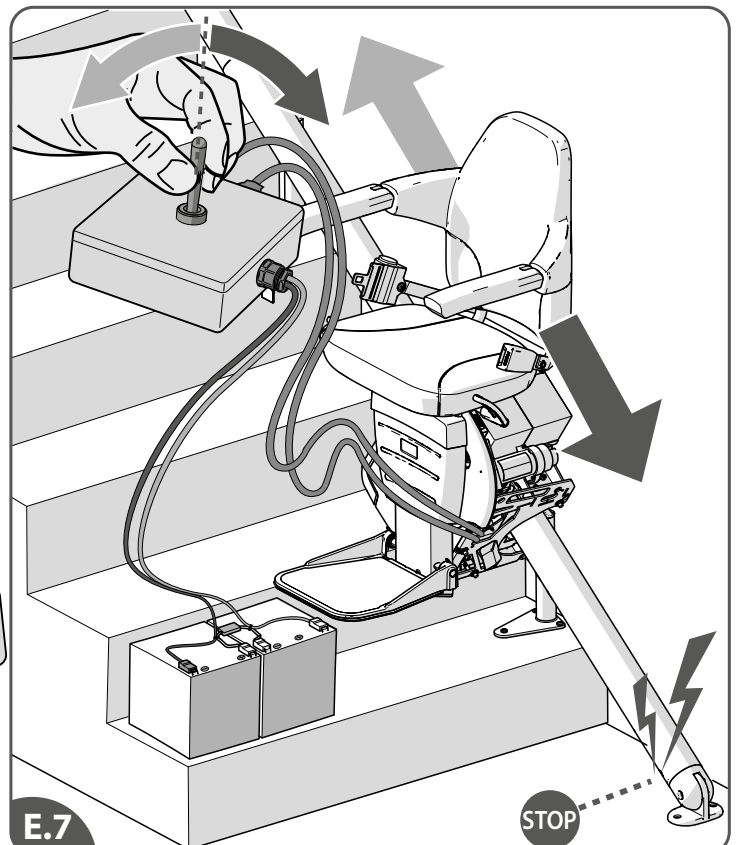E.7

- [E4]Only in case of no movement of the lift: connect the wires of the Tool unit traction/powered swivel with the extra batteries.
- [E5]Disconnect the battery loom{A}. Connect the Tool unit traction/powered swivel with the extra batteries{B}.
- [E6]Disconnect the motor looms{A}. Connect the motor looms with the tool{B}.
- [E7]Operate the joystick of the tool to move the lift to a safer loaction.

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

IT

E **EMERGENCY OPERATION** OPERATE MOVEMENT

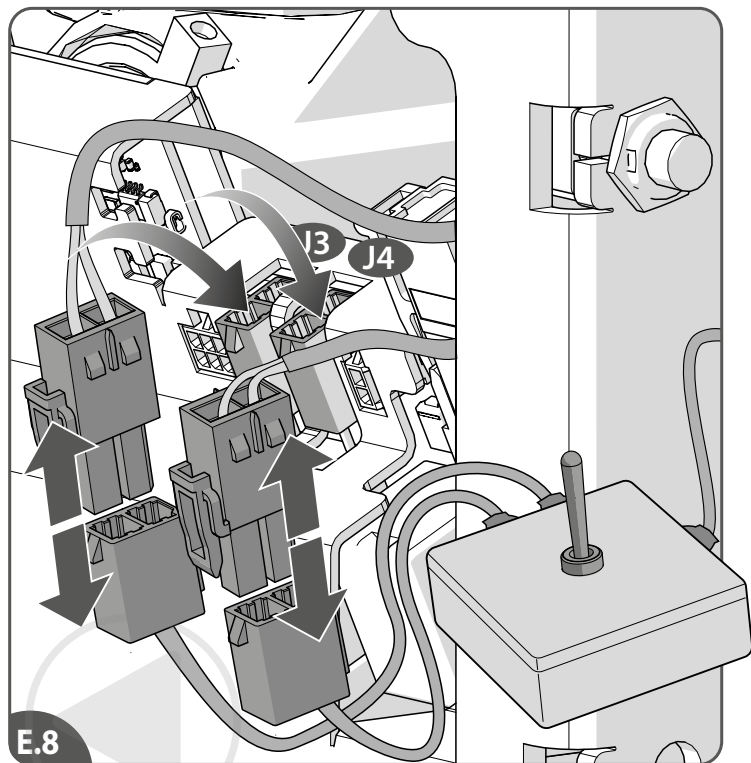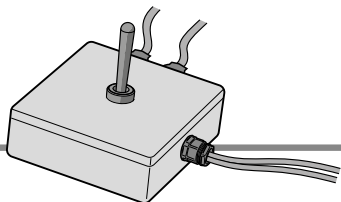REINSTALL LIFT > SEE MANUAL

[E8]Reconnect the motor looms the main PCB.

See this manual to reinstall the lift.

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




F FOOTREST ADJUSTMENT

UNSCREW BOLTS

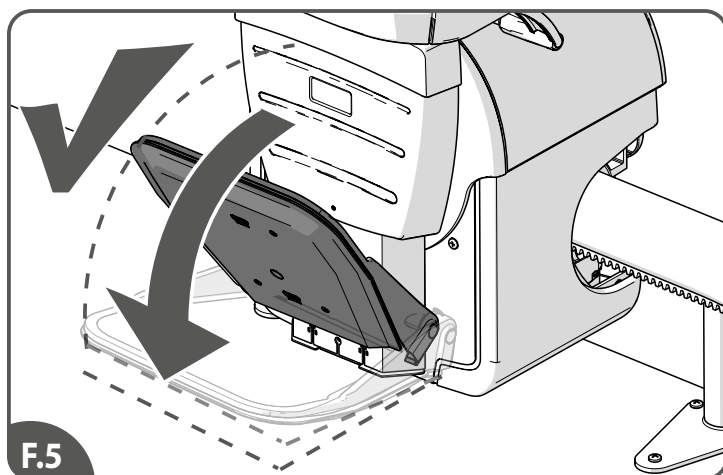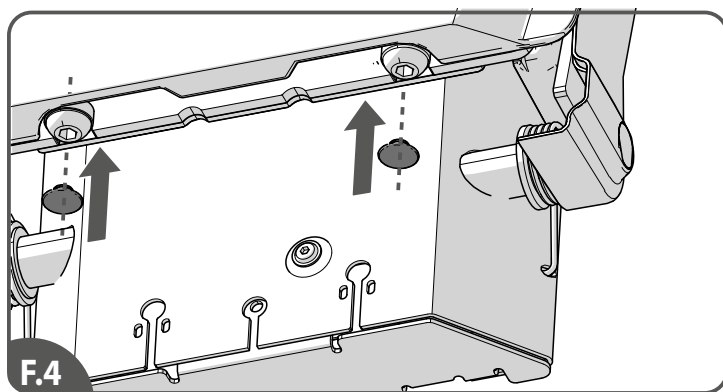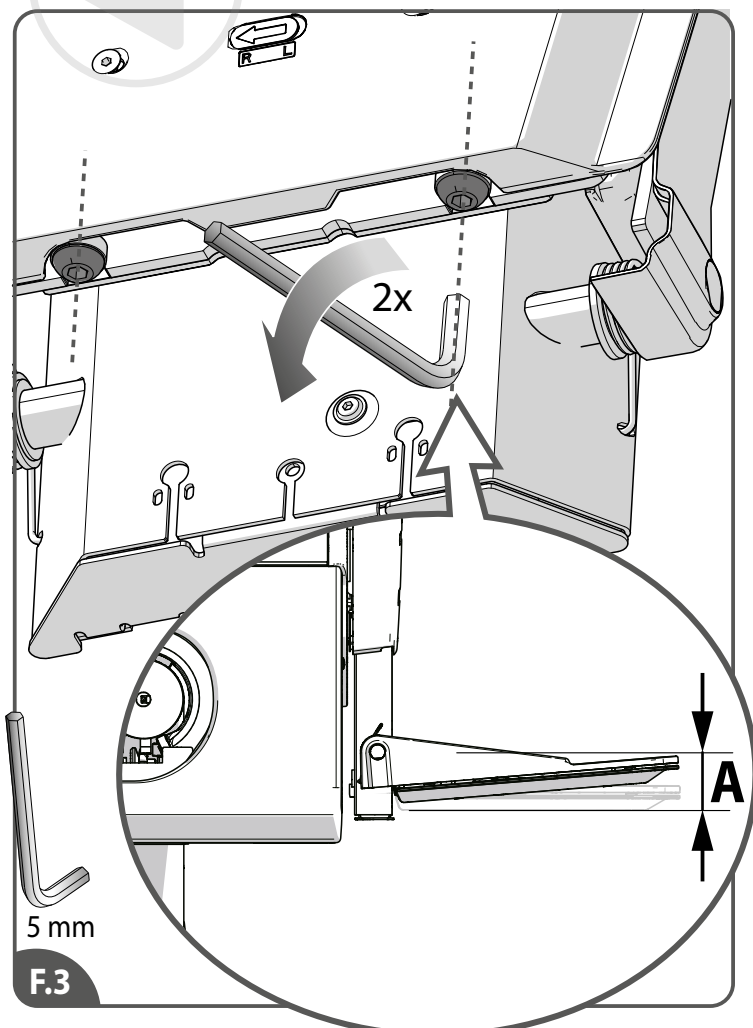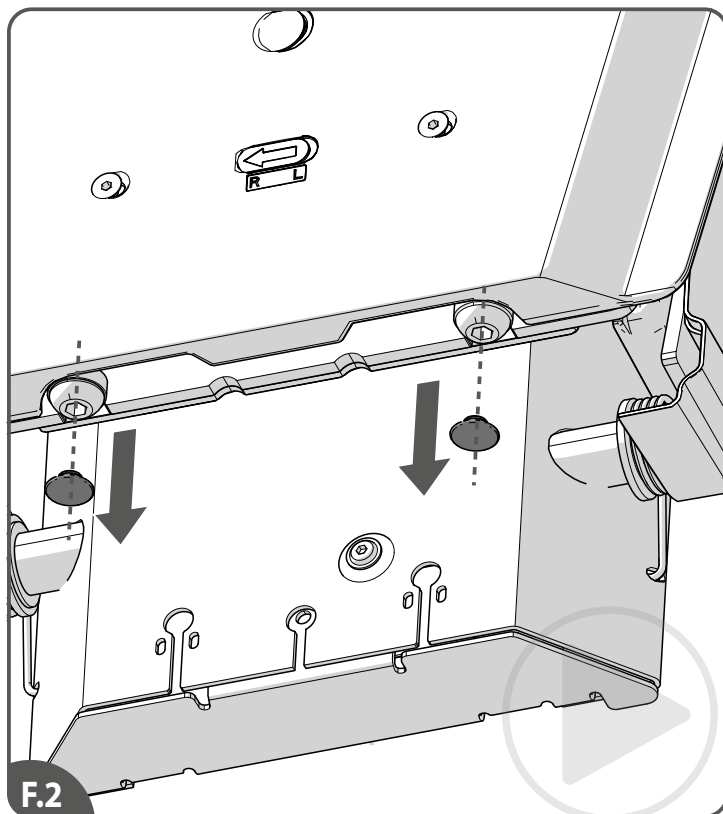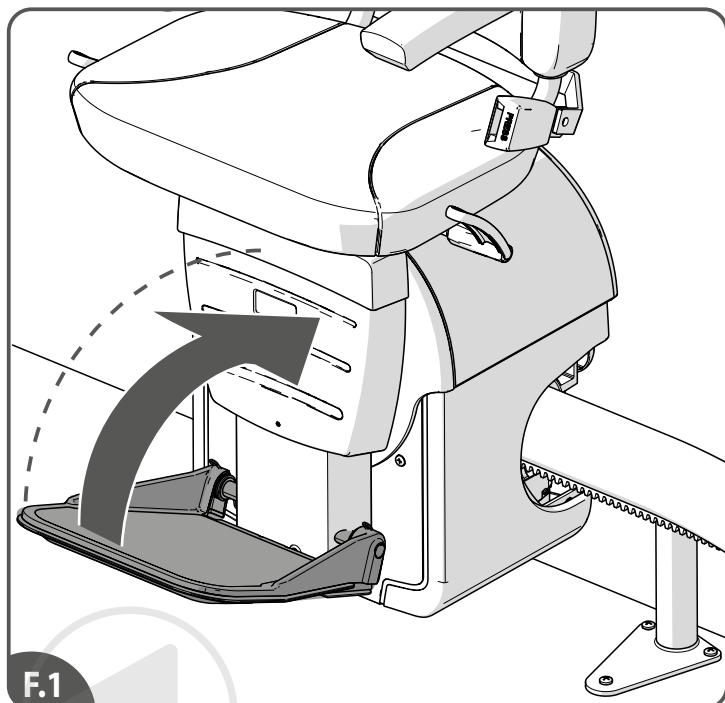- [F1]Fold up the footrest.
- [F2]Remove the cap (2x).
- [F3]Loosen the adjustment bolts (2x).
- [F4]Refit the cap (2x).
- [F5]Fold down the footrest and check the new position of the footrest.

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





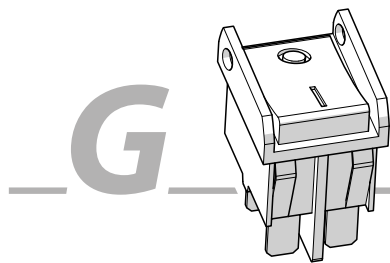
REMOVE COVERS

KIT: SP002000

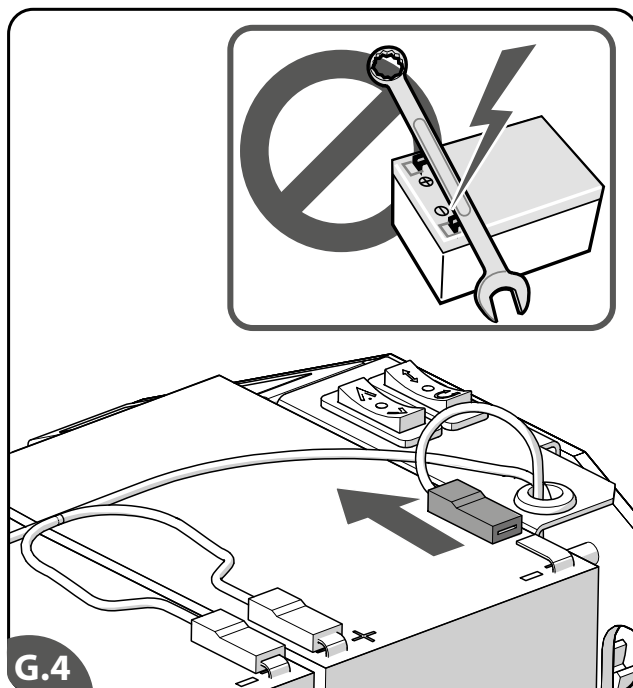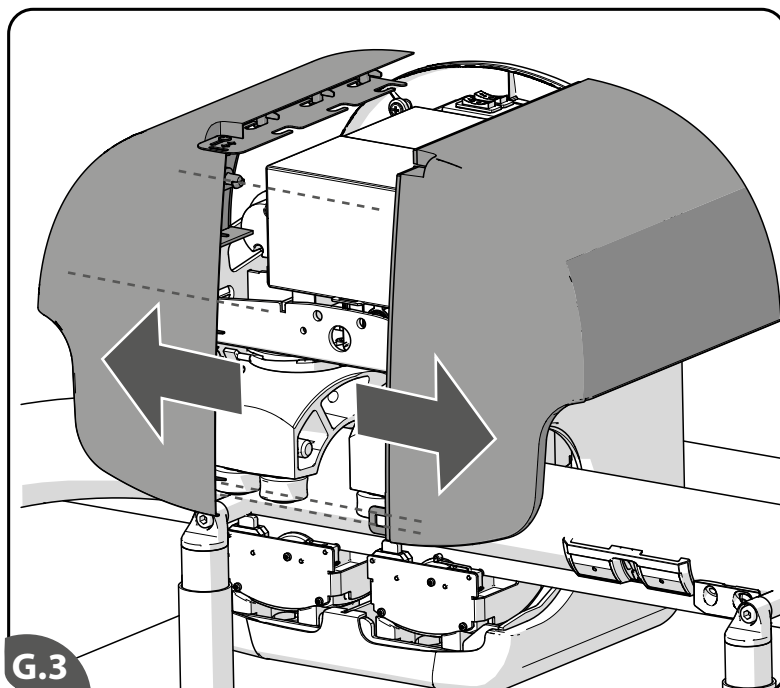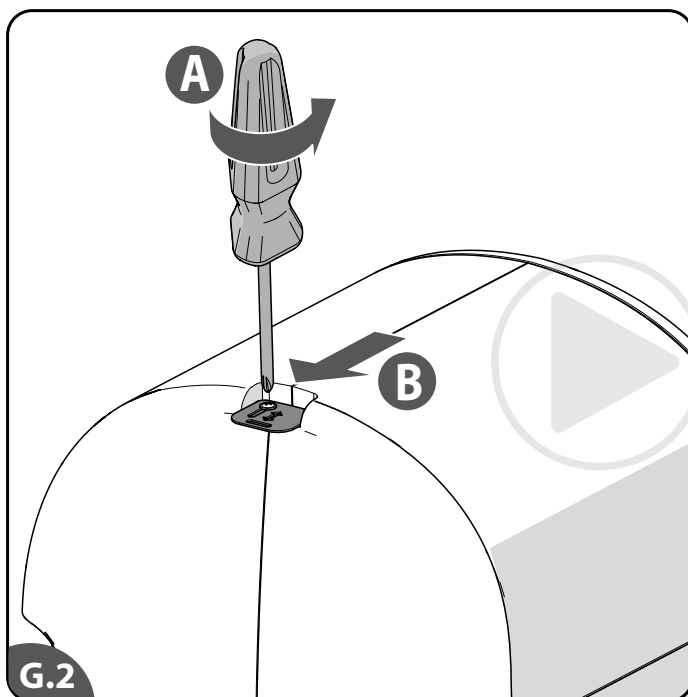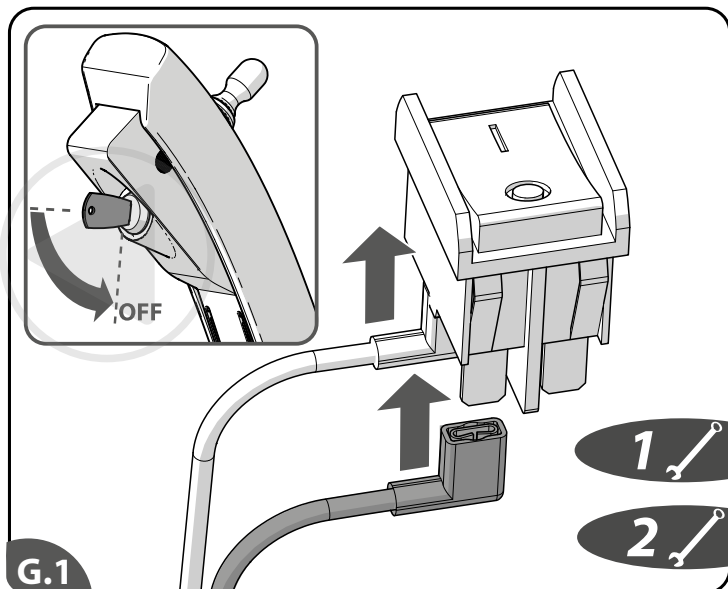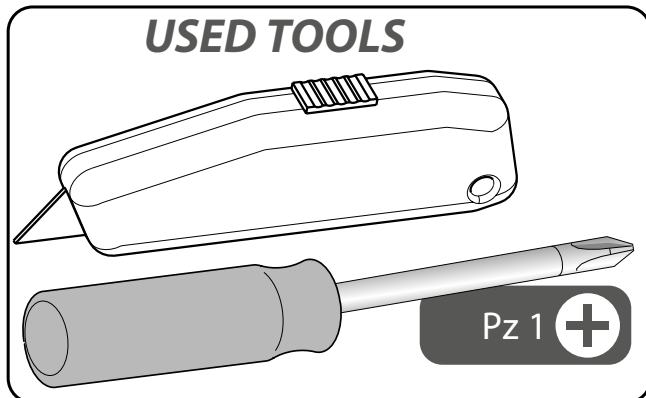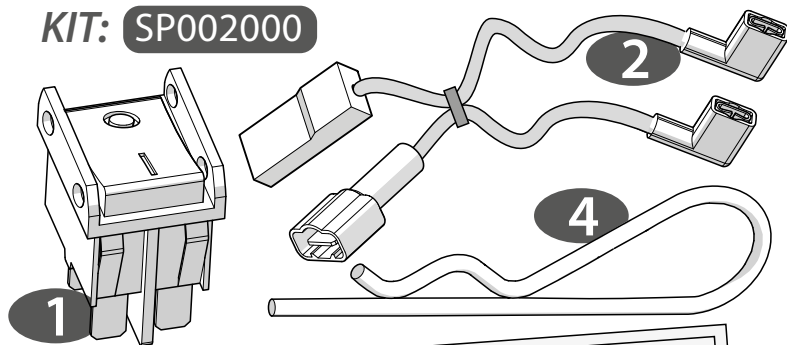- [G1] Connect the supplies wires with the on/off lockable switch.
- [G2] Loosen the top covers.
- [G3] Remove both covers.
- [G4] Disconnect the - wire of the battery.

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



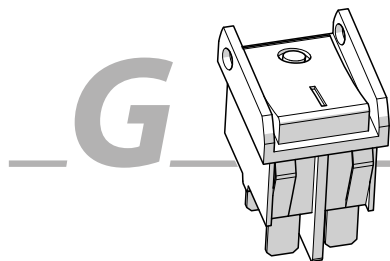
LOCKABLE SWITCH

CUT A HOLE

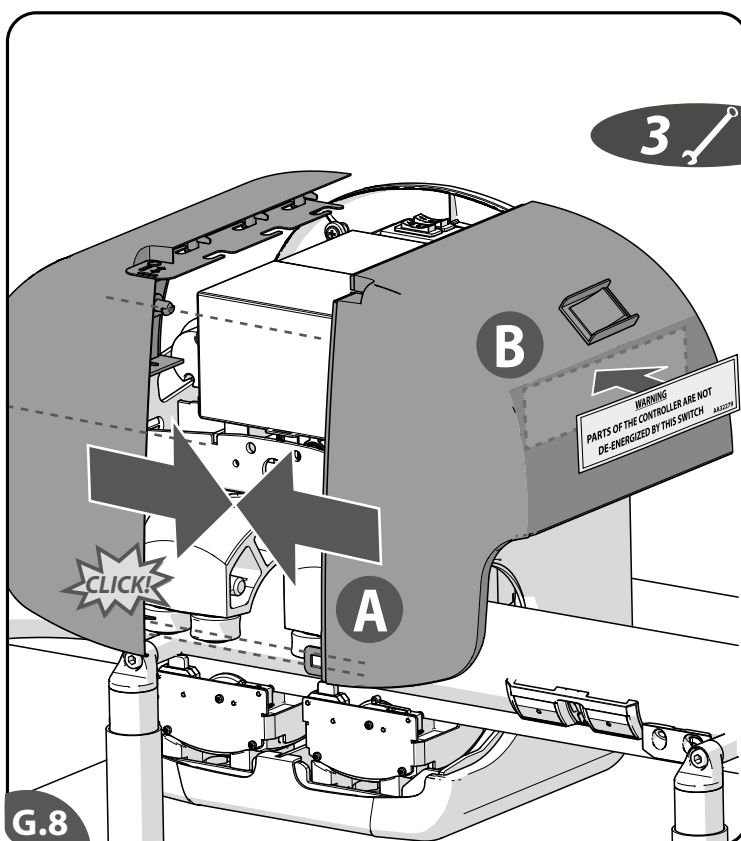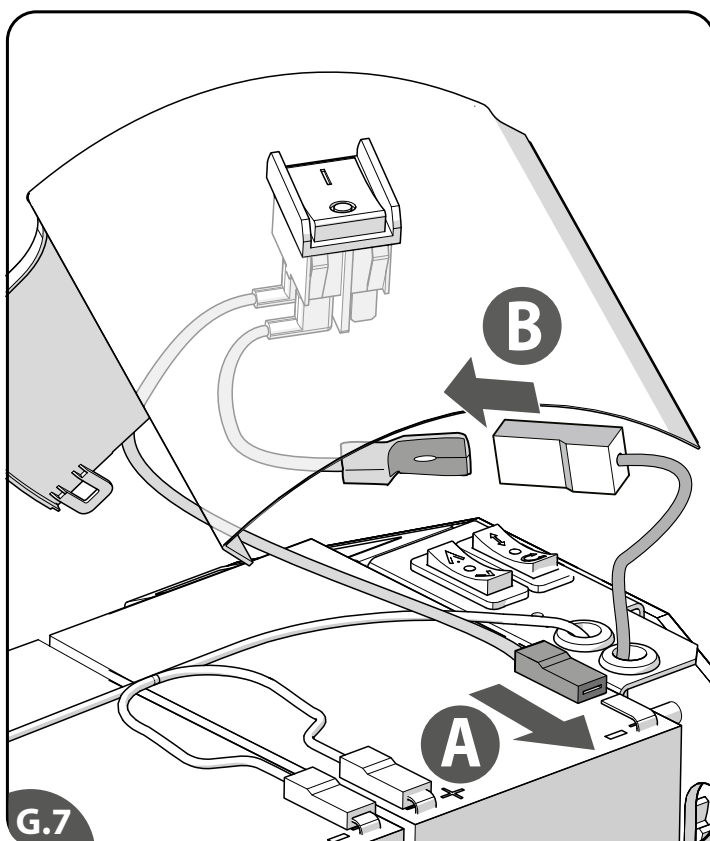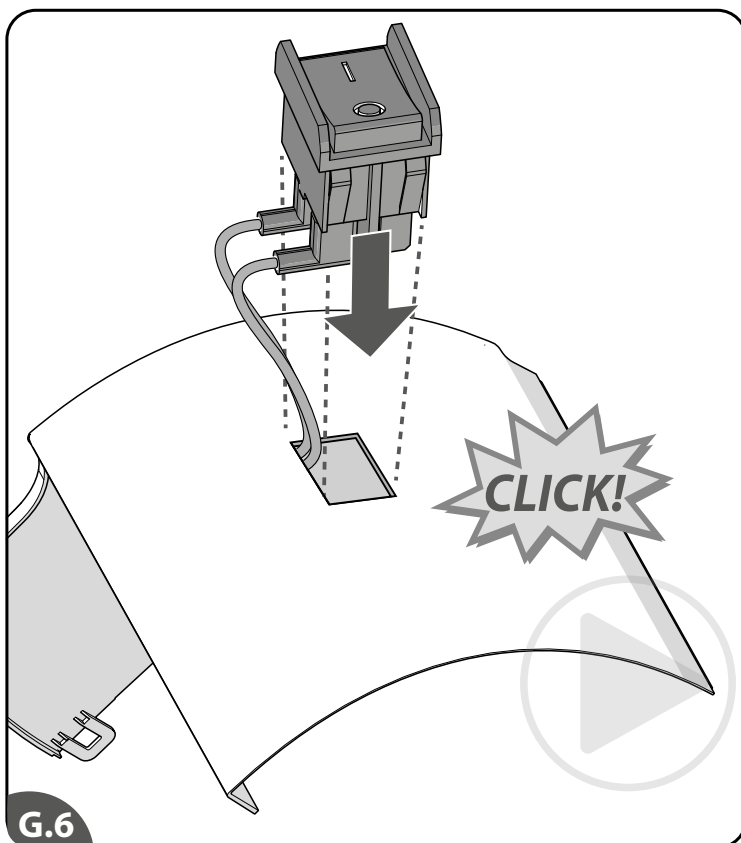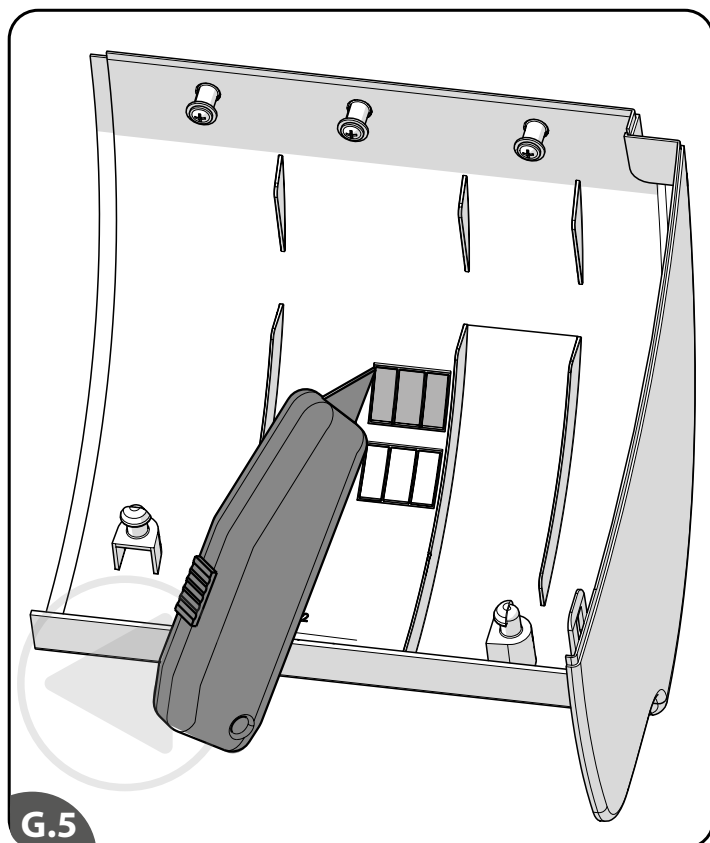- [G5]Cut, at the indicated location, a hole in the left top cover.
- [G6]Click the on/off lockable switch into the hole.
- [G7]Connect the on/off lockable switch with the battery{A}. Connect the on/off lockable switch with the “old” - wire{B}.
- [G8]Replace the covers{A}. Fit the warning label near the switch{C}.

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





IT

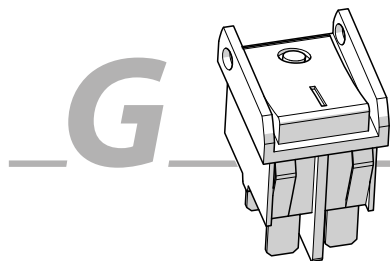LOCKABLE SWITCH TESTING

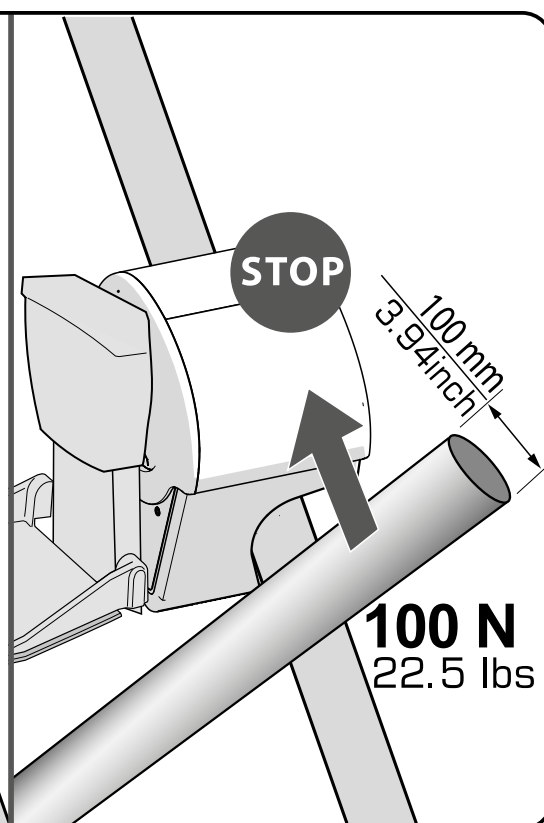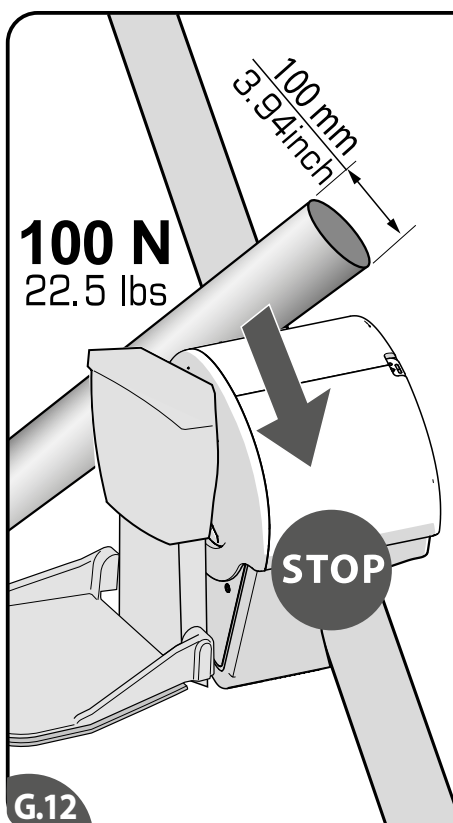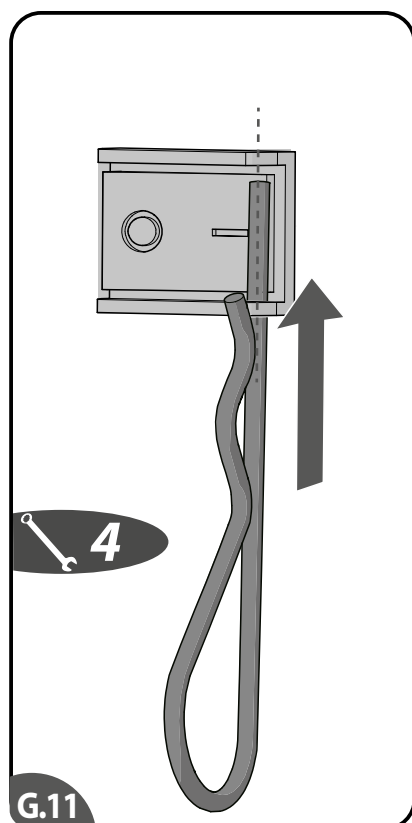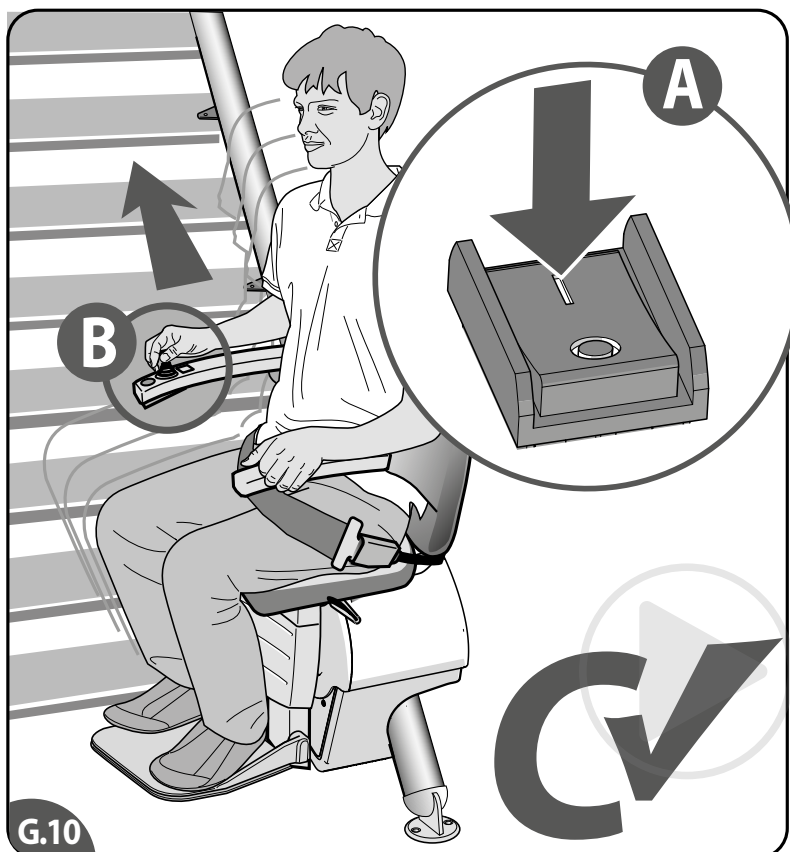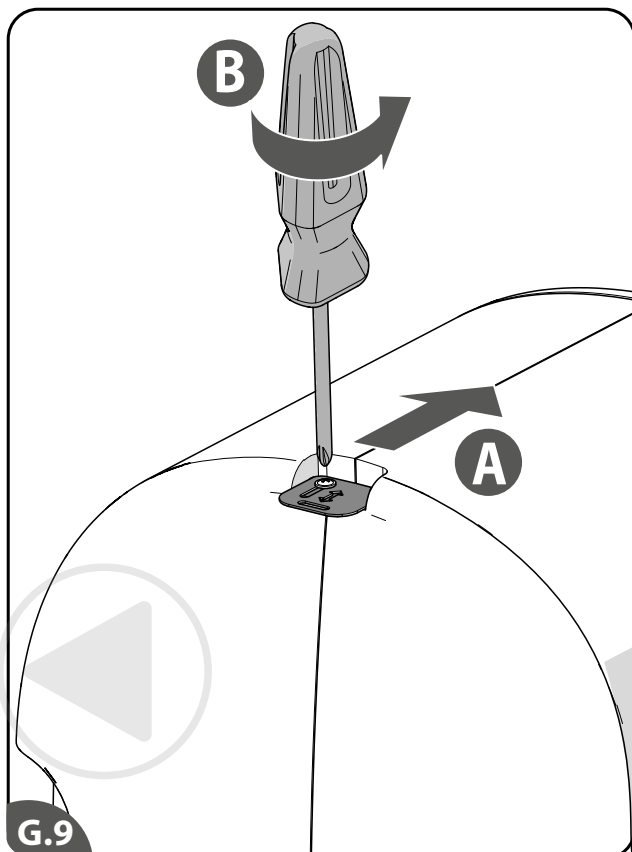- [G9] Refasten the covers.
- [G10] Switch on the on/off lockable switch{A} and check if the lift{B} is still good functioning.
- [G11] Install the locking pin.
- [G12] Check if the detection covers are still good functioning.

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





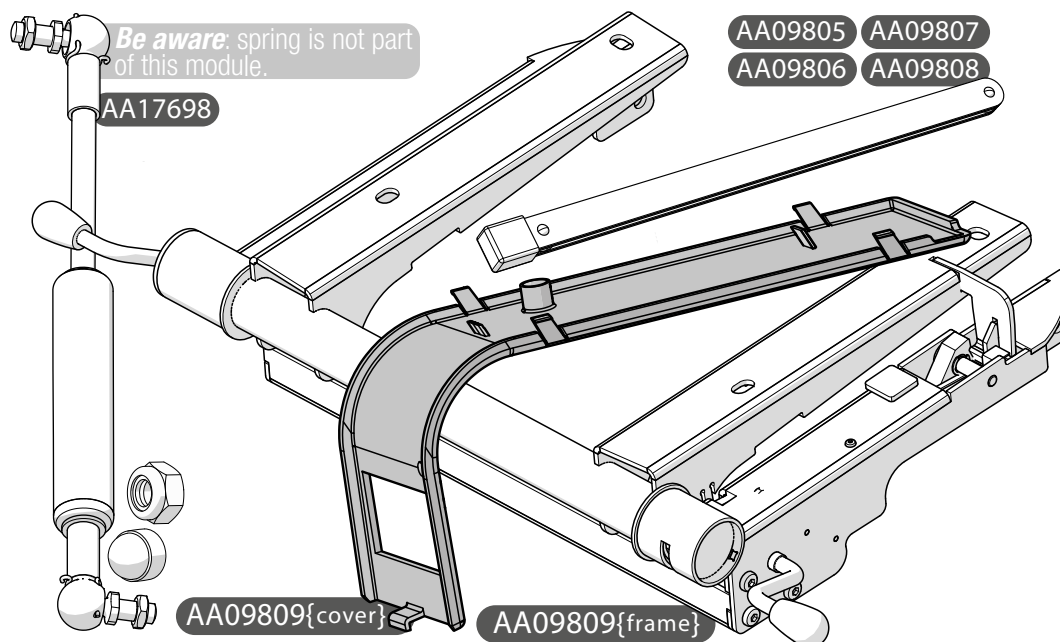
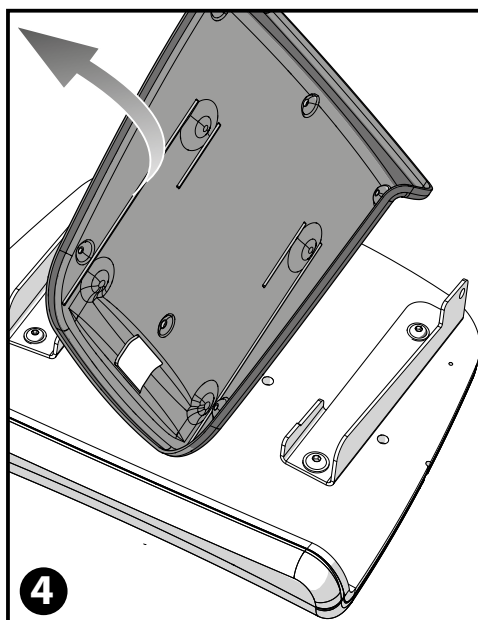
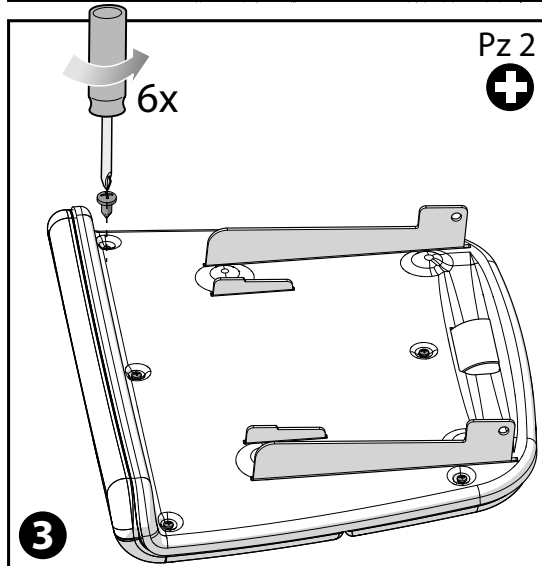
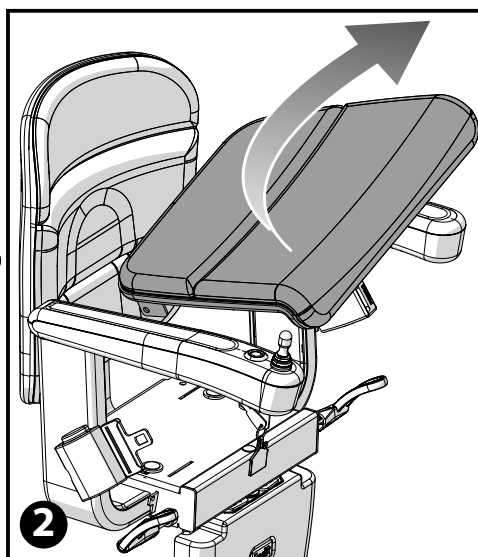
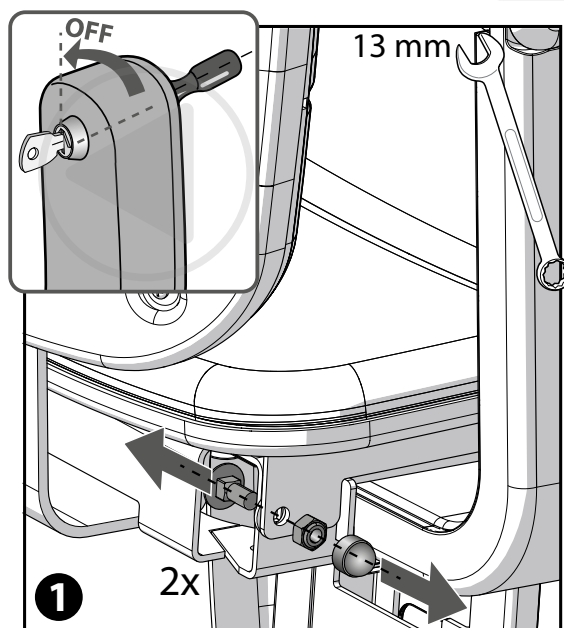
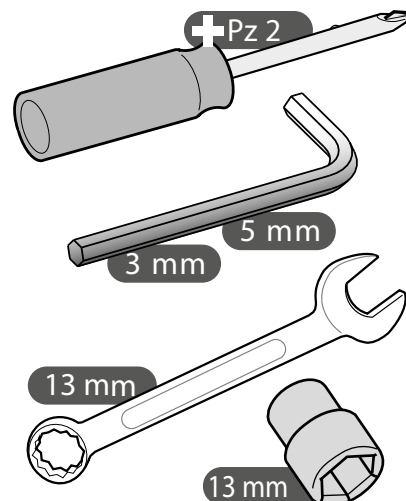
EN Check whether the lift-and-tilt module is complete and the appropriate tools are available. [1]Remove the bolts and nuts from the seat's swivel axle. [2]Remove the seat. [3]Loosen the screws in the bottom cover of the seat. [4]Remove the bottom cover.

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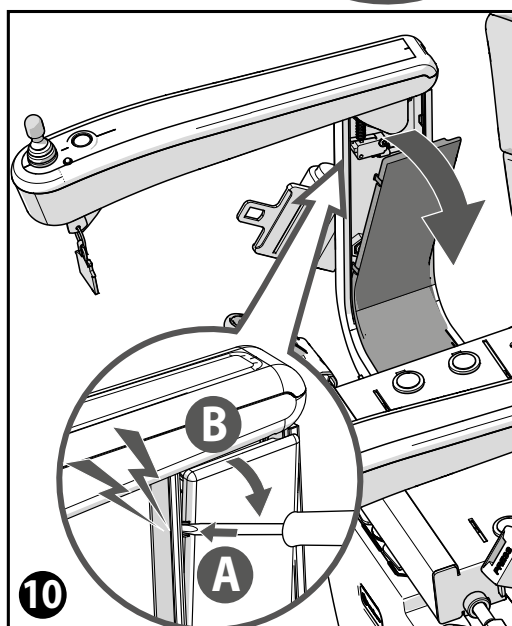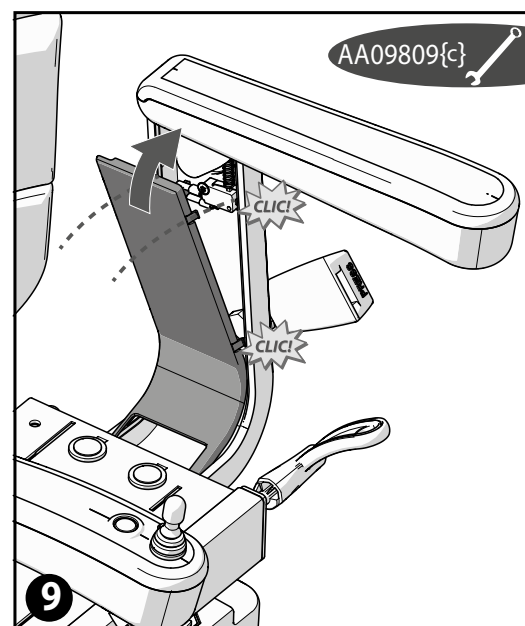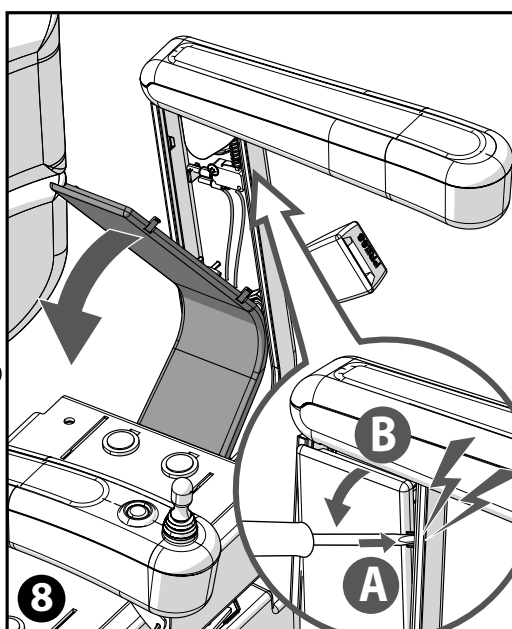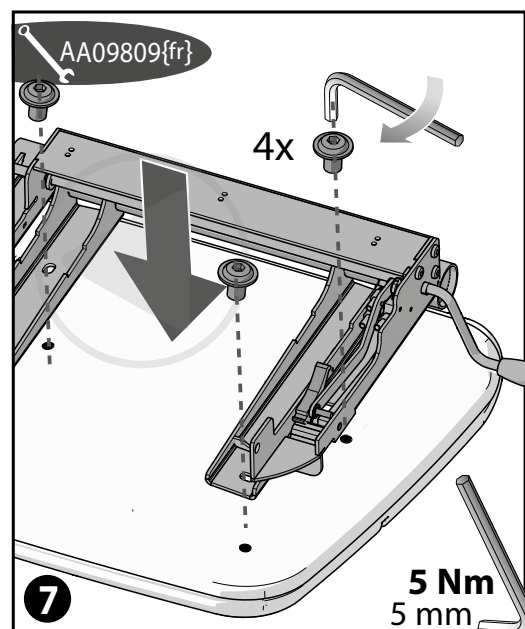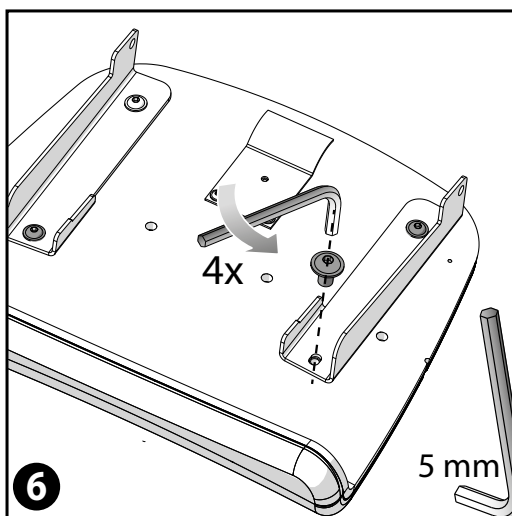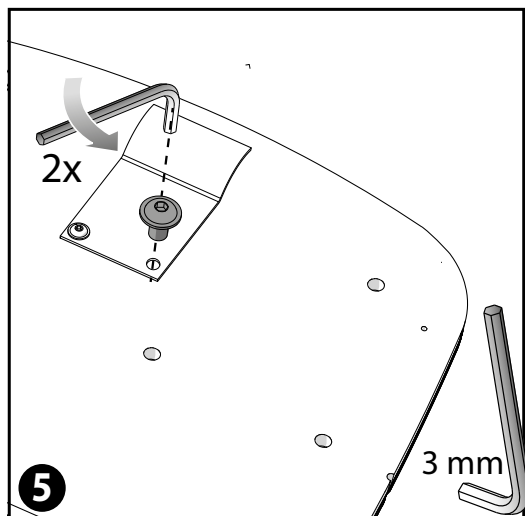EN [5]Remove the spring from below the seat. [6]Remove the brackets from the seat. [7]Fit the Active module to the seat. [8]Remove the cover from the armrest. Do so by carefully pushing the click latch on the side of the cover using a screwdriver. [9]Fit the new cover. [10]Remove the cover from the right armrest.

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





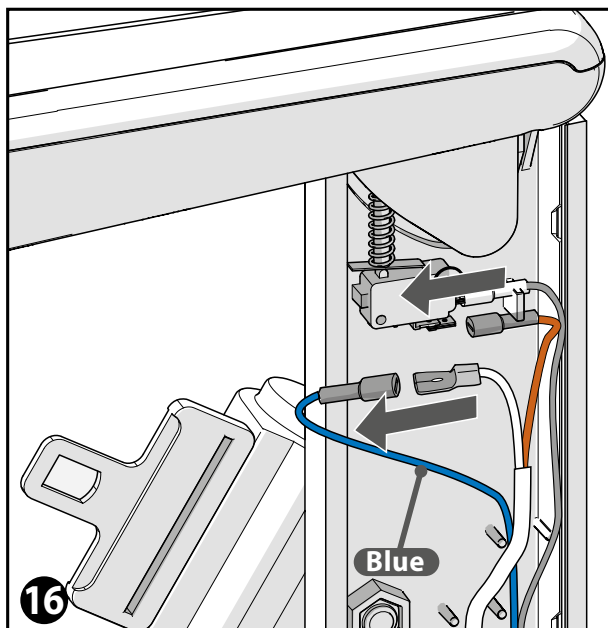
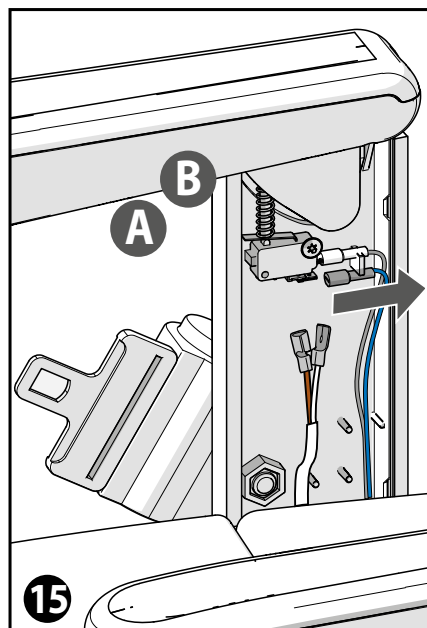
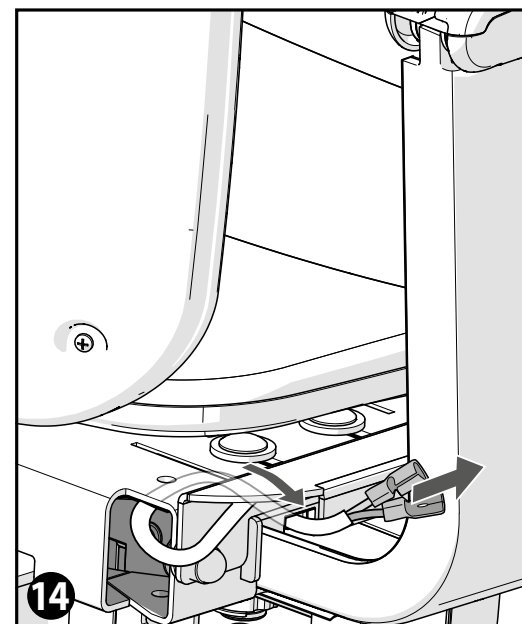
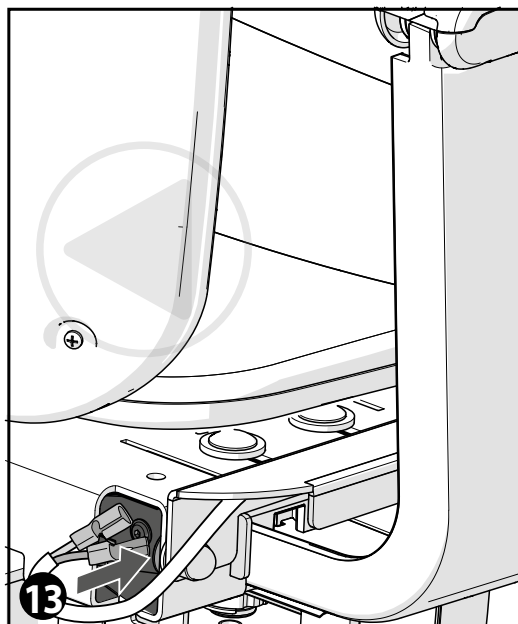
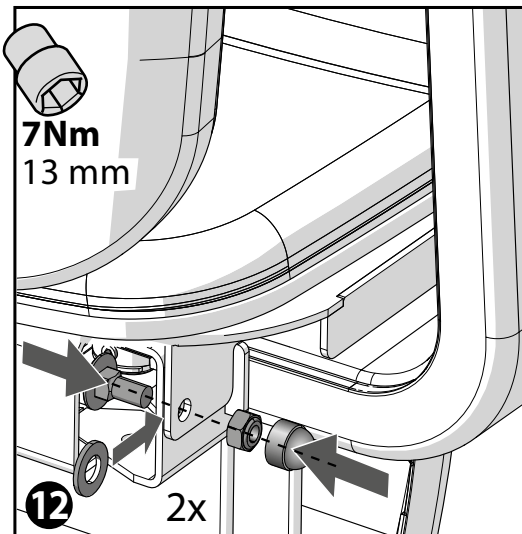
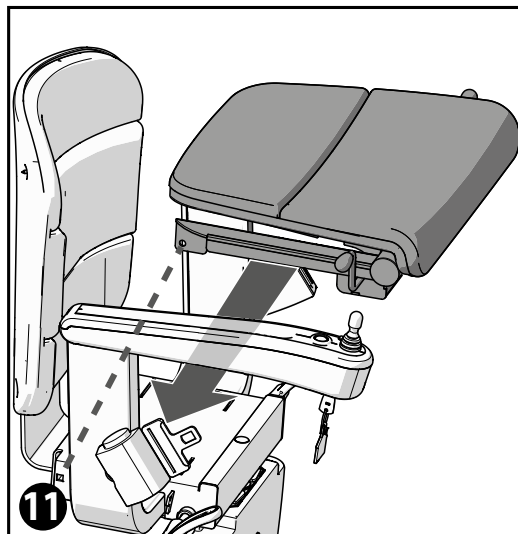
EN [11] Fit the complete Active module to the seat. [12] Reinstall the bolts and nuts of the seat's swivel axle. Remember to place the washer between the frame and the module. Hand tighten the nut. [13] Insert the module cable into to hole on the back of the module. [14] Guide the cable up through the square hole in the module near the armrest. [15] Disconnect the blue wire from the switch. [16] Connect this connector to the white wire of the module(A). Connect the other wire of the module to the available connector of the switch (B).

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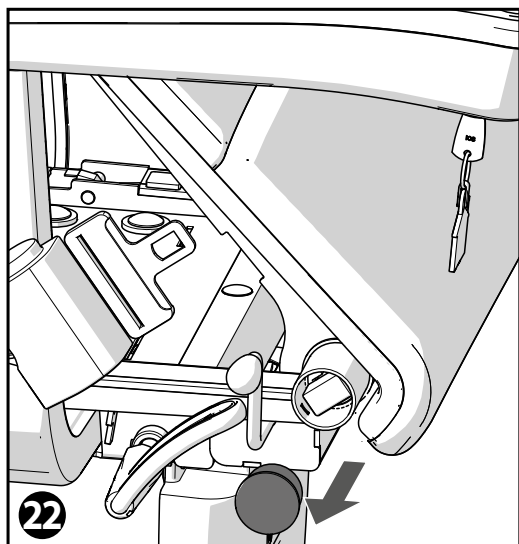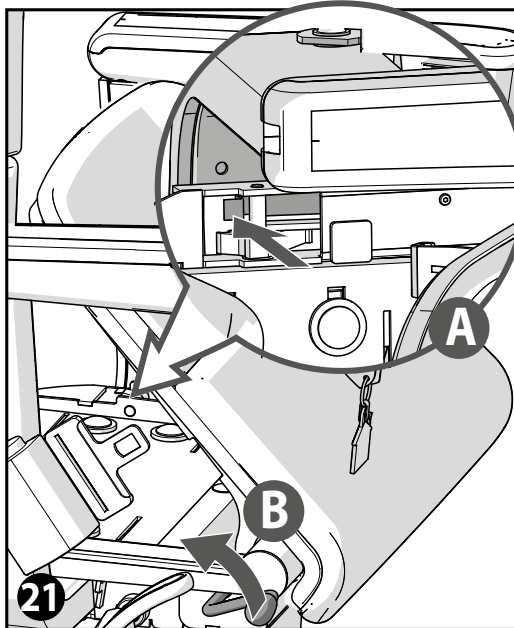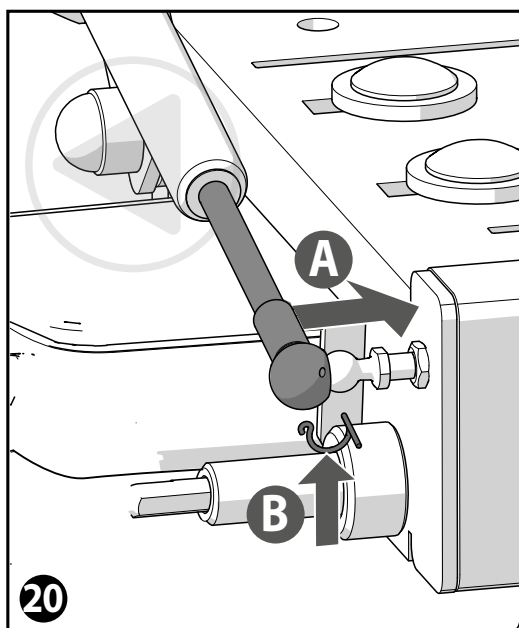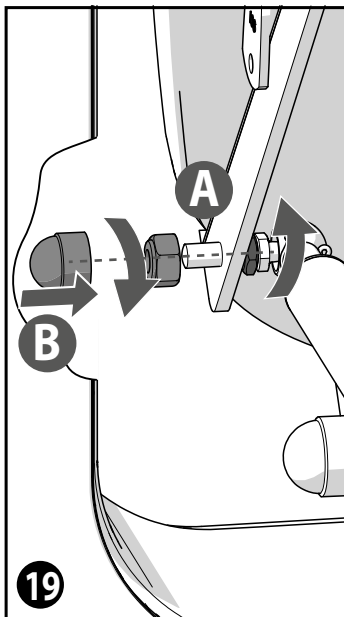
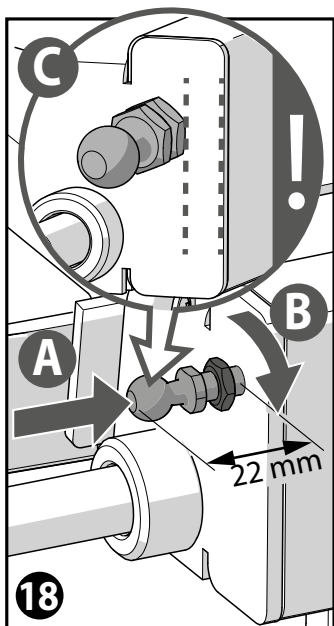
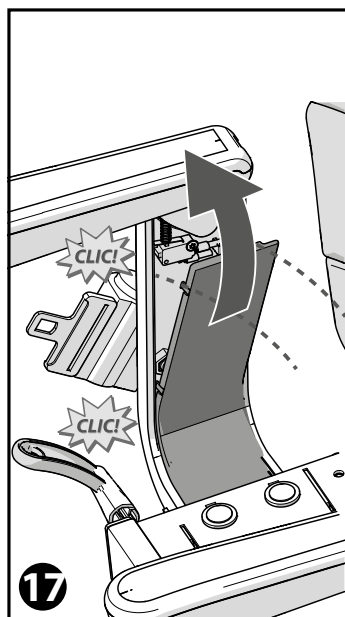
EN [17]Refit the cover. [18]Mount the fixing head of the gas spring in the threaded hole: it should stick out 22 mm[A]. Fasten the head with the locking bolt properly[B]. The position of the nuts must be parallel to the frame [C]. [19]Mount the gas spring to the seat bracket with the 2 nuts[A]. Place the nut cap[B], [20]Fix the other side of the gas spring to the fixing point[A] and lock the assembly with the locking pin[B]. [21]Lightly press the leaf spring at the indicated spot[A] and move the lever upward[B]. [22]Remove the revealed cap.

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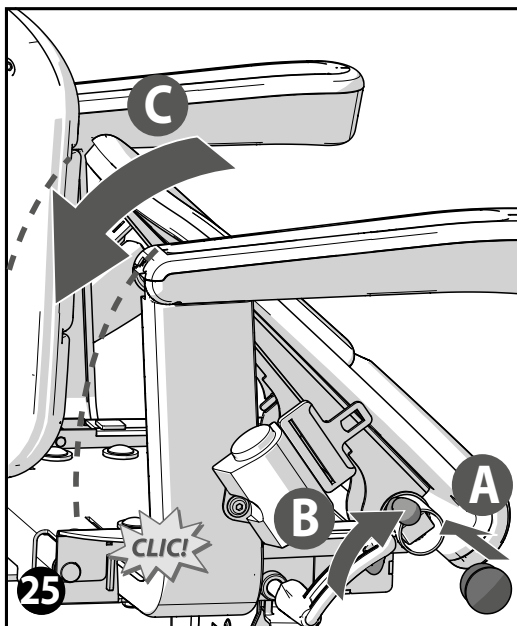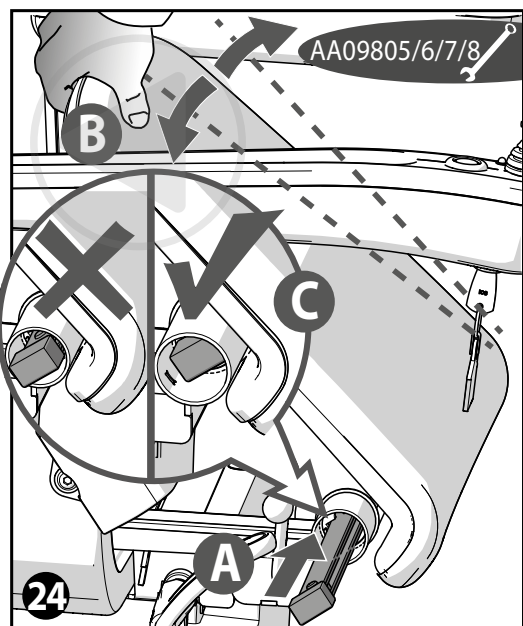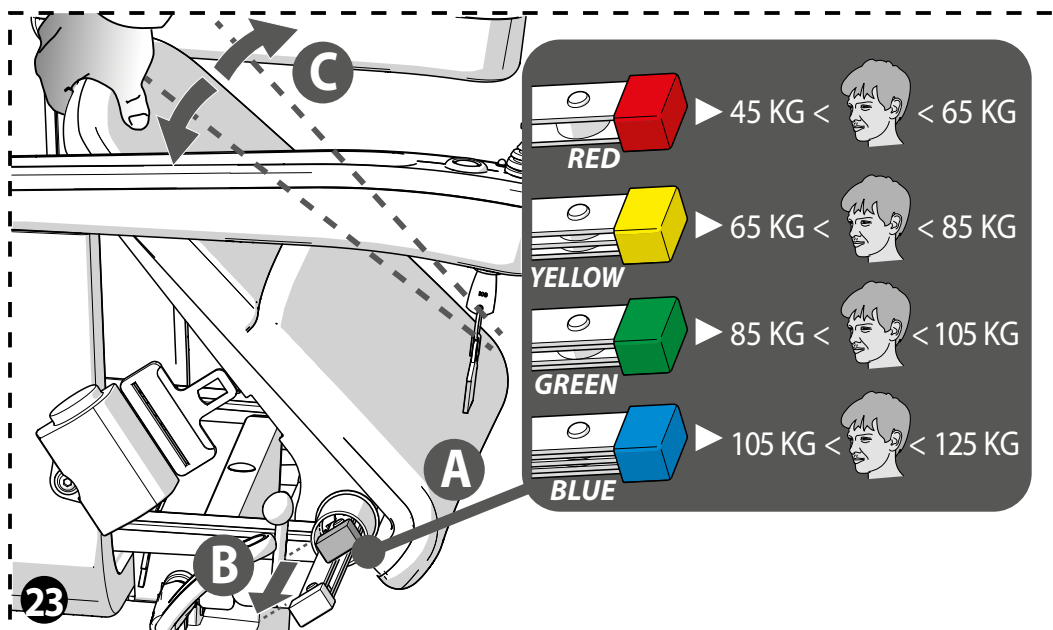
EN [23]Verify whether a spring has been preinstalled(A). Remove the spring if it is not the right one(B). Moving the seat back and forth will make it easier to remove the spring(C). [24]Fit the new spring. Make sure the spring is pushed far enough into the tube(C). [25]Refit the cap(A) and return the lever to its original position(B). Fold the seat down into its latch(C).

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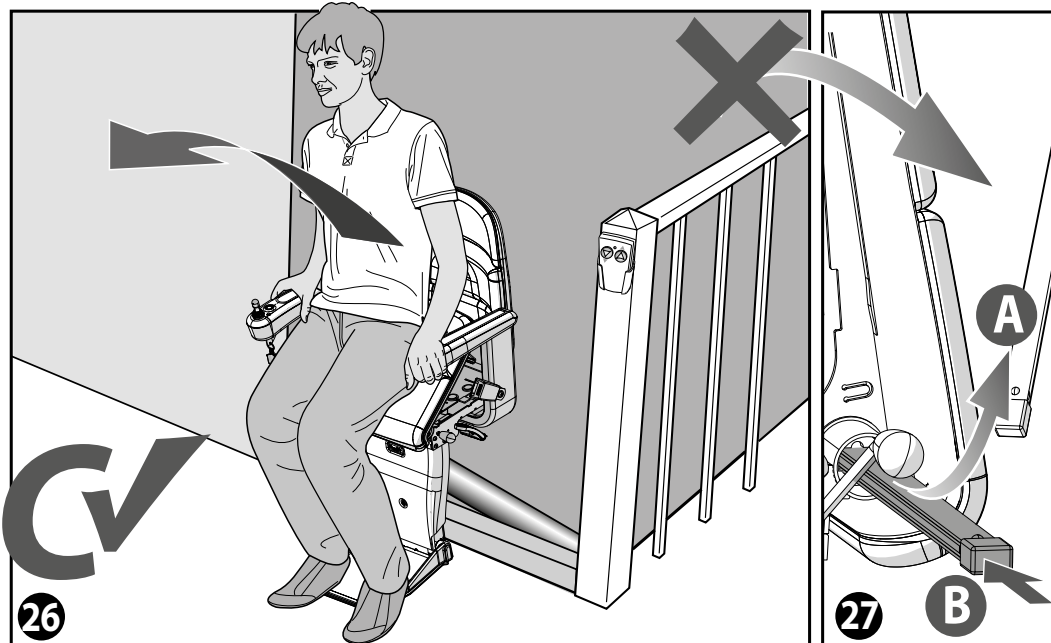EN [26]Check whether the user is correctly supported in getting up off the seat. [27]If necessary, fit a different leaf spring. [28]Check whether operation of the lever prevents the lift from moving.

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






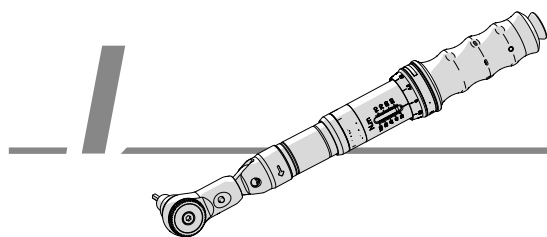
TORQUE WRENCH TEST


CHECK FASTENERS


After the test drive, check the torque of the load-carrying fasteners with a torque wrench

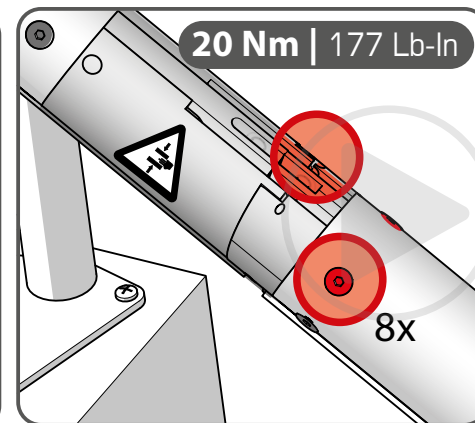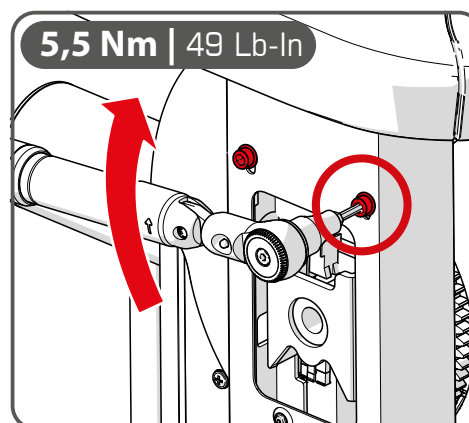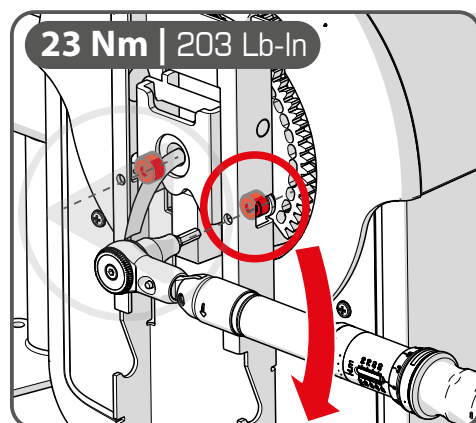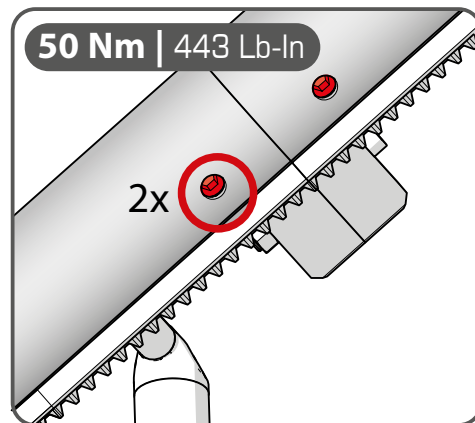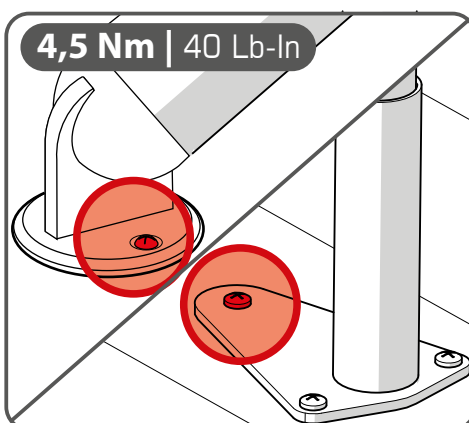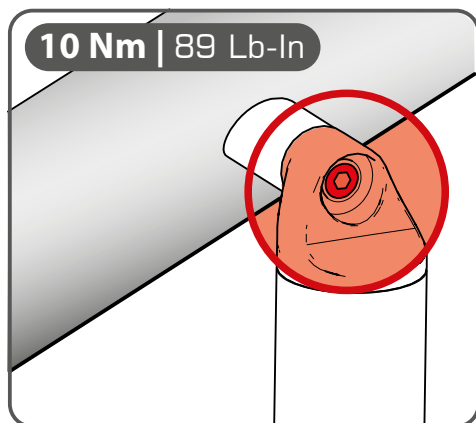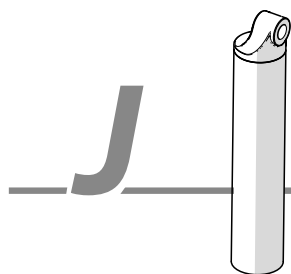UPPER SUPPORT TYPES

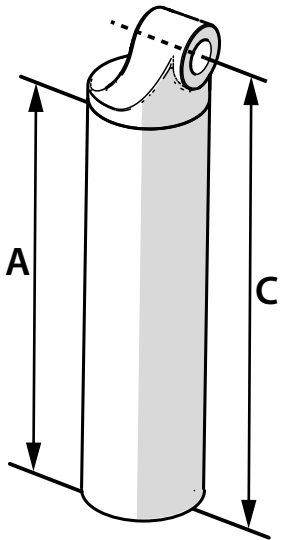	A	C	WHITE	GREY	BROWN	KIL/ZINC
	105	128	C 2051500	C 2051510	C 2051520	C 2051530
	135	158	C 2051380	C 2051390	C 2051400	C 2051470
	150	173	C 2051010	C 2051020	C 2051030	C 2051310
	175	198	C 2051410	C 2051420	C 2051430	C 2051480
	195	218	C 2051440	C 2051450	C 2051460	C 2051490
	250	273	C 2051040	C 2051050	C 2051060	C 2051320
	350	373	C 2051070	C 2051080	C 2051090	C 2051330
	450	473	C 2051100	C 2051110	C 2051120	C 2051340
	550	573	C 2051130	C 2051140	C 2051150	C 2051350

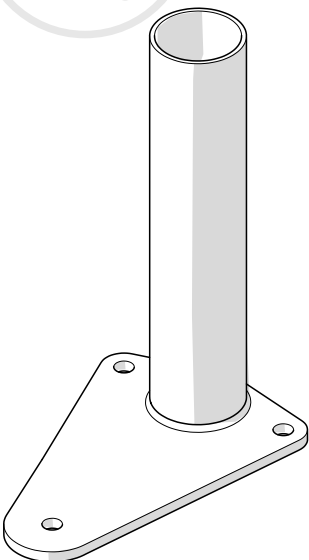	SIZE	WHITE	GREY	BROWN	KIL/ZINC
	SHORT	A0020070	A0020080	A0020090	
	MEDIUM	A0020040	A0020050	A0020060	
	LONG	A0020130	A0020140	A0020150	

K **PROTECTION EQUIPMENT** EXPLAINED

When performing any work on the stairlift, take the following precautions:

- Wear safety shoes type S3 or equivalent protection (S3 is an European classification; outside of Europe there are other classifications).
- Wear safety glasses for eye protection, i.e. to prevent injury because of parts suddenly striking into the direction of your face such as springs or when drilling in floors and walls.
- Use smooth fitting gloves to move and manage heavy / sharp parts of the unit and the rail.
- Wear a safety cap for installation with low bulkheads.
- Wear ear plugs when drilling and cutting. Exposure to sounds above 85 decibels can cause hearing loss.
- To prevent getting caught in moving parts of the stairlift: do not wear loose clothing or jewelry, such as necklace, bracelets and tie long hair.

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IT

**VERY IMPORTANT****EN**

This product must ONLY be installed by properly qualified and trained stairlift installation engineers, who have satisfactorily completed a Savaria-approved training course on how to install and/or commission the FC4-stairlift. Grave consequences regarding the safety of the stairlift in question and of its users may result from any violation of this requirement.

Intended Audience

This manual is intended for use by such qualified and trained stairlift installation engineers ONLY; it is NOT intended to show an unskilled person how to install and/or commission the stairlift, nor does it allow an unskilled person to install and/or commission the stairlift in an appropriate/safe manner.

This manual may NOT be used to replace the necessity for training.

The procedures described in this manual are 100% mandatory and MUST be followed.

Required Knowledge

All users of this manual MUST have the following knowledge:

Sufficient understanding of the key features, functions and capabilities of the FC4 stairlift;

Sufficient mechanical and electrical skills to be able to understand the function and purpose of each component of the stairlift, and sufficient mechanical and electrical skills as required for the correct installation and commissioning of the stairlift;

Sufficient Health and Safety awareness.

Disclaimer

It is the responsibility of the stairlift engineer to check that the manual is applicable to the type and version of the stairlift being installed; he/she MUST not fit a stairlift by reference to the manual if it is not the correct manual.

The stairlift engineer must read the manual to understand all the information to enable safe installation and/or commissioning of the stairlift.

In the absence of any negligence or breach of duty on the part of Handicare Accessibility Limited, no liability or responsibility will be accepted by it for any loss or damage arising from anyone's failure to follow and adhere to the guidance in the manual in connection with the installing and/or commissioning of the stairlift.

The latest and any previous versions of this manual can be obtained from www.handicarepartners.com

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Unit serial no:

Seat serial no:

- 1) Confirm that the installation dates have been completed on the EC Declaration of Conformity
- 2) Confirm that the remote controls are mounted on the wall and are functioning properly
- 3) Confirm that the safety belt is secure and functions correctly
- 4) Confirm that all warnings relating to the maximum working load on the lift have been affixed

2 Rail

- 1) Confirm that the rail joints are properly secured
- 2) Confirm that the rail supports are properly secured
- 3) Confirm that the folding rail is functioning correctly

3 Electrical drive unit tests

- 1) Confirm that the batteries are securely fitted in the carriage structure
- 2) Confirm that the following have been effectively indicated to the user:
 - a) location of the charge contacts.
 - b) "OFF CHARGE" alarm.
 - c) required action when stairlift is not charging

4 Clearance dimensions

- 1) Confirm the correct clearance dimensions as specified in the installation manual are maintained throughout the full travel of the stairlift
- 2) Confirm that there are no open bars/openings at the rear of the lift
- 3) Confirm that other users can safely use the staircase (minimum clearance of 400 mm between the folded seat and opposing obstacles). Install handrails if required.

5 Functional tests

- 1) Confirm that the final limit switches (upper and lower) operate correctly (position and operation)
- 2) Confirm that the control lever, key switch, emergency stop button, and the detection switches in the armrests operate correctly
- 3) Confirm that the safety detection devices operate correctly (footrest/unit/armrest)
- 4) Confirm that the swivel seat interlock and switch operate correctly
- 5) Confirm that the stopping distance of the lift is a maximum of 2 cm
- 6) Confirm that the restart time of the lift is at least 1 second
- 7) Confirm that the alarm device operates correctly when activated
- 8) Confirm that the lift has successfully passed the static (torque wrench test) and dynamic test ride in accordance with the instructions in the installation manual

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
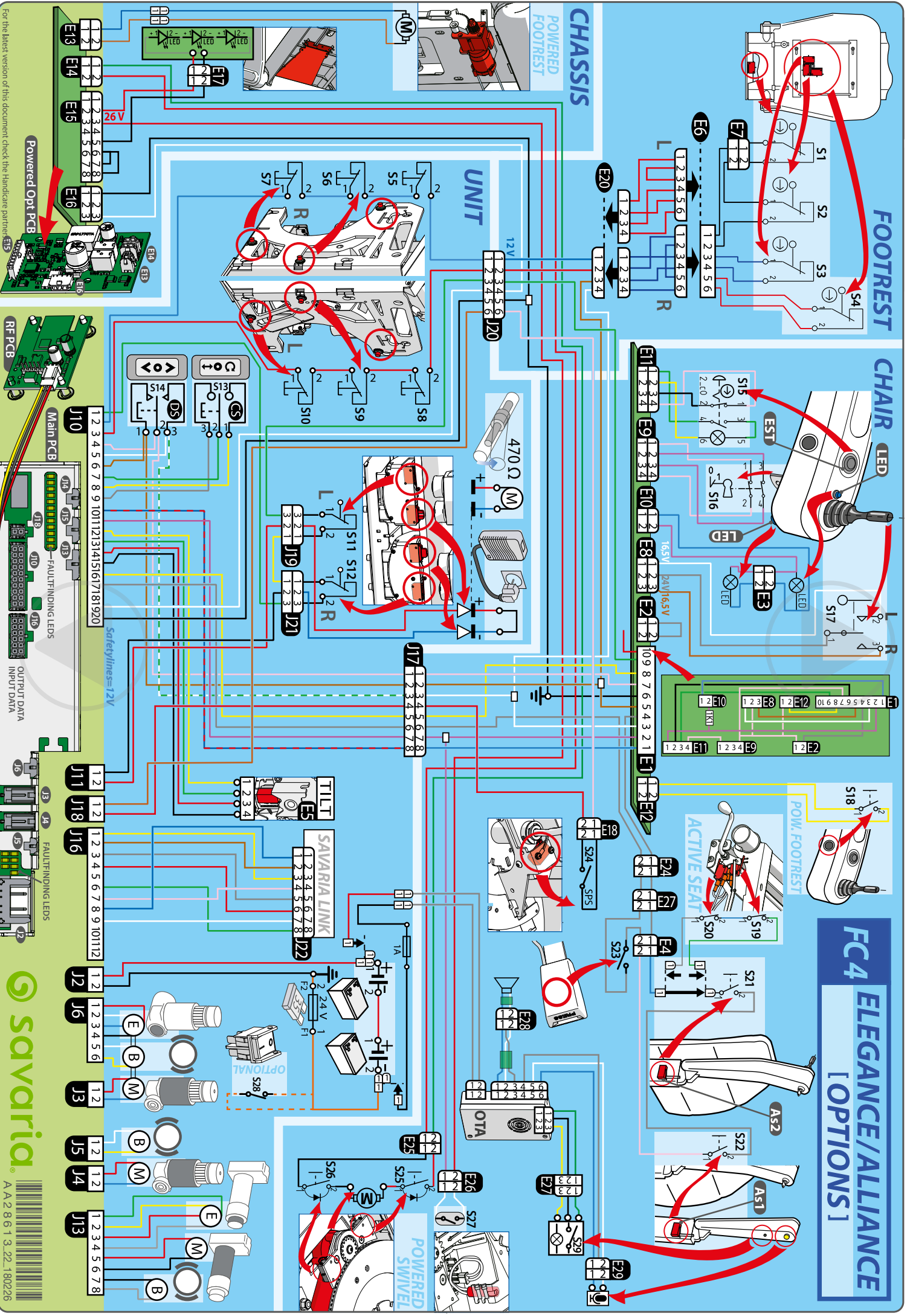
UNIT

POWERED SWIVEL

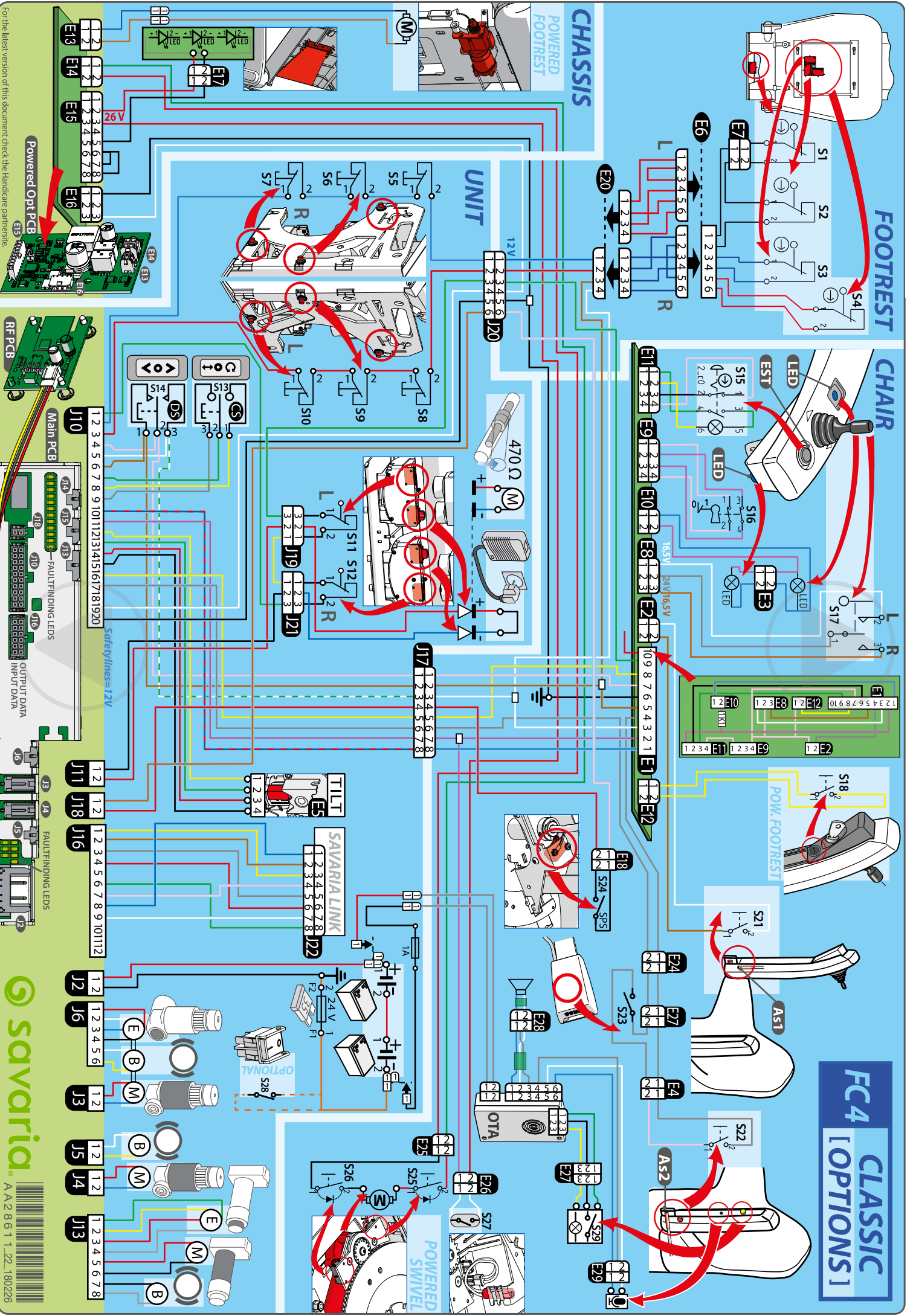For the latest version of this document check the Hardware partner site.

**Savaria Nederland B.V.**

Newtonstraat 35
1704 SB Heerhugowaard
The Netherlands
T +31 (0)72 576 88 88
traplift@handicare.com

Savaria Italia s.r.l.

Via Chiaviche 17
46020 Pegognaga (MN)
Italy T +39 (0)376 143 3516
Montascale@handicare.com

Savaria Lifts (UK) Ltd

82 First Avenue,
Pensnett Estate Kingswinford,
West Midlands DY6 7FJ,
England T +44 (0)1384 408700
enquiries@handicare.co.uk

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Savaria Canada

2 Walker Drive Brampton,
ON L6T 5E1 Canada
T 1.905.791.5555
T 1.855.728.2742
info@savaria.com

www.handicare .com

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Lübbecker Strasse 240
32429 Minden
Deutschland
tel+49 (0) 571 97 33 98 -55
tps.de@handicare.com

Date	Version	Change
08-05-26	11	Switching the installation method of linking remotes and calibrating unit
25-05-26	12	Changed bolts and bracket hinged rail cover and endstop size
01-06-26	13	Floorplates
13-07-26	14	Remove torque endstop bolts

INSTALLATION MANUAL

MANUAL D'INSTALLATION
INSTALLATIEHANDLEIDING
INSTALLATIONANLEITUNG
MANUALE PER L'INSTALLAZIONE
MANUAL DE INSTALACION
MANUAL DE INSTALACAO
INSTALLASJONSMANUALEN
INSTALLATIONS MANUAL
INSTRUKCJA

AA28609_14_130726

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Savaria Deutschland GmbH