H 4000

| INSTALLATION MANUAL


INSTALLATIEHANDLEIDING

MANUAL D'INSTALLATION

INSTALLATIONANLEITUNG

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MANUAL DE INSTALACION

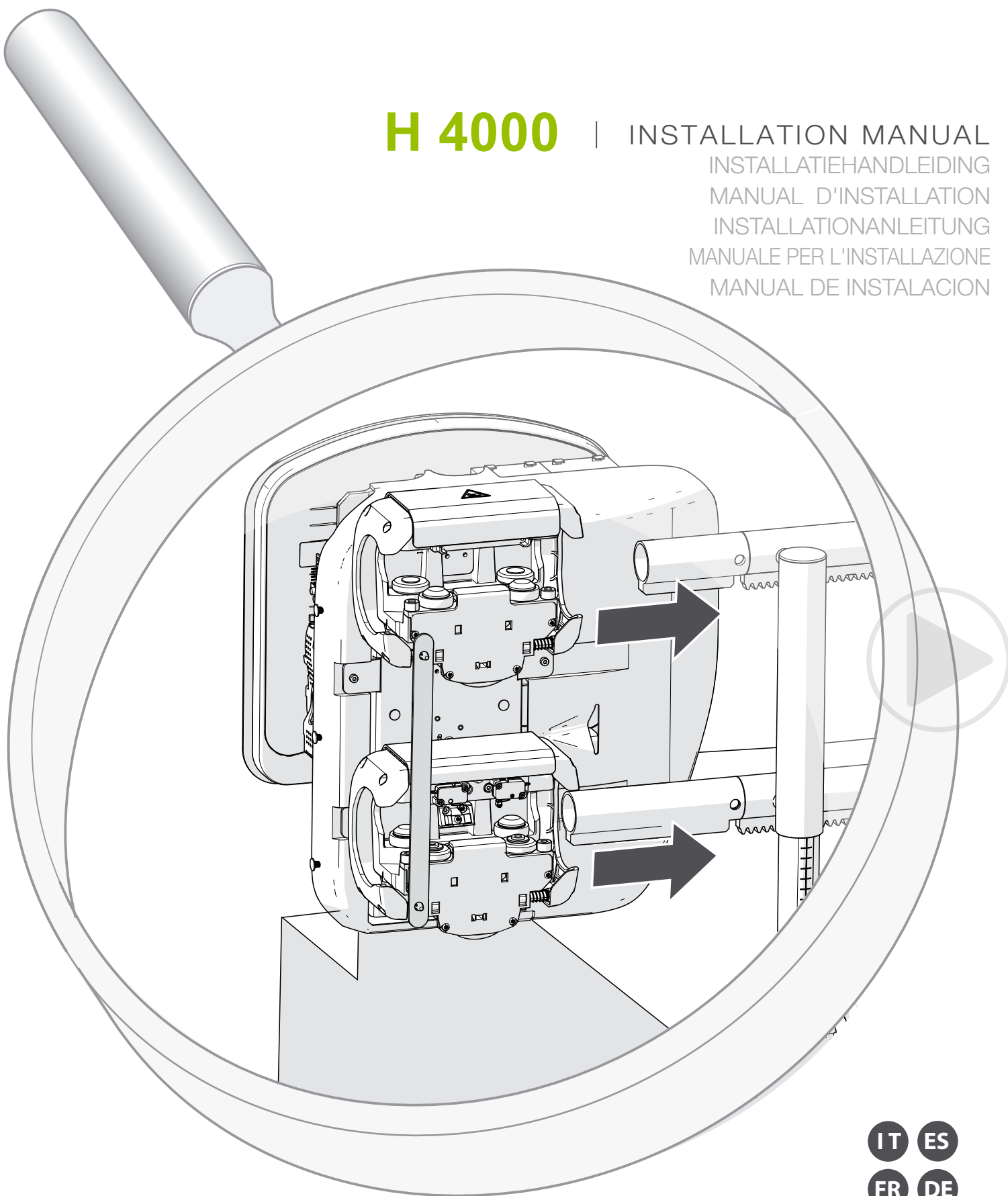Style
Simplicity
Smart

 **savaria®**

EN All rights reserved. No part of this publication may be reproduced and/or published by means of print, photocopy, microfilm or in any other way or form, without prior written permission from the manufacturer. This also applies for the accompanying drawings and diagrams. The manufacturer reserves the right to change parts at any time, without prior or direct notification to the customer. The contents of this manual may also be changed without prior warning. This manual is intended for the standard models of the lift. The manufacturer therefore cannot be held liable for any damage resulting from specifications that differ from those of the standard model of the lift supplied to you.

This manual has been compiled with the best intentions and every possible care but the manufacturer cannot accept responsibility for any errors in this manual or for the consequences thereof.

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EN
The English manual is the original language

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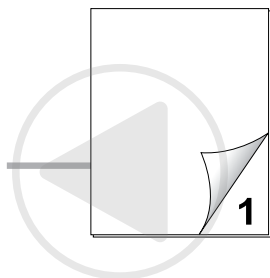
1 - USED TOOLS AND PARTS


2

Thank you for choosing a Handicare stairlift.
This manual is for the Handicare 4000 lift with a manual seat or powered seat type. The manual includes instructions for the standard and the optional features.

2 - TRACK ASSEMBLY


6

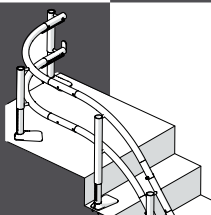Check stairs [p.06]
Caps [p.08]
First stanchion [p.10]
Second stanchion [p.12]
Position stanchion [p.14]
Top rail-first section [p.16]
Intermediate stanchion [p.18]
Last stanchion [p.20]
Top rail-sec. section [p.22]
Bottom rail-sec. section [p.24]
Repeat procedure [p.26]
Level stanchions [p.28]
Check track [p.30]

3 - POWER PACK


32

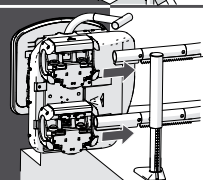Loading bars [p.32]
Batteries [p.34]
Preparing power pack [p.36]
Alignment trunnions [p.38]
Remove instal. parts [p.40]

4 - TRACK ELECTRICAL


42

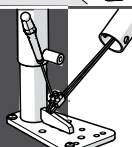Charge contacts [p.42]
Charge wires [p.44]
Remove the paint [p.46]
Earth wire and charger [p.48]
Terminal block [p.50]
End stop [p.52]
Tiewraps [p.54]
Power supply [p.56]

5a - STYLE SEAT


58

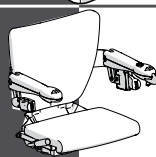Seat stem [p.58]
Seat stem heavy duty [p.60]
Locking bolt [p.62]
looms through stem [p.64]
Seat loom [p.66]
Dipswitch settings [p.68]
Closing front cover [p.70]

5b- SIMPLICITY SEAT


70

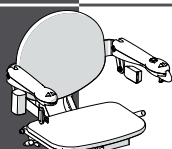Seat stem [p.72]
Seat stem heavy duty [p.74]
Locking bolt [p.76]
looms through stem [p.82]
Seat loom [p.80]
Closing front cover [p.82]

5c- SMART SEAT


82

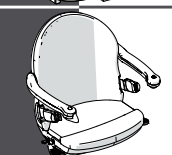Seat stem [p.82]
Seat stem heavy duty [p.84]
Fixation seat back [p.86]
Seat looms [p.88]
Seat base wiring [p.90]
Upholstery [p.92]
Closing front cover [p.94]

6 - TRACK FIXATION

96

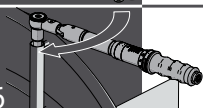Check [p.96]
Tighten bolts [p.98]
Wall brackets [p.100]
Apply grease [p.102]

7 - REMOTE CONTROL

104

Mounting [p.104]
Linking [p.106]
Repeater [p.108]
Changing channels [p.110]
Check [p.112]

8 - CHECK

114

Measurements [p.114]
Swivel - safetybelt [p.116]
Option functions [p.118]
Safety limits [p.120]
Armrests [p.122]
Keyswitch - stops [p.124]
Belt switch - hinge [p.126]
Instruct user [p.128]
Instruct user [p.130]
Declaration [p.132]
Placing label [p.134]

In this manual, the paperclip symbol refers to a relevant attachment.

ATTACHMENTS

137

A - Hinged rail [p.138]
B - Emergency operation{TR} [p.146]
C1 - Emergency operation{SM} [p.148]
C2 - Emergency operation{SL} [p.149]
D - OSG [p.152]
E - Intermediate stop [p.160]
F - Armrest conversion [p.162]
G - RF Repeater [p.174]
H - No charge bracket [p.186]
I - Dipswitch label [p.192]
J - Seat height [p.193]
K - Alternative fixing [p.194]
L - TorqueWrenchTest [p.195]
M - Protection equipment [p.196]
N - Competence [p.197]
O - Installation forms [p.200]

STAIRLIFT TIGHTENING TORQUE

UNLESS OTHERWISE STATED




Size (mm)	Material	Max. Torque (Nm)					
		Bolt Grade					
		4.6	4.8	5.8	8.8	10.9	12.9
M6	Self Colour	3	3.8	4.8	8.6	12	14.5
M6	Zinc	3.4	4.3	5.6	9.9	14	16.5
M8	Self Colour	7.2	9.1	11.7	21	29	35
M8	Zinc	8.3	10.5	13.6	24	34	40
M10	Self Colour	14.4	18	23.4	42	58	70
M10	Zinc	16.6	21	27	48	67	81

Nm

HAND TIGHT

Not loose, but not really tight either.
Somewhere around 1 or 2 Nm.

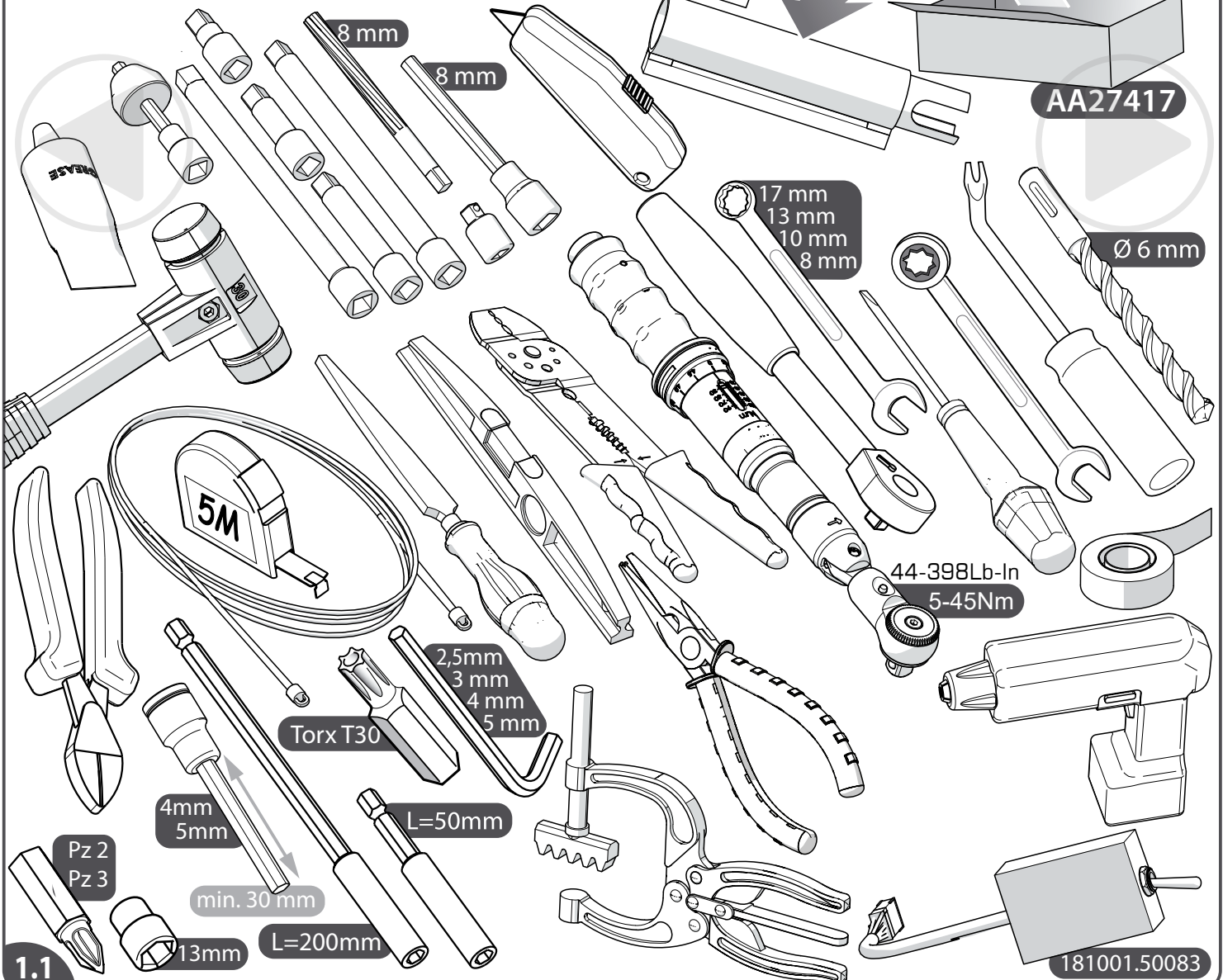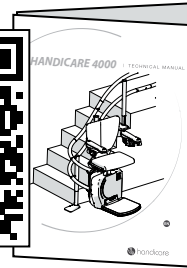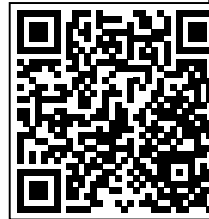Only a technician certified for this lift should install this lift (for more detailed information about the required competence see attachment of this manual). General required knowledge can be found in the 4000 Technical Manual. The latest instructions for the lift and instructions for options are available for download at www.handicarepartners.com.

Ensure good and sufficient tools; minimum required, in addition to installation kit, to install the stairlift:

File	Allen keys {2,5mm,3mm, 4mm, 5mm}	Tape	Crimper	Rail clamp
Loading toggle	Tape-measure {min. 5m}	1/2" racheting socket wrench	Pliers	Pull through
Loading bars	Spanners {8mm,13mm,17mm}	Shorting connectors	Track lubricant(HTEP)	Torx bit {T30}
Drill, Drills {8mm}	Screwdrivers with fine flat head	Stanley knife	1/2" in hex socket driver	Spirit level
Phillips bit {Pz 2, Pz 3}	Bit extenders {2x250mm, 125mm}	Allen key(ball ended) {6mm}	Torque wrench	Reamer {8 mm}
Adaptor cable	Rivet Removal Tool	Hamer		

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
IT

The diagram shows two rods. On the left is a straight rod with a rectangular cross-section. On the right is a rod that is bent into a curved shape, also with a rectangular cross-section. This illustrates the concept of bending or deformation.

TOOLS & PARTS
USED PARTS

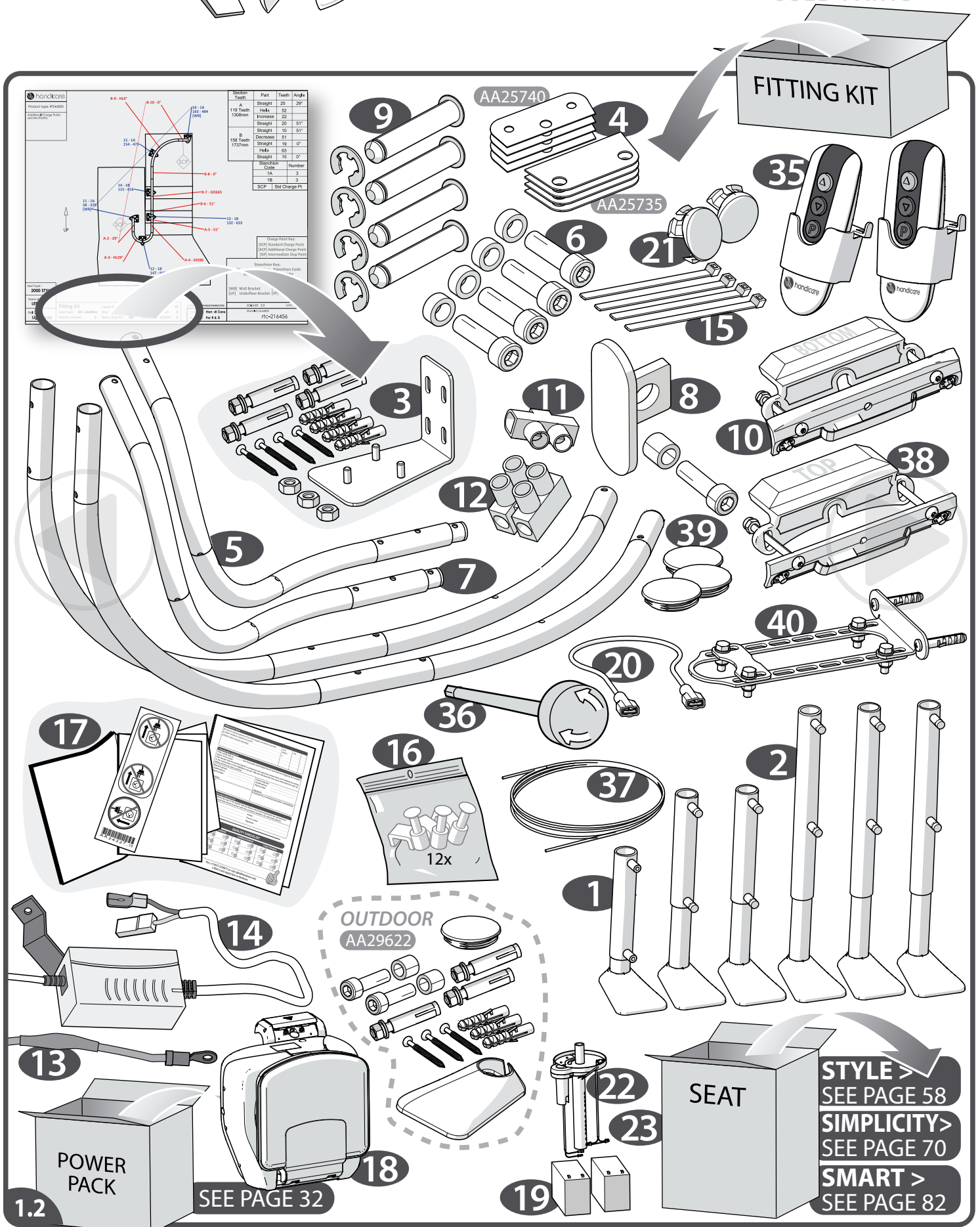For the best digital display, use the settings below (display>page display>display two pages+display frontpage):

[1.2] Check that all the parts are present. The part numbers shown next to the parts are given in the installation manual at the time of installation.

When replacing parts on an existing stairlift, always use parts with identical part numbers.

EN **Hand tight** = not loose, but not really tight either. Somewhere around 1 or 2 Nm.

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
2.1

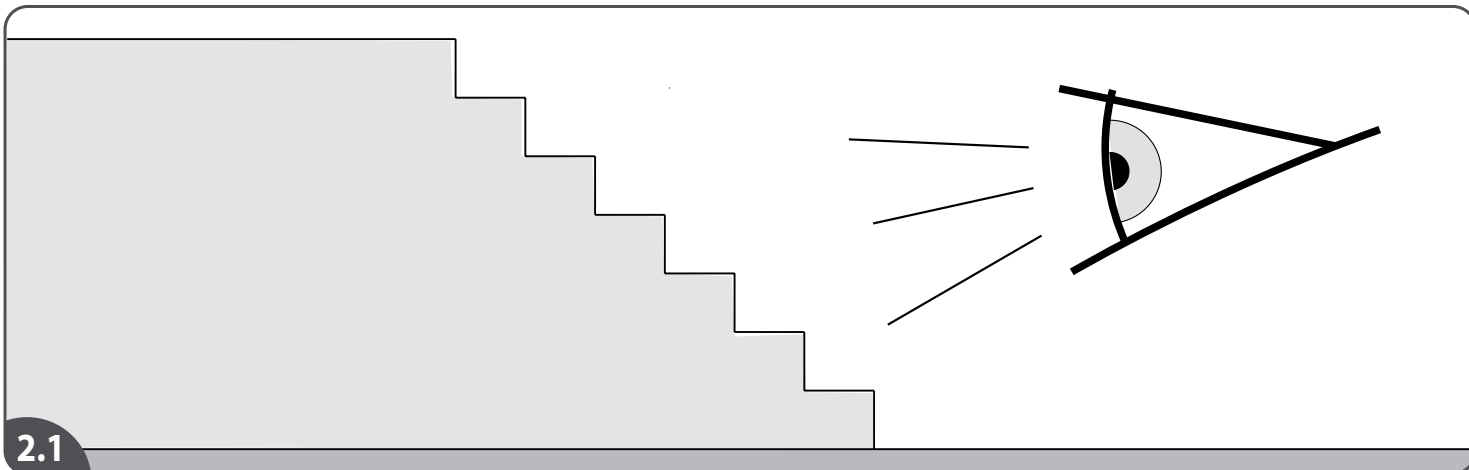Product type: RTC-4000

Additional Charge Points
(section/Tooth)

B-9 - HL0°

B-10 - 0°

16 - A
163 - 484
[UF]

15 - A

14 - B
115 - 616

11 - A
50 - 379
[UF]

UP

A-2 - 29°

A-3 - HL29°

12 - B

A-5 - 51°

A-4 - GI596

B-6 - 51°

B-7 - GD665

B-8 - 0°

Section Teeth	Part	Teeth	Angle
A 119 Teeth 1308mm	Straight	25	29°
	Helix	52	
	Increase	22	
	Straight	20	51°
B 158 Teeth 1737mm	Straight	10	51°
	Decrease	51	
	Straight	19	0°
	Helix	63	
	Straight	15	0°
	Stanchion Code		Number
	A		3
	B		3
	SCP	Std Charge Pt	

Charge Point Key:

[SCP] Standard Charge Point
[ACP] Additional Charge Point
[ISP] Intermediate Stop Point

Stanchion Key:

Stanchion Number - Stanchion Code
Setting Height - Total Height

[UF] Underfloor Bracket

Seat Type :
4000 STYLE

Hand of Rail :
LEFT

Rail Colour :
Light Grey

Fitting Kit

Hand Sets : **RF - 868MHz**
Number needed : **2**

Carpet Spacers : **0**

Wire (m) : **5**

Heavy Duty Kit : **NO**

Spread Plate Std : **16**

Spread Plate Fan : **0**

UnderFloor Standard : **0**

ALL DIMENSIONS IN MM UNLESS OTHERWISE STATED

CUSTOMER **Han di Care**
REFERENCE **For R & D**

DESIGNER: S.S
DRAWING NUMBER:

DATE: 20/11/18

rtc-216456

2.2

[2.1] Before you start, inspect the stairs. Clarify with your customer regarding any floor heating, alarm wires, pads or pipework that may affect/hinder installation prior to starting.

[2.2] The stanchion info is given in blue, near the stanchions{A}. The track parts info is given in red, near the parts{B}. If extra additions for the stanchions are needed, its given and explained on the drawing{C}. The locations of the charge points are given{D}. And as mentioned the fitting kit parts (kind and amount) are on the drawing {E}.

Note: Before assembly you must ensure that all of the required track is to the correct specification, the stanchions supplied are the correct length, the correct number have been supplied and the corresponding number of feet are available. In order to do this you need to cross reference the items supplied to the drawing supplied.

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
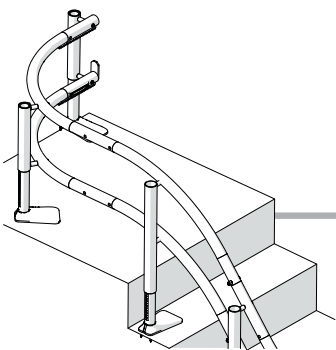
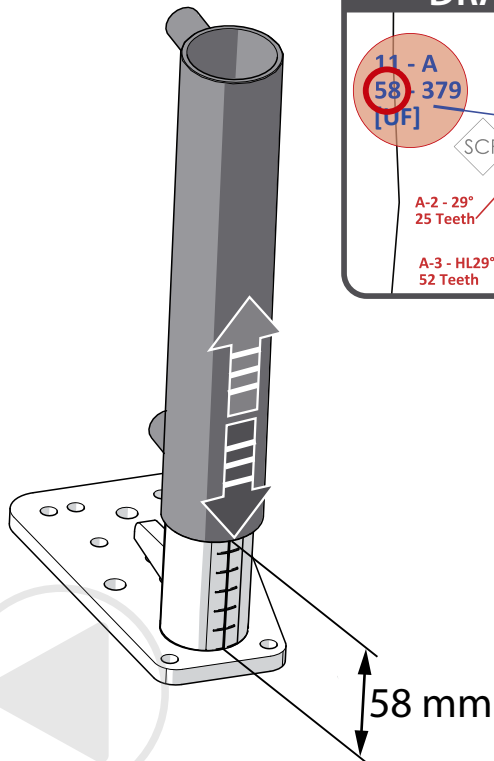
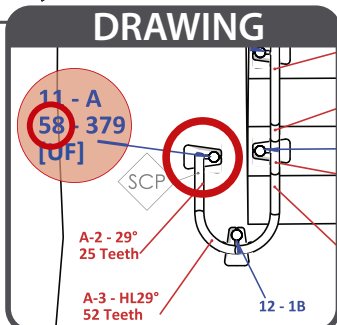
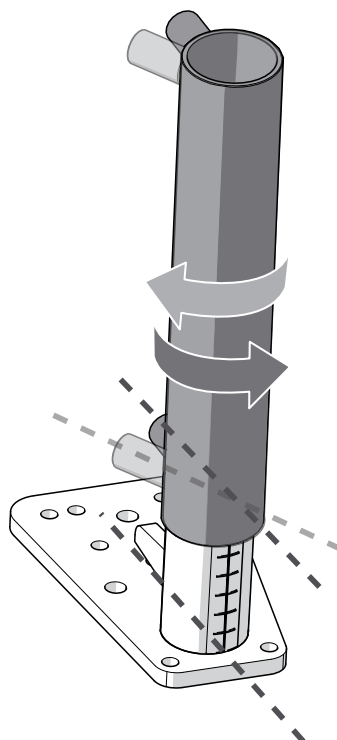
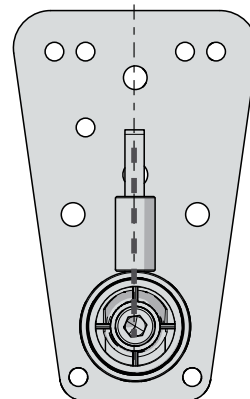
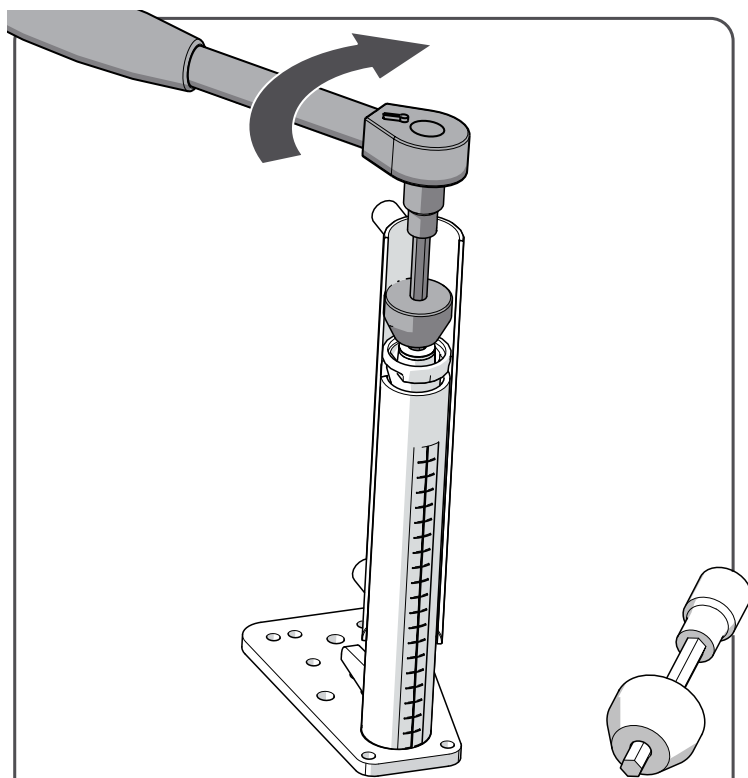
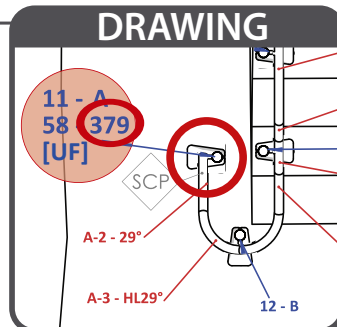
8 Nm

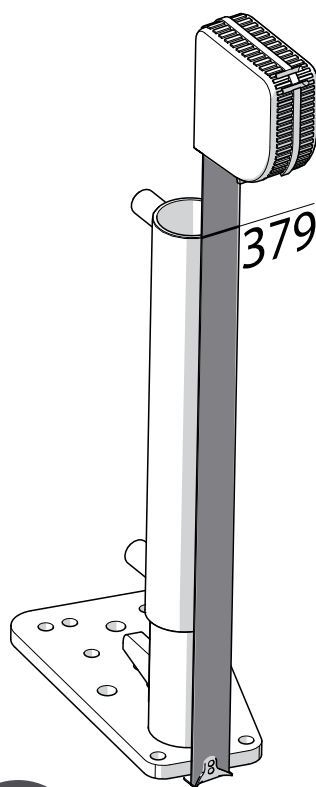2.9

[2.6] Prepare the first stanchion: *set the height of the first stanchion.*

[2.7] Prepare the first stanchion: *set the lugs in line with the foot.*

[2.8] Prepare the first stanchion: *tighten the lock bolt [hand tight]. **After the complete track installation the torque of the lock bolt need to be set on 45 Nm[398 Lb-In].***

[2.9] Prepare the first stanchion: *always check the total length of the first stanchion on a rail section.*

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




TRACK ASSEMBLY SECOND STANCHION

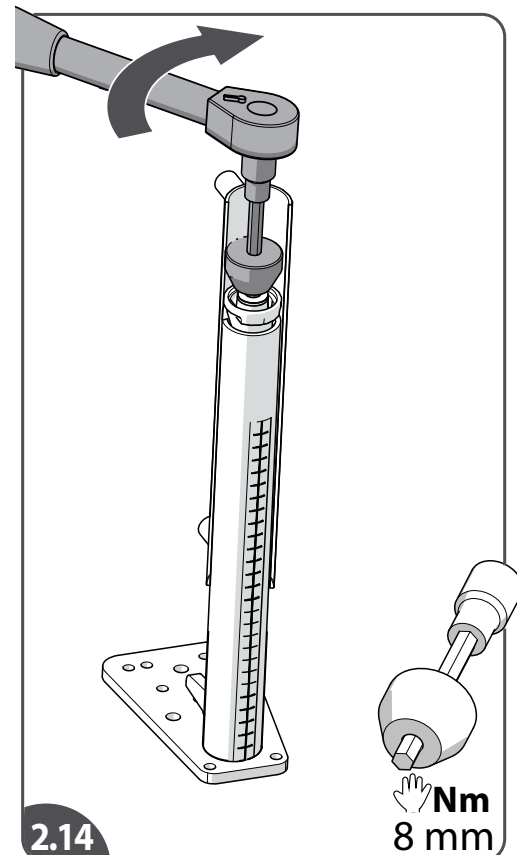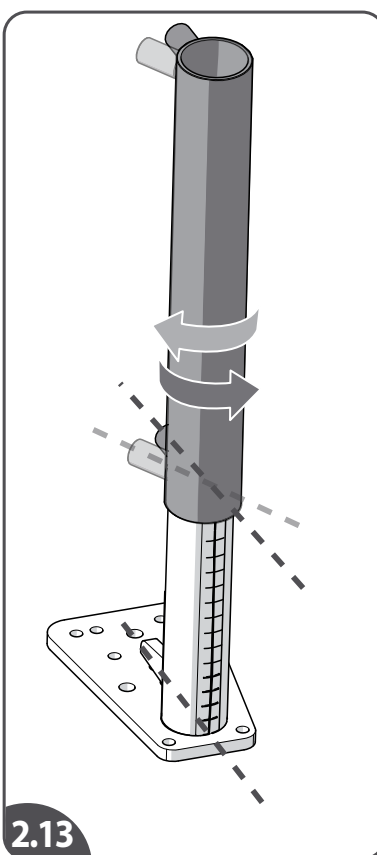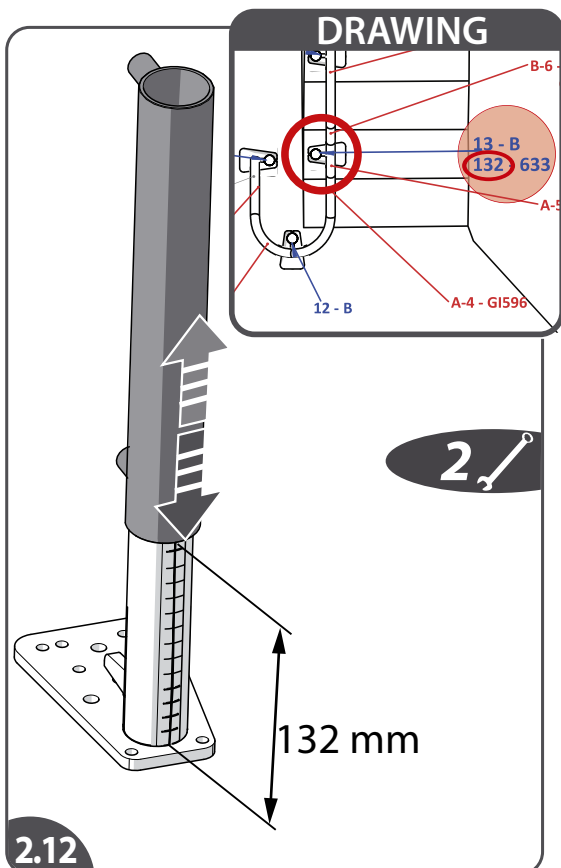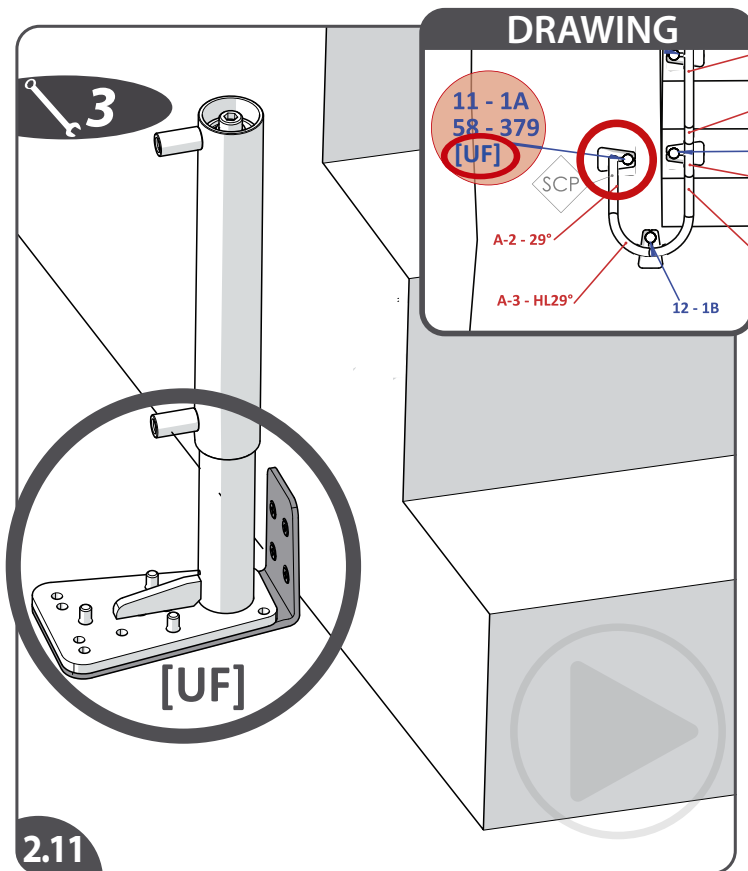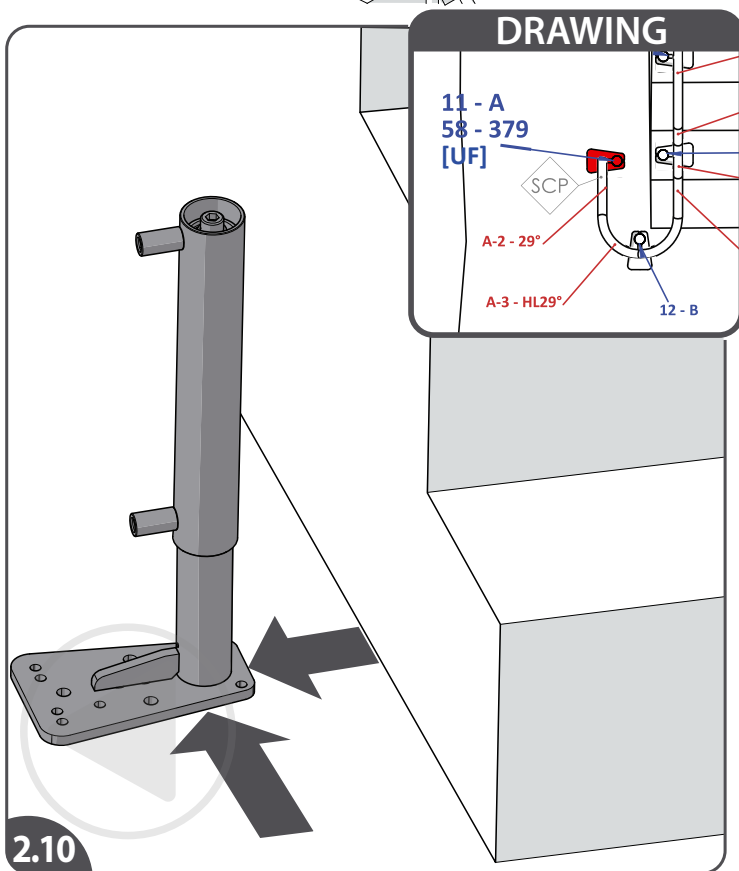[2.10] Position the stanchion.

[2.11] Check the drawing for additional fixations for the stanchion, e.g. for a Wall Mounted Stanchion Plate. A Wall Mounted Stanchion Plate is intended for when a stanchion cannot be fixed directly into the floor. They are designed to withstand the forces of the loaded stairlift. Do not install the additions at this stage.

[2.12] Set the last stanchion of the first section: *set the height of this stanchion*.

[2.13] Set the last stanchion of the first section: *set the lugs in line with the foot*.

[2.14] Set the last stanchion of the first section: *tighten the lock bolt [hand tight]*. Hand ratchet all the bolts of the track: do not use power tools.

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



IT

2

TRACK ASSEMBLY TOP RAIL- FIRST SECTION

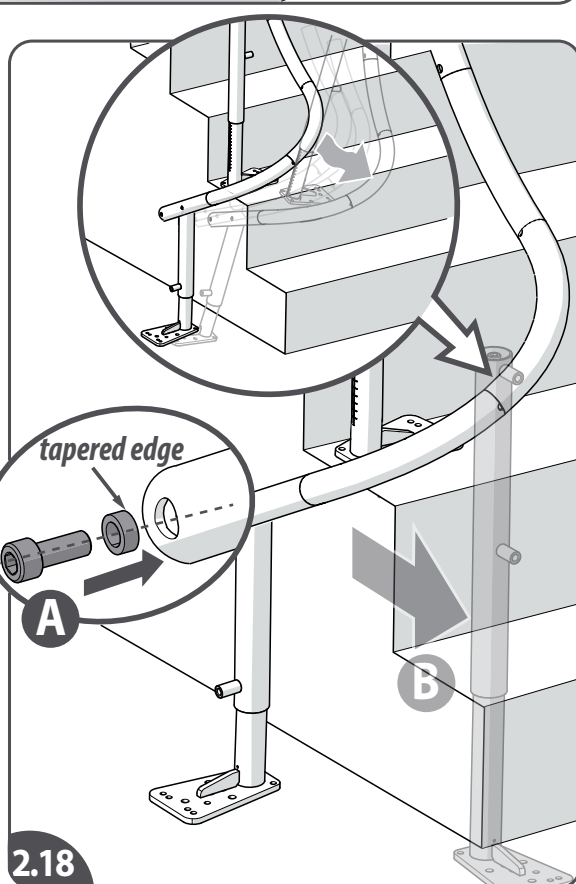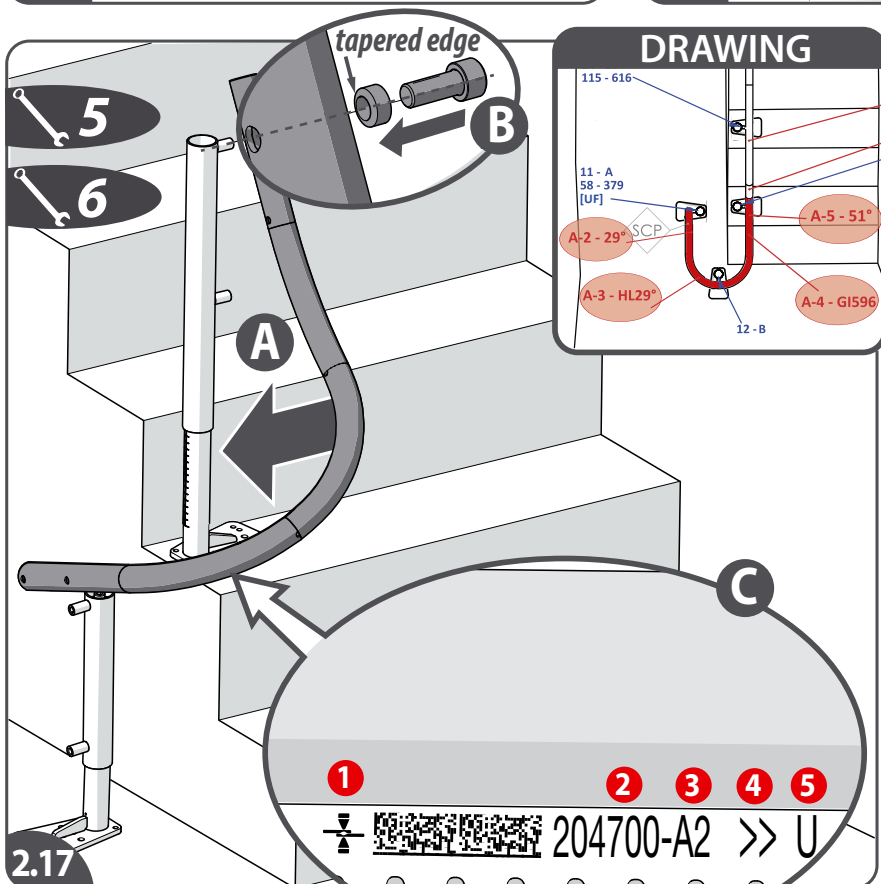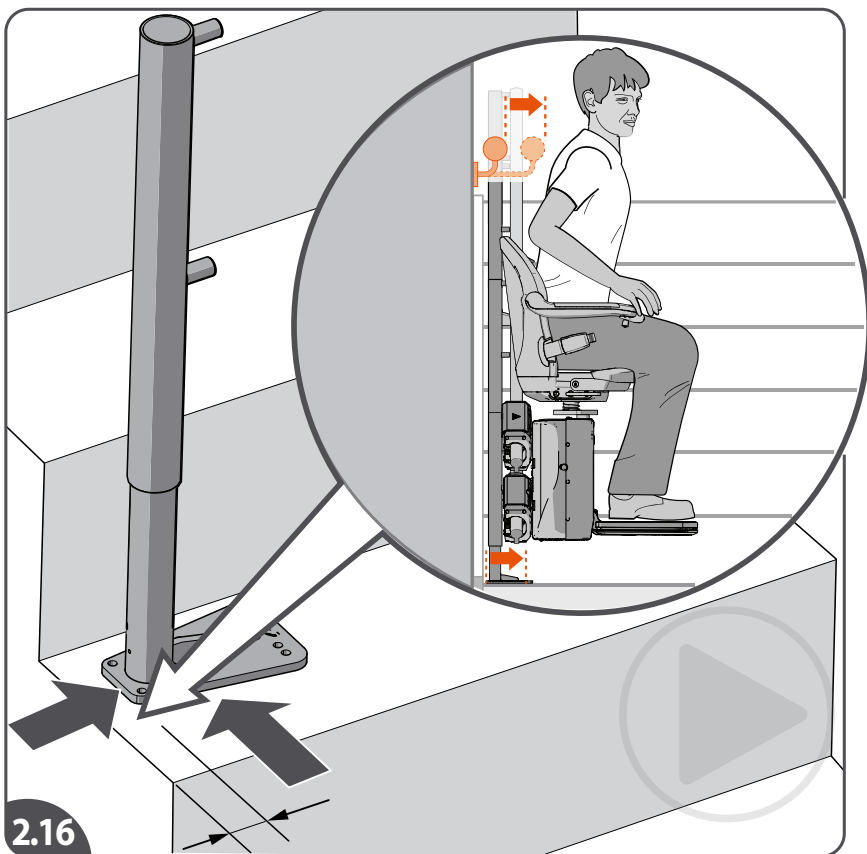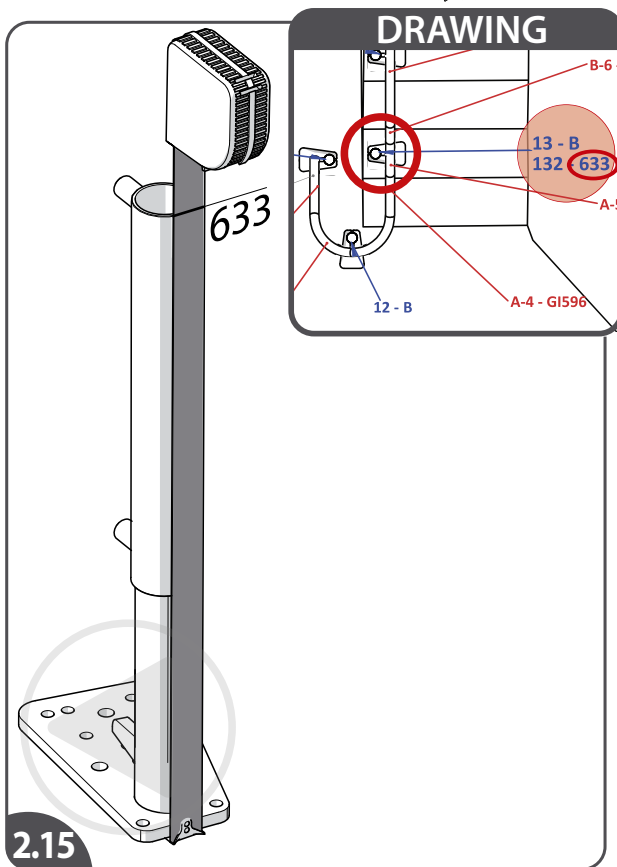[2.15] Set the last stanchion of the first section: *always check the total length of the last stanchion on a rail section.*
 [2.16] Position this stanchion in accordance with the track offsets stated on the drawing.
 [2.17] Install the top rail of the first section.
 {C} You can find specific information about a rail part on the label on the toothrack of this part: 1-Joint type (crimped, bolted or welded), 2-Project number, 3-Rail part number seen from below (indicates the rail section i.e A,B,C – A is always the lowest), 4-Direction (arrow point up the stairs), 5-Upper or lower railpart.
 [2.18] Use the bolts+spacers to install the top rail. For more stability additional stanchions in the rail section can be installed.

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



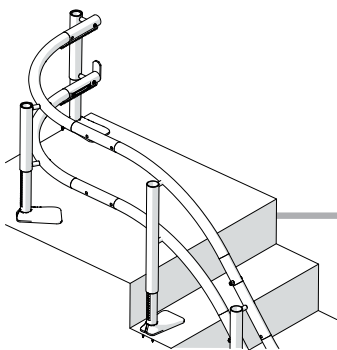
TRACK ASSEMBLY

BOTTOM RAIL- FIRST SECTION

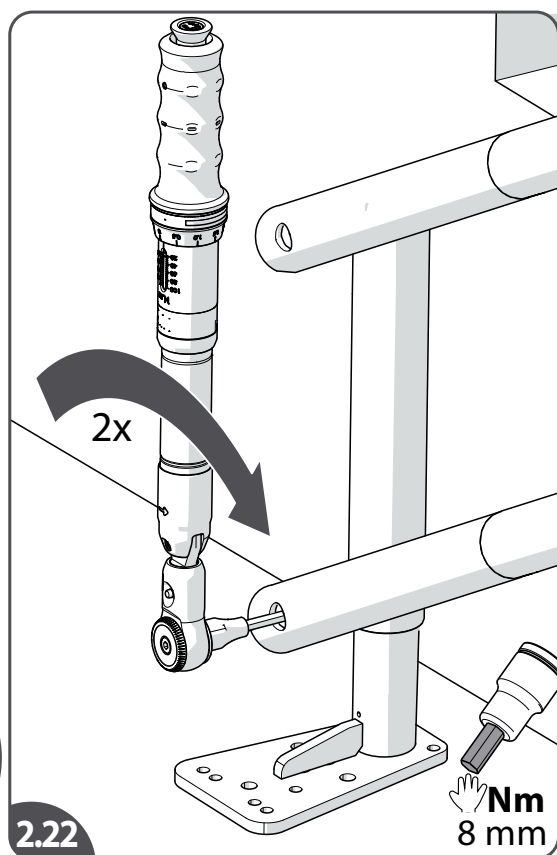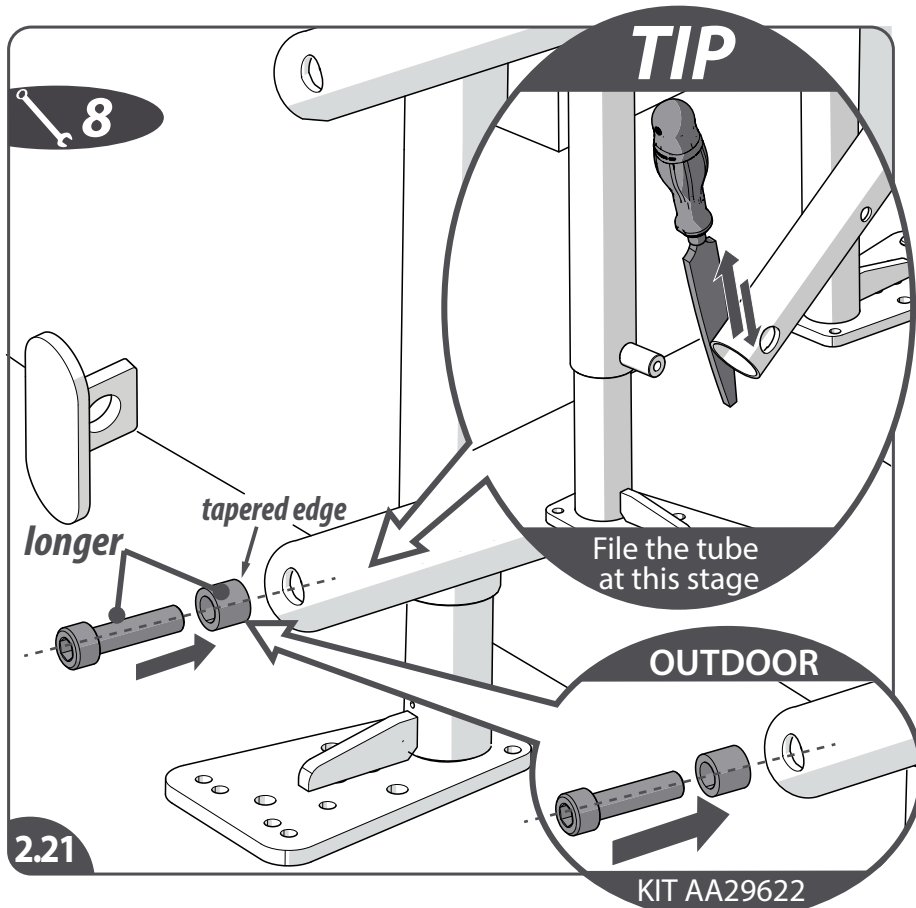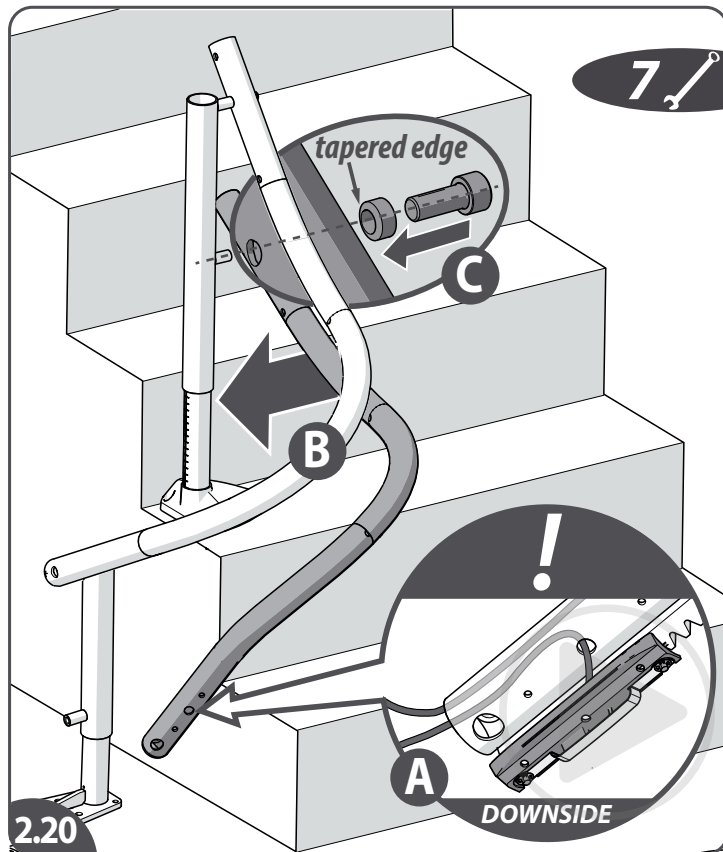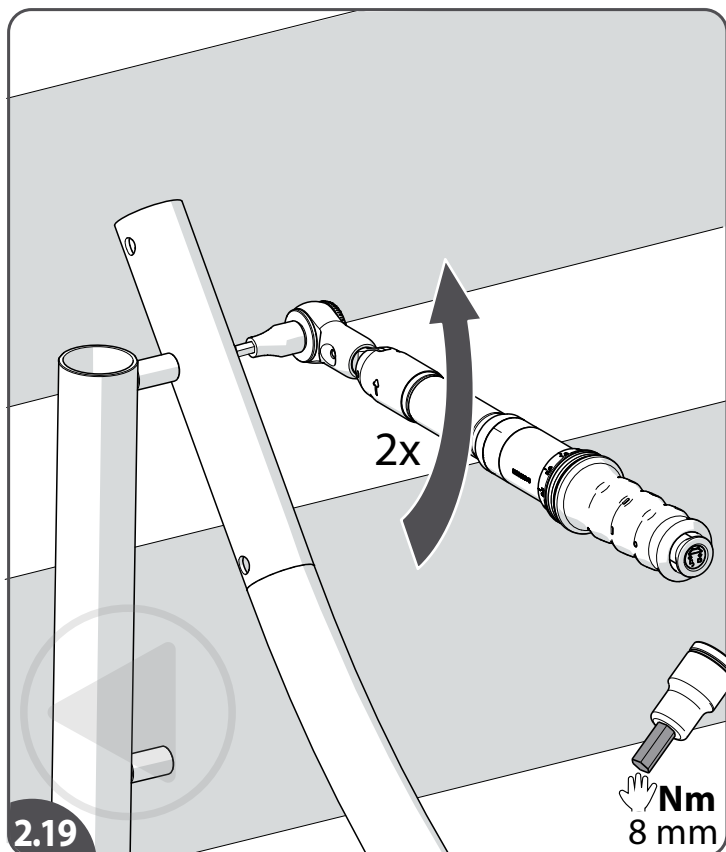[2.19] Tighten the stanchion-rail bolts [hand tight]. **After the complete track installation the torque of the stanchion-rail bolt need to be set on 45 Nm[398 Lb-In].**

[2.20] Install the bottom rail of the first section. Note the mounting holes for the charge contacts; all charge contacts are mounted in the bottom rail.

[2.21] Use the bolts+spacers to install the bottom rail. Use a longer spacer where an endstop is used. In case of an OUTDOOR installation, use kit AA29622 to mount the stanchions.

[2.22] Tighten the stanchion-rail bolts [hand tight].

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


TRACK ASSEMBLY

INTERMEDIATE STANCHION

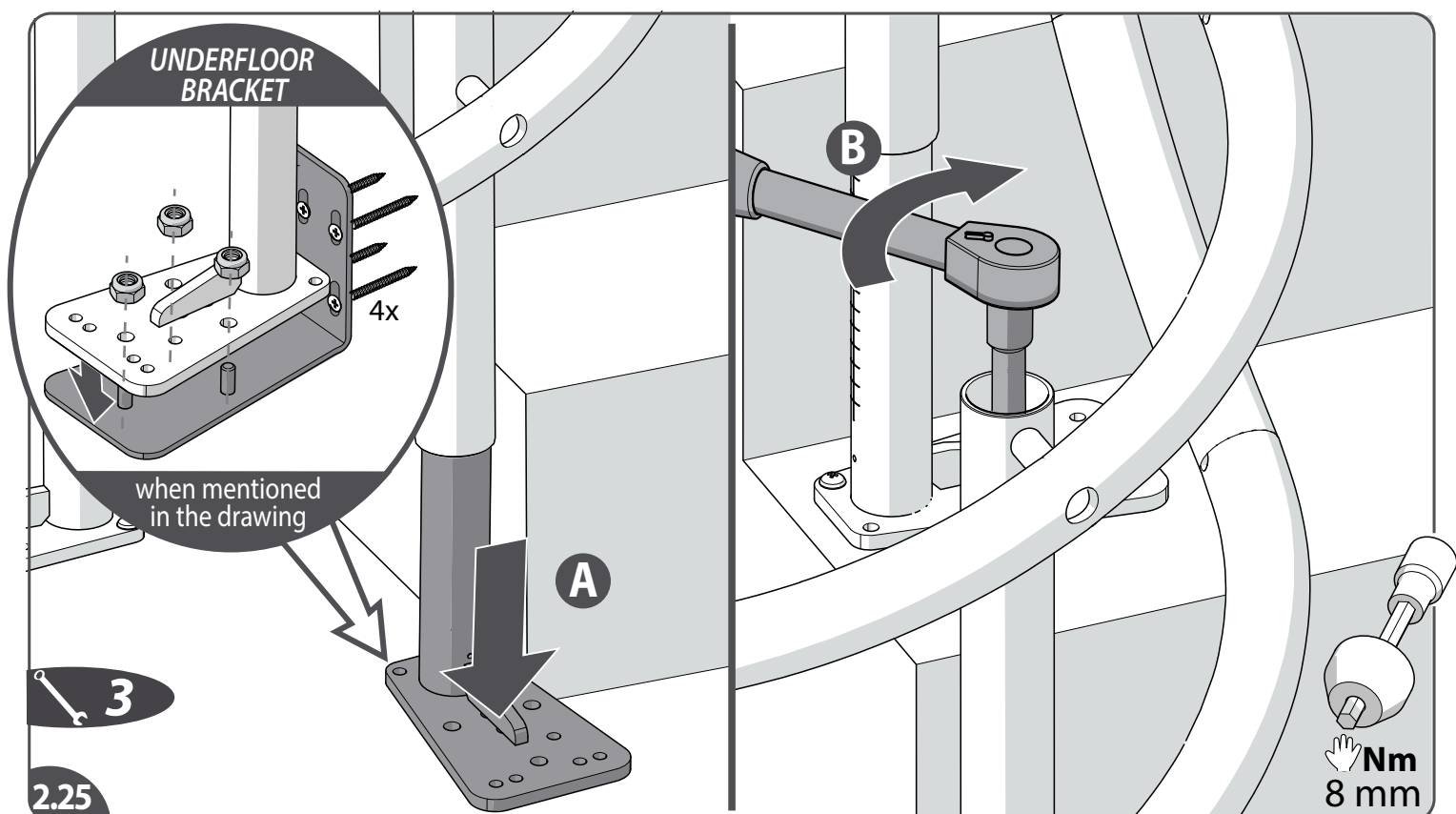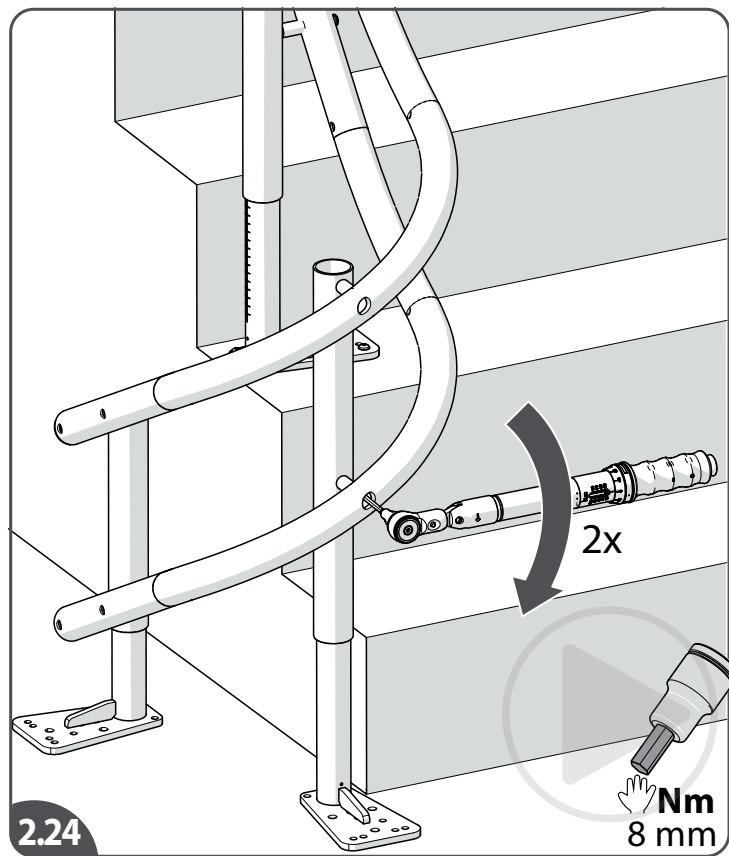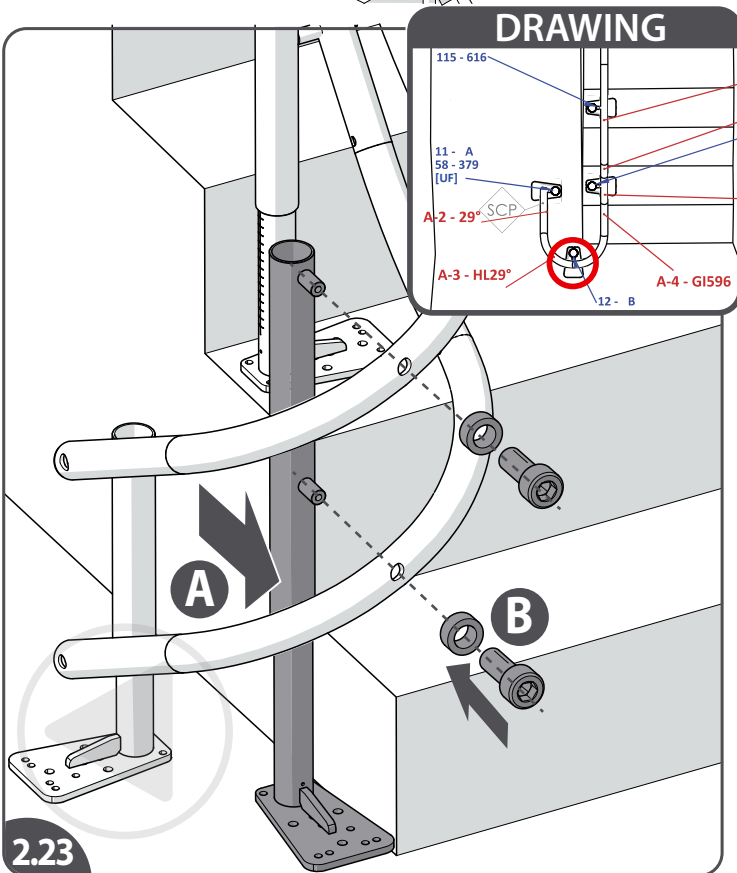[2.23] Install more stanchions of first section to get the track enough stable. **The first and last stanchion of each section must be set to the height stated on the drawing.** Specially in case of a stone/marble staircase it is important to place all stanchions of the first section to get the track enough stable.

[2.24] Tighten the stanchion-rail bolts [hand tight].

[2.25] Before fastening the stanchion lock bolt of a stanchion{B}, take away the space underneath the foot{A}. When mentioned in the drawing use **underfloor brackets**.

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

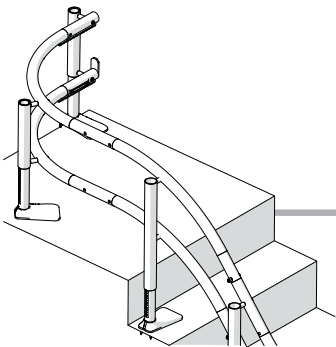
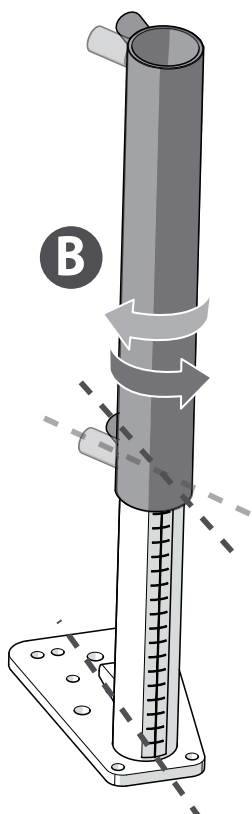
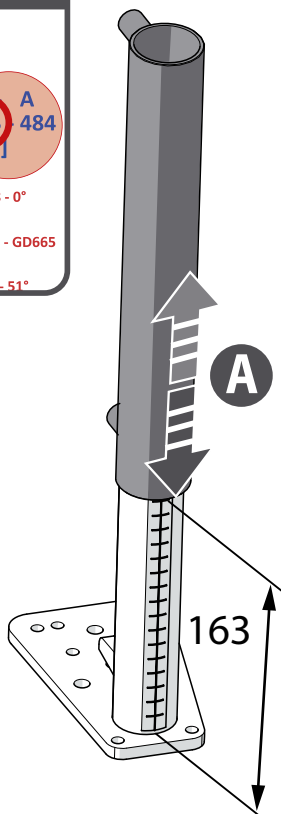
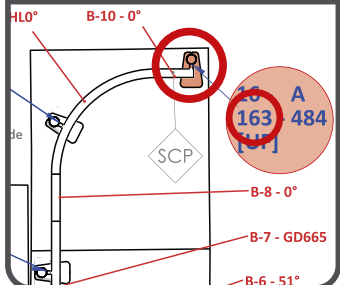
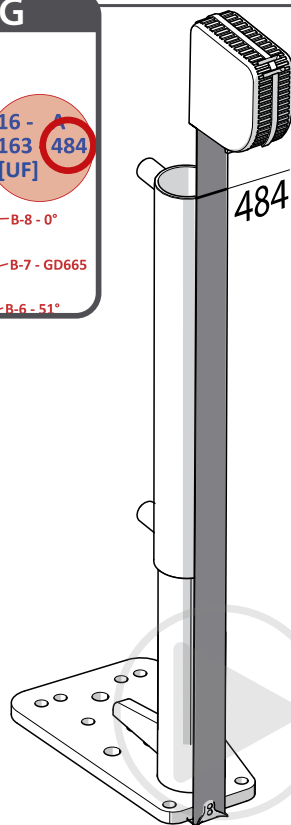
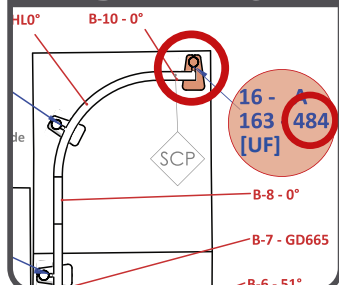
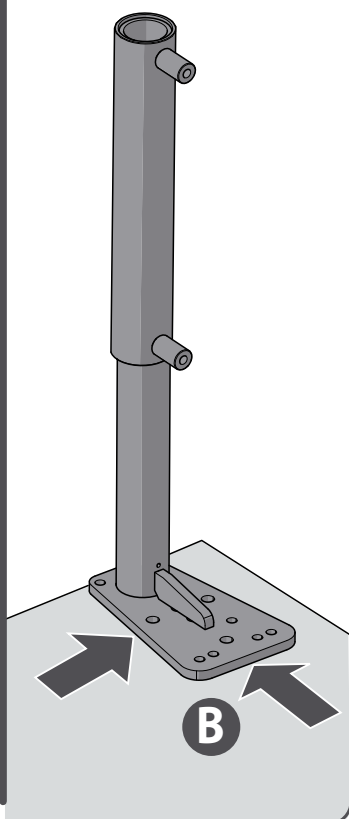
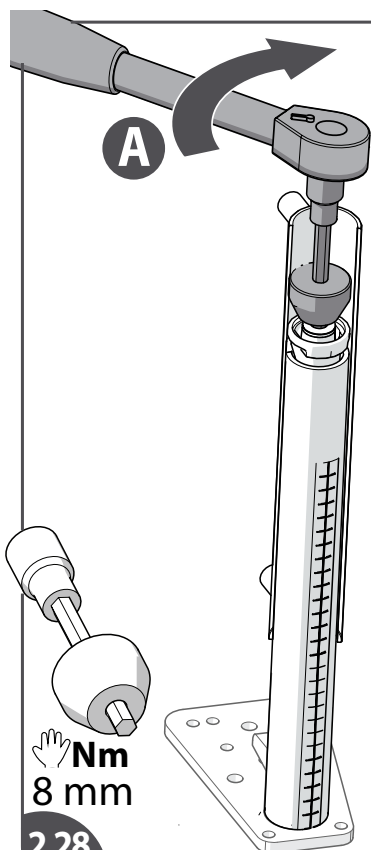
DRAWING

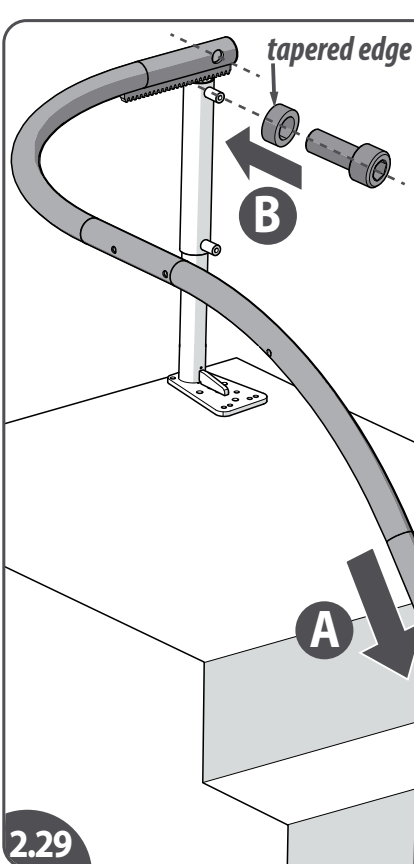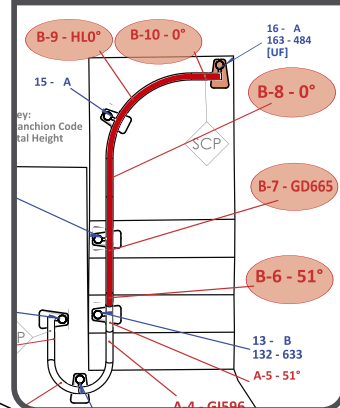2.29

- [2.26] Set the last stanchion of the second section: *set the height of this stanchion and set the lugs in line with the foot.*
 [2.27] Always check the total length of this stanchion.
 [2.28] Tighten the lock bolt [hand tight]{A}. Position this stanchion in accordance with the track offsets stated on the drawing{B}.
 [2.29] Install the top rail of the second section. For more stability additional stanchions in the rail section can be installed.

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



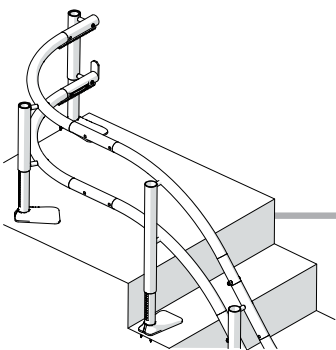
TRACK ASSEMBLY

TOP RAIL - SECOND SECTION

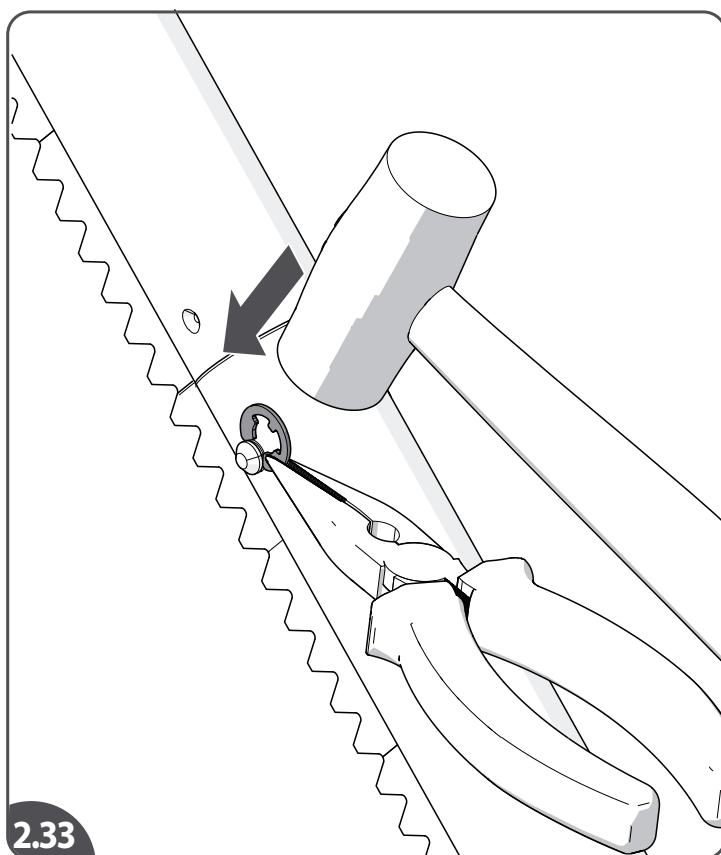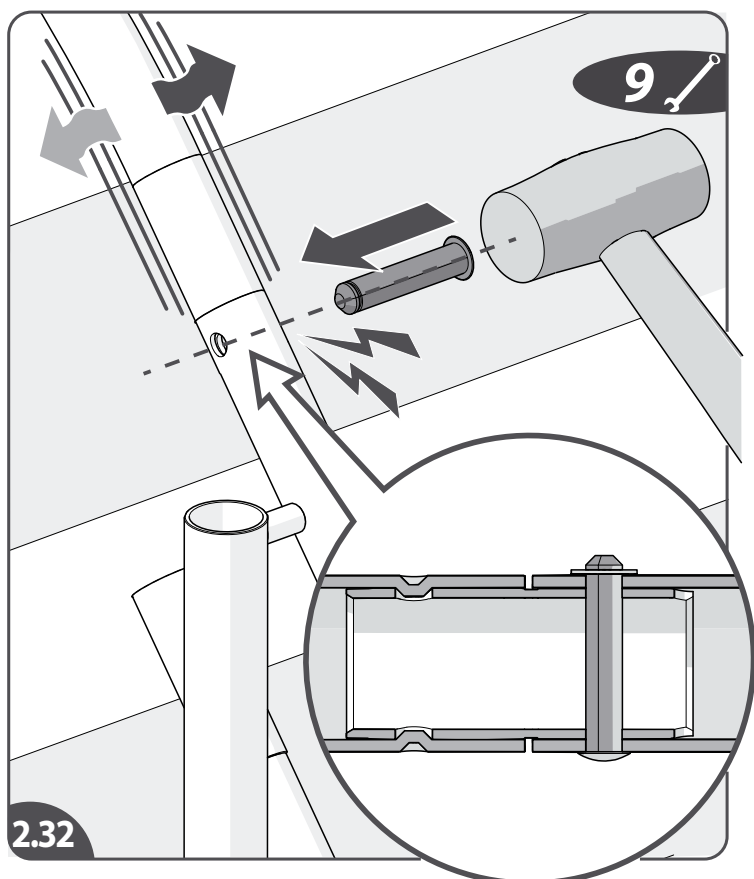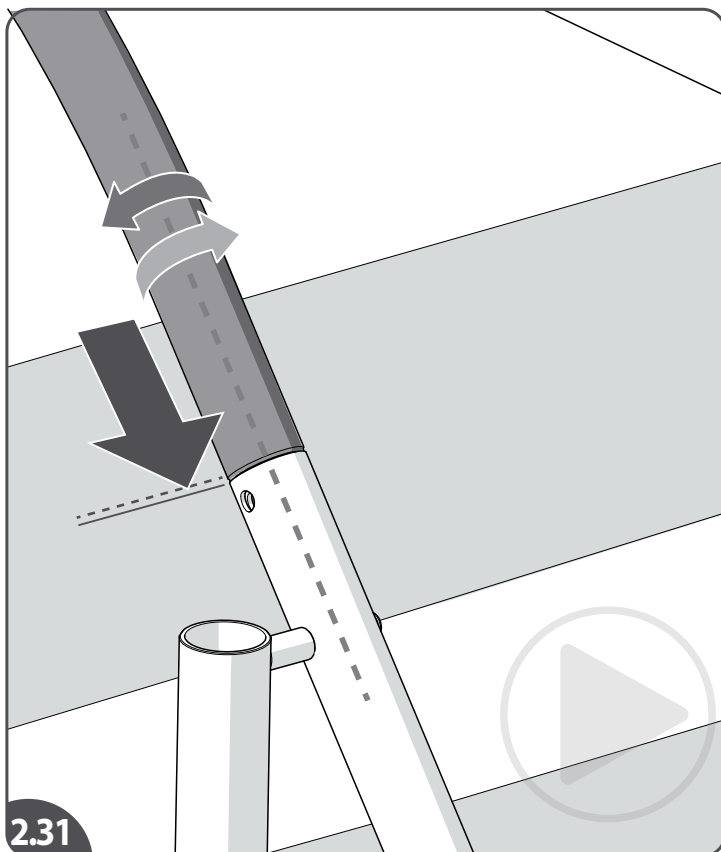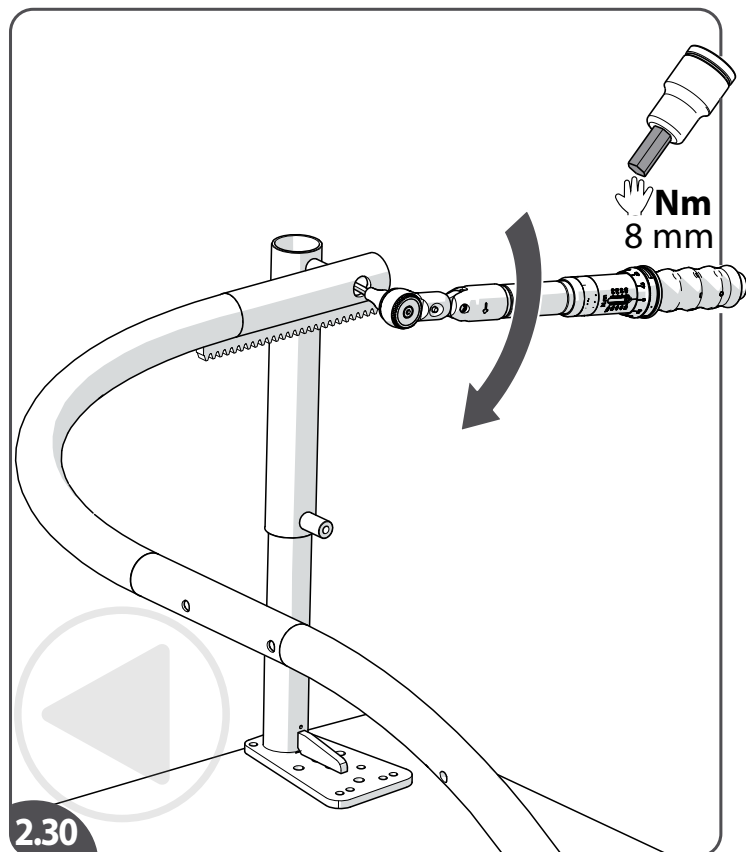- [2.30] Tighten the stanchion bolts [hand tight].
- [2.31] Carefully close the two rail sections{A}.
- [2.32] Carefully attach the rail joint pin. Moving the upper section a little bit can be helpful.
- [2.33] Install the E-clip.

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



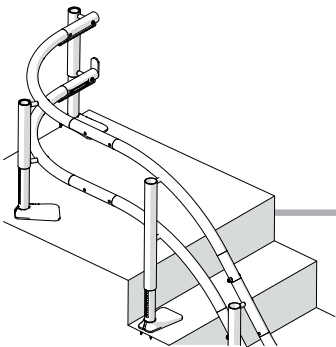
TRACK ASSEMBLY

BOTTOM RAIL - SECOND SECTION

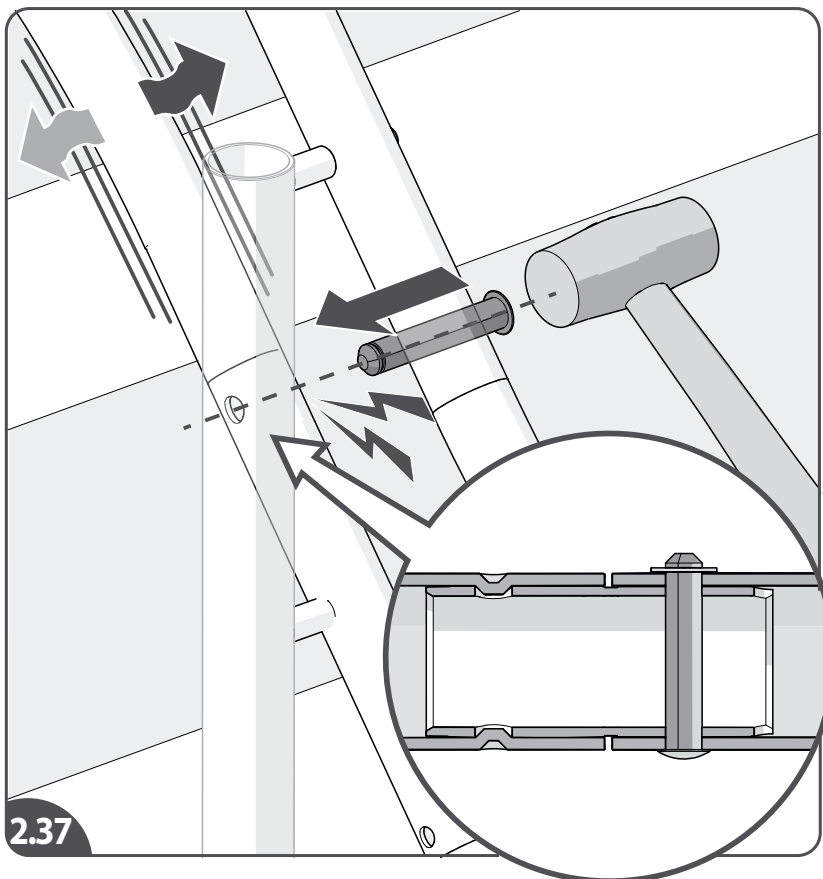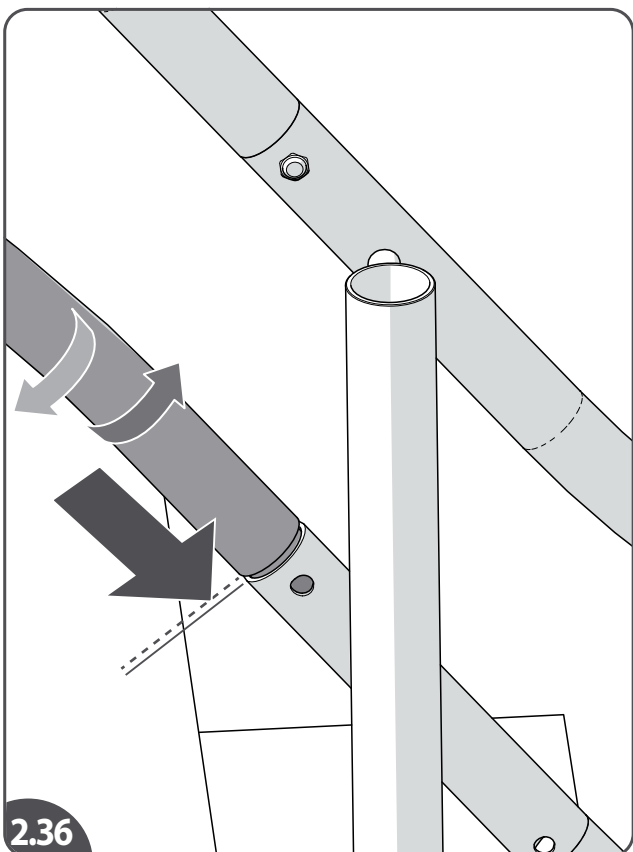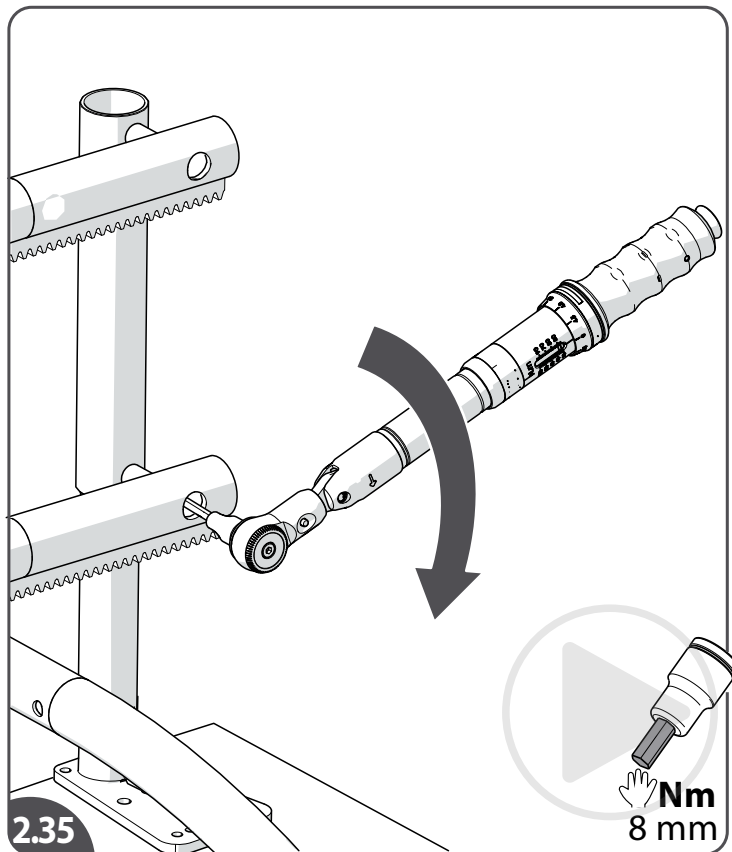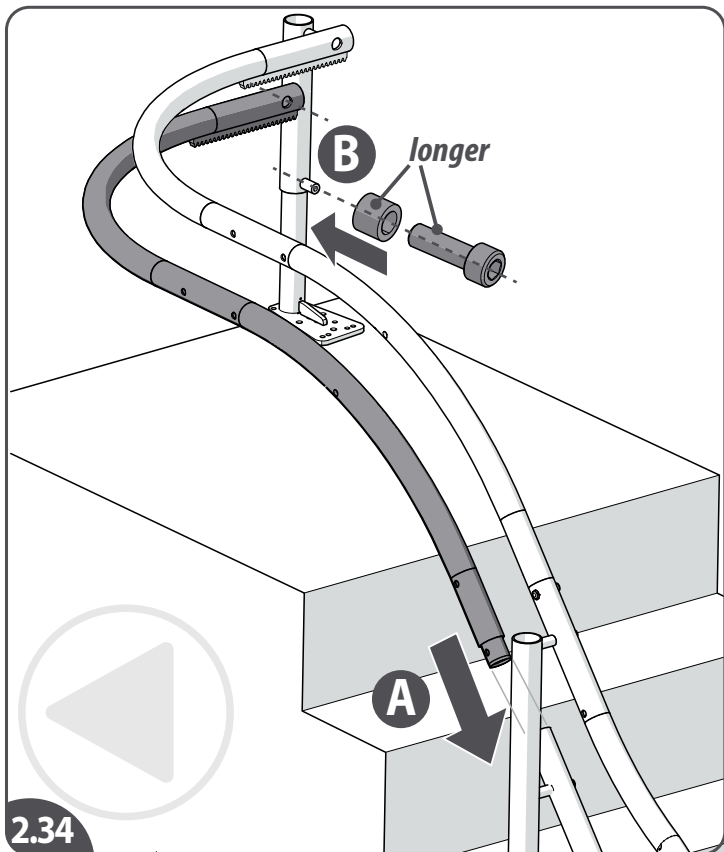- [2.34] Attach the bottom rail of the second section.
- [2.35] Fasten this rail section to the stanchion [hand tight].
- [2.36] Carefully close the two rail sections{A}. Use a railclamp (art.nr.101001.54002) to align the toothrack{B}.
- [2.37] Carefully attach the rail joint pin. Moving the upper section a little bit can be helpfull.

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



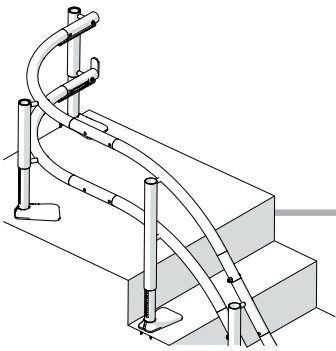
REPEAT PROCEDURE

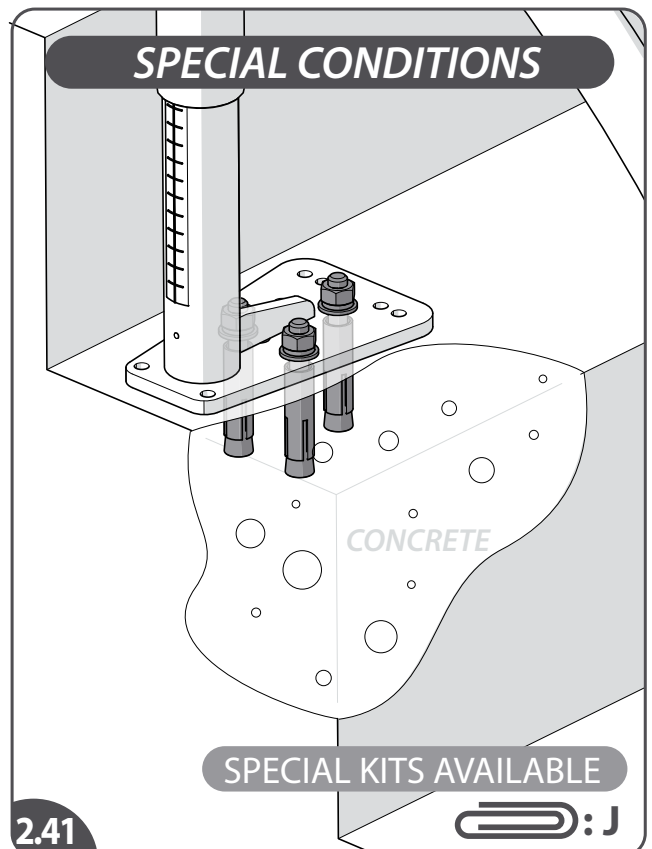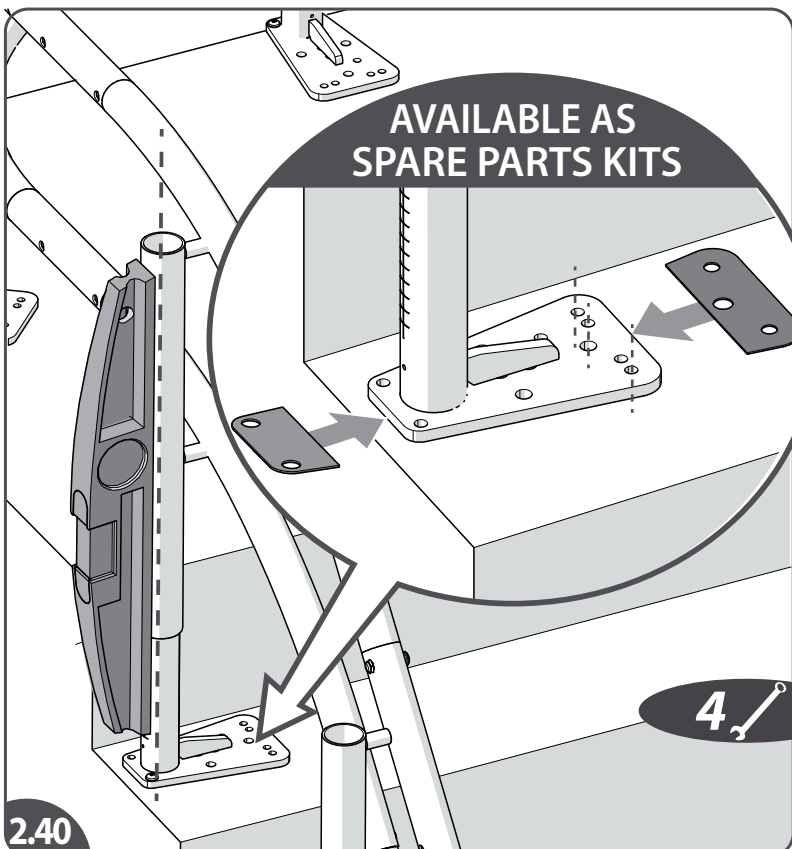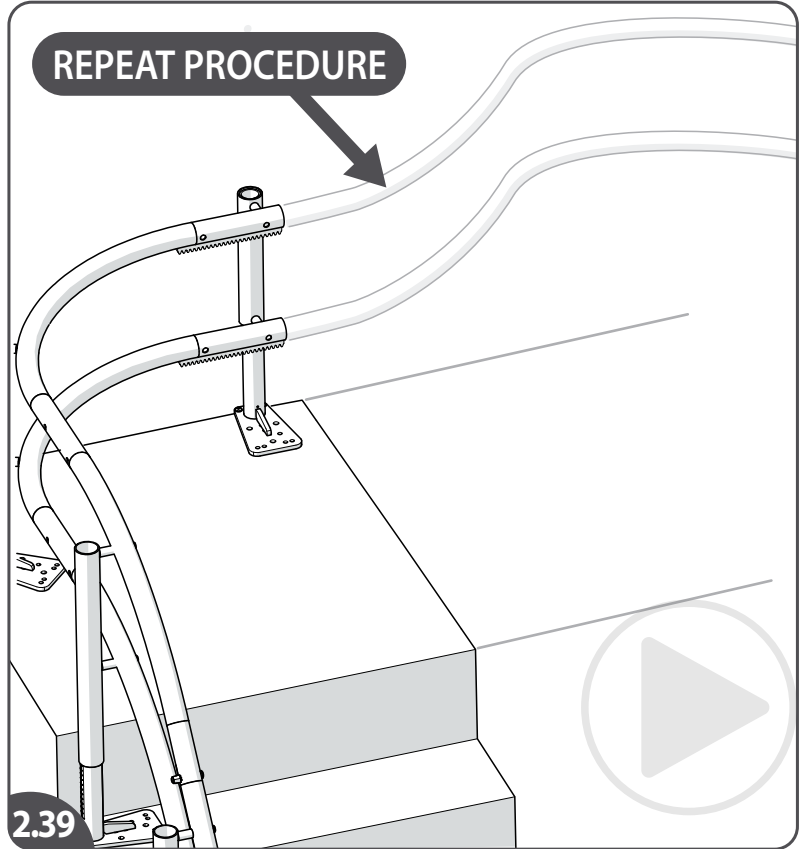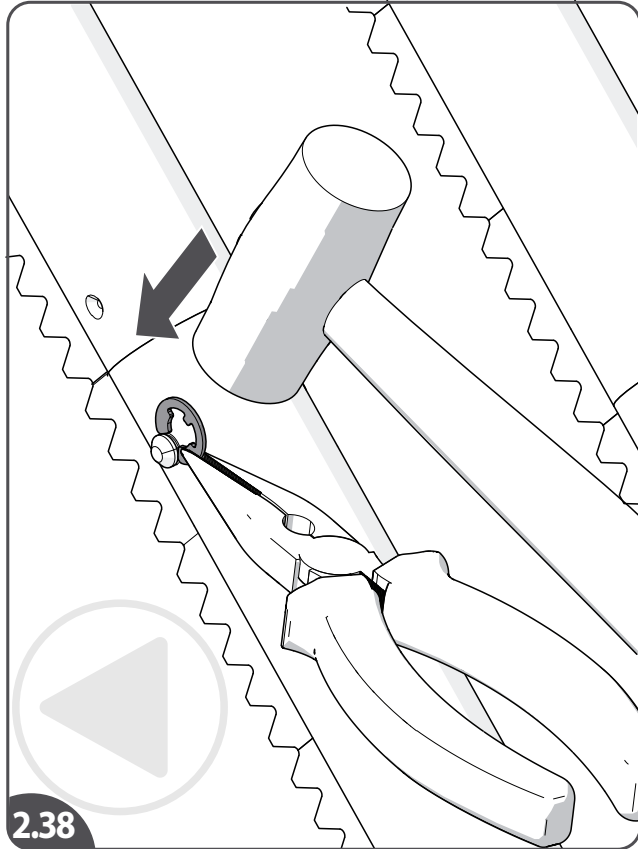[2.38] Install the E-clip.

[2.39] If there are more sections: repeat procedure for the next sections.

[2.40] Check if the stanchions are vertical and flat on the steps. If needed: add the stanchion spacers. These are available as spare parts and can be used as spacers under the front (front spacer kit, supplied in packs of 30 as **AA25740**) or rear (rear spacer kit, supplied in packs of 30 as **AA25735**) of the stanchion to help make it level – e.g. if a step is worn or has a very thick carpet.

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


TRACK ASSEMBLY.

LEVEL STANCHIONS

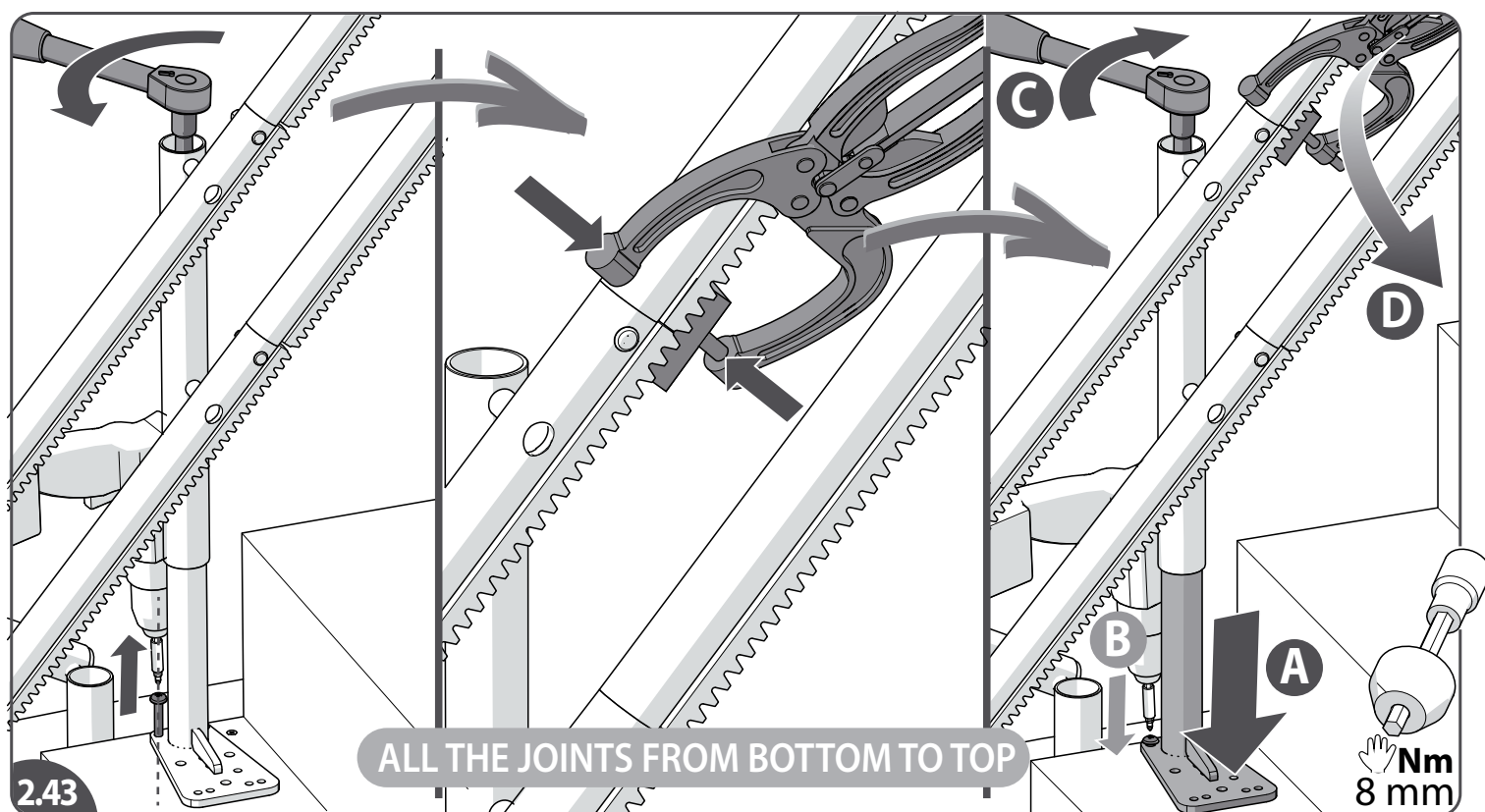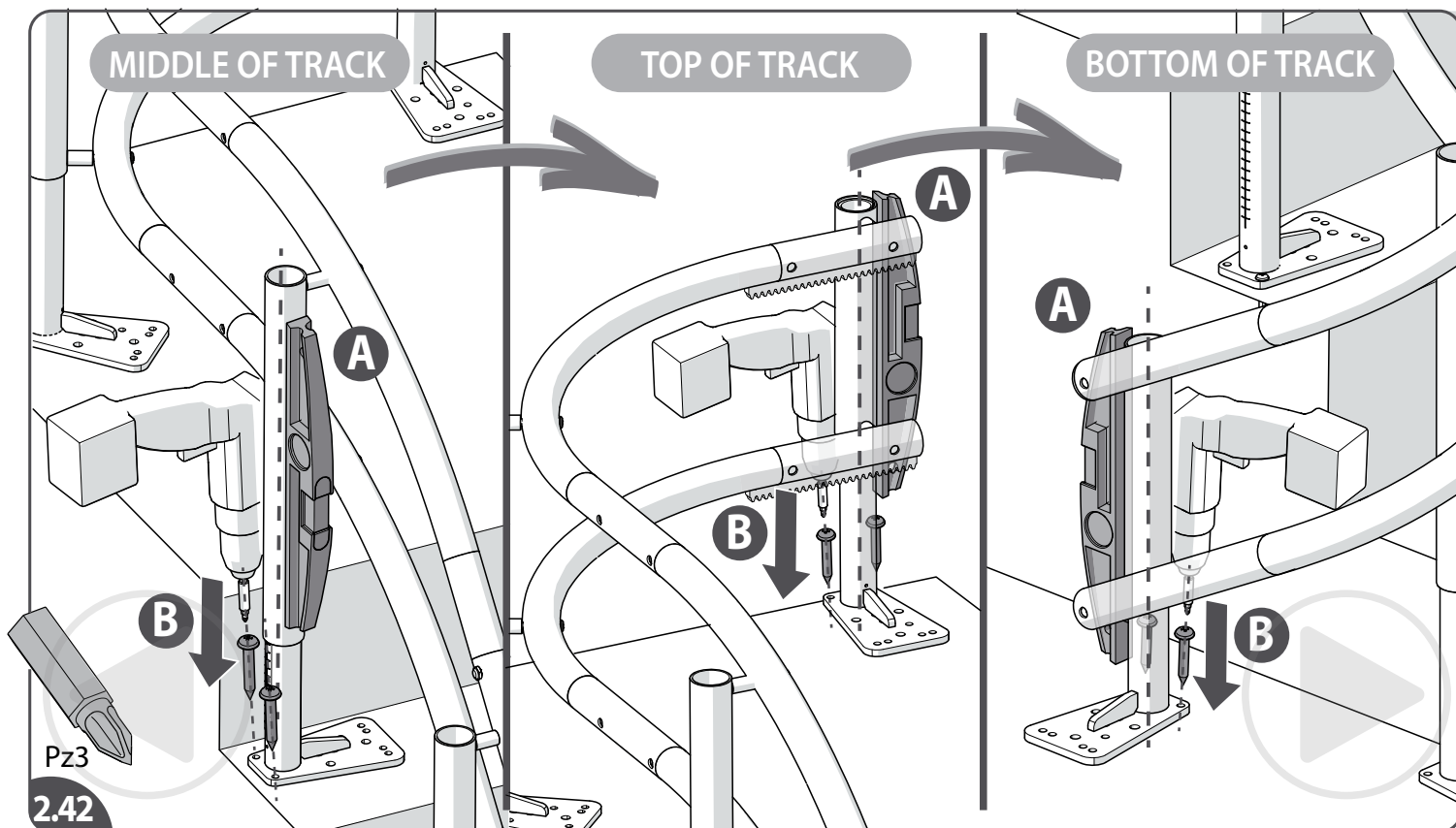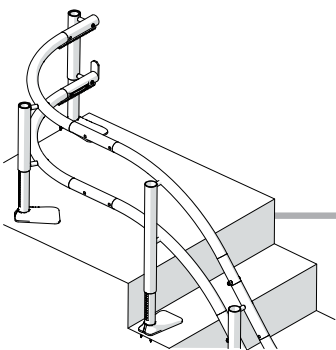[2.42]Align the stanchions vertical and fasten the feet with two screws (both sides). First the middle stanchion of the track, then the top stanchion and finally the bottom stanchion.

Note: in case of a marble staircase fix the stanchions at the end of the installation (chapter 6). Still take care that the lift is stable enough when making the test rides.

[2.43]Check all the joints from the bottom to the top: *loosen the screws and the lock bolt. Use a railclamp (art.nr.101001.54002) to align the toothrack. This is essential for a good joint. Take away the space underneath the foot{A} and fix it to the steps{B}. Tighten the lock bolt [hand tight]{C}.*

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
IT

3

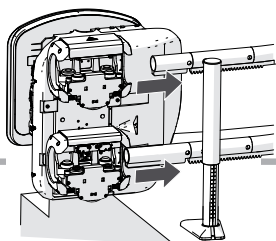
YOU MUST BE LOGGED IN TO
THE HANDICARE PARTNERSITE!

POWER PACK LOADING BARS

POWER
PACK


3.0

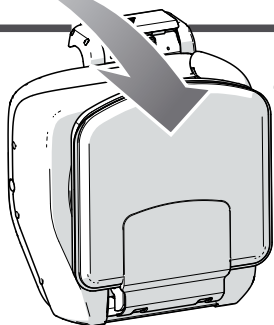
18

TIP


3.2

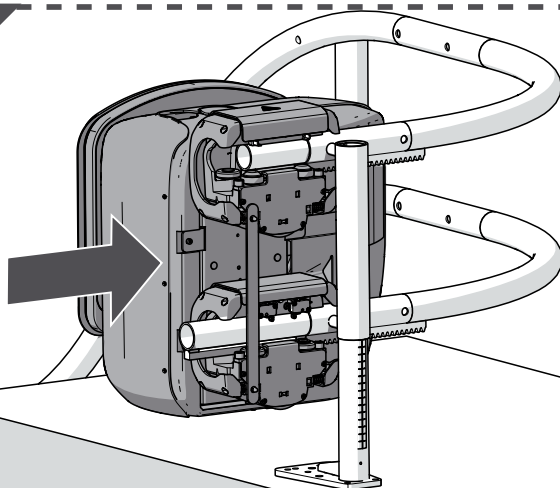
3.1

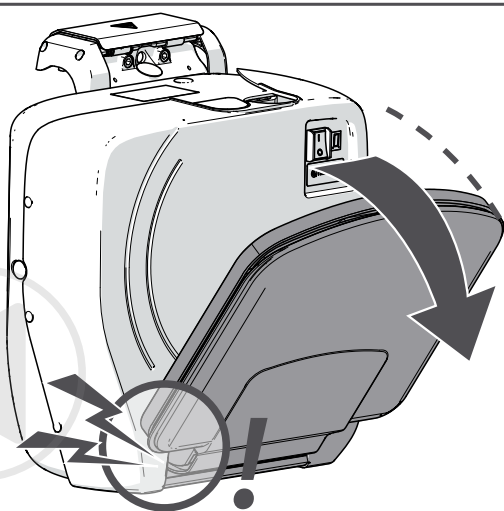
3.3

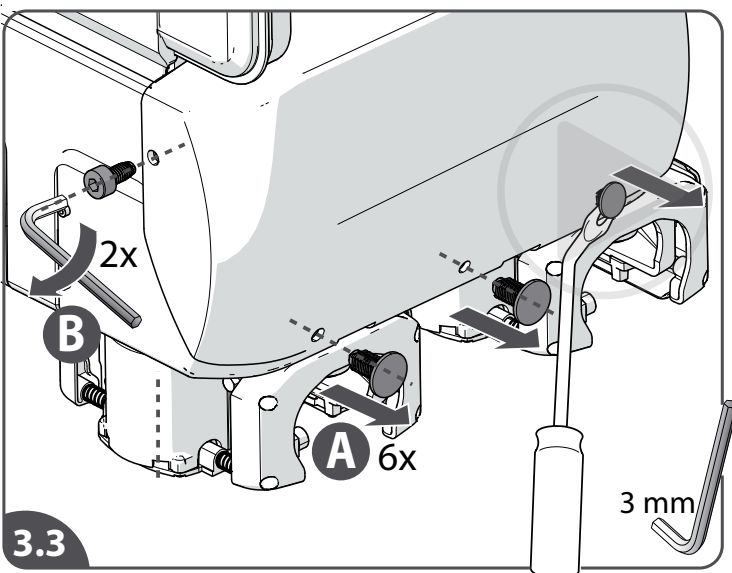3.4

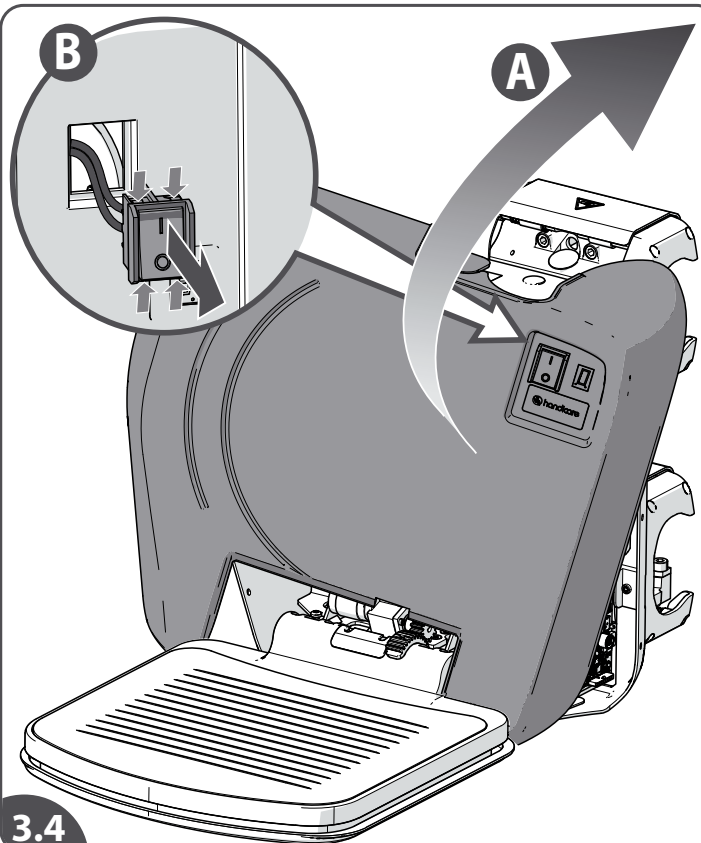3.5

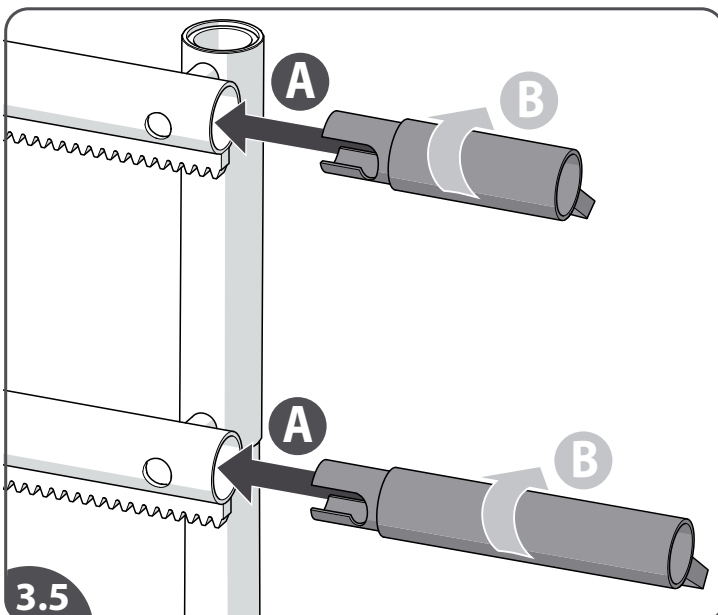- [3.0] Take the power pack out of the box.
 [3.1] Fold down the footplate. **Note:** careful, do not damage the footplate.
 [3.2] **TIP:** place the power pack, from the box, directly on the rails to avoid damage to footplate/power pack.
 [3.3] Loosen the rivets(6x){A} and screws(2x){B} and remove the front cover{C}.
 [3.4] Remove the front cover{A}. Remove the on/off-switch from the front cover and disconnect the loom of the display{B}.
 [3.5] Fit the loading bars to the top of the track: *mount the top loading bar (short).*

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




IT

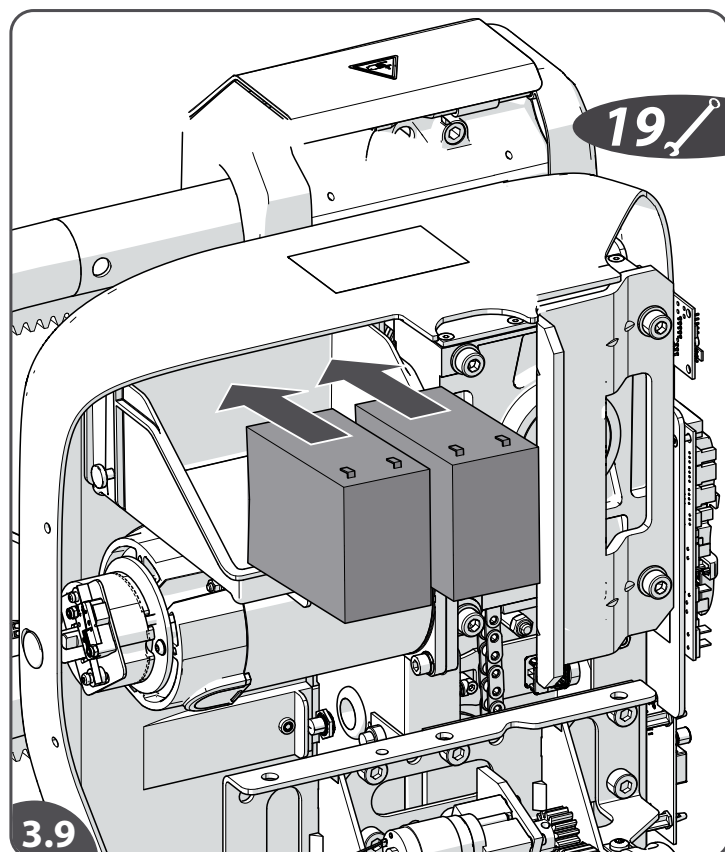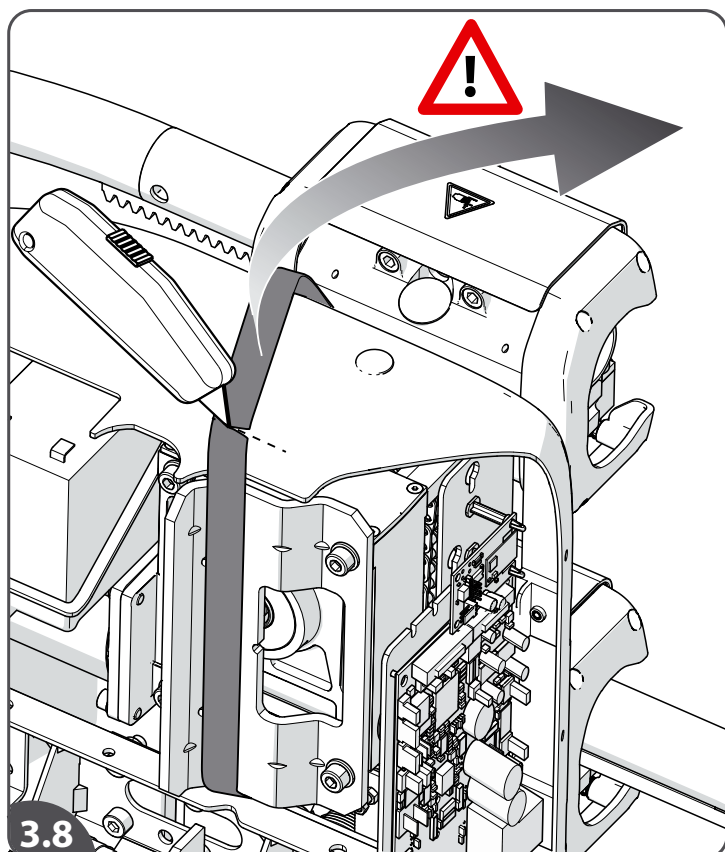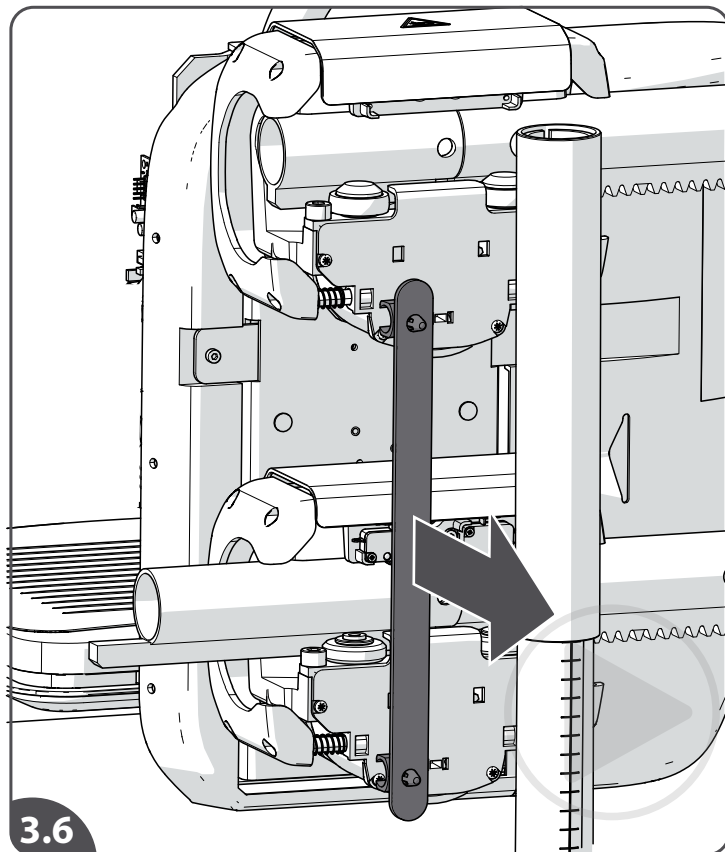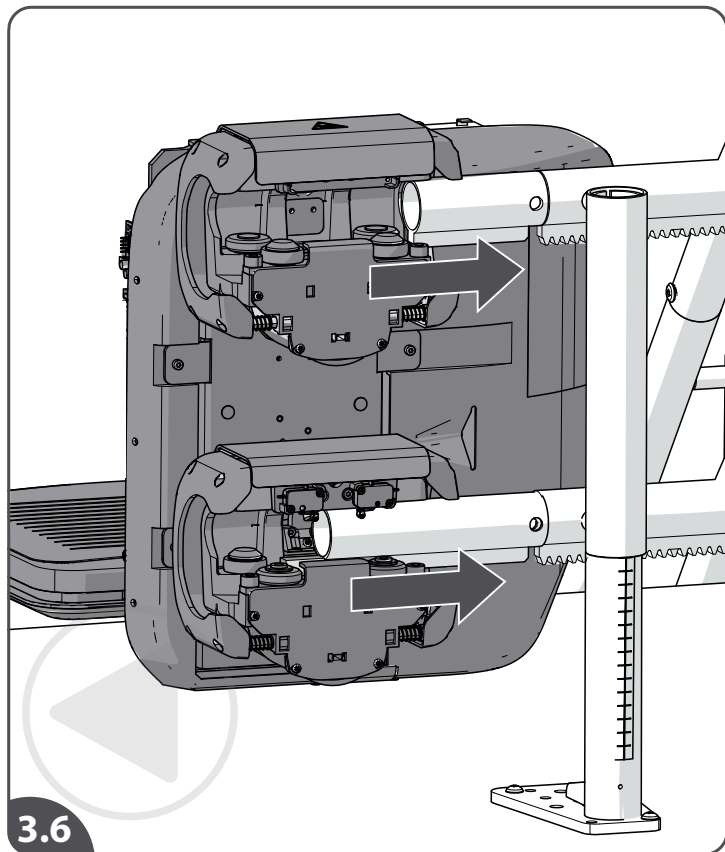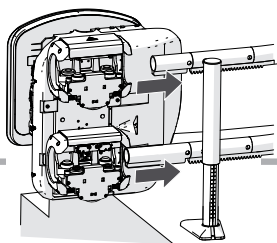- [3.6] Feed the power pack on to the loading bars.
- [3.7] Remove the trunnion loading strip.
- [3.8] Remove the carry strip !
- [3.9] Install the batteries.

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



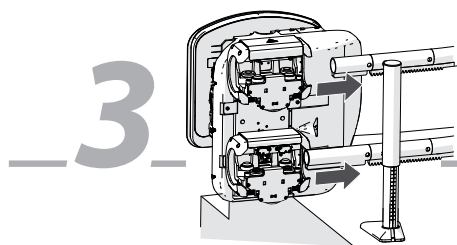
POWER PACK

PREPARING POWER PACK

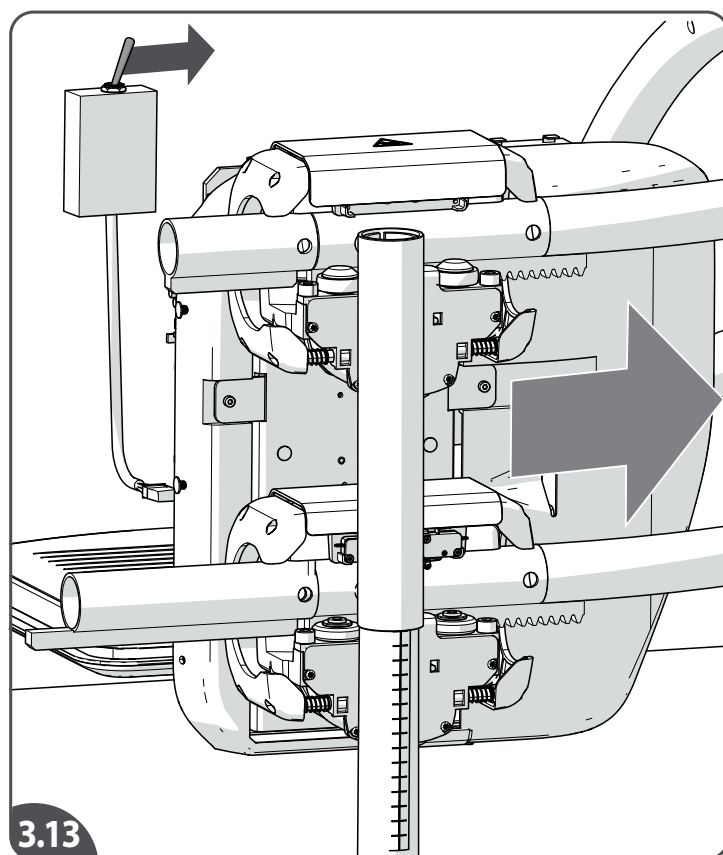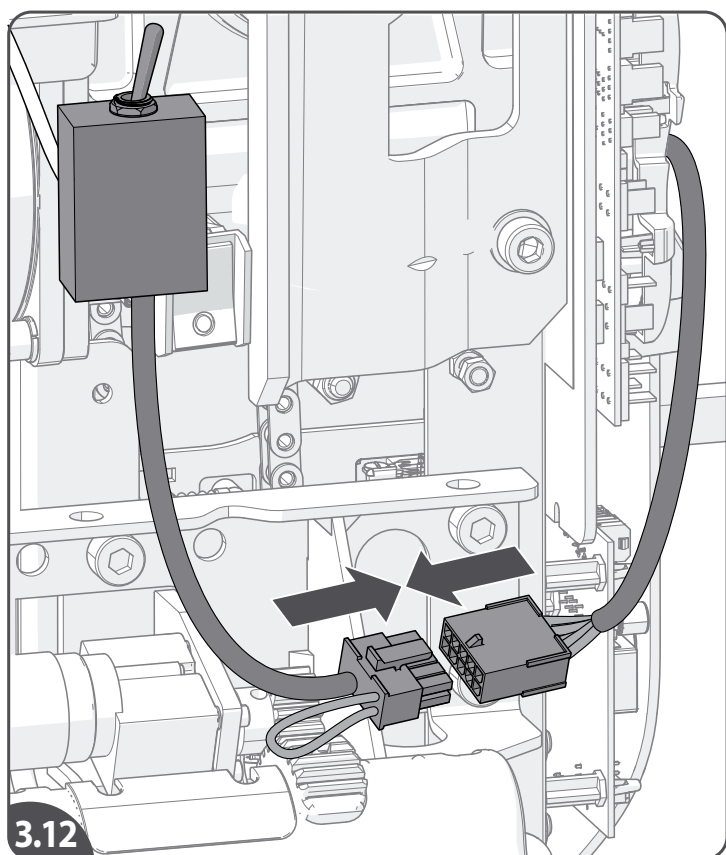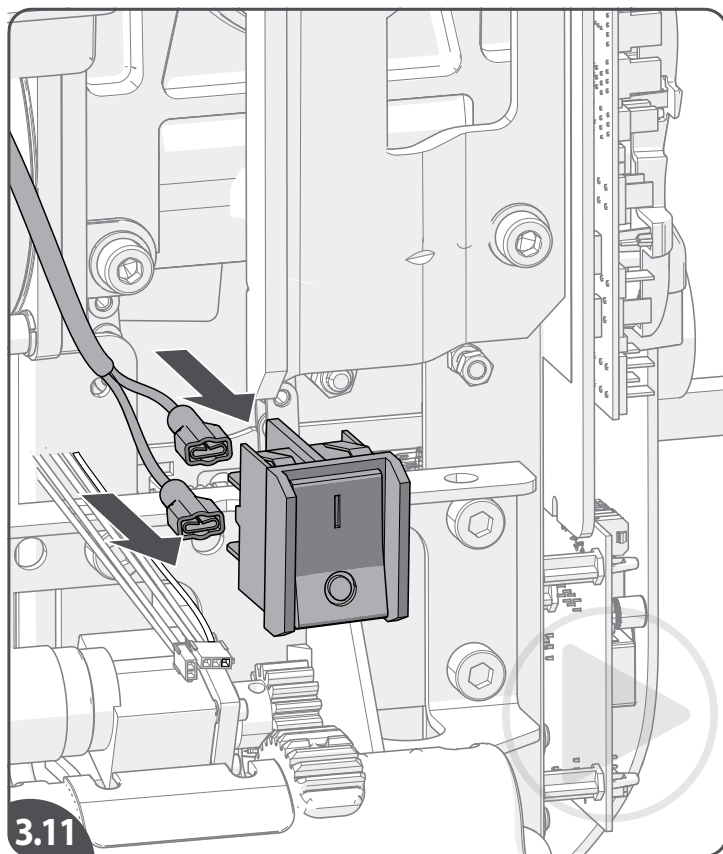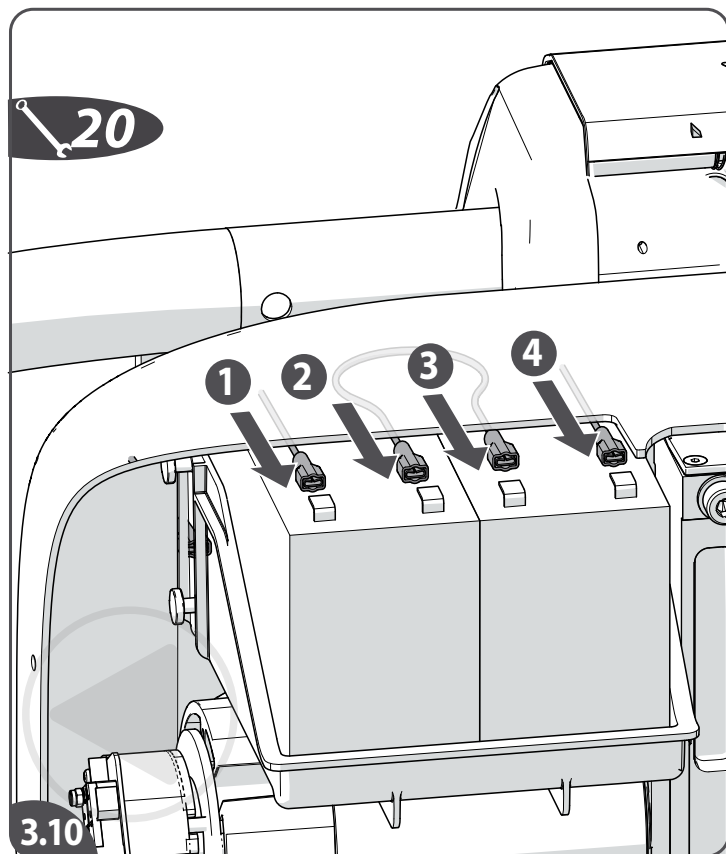- [3.10] Connect the wires of the batteries.
- [3.11] Connect the on/off-switch.
- [3.12] Connect the loading toggle.
- [3.13] Carefully drive the power pack on to the track.

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
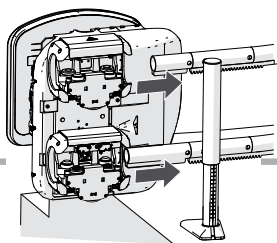
**A $\approx$ B, maximum
5 mm difference**

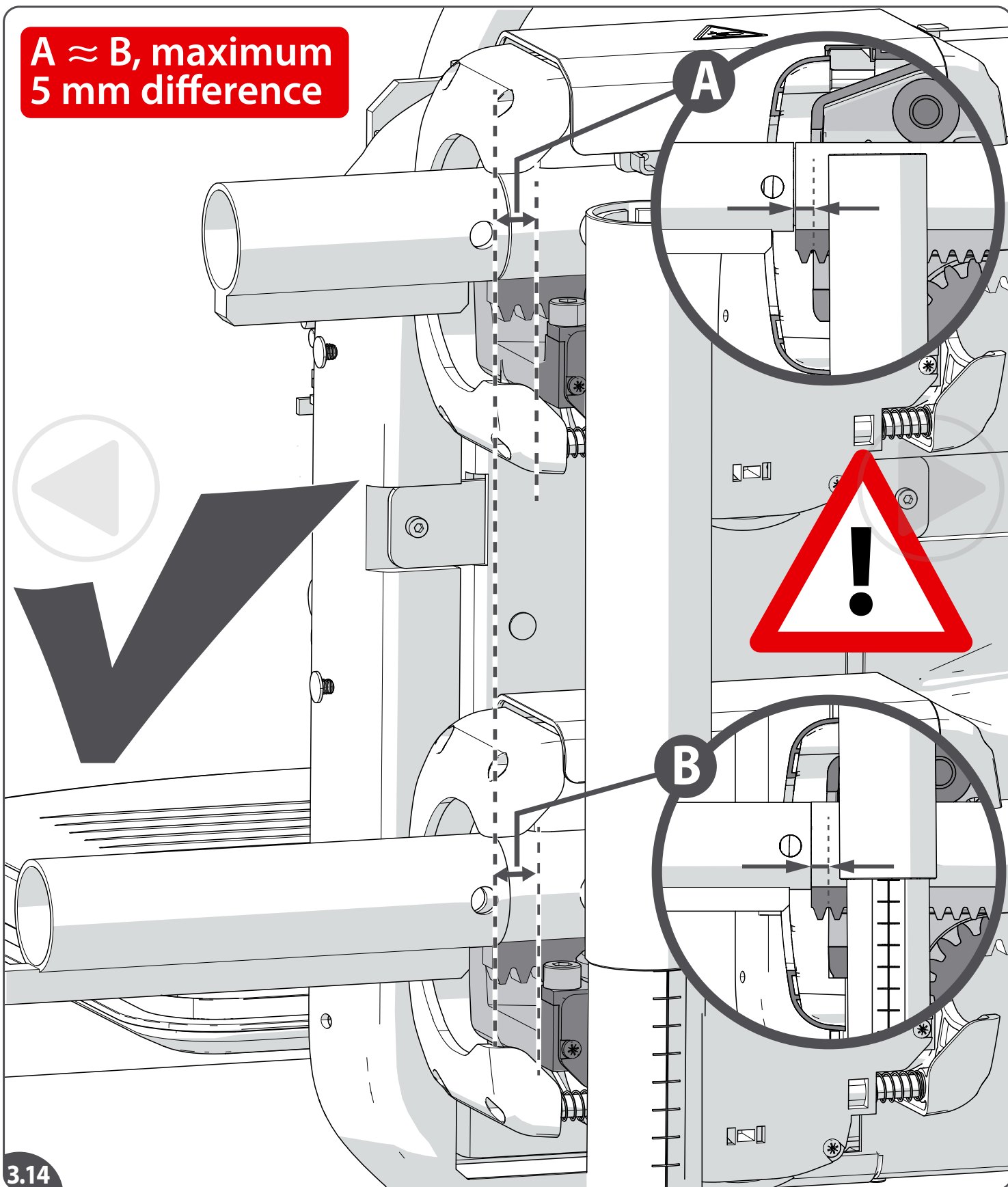3.14

[3.14] Make sure the trunnions are vertically aligned ($A \approx B$, maximum 5 mm difference). ***If the power pack is loaded without the trunnions correctly aligned, it will cause permanent damage to the power pack. If this occurs a new power pack MUST be ordered as a matter of safety.***

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



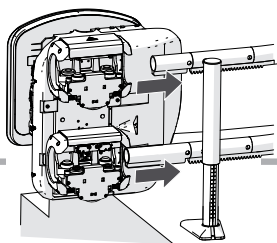
POWER PACK

REMOVE INSTALLATION PARTS

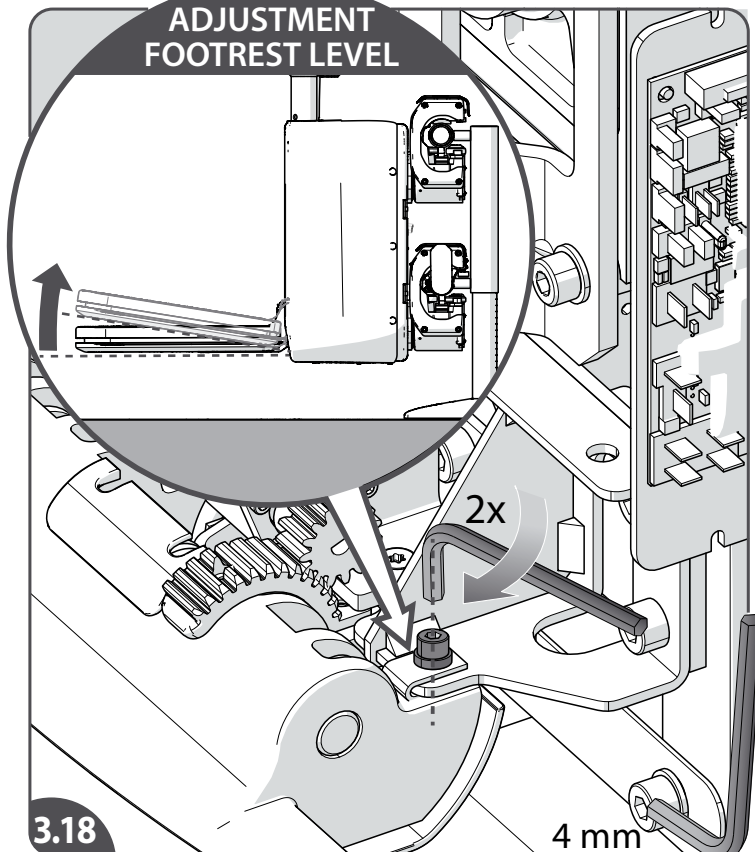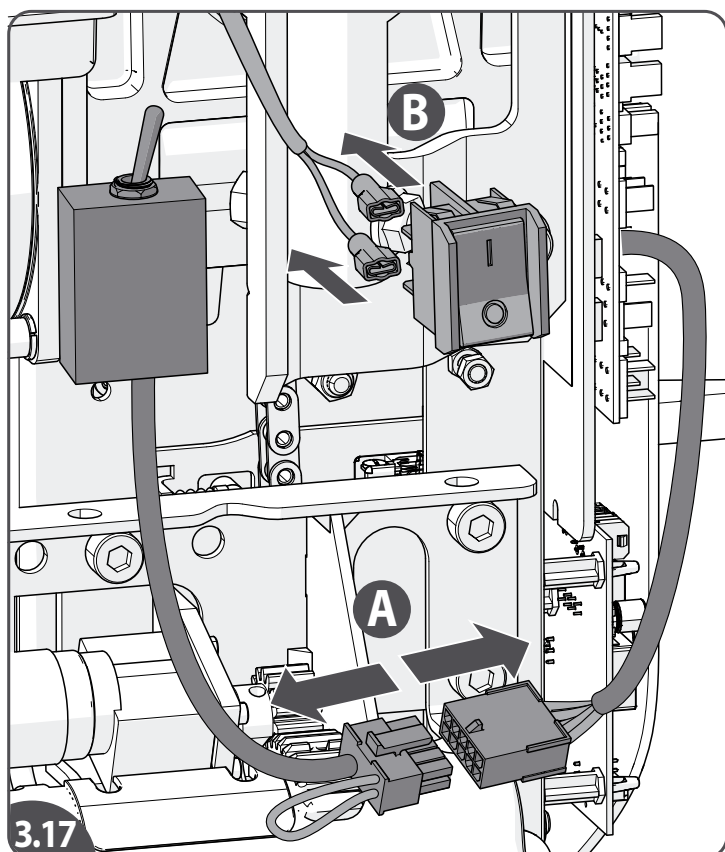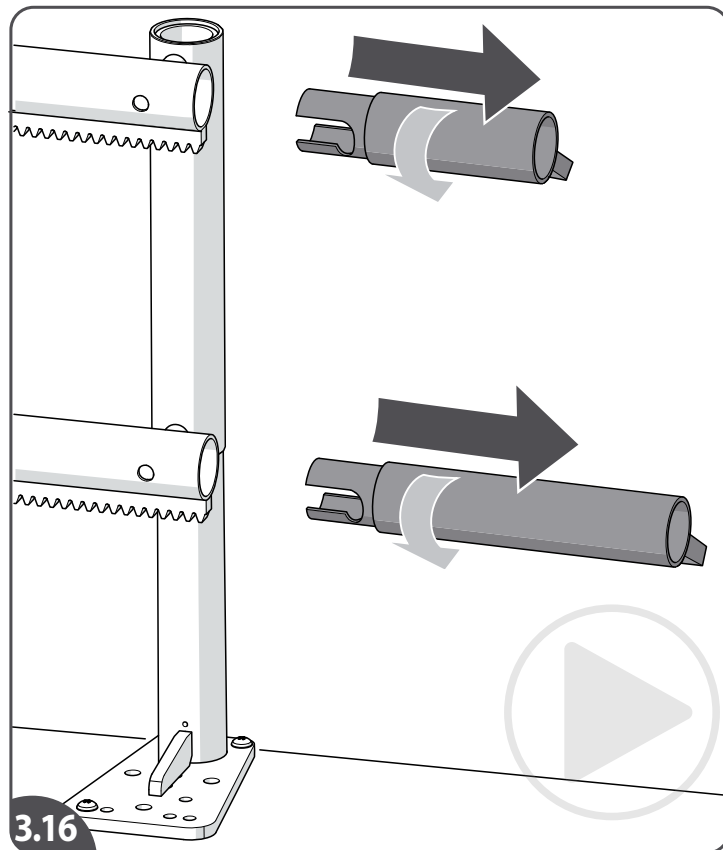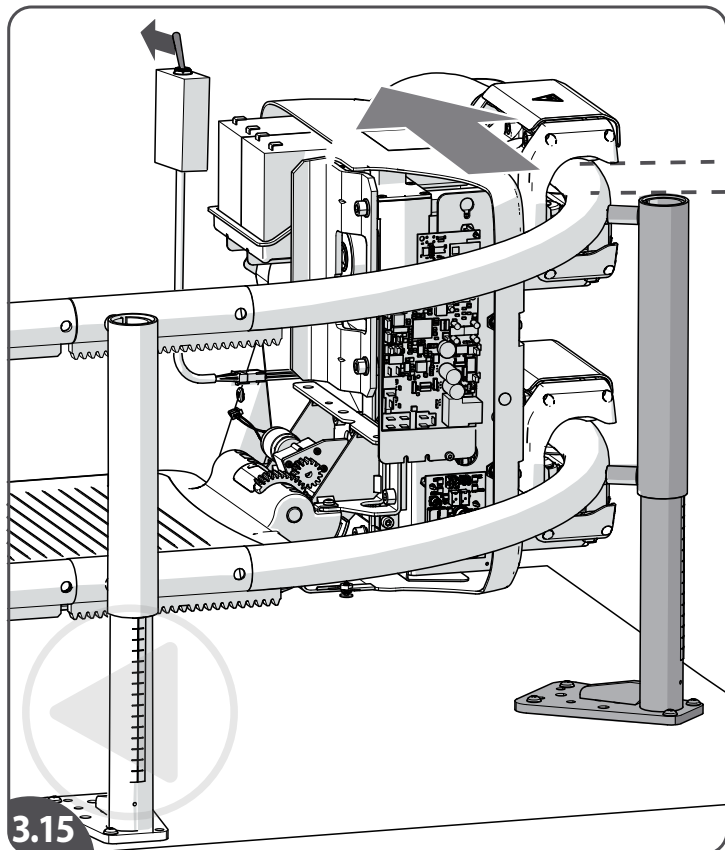- [3.15] Drive the power pack beyond the next stanchion.
- [3.16] Remove the long loading and the short loading bar.
- [3.17] Remove the dummy seat toggle looms{A}. Remove the on/off-switch{B}.
- [3.18] If needed, adjust the footrest level/angle.

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




IT

4

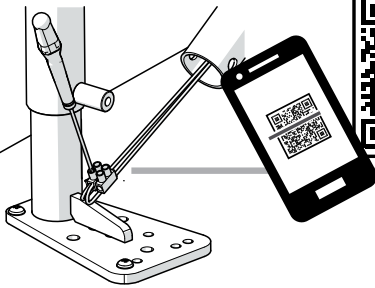YOU MUST BE LOGGED IN TO
THE HANDICARE PARTNERSITE!

TRACK ELECTRICAL CHARGE CONTACTS

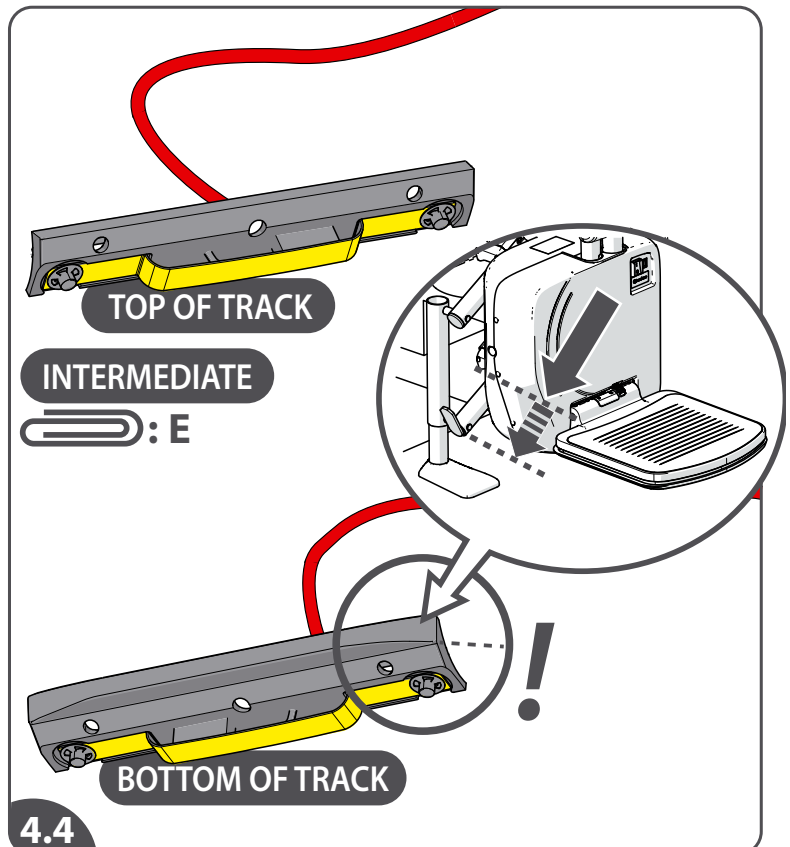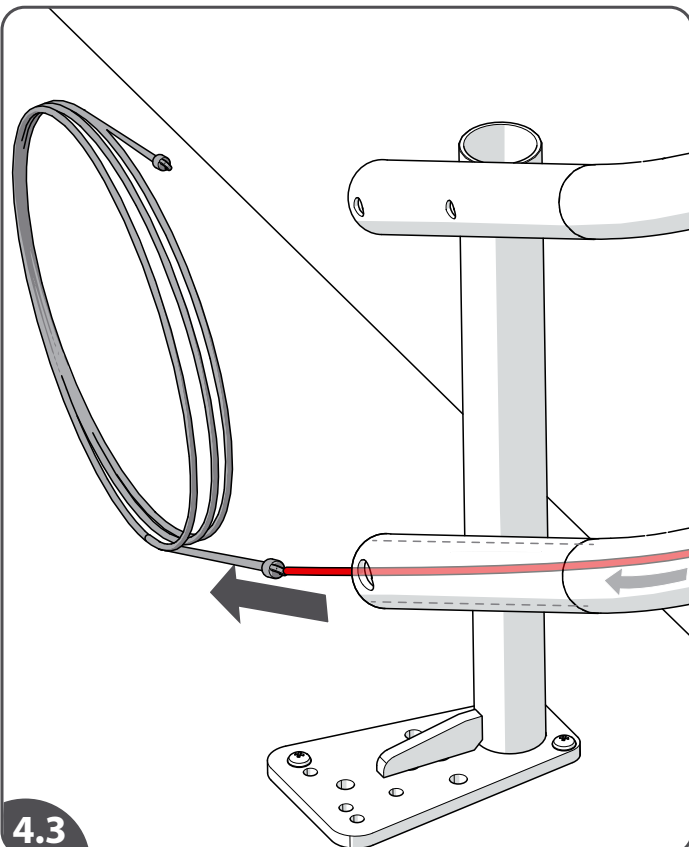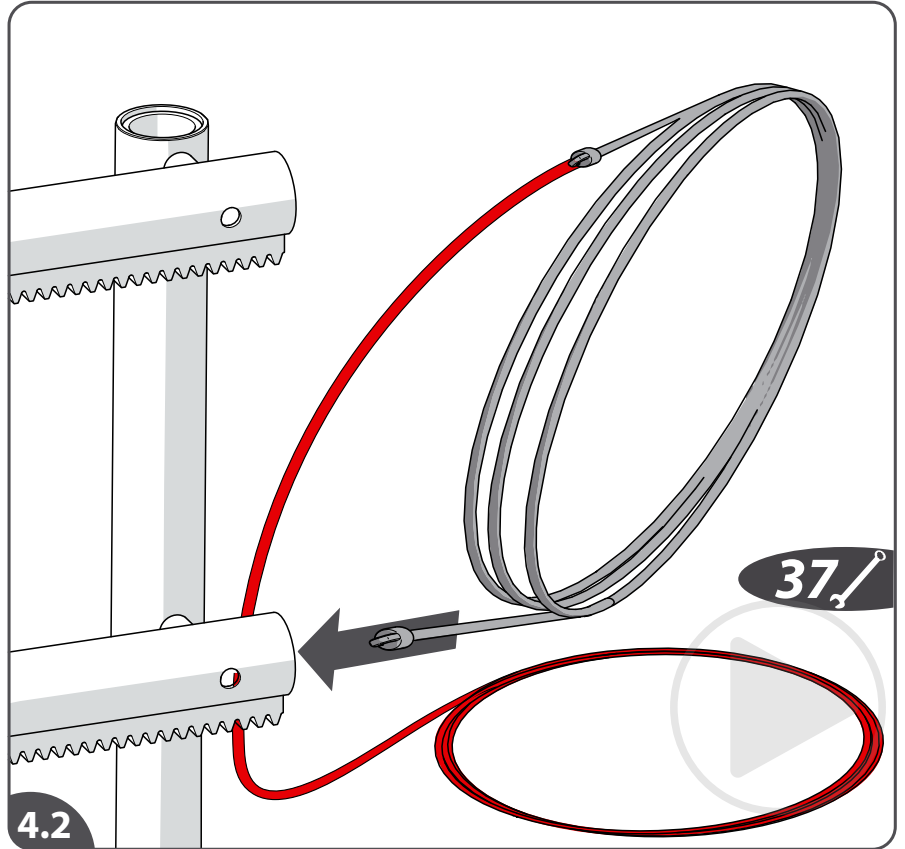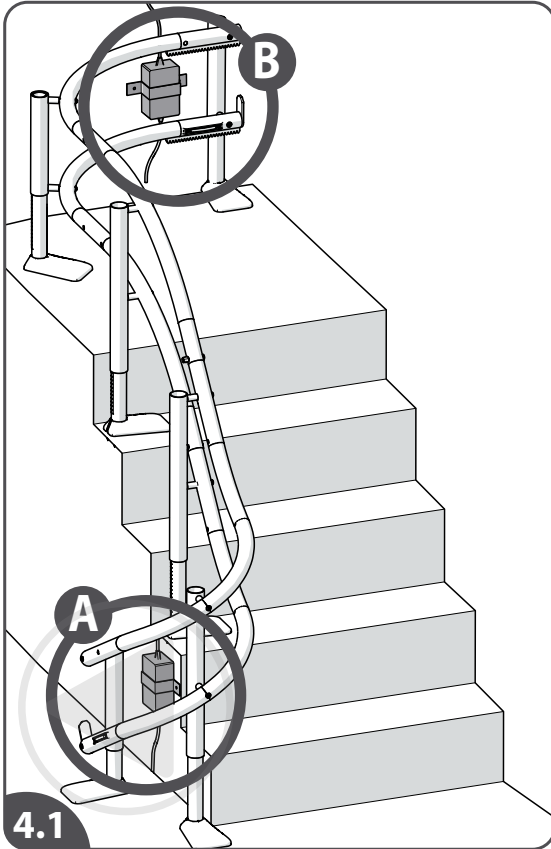- [4.1] Decide where to install the power supply: depending on application this could be at the top{B} or bottom{A} of the stairs.
- [4.2] Guide the pull through spring through the bottom rail.
- [4.3] Pull the charger wires completely through the bottom rail.
- [4.4] The upper charge contact is different from bottom charging contact: the bottom charging contact has a sloping surface to reduce the speed of the lift. See attachment **E** in case also an intermediate stop needs to be installed.

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



TRACK ELECTRICAL

CHARGE WIRES

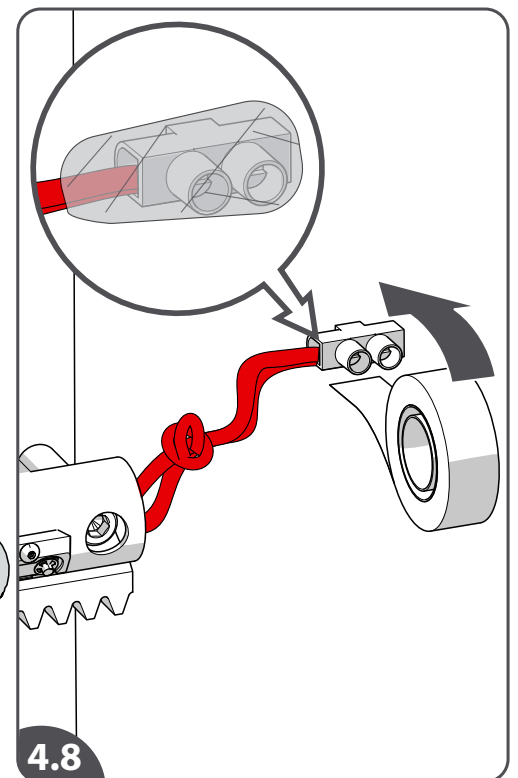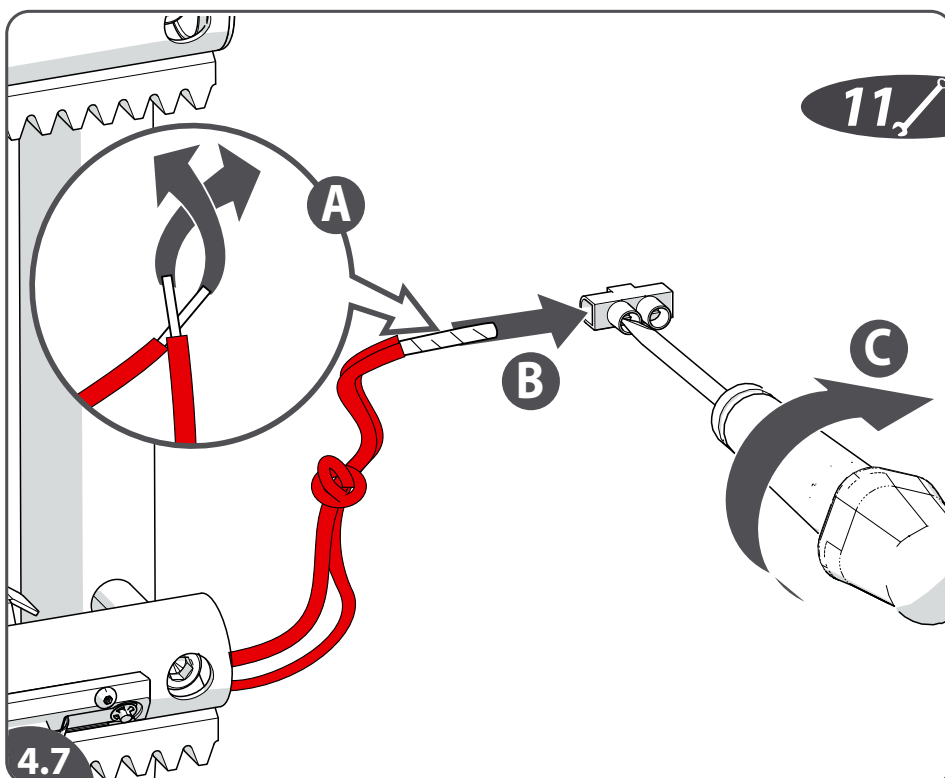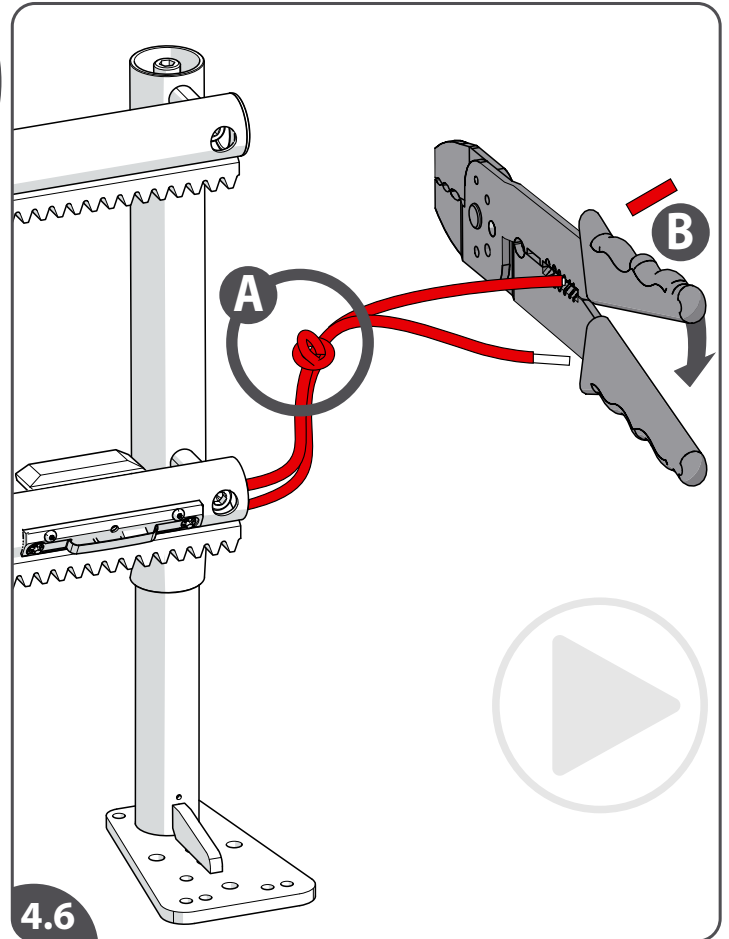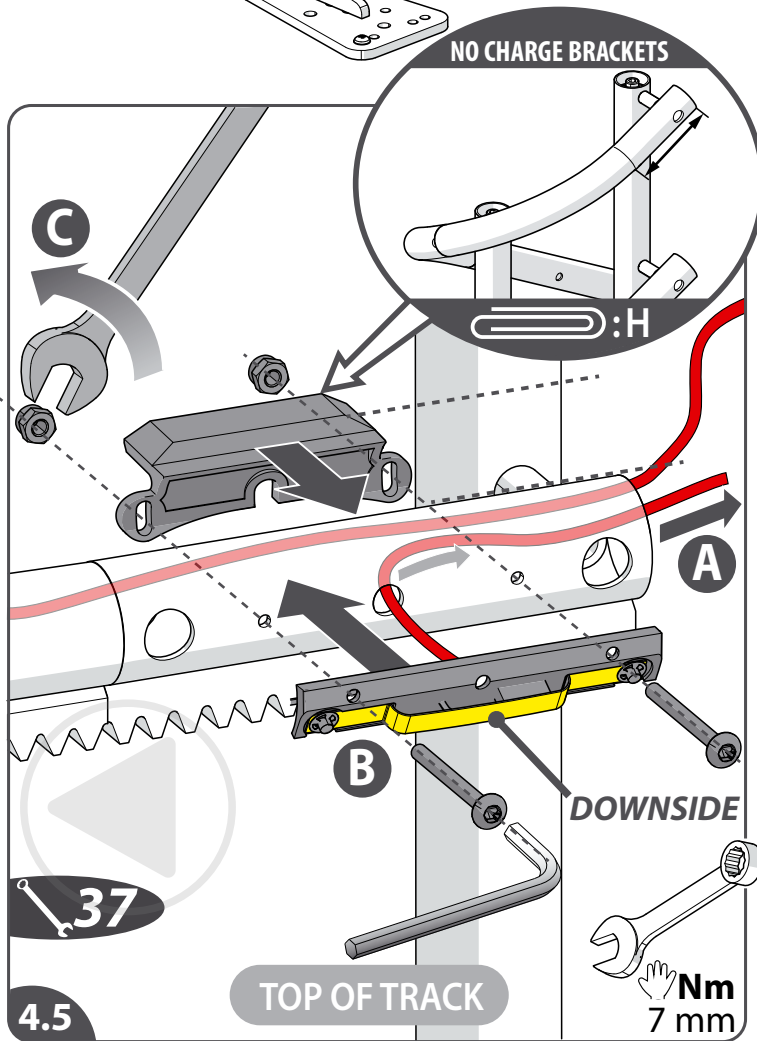[4.5] Feed the wire from the upper charge contact in the hole and guide it through the rail tube{A}. Install the upper charge contact and the limit actuator{B}. Tighten the bolts{C}.

Note: the limit actuator needs to be correctly aligned on the rail so that both switches can contact the limit actuator (see check 8.11 on page 118).

In case the last straight railpart is very short, mount special brackets (see attachment H).

[4.6] Tie the extra wiring to the contact wire{A}. Remove the insulation of the two wires (15 mm){B}.

[4.7] Twist the two ends bare wires together{A}. Fasten the wires in a connector block{B+C}.

EN [4.8] Secure and insulate the connector block.

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



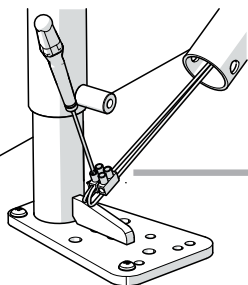
TRACK ELECTRICAL

REMOVE THE PAINT

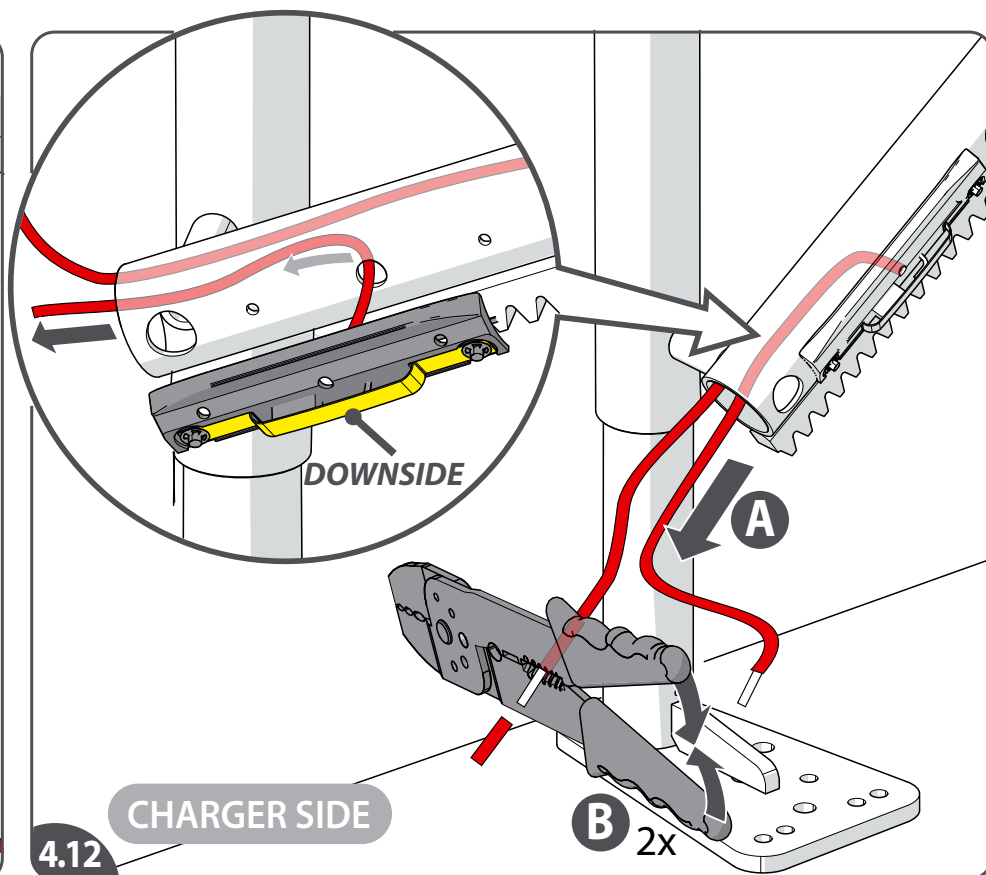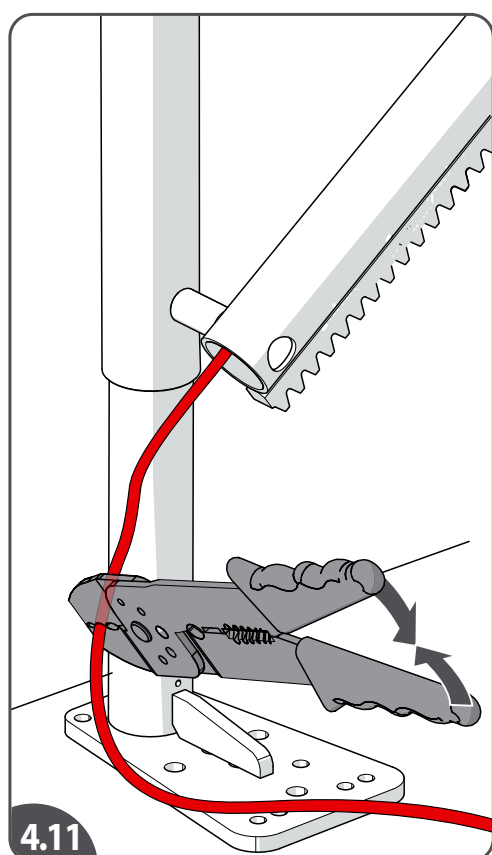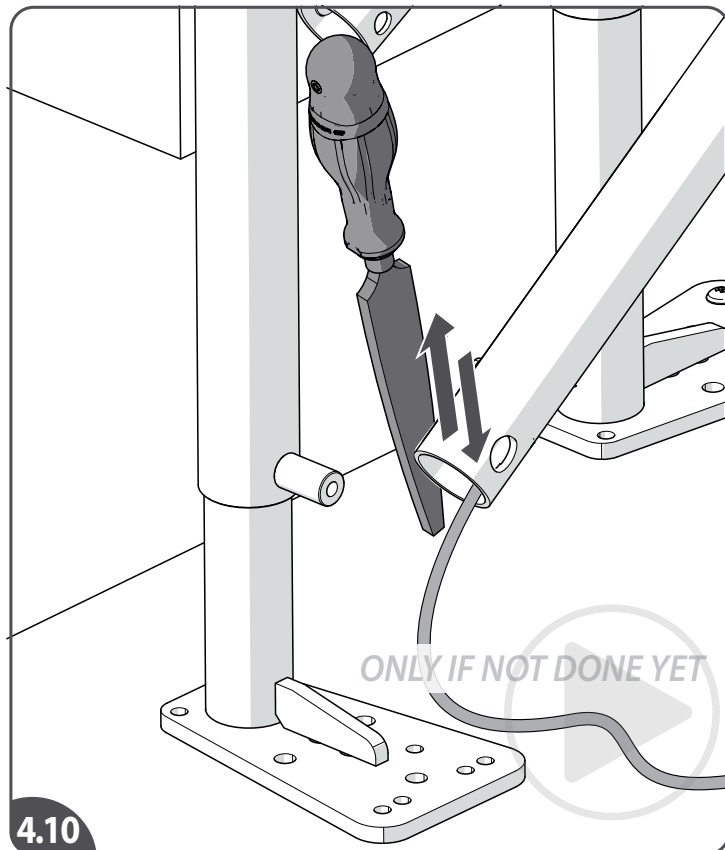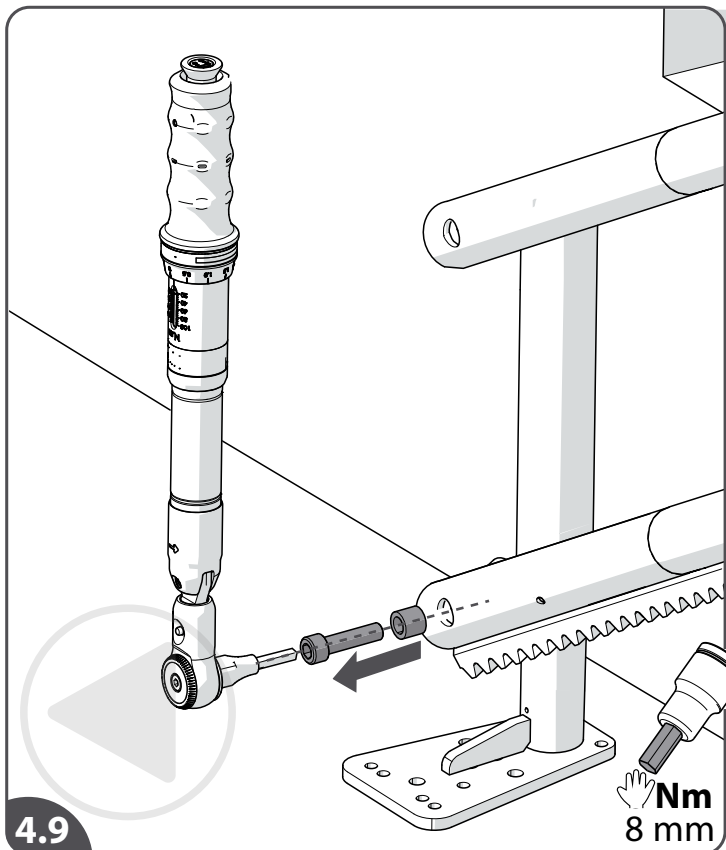- [4.9] Remove the bolt from the bottom stanchion.
- [4.10] Remove the paint from the rail.
- [4.11] Cut the excess wire off the charge cable.
- [4.12] Feed the wire from the bottom charge contact (different from upper charging contact !) in the hole and guide it through the rail tube{A}. Remove the insulation of the two wires (15 mm){B}.

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




IT

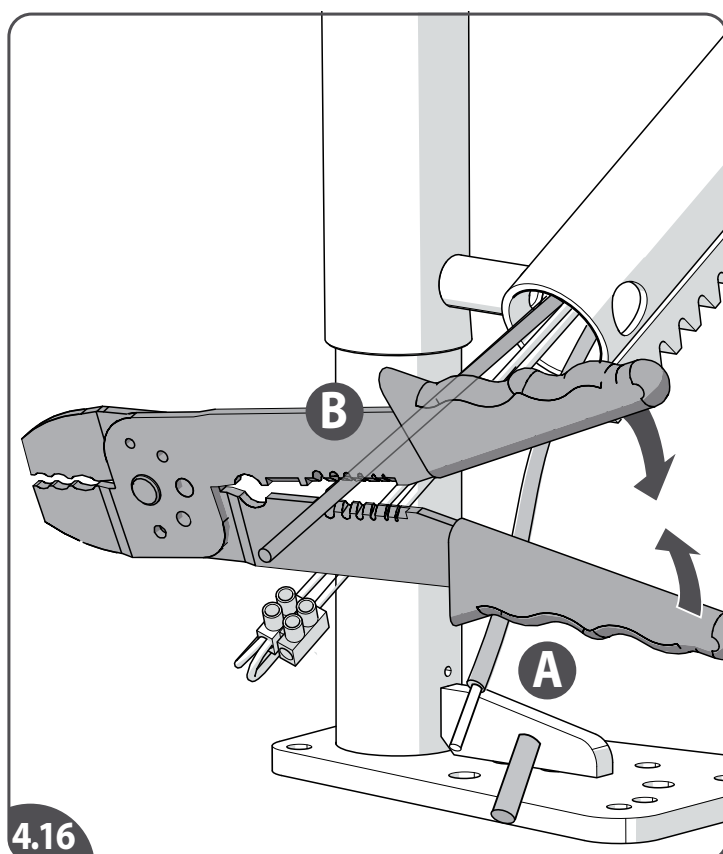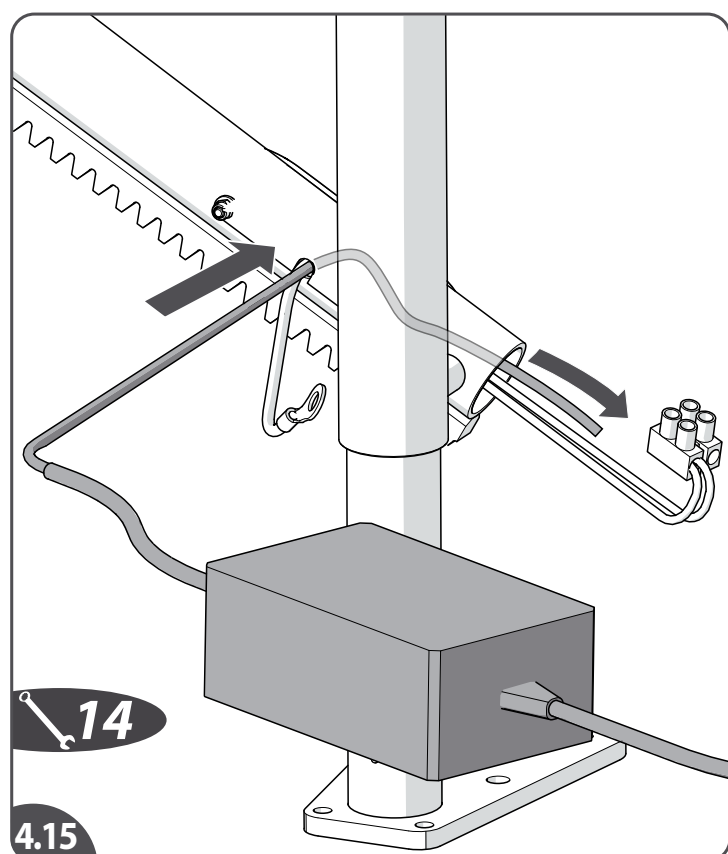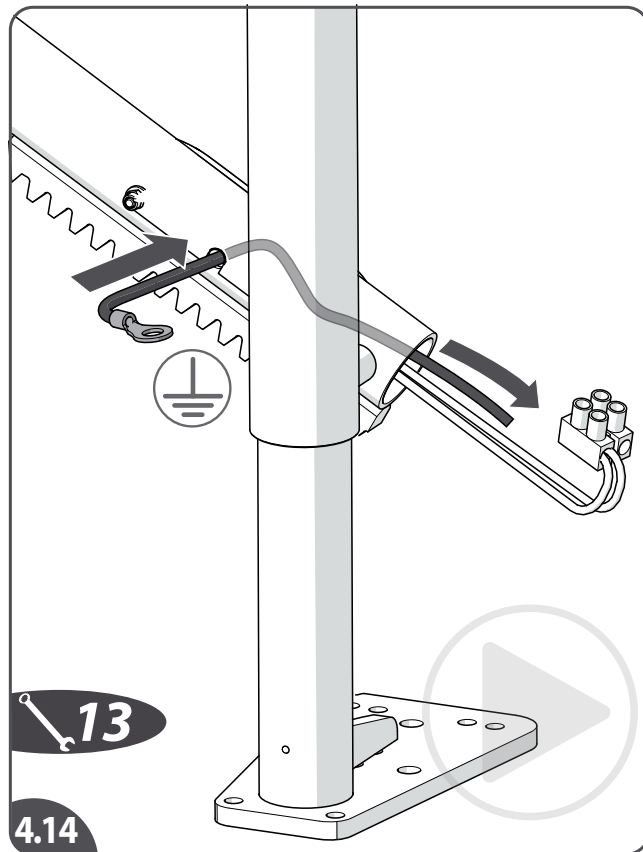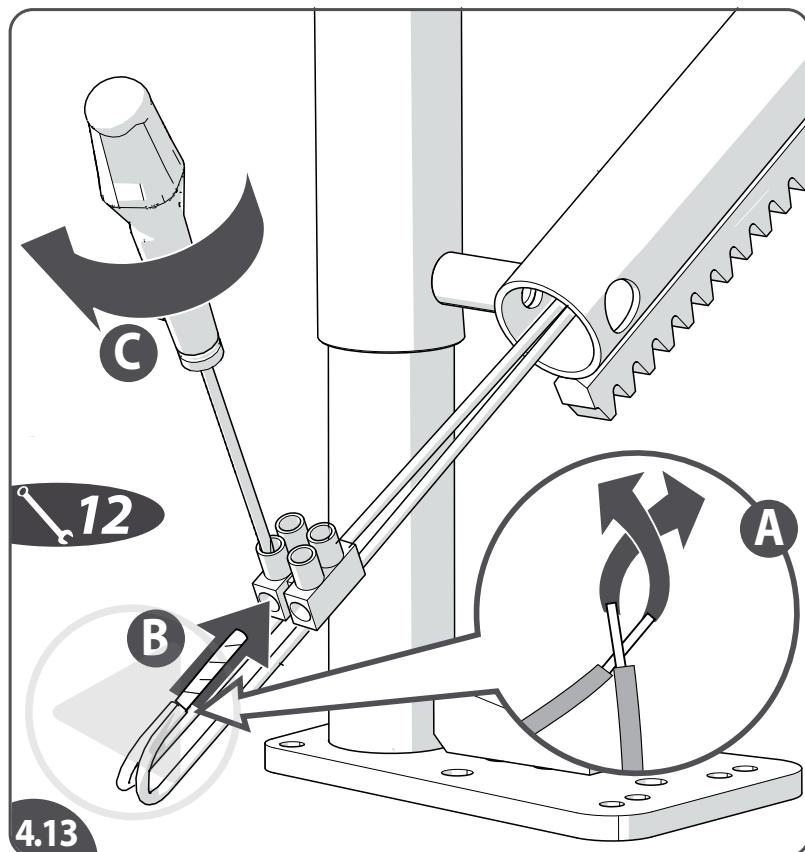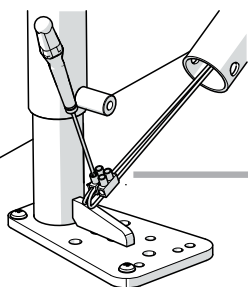- [4.13] Wrap the two red wires together{A}. Fasten the wires in the connector block{B+C}.
- [4.14] Feed the black earth wire in the hole at the rear and guide it through the rail tube.
- [4.15] Feed the charger cable through the same hole and guide the cable through the rail tube.
- [4.16] Remove the insulation of the earth wire (15 mm){A}. Remove the insulation of the charger cable (25 mm){B}.

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




IT

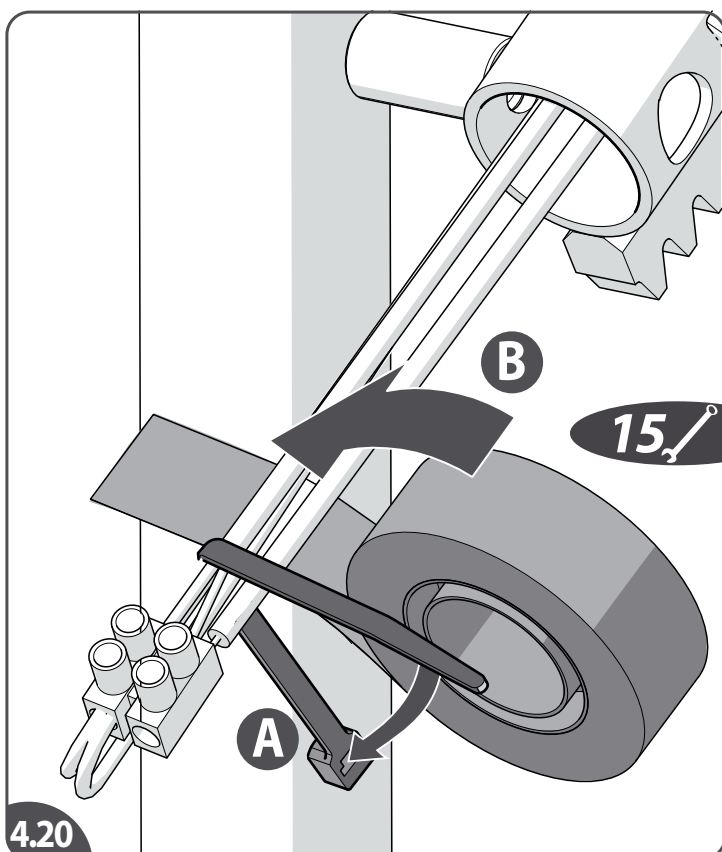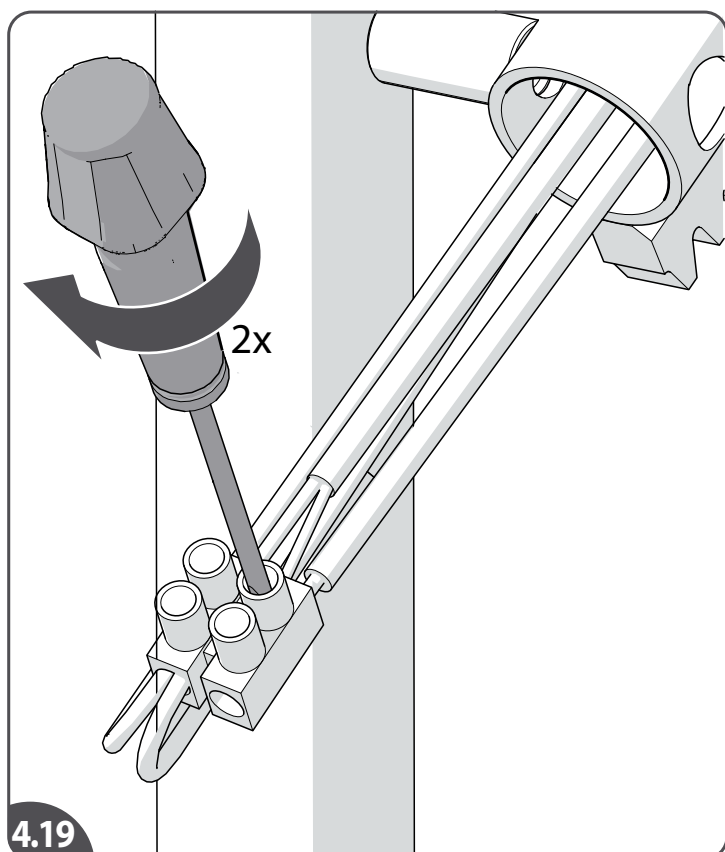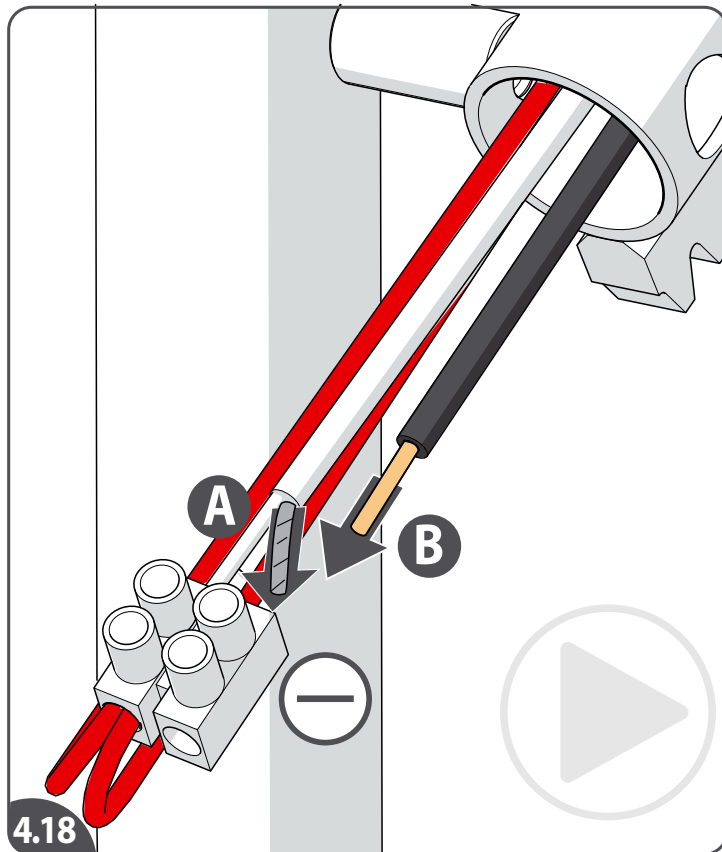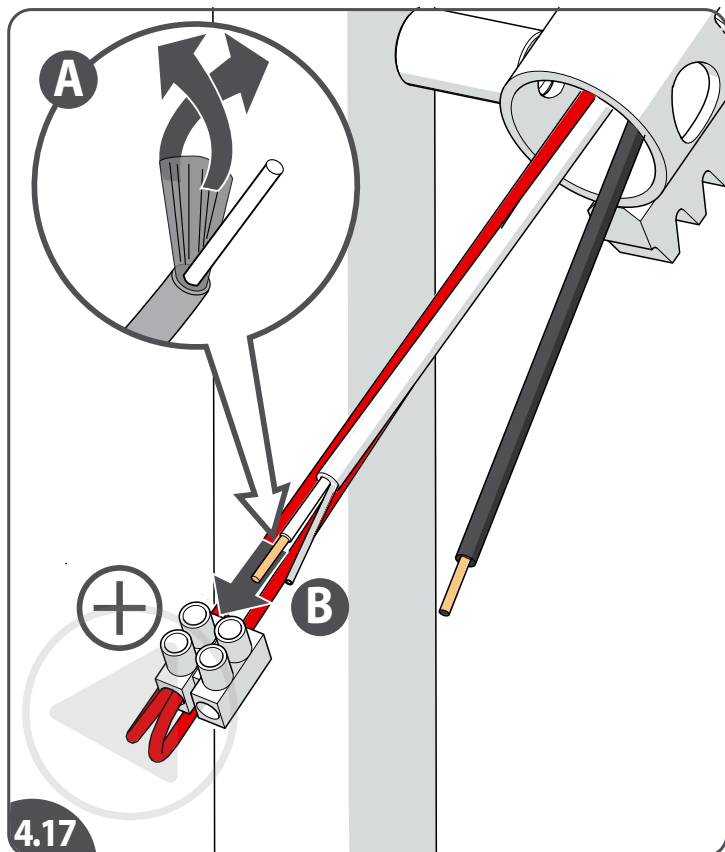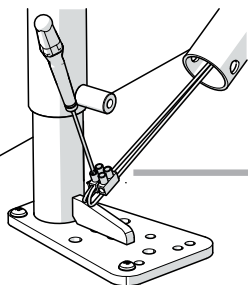- [4.17] Twist the outer wires together{A}. Feed the charger cable in the connector block: the positive wire against the red wires{B}.
- [4.18] The negative wire of the charger cable in the other terminal of the block{A}. The earth wire also in this terminal of the block{B}.
- [4.19] Fasten the wires.
- [4.20] Secure the wires with a cable tie{A}. Secure and insulate the connector block{B}.

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




IT

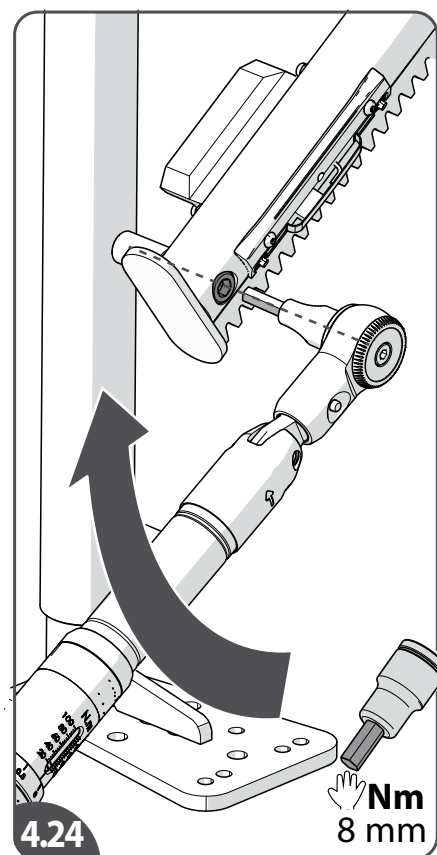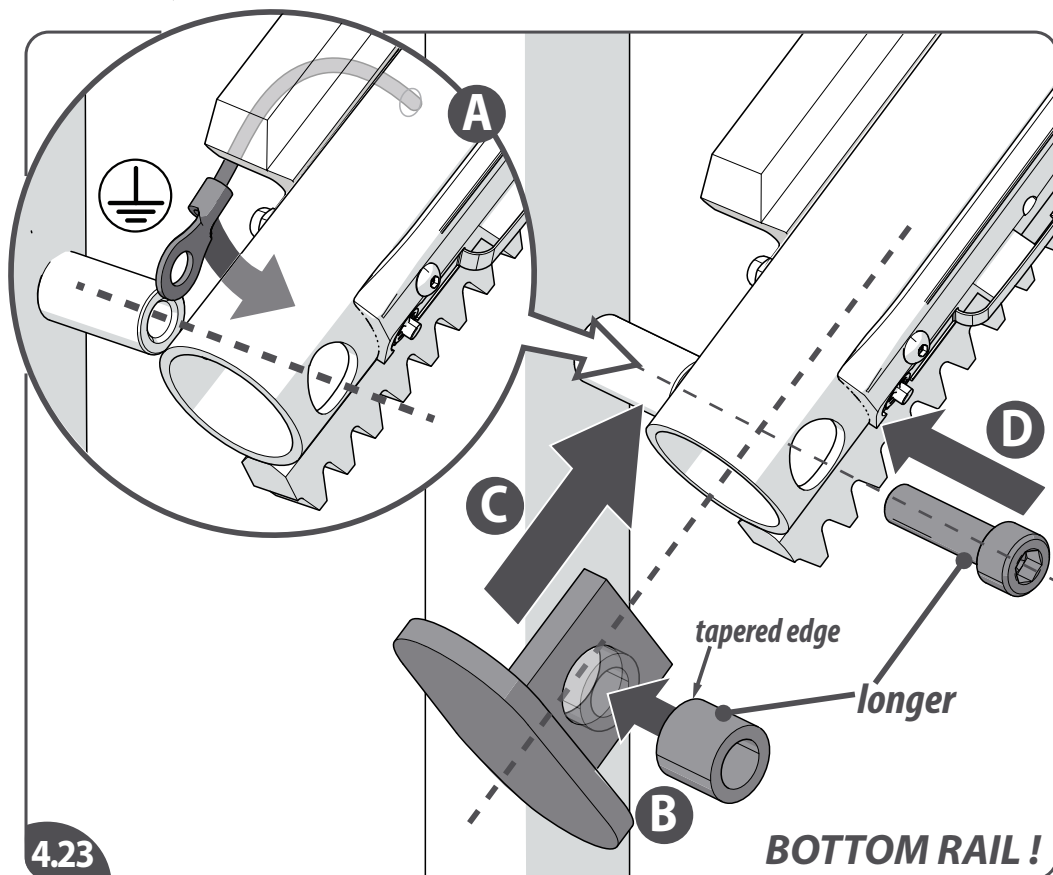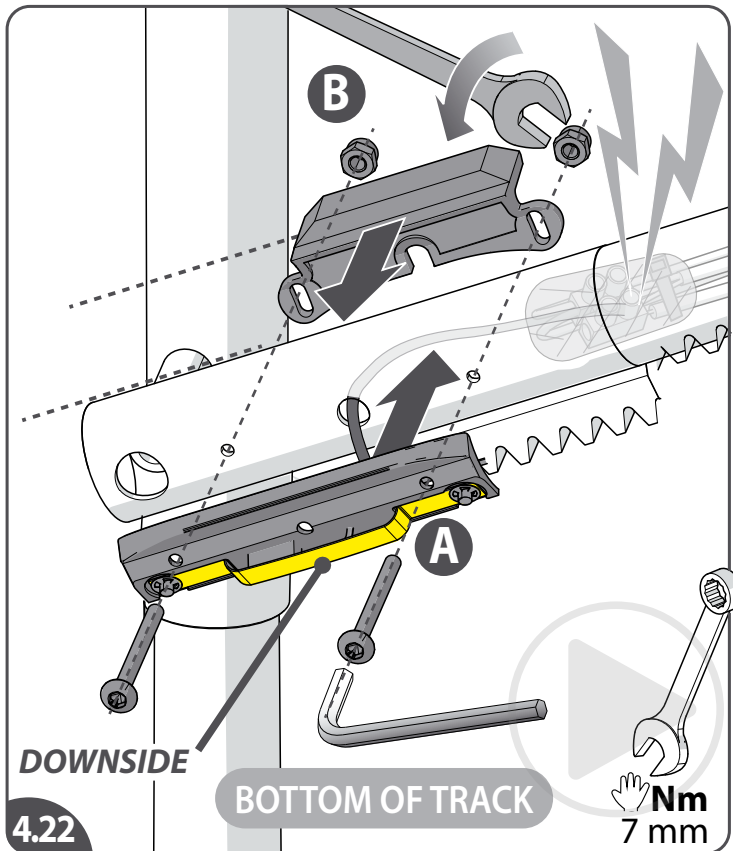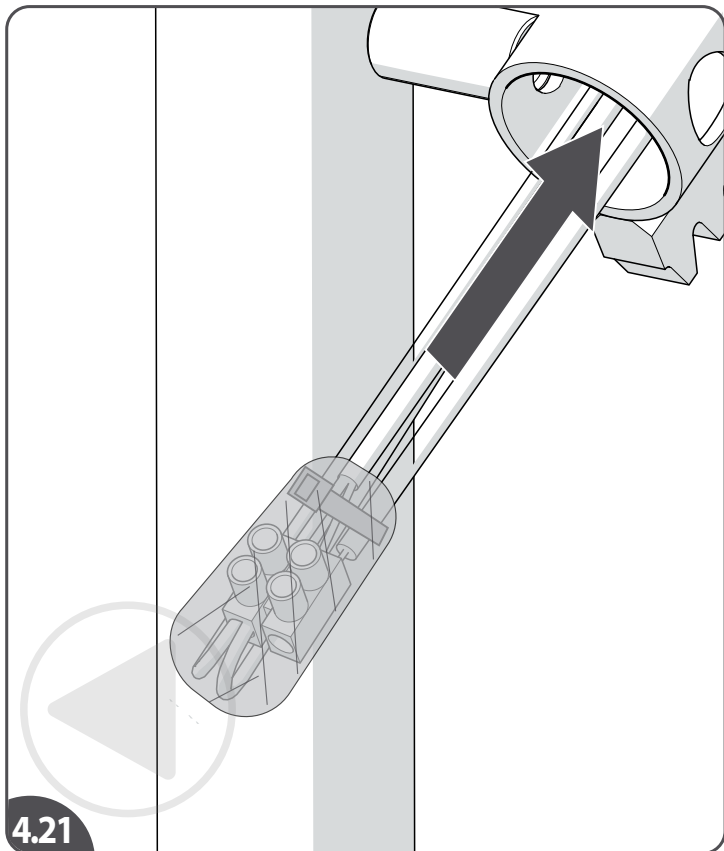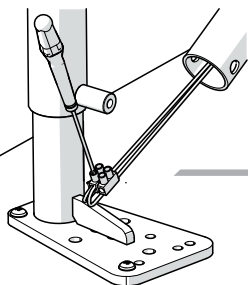[4.21] Feed all the wires into the end of the rail.

[4.22] Install the bottom charge contact and the limit actuator{A}. Tighten the bolts{B} (Careful: do not damage wiring inside tube!). **Note:** the limit actuator needs to be correctly aligned on the rail so that both switches can contact the limit actuator (see check 8.11 on page 118).

[4.23] Position the earth wire between lug and rail{A}. Position the spacer in the endstop{B}. Install the end stop{C}. Place stanchion bolt{D}.

[4.24] Tighten the stanchion-rail bolt [hand tight].

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




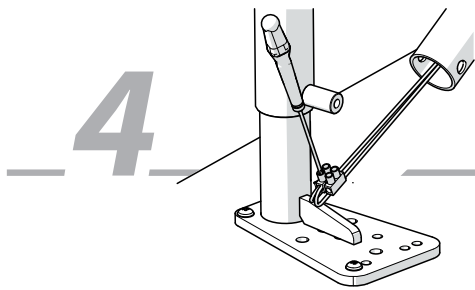
TRACK ELECTRICAL

TIEWRAPS

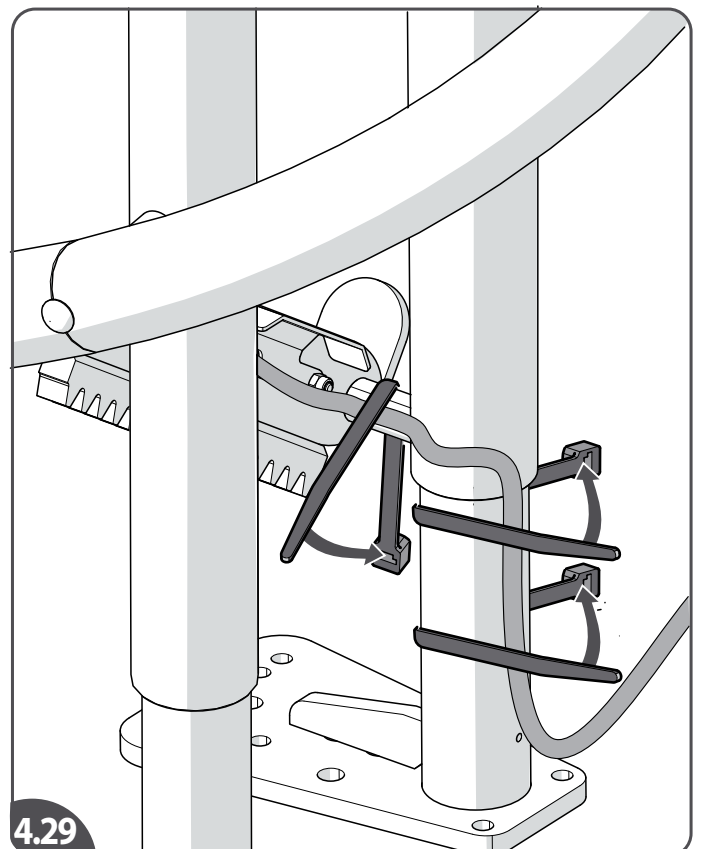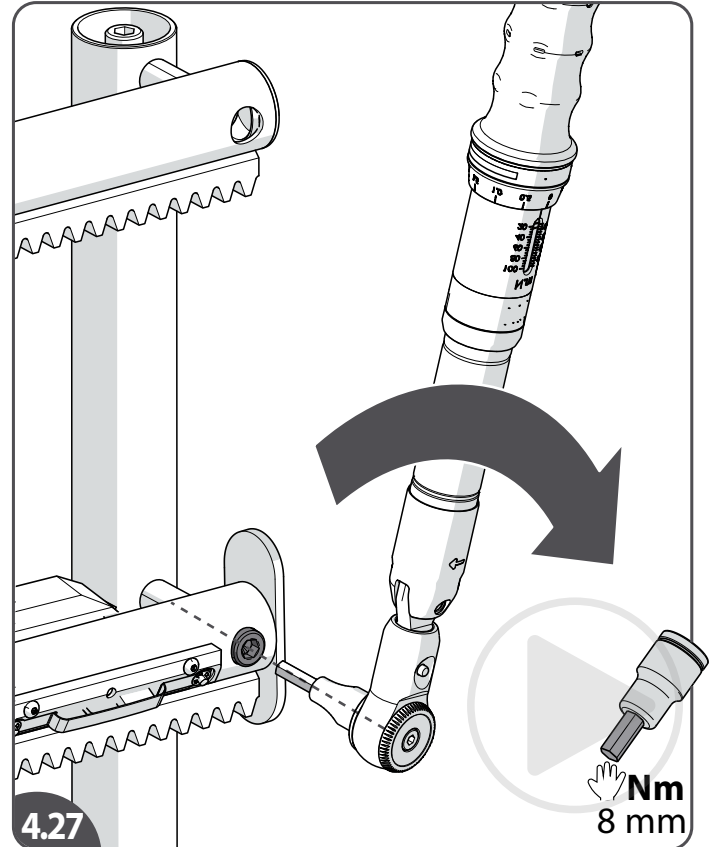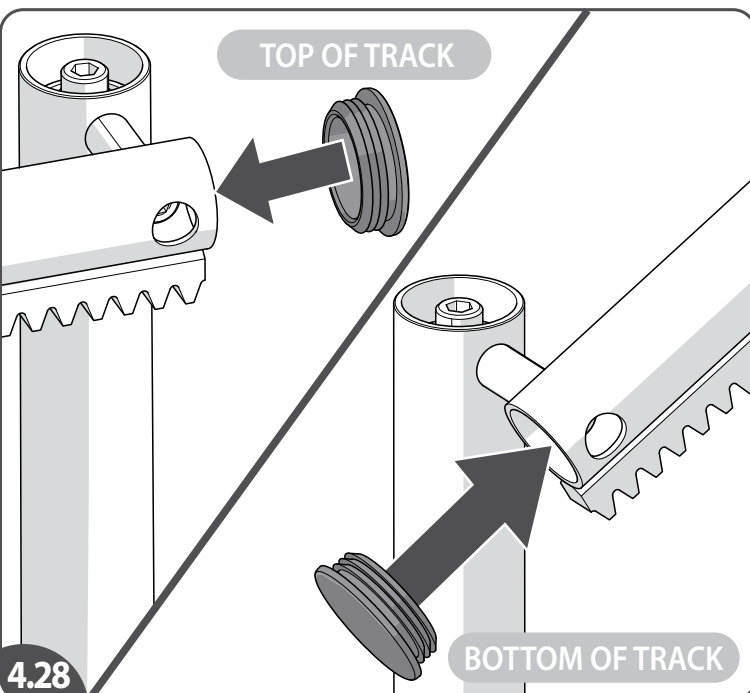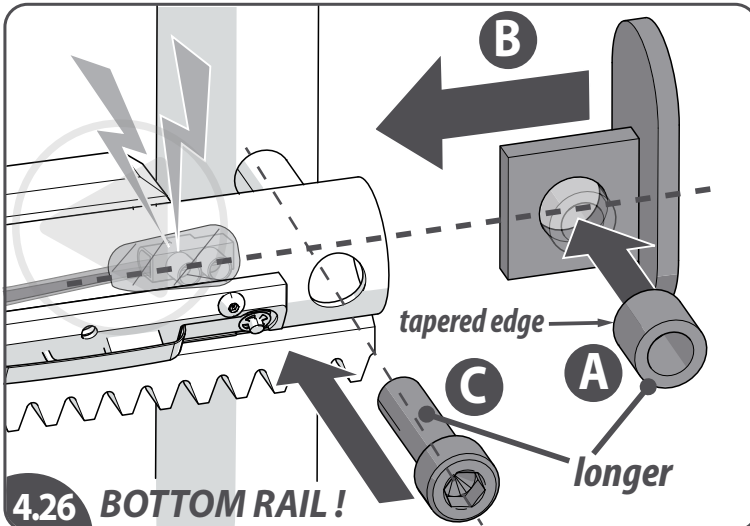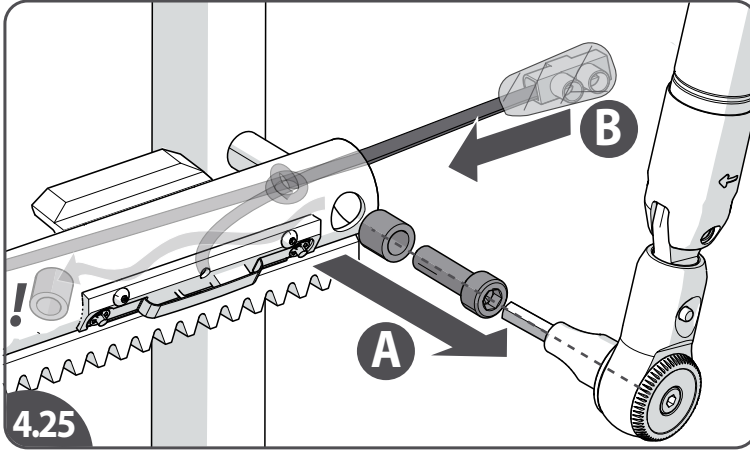- [4.25] Remove the stanchion-rail bolt+spacer (do not let the spacer fall into the tube!){A}. Feed all the wires into the end of the rail{B}.
- [4.26] Position the spacer in the endstop{A}. Install the end stop{B}. Place stanchion bolt (Careful: do not damage wiring inside tube !){C}.
- [4.27] Tighten the stanchion-rail bolt [hand tight].
- [4.28] Fit the caps{A}. Replace stanchion bolt (Careful: do not damage wiring inside tube !){B}.
- [4.29] Use cable ties to fasten the charge cable.

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




IT

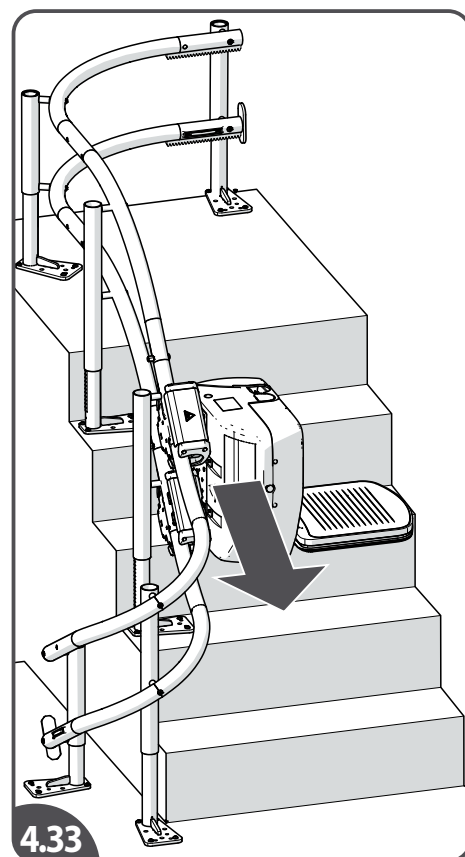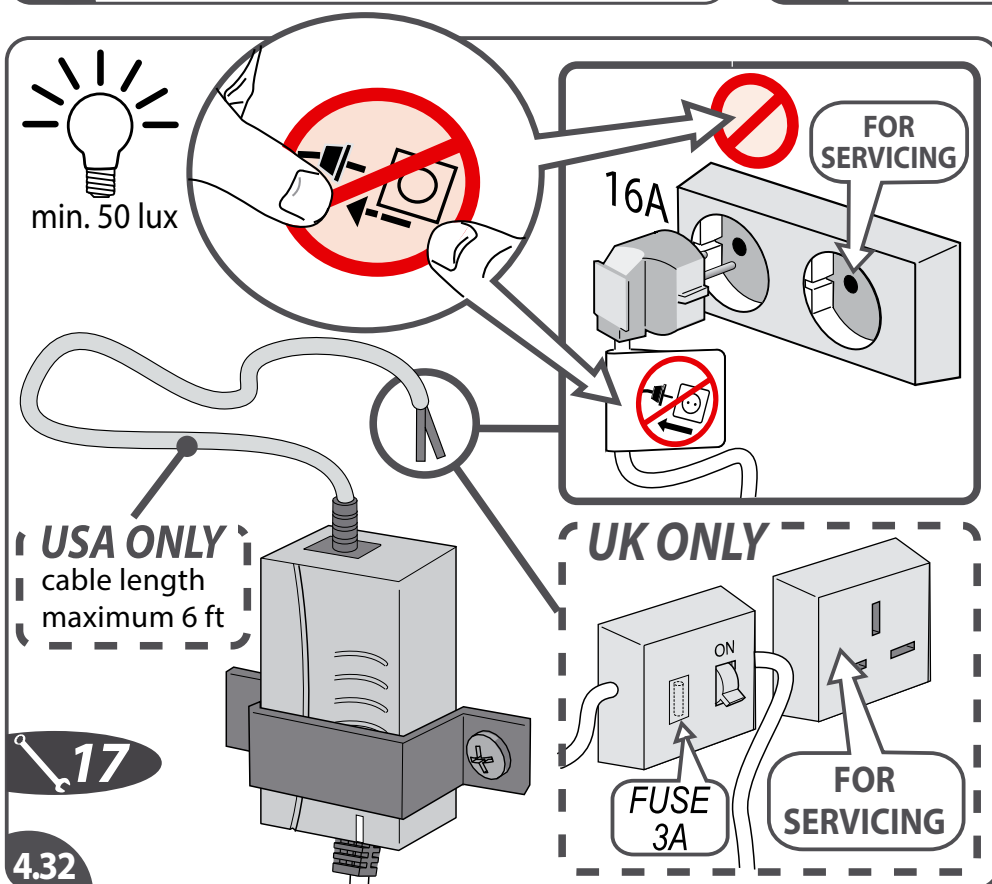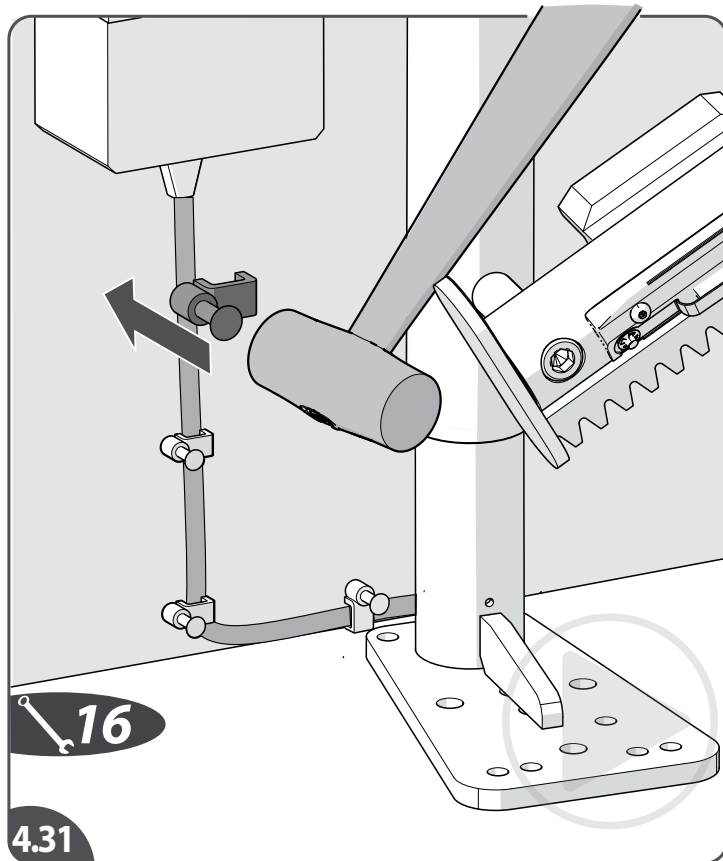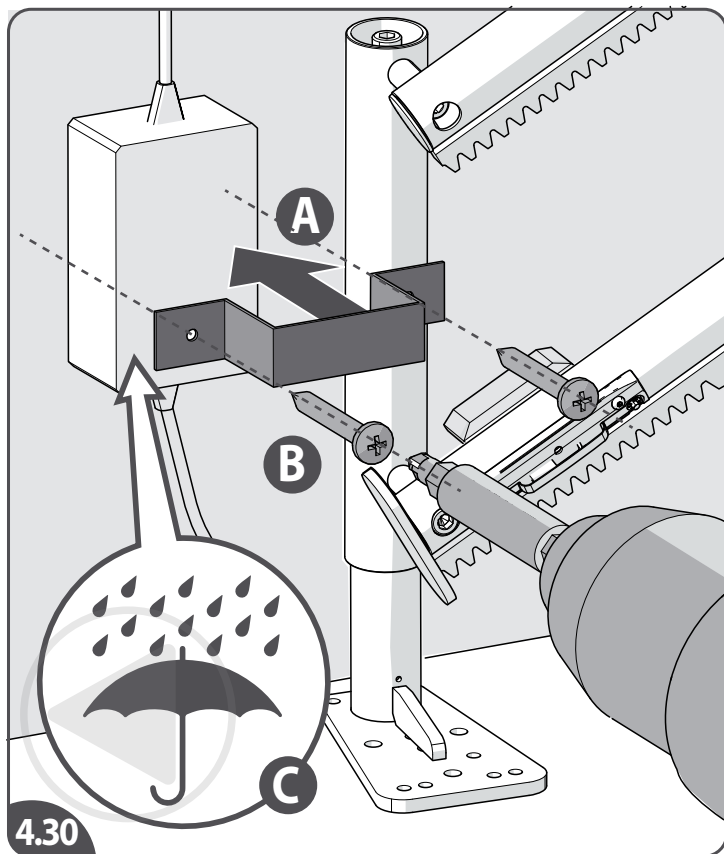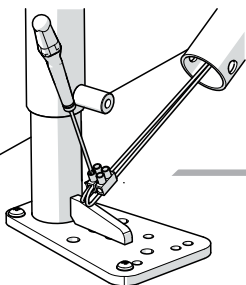[4.30] Fit the power supply to the wall {A+B}. For outdoor cases: the charger must to be placed **indoor** or in a **waterproof** casing {C}.

[4.31] Fasten the charge cable to the wall.

[4.32] Terminate the power supply to the ac supply. All mains wiring must comply with the local standards in force at time of installation and with section 4.3 of EN60204/1 and also with the requirements of NEC and ANSI/NFPA70. 3-pin plug not supplied.

[4.33] Drive the power pack all the way down.

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




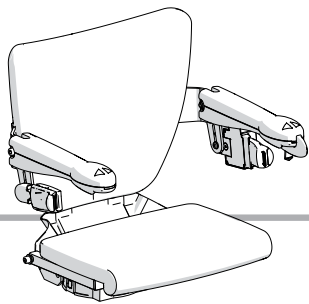
STYLE SEAT

LOOMS THROUGH STEM

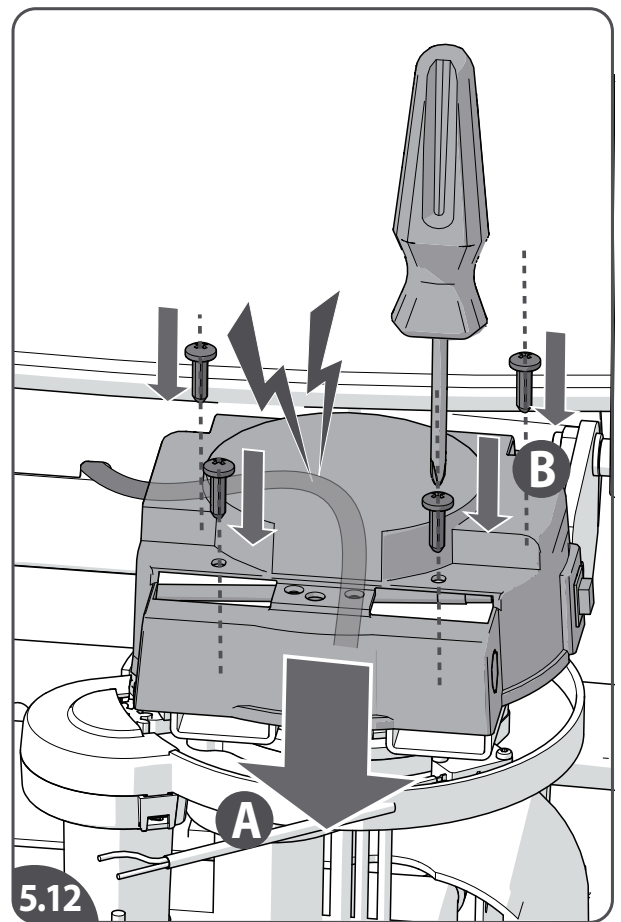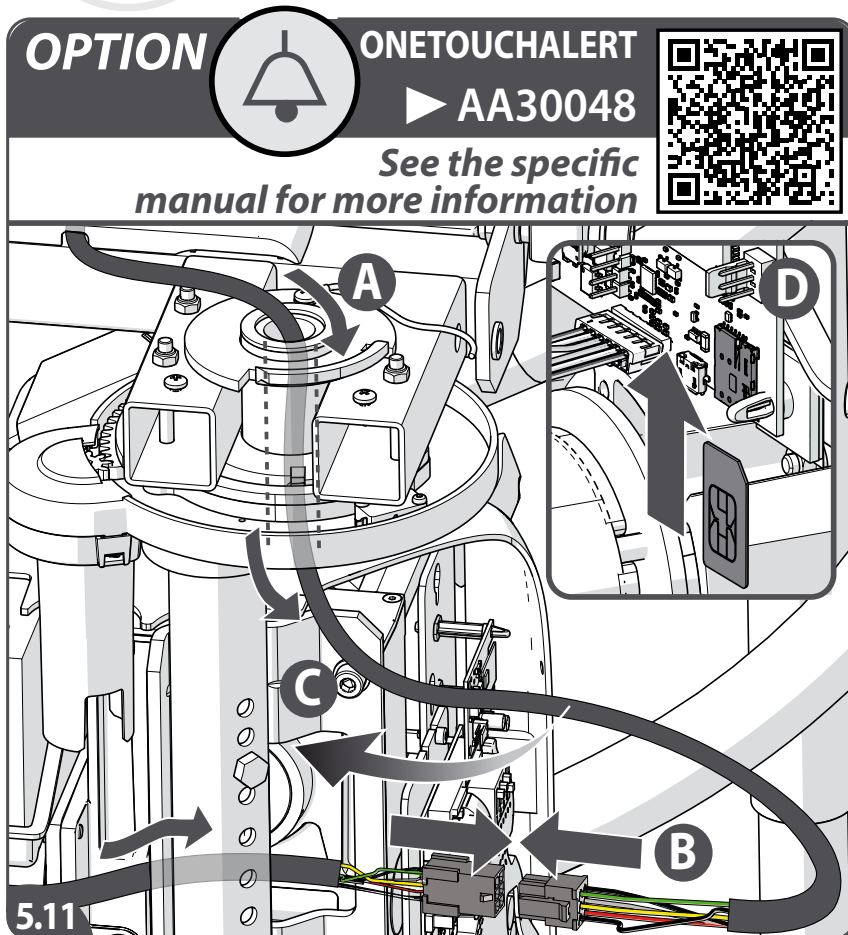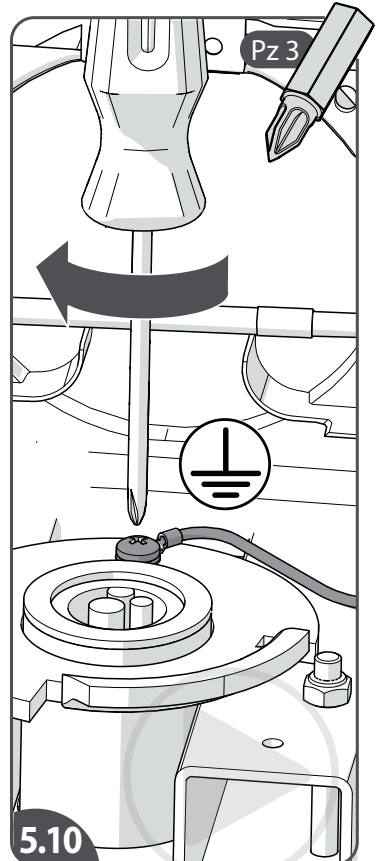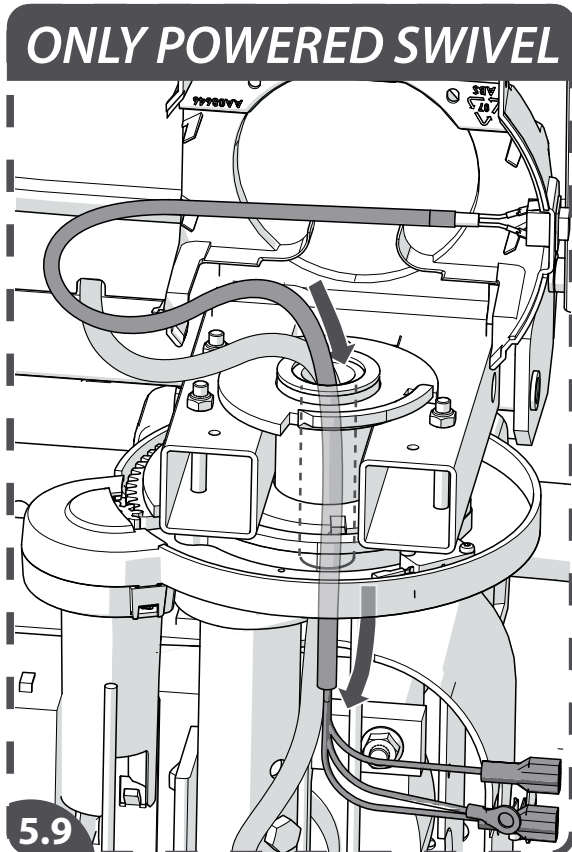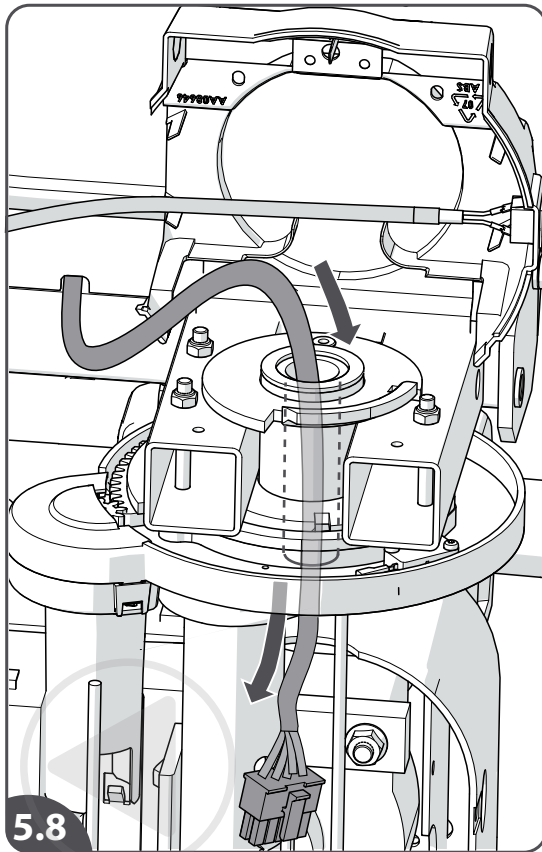- [5.8] Feed the seat loom through the stem.
 [5.9] If present, feed the swivel release loom through the stem.
 [5.10] If not fastened yet: fasten the earth wire.
 [5.11] If present on the lift, guide the OTA loom {armrest part} through the chair and chair axle, guide the OTA loom {print part} behind the seat stem {A}. Connect the OTA loom parts (U28) {B}. Tuck the excess cable + connectors behind the seat stem {C}. Please refer to instruction AA30048 for more information. Install the SIM card {D}.
 [5.12] Reinstall the stem cap.

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
5a

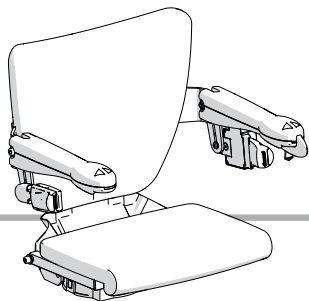STYLE SEAT DIPSWITCH SETTINGS

IN THE COLOURED SCHEME ON THE BACKCOVER YOU CAN FIND THE MARKED CONNECTIONS.

U6

U16

5.13

ONLY POWERED SWIVEL

U27

5.14

ONLY POWERED SWIVEL

7 mm

A

B

5.15

D1

- 1 intermediate stop
- 2 always deactivated
- 3 always deactivated
- 4 activated

D2

- 1 foldable track
- 2 right/left
- 3 activated
- 4 always deactivated

* Not during installation

D1

D2

5.16

- [5.13] Connect the seat loom (**U6**) and connect the swivel cable (**U16**) .
- [5.14] Connect the swivel release loom to the motor loom (**U27**).
- [5.15] Connect the earth wire to the chassis: fit the earth wire of the swivel release loom{A} and fasten these earth wire{C}.
- [5.16] If needed, change the settings on the dip switches on the main PCB (power supply/power switch should be off before changing): e.g. in case of a hinged rail(D2-1) or in case of a right hand configuration(D2-2). D2-3 should be always on. Look at the inside of the front cover for the settings (or attachment D).

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








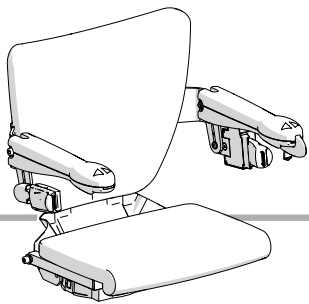
STYLE SEAT

CLOSING FRONT COVER

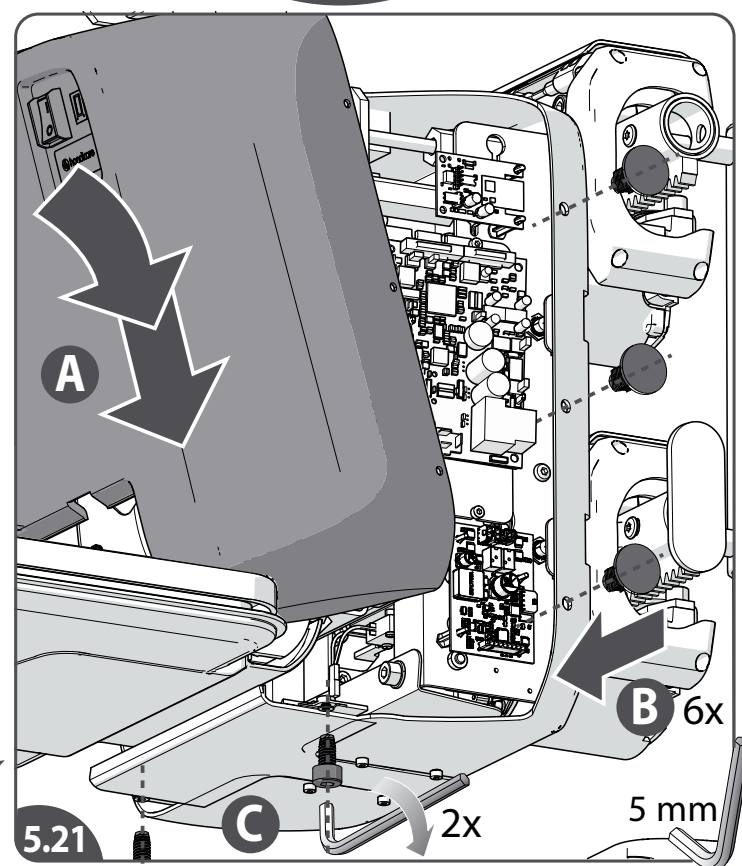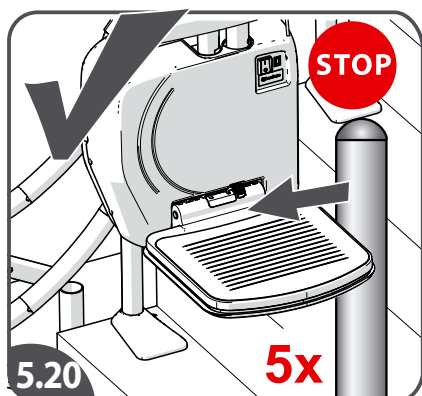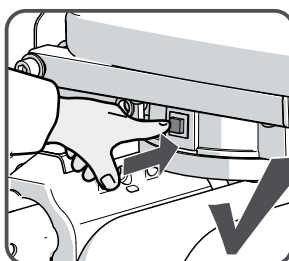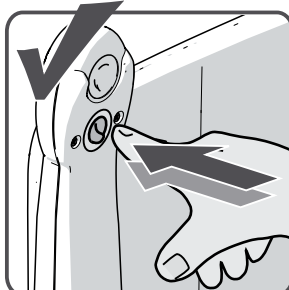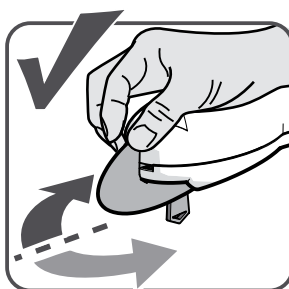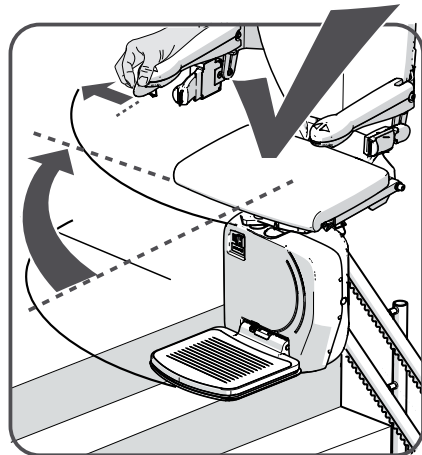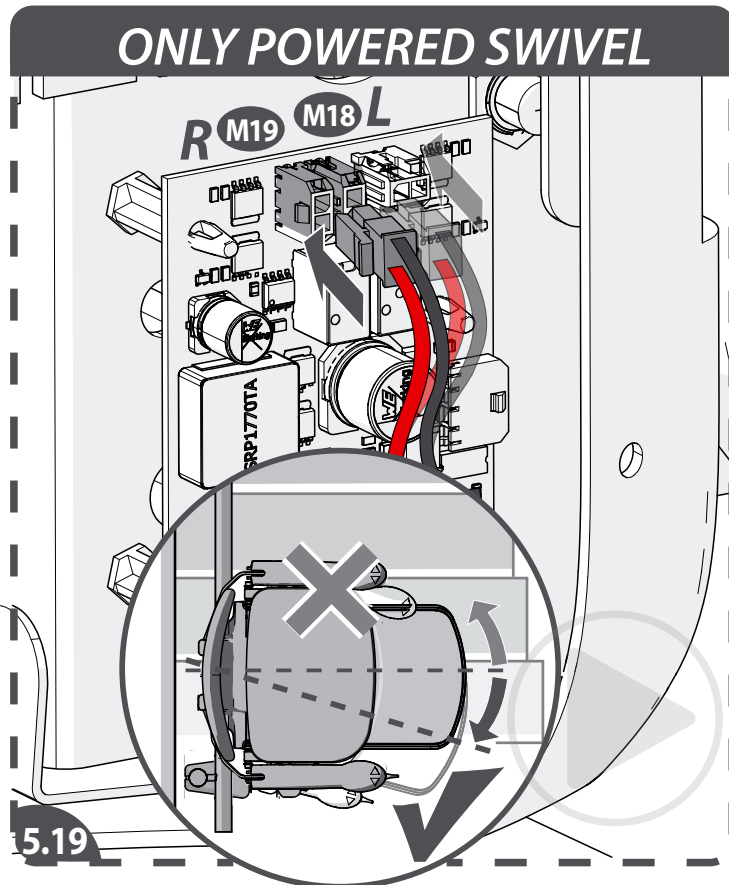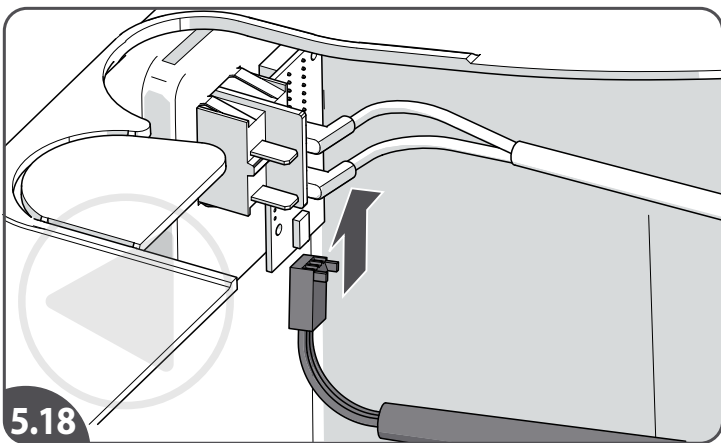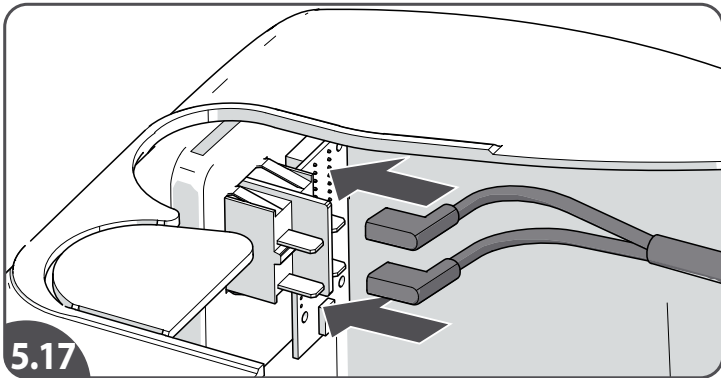- [5.17] Connect the power switch cables.
- [5.18] Connect the diagnostic display loom.
- [5.19] Connect the wires of the swivel motor.
- [5.20] Check the functions and safeties before installing the cover.
- [5.21] Reinstall the front cover{A} and secure with the rivets(6x){B} and screws(2x){C}.

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




IT

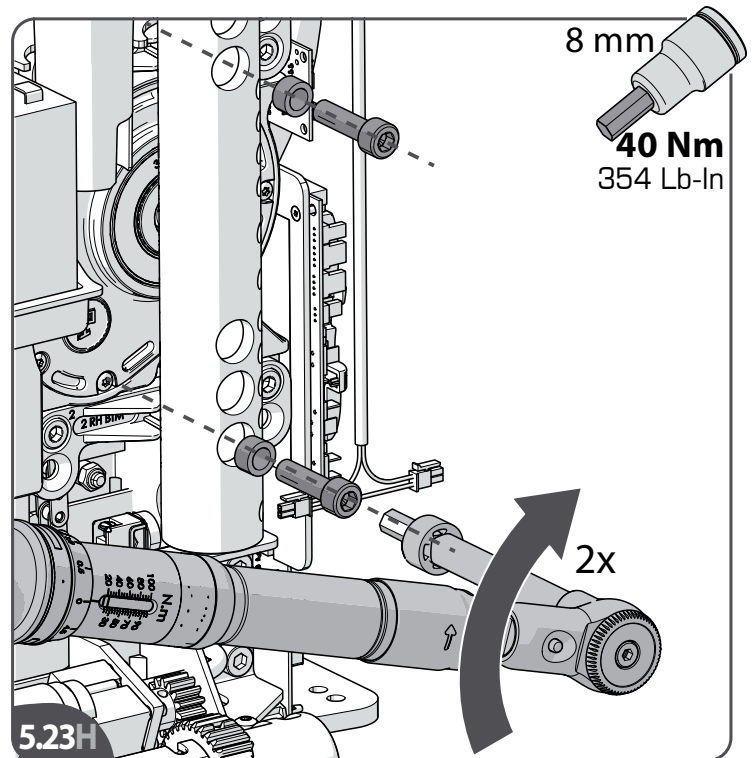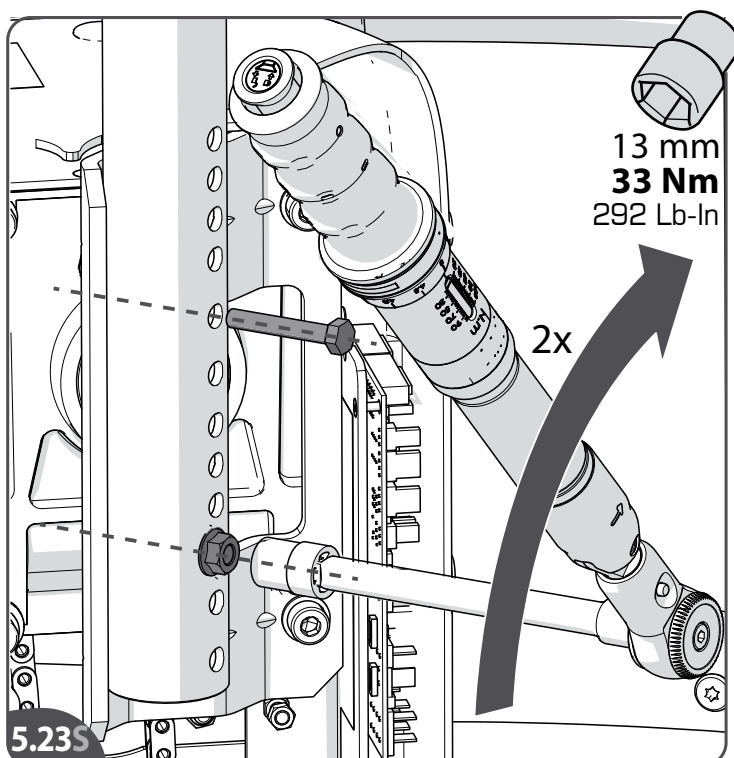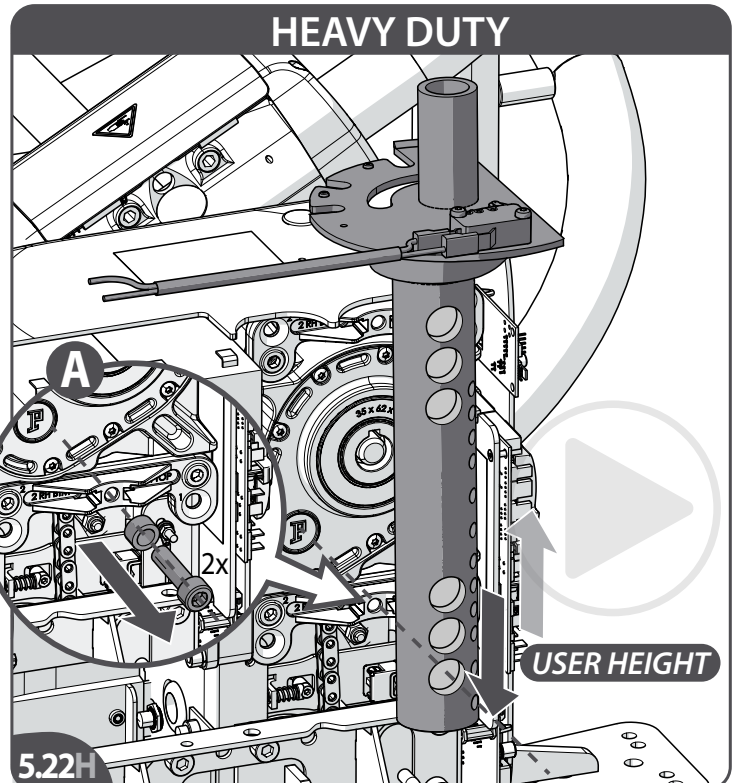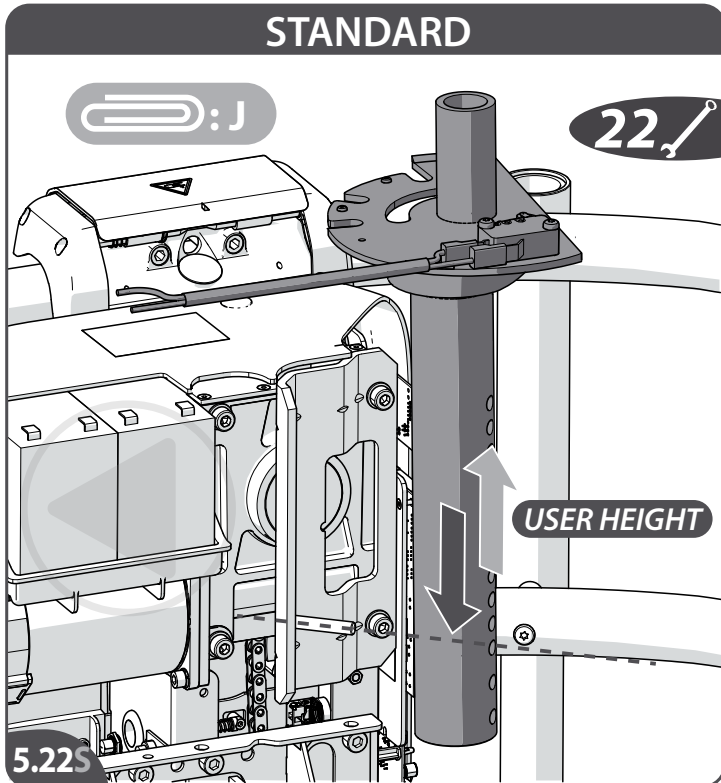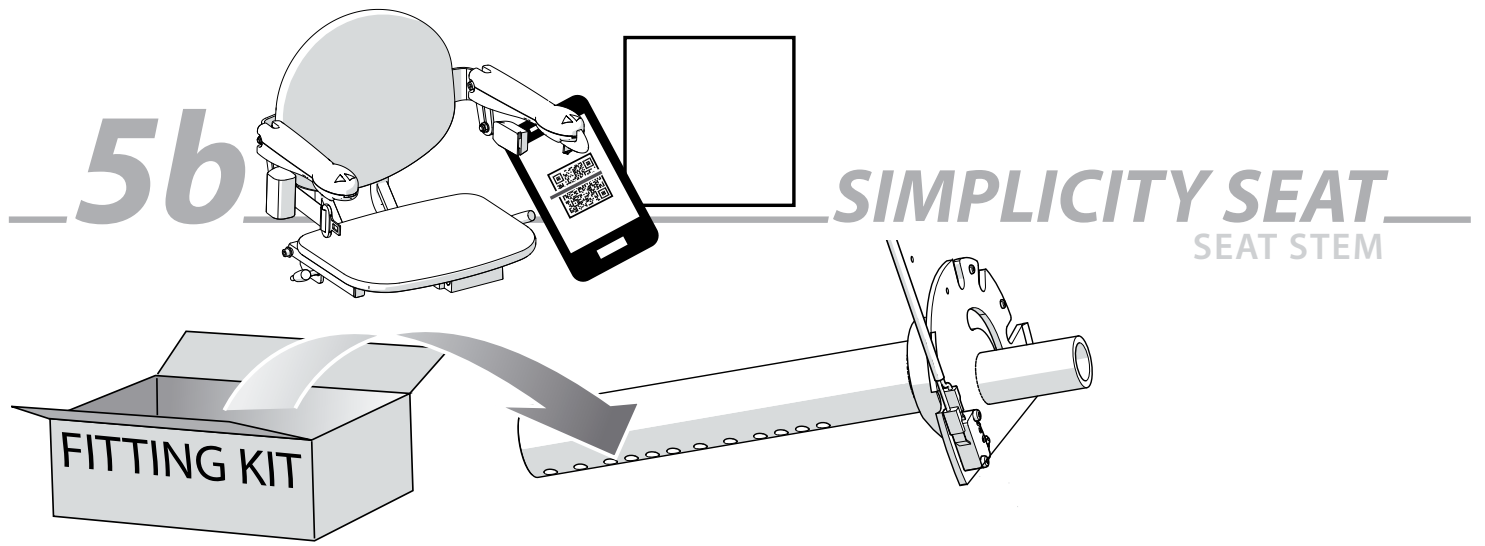Take the seat post out of the box. Carry the post to the power pack.

[5.22S]Position the seat post: fix the seat to the seat height required for the primary user of the stairlift.

[5.22H]Remove the bolts (reuse them){A}. Position the seat post: fix the seat to the seat height required for the primary user of the stairlift.

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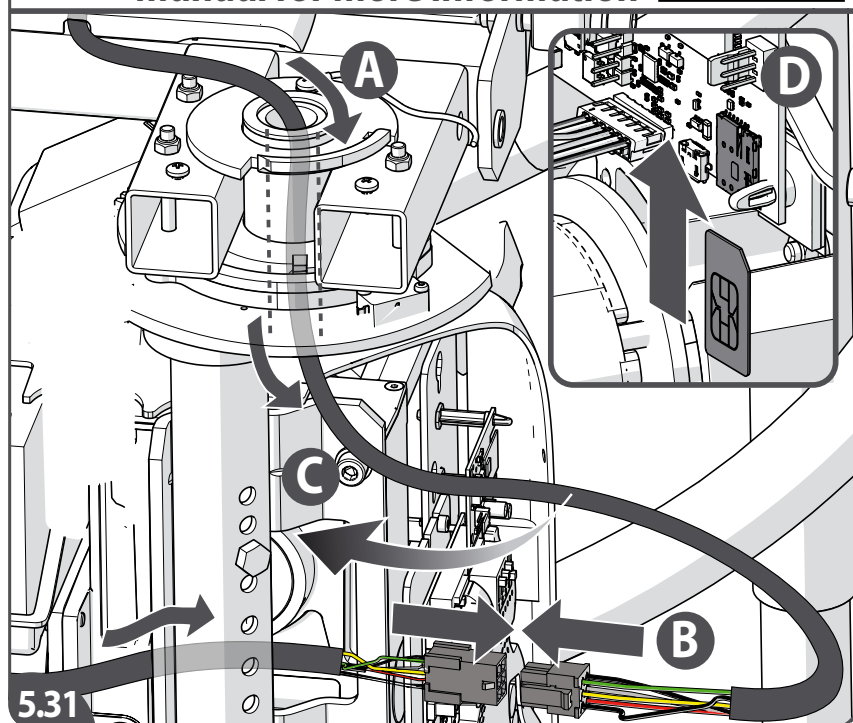
[5.29] Feed the seat looms through the stem.

[5.30] If not fastened yet: fasten the earth wire.

[5.31] If present on the lift, guide the OTA loom {armrest part} through the chair and chair axle, guide the OTA loom {print part} behind the seat stem {A}. Connect the OTA loom parts (U28) {B}. Tuck the excess cable + connectors behind the seat stem {C}. Please refer to instruction AA30048 for more information. Install the SIM card {D}.

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


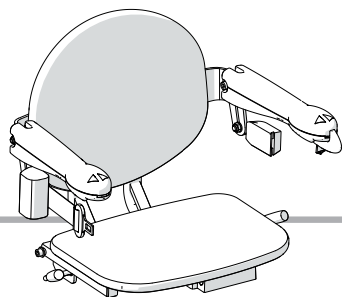
SIMPLICITY SEAT

CONNECTING SEAT LOOM

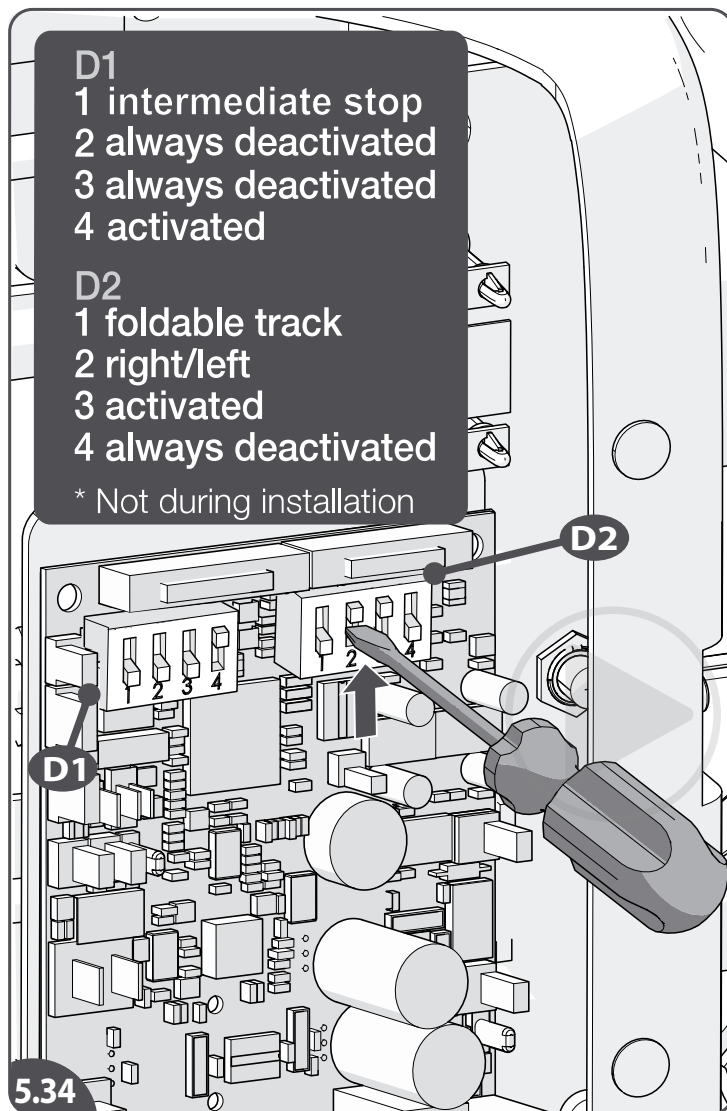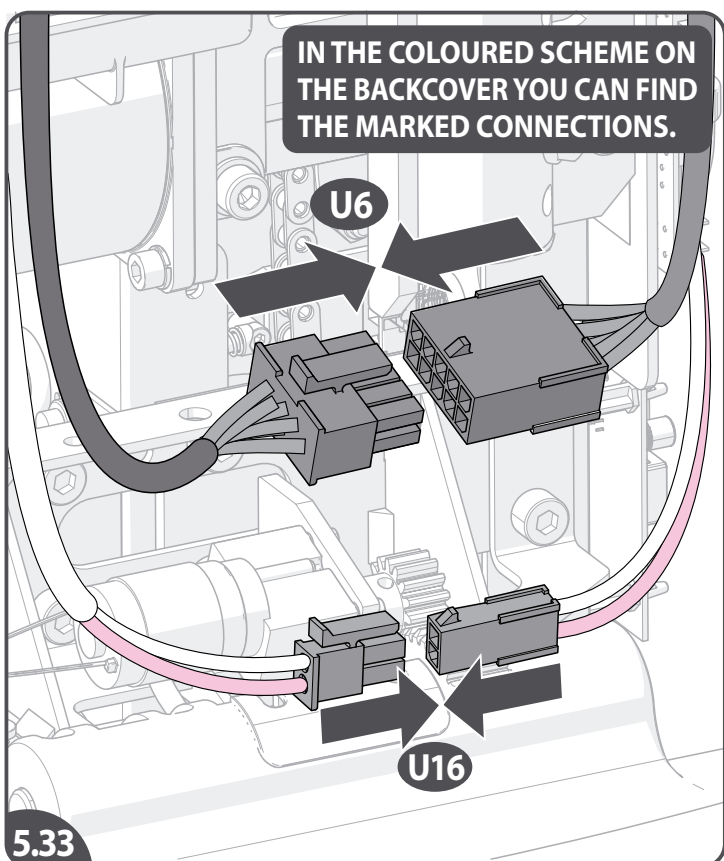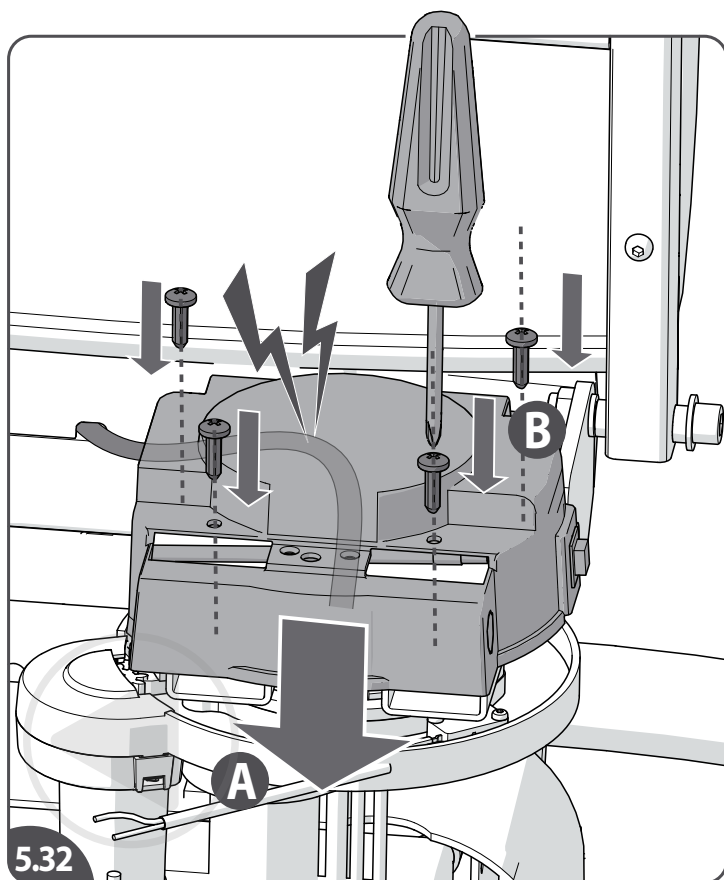D1
 1 intermediate stop
 2 always deactivated
 3 always deactivated
 4 activated

D2
 1 foldable track
 2 right/left
 3 activated
 4 always deactivated

* Not during installation

[5.32]Reinstall the stem cap.

[5.33]Connect the seat loom (**U6**) and connect the swivel cable (**U16**) .

[5.34]If needed, change the settings on the dip switches on the main PCB (power supply/power switch should be off before changing): e.g. in case of a hinged rail(D2-1) or in case of a right hand configuration(D2-2). D2-3 should be always on. Look at the inside of the front cover for the settings (or attachment D).

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



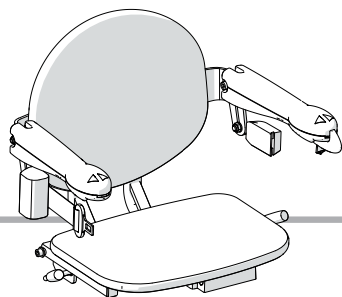
SIMPLICITY SEAT

CLOSING FRONT COVER

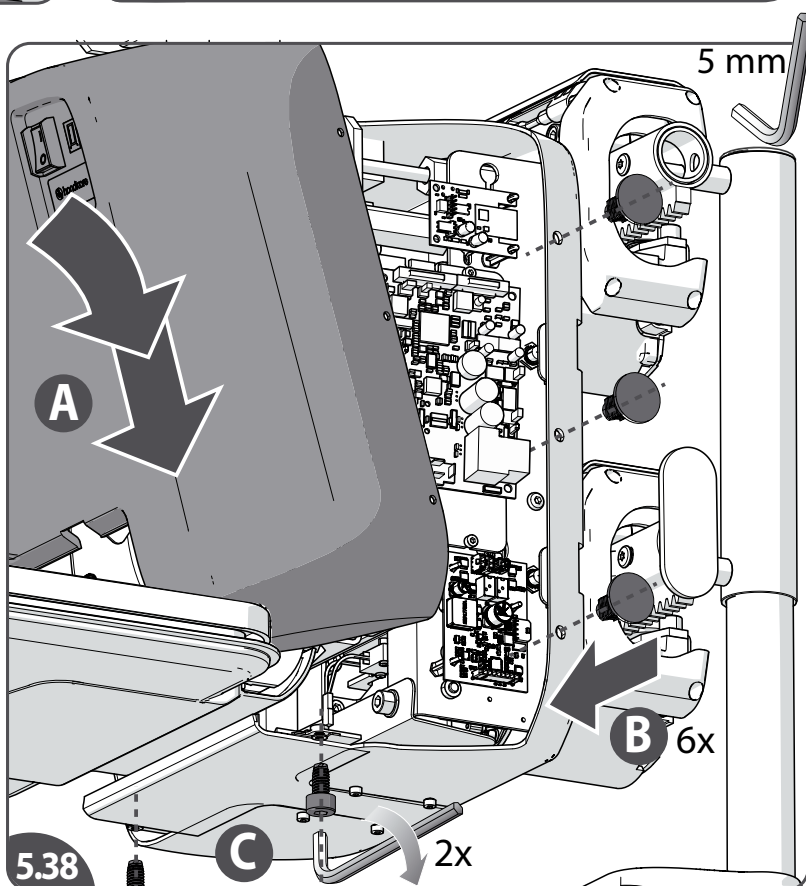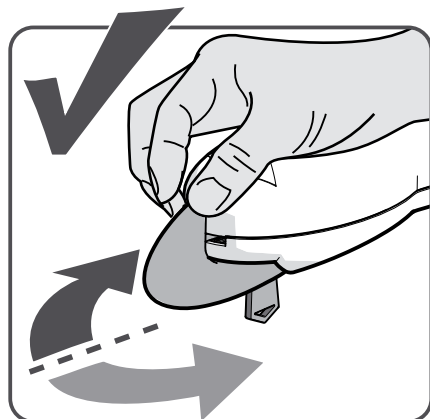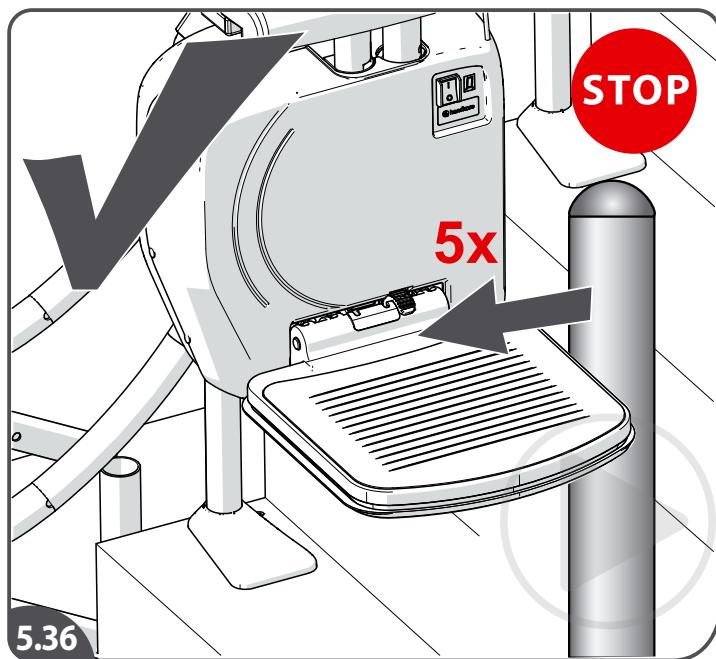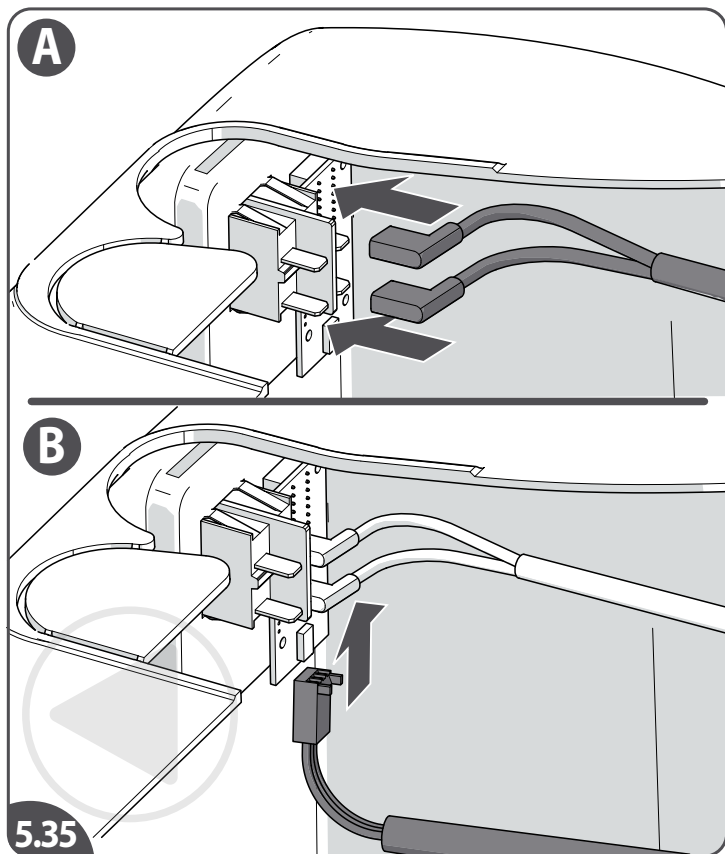- [5.35] Connect the power switch cables{A}. Connect the diagnostic display loom{B}.
- [5.36] Check the safeties before installing the cover.
- [5.37] Check the functions before installing the cover.
- [5.38] Reinstall the front cover{A} and secure with the rivets(6x){B} and screws(2x){C}.

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



5c

SMART SEAT SEAT STEM

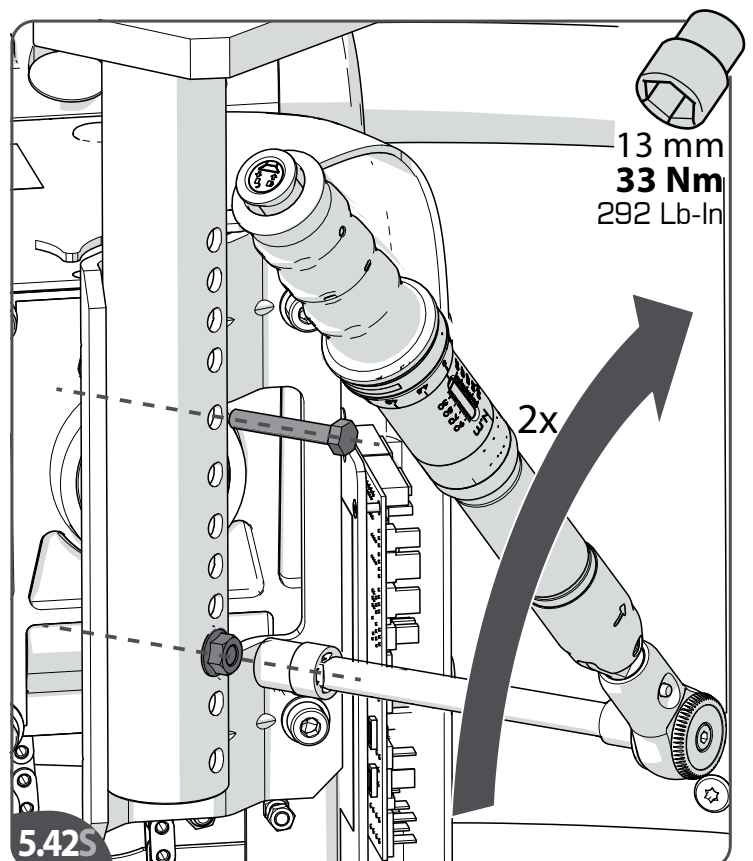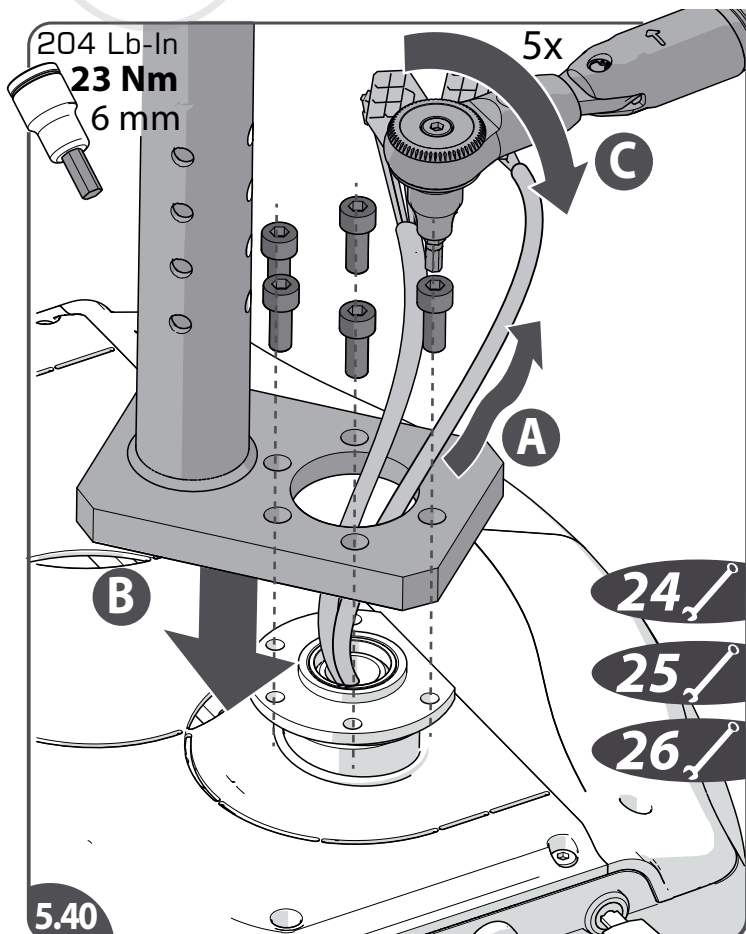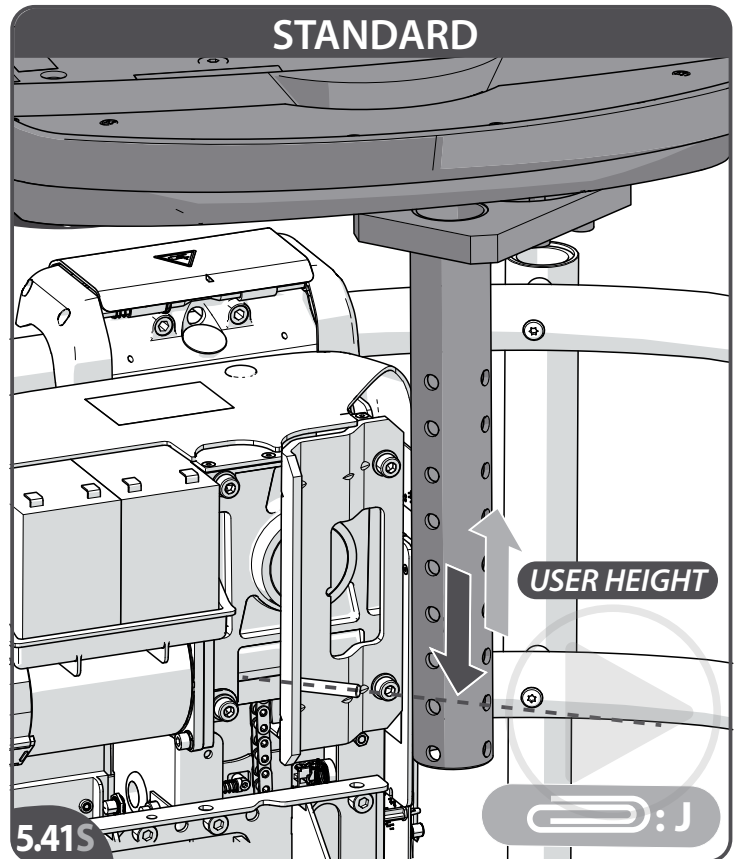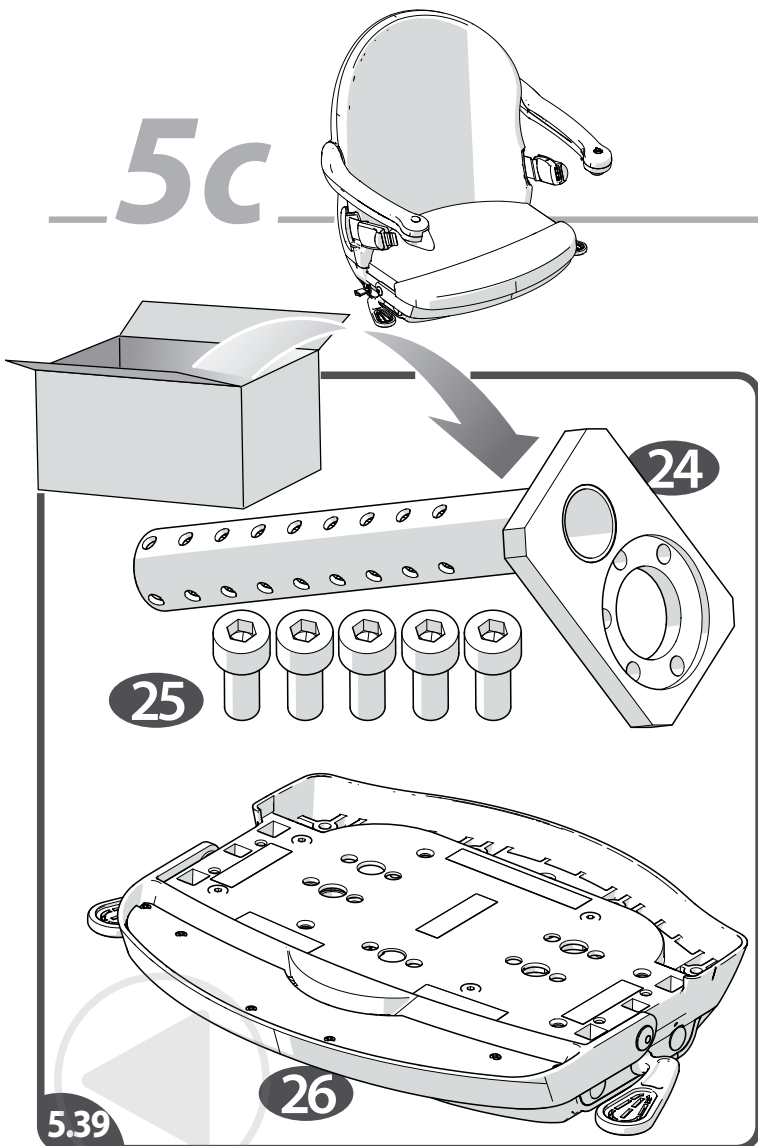- [5.39] Take the seat post out of the footplate box. Carry the post to the power pack. Check the configuration.
 [5.40] Guide the looms through the hole of the stem{A}. Position the stem{B}. Fasten the stem with the supplied bolts (23 Nm | 204 Nm){C}.
 [5.41S] Position the seat post: fix the seat to the seat height required for the primary user of the stairlift.
 [5.42S] Fasten the seat post [33 Nm | 292 Lb-In].

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

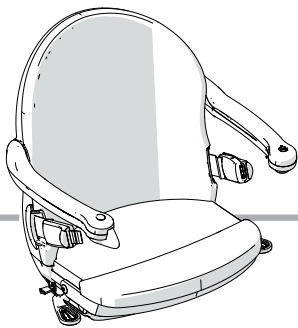
SEAT STEM HEAVY DUTY

HEAVY DUTY

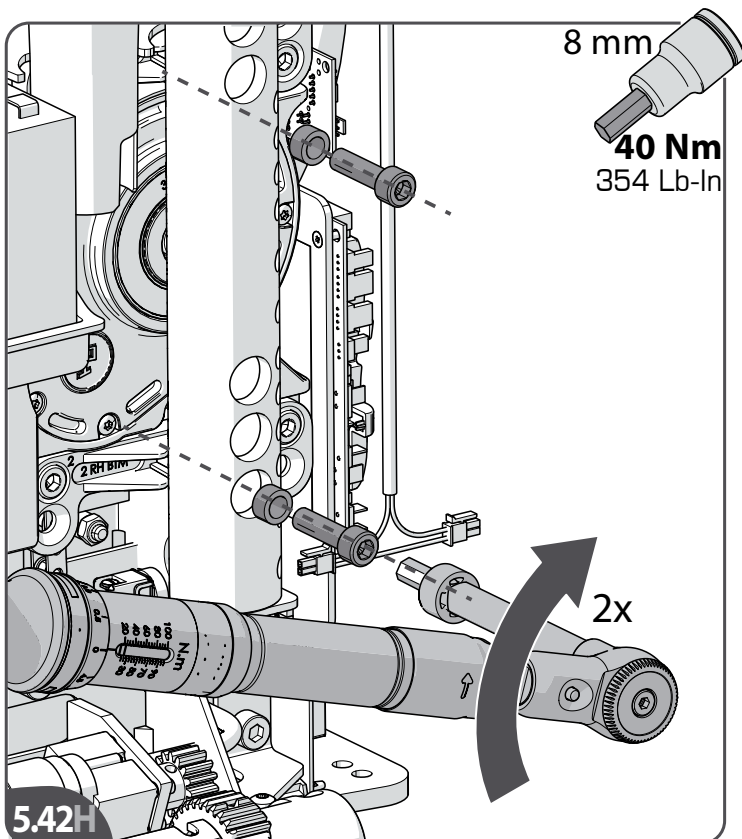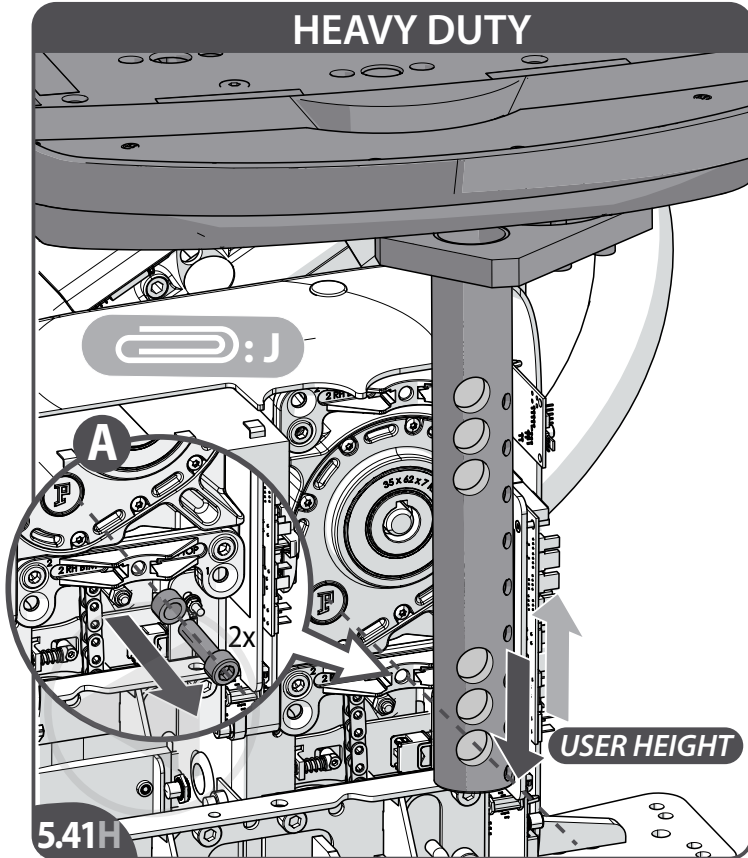[5.41H]Remove the bolts (reuse them){A}. Position the seat post: fix the seat to the seat height required for the primary user of the stairlift.

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



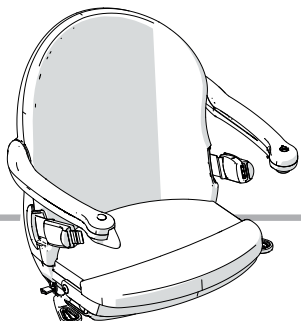
SMART SEAT

FIXATION SEAT BACK

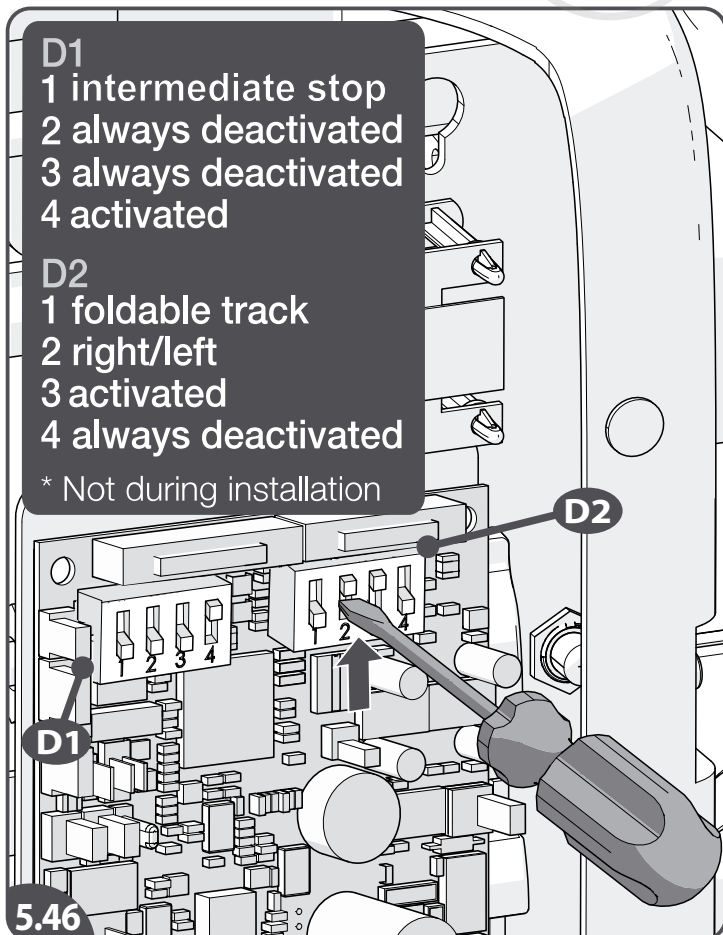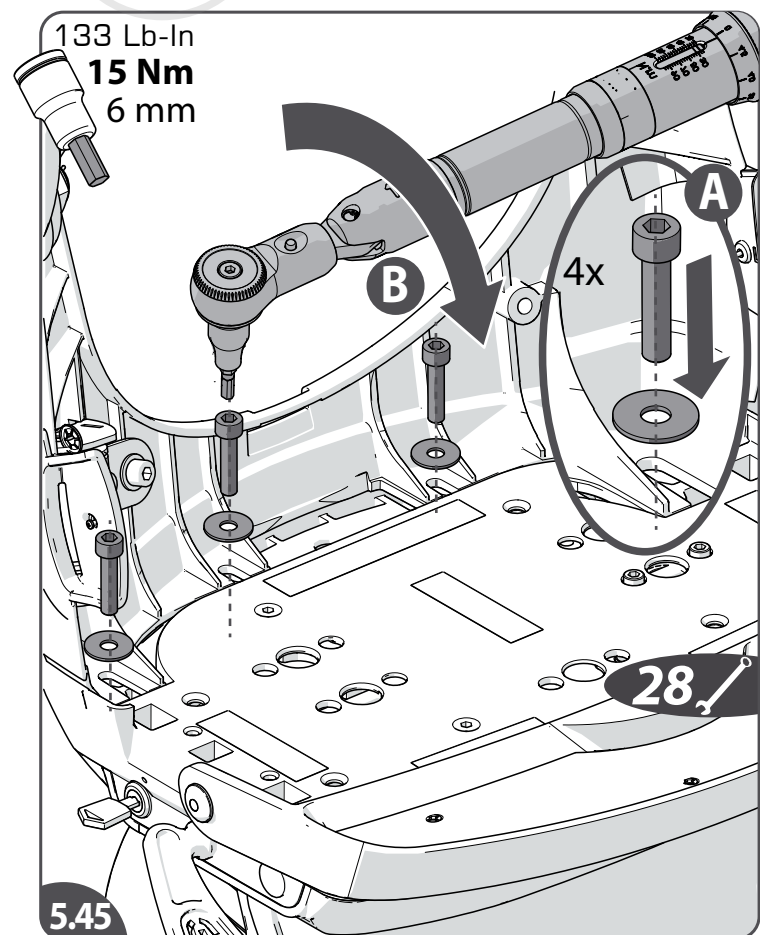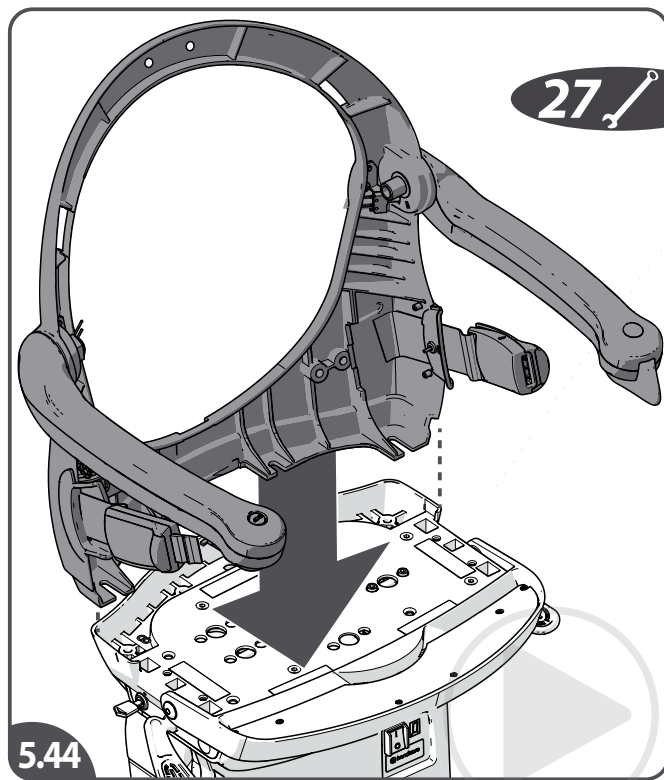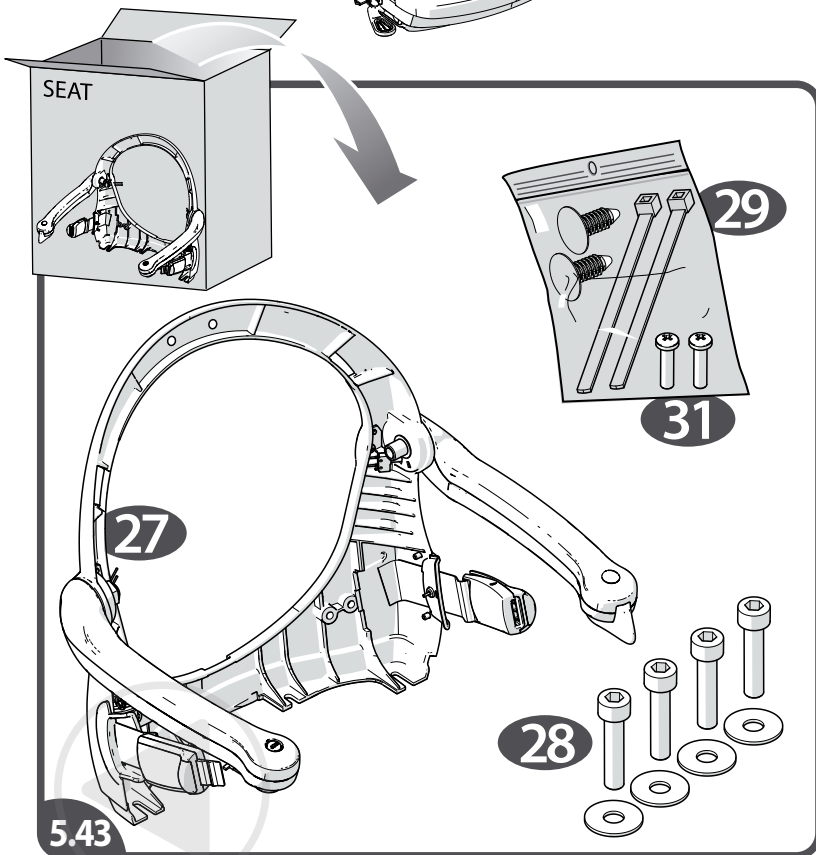[5.43] Remove the seat back components from the box.

[5.44] Position the seat back on the seat base.

[5.45] Install the washers+bolts(4x) {A} and fasten the bolts [15 Nm | 233 Lb-In]{B}.

[5.46] If needed, change the settings on the dip switches on the main PCB (power supply/power switch should be off before changing): e.g. in case of a hinged rail(D2-1) or in case of a right hand configuration(D2-2). D2-3 should be always on. Look at the inside of the front cover for the settings (or attachment G).

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

OPTION

ONETOUCHALERT
▶ **AA30047**

See the specific manual for more information

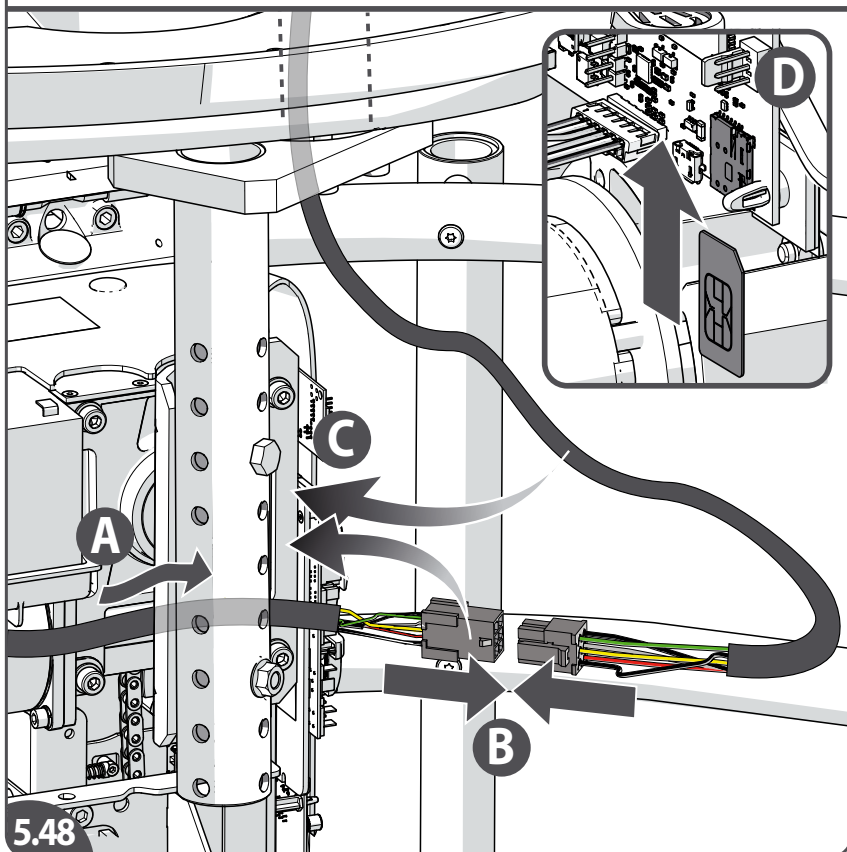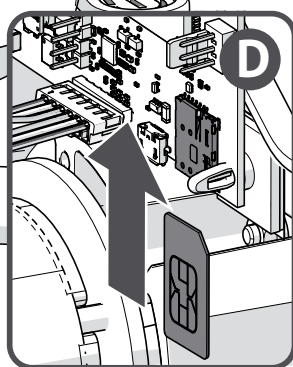5.48

[5.47] Connect the seat loom (**U6**).

[5.48] If present on the lift, guide the OTA loom {print part} behind the seat stem{A}. Connect the OTA loom parts(**U28**){B}. Tuck the excess cable+connectors behind the seat stem{C}. Please refer to instruction AA30047 for more information. Install the SIM card{D}.

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



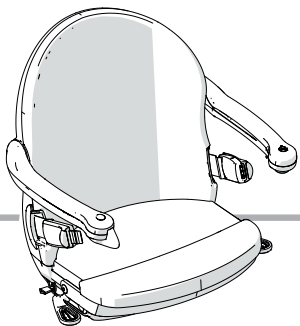
SEAT BASE WIRING

ONLY POWERED FOOTREST

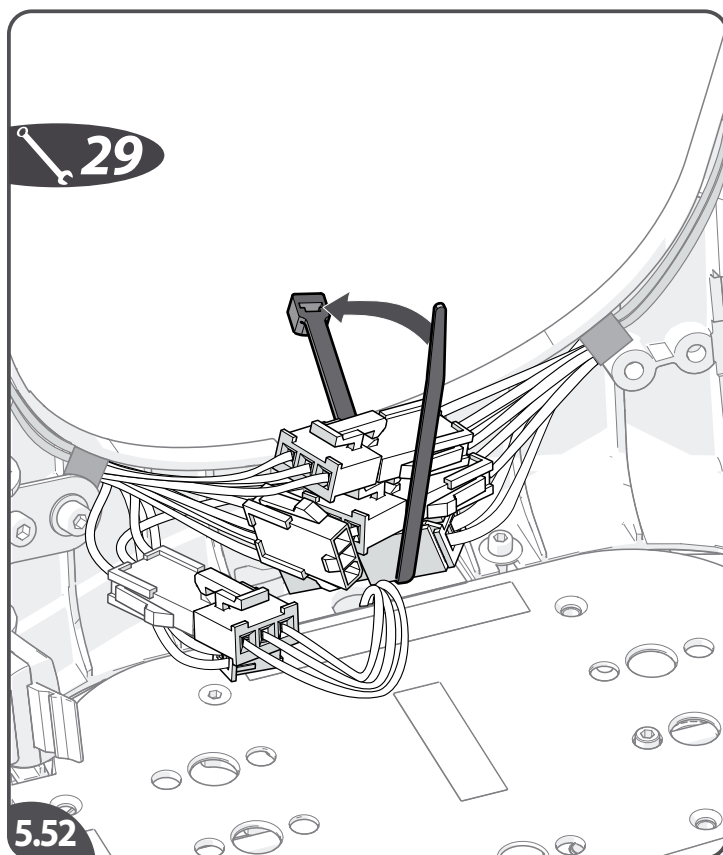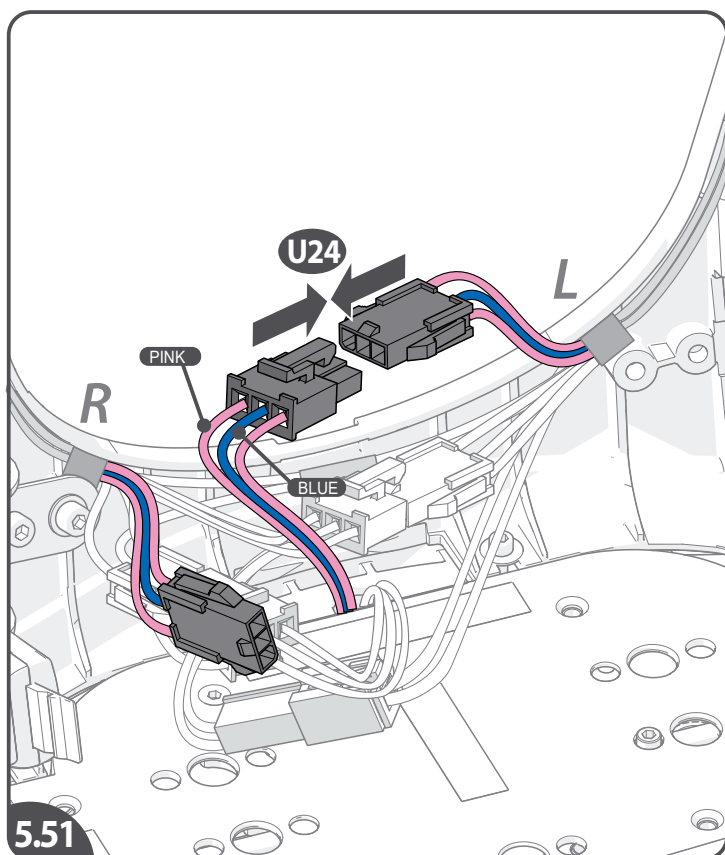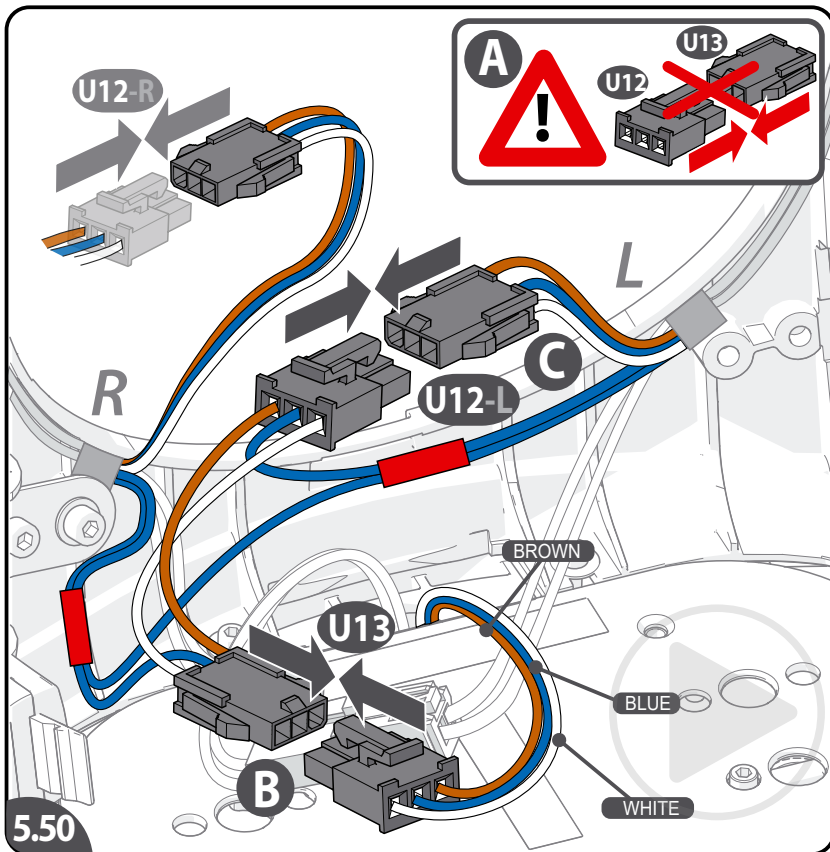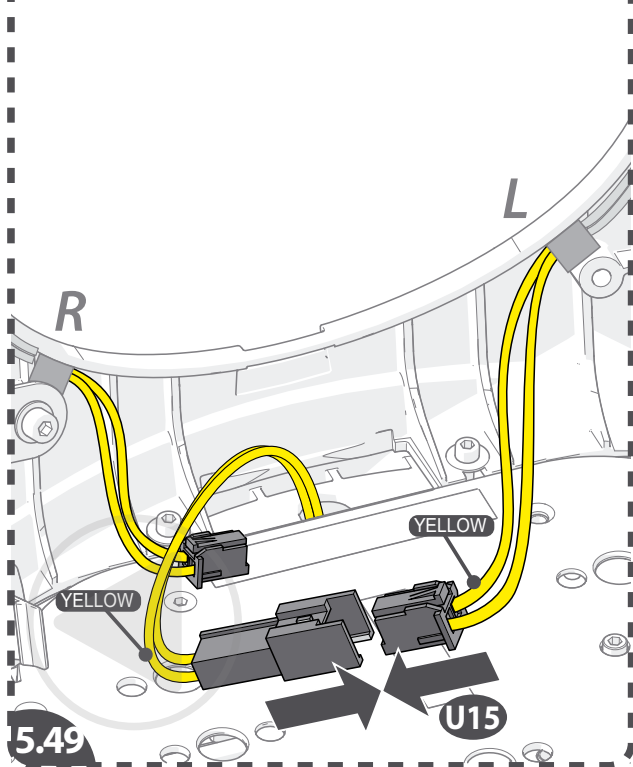- [5.49] *Only the powered footrest:* connect the footrest loom to the footrest switch {U15} (pay attention to the correct armrest!).
- [5.50] Connect the toggle loom to the armrest safety loom {U13} {B}. Do not bypass the armrest safety loom by connecting the toggle loom directly to the toggle {A}. Connect the armrest safety loom to the toggle {U12} (pay attention to the correct armrest!) {C}.
- [5.51] Connect the e-stop loom to the stopbutton {U24} (look for the correct arm rest!).
- [5.52] Tie the looms/connectors together with a tie-wrap.

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






5c

SMART SEAT UPHOLSTERY

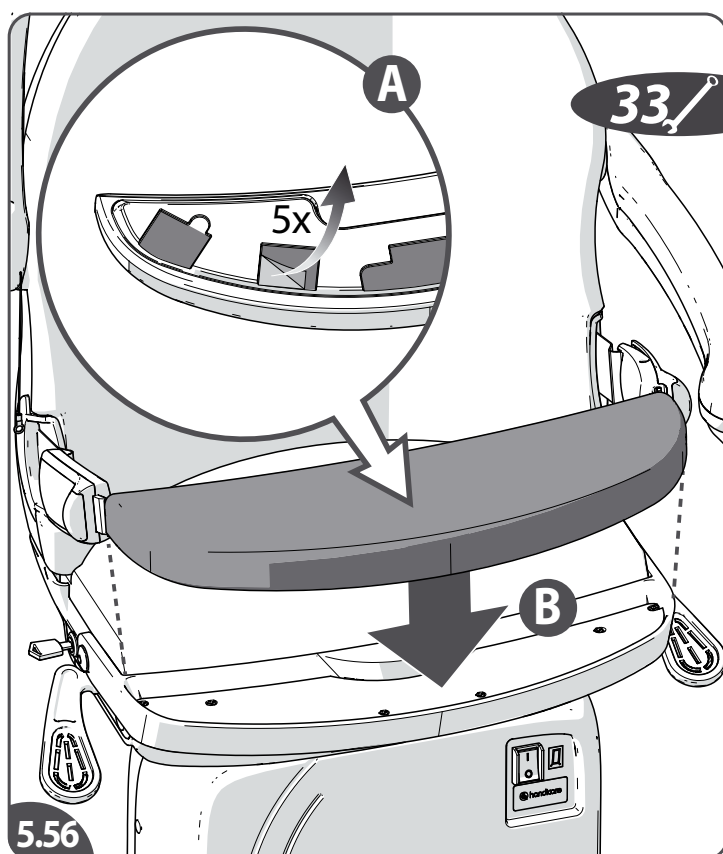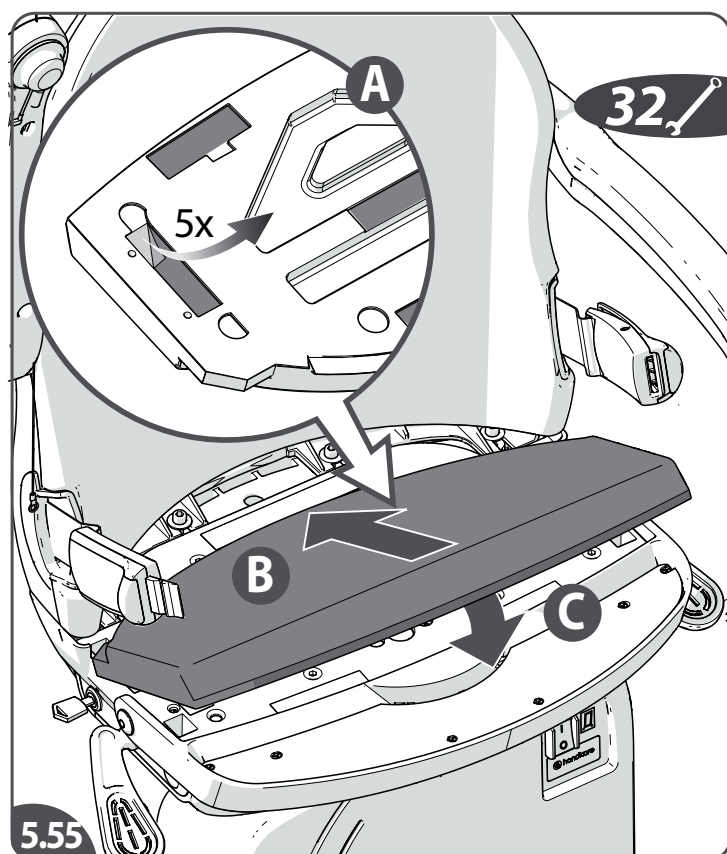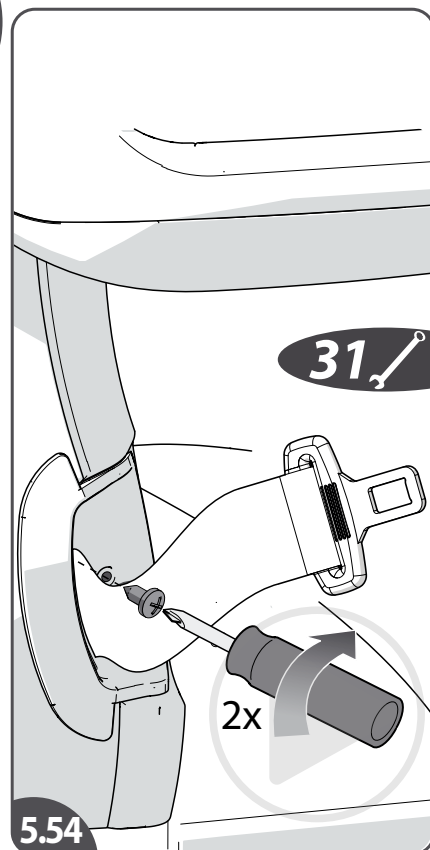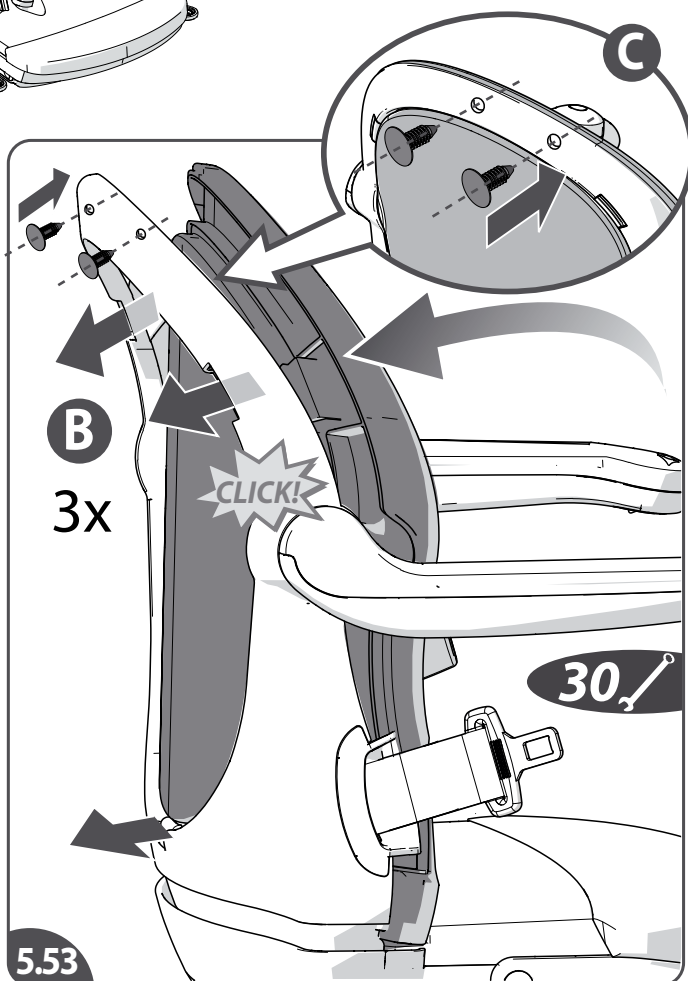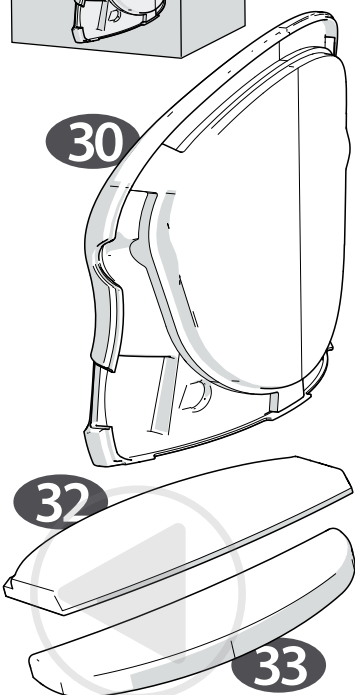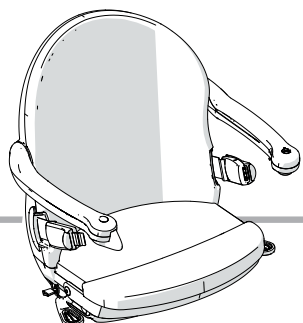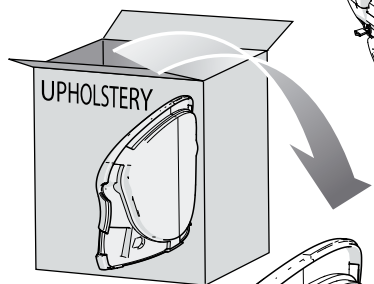Remove the upholstery components from the box.

[5.53] Install the back seat upholstery {A}. Press firmly on the locks to fix the upholstery {B}.

[5.54] Fasten the upholstery with the screws (2x).

[5.55] Remove the backings of the velcro pads (6x) {A}. Fit the main seat cushion {B+C}.

[5.56] Remove the backings of the velcro pads (5x) {A}. Fit the seat flap cushion {B}.

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







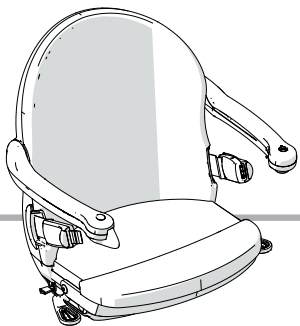
SMART SEAT

CLOSING FRONT COVER

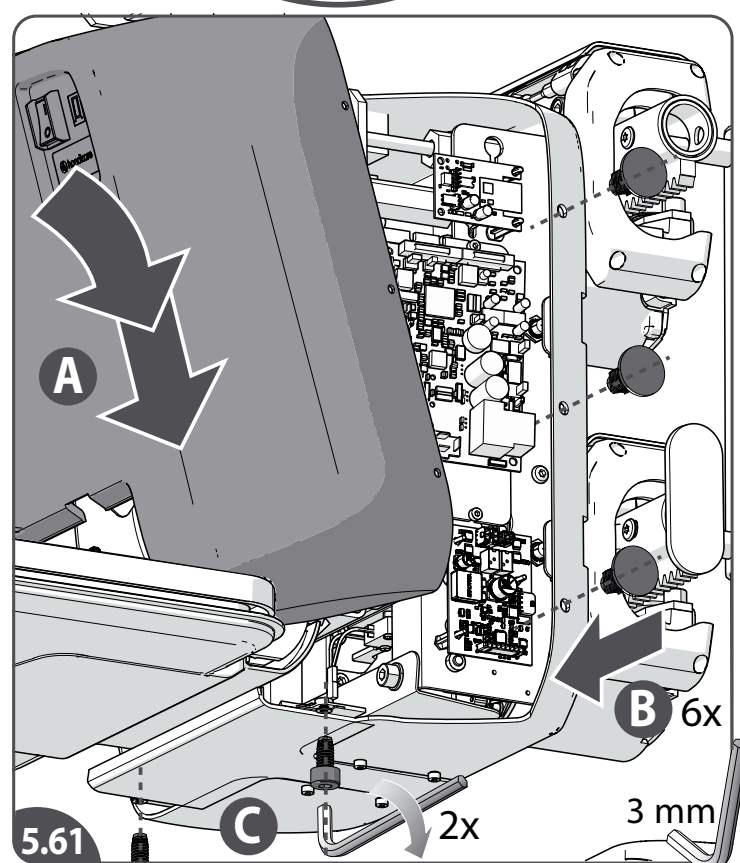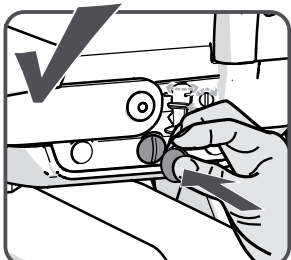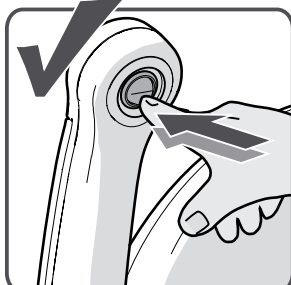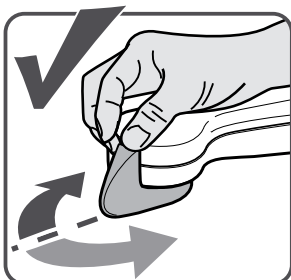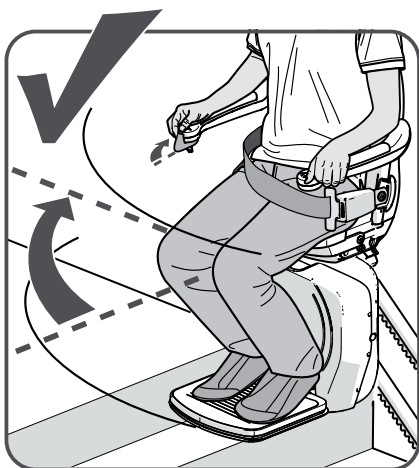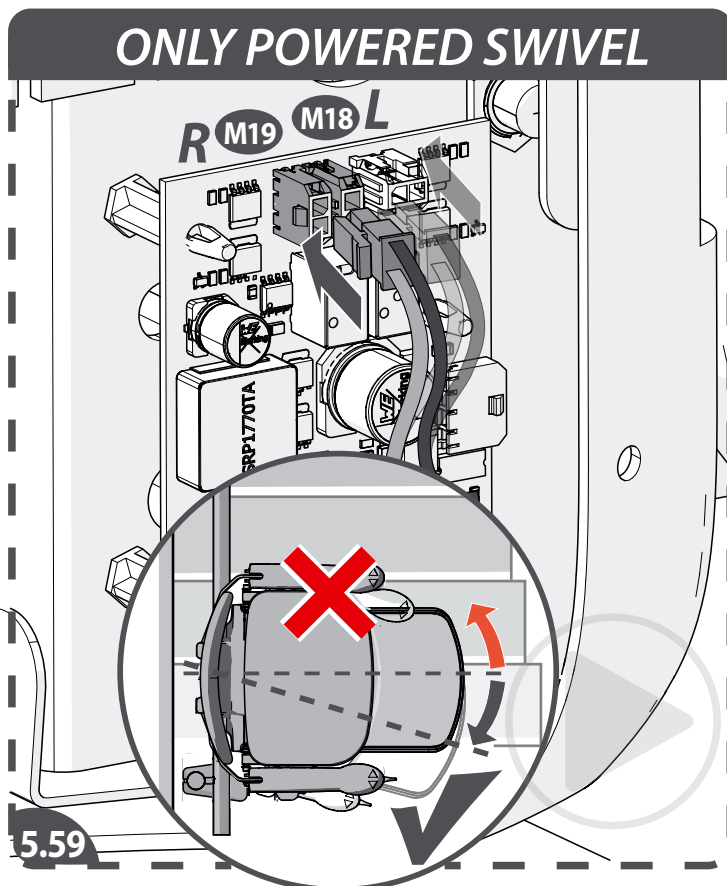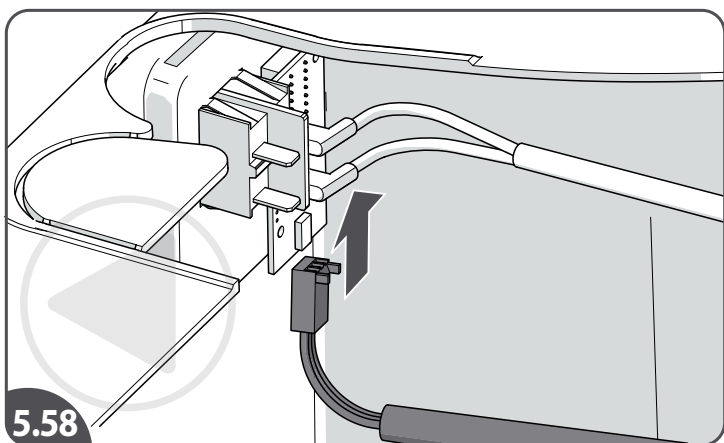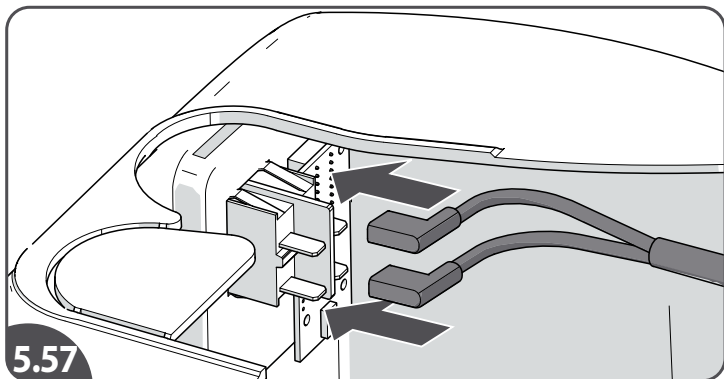- [5.57] Connect the power switch cables.
- [5.58] Connect the diagnostic display loom.
- [5.59] Connect the wires of the swivel motor.
- [5.60] Check the functions and safeties before installing the cover.
- [5.61] Reinstall the front cover{A} and secure with the rivets(6x){B} and screws(2x){C}.

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



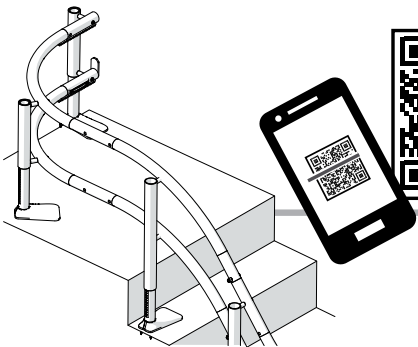
TRACK FIXATION CHECK

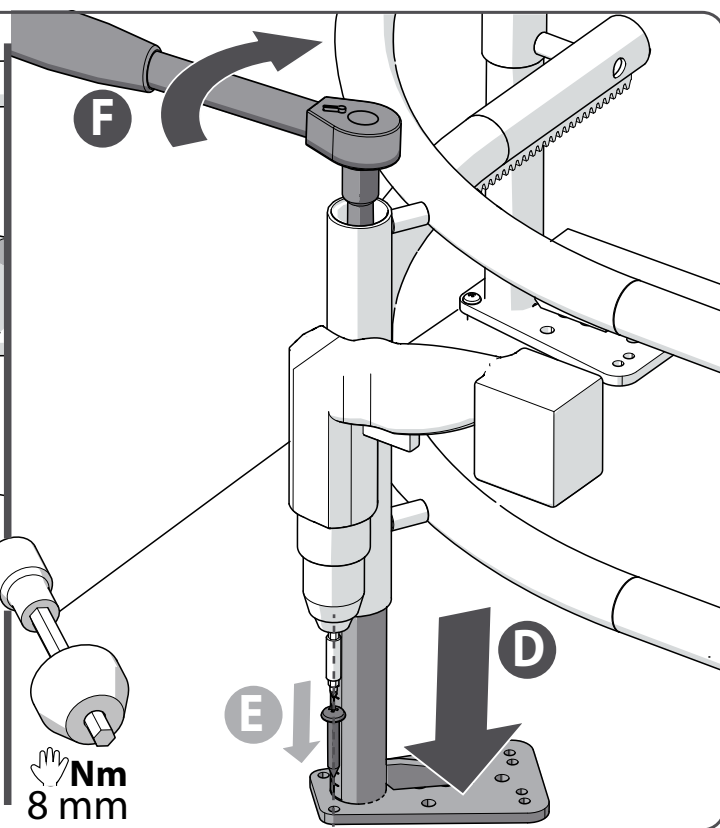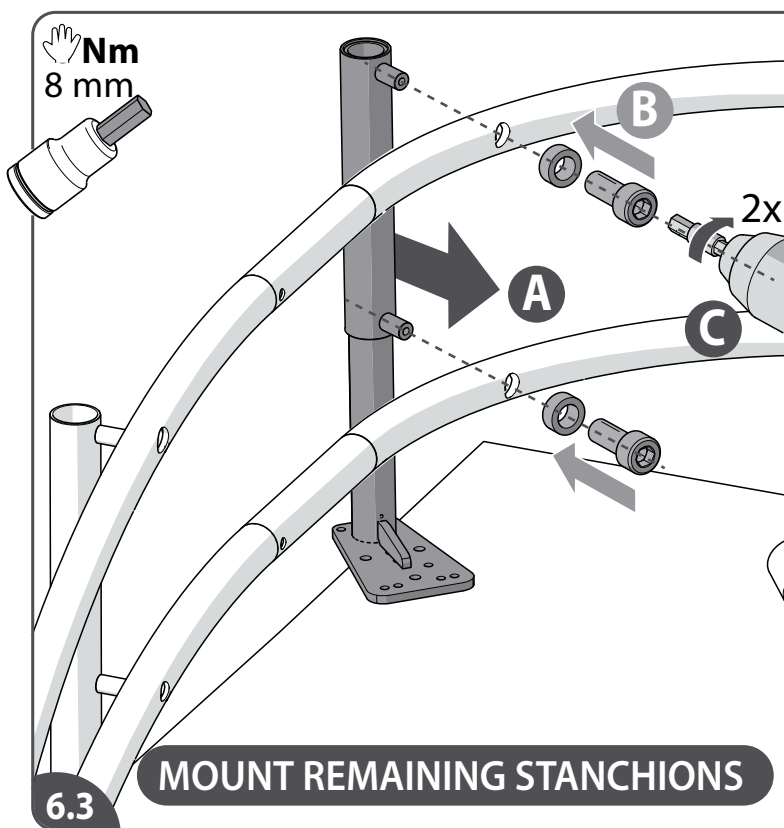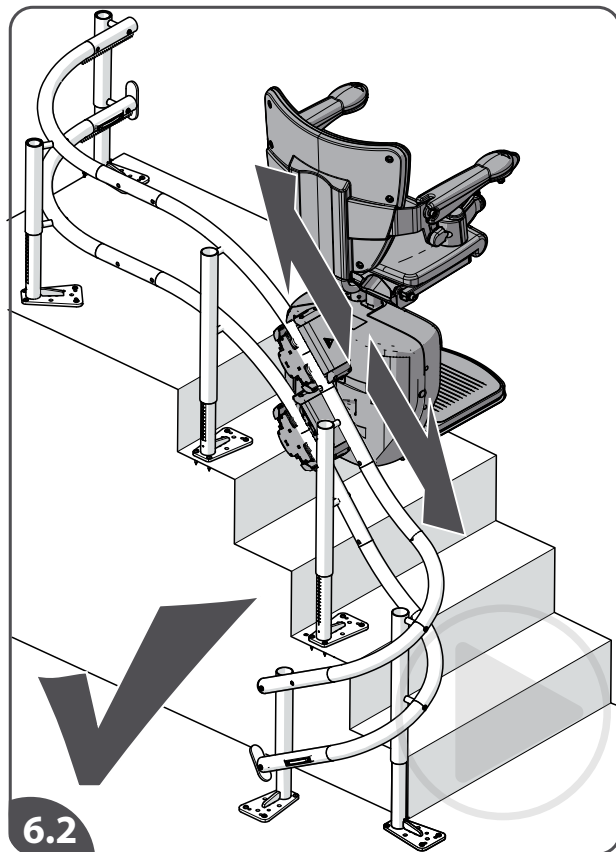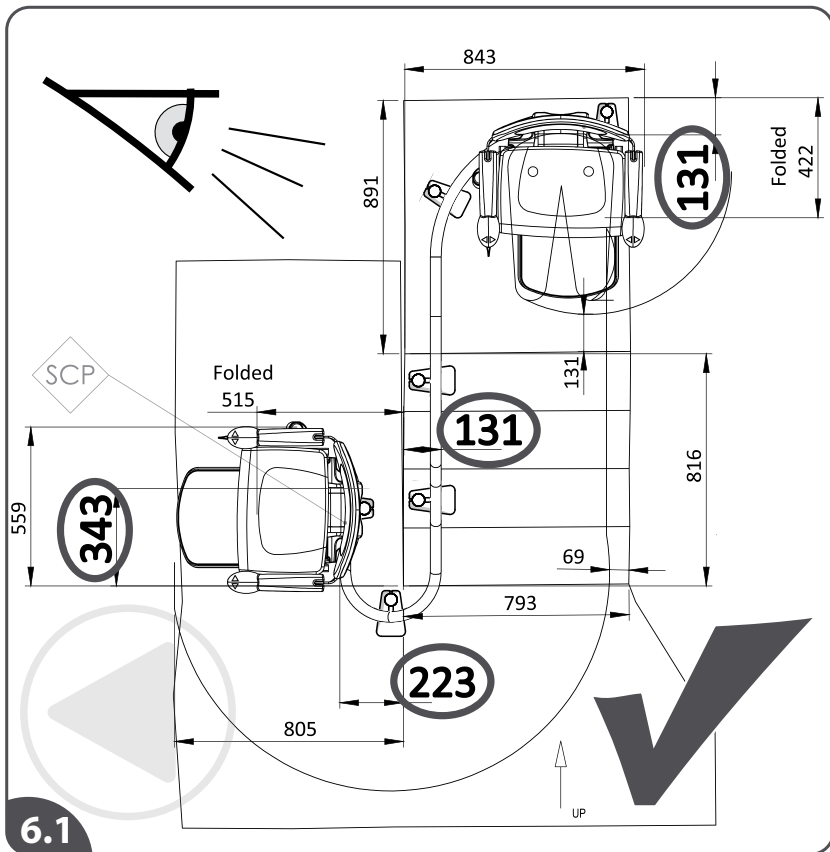- [6.1] Final check for obstacles along the track and a final check for the track offsets (mentioned on the drawing).
- [6.2] Final check track position: drive the lift down/up.
- [6.3] Install the remaining stanchions.

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



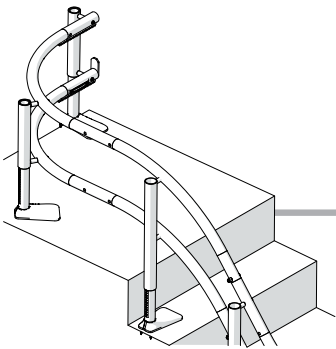
TRACK FIXATION

TIGHTEN BOLTS

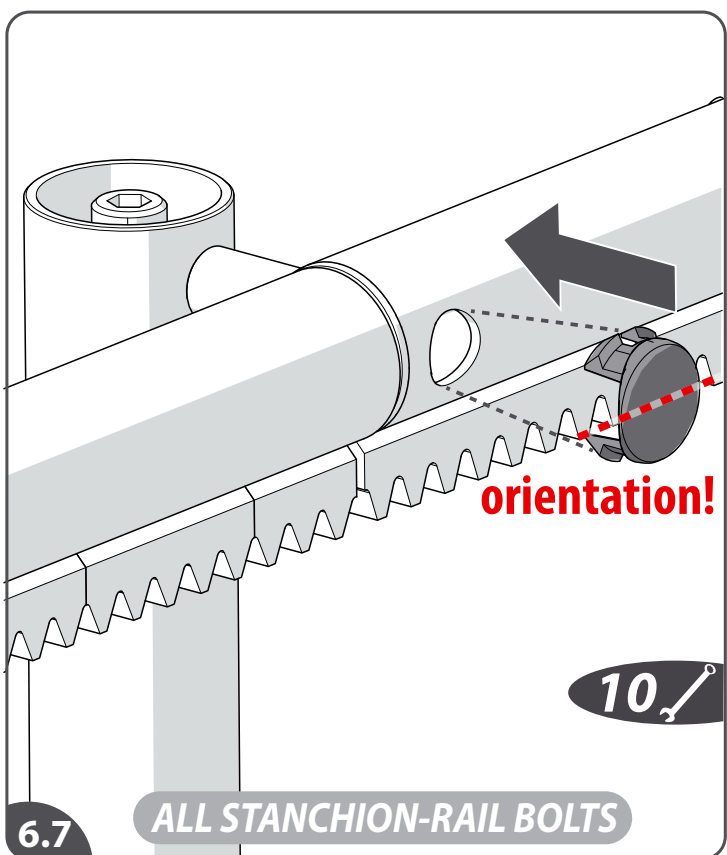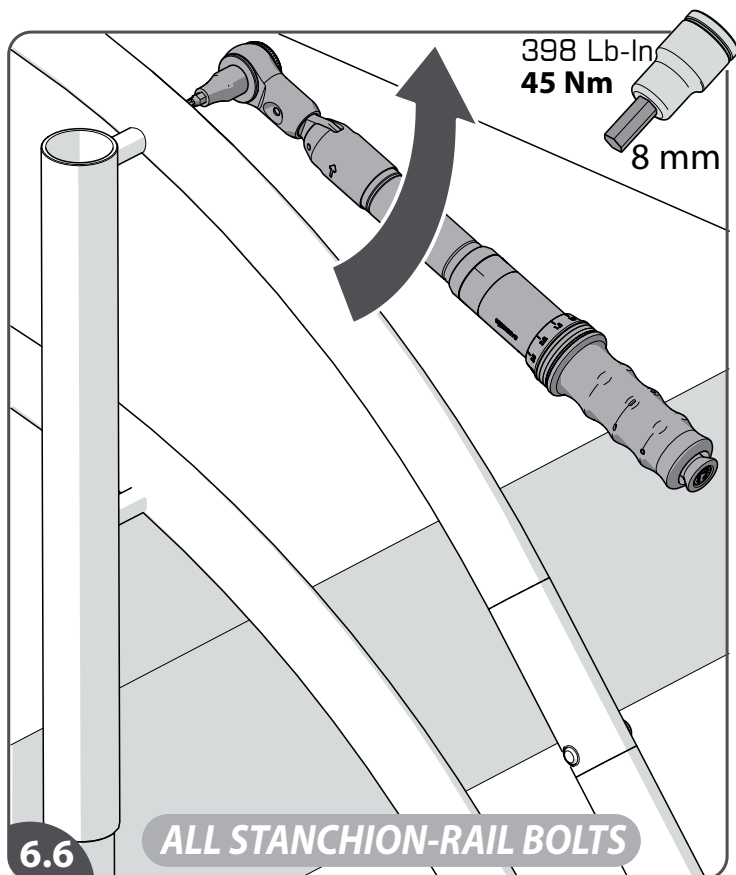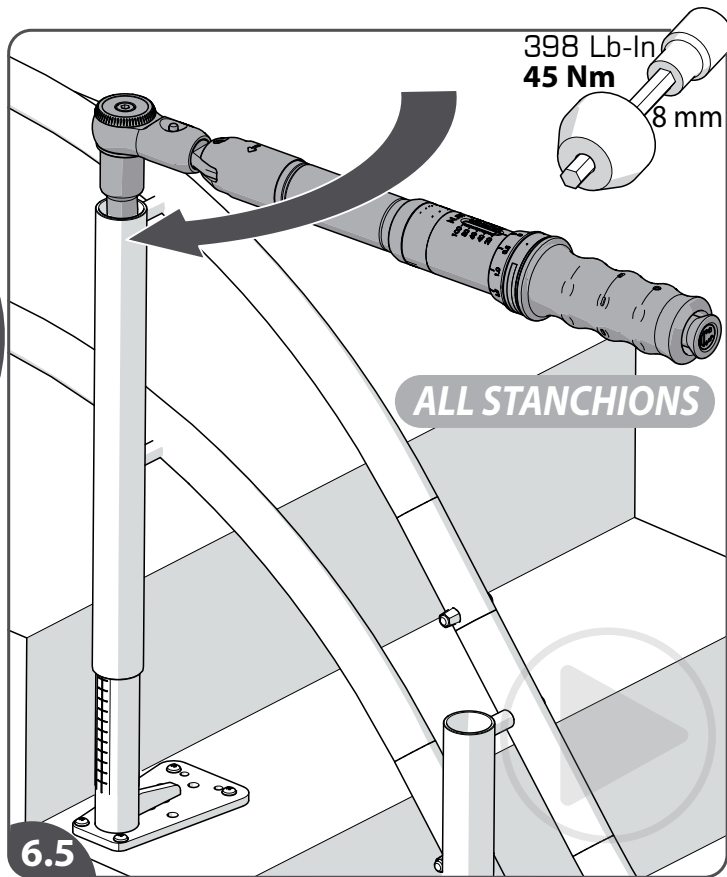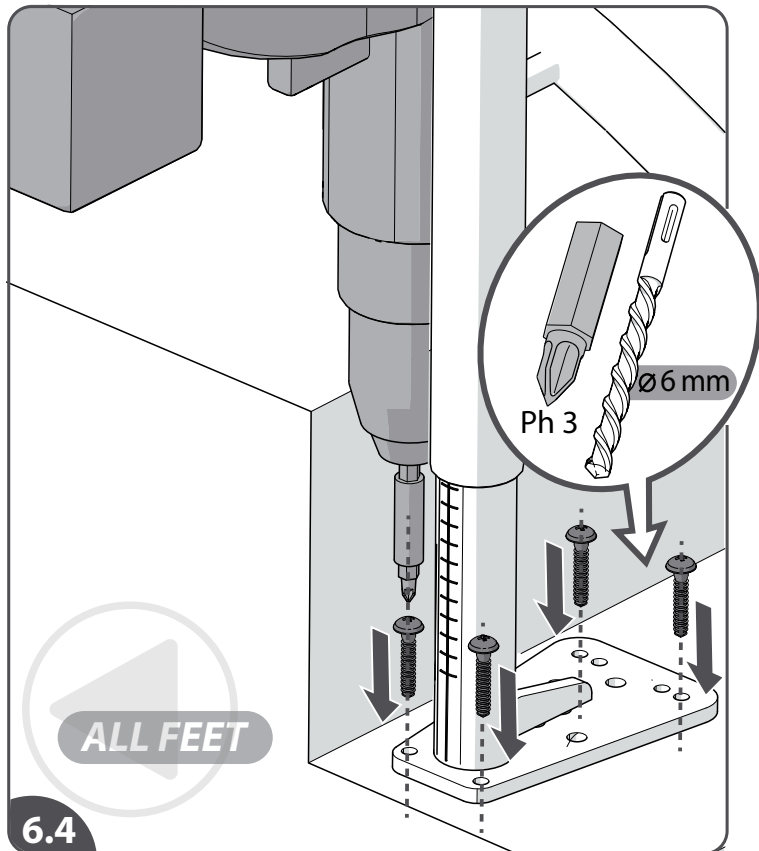- [6.4] Fasten all the feet to the steps/floor using screws (for wooden stairs) or rawl plugs/screws (for concrete/stone stairs). Make sure there is no space between the foot and the step before fixing the rail.
- [6.5] Torque up all the lock bolts [45 Nm | 398 Lb-In].
- [6.6] Torque up the stanchion-rail joint bolts [45 Nm | 398 Lb-In].
- [6.7] Install the caps of the rail.

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



6

TRACK FIXATION WALL BRACKETS

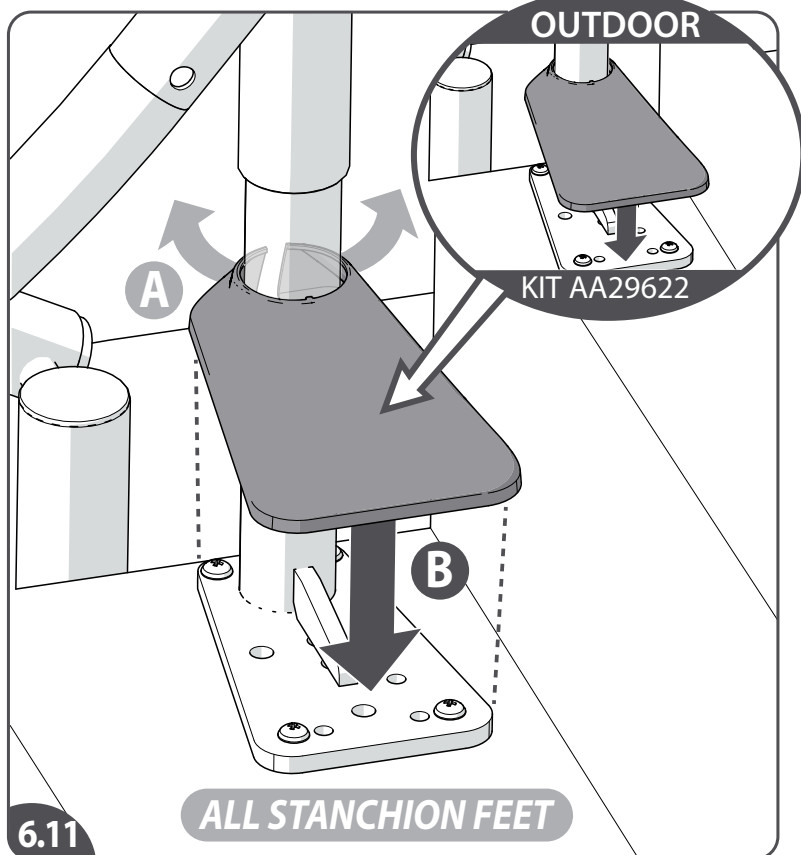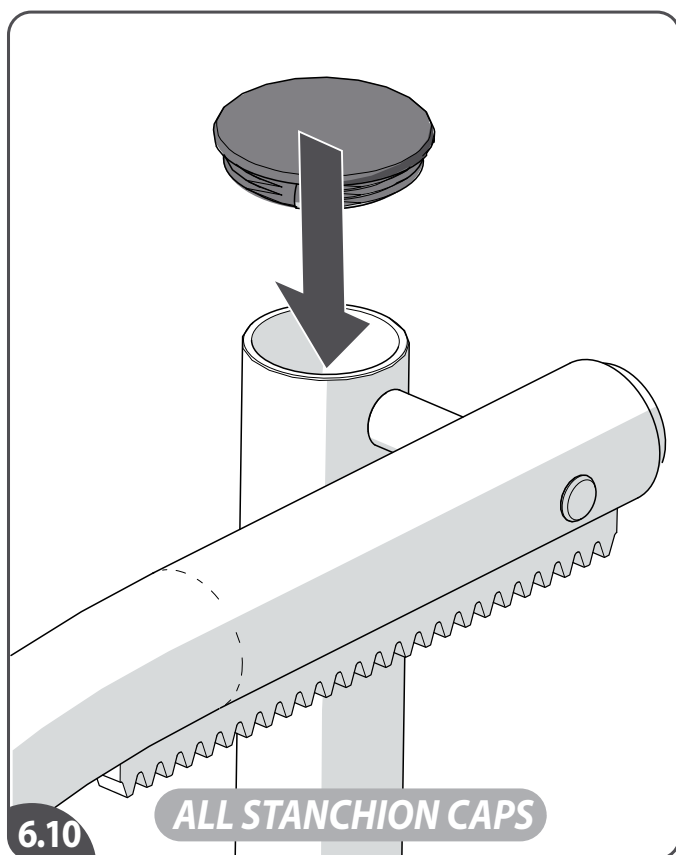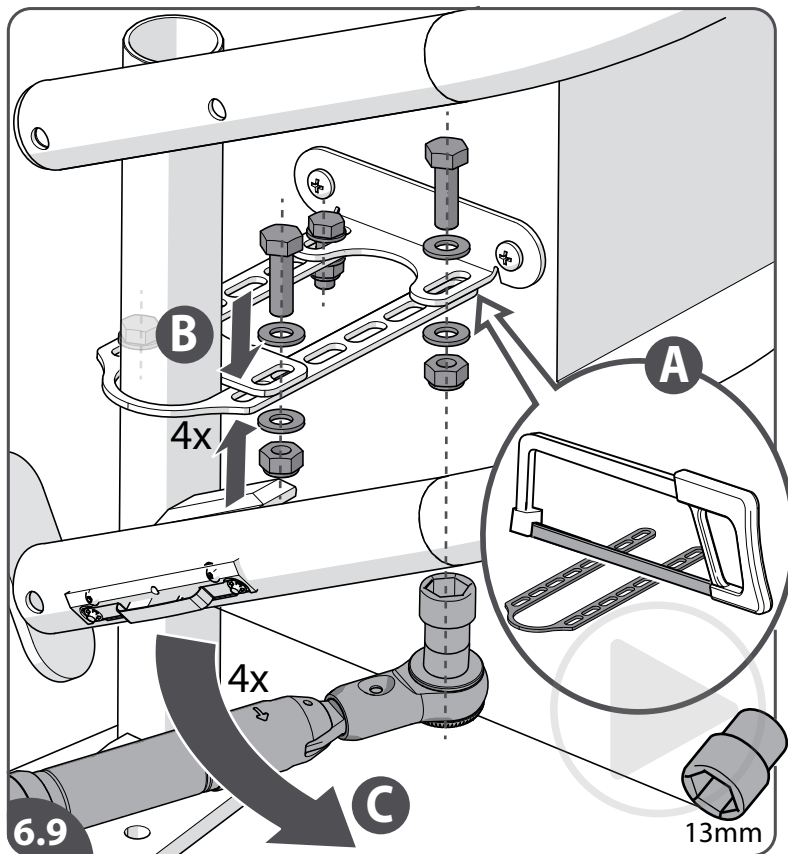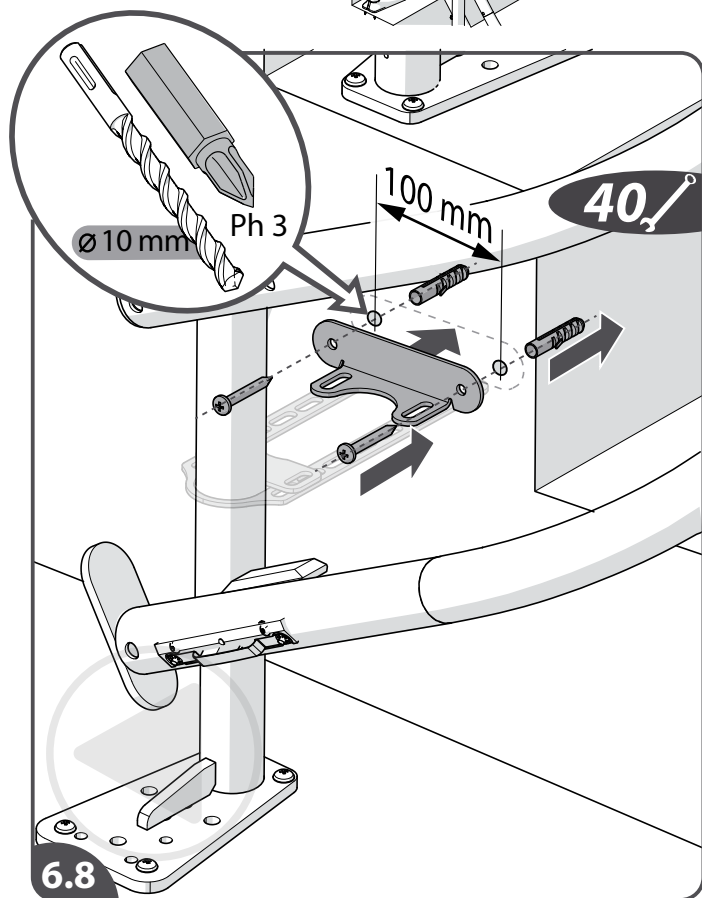[6.8] Wall Brackets are available to use for additional support for a stanchion that is fixed directly into the floor. They are not designed to withstand the forces of the loaded stairlift on their own and so can only be used as secondary fixings. They should be attached to a structural wall.

To install the wall brackets; position the bracket in-between the two stanchion lugs. position the bracket in-between the two stanchion lugs. Secure the bracket to a structural wall using screws or rawl plugs/screws (for concrete/stone walls).

[6.9] Shorten the brackets if necessary{A}. Fasten the bracket parts securely with the supplied bolts/nuts{B}.

[6.10] Install the caps of the stanchions.

EN [6.11] Install the covers of the stanchions. In case of an OUTDOOR installation use the UV-stabilized covers.

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

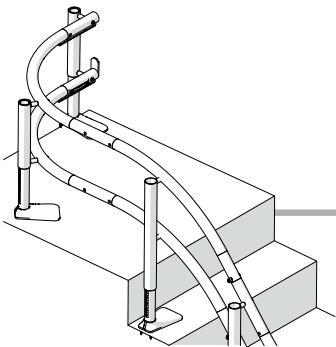
TRACK FIXATION

APPLY GREASE

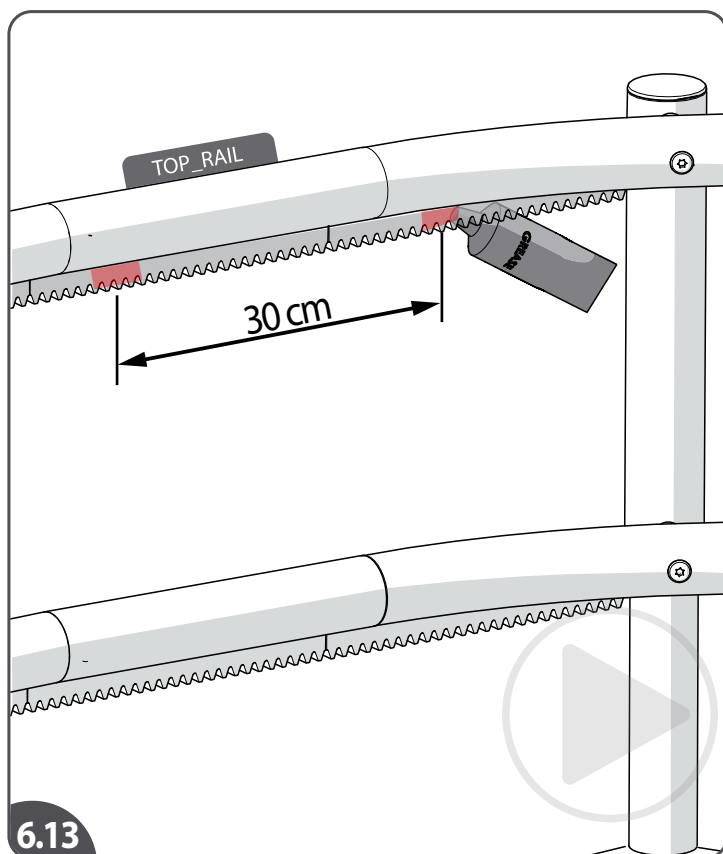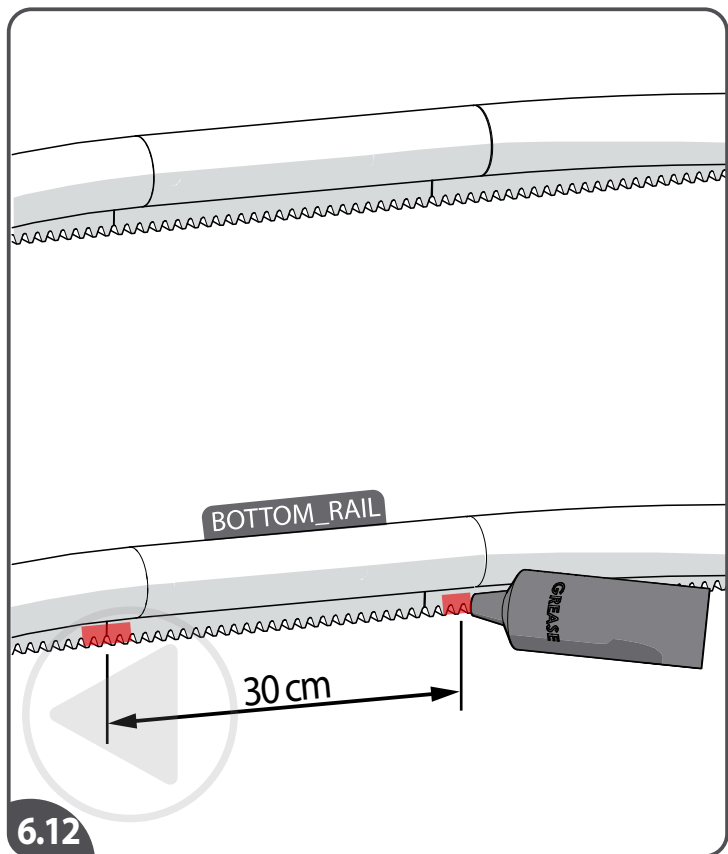[6.12]Apply the grease to the bottom track racks; apply to the front face, intermittently on every 30 cm.
[6.13]Apply the grease to the top track racks; apply to the back face, intermittently on every 30 cm.

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




REMOTE CONTROL MOUNTING

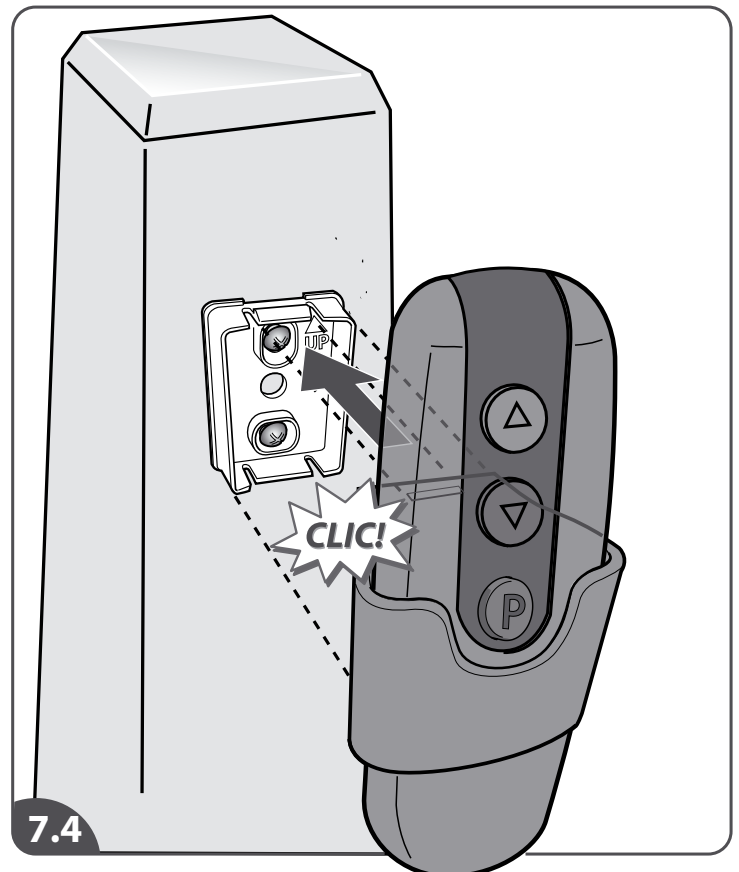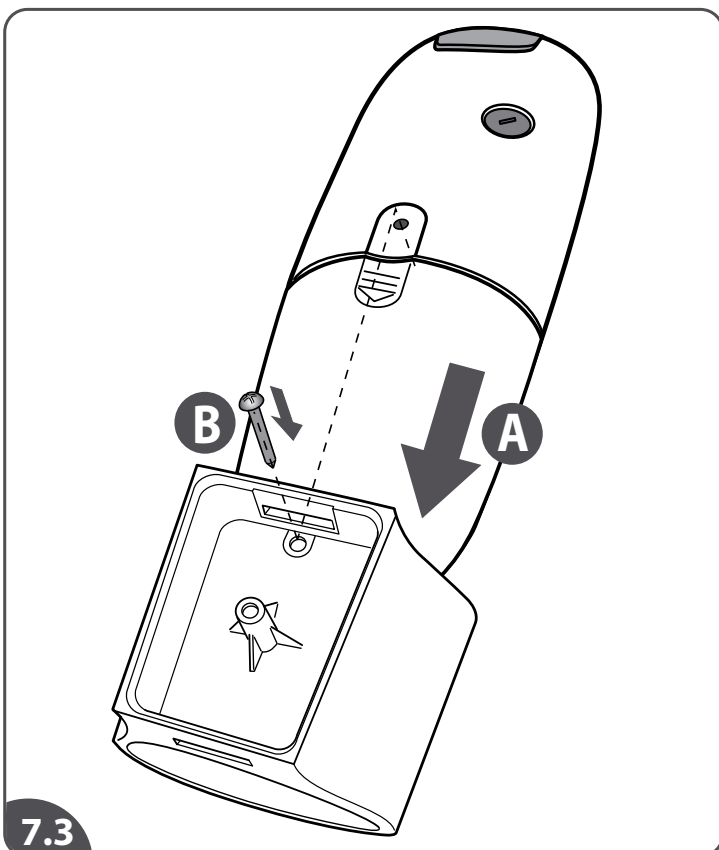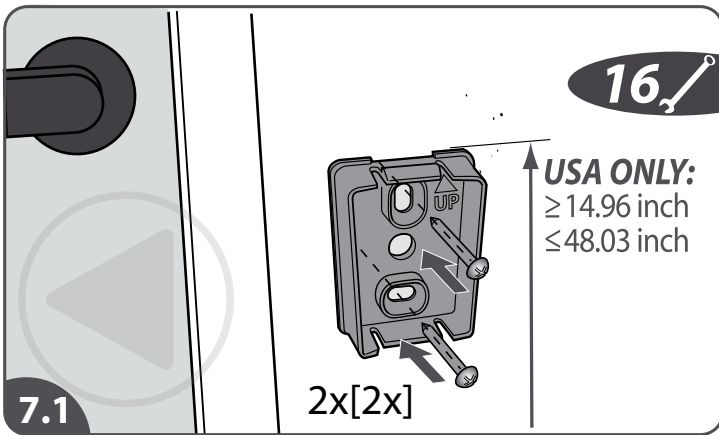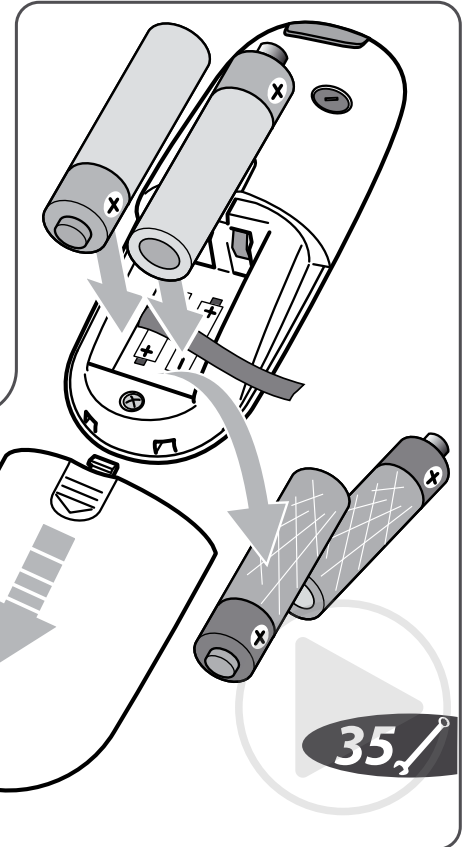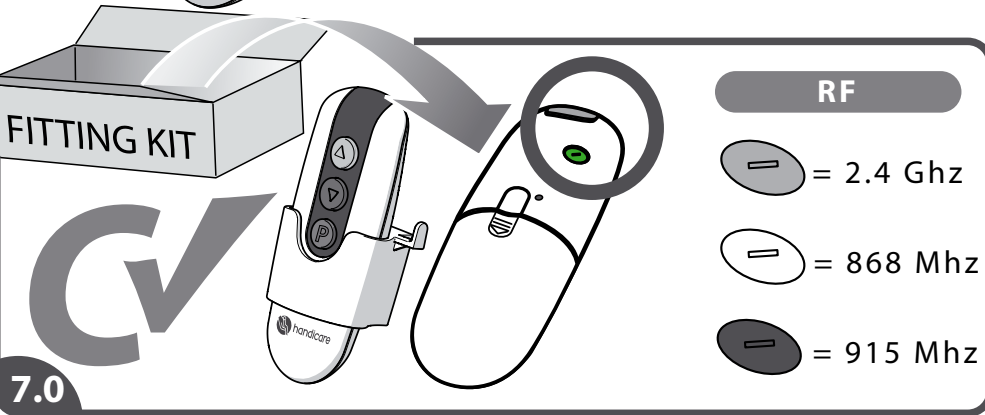[7.0] Check the amount and type of remote controls. The cap on the back of the remote control indicates the type: Grey=RF{**2.4Ghz**}, White=RF{**868Mhz**}, Black=RF{**915Mhz**}. The supplied remote controls should work out of the box. If not or in case of multiple lift, they can be reprogrammed (see manual, page 110).

[7.1] Measure the position for the remote control bracket. Mount the bracket of the remote control with two screws on a solid background. There must be a minimum of 2 remote controls, with at least one at the top and bottom of the installation.

[7.2] Insert the batteries into the remote control.

[7.3] Mount the remote control into the holster{A} and fix it with a screw{B}.

EN [7.4] Click the holster/remote control over the wall plate. Leave the remote control in the holster so that it does not get lost.

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


REMOTE CONTROL LINKING

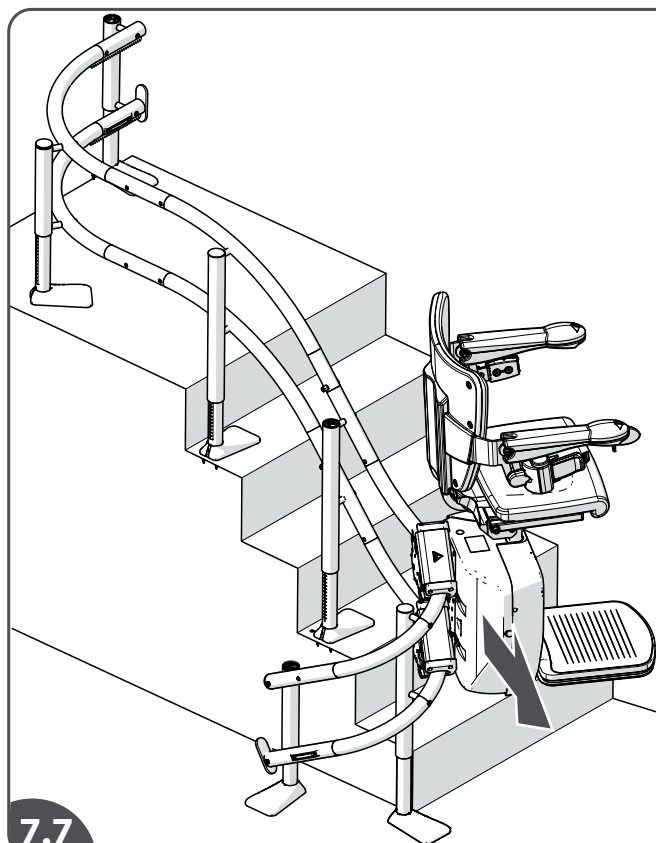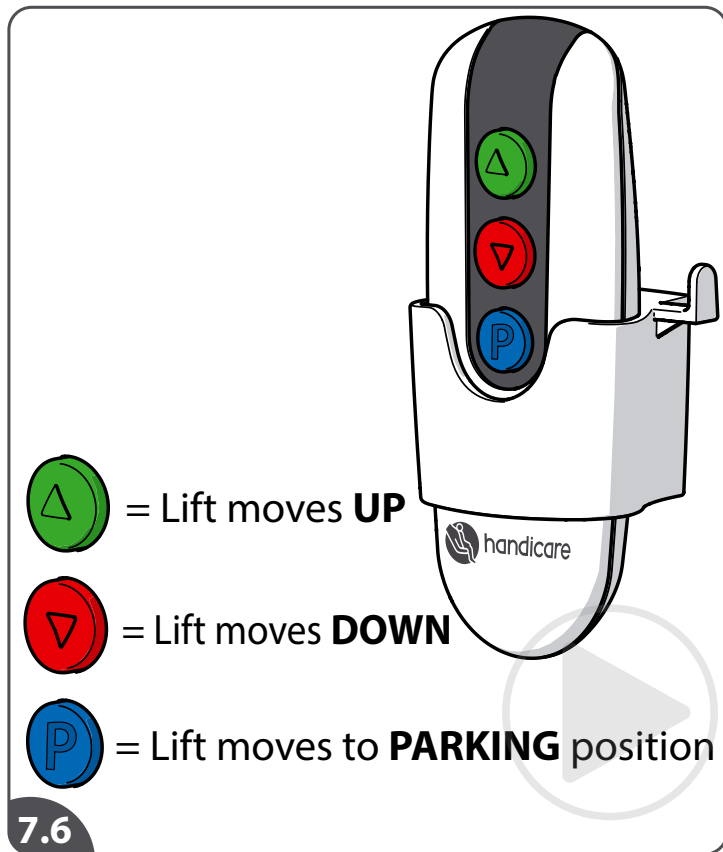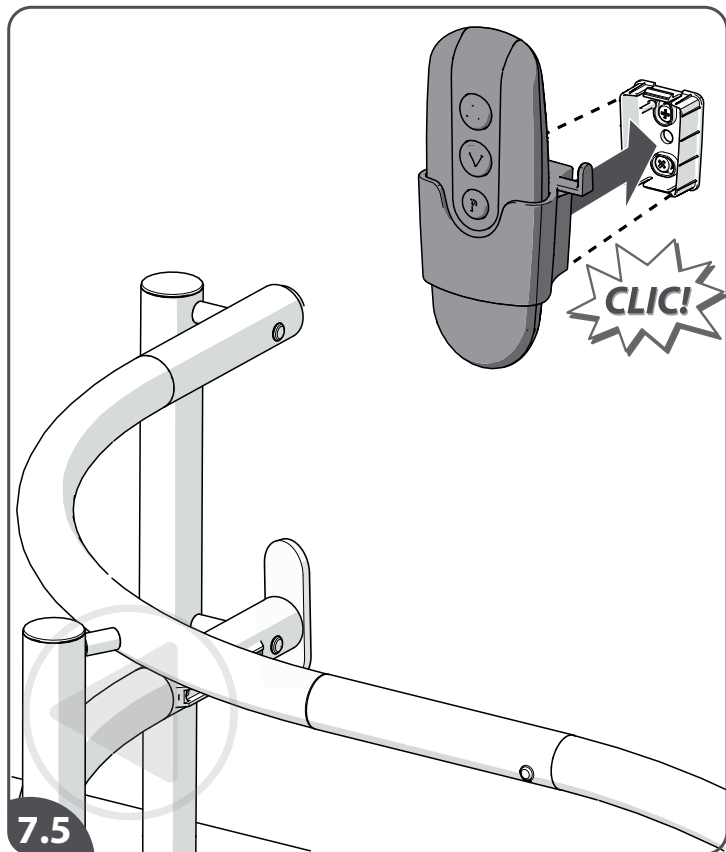[7.5] Repeat this procedure on other floor(s).

[7.6] Function of the buttons: green = **UP**, red = **DOWN**, blue = **PARK**.

[7.7] Press and hold the button (only press one button at a time) of the remote control to check if the remote control is linked with the lift. The remote control will show a flashing green light (searching). Once this signal has connected with your lift the green light will stop flashing and stay lit while your lift is moving.

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



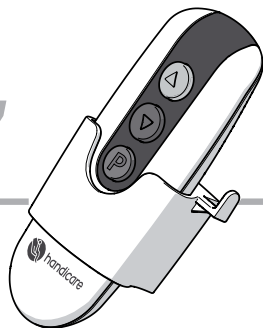
REMOTE CONTROL REPEATER

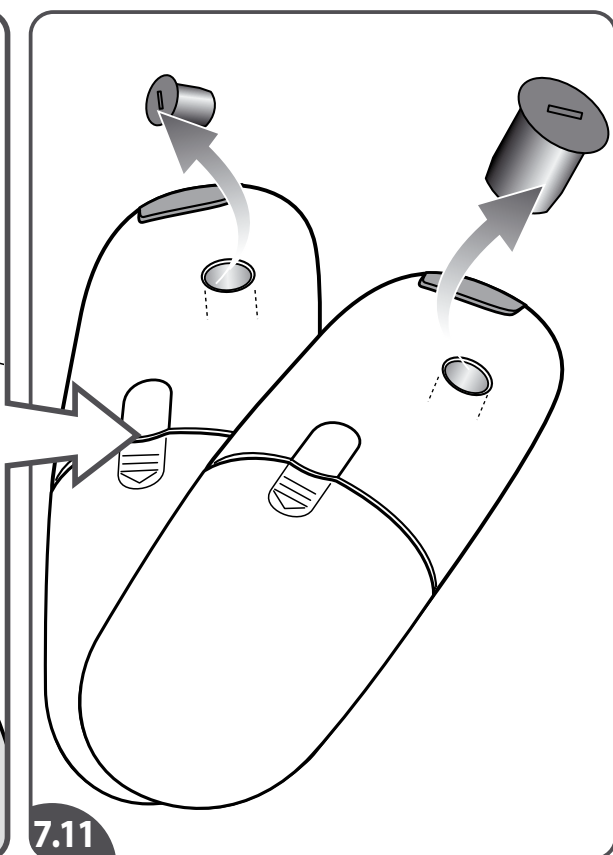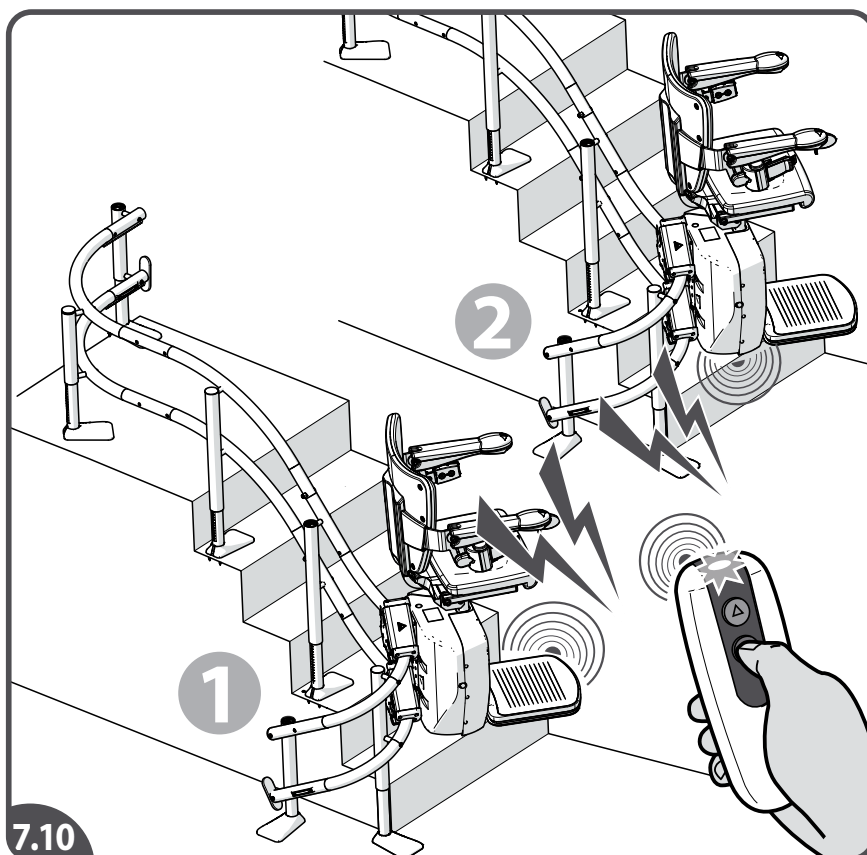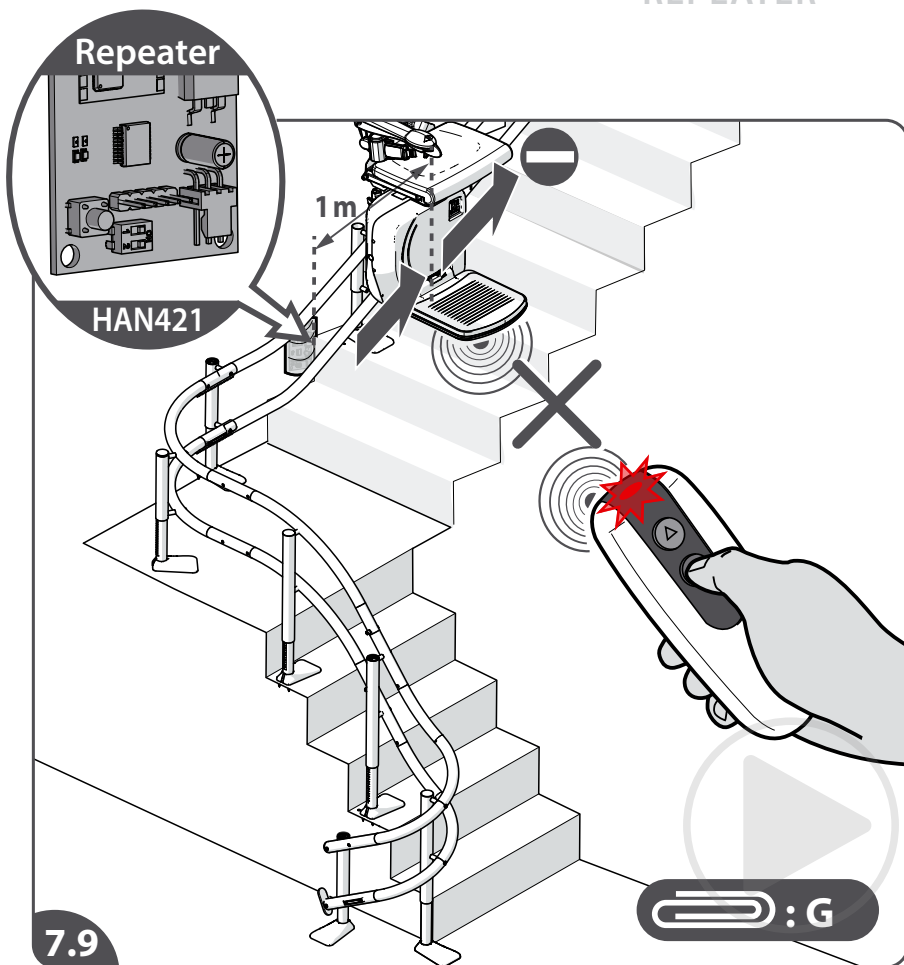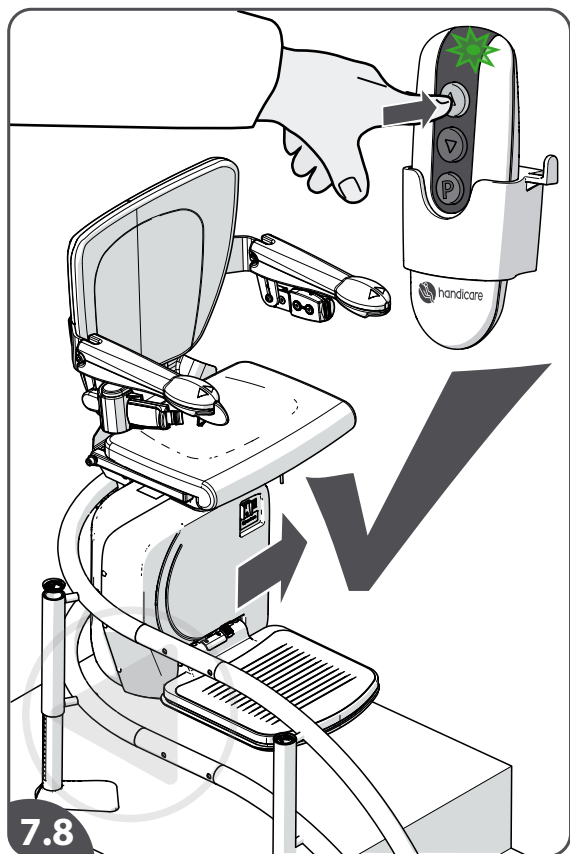[7.8] Check the other remote controls.

[7.9] In case the remote controls are out of reach install a repeater (**HAN421**) in between. See attachment G.

[7.10] In case there is interference in the radiosignal, change the channel of the radiosignal of the lift.

Note: another lift nearby can interfere with the signal.

[7.11] **IF NEEDED FOLLOW STEPS 7.11 – 7.18** > To change the channel of the radio signal it is necessary to change both the dip switches in the unit and the remote controls. Remove the caps on the backside of the remote controls.

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



REMOTE CONTROL

CHANGING CHANNELS

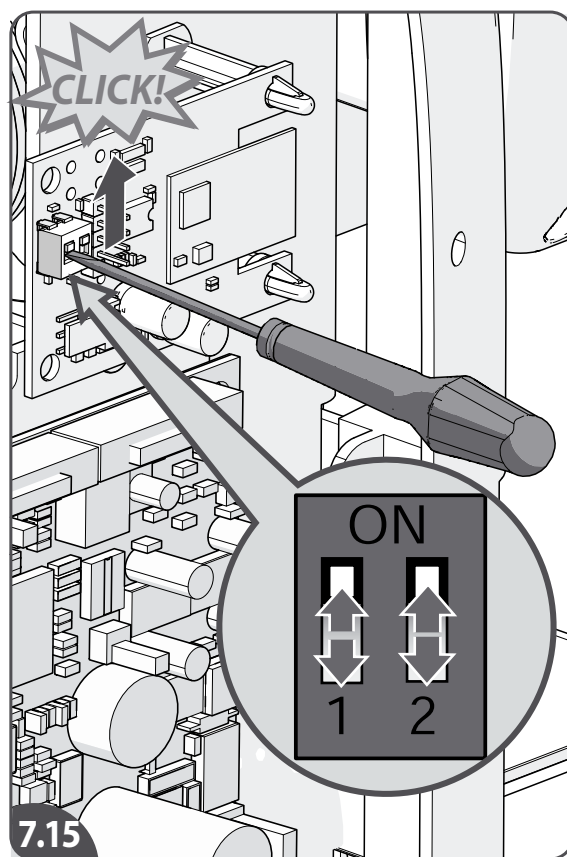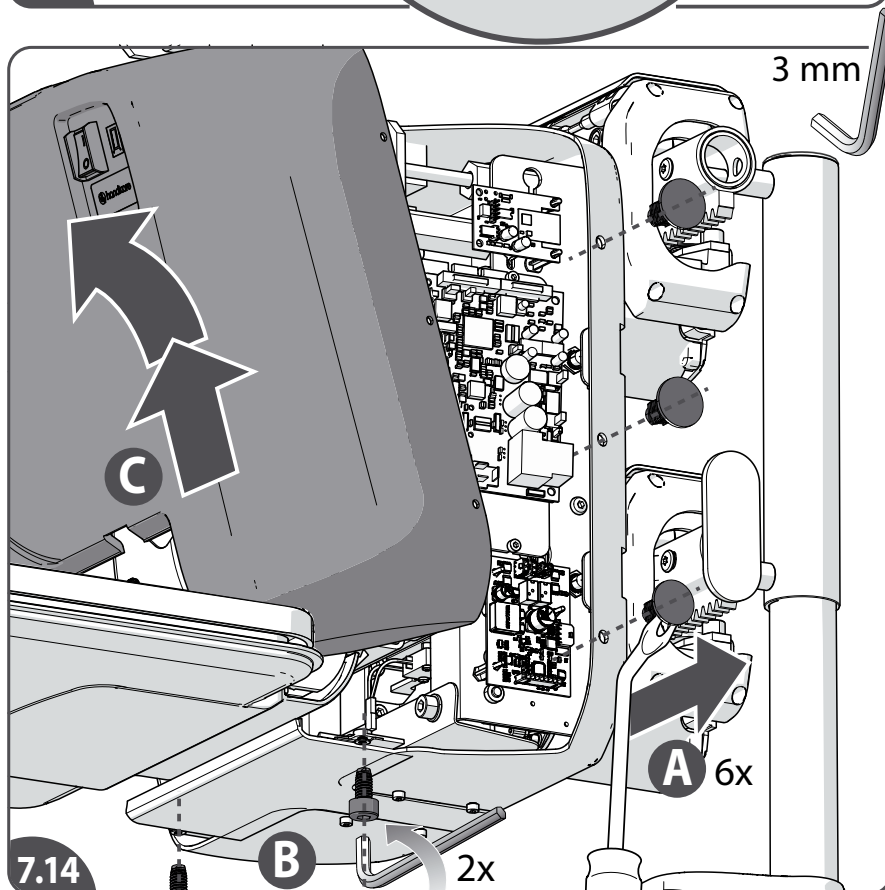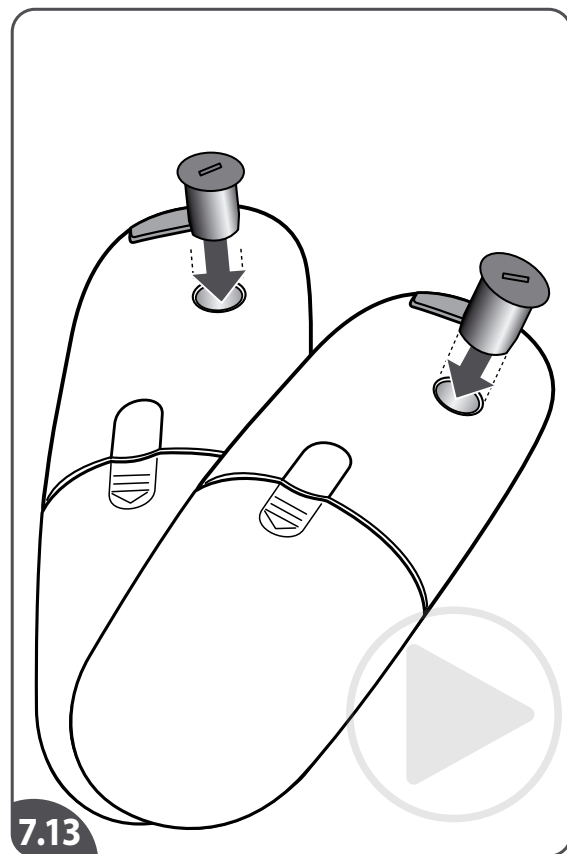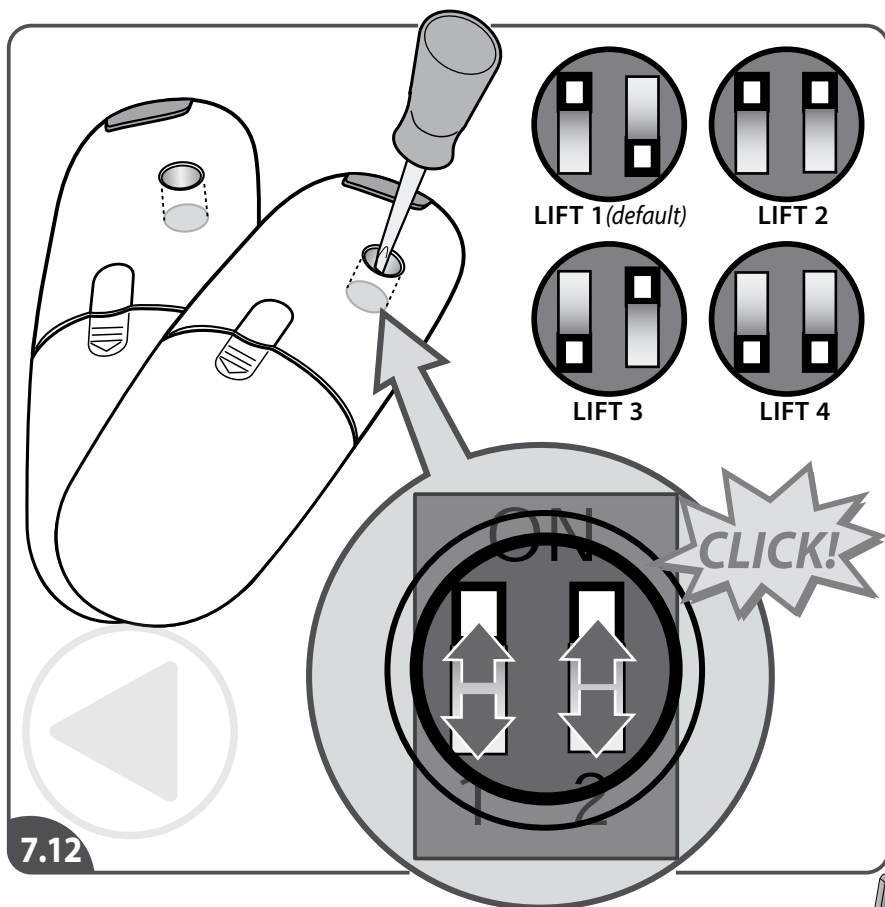- [7.12]Change the setting on the dipswitches on all remote controls.
- [7.13]Refit the caps on the backside of the remote controls.
- [7.14]Loosen the rivets(6x){A} and screws(2x){B} and remove the front cover{C}.
- [7.15]Change the dipswitches on the RF board to the same settings.

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


IT

7 REMOTE CONTROL CHECK

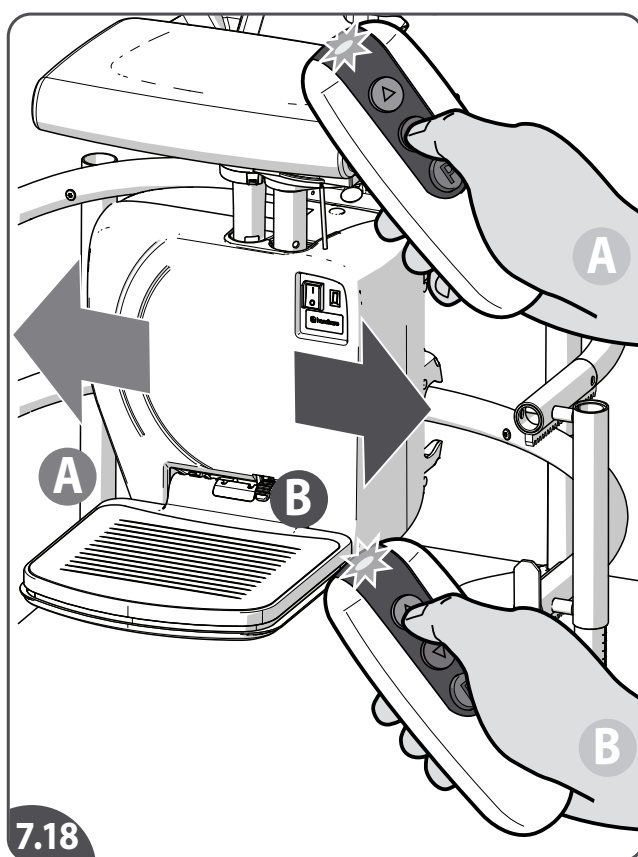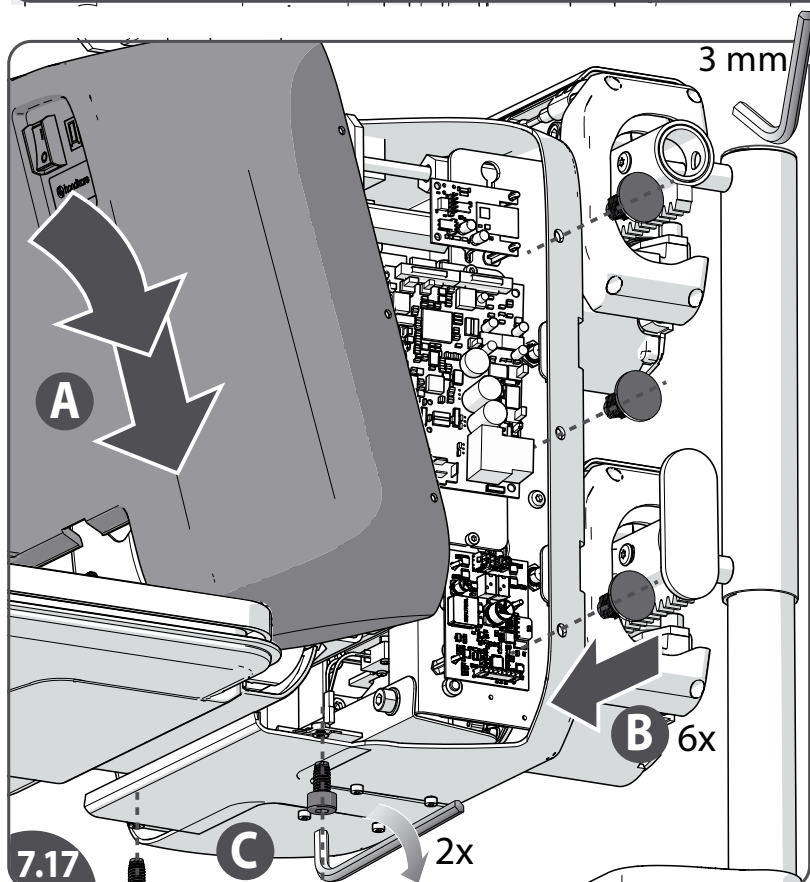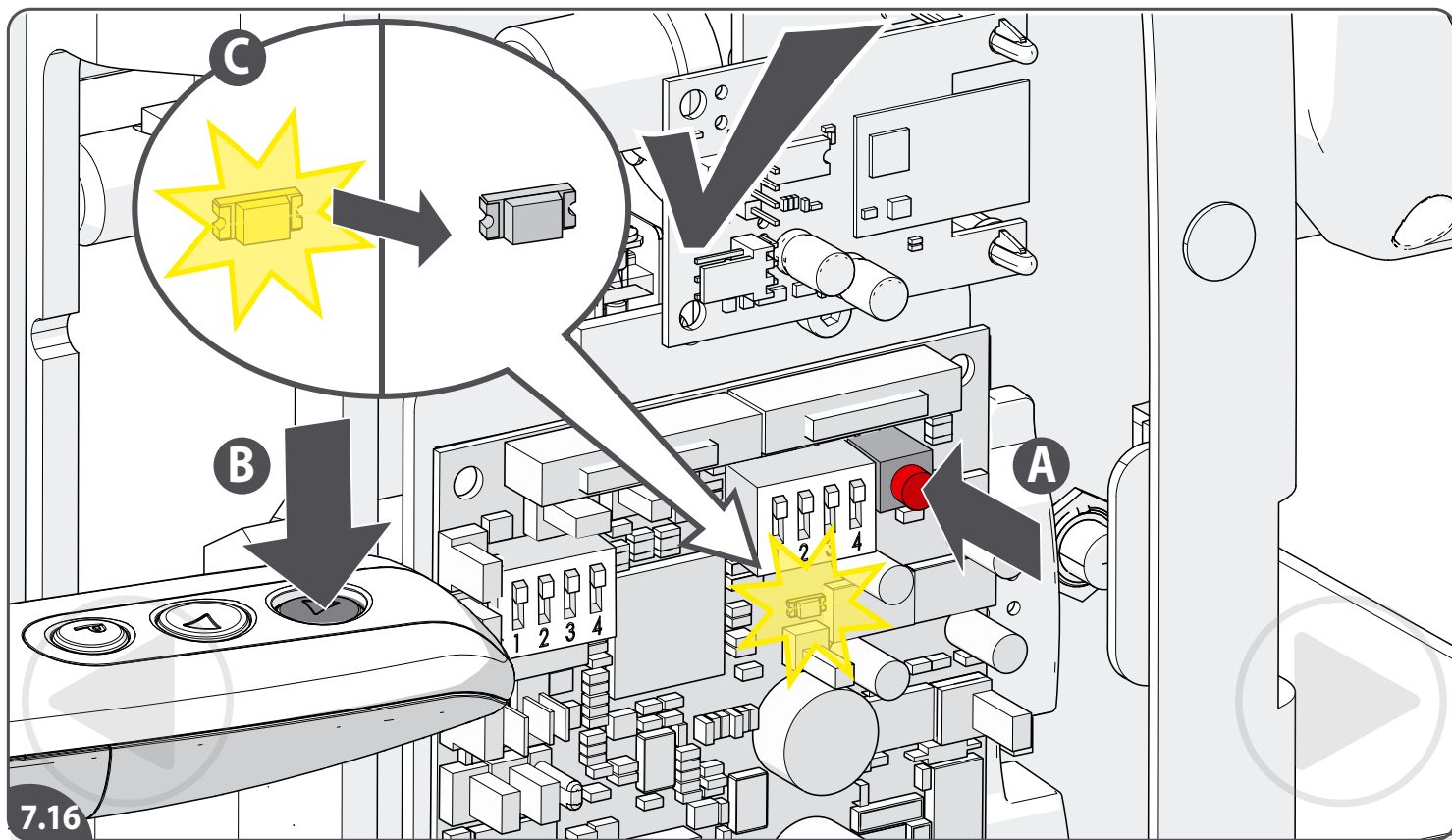[7.16] Press and hold the red button on the main PCB{A}. Press any button on the remote control{B}. When the light will go out, the programming is complete{C}.

Remark: If the red button is depressed for more than 10 seconds then the remote control will not be programmed. You must repeat the programming again.

Note: Always remember after any changes to main PCB to switch off for 15 seconds and switch back on again.

[7.17] Reinstall the front cover{A} and secure with the rivets(6x){B} and screws(2x){C}.

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




[7.18] Provare i telecomandi.

IT

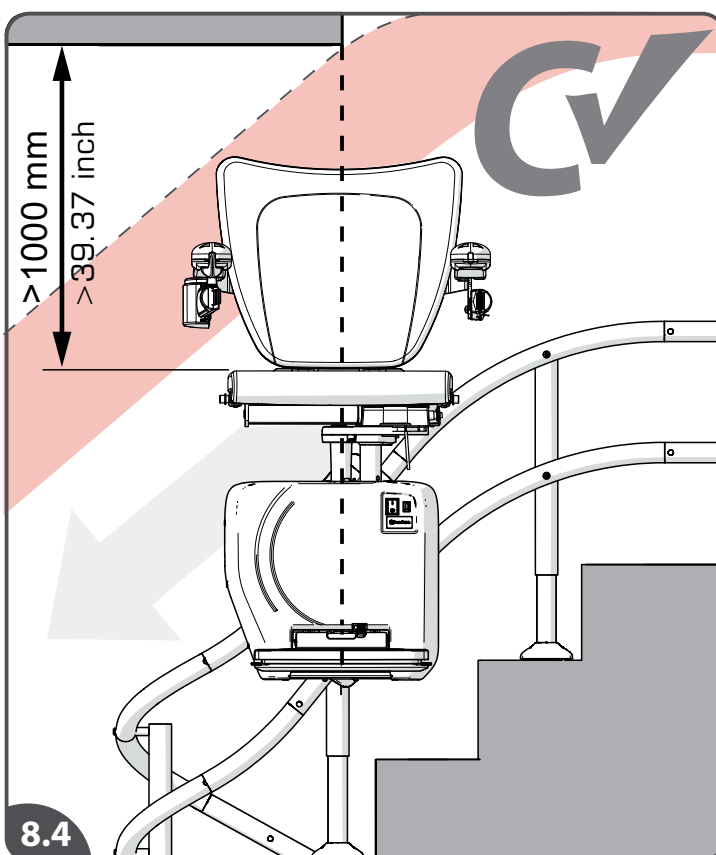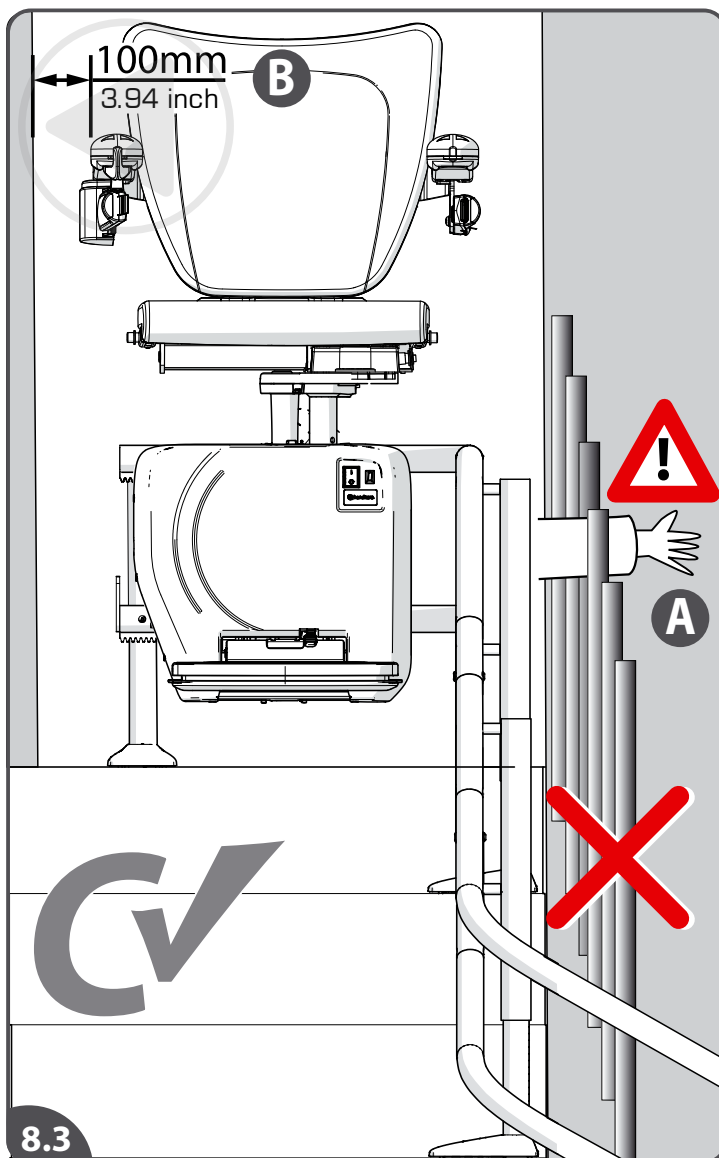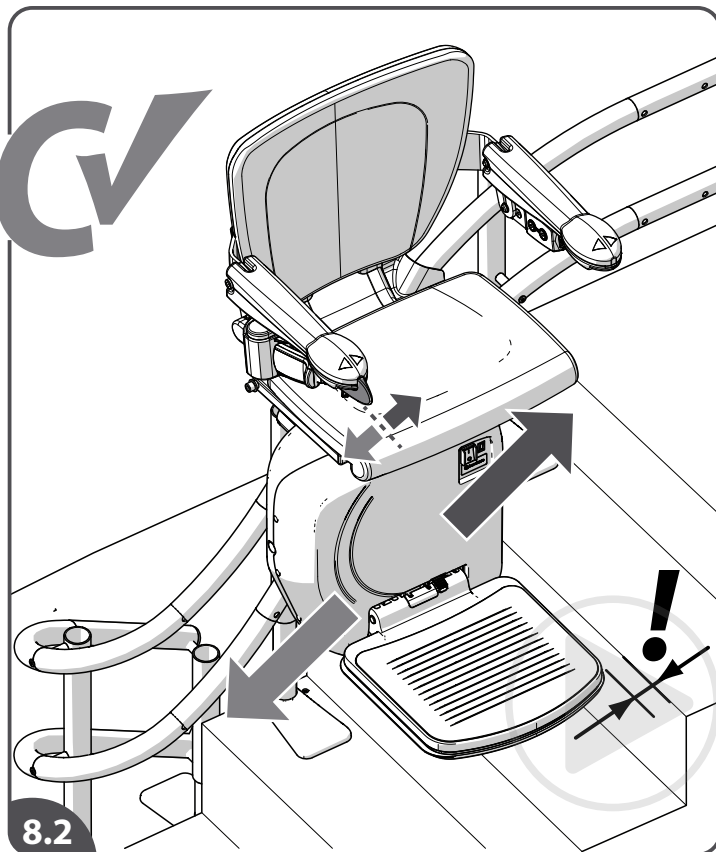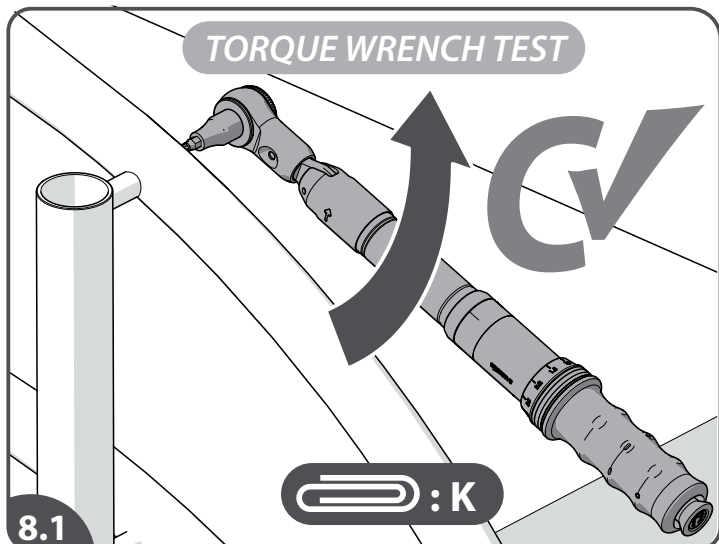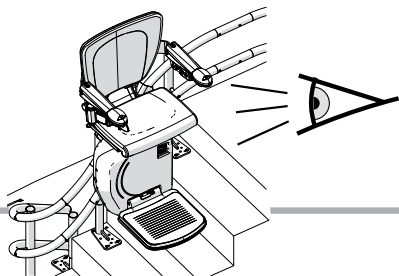- [8.1] Check the torque of the load-carrying fasteners. See attachment **K** how to do this.
- [8.2] Move the toggle/joystick up and down. Check whether the stairlift responds correctly.
- [8.3] Make sure that there is nothing behind the stairlift that body parts could get trapped in or on.
- [8.4] Check the position of the bulkhead.

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




IT

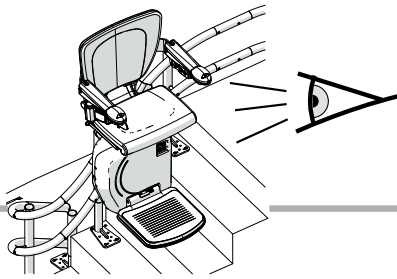CHECK SWIVEL-SAFETY BELT

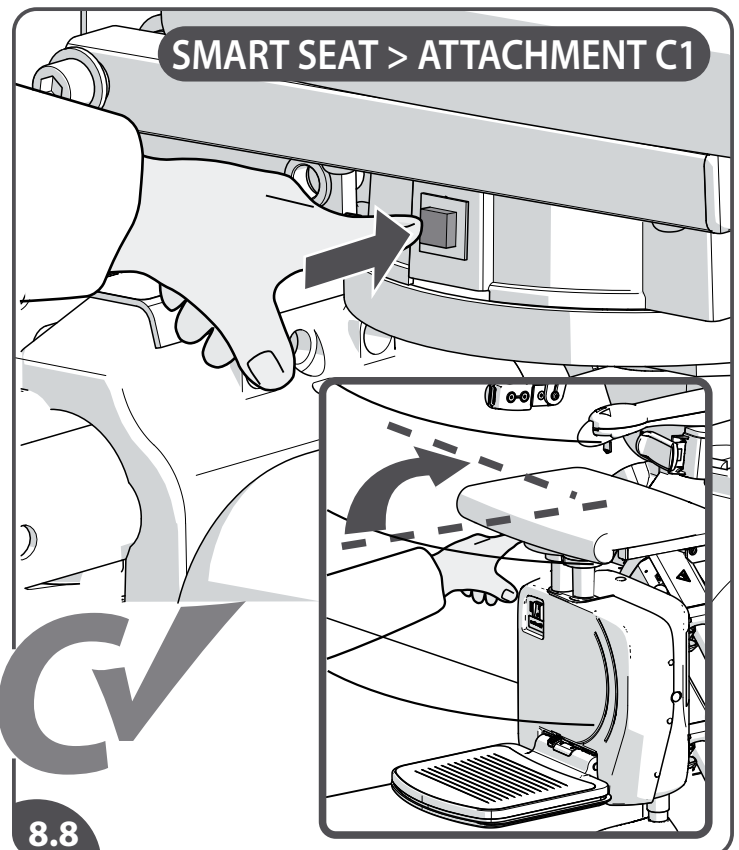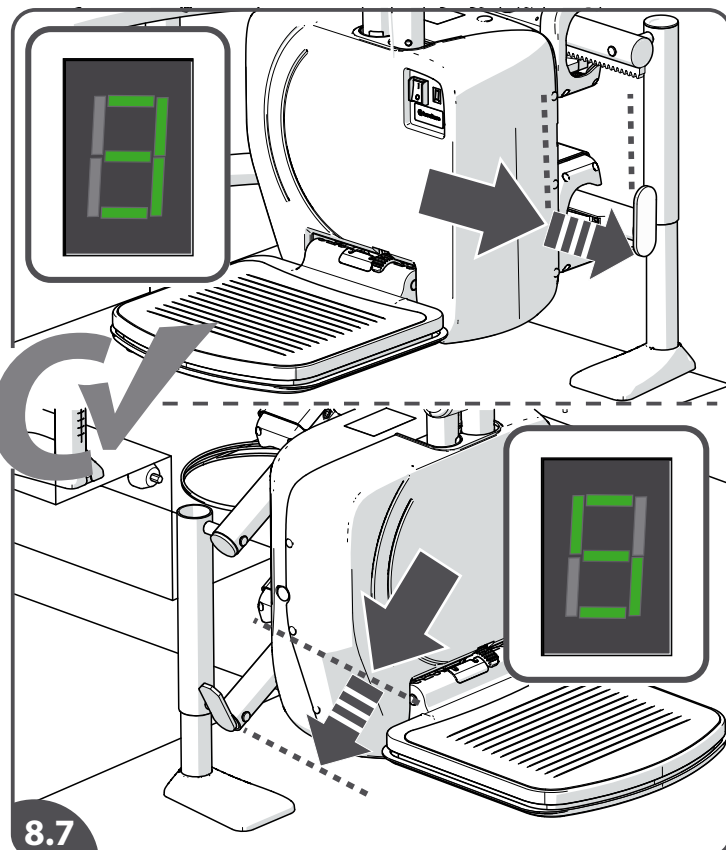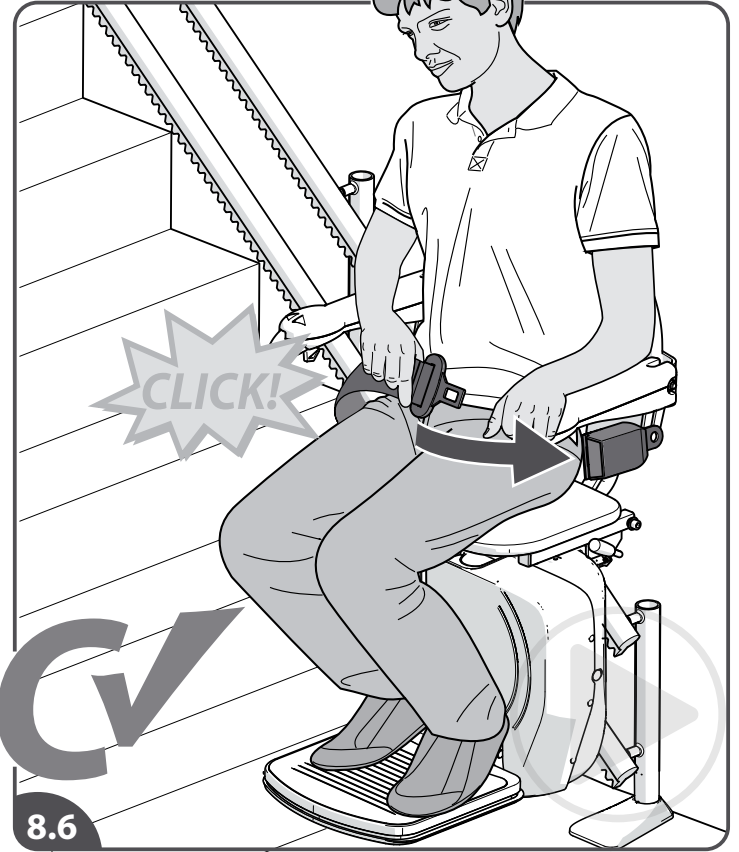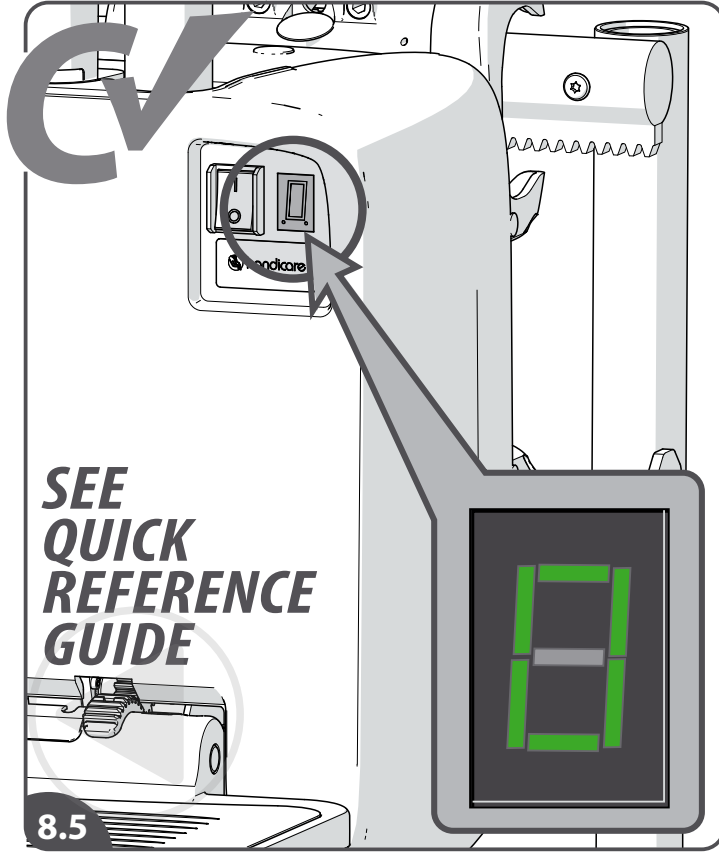- [8.5] Check the diagnostic display; make sure there are no fault messages. See the user manual for the error codes.
- [8.6] Check the seat belt.
- [8.7] Check that the stairlift stops correctly at the top and bottom of the rail.
- [8.8] Powered swivel only > check the release button: press the release button and swivel the seat manually.

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




IT

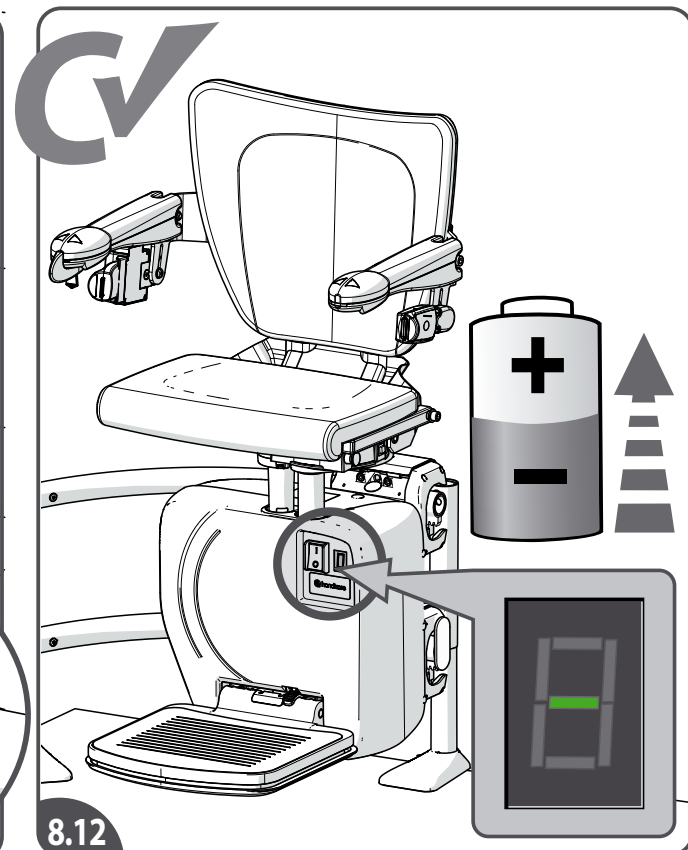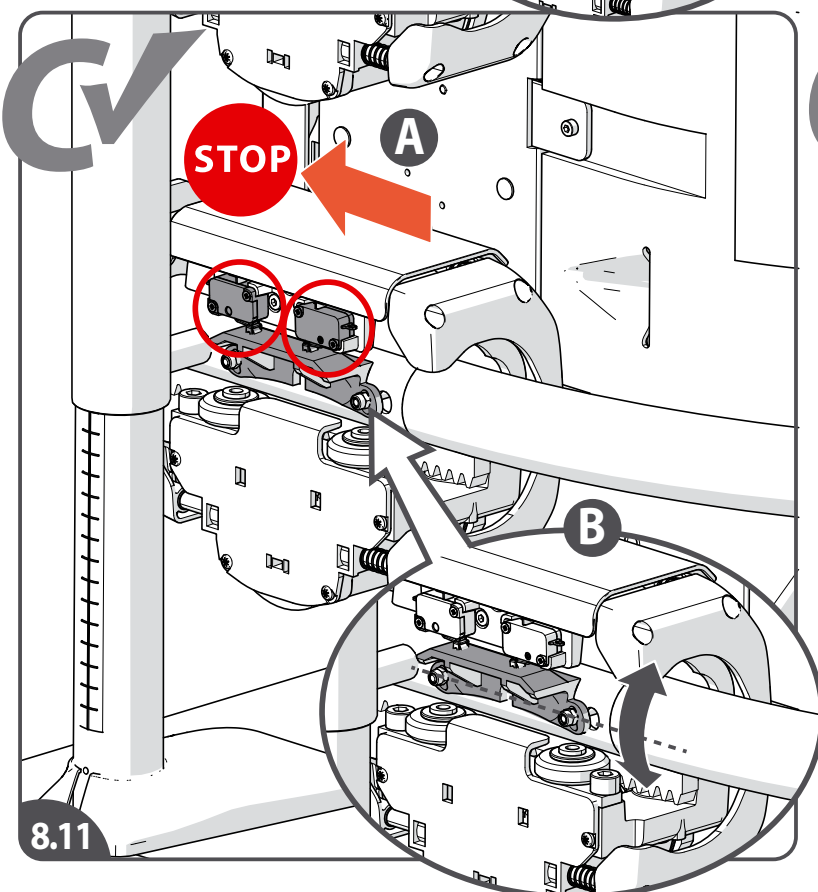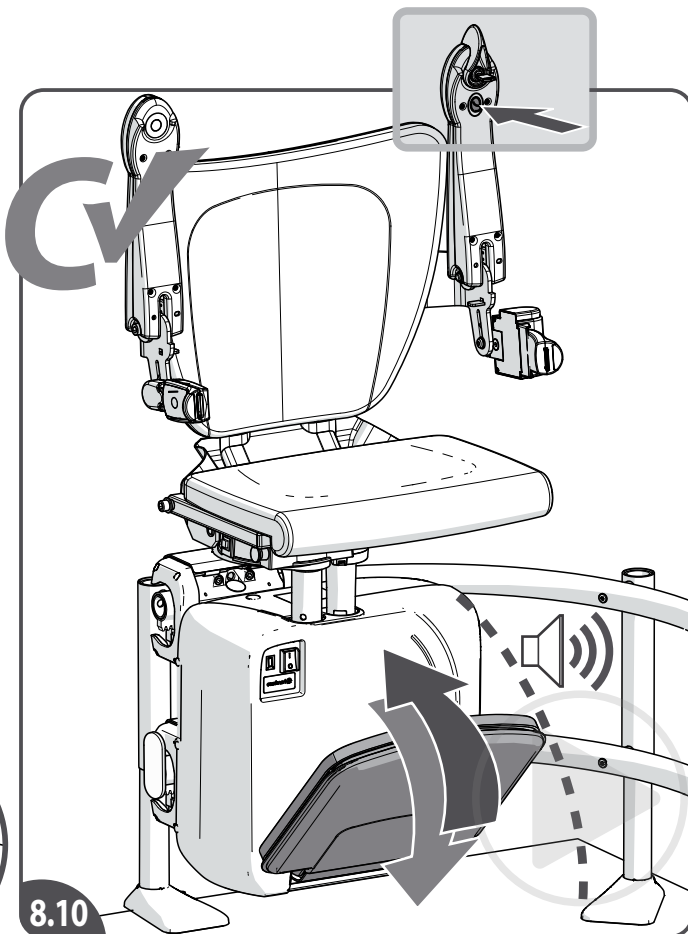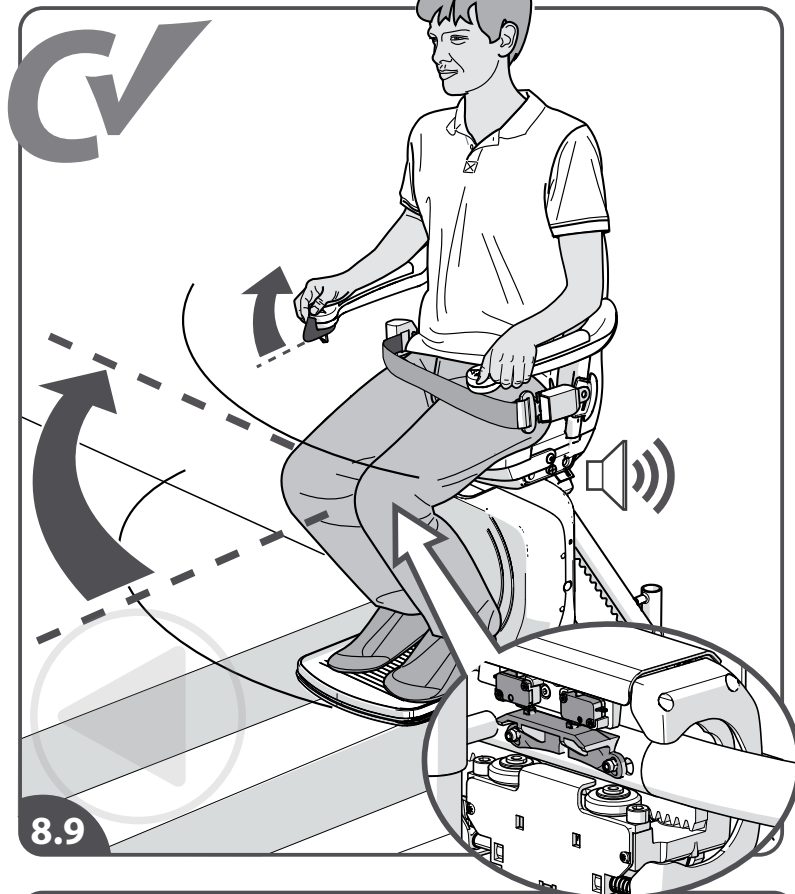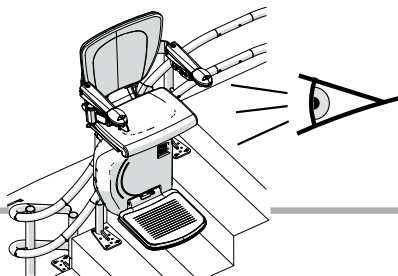- [8.9] Check the powered swivel function. A lift (with an automatic swivel configuration) is able to swivel when:
1. The stop switch hits the limit actuator at the top of the staircase.
 2. There is 33 volt on the charge point.
- [8.10] Check the powered footrest function.
- [8.11] Check the limit actuator (both switches activated){A}. If needed, adjust{B}.
- [8.12] If the lift has an intermediate charge point, test the park function. Ensure that the diagnostic display shows a “-” symbol at the charging point.

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

IT

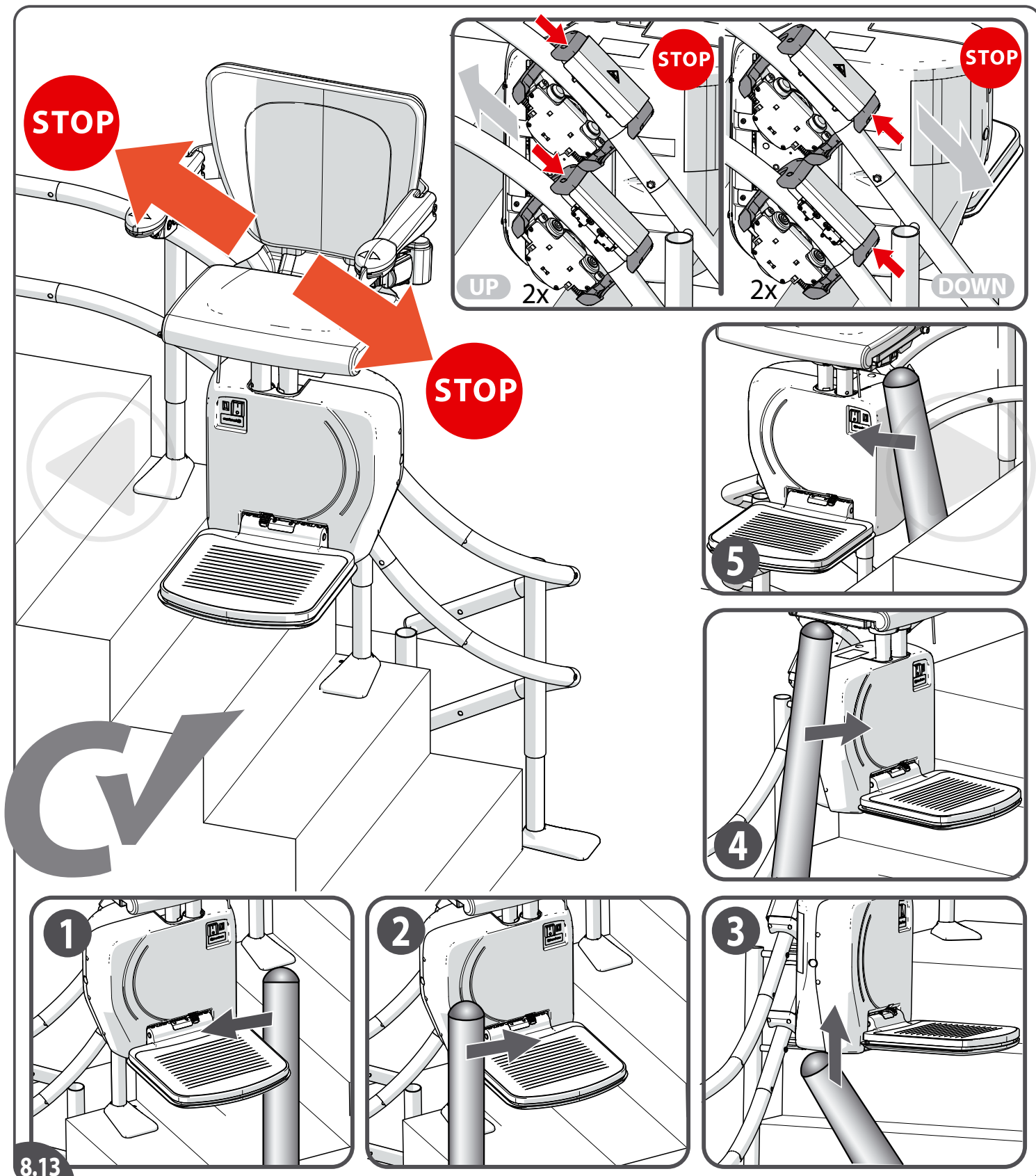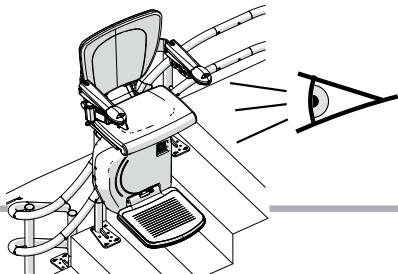[8.13] Check the safety stops.

Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to the requirements of the latest adopted ASME A18.1 code.

Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020

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



IT

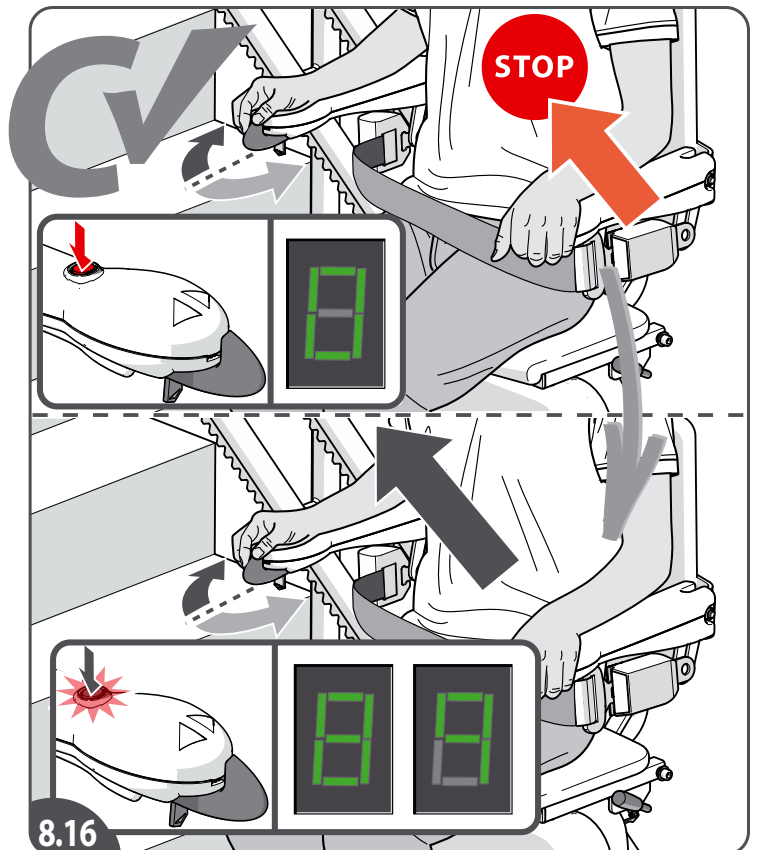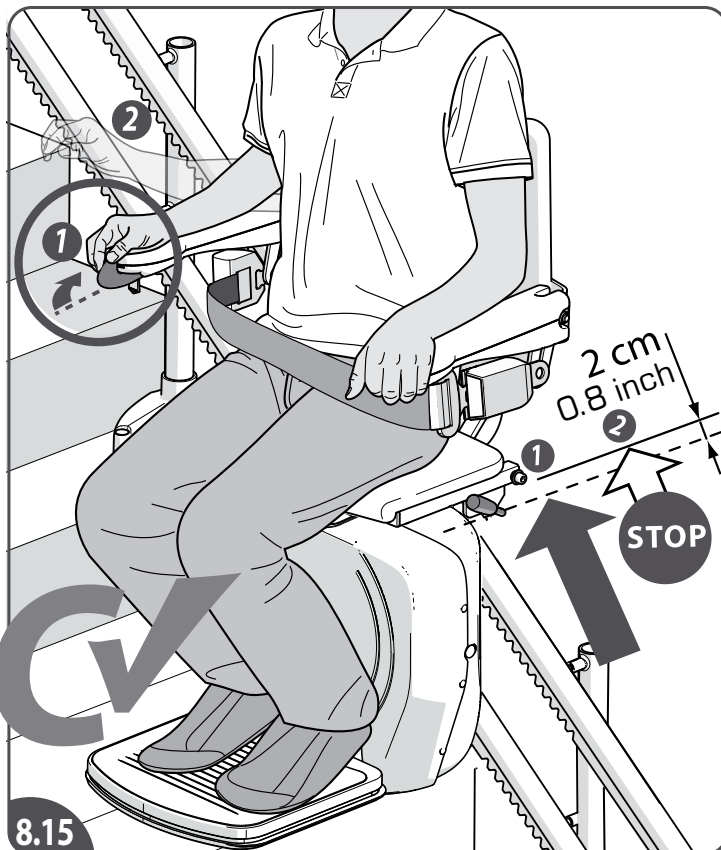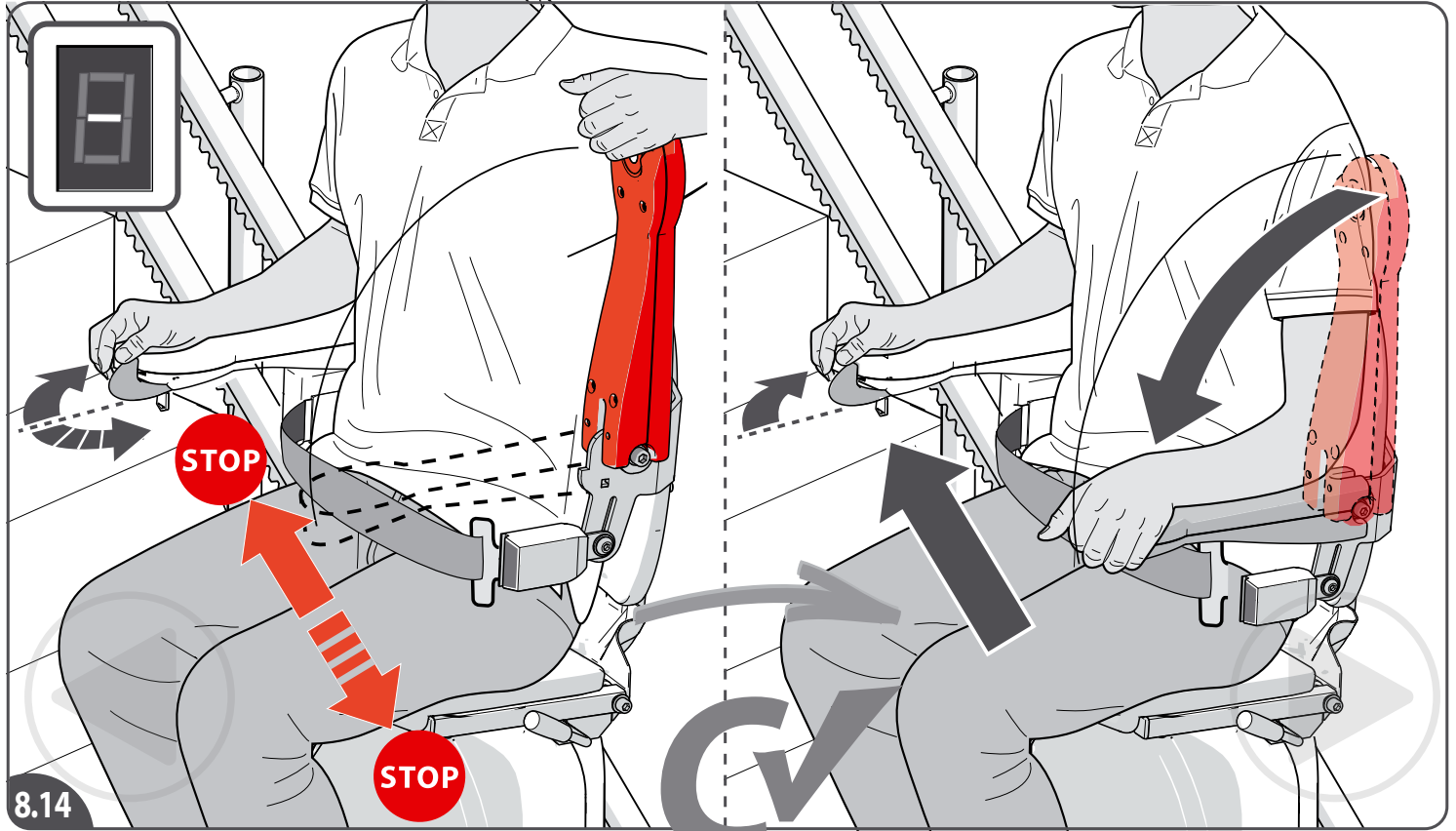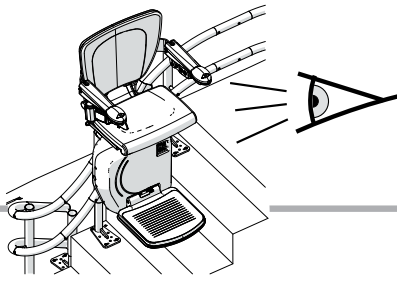- [8.14] Fold the downside armrest up. The stairlift must not operate (lift will still operate from remote controles). Put the armrest back down and check that the stairlift operates normally again.
- [8.15] Release the joystick toggle switch. The stairlift should stop within 2 cm (0.8 inch).
- [8.16] Press the stop button. The lift must stop immediately. Press the stop button again; check that the stairlift operates normally.

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




IT

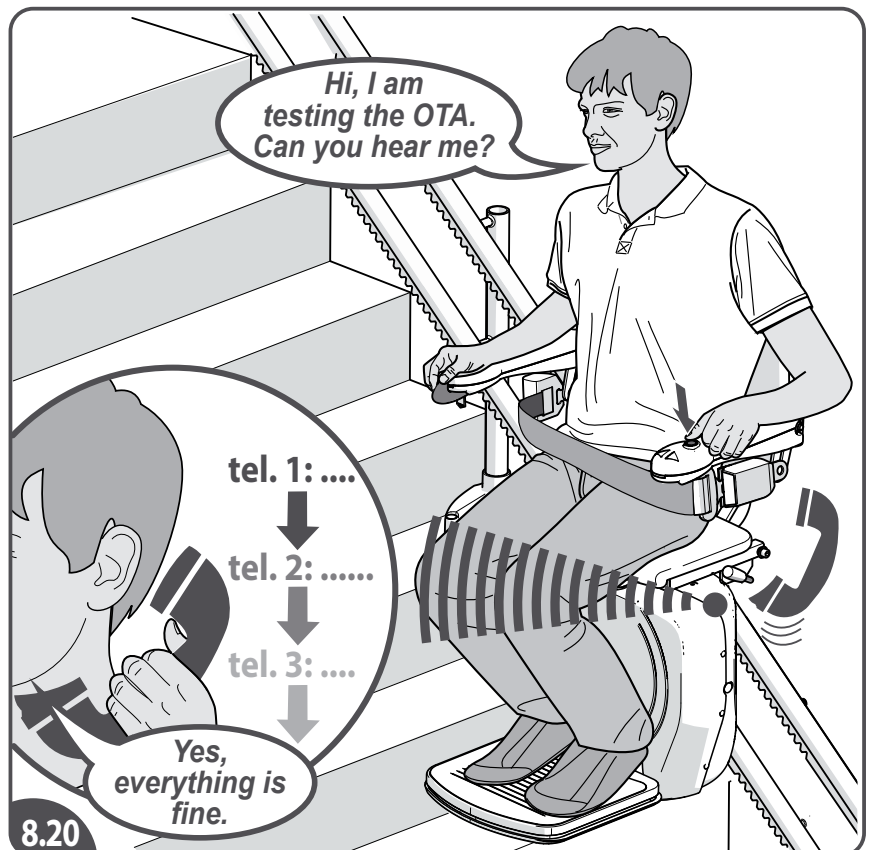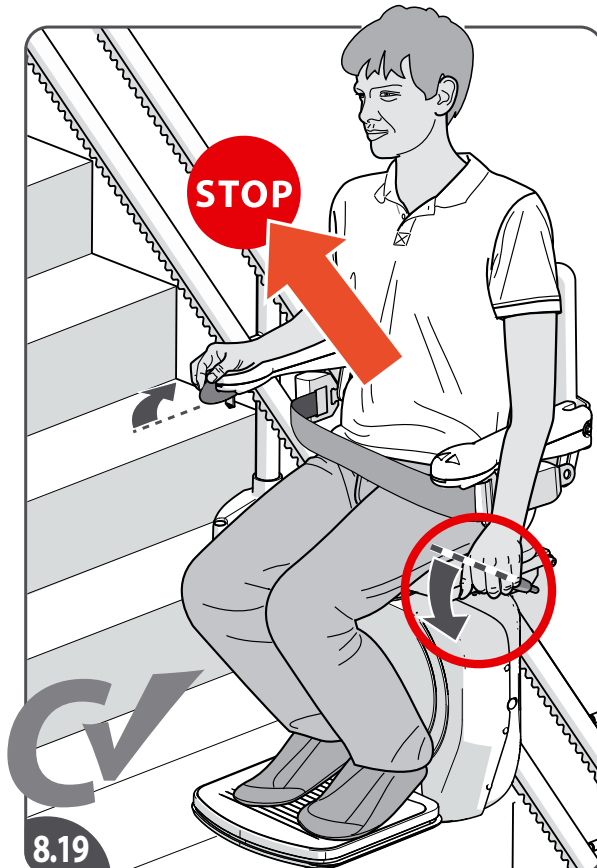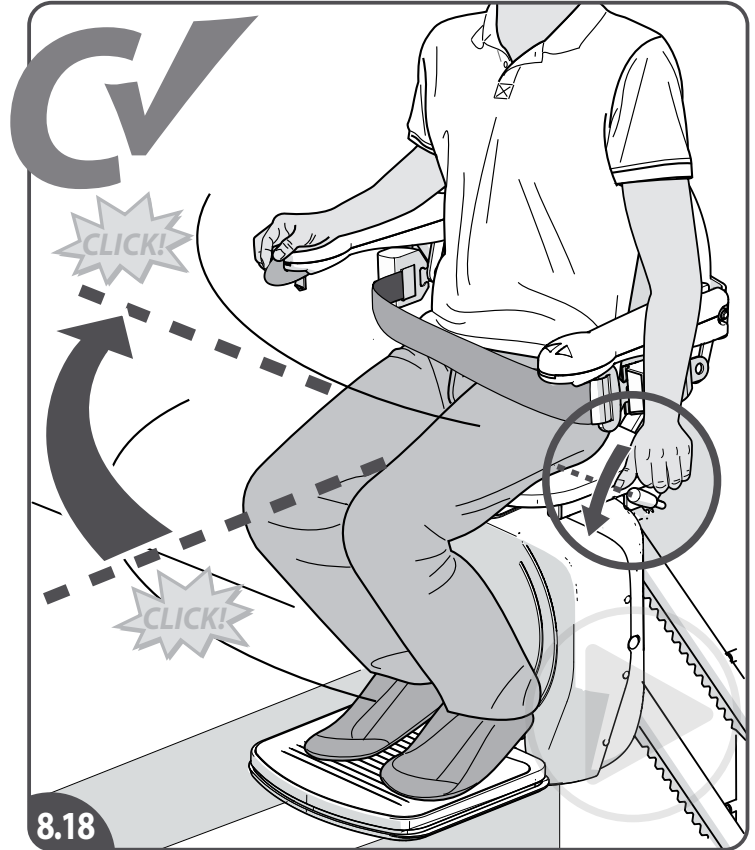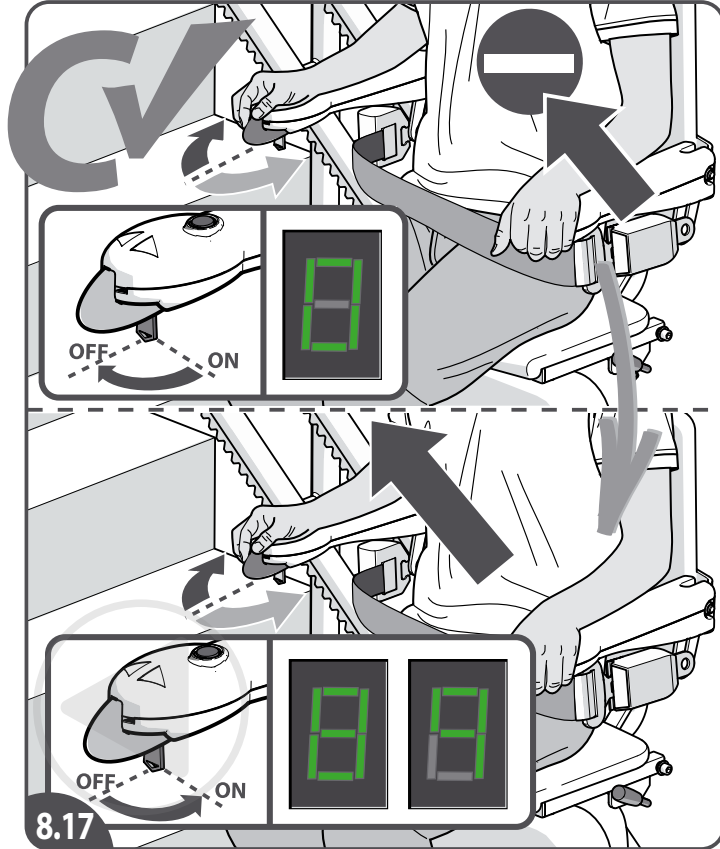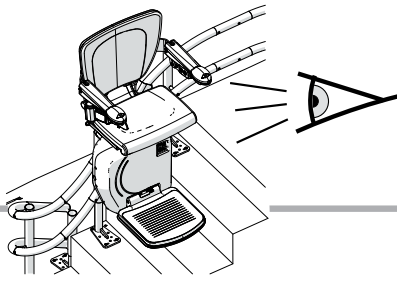- [8.17] Turn the key to the OFF position and ensure that the stairlift will not move. The diagnostic display should show a "O". Turn the key to the ON position. The stairlift should be able to move again.
- [8.18] Check the manual swivel function and that the end stops work satisfactorily (it clicks at the end).
- [8.19] While moving the stairlift, swivel the seat; the stairlift must stop immediately.
- [8.20] It is the responsibility of the Dealer to ensure that the Integrated TWVC System is checked and test calls are performed at installation and normal stairlift service intervals. As it is the responsibility of the customer (end-user) to provide and manage the Sim Card used within the Integrated TWVC System.

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



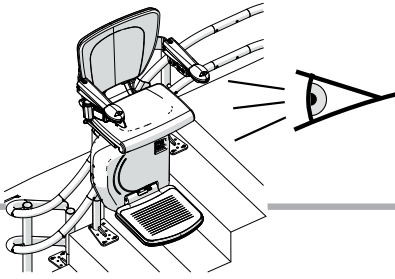
CHECK

BELT SWITCH - HINGE

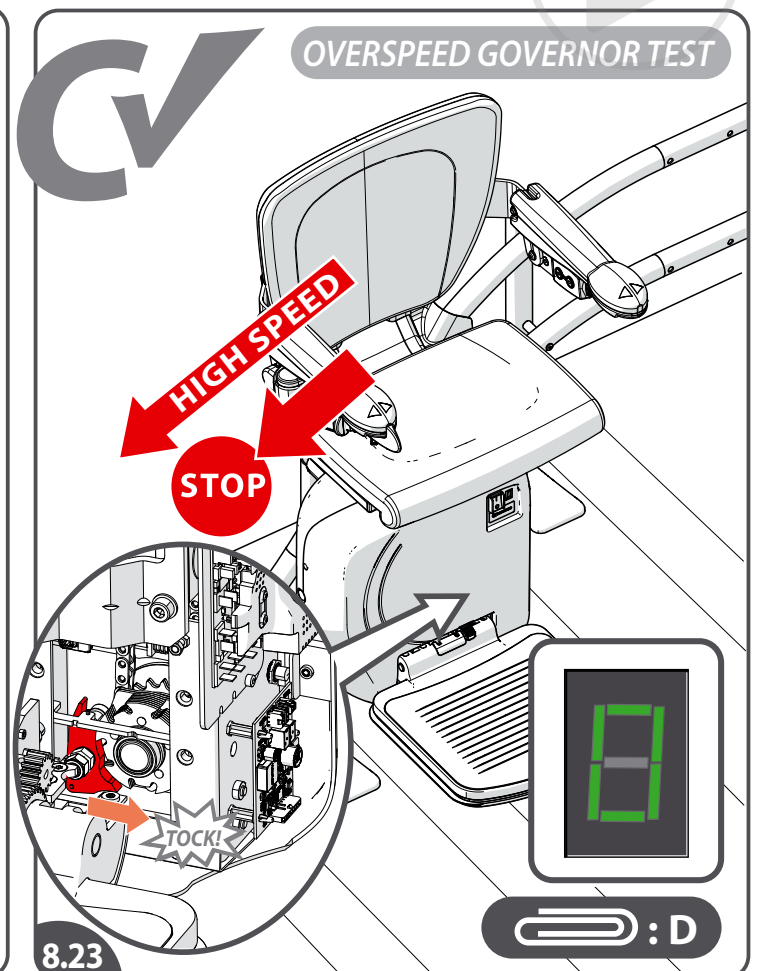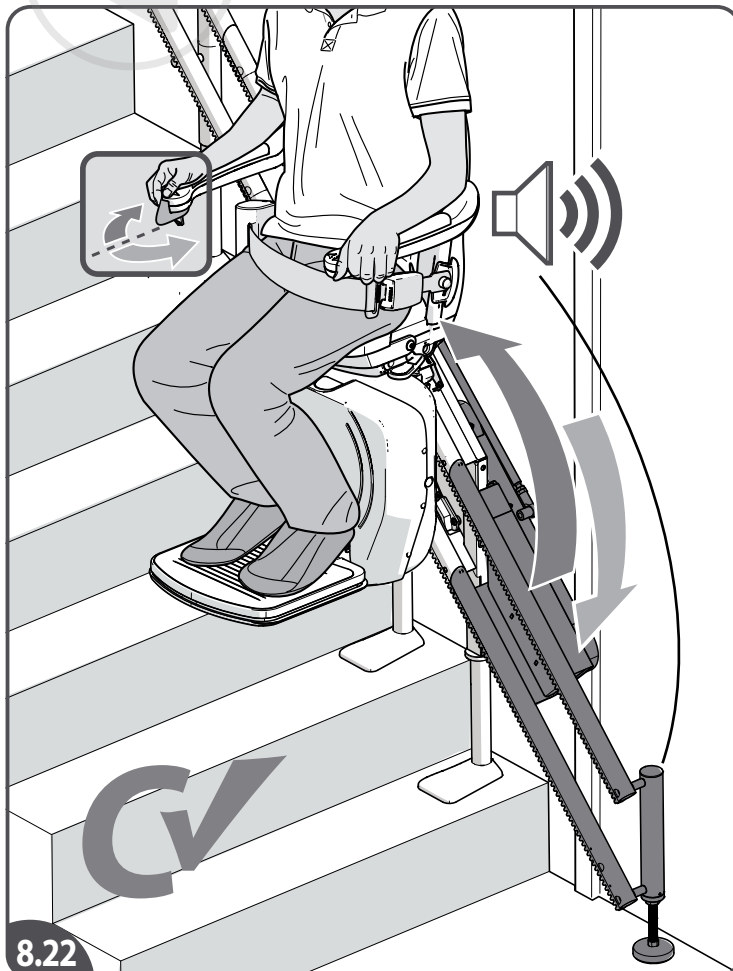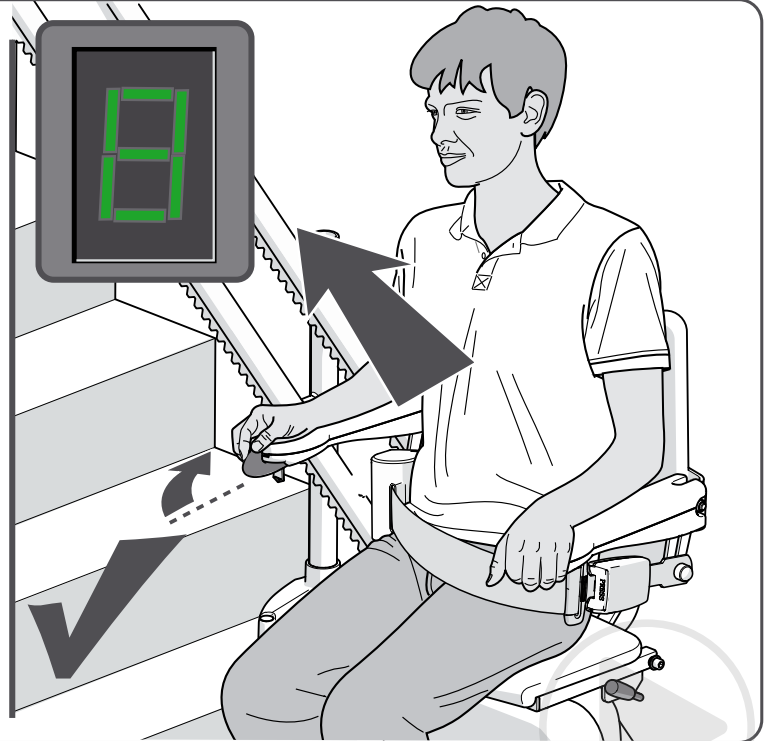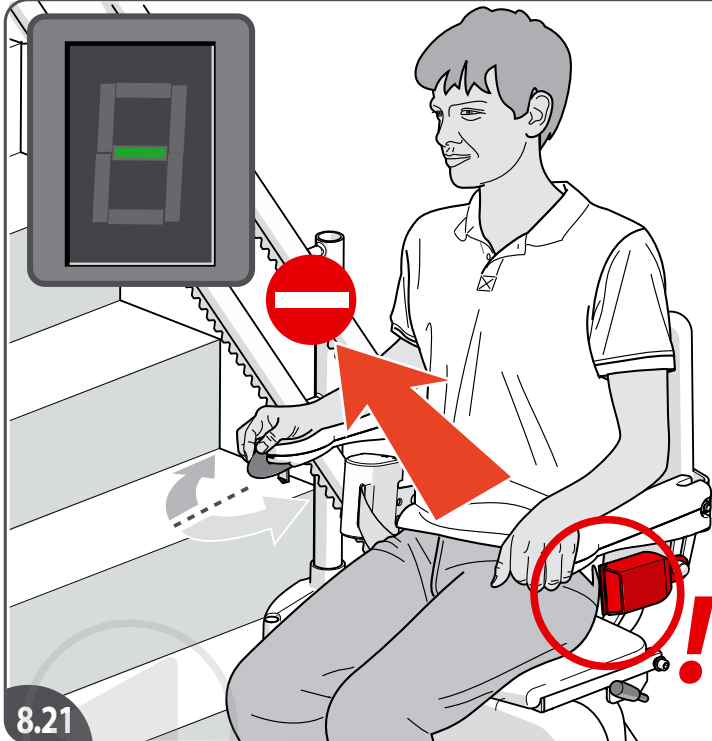- [8.21] Check the seat belt safety.
 [8.22] Check the powered hinge track function.
 [8.23] Check the OverSpeed Governor(OSG) safety (see attachment D).

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

IT

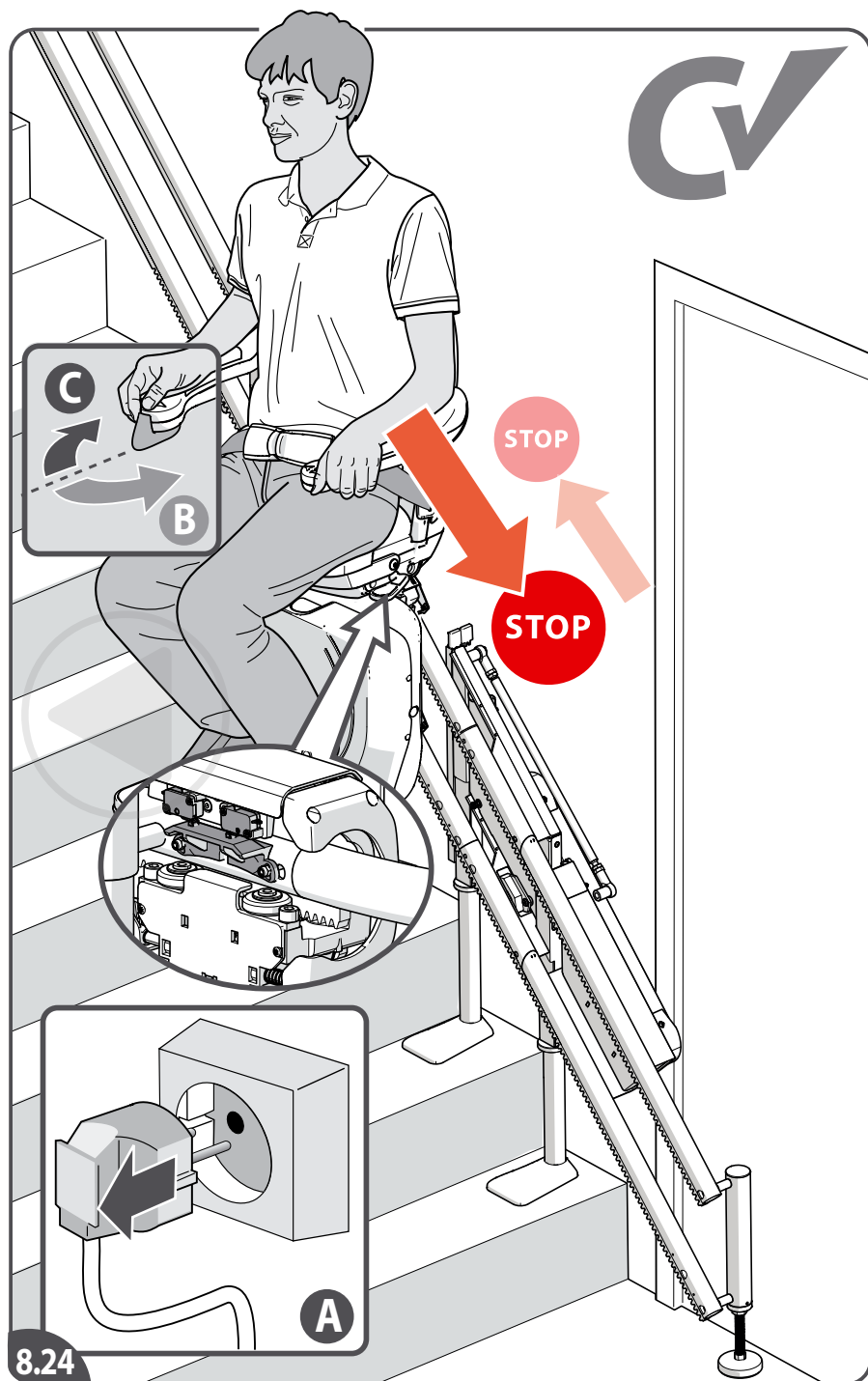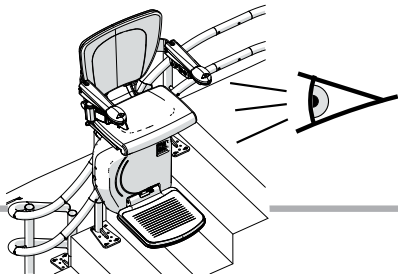8.24

[8.24]The following test must be carried out to check for the secondary interlock operation:

- 1-Ensure lift is above intermediate charge contacts.
- 2-Turn the 33VDC power supply OFF.
- 3-Drive the lift in the DOWN direction.
- 4-The safety interlock switch (on the top trunnion) should strike the bracket on the top rail and stop. The hinge will lower, and then the lift can drive down.
- 5-Drive the lift back UP to raise the hinge.
- 6-The switch should also stop the lift in the UP direction to raise the hinge.

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IT

[8.26] A full ride is made along the complete travel, include hinged rail guide sections and swivelling of the seat where applicable; fill in the installation checklist. Corrective actions required for boxes that are ticked as "No".
 [8.27/28/29] The declaration of conformity is legal after the declaration engineer and the declaration instructed person have been signed off by the engineer and the user. The declaration engineer and the declaration instructed person shall be saved by the dealer and made available to Handicare on request.

EN *Routine, Periodic, and Acceptance inspections and test in the USA shall be done according to the requirements of the latest adopted ASME A18.1 code. Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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NL *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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ES *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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
DE *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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FR *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

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
IT *Routine, Periodic, and Acceptance inspections and test in Europe shall be done according to EN 81-40:2020.*

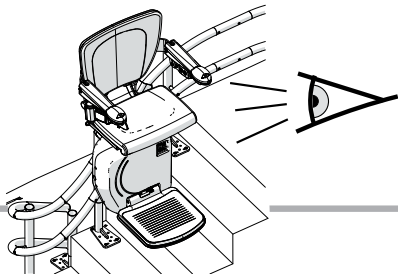CHECK

PLACING LABEL

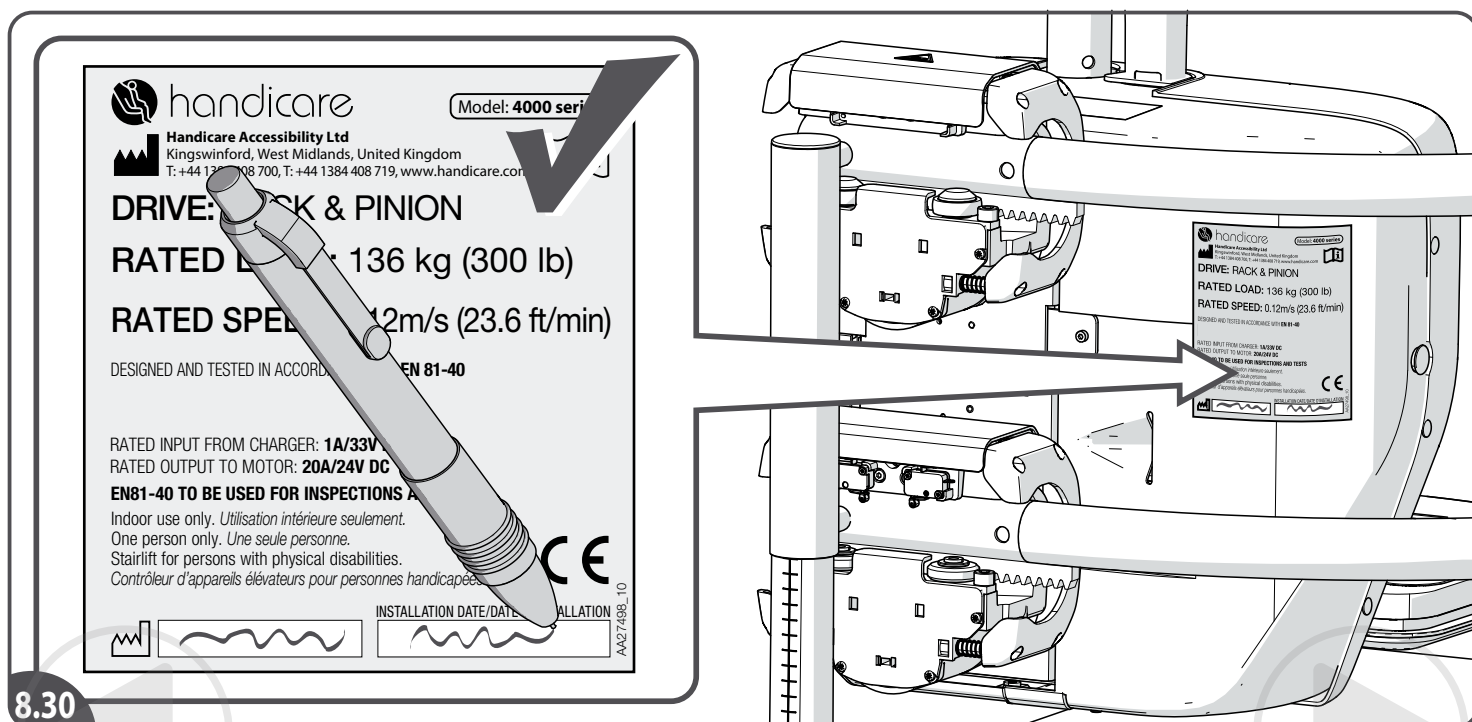8.30

[8.30] Fill in the installation dates on the label. The label is located on the rear of the power pack cover.

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




A

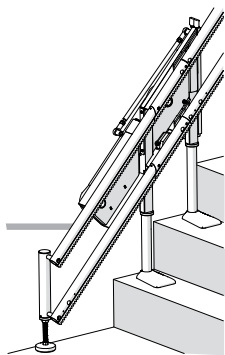HINGED RAIL POSITION

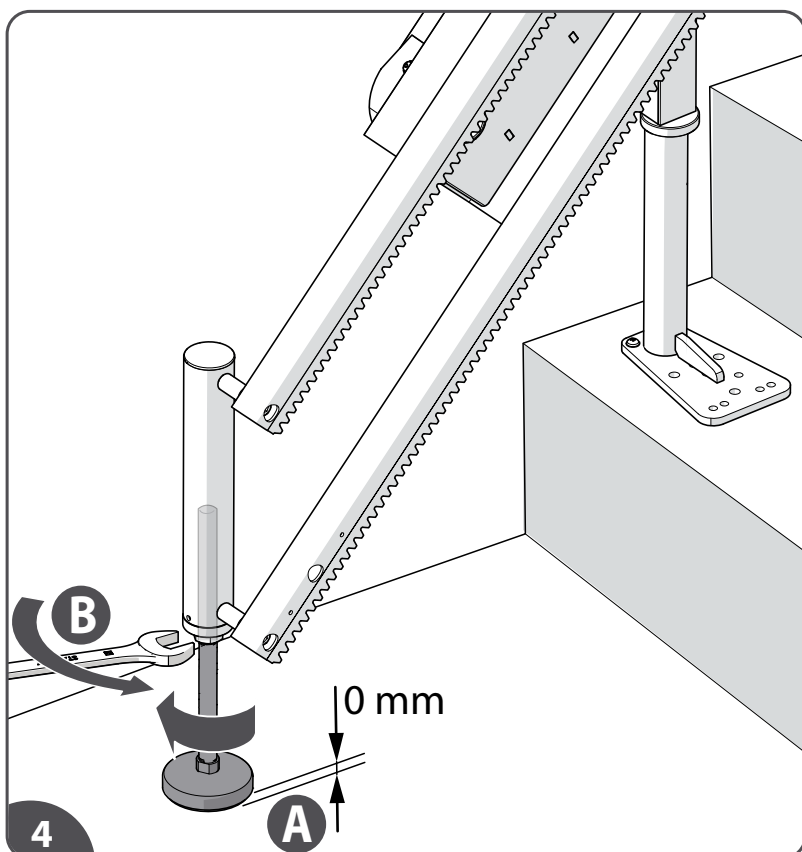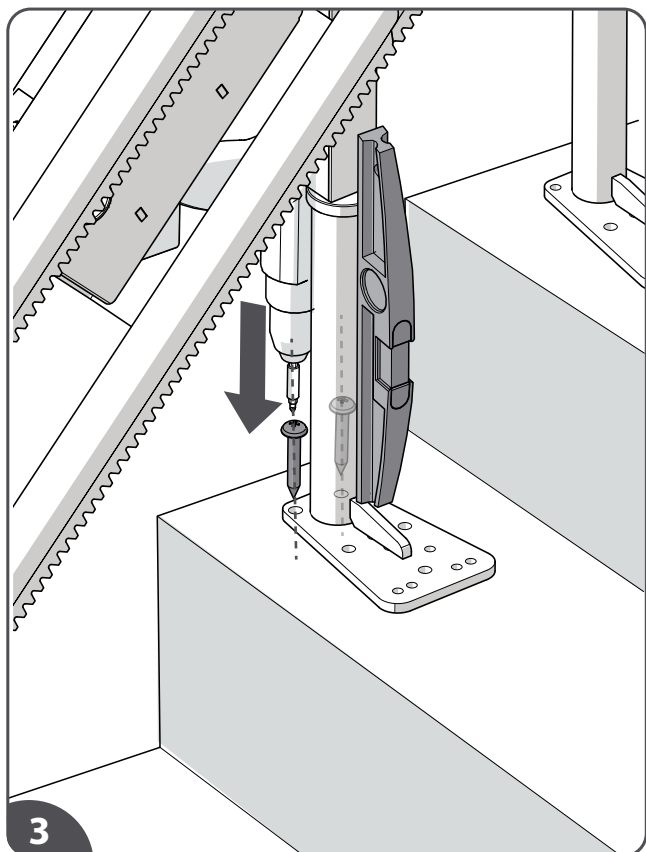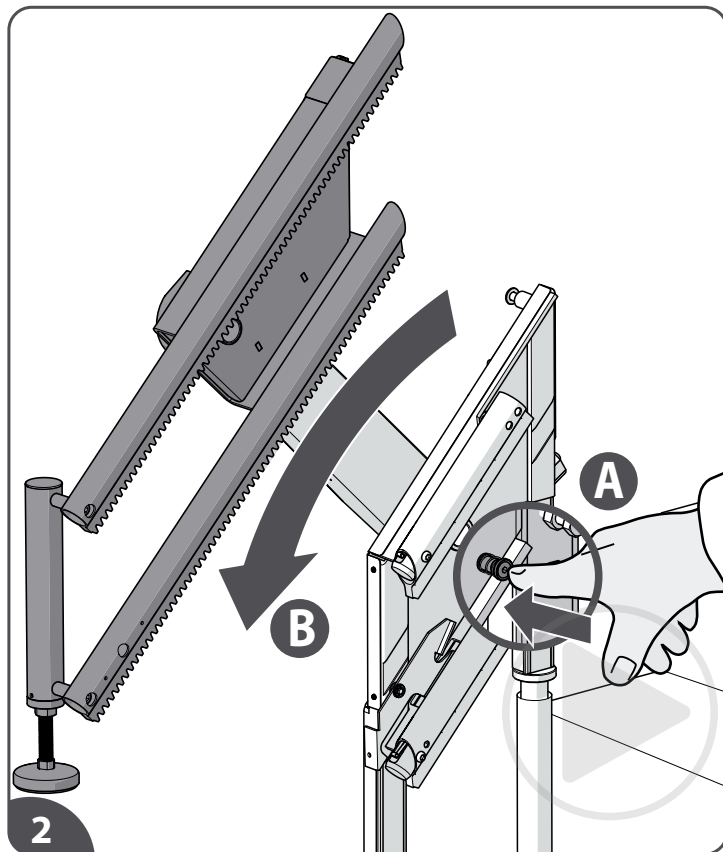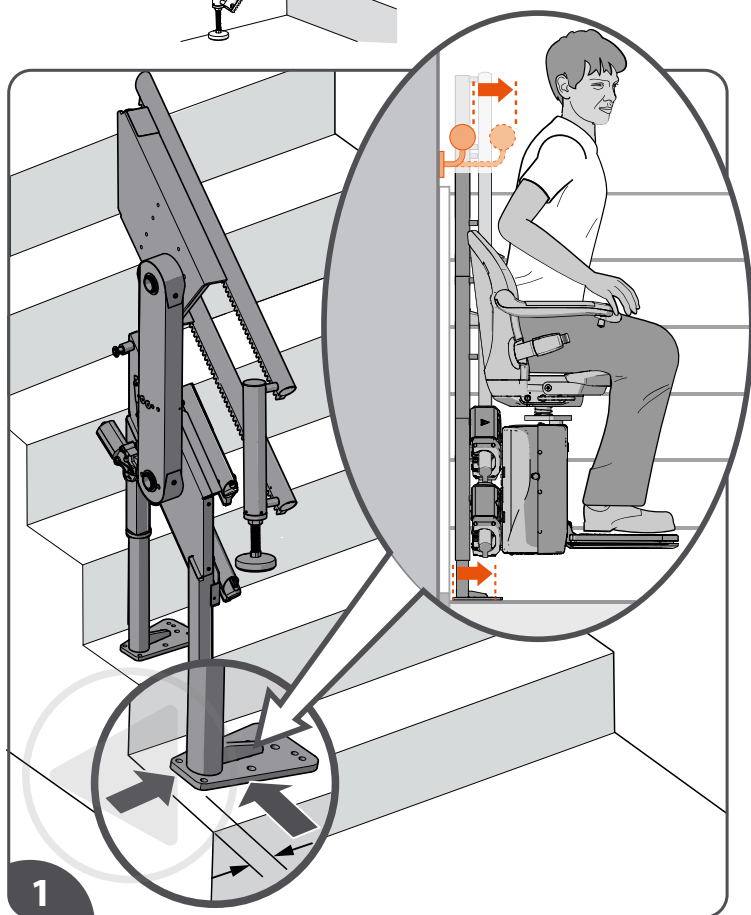- [1]Position the hinge track section in accordance with the track offsets stated on the drawing.
- [2]Release the motor by pressing on the motor release to allow the arm of the hinge to be moved.
- [3]Align the stanchions vertical and fasten the feet with two screws (both sides).
- [4]Lower the arm and position the hinge on the treads correctly so that the track reaches the floor{A}. Fasten the support bolt{B}.

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



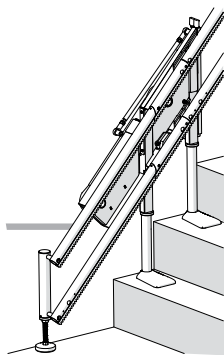
HINGED RAIL

INTERMEDIATE AND END STOPS

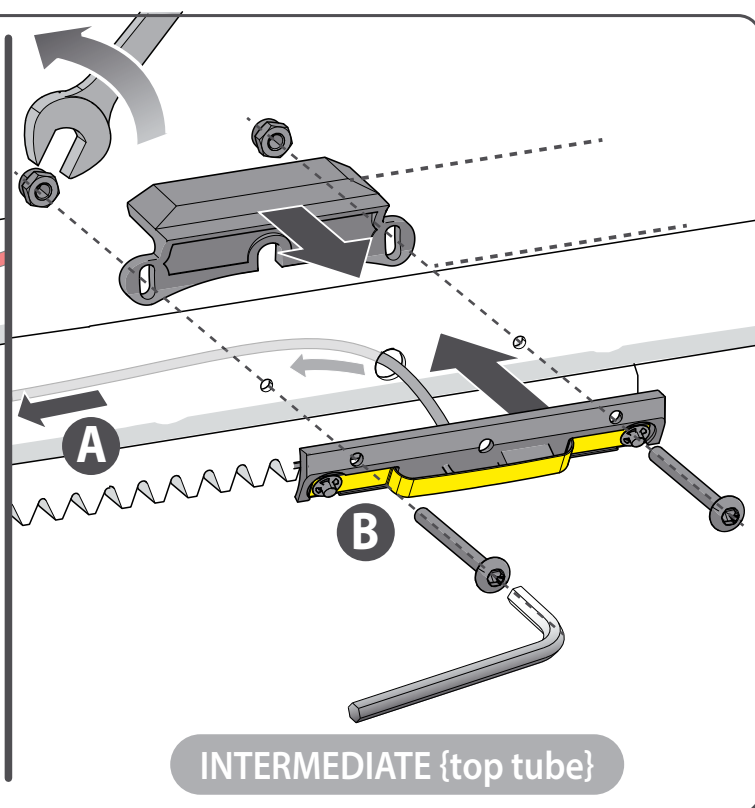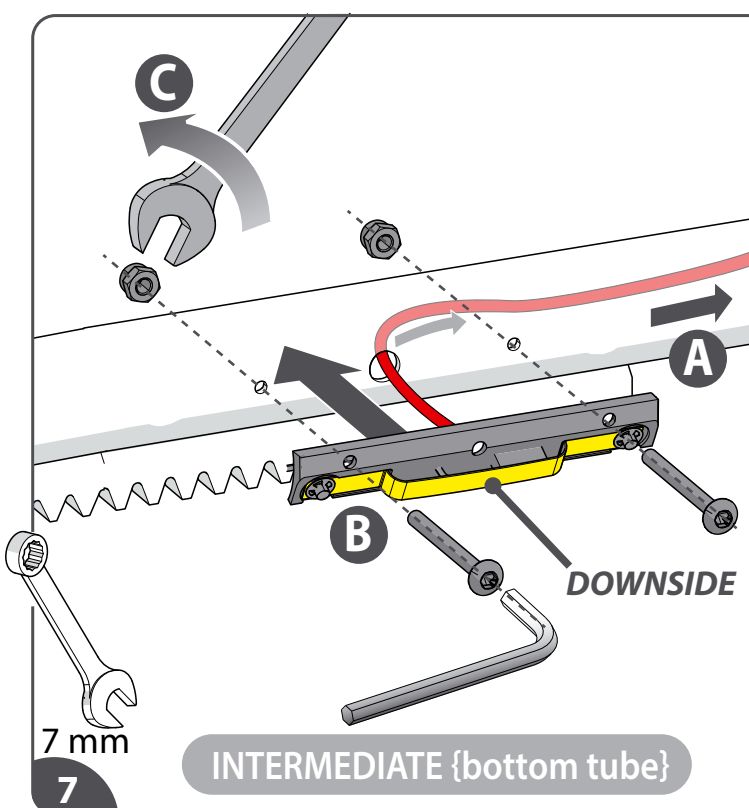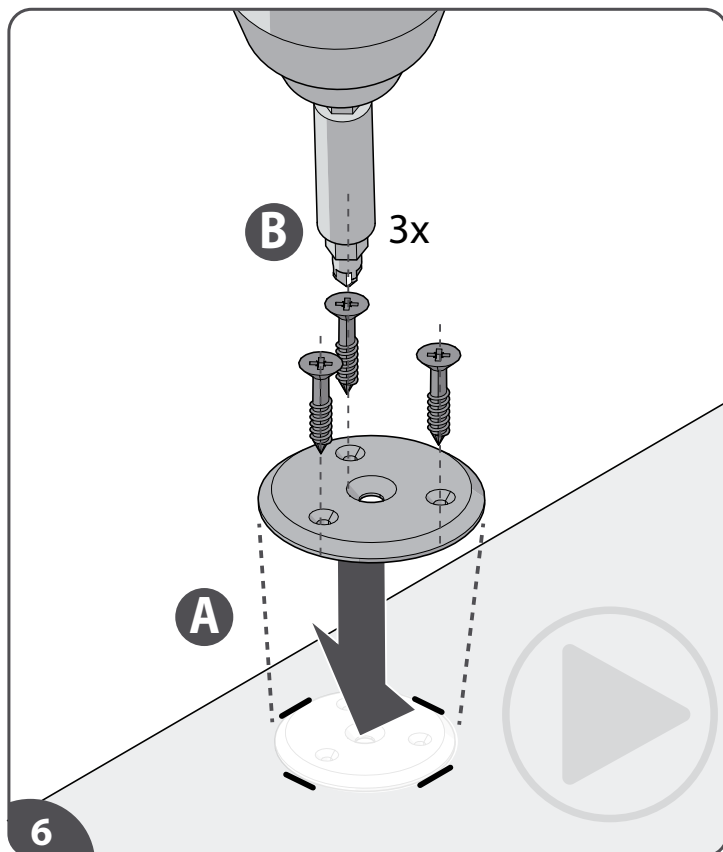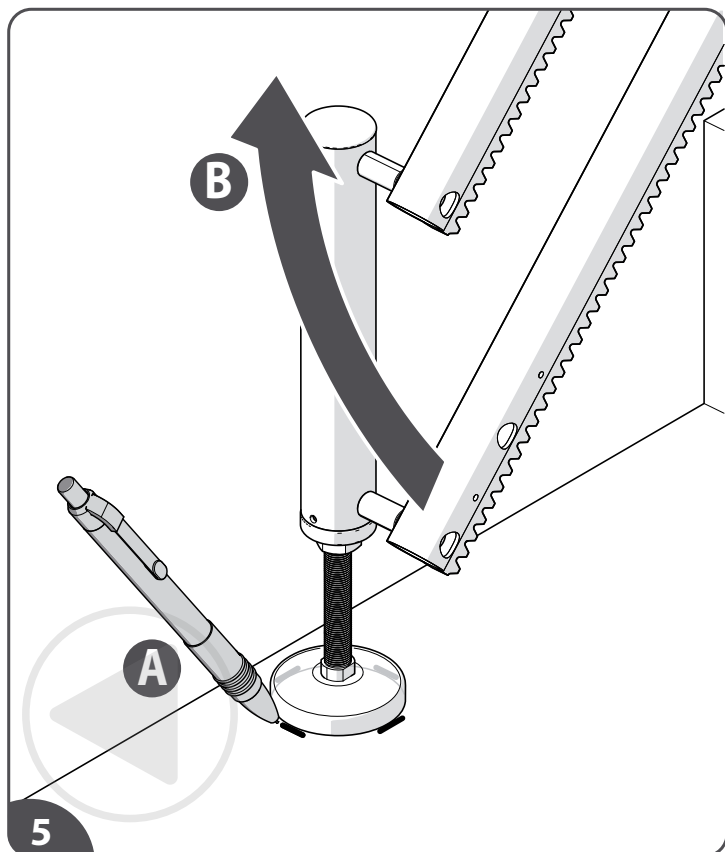- [5] Mark the position of the adjustable foot of the hinged rail {A}. Fold up the hinge {B}.
- [6] Position the mating plate {A}. Fasten the plate {B}.
- [7] Install the intermediate stop/chargepoints on the bottom and top tube on the indicated position.

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


CONNECT INTERMEDIATE STOP/CHARGE

IF NOT INSTALLED YET

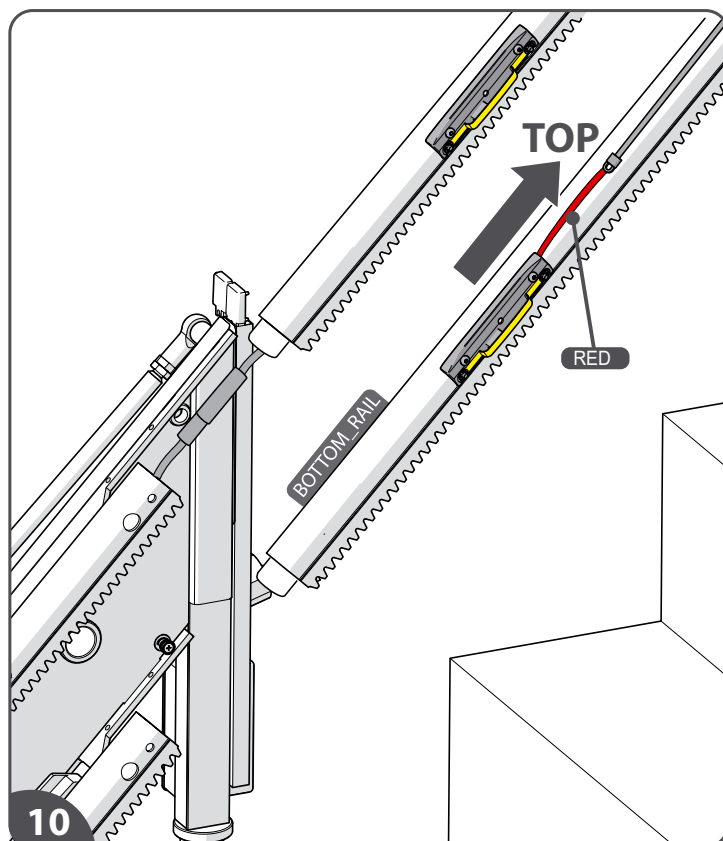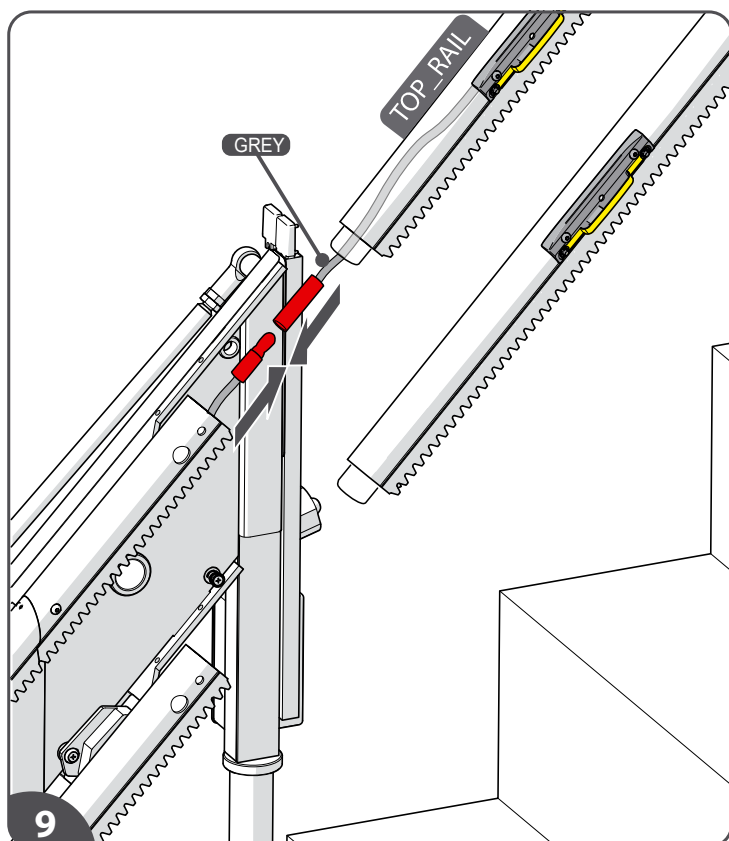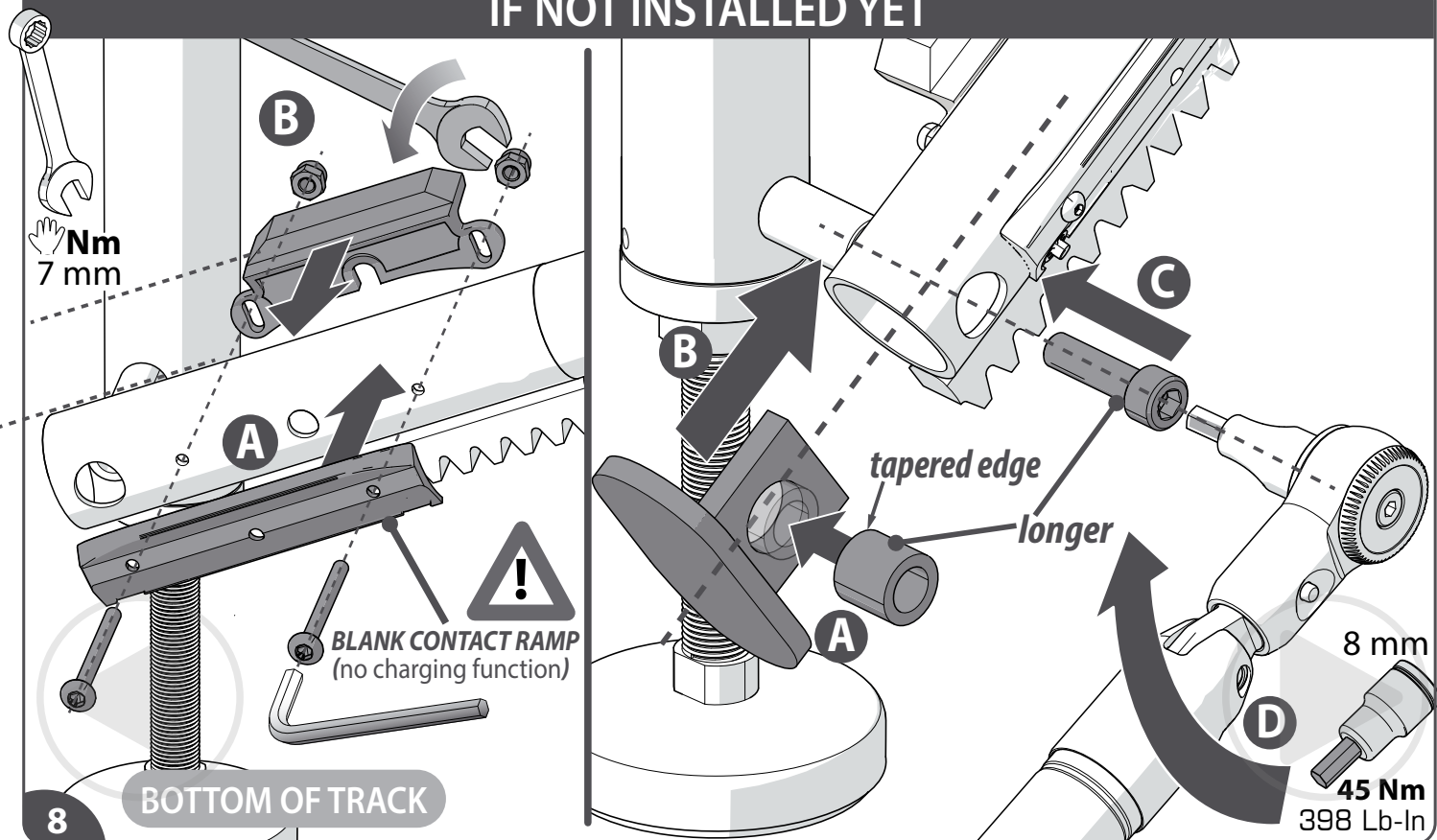[8]If not done yet: install the limit bracket, a blank contact ramp (no charging function) and the end stop on the lower tube at the bottom of the track.

[9]The grey wire with a bullet connector must be connected to the communication contact on the top track.

[10]In the **BOTTOM** tube connect a red cable from the +ve charge contact, pull up to the top of the stairs and connect to the upper charge +ve contact.

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
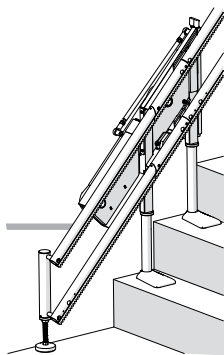
DIPSWITCH SETTINGS

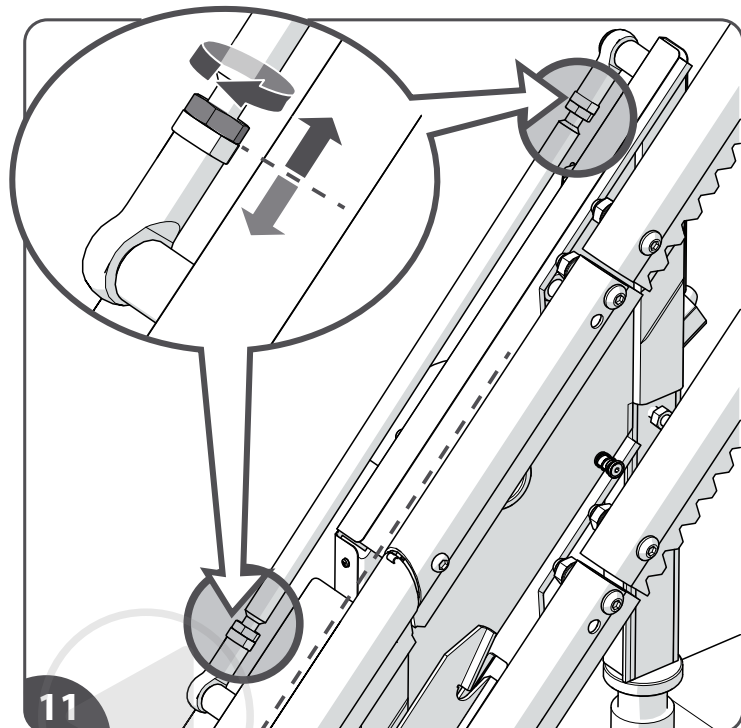
11

**CONTINUE INSTALLATION ►
GO TO PAGE 20 (fig. 2.26) OF
THE INSTALLATION MANUAL**

POWER PACK INSTALLATION

D2
1 foldable track
2 right/left
3 activated
4 always deactivated
* Not during installation

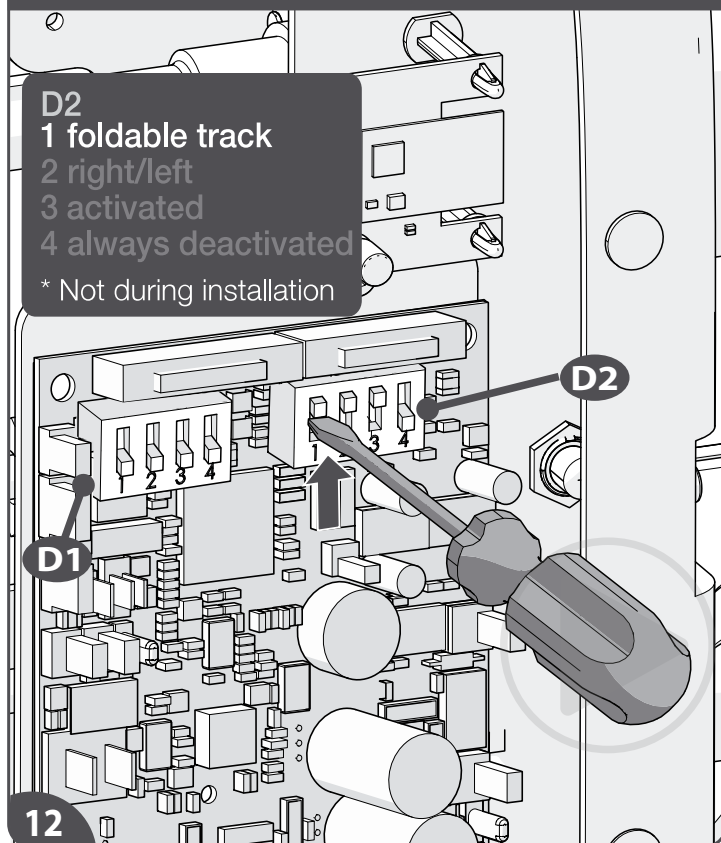12

[11]Adjust the hinge top links to optimize the hinge track joint alignment.
Follow the normal track installation procedure.
[12]Change the setting on the dip switches on the main PCB: a hinged rail(D2-1).

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



IT

B EMERGENCY OPERATION {TRACK}

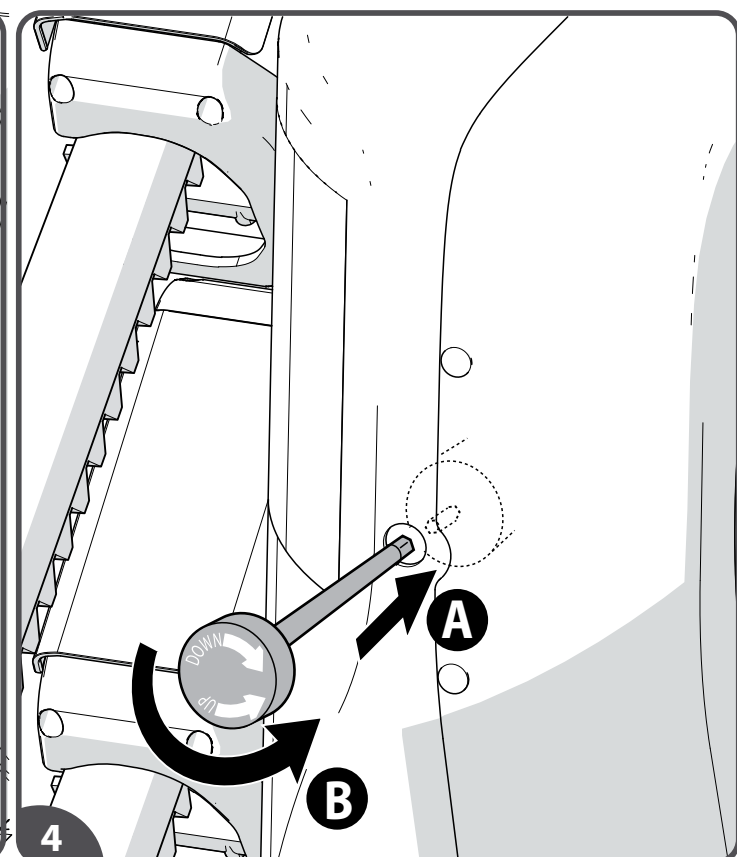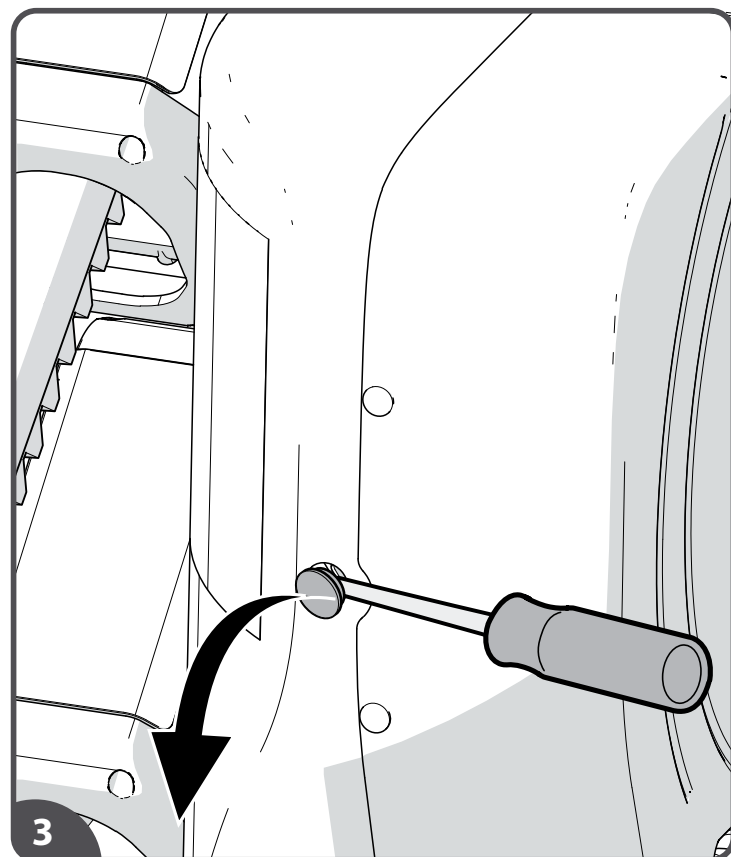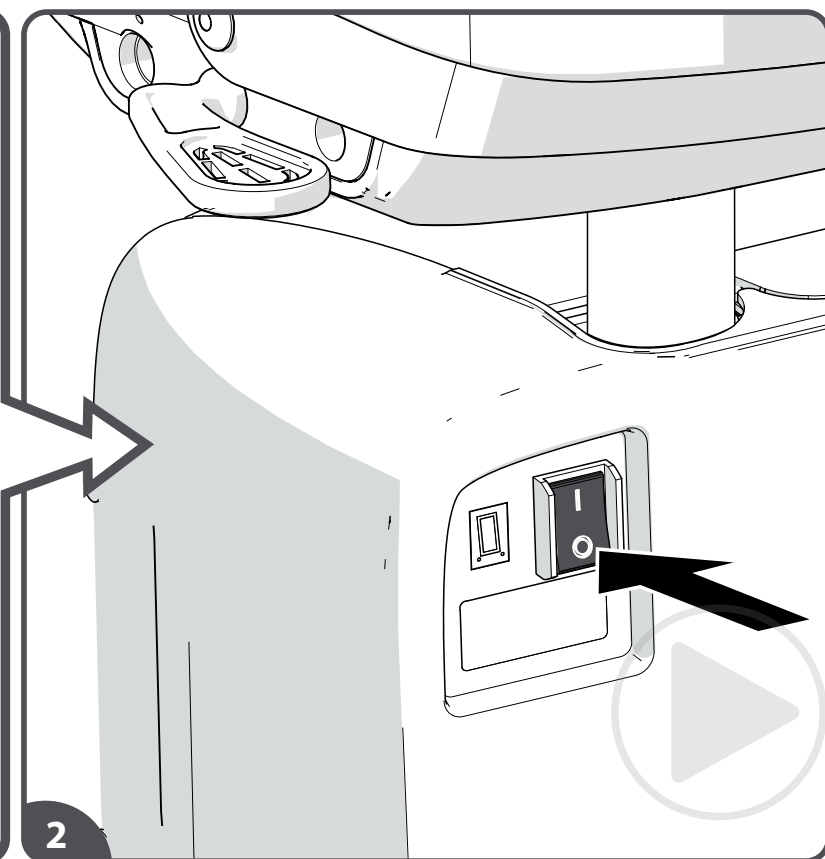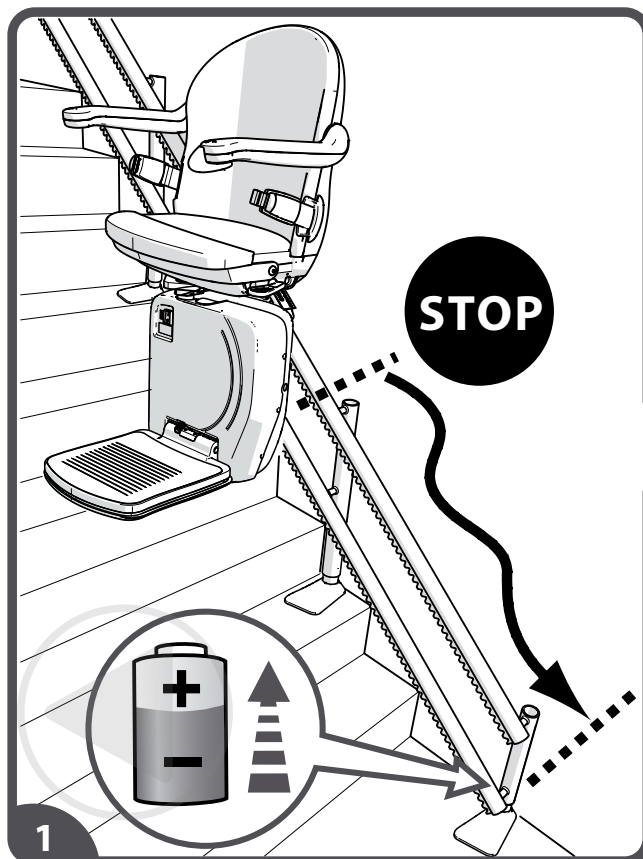- [1] If the stairlift stops unexpectedly and no longer moves any further, the hand winder can be used to move the lift by an authorized and trained person.
- [2] Turn off the lift.
- [3] Remove the cap.
- [4] Turn the hand winder as indicated on the winder. Then wind the lift in the required direction.

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


IT

C1 EMERGENCY OPERATION {SMART SEAT}

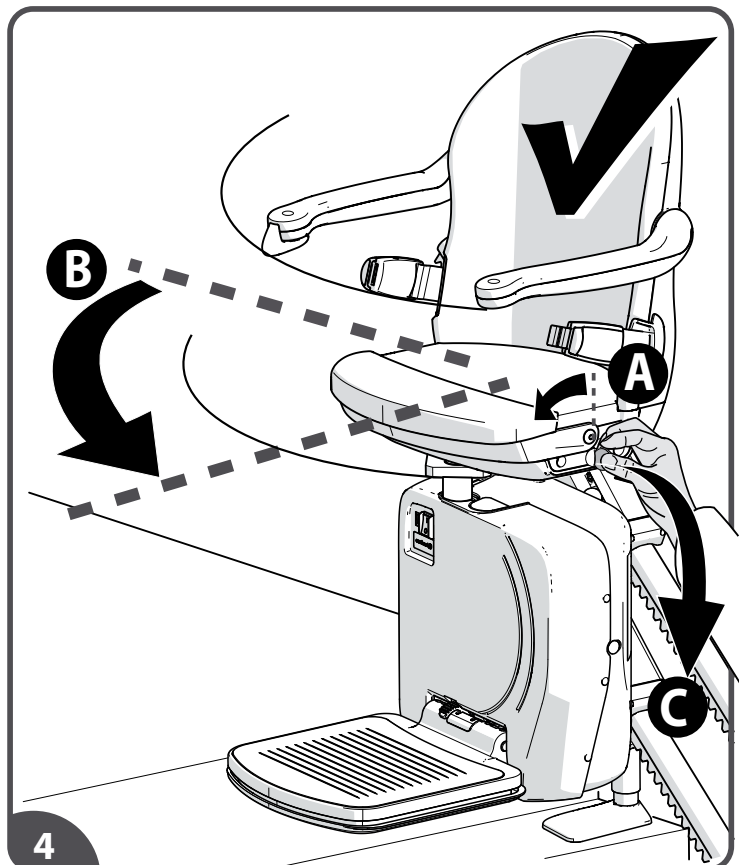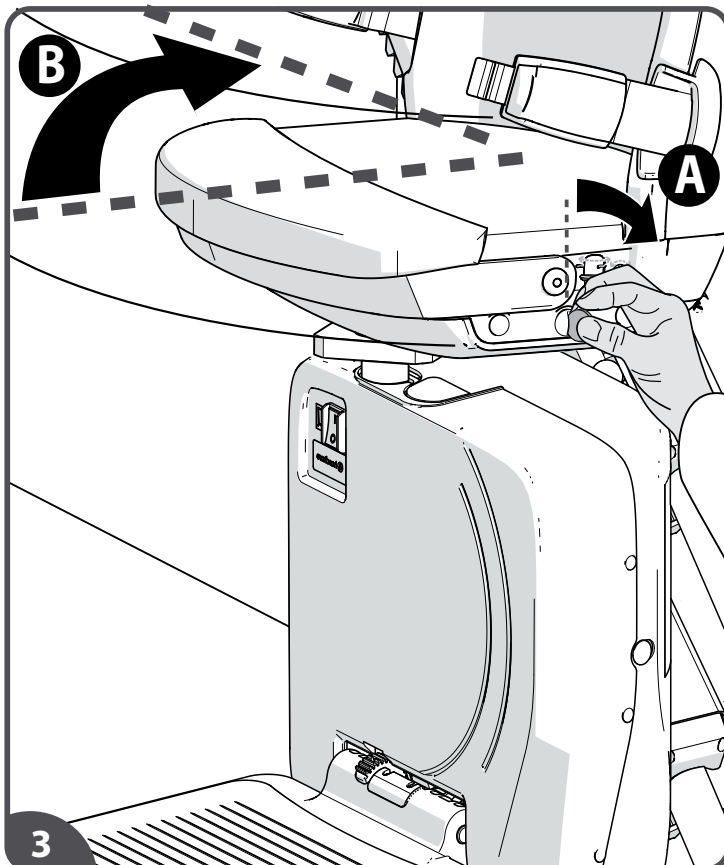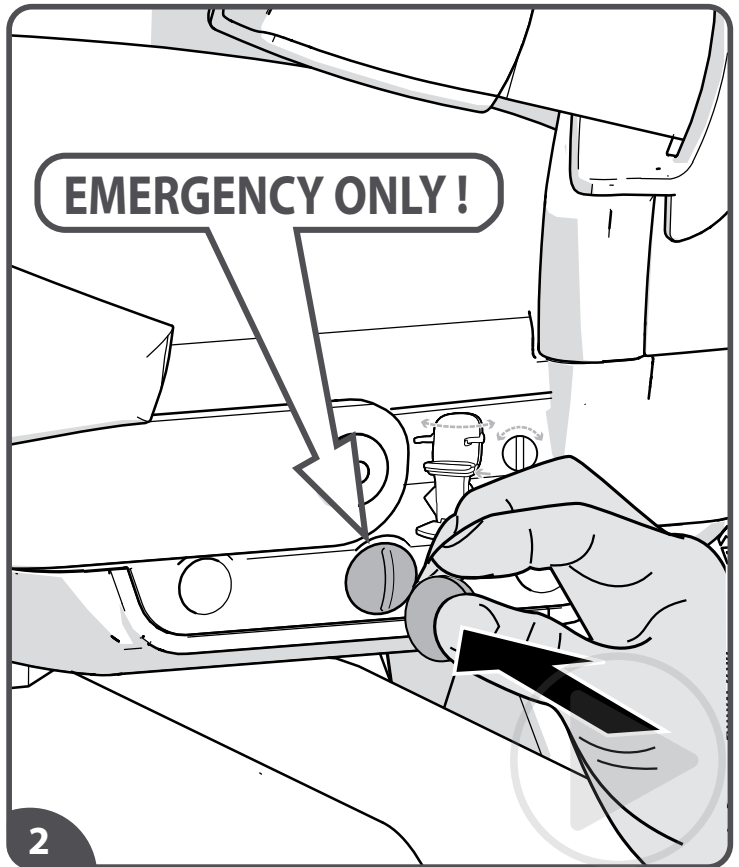- [1]If the power swivel is not functioning it can be swivelled manually. Do this only in an emergency.
- [2]Insert a coin in the button on the side of the seat.
- [3]Swivel the seat manually.
- [4]The seat must be returned to the start position. Remove the coin. Do this before you use the lift again.

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



IT

C2 EMERGENCY OPERATION {STYLE SEAT}

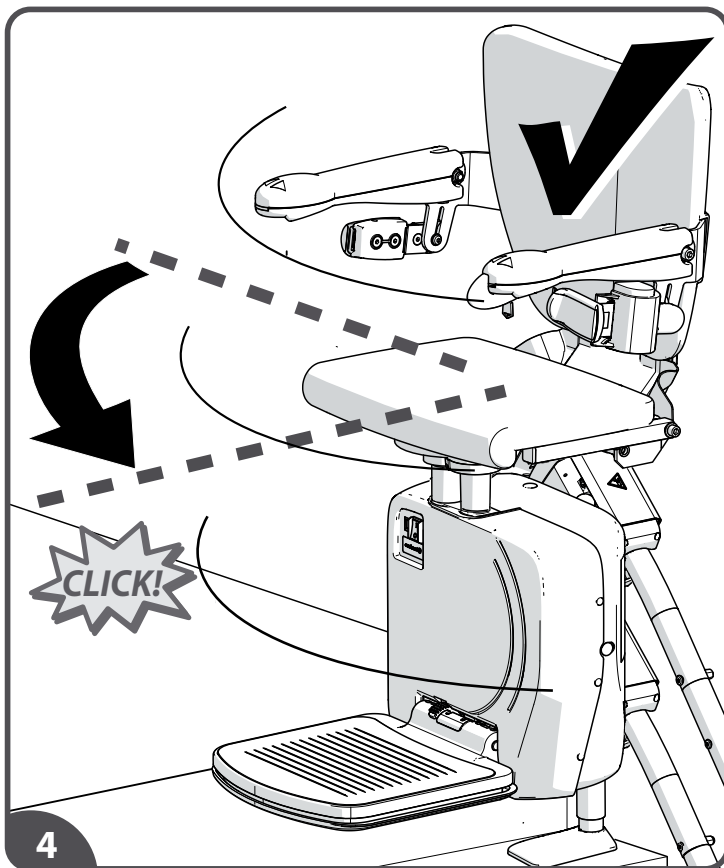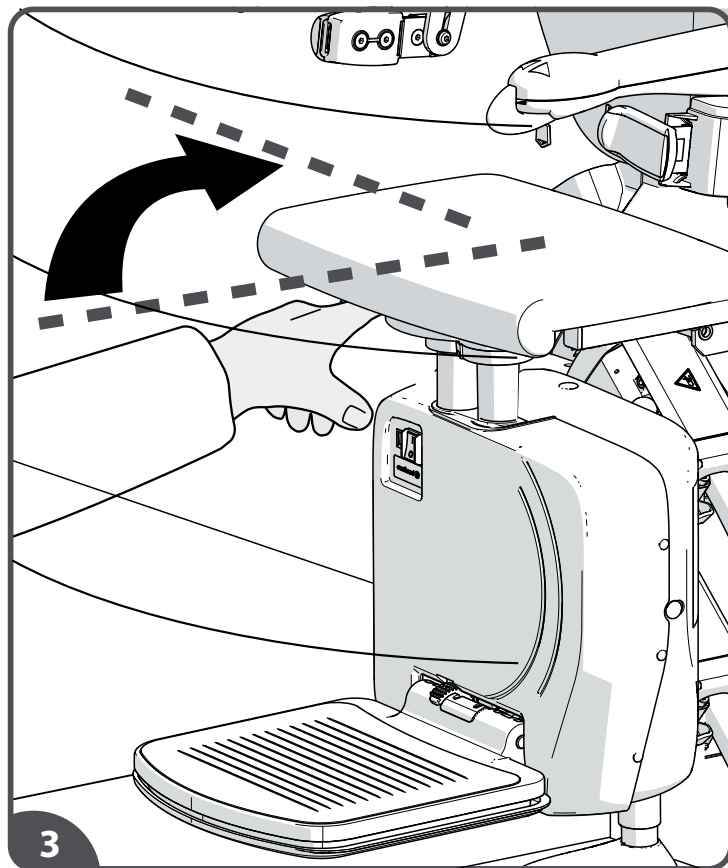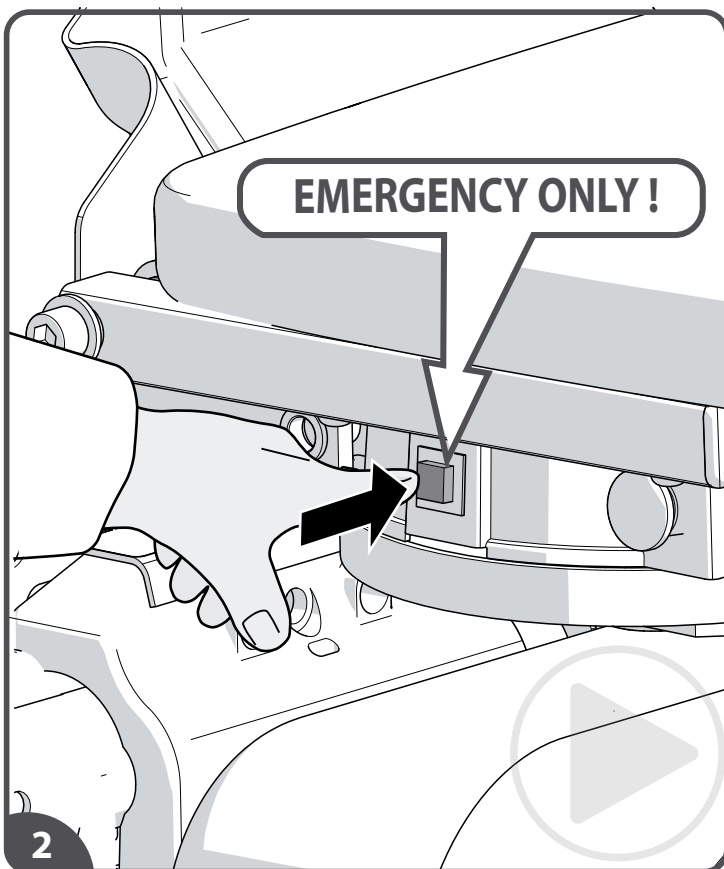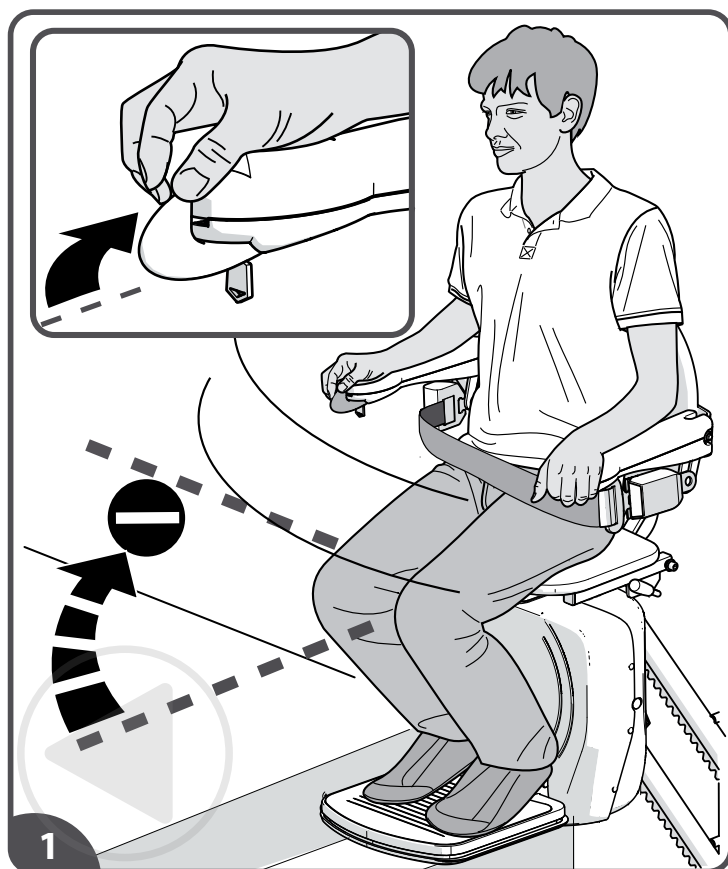- [1] If the power swivel is not functioning it can be swivelled manually. Do this only in an emergency.
- [2] Press the release button.
- [3] Swivel the seat manually.
- [4] The seat must be returned to the start position. Do this before you use the lift again.

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




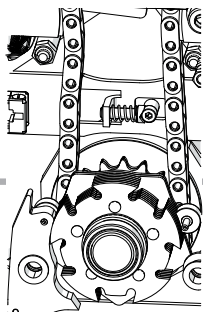
OSG TEST

TOOLS REQUIRED

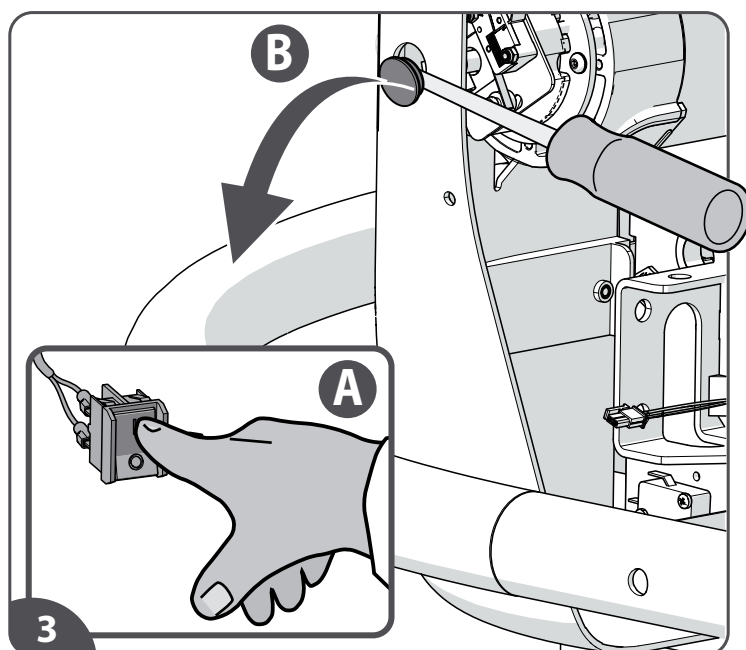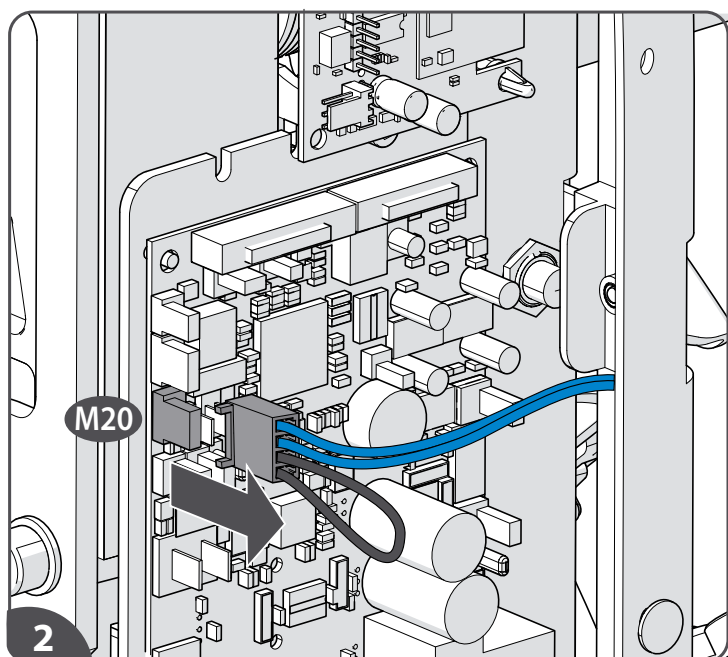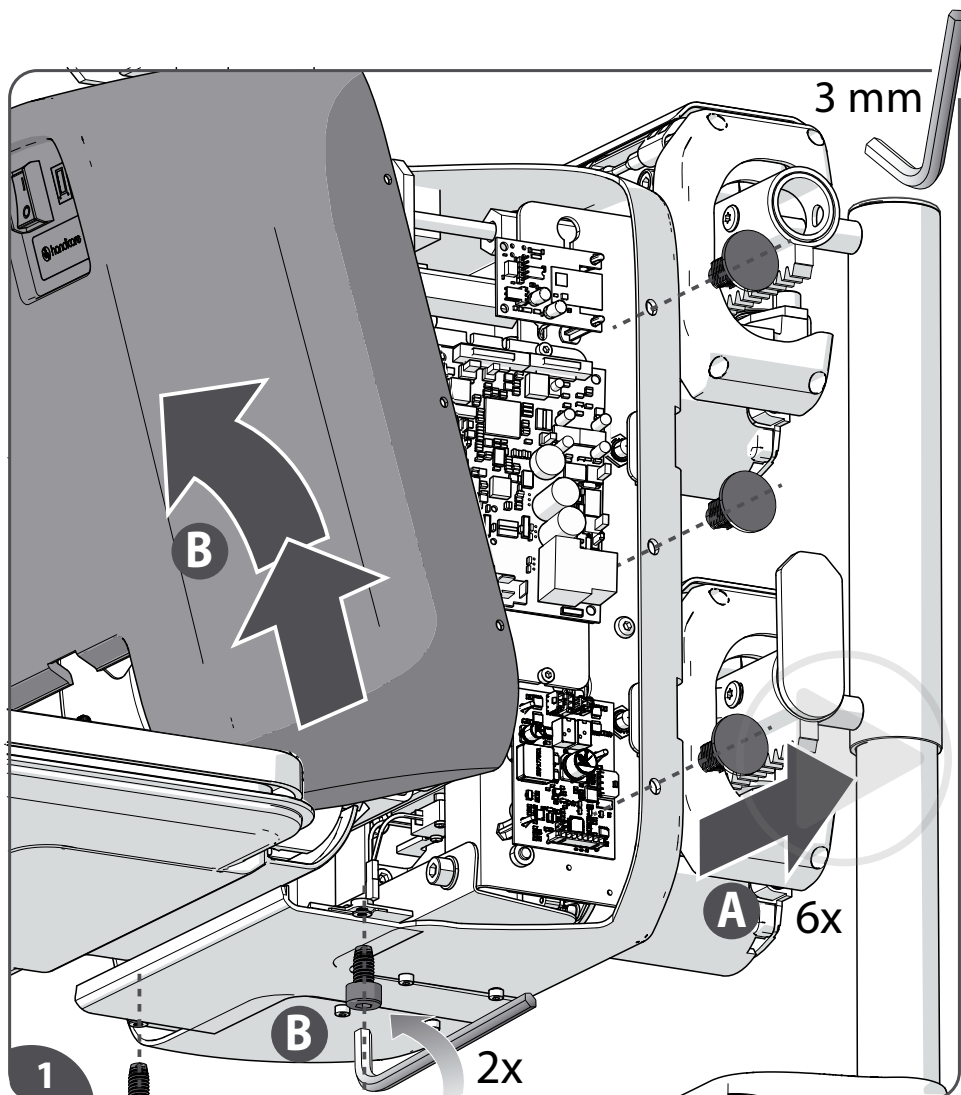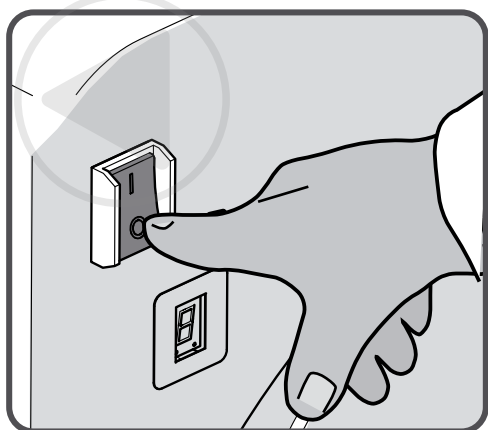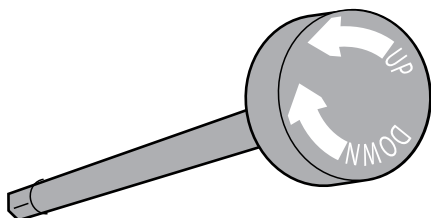- [1] Turn off the lift. Loosen the rivets(6x){A} and screws(2x){B} and remove the front cover{C}.
- [2] Disconnect the 4-way connector (**M20**) on the HAN125 PCB.
- [3] Turn on the lift{A}. Remove the hand winder access hole blanking plug{B}.

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
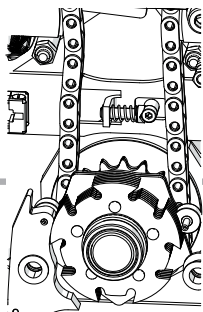
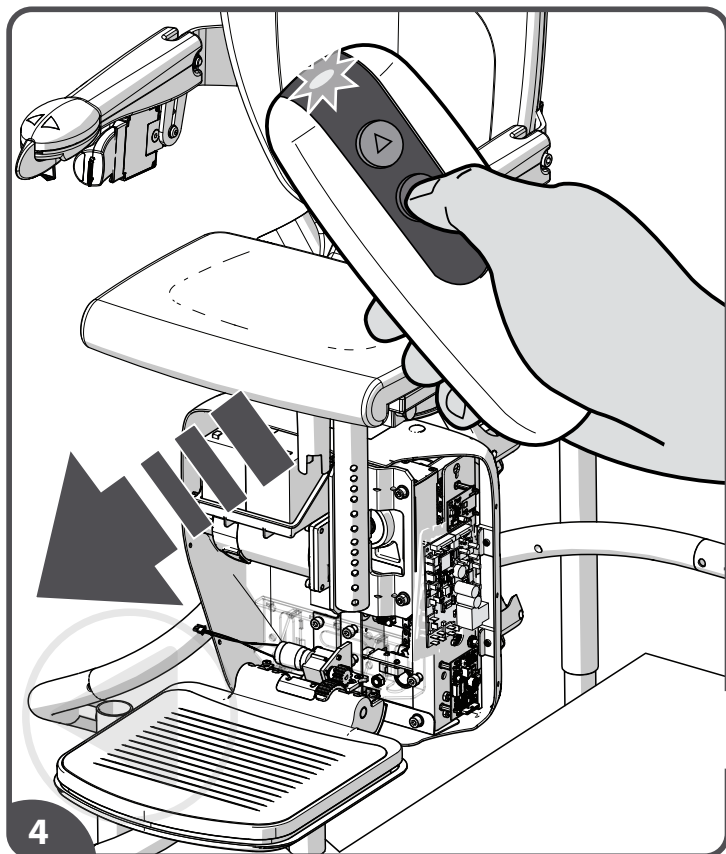
4

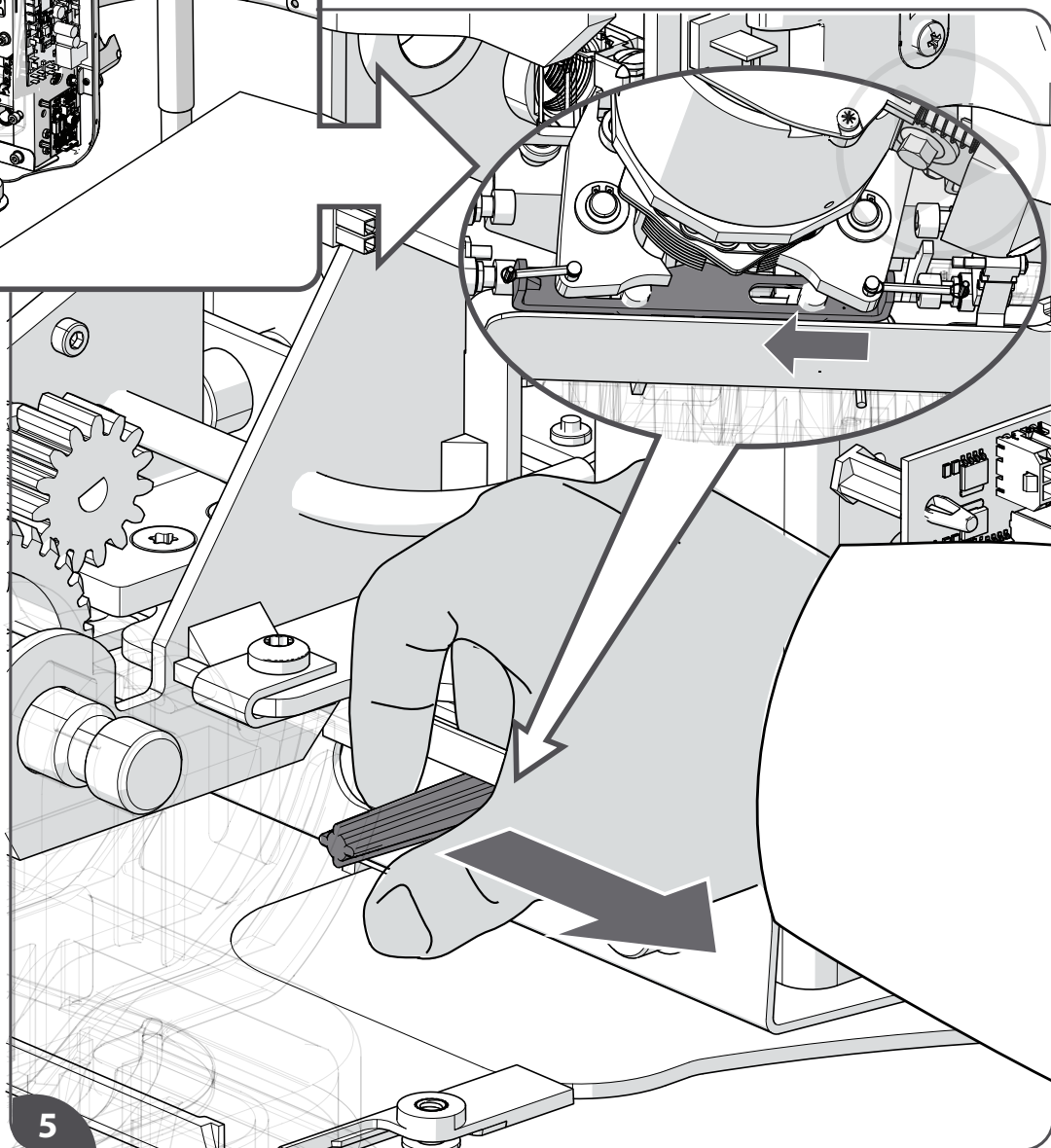5

[4] Firstly, drive the lift down with the remote control. The lift will be travelling at half speed.

[5] Then, using two fingers either side of the underpan plate on the OSG activator, and a thumb on the side of the chassis, pull the activator to engage the OSG. The lift will stop.

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




IT

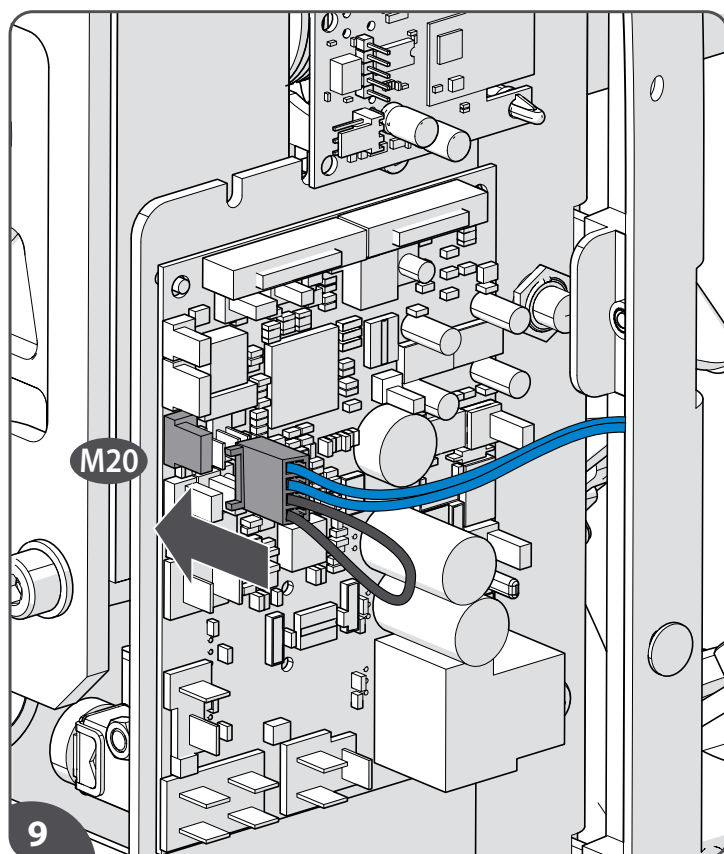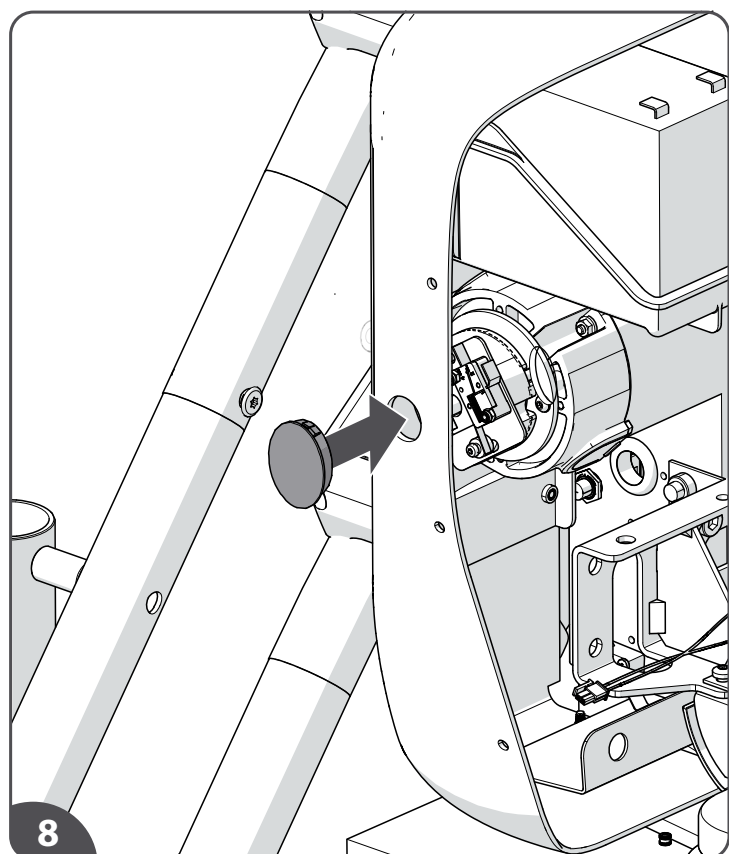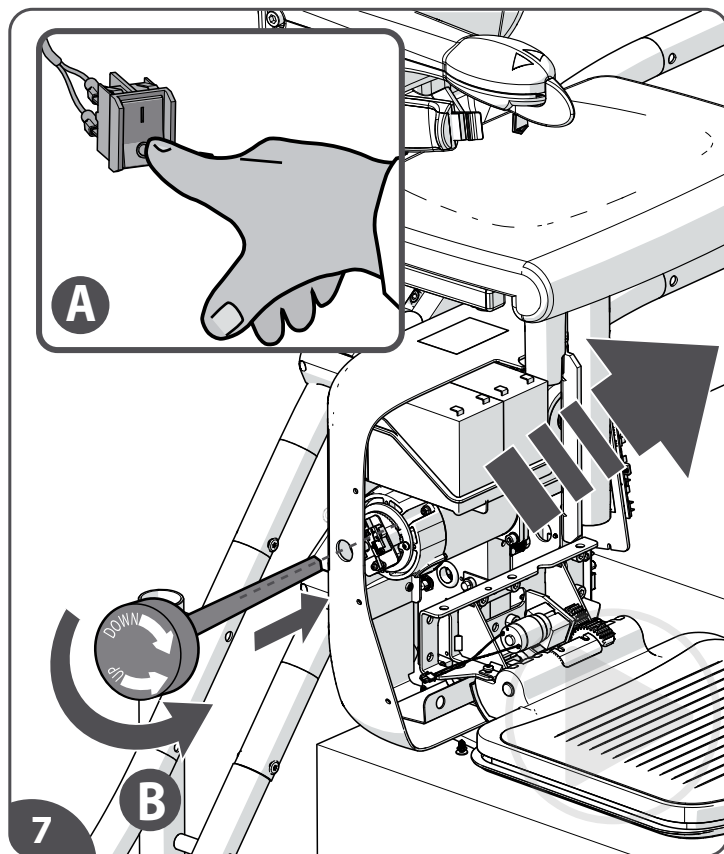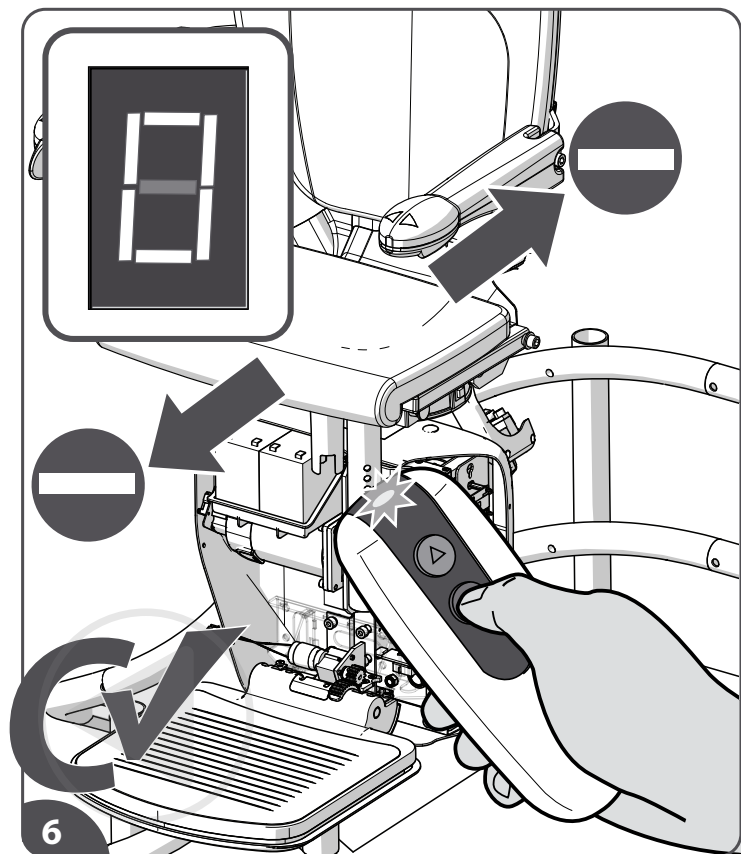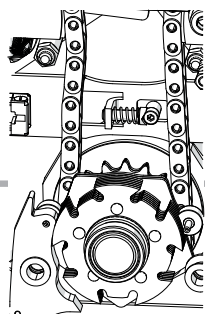- [6] Attempt to drive the lift. Operation of the OSG is confirmed by the lift not driving and Code **"0"** on the display on the frontcover (Safeties / OSG active).
- [7] Turn off the lift{A}. Use the hand winder to move the lift up a short distance to release the OSG{B}.
- [8] Refit the blanking plug.
- [9] Refit the 4-way connector on the HAN125 PCB.

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




.D.

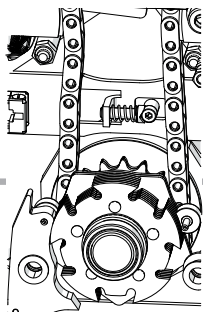OSG TEST

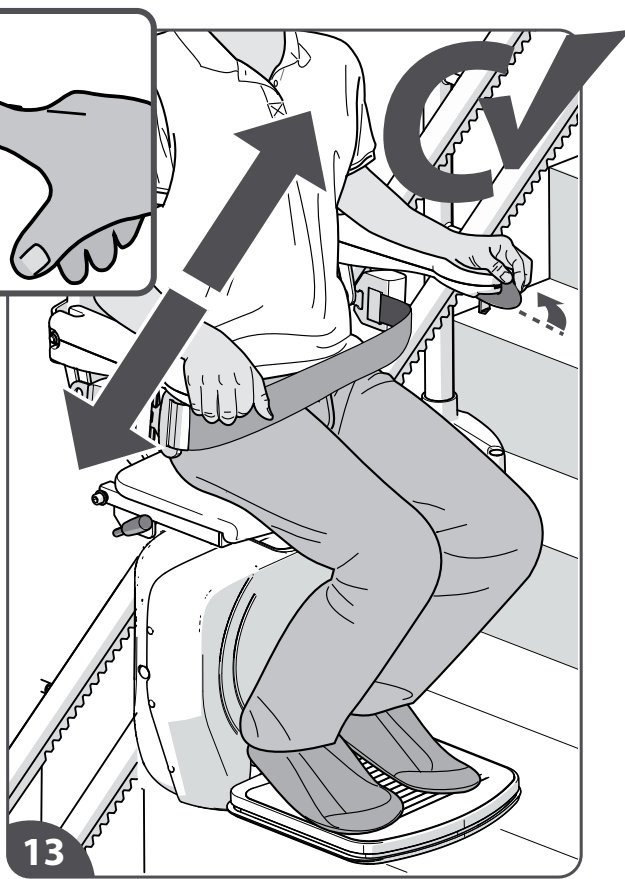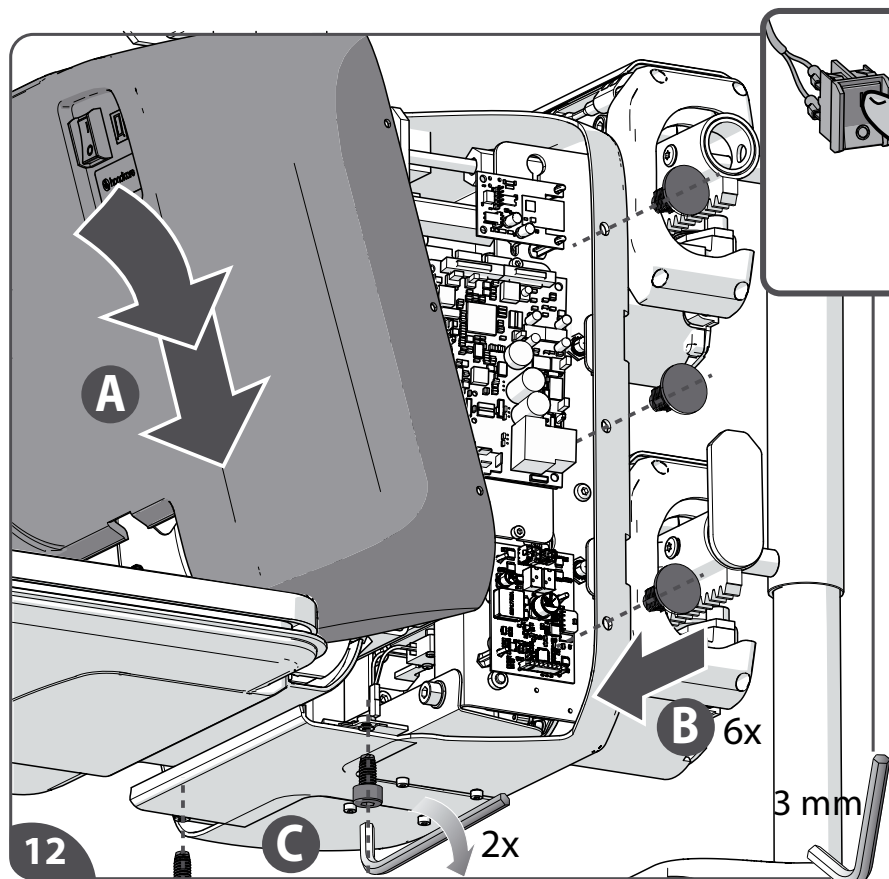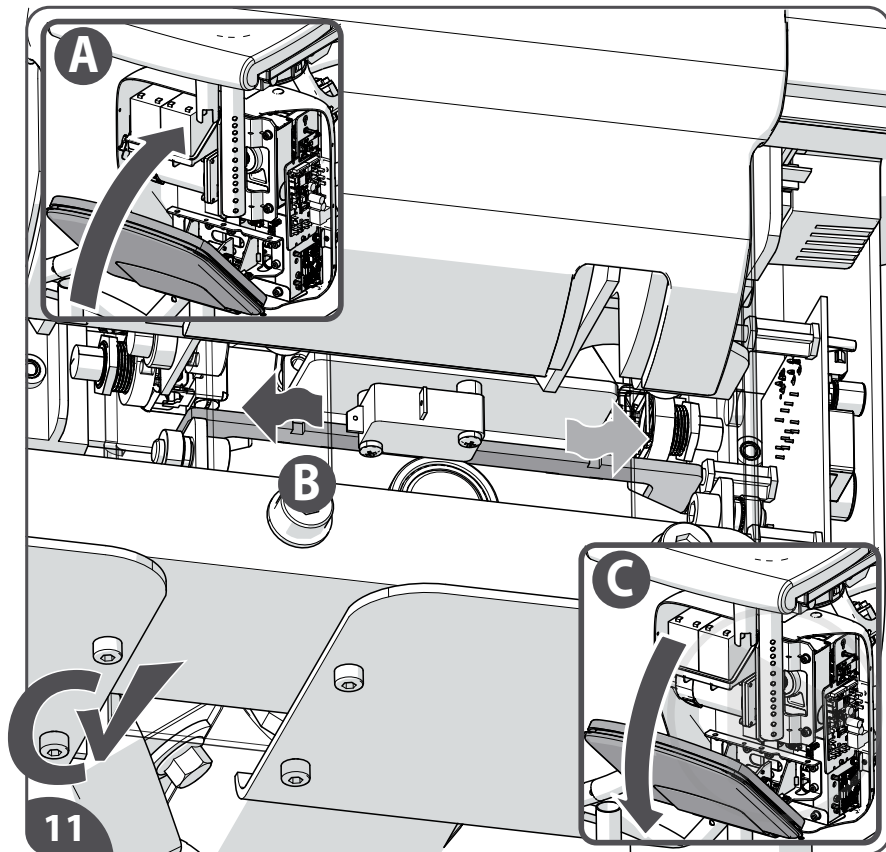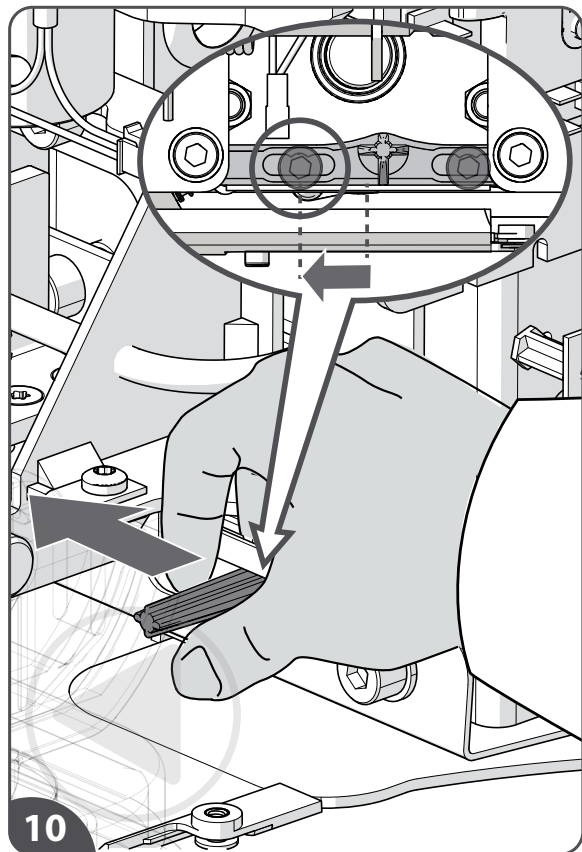- [10]Set the OSG activator back to its neutral position.
- [11]Raise the footplate{A} and check that the OSG Switch Plate is free to move left and right{B}. Lower the footplate{C}.
- [12]Refit{A} and fasten{B} the front cover from the power pack.
- [13]Turn on the lift. Drive the lift UP / DOWN to confirm that the OSG has been released and the lift drives normally.

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




E _____ *INTERMEDIATE STOP* _____

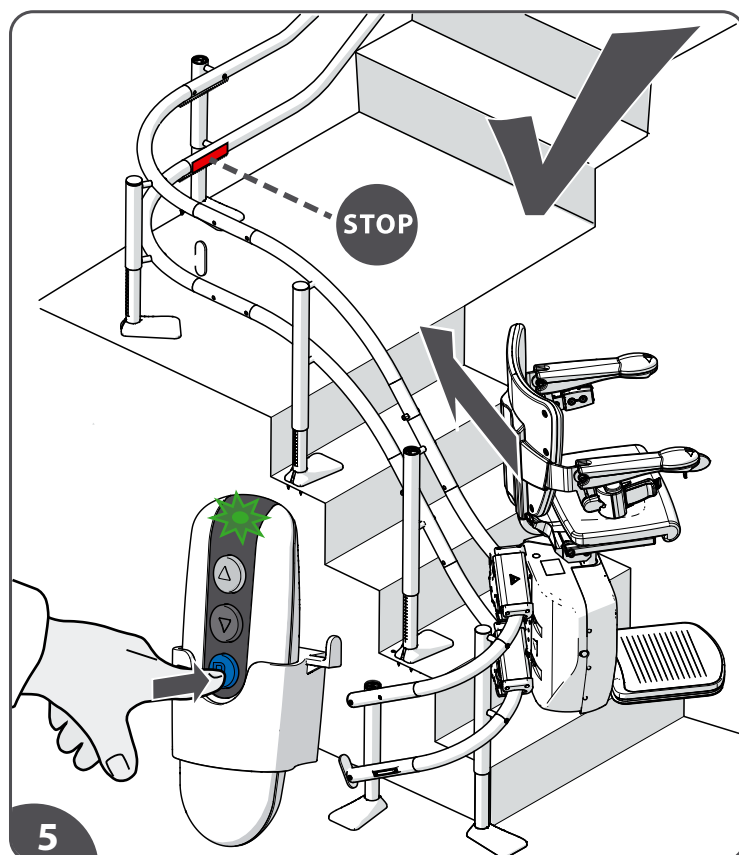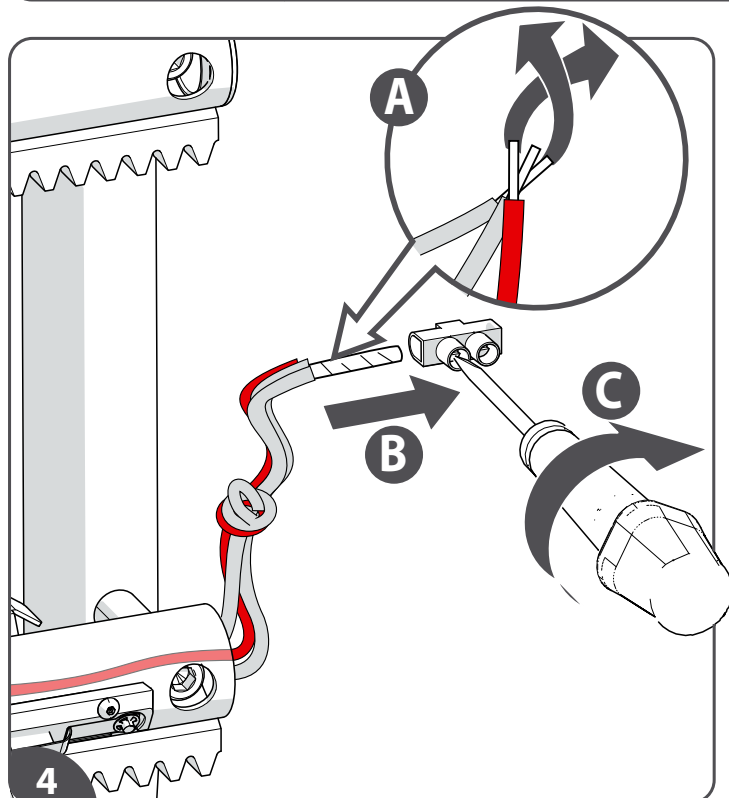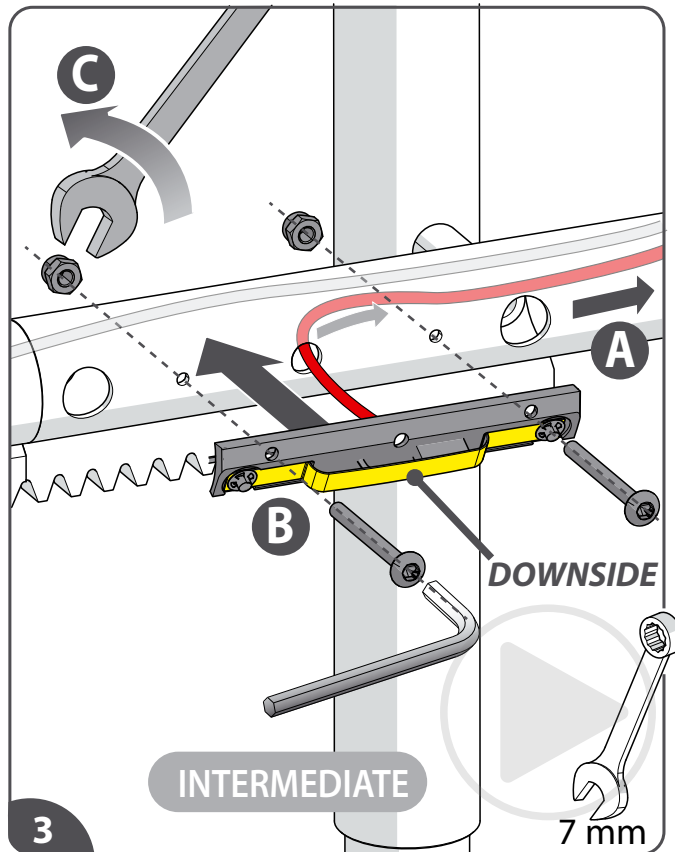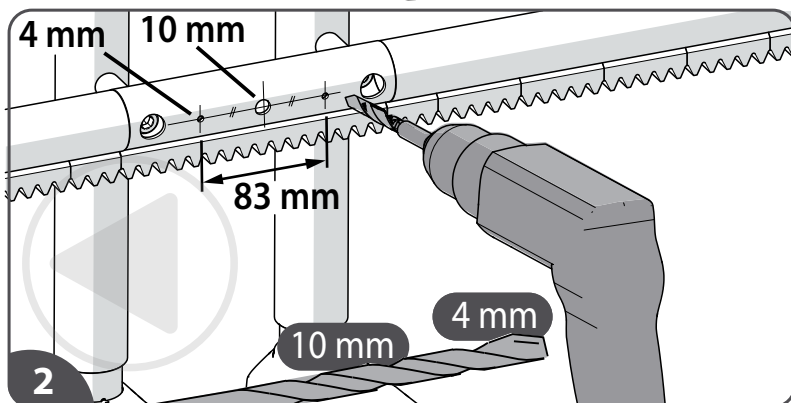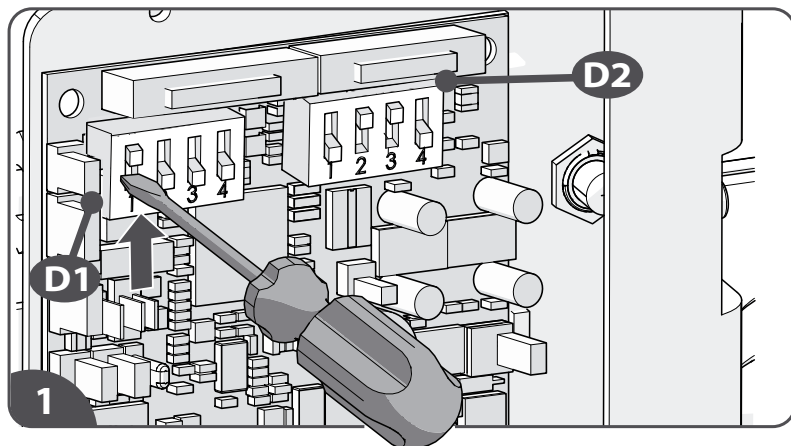CONTINUE NORMAL INSTALLATION

- [1] If not preinstalled, enable the dipswitch setting for the intermediate stop (D1.1).
 - [2] If not preinstalled, drill the holes for the stoppoint on the desired position.
 - [3] Install the stoppoint.
 - [4] Twist the wire of the stoppoint with the two other ends bare wires together {A}. Fasten the wires in a connector block {B+C}.
- Continue with the normal installation.*
- [5] After the installation test the intermediate stop.

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




TOOLS REQUIRED

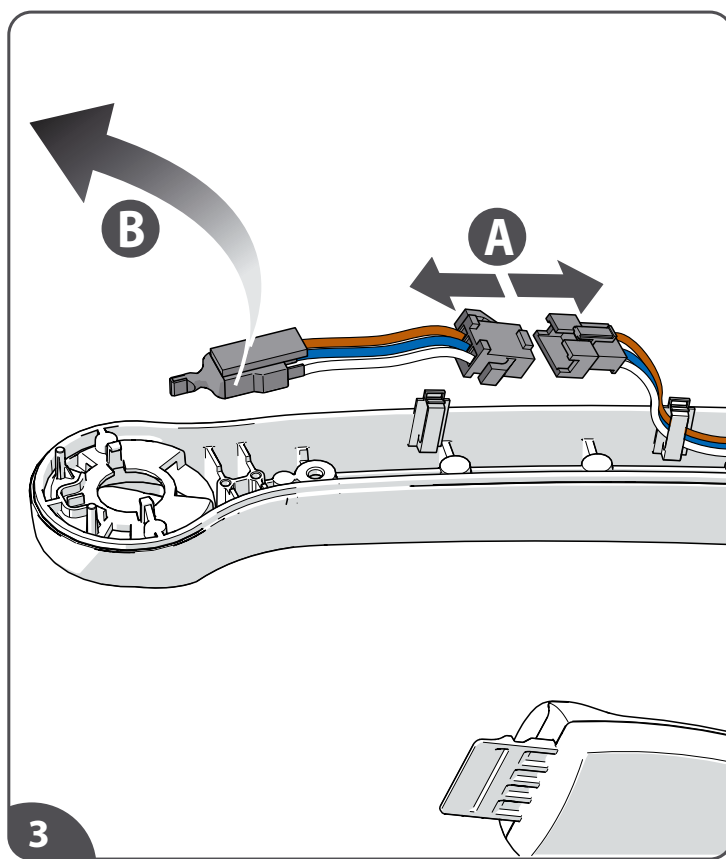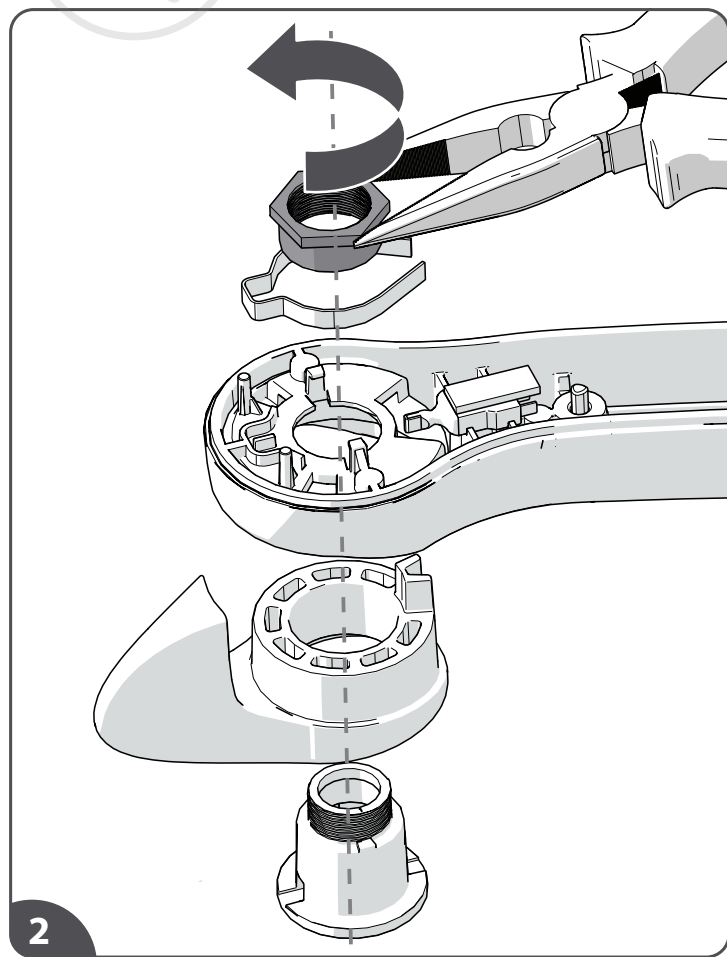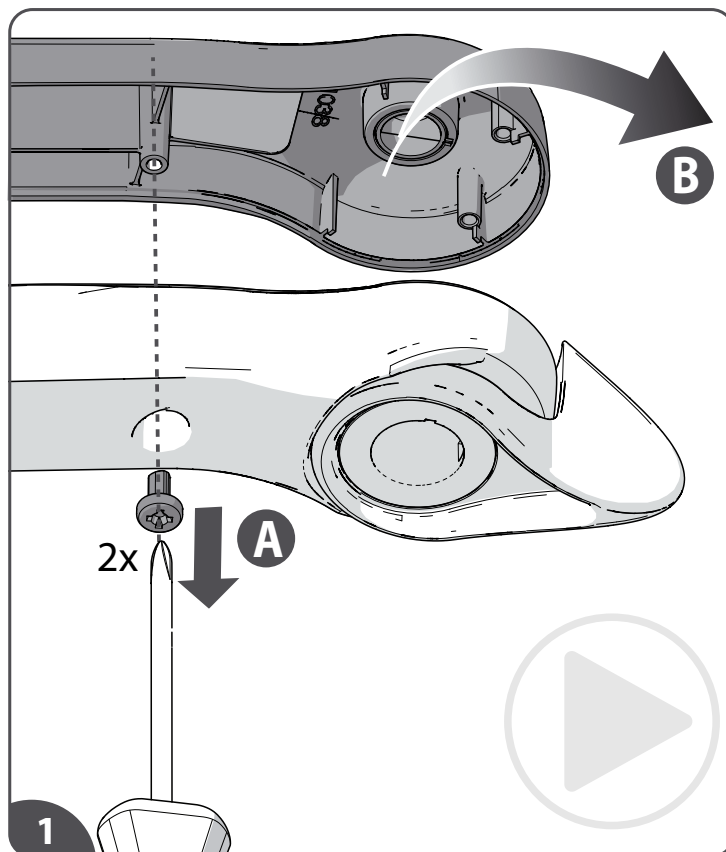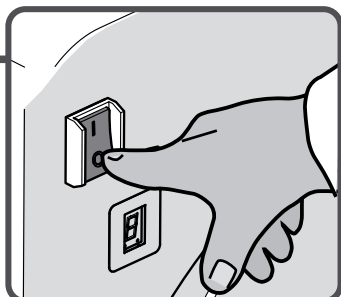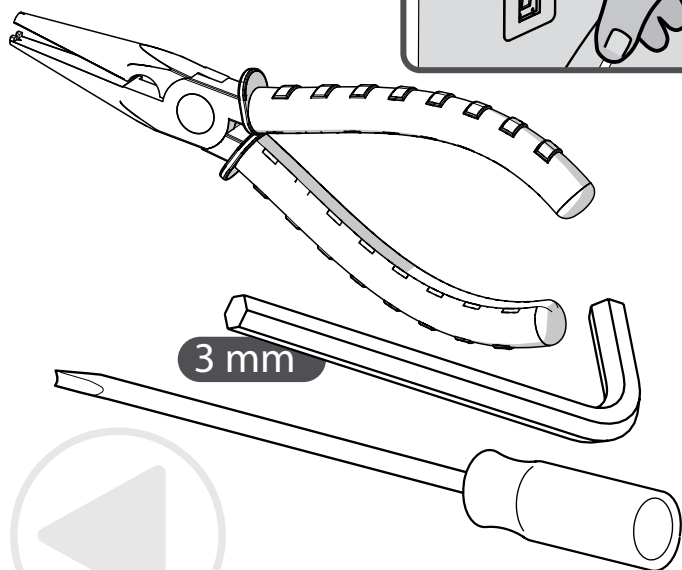Turn off the lift.

[1]Remove top mouldings from current control arm.

[2]Using pliers to undo top-hat nut. Remove toggle, leaf spring and a toggle fastener.

[3]Disconnect toggle loom (Brown, Blue, White){A}. Remove the toggle switch{B}.

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



IT

F ARMREST CONVERSION

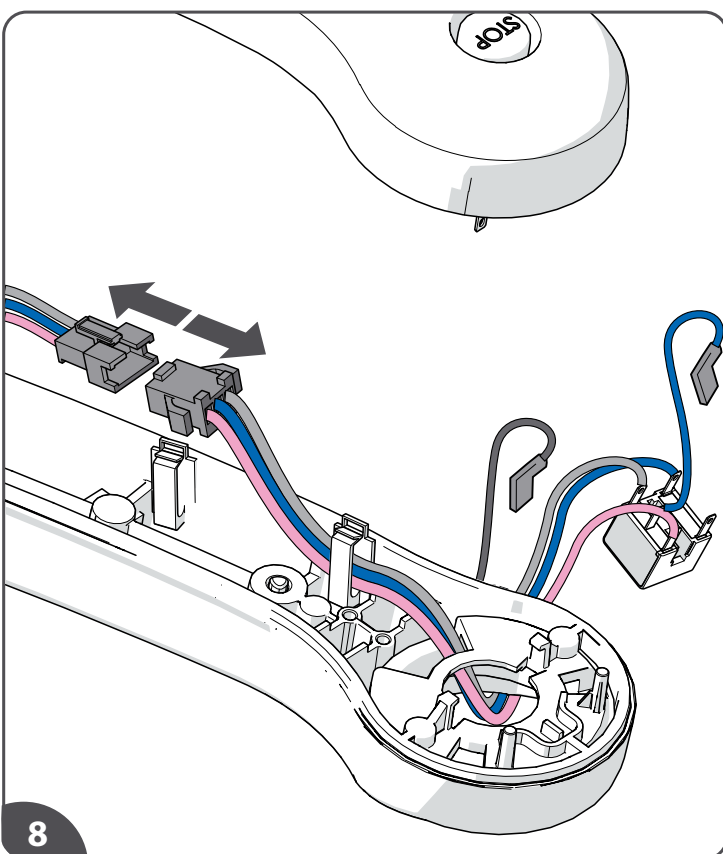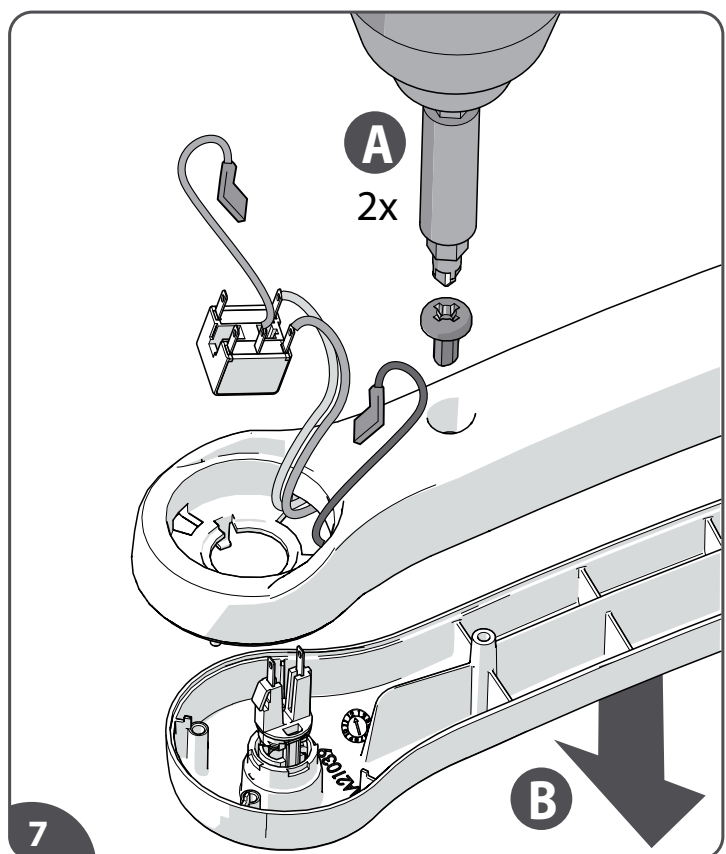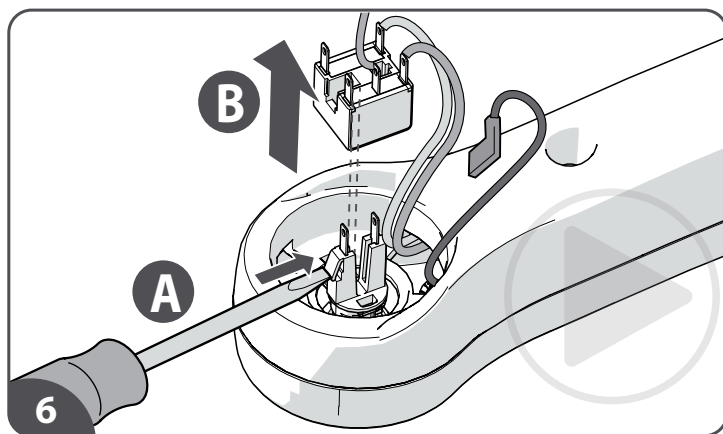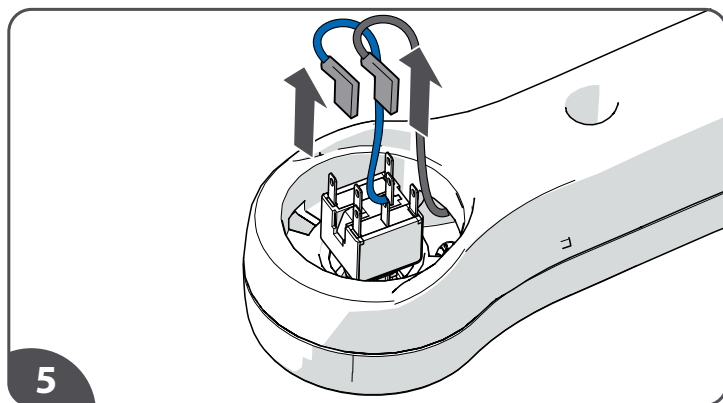- [4] Carefully remove plastic bung from non-control arm by inserting small flat-head screw-driver.
- [5] Disconnect earth loom from pole-block.
- [6] Using small flat-head screwdriver, insert into the notch on stop switch pole block{A} and separate block from push-button{B}.
- [7] Remove screws(2x) which hold top and bottom mouldings together{A}. Remove the top moulding{B}
- [8] Disconnect stop switch pole block loom (pink, blue, grey).

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




IT

F ARMREST CONVERSION

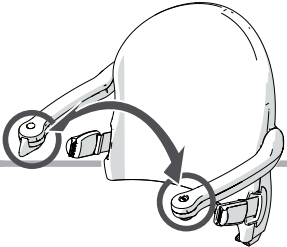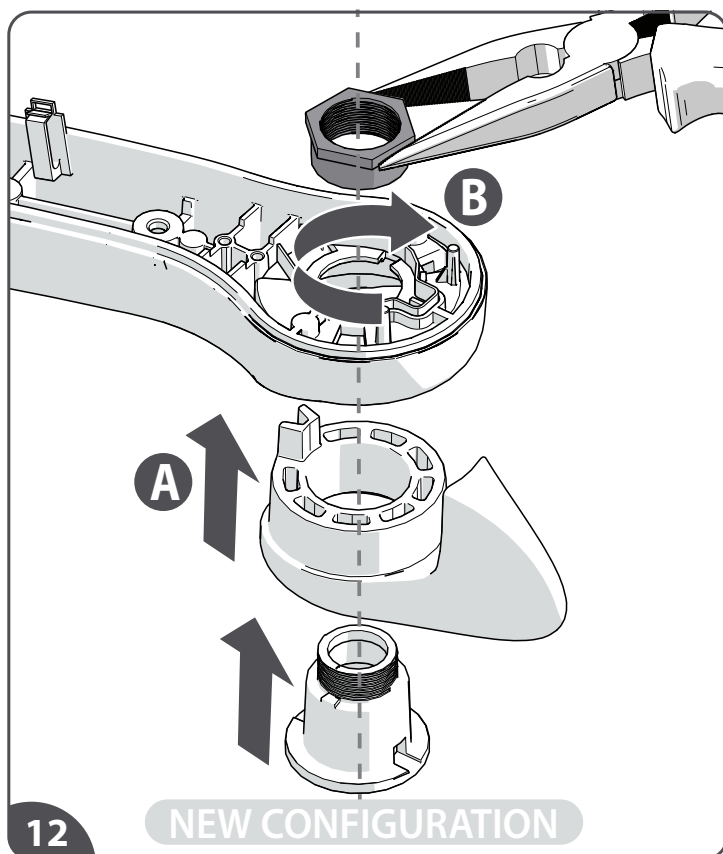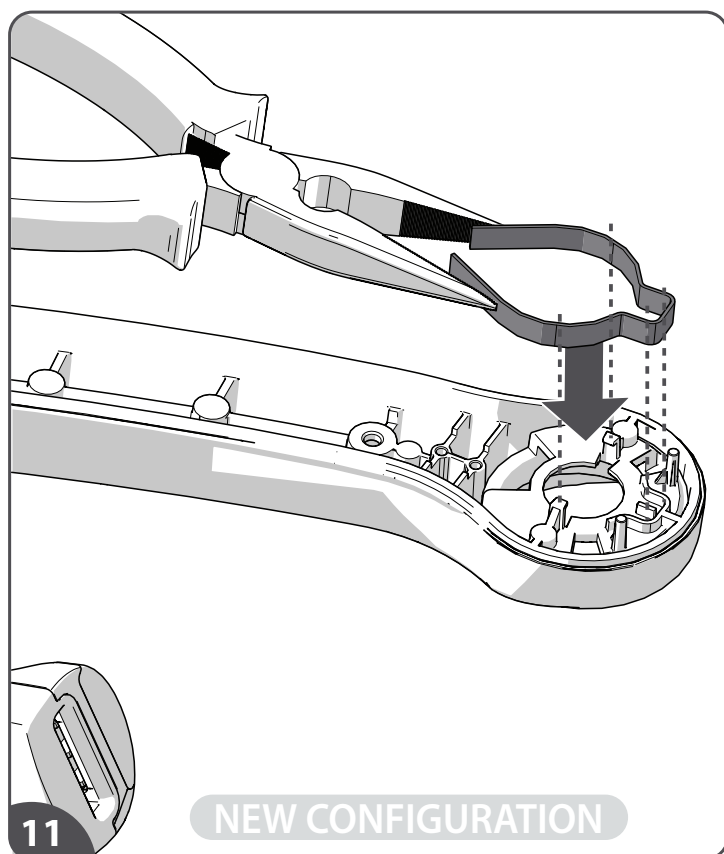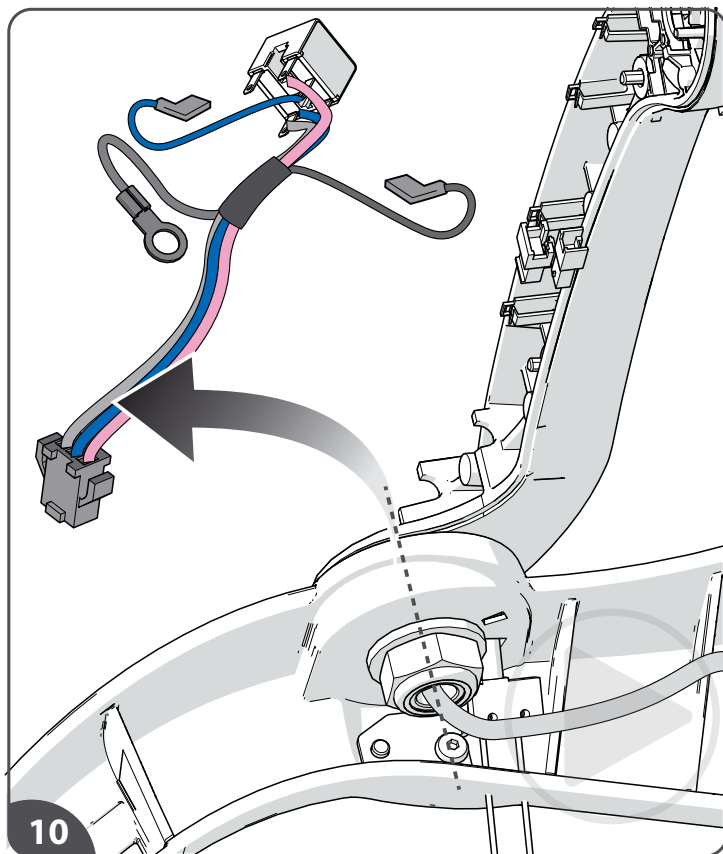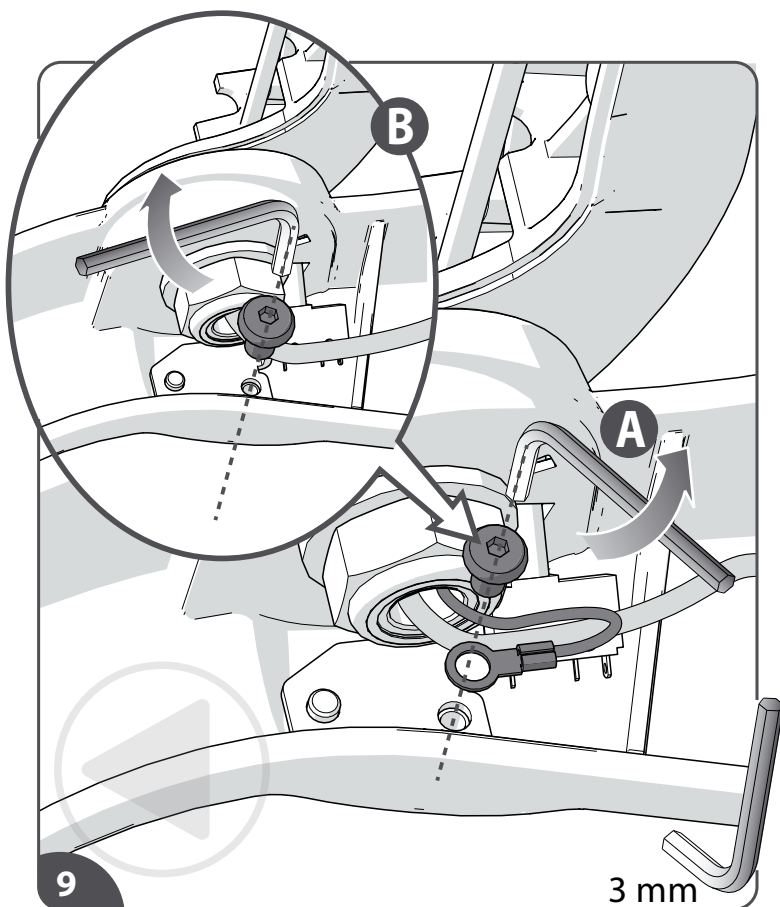- [9] Remove allen key screw which is currently securing earth wire {A}. Retighten the screw once the earth wire has been removed {B}.
- [10] Remove earth wire and pole block from non-control arm.
- [11] Refit the leaf spring.
- [12] Secure using top-hat nut, the toggle and arm lever switch fastener inside new control arm.

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




IT

F ARMREST CONVERSION

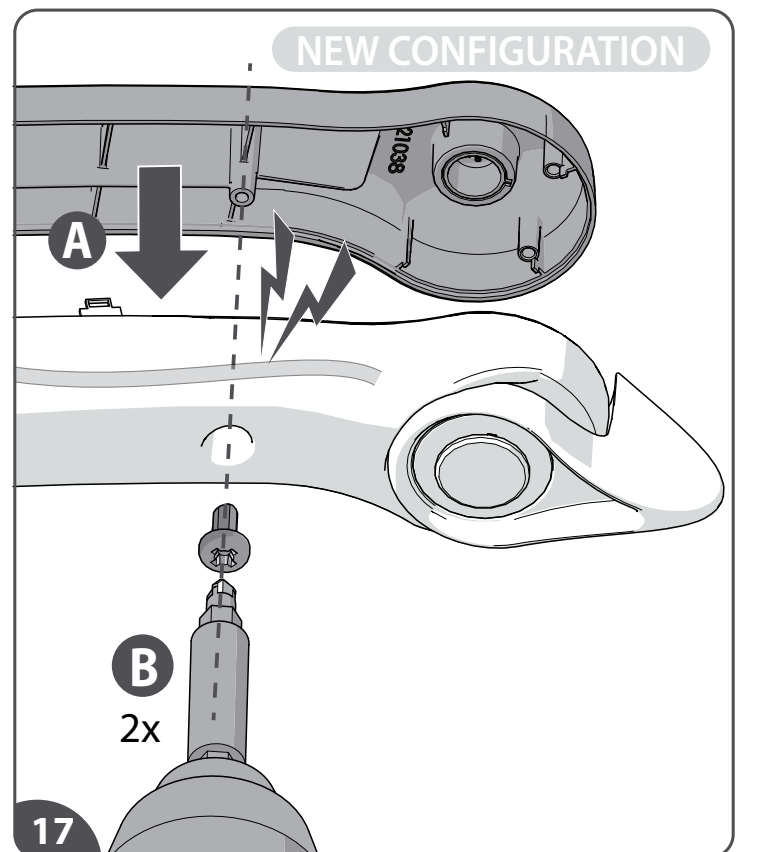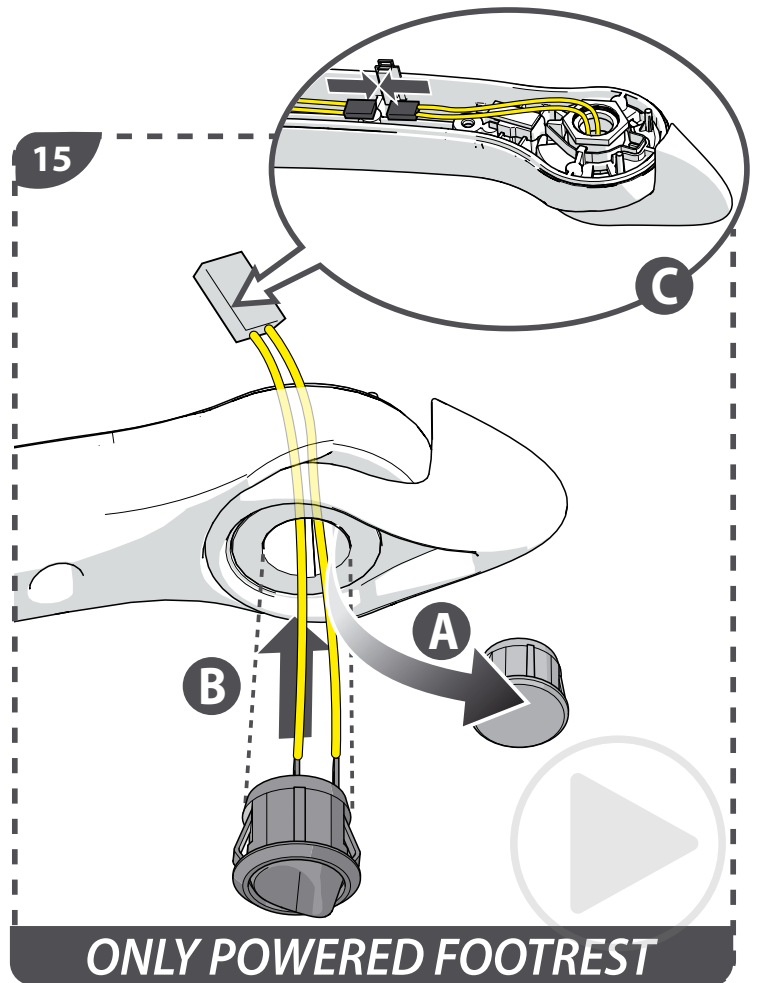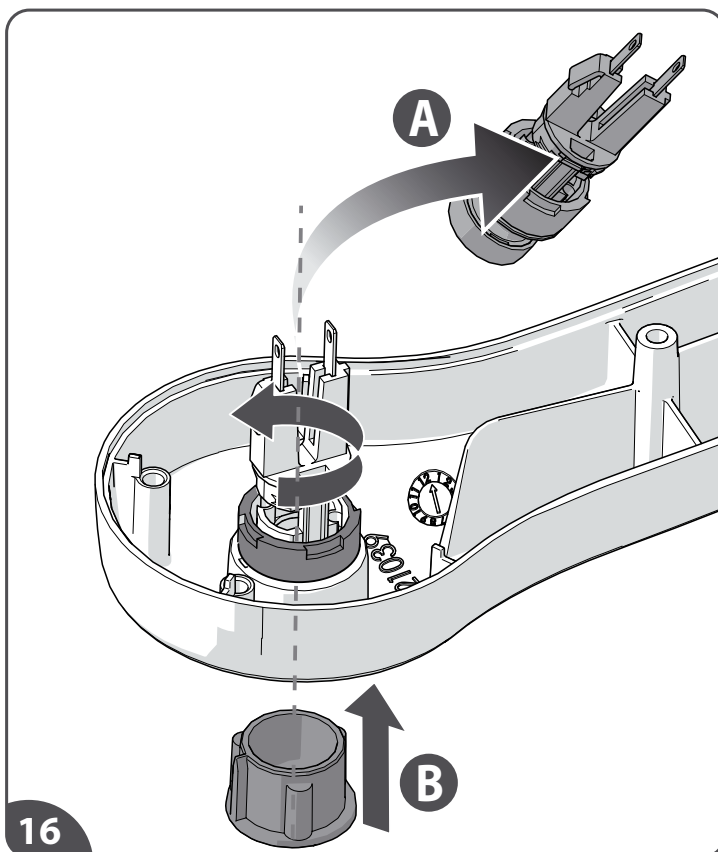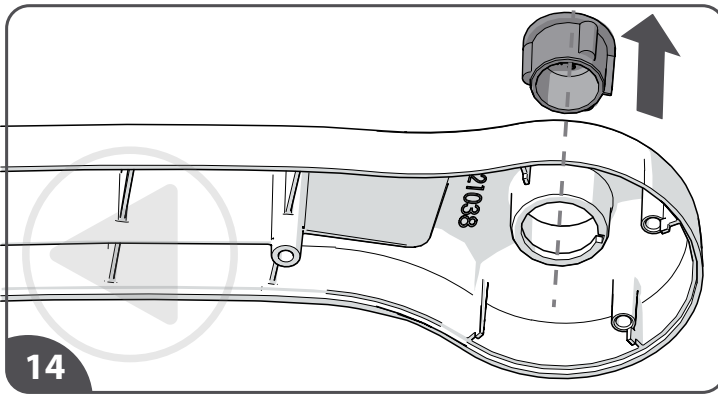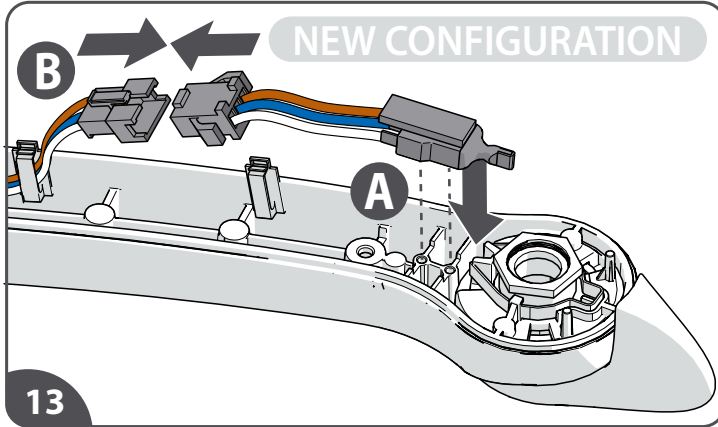- [13] Insert the toggle switch{A} and make the loom connection{B}.
- [14] Remove the dop from the top moulding of the old control arm.
- [15] *Only powered footrest*: remove the cap and install the powered footrest switch{B}
- [16] Remove stop switch from top moulding by undoing plastic washer/nut{A}. Press-in plastic component into new control arm{B}.
- [17] The new control arm is now finished and the top moulding can be secured to the bottom.
Be careful not to pinch the wiring.

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




IT

F ARMREST CONVERSION

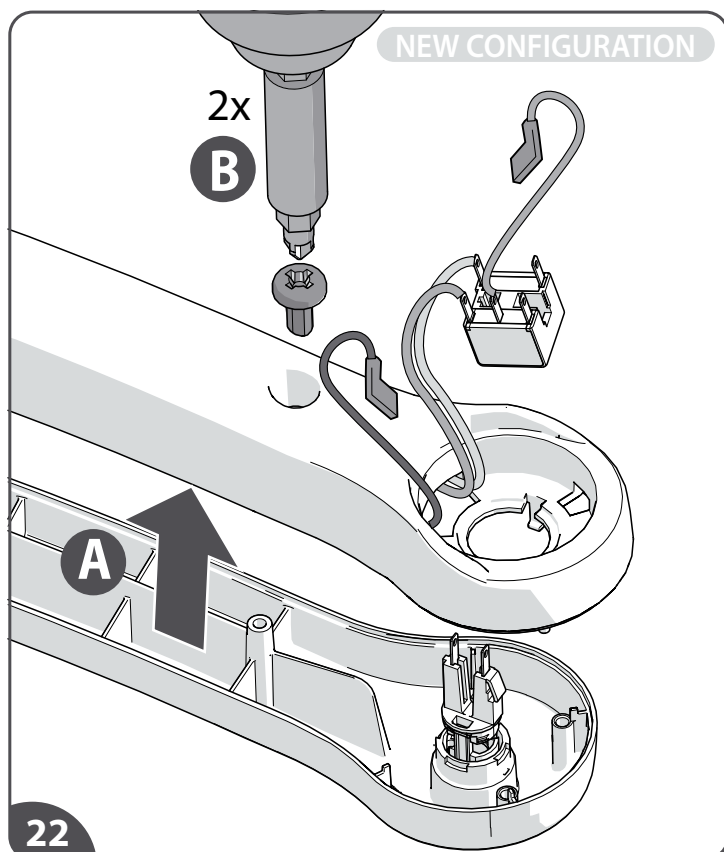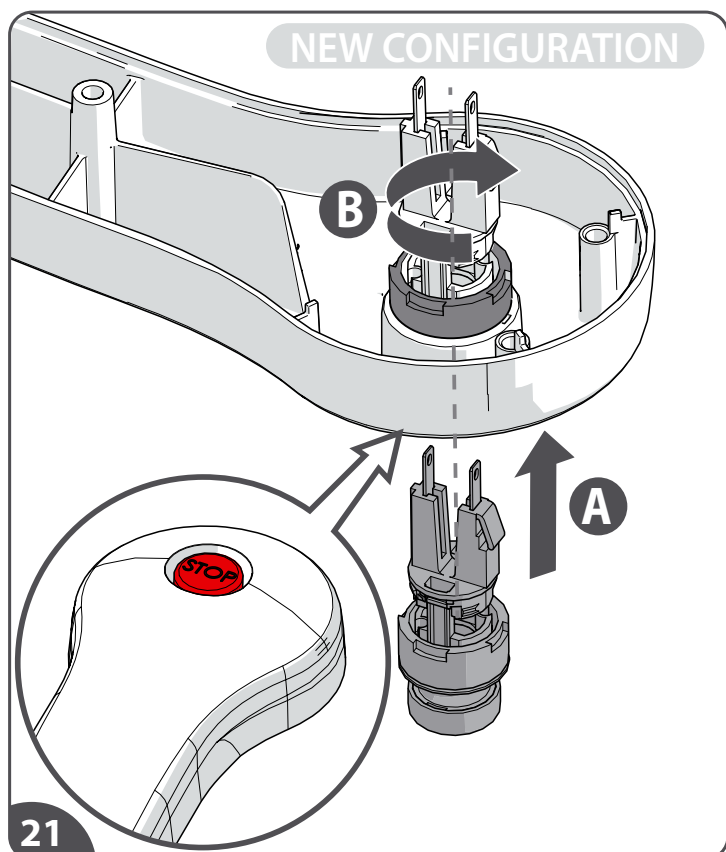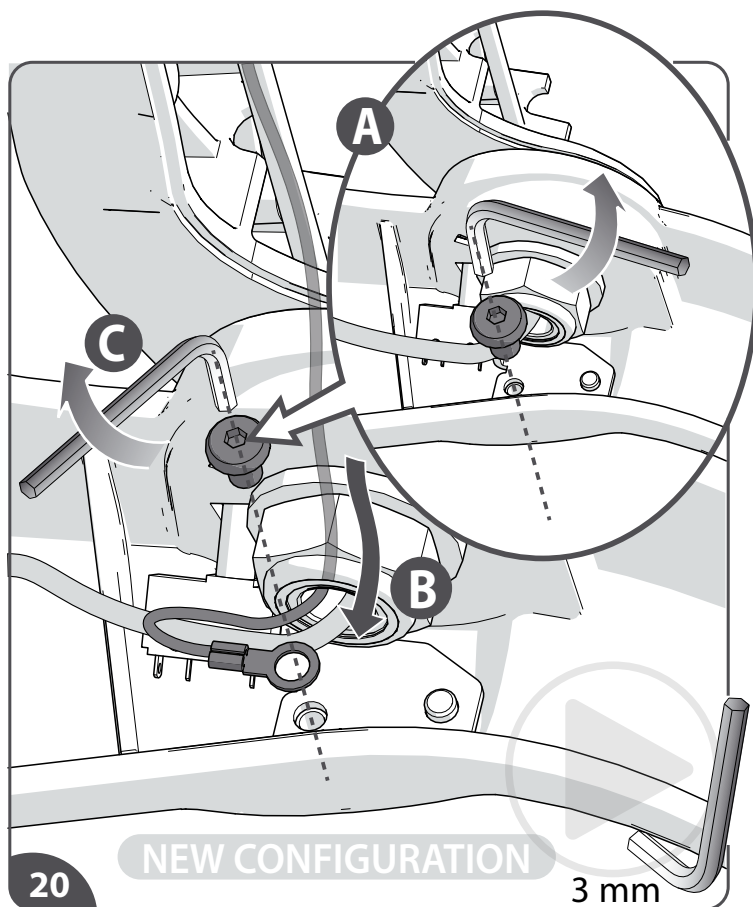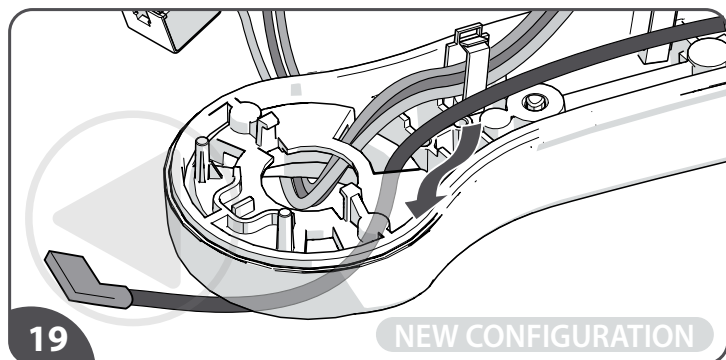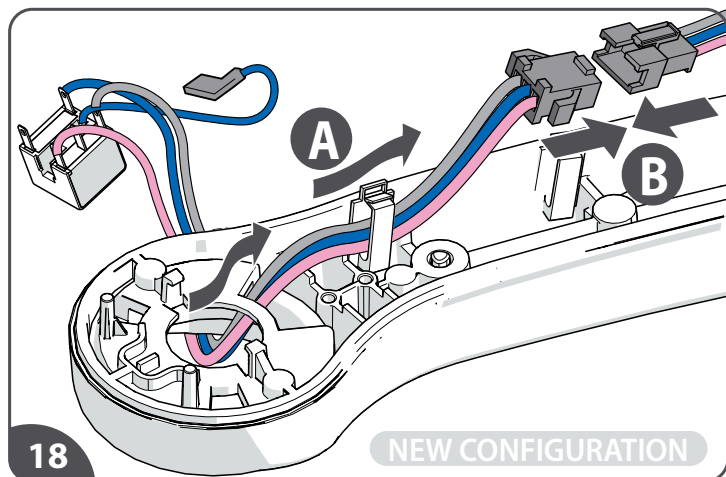- [18]Feed pole block loom through bottom moulding{A}. Make the stop switch loom connection{B}.
- [19]Feed ground wire through insert in bottom moulding.
- [20]Loosen allen key screw{A}. Feed ground wire through new non-control arm pivot{B} and secure with allen key screw{C}.
- [21]Place stop switch into new non-control arm, ensure 'STOP' is in the correct orientation{A}. Secure stop switch with plastic washer/nut{B}.
- [22]The top moulding can be secured to the bottom.

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



IT

F *ARMREST CONVERSION*

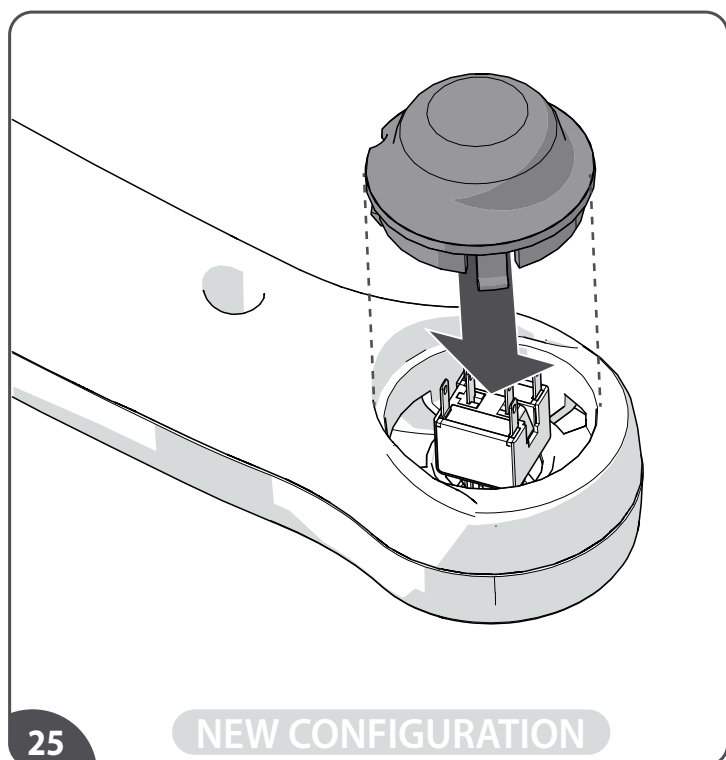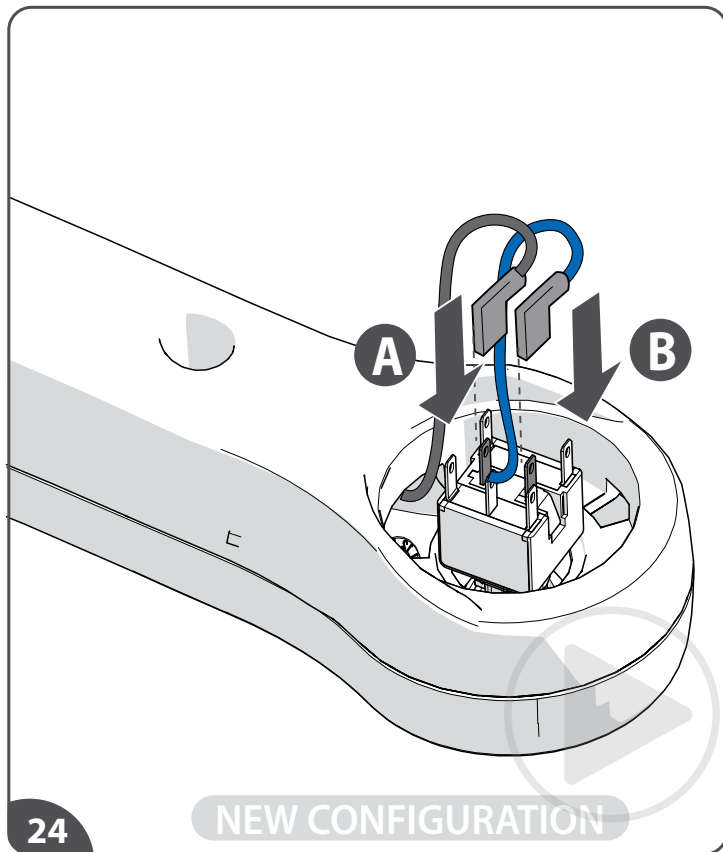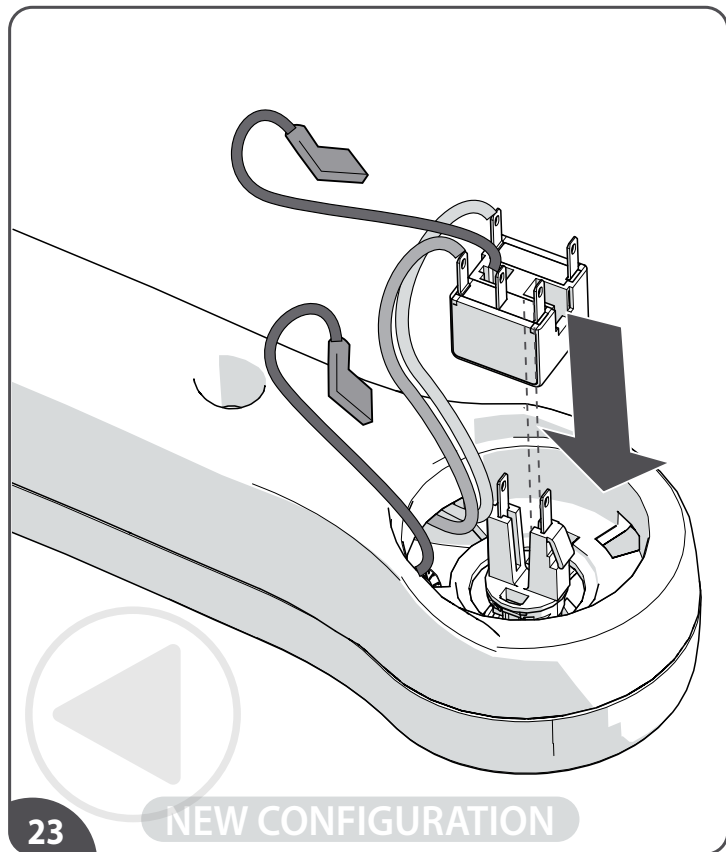CONTINUE WITH A NORMAL INSTALLATION

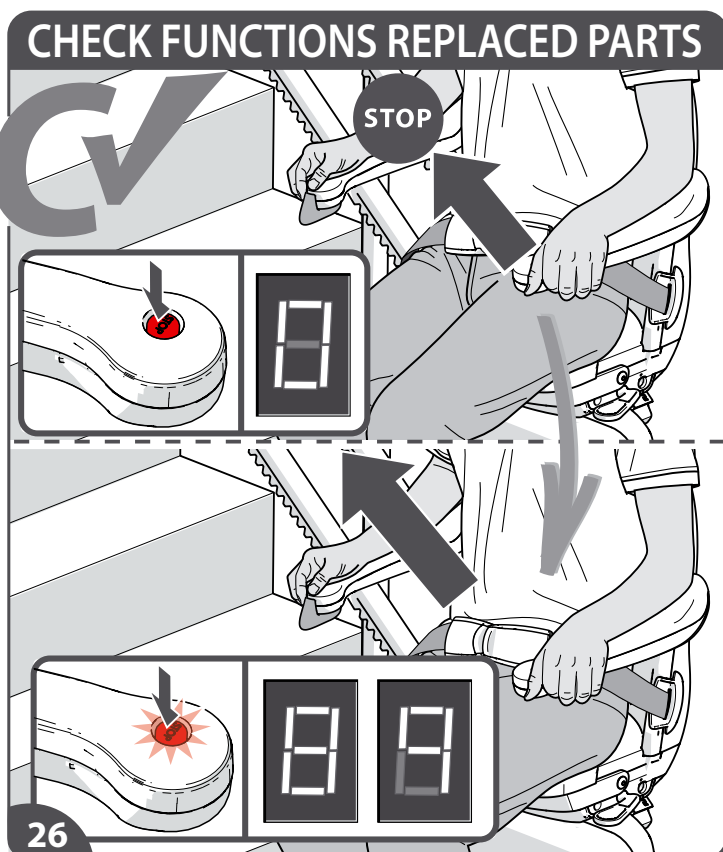- [23]Connect pole block to E-Stop.
 [24]Connect the blue wire to the pole block{A} and connect the earth wire to the pole block{B}.
 [25]Place plastic bung, to cover pole block. Ensure bung is in correct orientation.
Continue with a normal lift installation
 [26]After installation check the function of the replaced parts.

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



AA33961 2.4 GHz

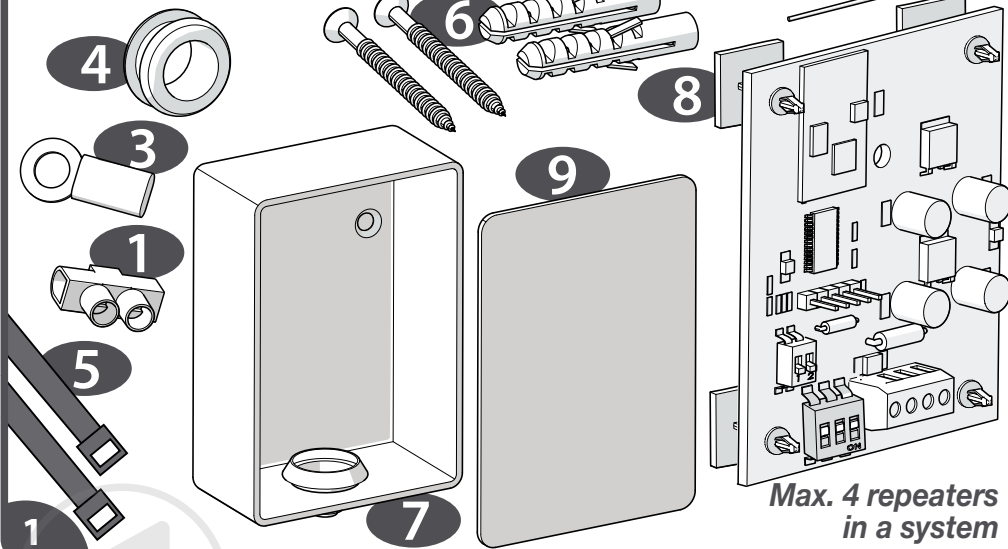TOOLS REQUIRED

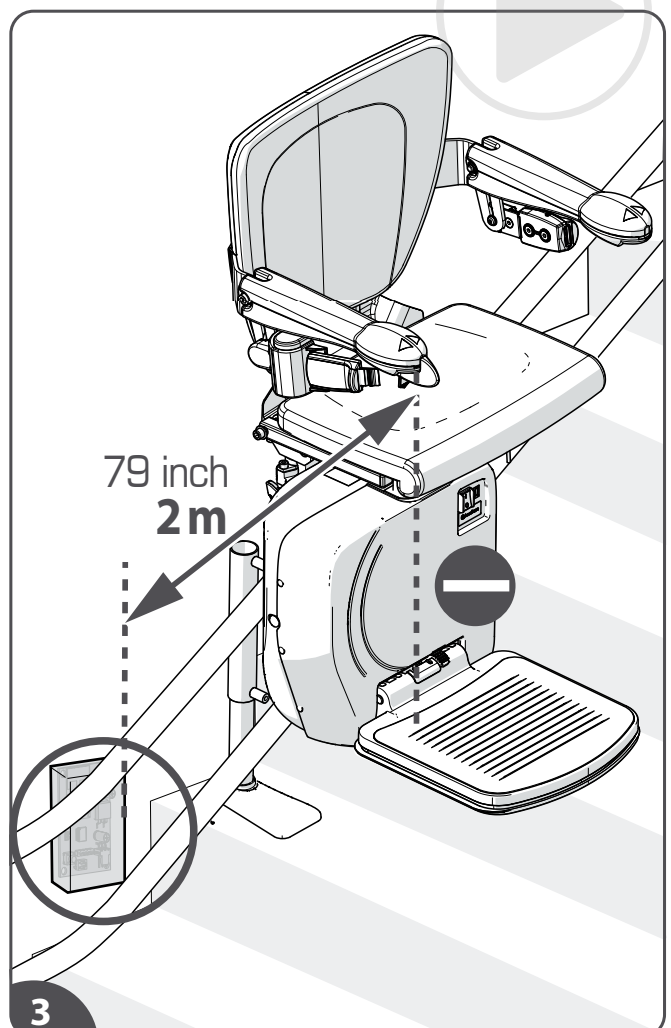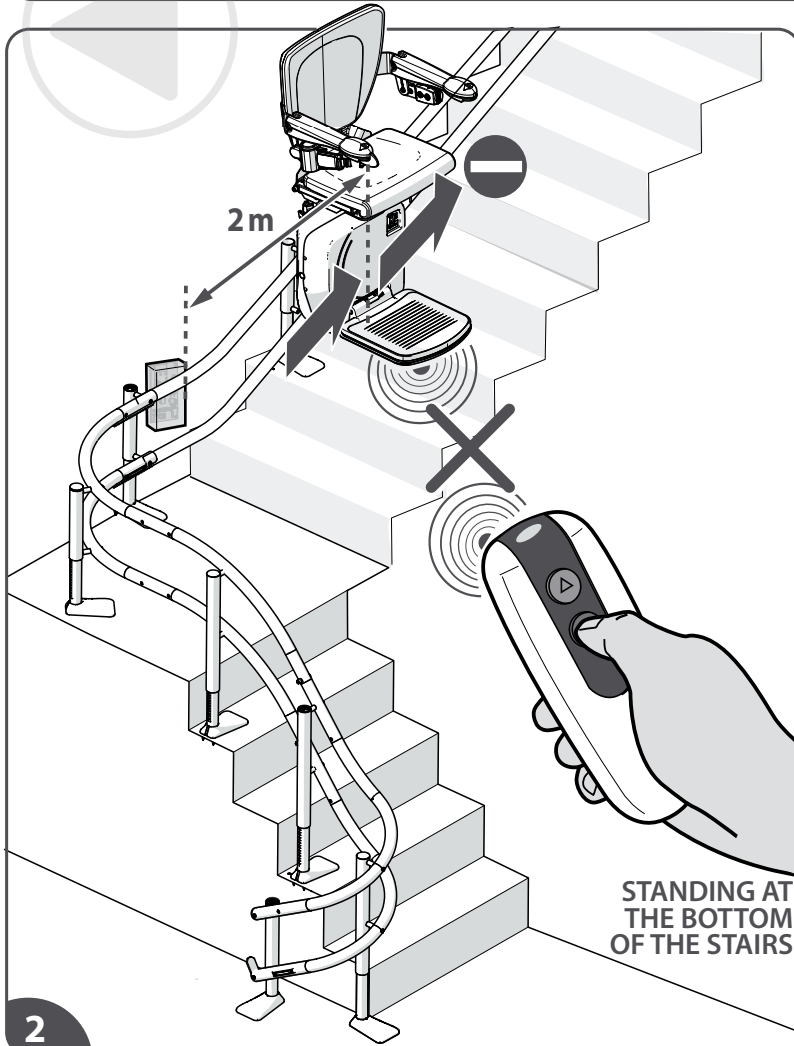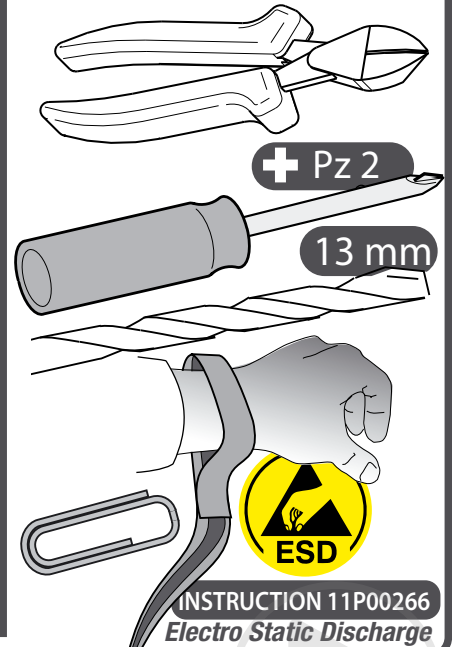- [1] Check that all the components of the concerning kit are present. Max. 4 repeaters in a system.
- [2] While standing at the bottom of the stairs, use the remote handset to drive the lift towards the other end of the track. When the lift loses the radio signal it will come to a stop.
- [3] The RF repeater should be fitted approximately **2 m** back from where the lift has stopped. The repeater does not require line of sight so it can be placed discreetly, for example behind the track. Where possible, the repeater should be mounted at least 150 mm [6 inches] away from any metalwork.

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



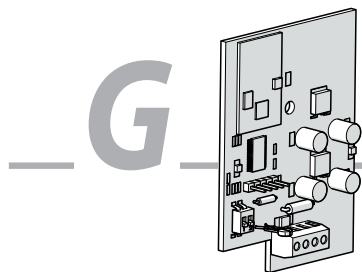
G

RF REPEATER

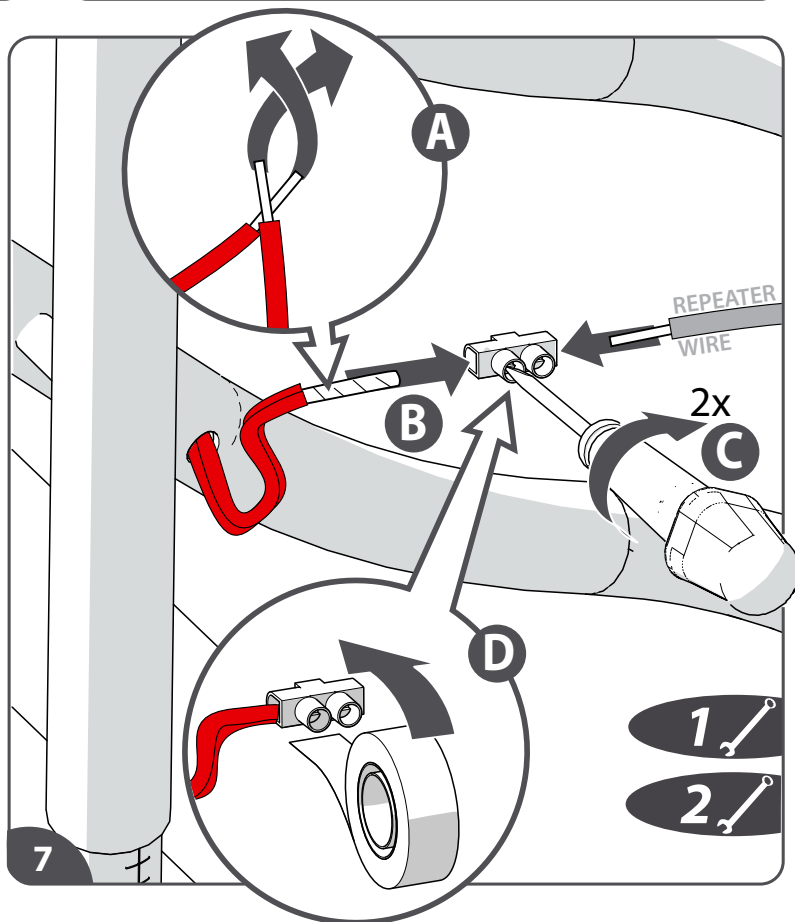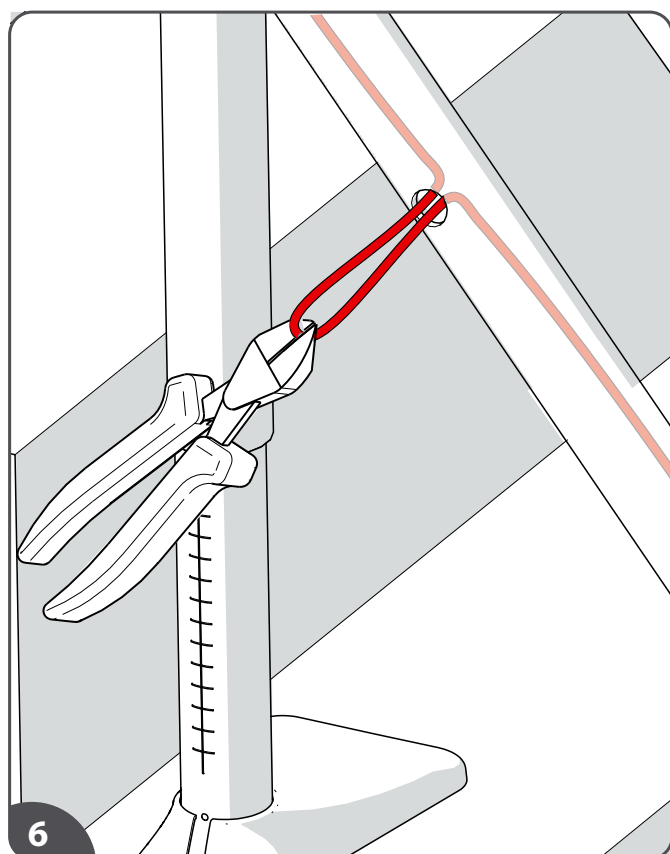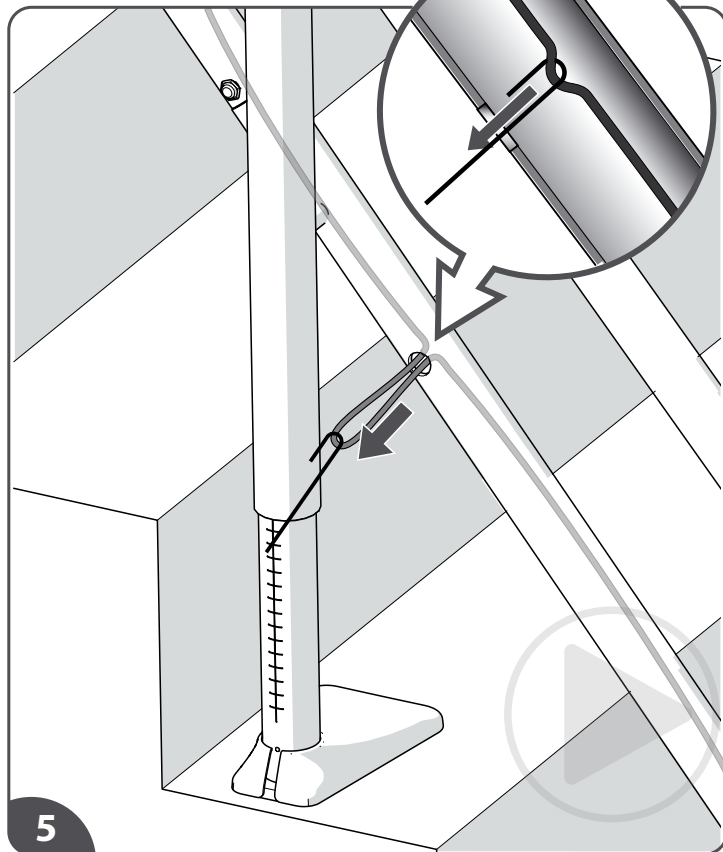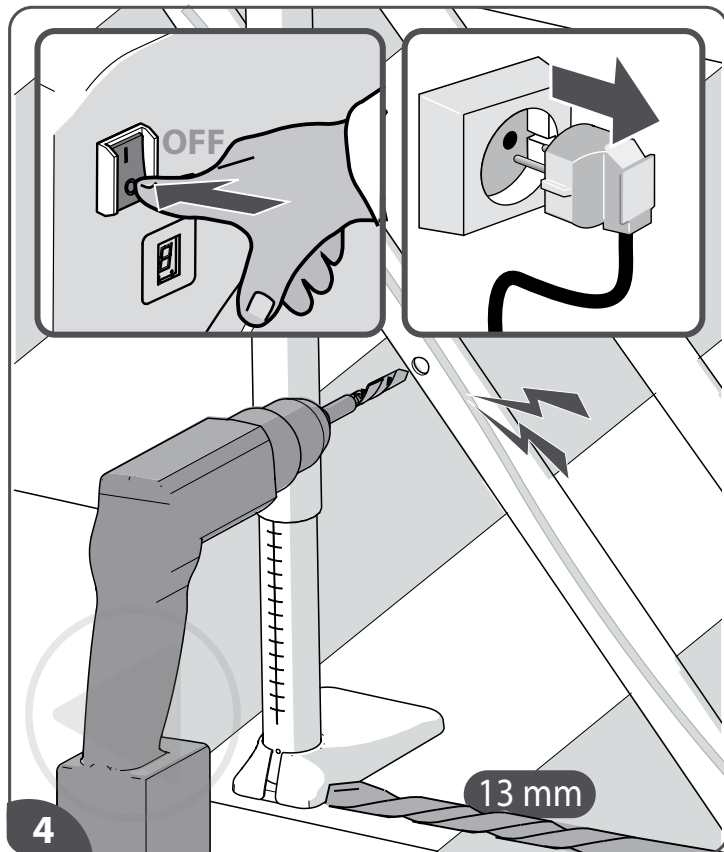- [4] Turn off/disconnect the power supply and switch off the powerpack. Drill a 13 mm hole in the bottom rail. Be careful not to damage the wire inside.
- [5] Use a bent piece of iron wire to pull the cable from inside the rail through the hole.
- [6] Cut through the cable. Remove the insulation of the two wires (15 mm).
- [7] Twist the two ends bare wires together {A}. Fasten the wires in a connector block {B+C}. Secure and insulate the connector block {D}.

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




IT

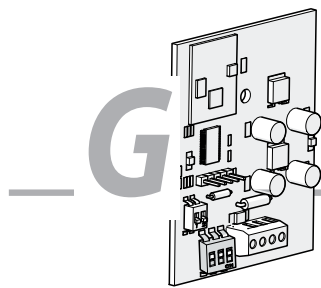RF REPEATER

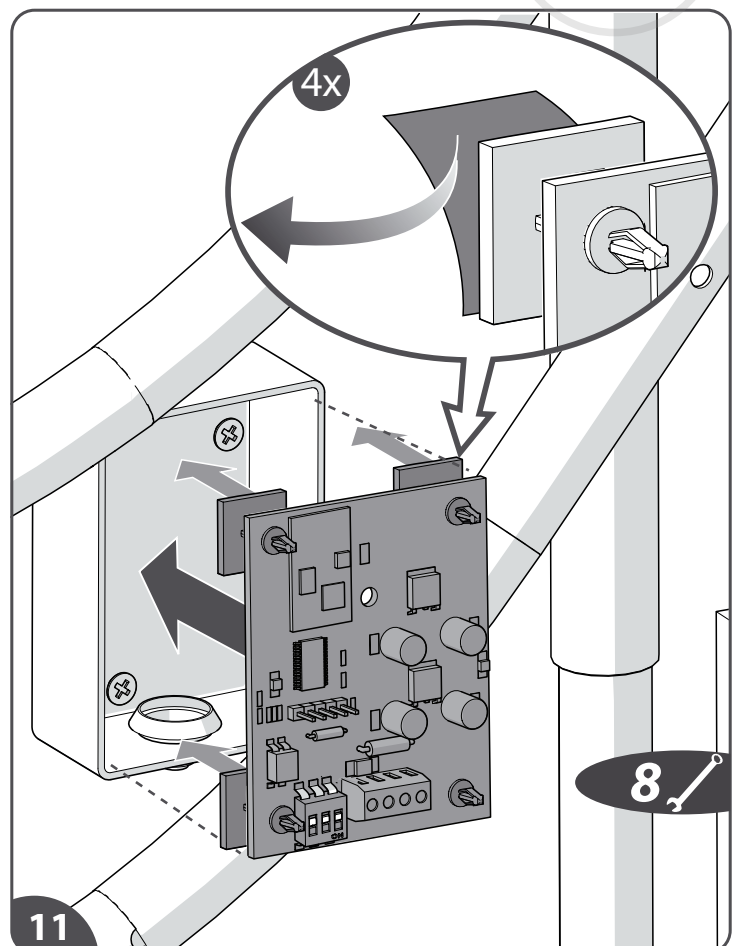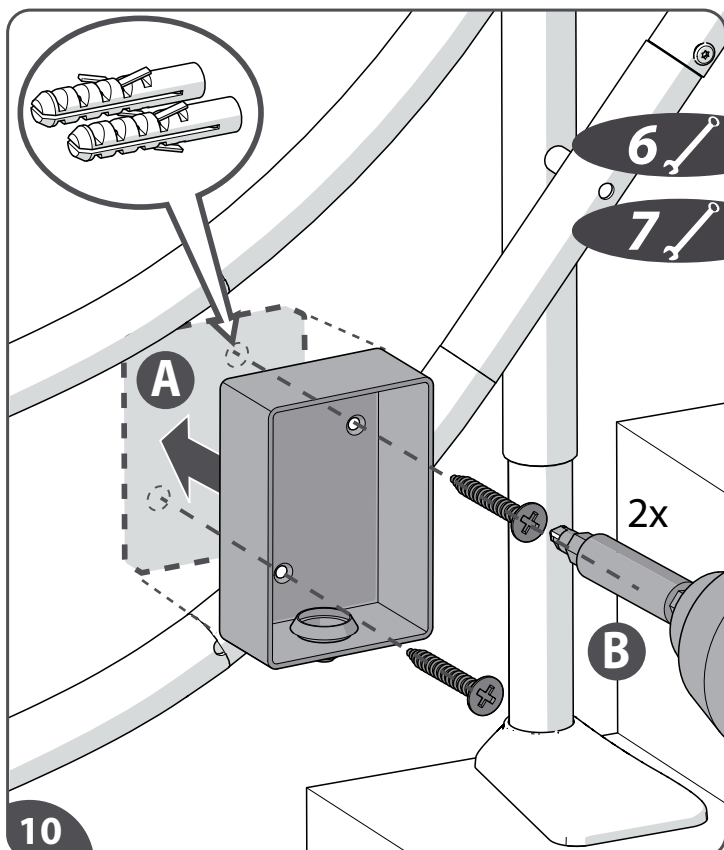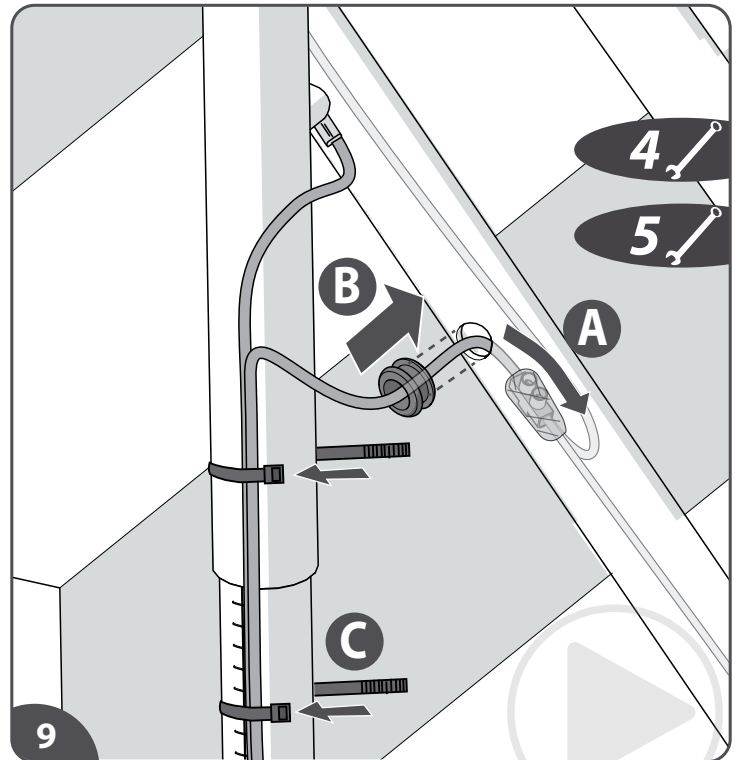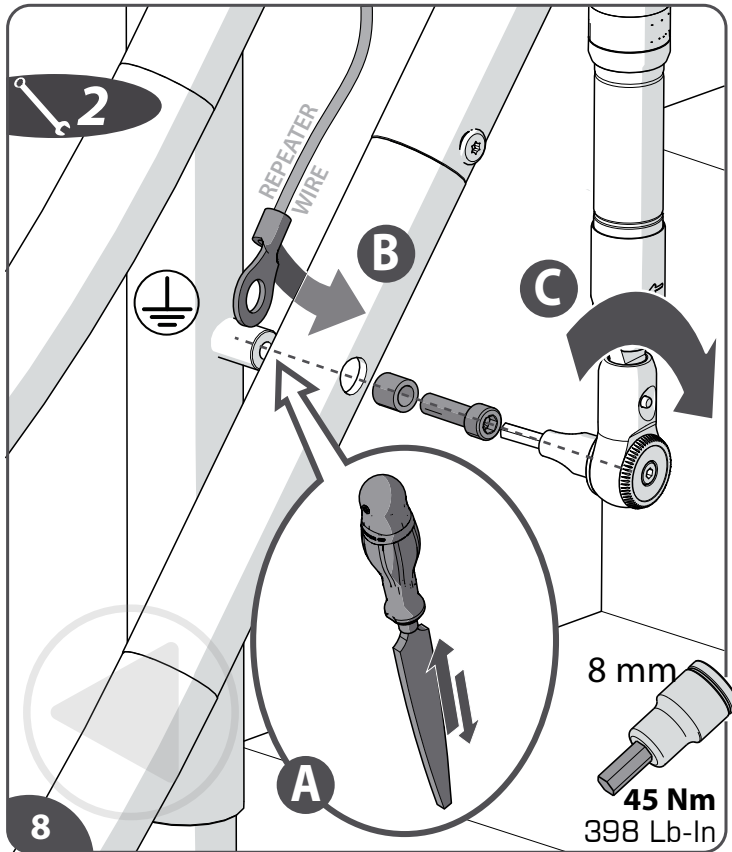- [8] Remove the paint from the rail{A}. Position the earth wire between lug and rail{B}. Place stanchion bolt and spacer. Tighten the stanchion-rail bolt [45 Nm | 398 Lb-In]{C}.
- [9] Feed all the wires into the hole of the rail{A}. Insert the grommet into the hole{B}. Use cable ties to fasten the excessive wires to the stanchion{C}.
- [10] Position the holder at the desired location and fix it to the wall.
- [11] Press the repeater PCB in the holder.

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




IT

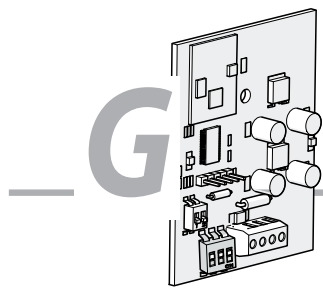RF REPEATER

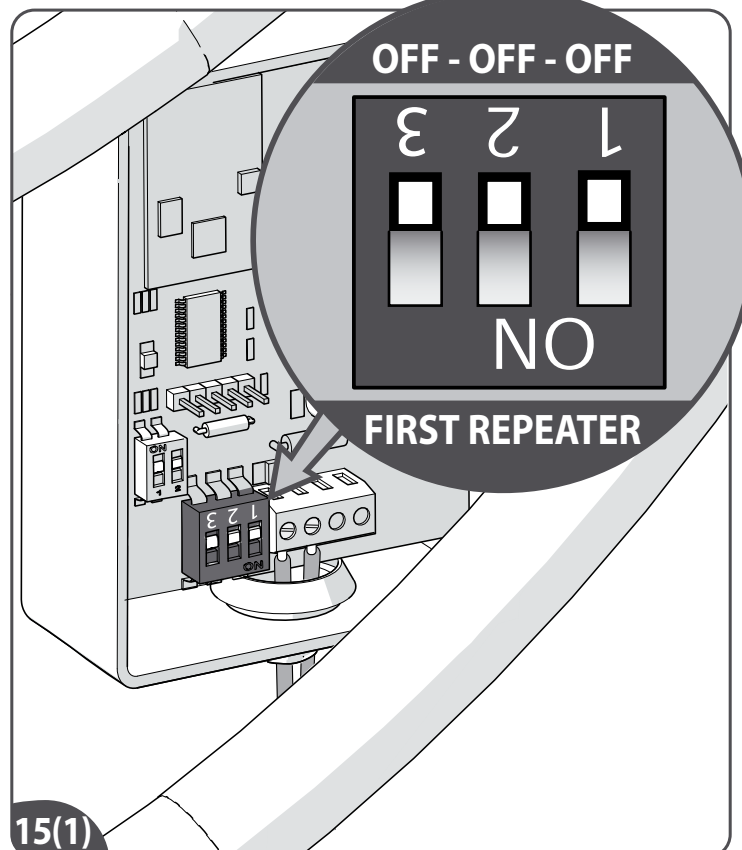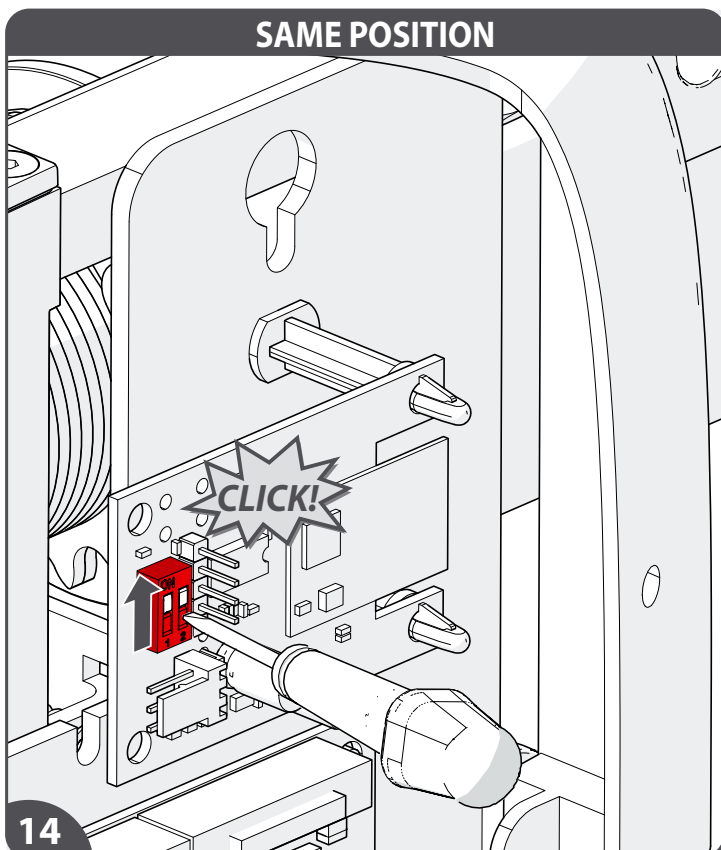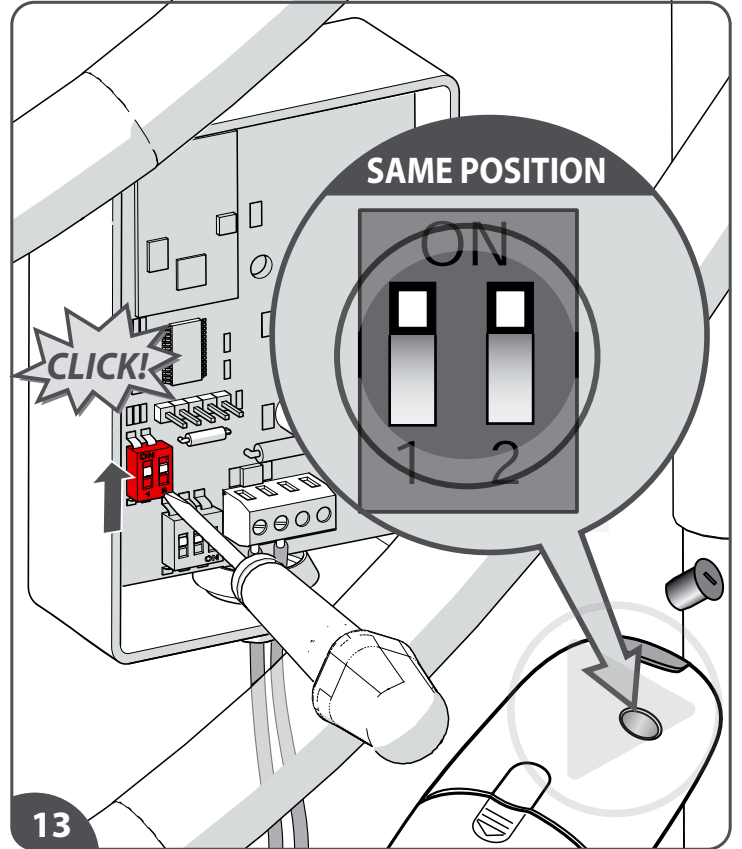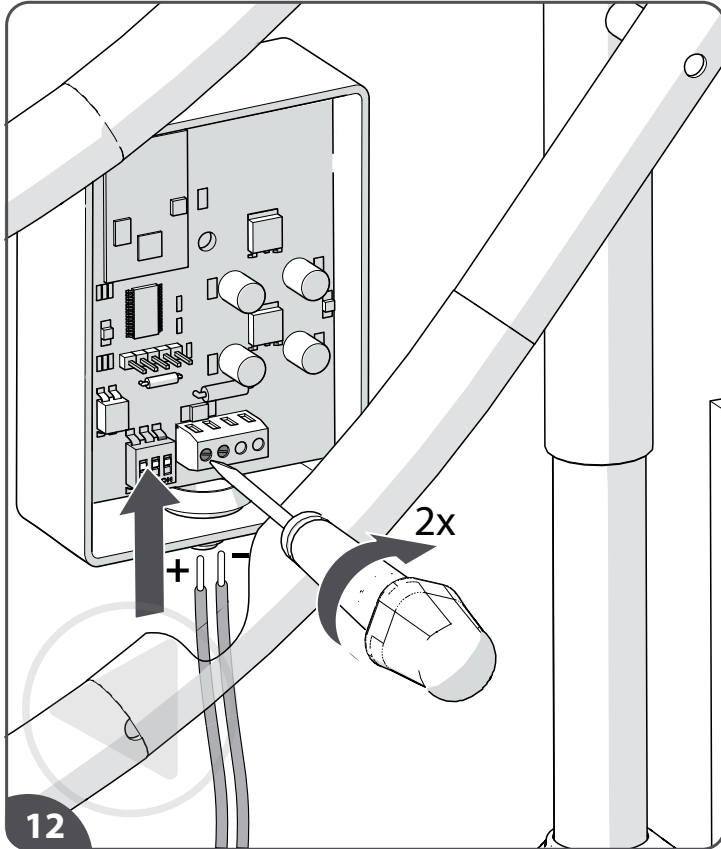- [12]Connect the wires with the repeater PCB. The **+** and **-** connections are labelled on the RF Repeater.
 [13]Set the dipswitch settings on the PCB: each RF remote should have the same dipswitch settings as the RF Repeater.
 [14]Also the RF receiver should have the same dipswitch settings as the RF Repeater.
 [15_1]Settings 3 way DIP switch on the repeater: 3 x OFF.

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



IT

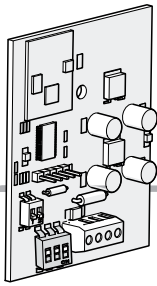RF REPEATER

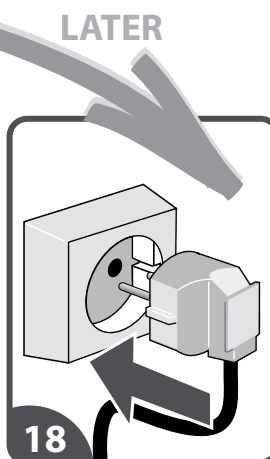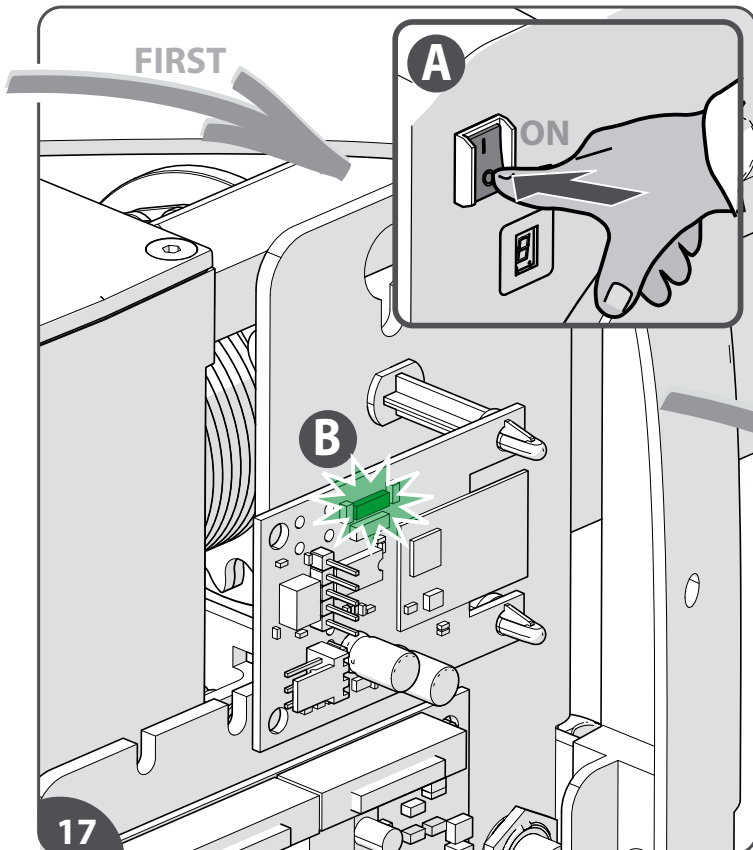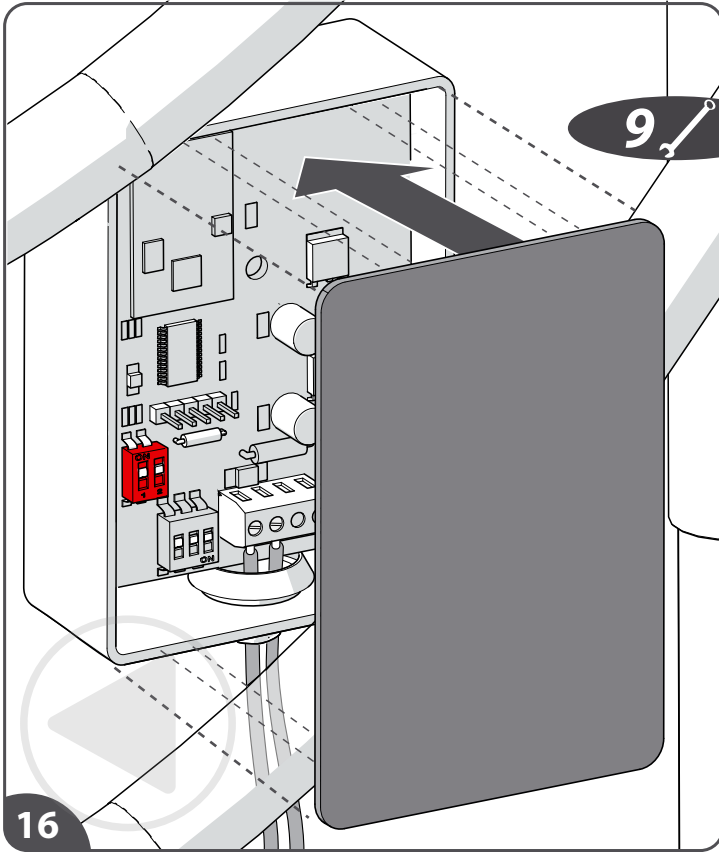MORE REPEATER CAN BE PLACED IF SIGNAL WILL BE LOST FURTHER ON THE TRACK ► FOLLOW THE SAME PROCEDURE EXCEPT STEP 15 (see next page for the step that takes its place)

[16]Close the holder with the screen in between{A}. Fasten the cover of the holder{B}.

[17+18]Turn on the powerpack{A} (and wait until receiver LED is solid green{B}) before restoring power to the power supply.

In the event that the stair configuration results in the signal being lost, further repeaters can be fitted by returning to the start of the rail and repeating the steps above, except STEP 15 (see next page for the step that takes its place)

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




IT

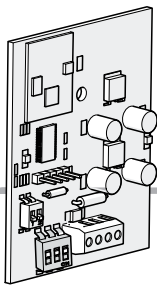RF REPEATER

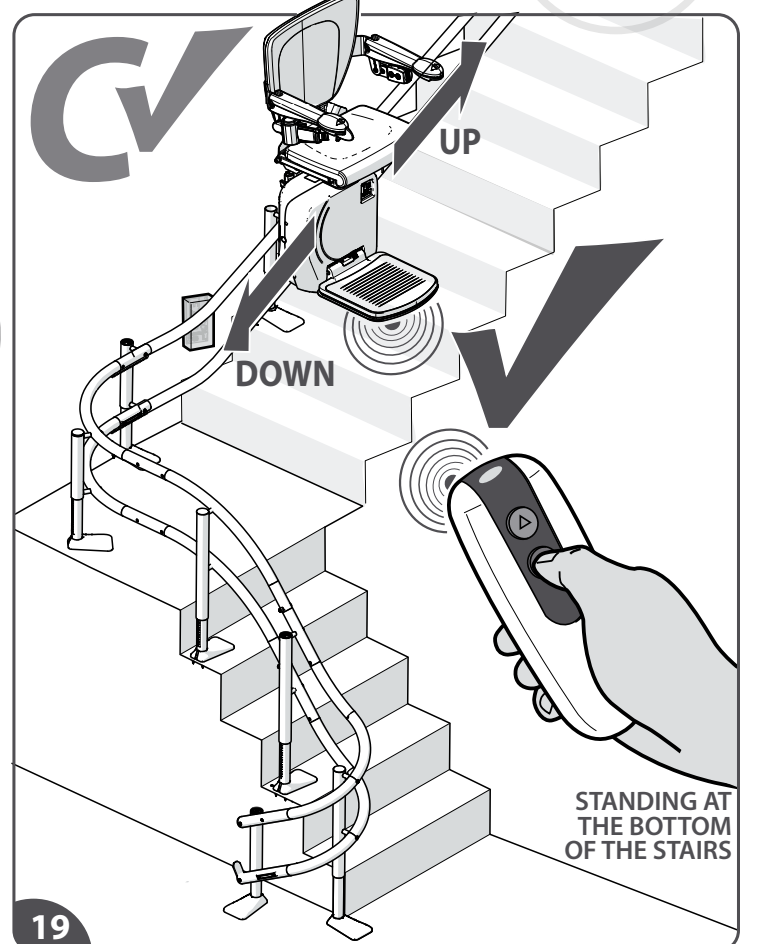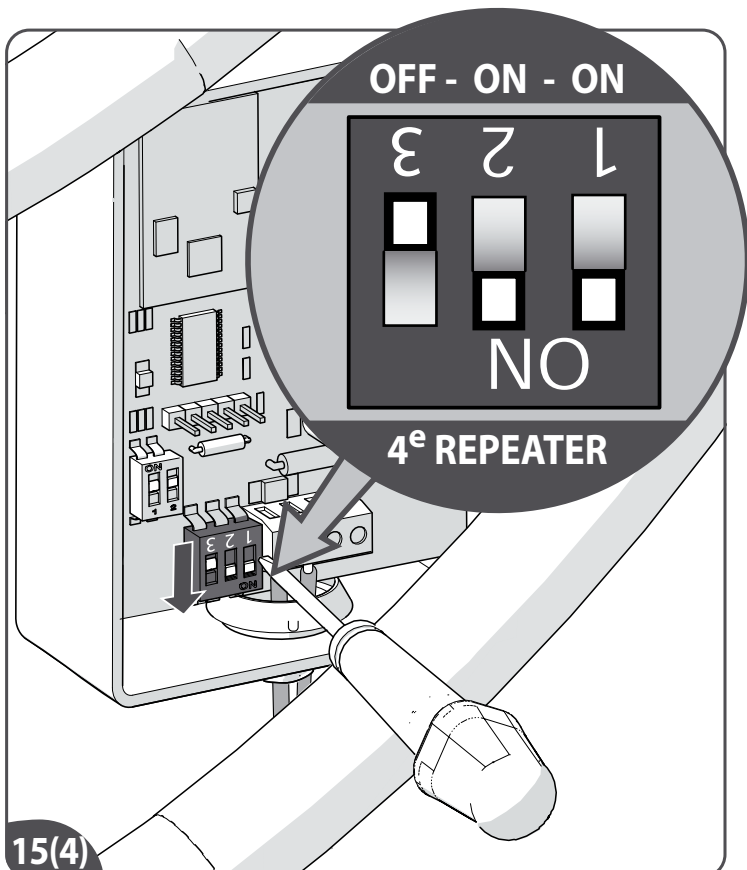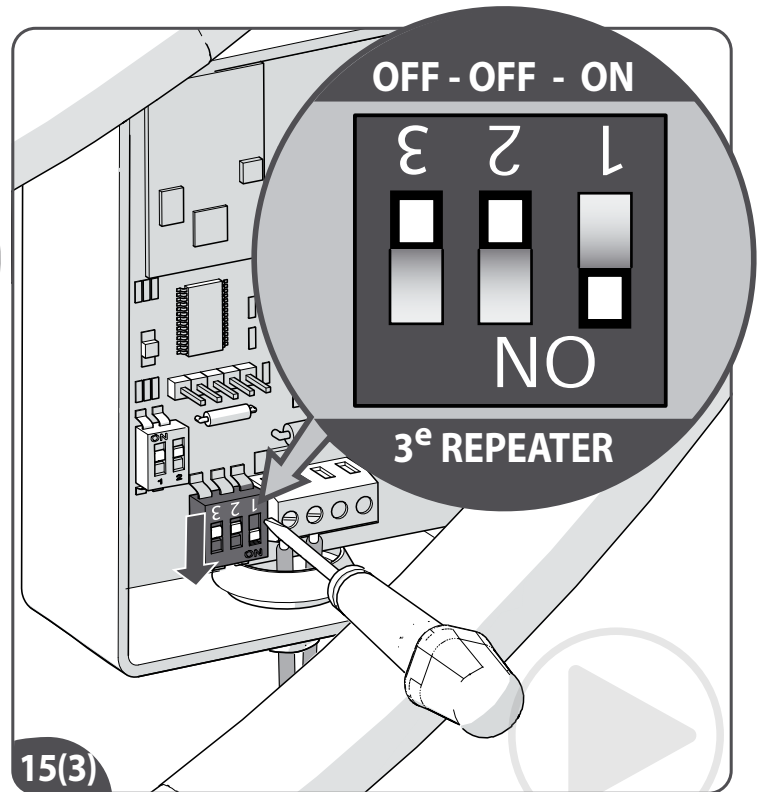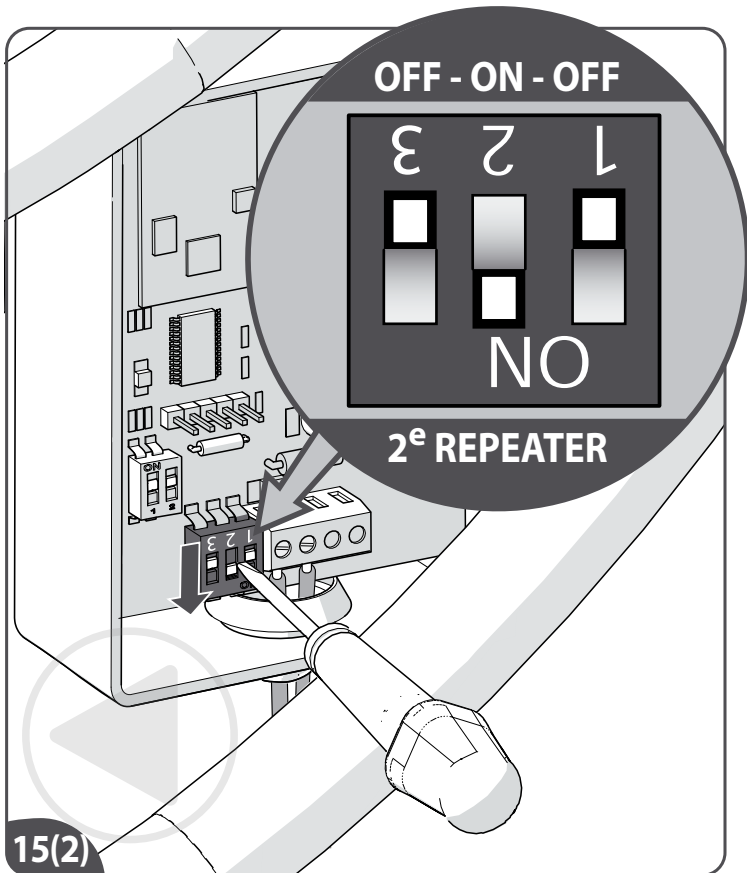[15_2]In case 2 repeaters are needed > Settings 3 way DIP switch on the repeater: OFF - ON - OFF.
 [15_3]In case 3 repeaters are needed > Settings 3 way DIP switch on the repeater: ON - OFF - OFF.
 [15_4]In case 4 repeaters are needed > Settings 3 way DIP switch on the repeater: ON - ON - OFF.
 [19]When all repeaters have been installed, drive the lift from bottom landing to top landing and from top landing to bottom landing to make sure a reliable signal is present throughout the journey of the lift.

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


IT

—H— **NO CHARGE BRACKETS** — FIT THE TOP RAIL BRACKET

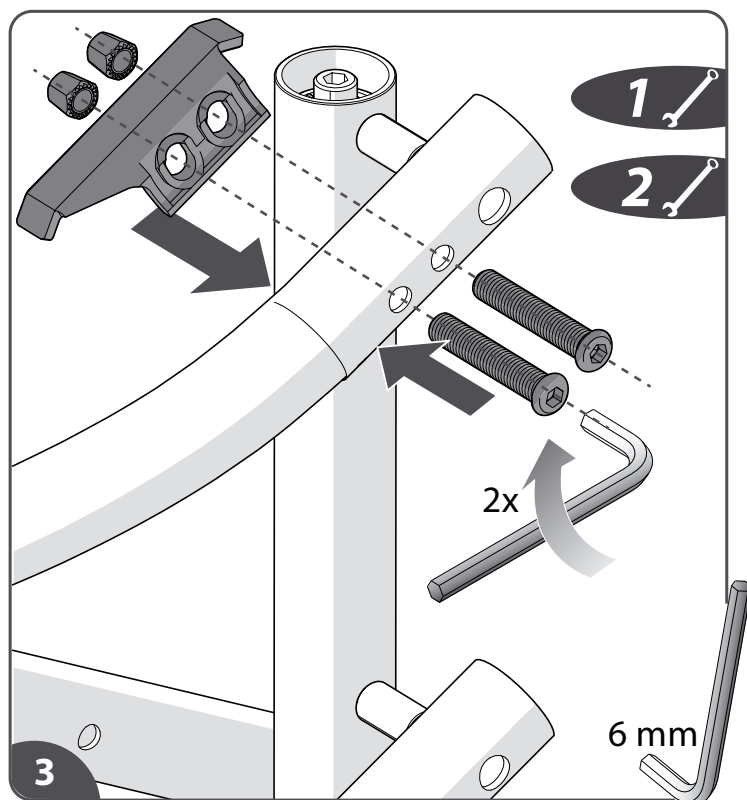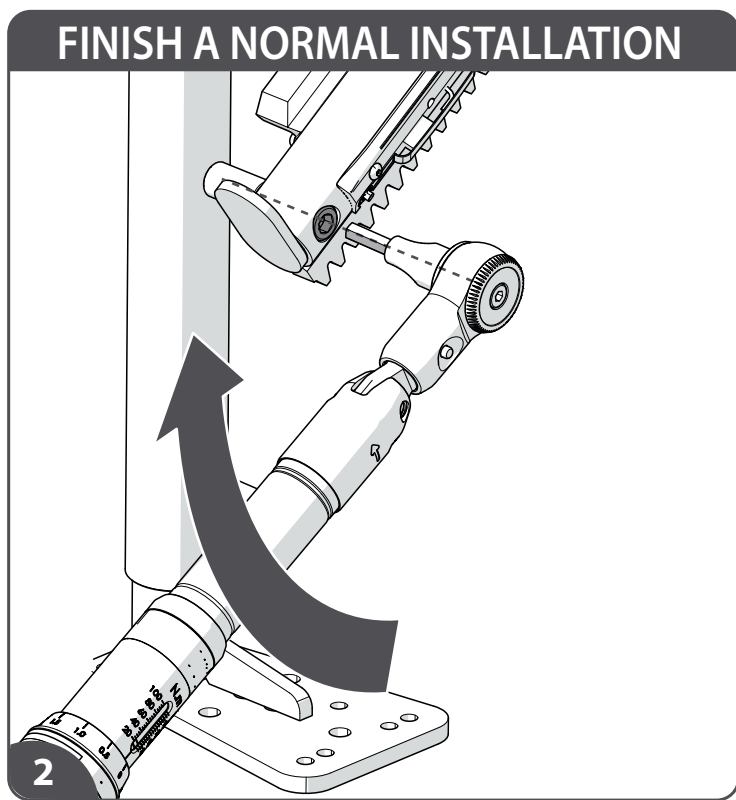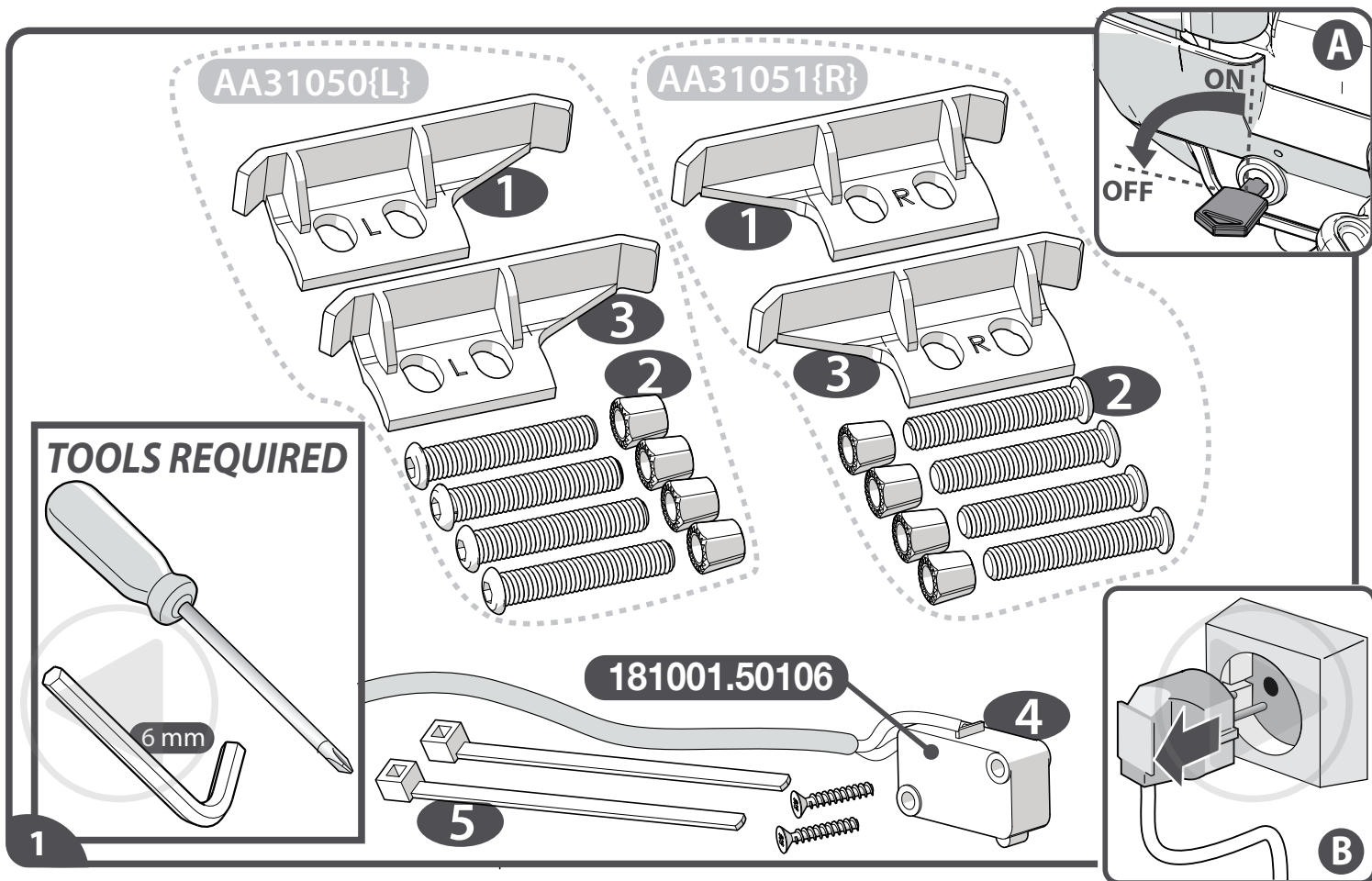Only a technician certified for this lift should install this lift (for more detailed information about the required competence see attachment of the installation manual). **General required knowledge can be found in the 4000 Technical Manual.**

The latest instructions for the lift and instructions for options are available for download at www.handicarepartners.com.

[1] Check that all the components of the concerning kits are present. Remove the key from the switch{A}.

Disconnect the charger{B}. Ensure good and sufficient tools.

[2] Do a normal 4000 rail installation, except the instructions mentioned below.

[3] Install the bracket on the upper rail.

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NL

- [4] Install the bracket on the lower rail.
- [5] Guide the cable through the trunnion.
- [6] Install the swivel switch.
- [7] Guide the cable into the power pack{A}. Fasten, without tension, the cable to the other trunnion cables{B}.

ES

- [4] Install the bracket on the lower rail.
- [5] Guide the cable through the trunnion.
- [6] Install the swivel switch.
- [7] Guide the cable into the power pack{A}. Fasten, without tension, the cable to the other trunnion cables{B}.

DE





- [4] Install the bracket on the lower rail.
- [5] Guide the cable through the trunnion.
- [6] Install the swivel switch.
- [7] Guide the cable into the power pack{A}. Fasten, without tension, the cable to the other trunnion cables{B}.

FR

- [4] Install the bracket on the lower rail.
- [5] Guide the cable through the trunnion.
- [6] Install the swivel switch.
- [7] Guide the cable into the power pack{A}. Fasten, without tension, the cable to the other trunnion cables{B}.

IT

H **NO CHARGE BRACKETS** CHECK THE BRACKETS

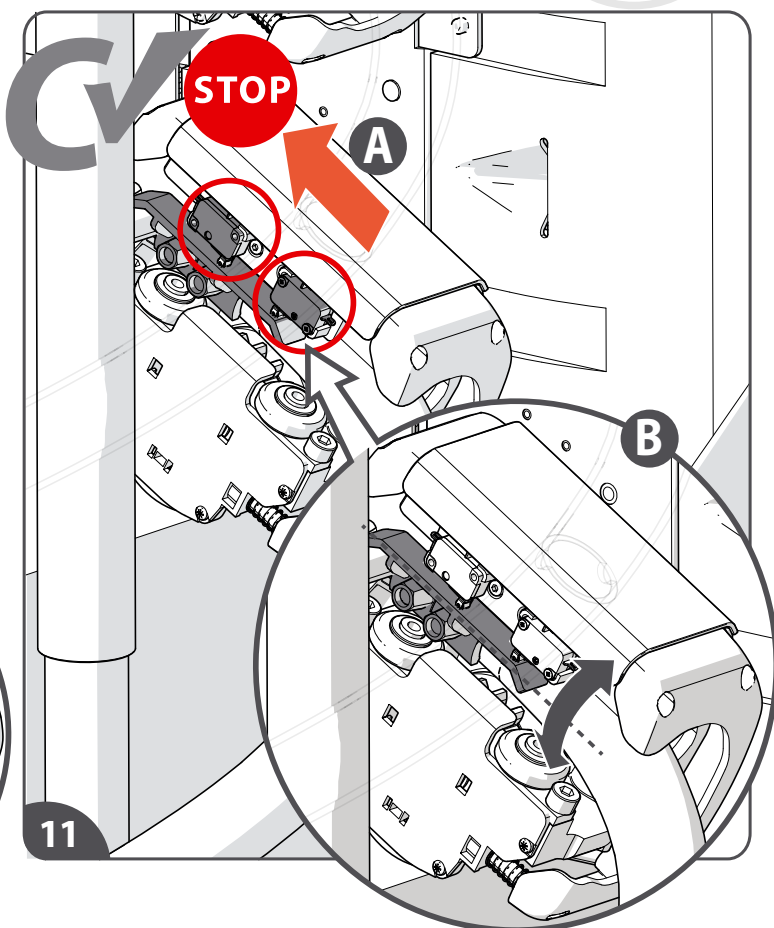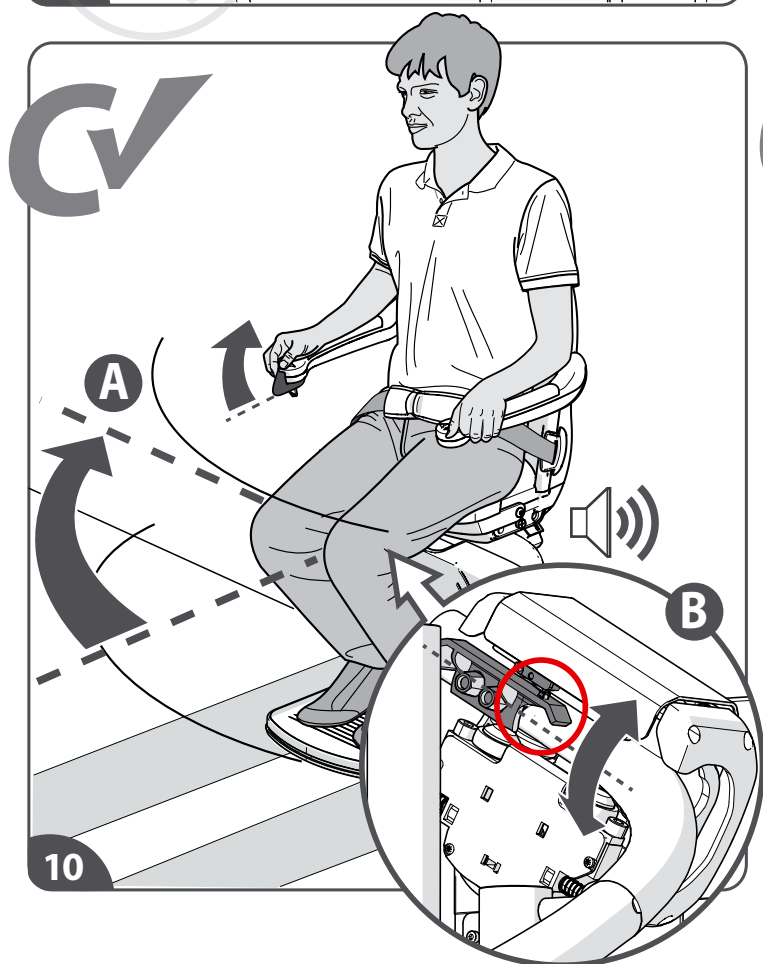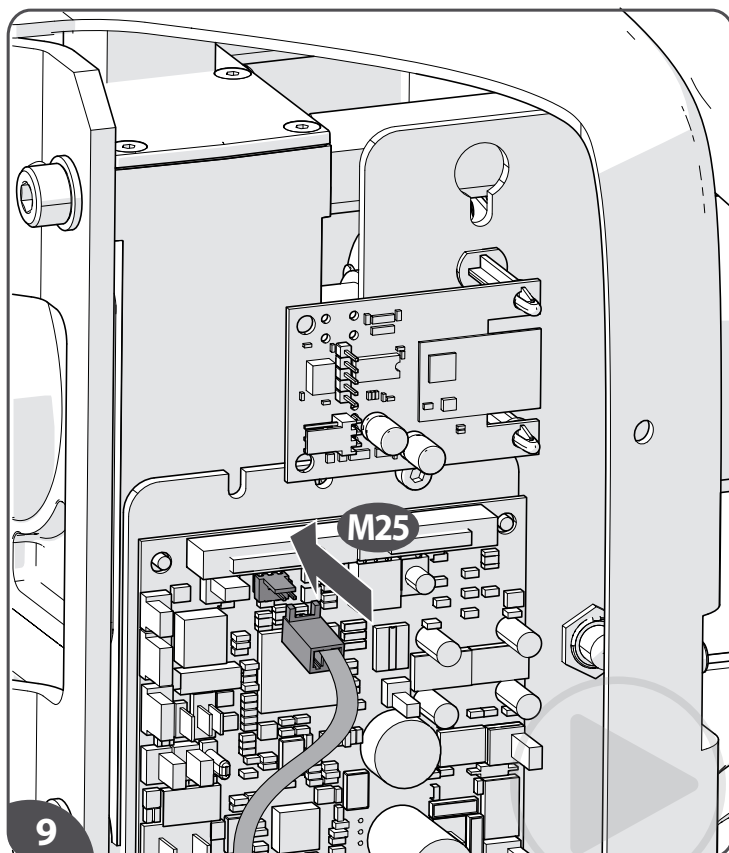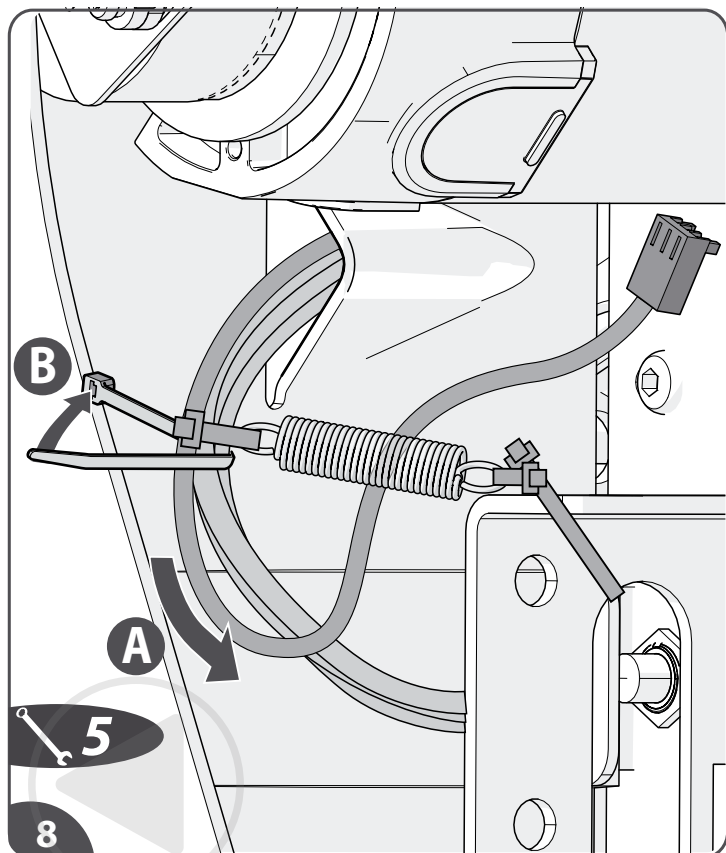- [8]Apply a cable tie to the indicated location.
 [9]Connect the swivel cable (**M25**) .
 [10]Check the powered swivel function. A lift (with an automatic swivel configuration) is able to swivel when:
 1.The stop switch hits the limit actuator at the top of the staircase.
 2.There is 33 volt on the charge point.
 If needed, adjust{B}.
 [11]Check the limit actuator (both switches activated){A}. If needed, adjust{B}.

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FUSE DETAILS*

FUSE REFERENCE	AMPERAGE RATING	VOLTAGE RATING	MANUFACTURER	TYPE	MANUFACTURER'S PART NUMBER
F1	25A	32V	MTA	ATO	0700360/20
F2	5A	32V	MTA	ATO	0700310/20

POWEROPTS

☒

 FOOTREST

☒

 SWIVEL LH

☒

 SWIVEL RH


*OR DIRECT EQUIVALENTS


192 INSTALLATION MANUAL


J SEAT HEIGHT IN RELATION TO MOUNTED STEM BOLT

SMART - STYLE - SIMPLICITY

Seat base height in mm in relation to the lowest mounted stem bolt.

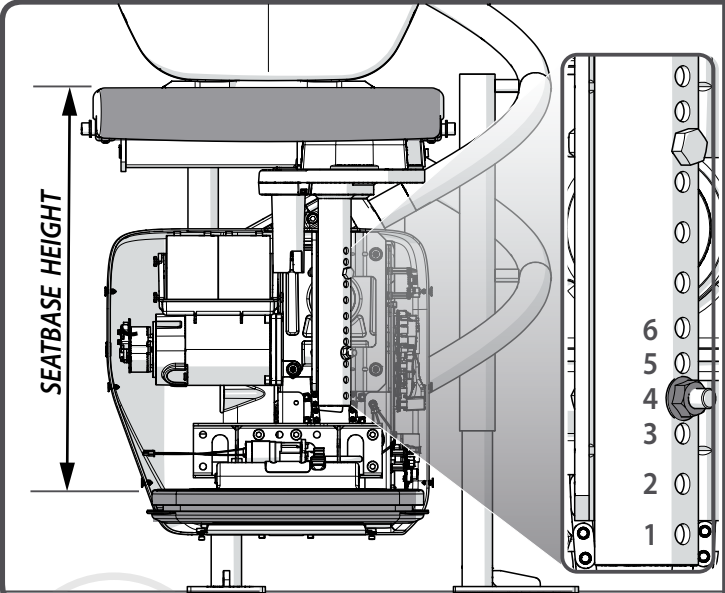	BOLT POS.	STYLE	SIMPLICITY	SMART
	6	528	528	
	5	514	514	525
	4	500	500	500
	3	486	486	475
	2	466	466	450
	1	446	446	425

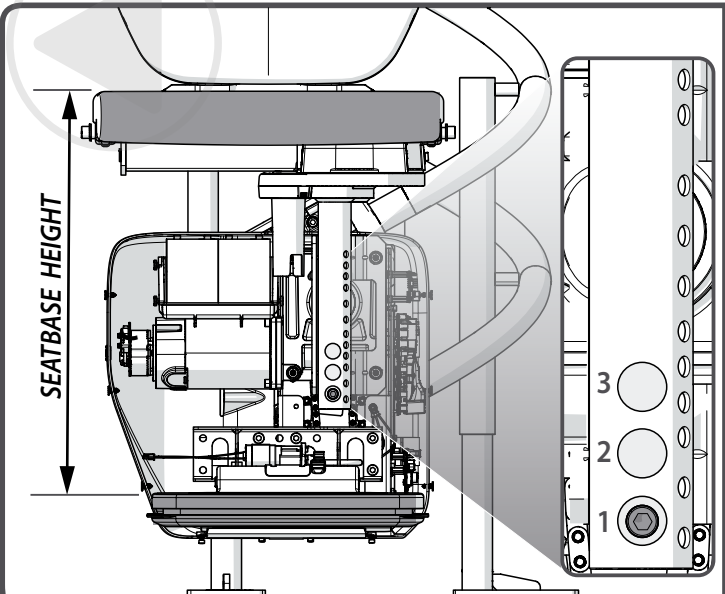
	BOLT POS.	STYLE HD	SIMPLICITY HD	SMART HD
	3	502	502	481
	2	474	474	453
	1	446	446	425

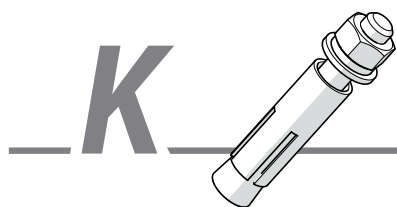ALTERNATIVE FIXING KITS


CONCRETE/STEEL/OUTDOOR

H4000 STEEL FIT KIT

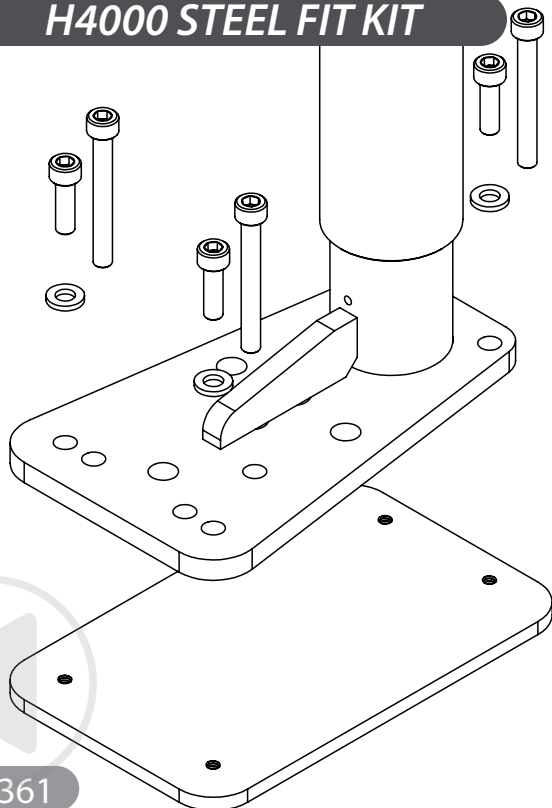
AA29361

H4000 CONCRETE FIT KIT

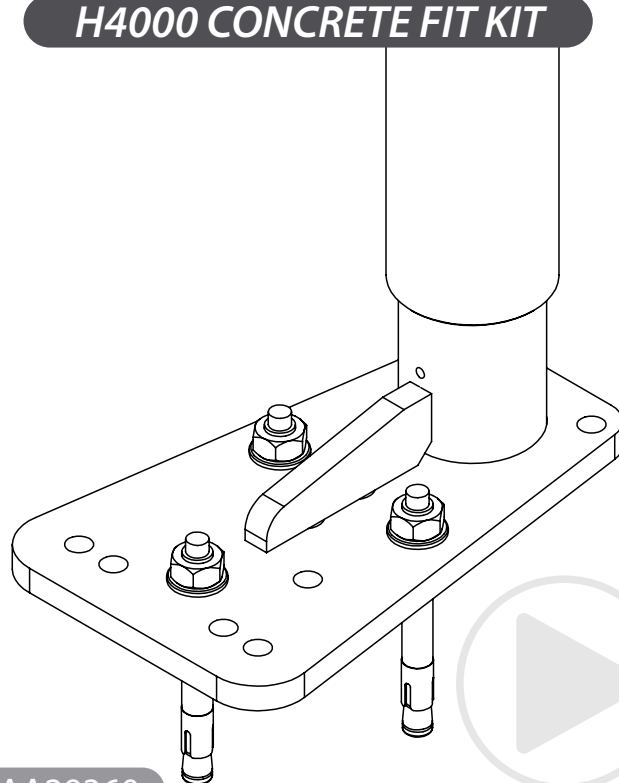
AA29360

H4000 SLOPED FLOOR FIXING OUTDOOR HEAVY DUTY KIT






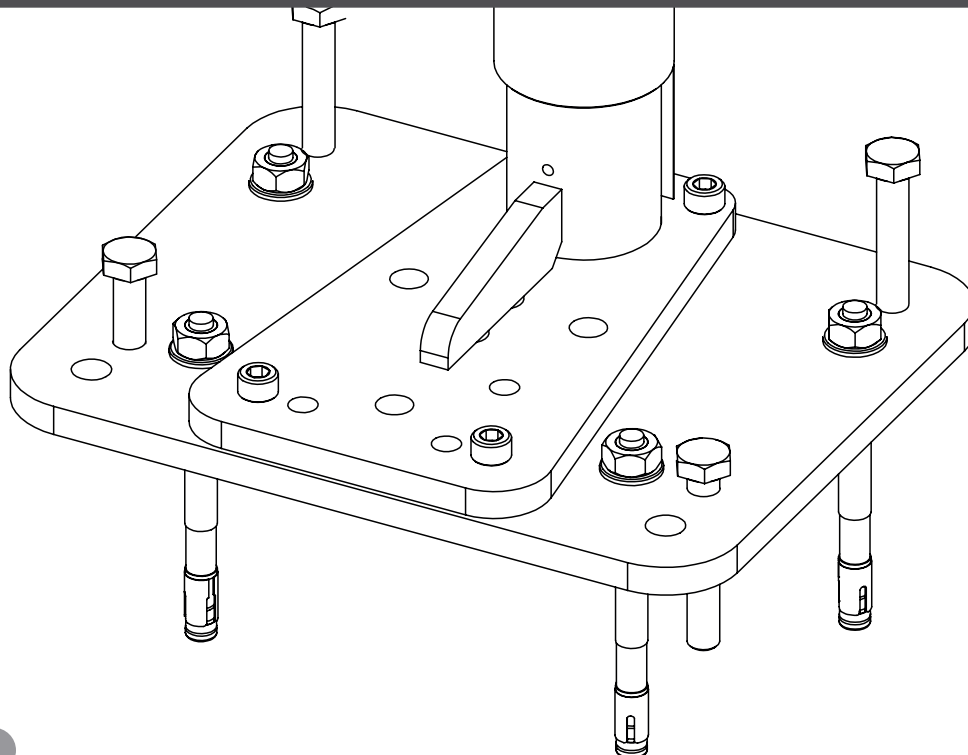
AA23321

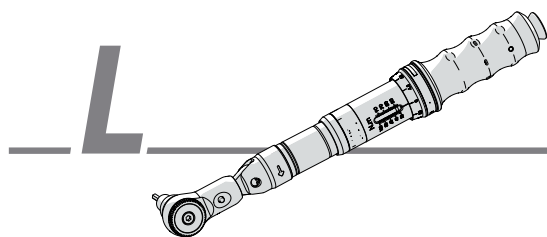TORQUE WRENCH TEST

CHECK FASTENERS

After the test drive, check the torque of the load-carrying fasteners with a torque wrench

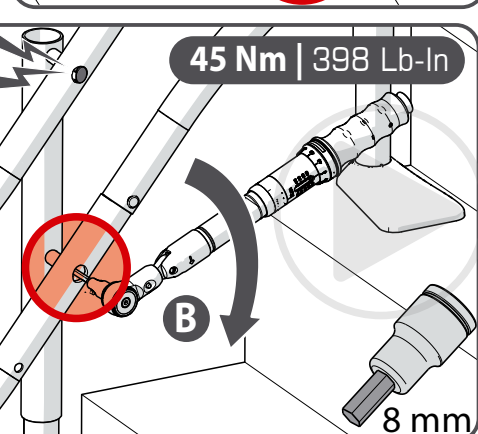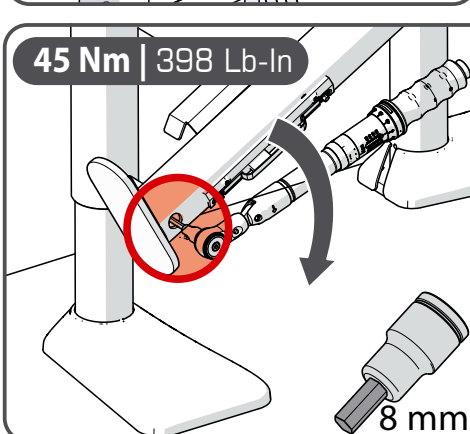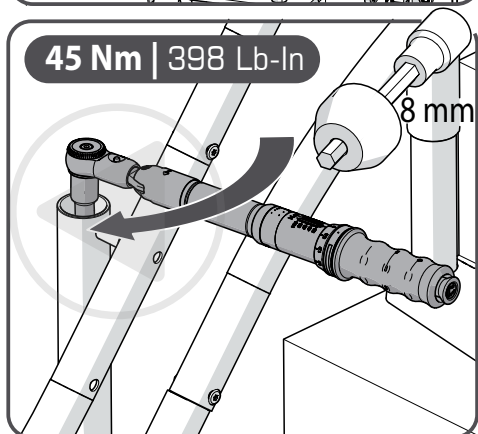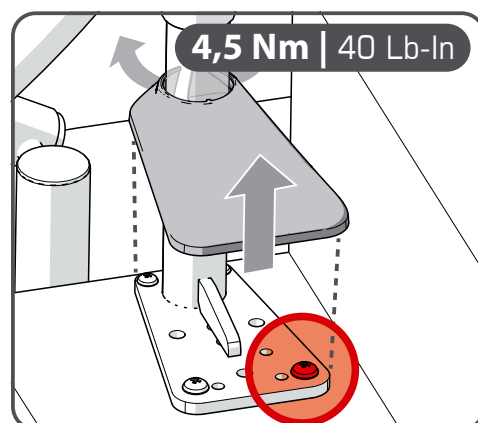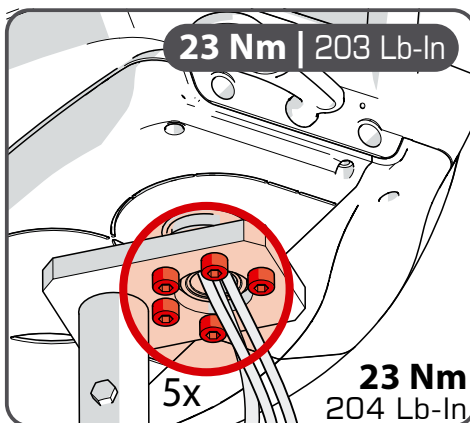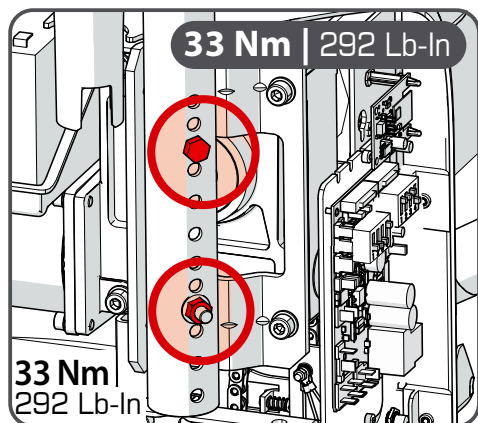M PROTECTION EQUIPMENT EXPLAINED

When performing any work on the stairlift, take the following precautions:

- Wear safety shoes type S3 or equivalent protection (S3 is an European classification; outside of Europe there are other classifications).
- Wear safety glasses for eye protection, i.e. to prevent injury because of parts suddenly striking into the direction of your face such as springs or when drilling in floors and walls.
- Use smooth fitting gloves to move and manage heavy / sharp parts of the unit and the rail.
- Wear a safety cap for installation with low bulkheads.
- Wear ear plugs when drilling and cutting. Exposure to sounds above 85 decibels can cause hearing loss.
- To prevent getting caught in moving parts of the stairlift: do not wear loose clothing or jewelry, such as necklace, bracelets and tie long hair.

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IT

**VERY IMPORTANT****EN**

This product must ONLY be installed by properly qualified and trained stairlift installation engineers, who have satisfactorily completed a Handicare-approved training course on how to install and/or commission the Handicare 4000-stairlift. Grave consequences regarding the safety of the stairlift in question and of its users may result from any violation of this requirement.

Intended Audience

This manual is intended for use by such qualified and trained stairlift installation engineers ONLY; it is NOT intended to show an unskilled person how to install and/or commission the stairlift, nor does it allow an unskilled person to install and/or commission the stairlift in an appropriate/safe manner.

This manual may NOT be used to replace the necessity for training.

The procedures described in this manual are 100% mandatory and MUST be followed.

Required Knowledge

All users of this manual MUST have the following knowledge:

Sufficient understanding of the key features, functions and capabilities of the Handicare 2000 stairlift;

Sufficient mechanical and electrical skills to be able to understand the function and purpose of each component of the stairlift, and sufficient mechanical and electrical skills as required for the correct installation and commissioning of the stairlift;

Sufficient Health and Safety awareness.

Disclaimer

It is the responsibility of the stairlift engineer to check that the manual is applicable to the type and version of the stairlift being installed; he/she MUST not fit a stairlift by reference to the manual if it is not the correct manual.

The stairlift engineer must read the manual to understand all the information to enable safe installation and/or commissioning of the stairlift.

In the absence of any negligence or breach of duty on the part of Handicare Accessibility Limited, no liability or responsibility will be accepted by it for any loss or damage arising from anyone's failure to follow and adhere to the guidance in the manual in connection with the installing and/or commissioning of the stairlift.

The latest and any previous versions of this manual can be obtained from www.handicarepartners.com

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Projectnr.....

INSTALLATION FORM**H4000****1 Description of installation**

Location:

Manufacturer:

Type:

Unit serial no:

Seat serial no:

- 1) Confirm that the installation dates have been completed on the EC Declaration of Conformity
- 2) Confirm that the remote controls are mounted on the wall and are functioning properly
- 3) Confirm that the safety belt is secure and functions correctly
- 4) Confirm that all warnings relating to the maximum working load on the lift have been affixed

2 Rail

- 1) Confirm that the rail joints are properly secured
- 2) Confirm that the rail supports are properly secured
- 3) Confirm that the folding rail is functioning correctly

3 Electrical drive unit tests

- 1) Confirm that the batteries are securely fitted in the carriage structure
- 2) Confirm that the following have been effectively indicated to the user:
 - a) location of the charge contacts.
 - b) "OFF CHARGE" alarm.
 - c) required action when stairlift is not charging

4 Clearance dimensions

- 1) Confirm the correct clearance dimensions as specified in the installation manual are maintained throughout the full travel of the stairlift
- 2) Confirm that there are no open bars/openings at the rear of the lift
- 3) Confirm that other users can safely use the staircase (minimum clearance of 400 mm between the folded seat and opposing obstacles). Install handrails if required.

5 Functional tests

- 1) Confirm that the final limit switches (upper and lower) operate correctly (position and operation)
- 2) Confirm that the control lever, key switch, emergency stop button, and the detection switches in the armrests operate correctly
- 3) Confirm that the safety detection devices operate correctly (footrest/unit/armrest)
- 4) Confirm that the swivel seat interlock and switch operate correctly
- 5) Confirm that the stopping distance of the lift is a maximum of 2 cm
- 6) Confirm that the restart time of the lift is at least 1 second
- 7) Confirm that the alarm device operates correctly when activated
- 8) Confirm that the lift has successfully passed the static (torque wrench test) and dynamic test ride in accordance with the instructions in the installation manual

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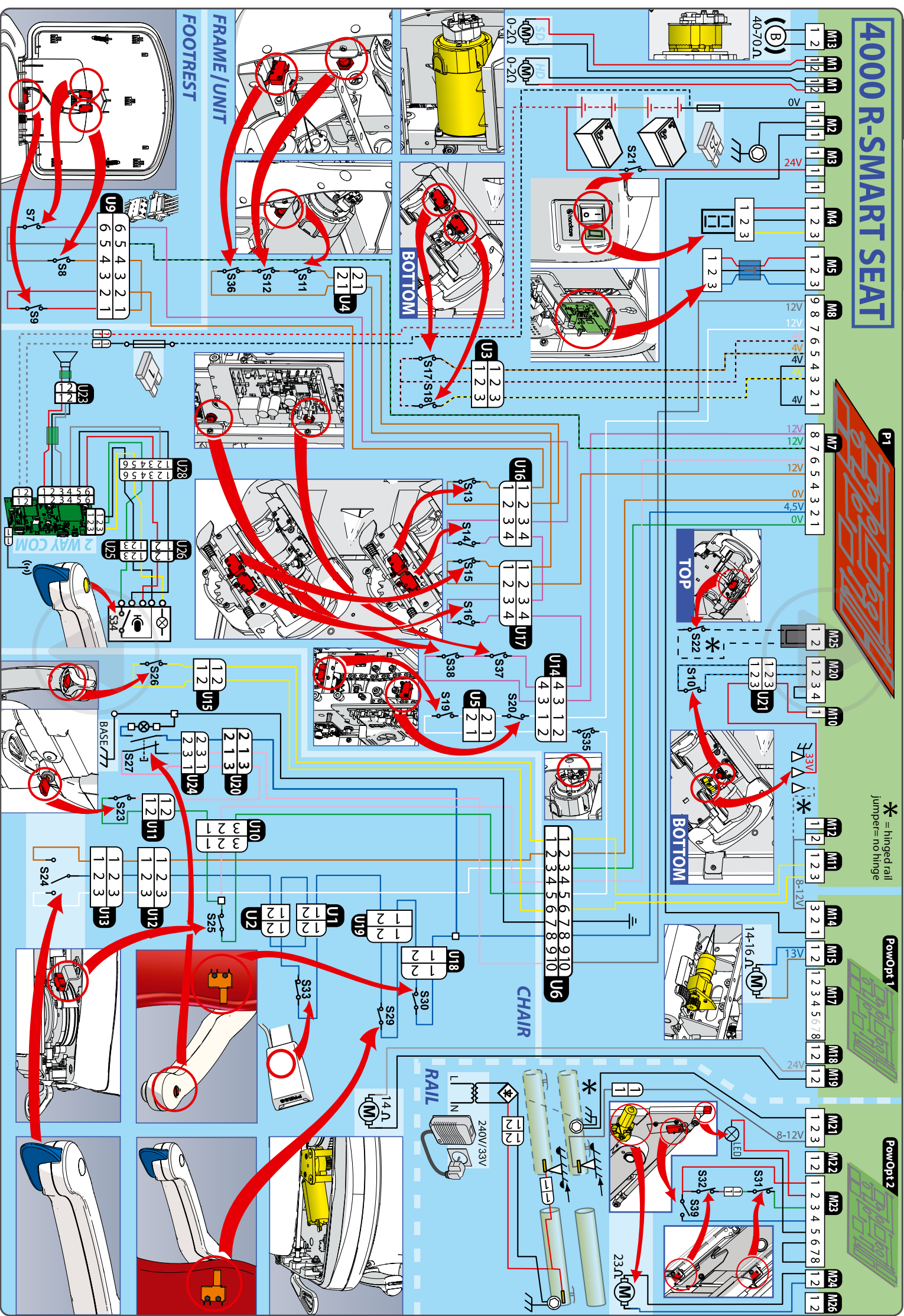
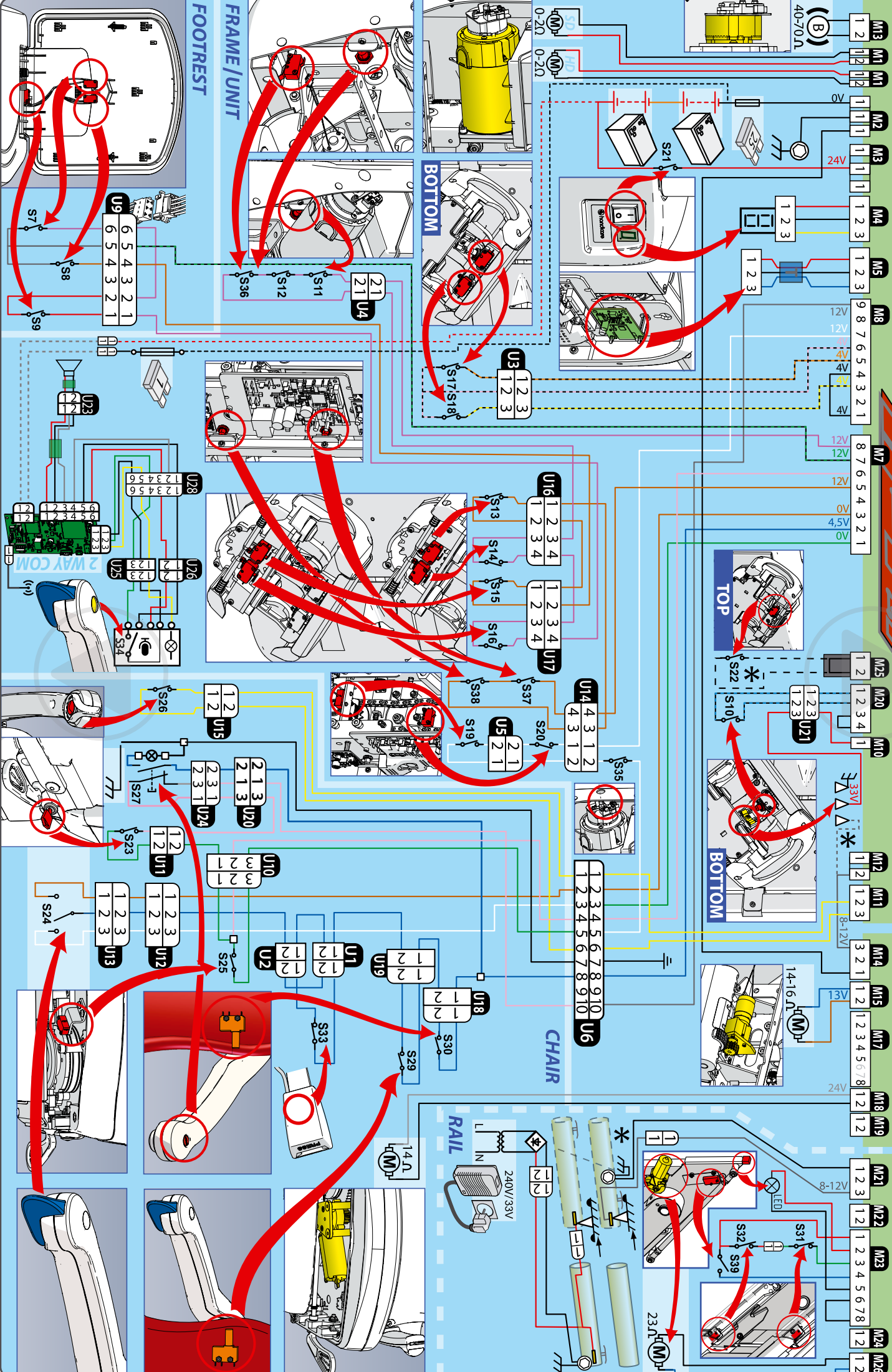
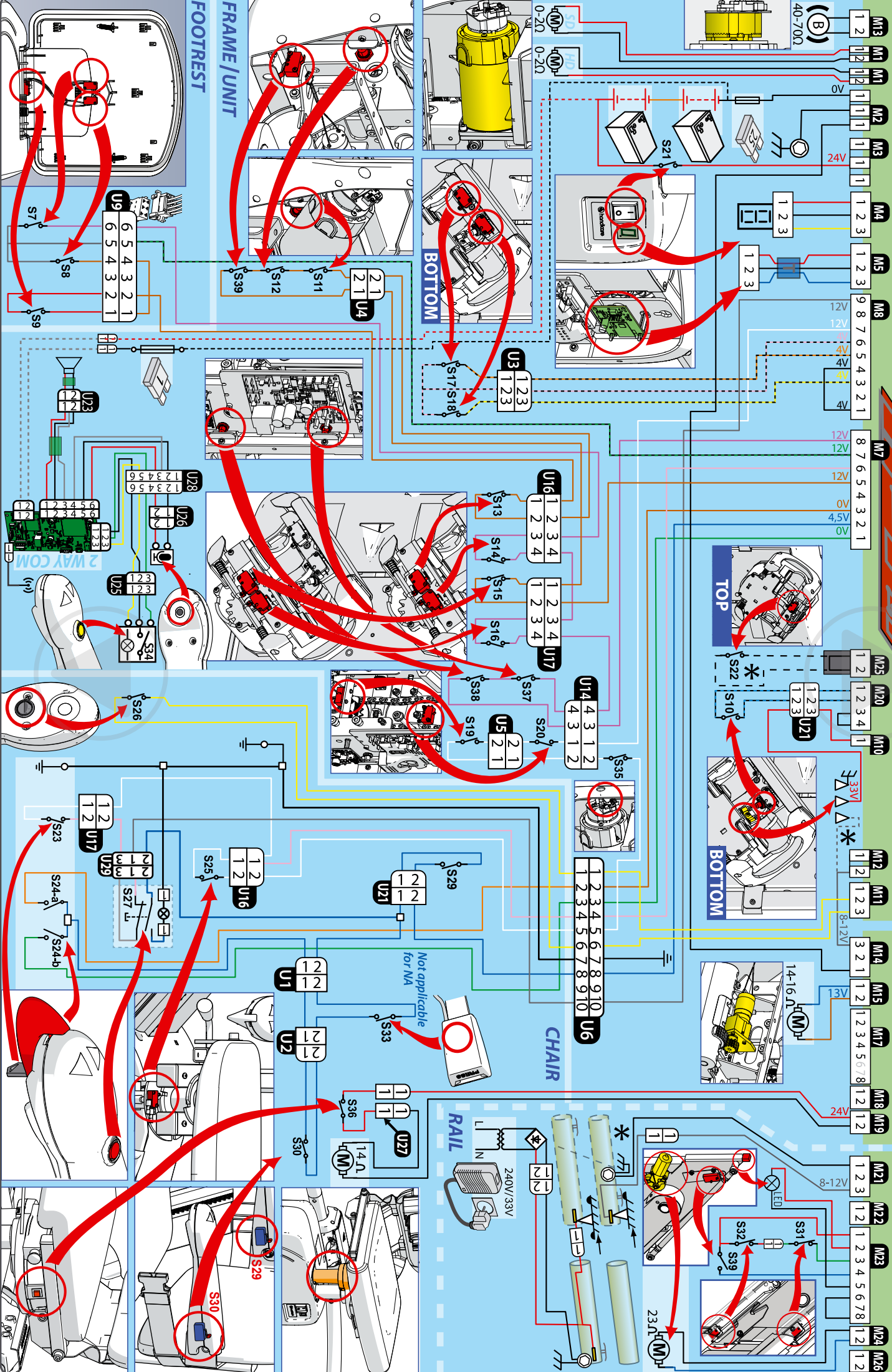
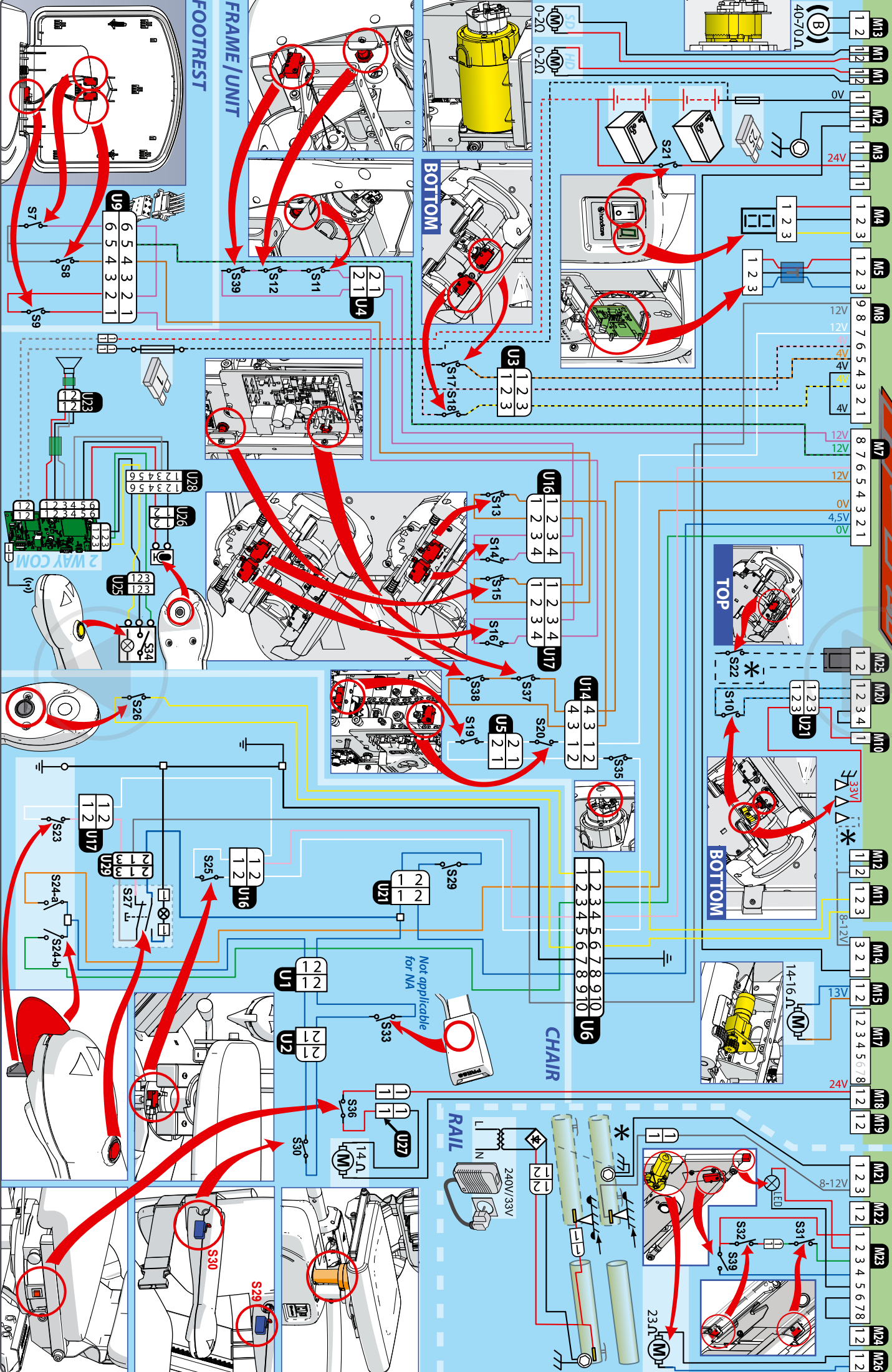
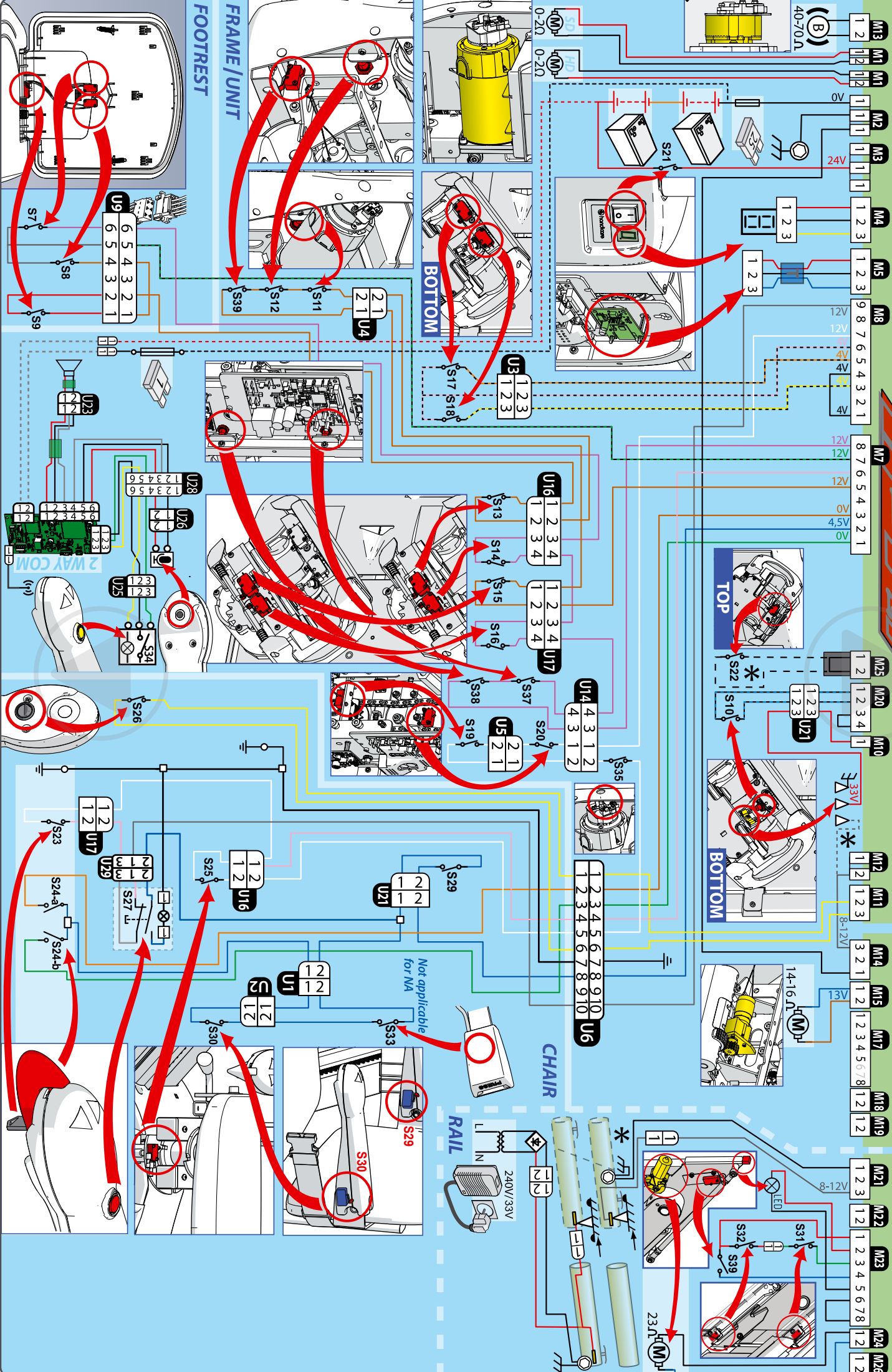
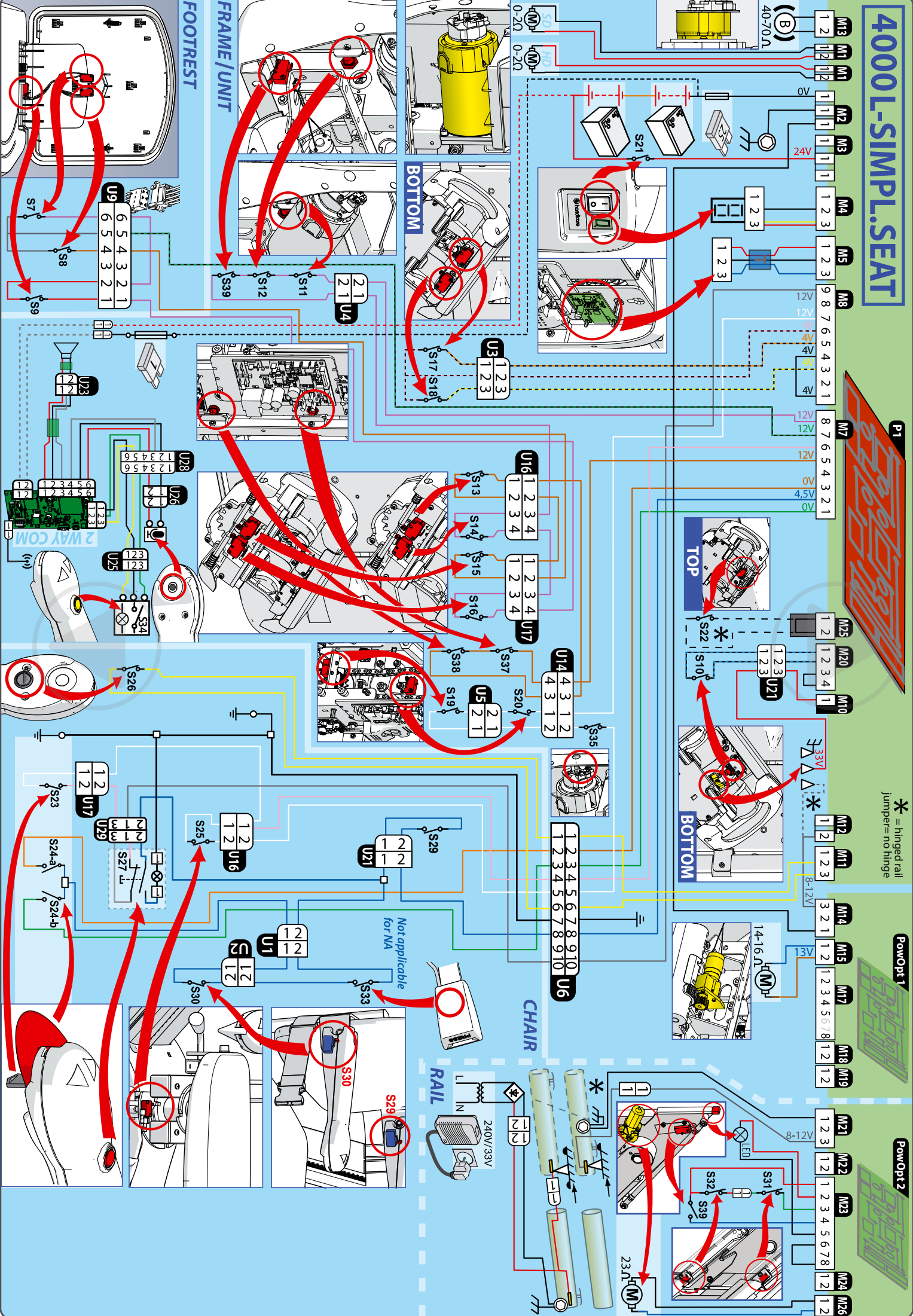
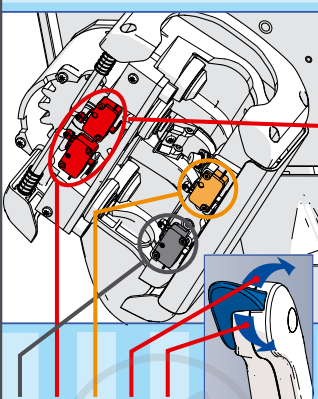
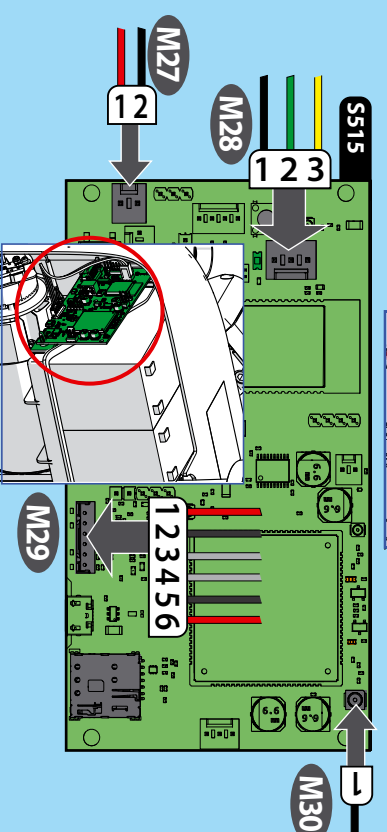
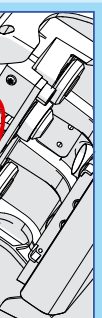
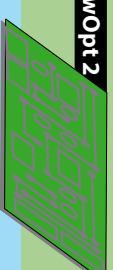
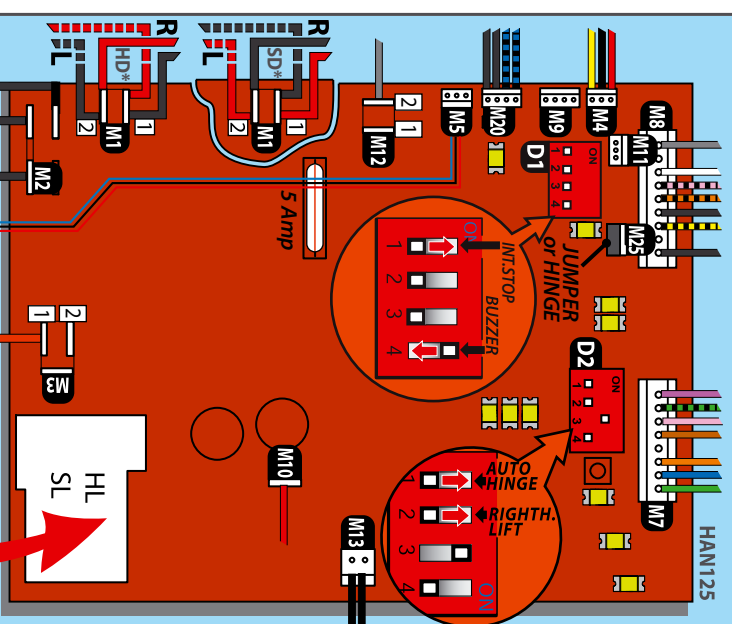
Savaria Lifts (UK) Ltd

82 First Avenue,
Pensnett Estate Kingswinford,
West Midlands DY6 7FJ,
England T +44 (0)1384 408700
enquiries@handicare.co.uk

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Savaria Canada

2 Walker Drive Brampton,
ON L6T 5E1 Canada
T 1.905.791.5555
T 1.855.728.2742
info@savaria.com

www.handicare .com

AA27496_35_290126

Date	Version	Change
10-11-24	27	Website links and QR codes
29-11-24	28	Add installation automatique hinge
23-01-25	29	Extra notification that charge contact are bottom rail mounted
27-02-25	30	Change PCB H125 and hinge illustrations
11-05-25	31	Add Outdoor stanchion fitting kit and HD seat height
03-07-25	32	Add OSG test in check chapter and reference to appendix
28-08-25	33	hinge stop bracket has no charging function
27-10-25	34	Extended endstop
28-01-26	35	Hinge adjustable foot with a mating plate

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Savaria Deutschland GmbH
Inc. 10888 Metro Court